




order on
VWR.COM

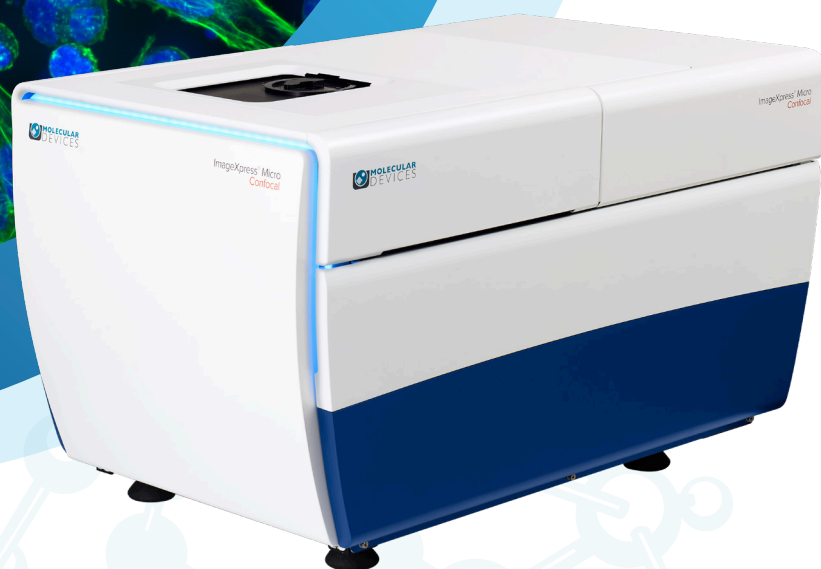


ImageXpress Micro Confocal

High-Content Imaging System


The confocal solution for your complex biology

Deeper insights into complex biology

Higher quality images, faster throughput and more powerful analysis

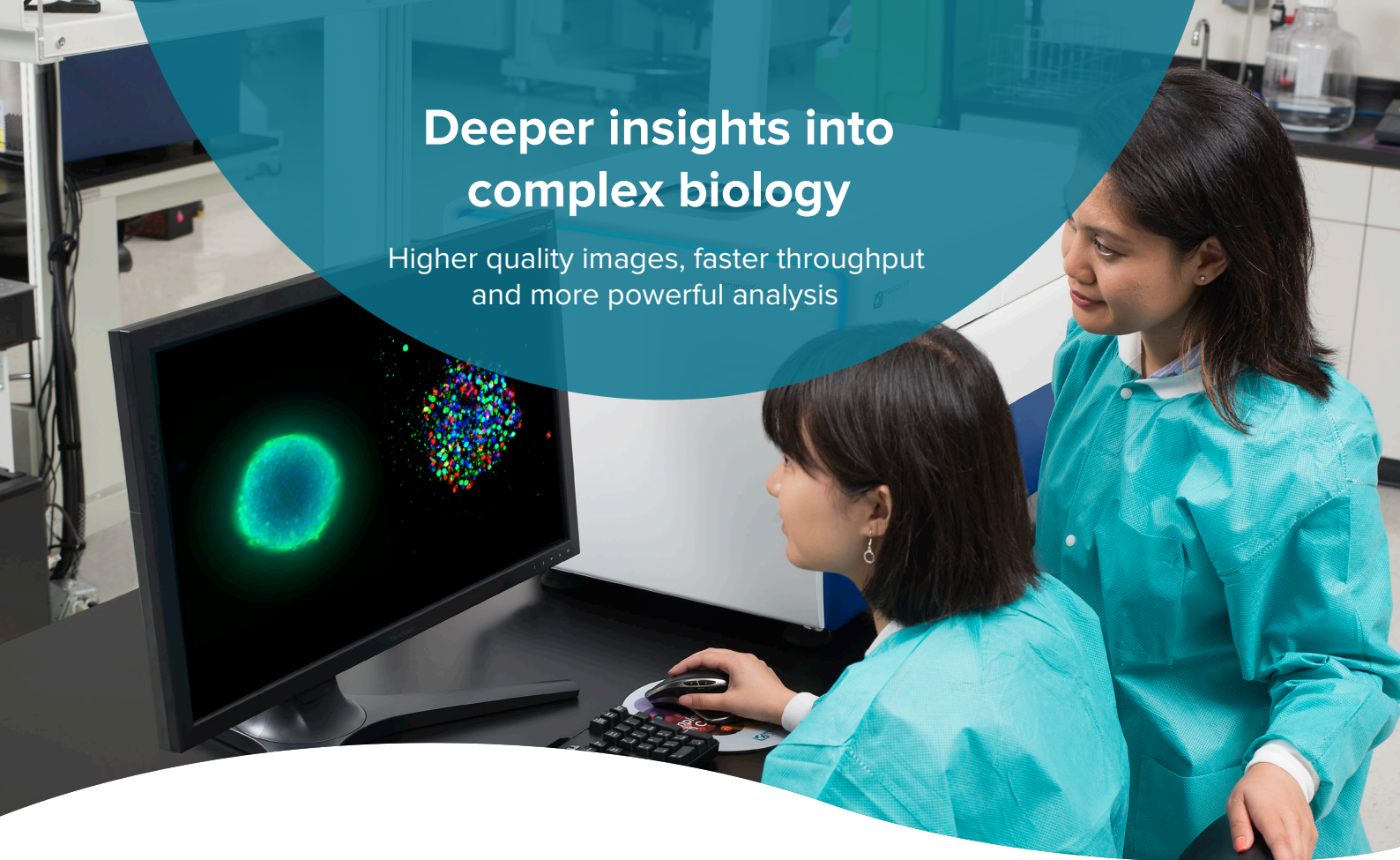

The ImageXpress Micro Confocal High-Content Imaging System provides improved quantification for live or fixed cell assays. This versatile imaging system features a unique confocal technology which allows you to explore more physiologically relevant, complex three dimensional models including spheroids, tissues, and whole organisms and to generate publication quality images at high throughput for samples in slides or one to 1536-well microplates.

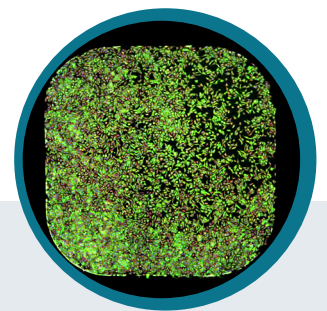

Key capabilities

- Select a confocal geometry optimal for your assay and throughput needs
- Achieve excellent image quality without sacrificing throughput via our unique optical path technology
- Acquire statistically relevant data quickly with an advanced scientific CMOS detector, enabling >3 log dynamic range
- Large field of view enables whole-well imaging
- Expand your research capabilities with water immersion objectives, transmitted light, phase contrast optics, on-board liquid handling, and environmental control options


Confocal technology at the speed of widefield imaging

- Capture an entire well of a 384-well plate with a single image at 4X magnification
- Capture four wells of a 1536-well plate in a single image at 4X magnification
- Throughput of >160K wells/day confocal, >200K wells/day widefield

AgileOptix technology at the heart of the ImageXpress Micro Confocal system

Software-selectable configurations right for your research

The ImageXpress Micro Confocal system features AgileOptix Spinning Disk Technology. Our optical options make it easy to select and configure your system to ensure the best read for your assay.

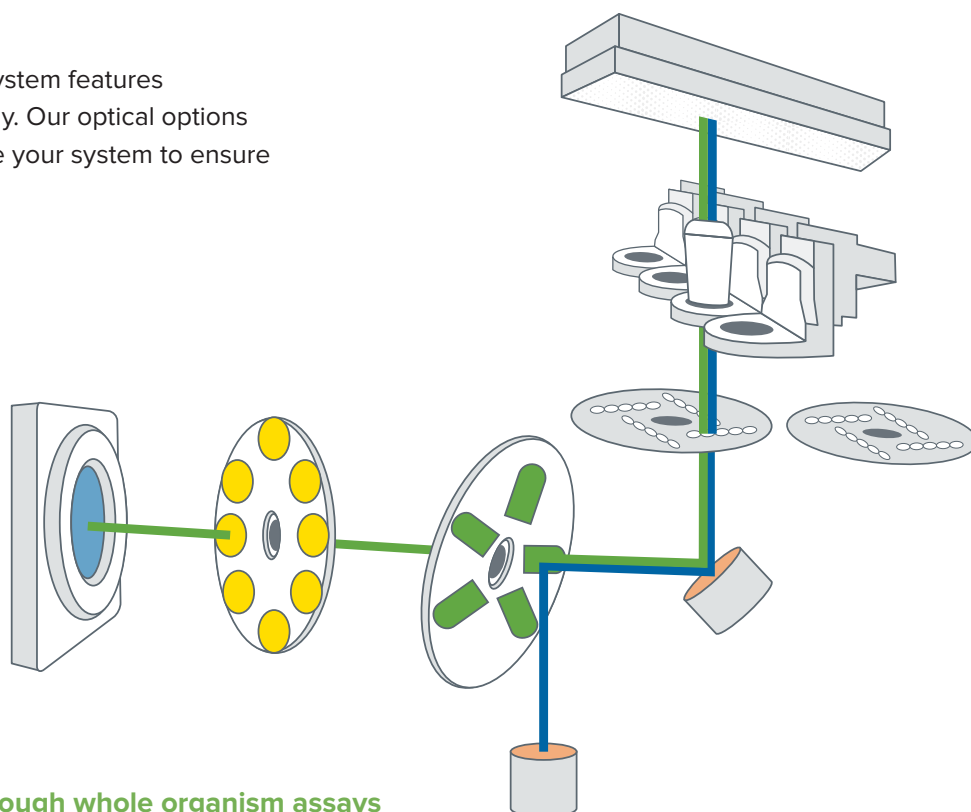


Supports subcellular assays through whole organism assays


- Widest selection (> 25) of objectives
- Oil objectives with up to 1.4 NA
- Air objectives with 0.05 to 0.95 NA
- Water objectives with up to 1.2 NA

Select a spinning disk confocal geometry matched to your assay requirements.

Spinning disk geometry	60 μm pinhole (single disk)	60 μm dual disk with 50 μm slit	60 μm dual disk with 42 μm pinhole
High-sensitivity detection	•	•	•
Fast acquisition	•	•	•
>3 log dynamic range*	•	•	•
Widfield mode for flat biology	•	•	•
Most confocal applications	•	•	•
Highest resolution imaging			•
High throughput applications		•	

*Powered by our highly responsive sCMOS sensor and advanced solid state light source.

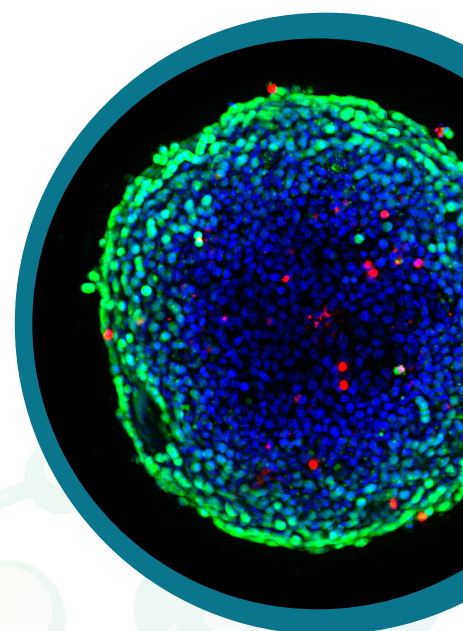


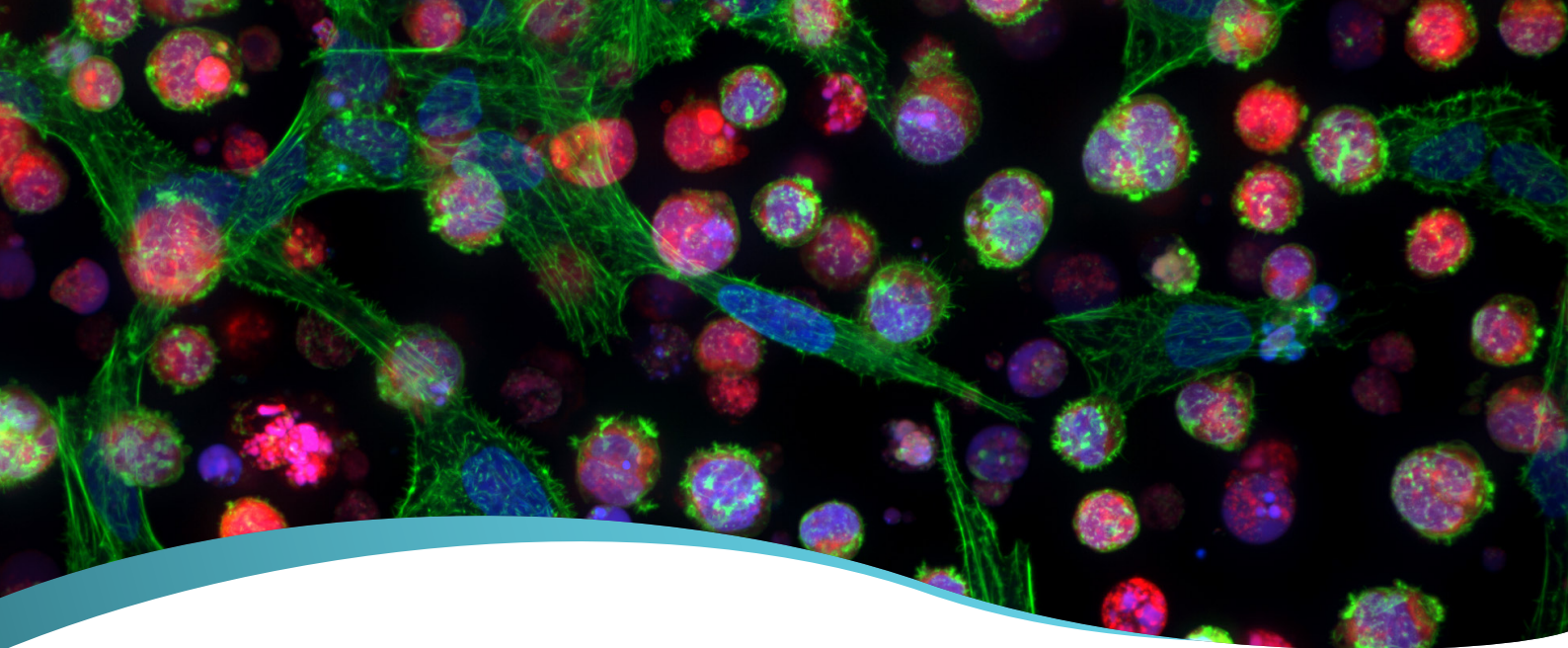

Expand your research into a new dimension

More relevant results with 3D assay models into complex biology



Complex, 3D cellular models yield more predictive, physiologically relevant results versus monocultures or other 2D cellular models. Explore the complexities of these models faster and gain better results even in samples grown in a thick extracellular matrices using the ImageXpress Micro Confocal system. It offers flexible options to meet your specific research needs and ensure the best images for your assay. Get more insight into 3D cells with options like water immersion objectives to improve signal up to 4 times at depth and decreasing optical aberrations for better data quality or monitor cell health kinetics over time using environmental control.

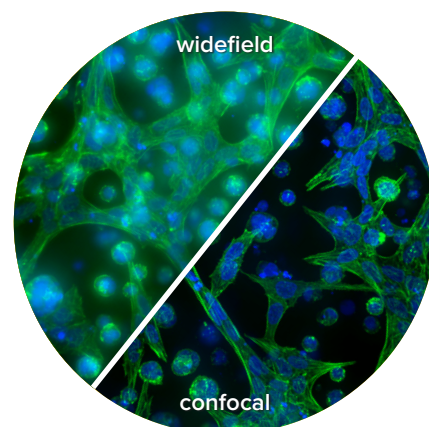

Clearer images and improved quantitative screening for:

- Spheroids
- Zebrafish and *C. elegans*
- Thick tissue samples
- Homogenous no-wash assays
- Organoids

Cell stained for nuclei and actin growing in a 3D gel. Projection image of seven planes acquired with a 40X Plan Apo objective.

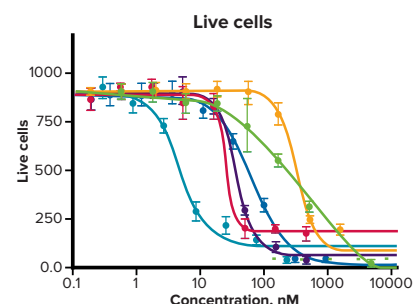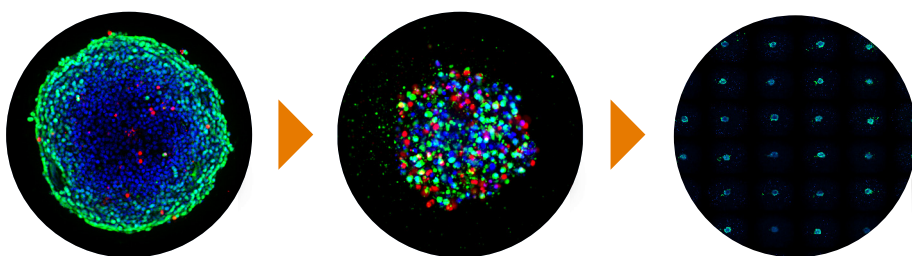

Spheroids in round bottom plates. Dose-dependent effects of selected compounds. Image montage of the 384 plate. Entire spheroid is captured with one image, example segmentation of spheroid and 4-parametric curve fits for the number of live cells in spheroids.


Improve visualization and quantitation with 3D assay models

Confocal capability improves image clarity and data quality

Compared with widefield imaging, images of thick samples captured with confocal have reduced background and improved sharpness, resulting in improved image segmentation. Further improve signal and sharpness by using water immersion objectives to capture brighter images at lower exposure times.

Enhanced imaging of tissues and 3D matrices


- Select specific cells of interest in 3D matrices such as neurons and stem cells
- Reject high-background fluorescence in thick tissue samples
- Acquire Z-stacks easily with 3D reconstruction capability

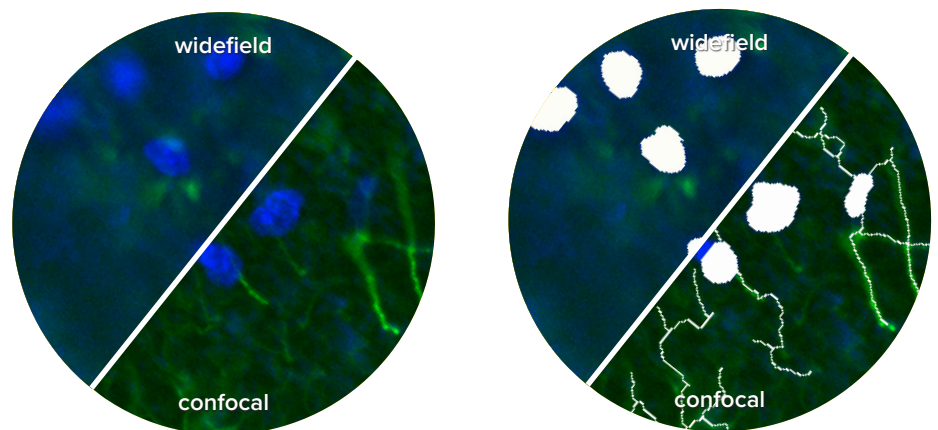

Rat brain section, stained for nuclei and neural outgrowths. Images taken with a 20X Plan Fluor ELWD objective, confocal with 60 μm pinhole. **Left** Without segmentation and **right** with image segmentation with Neurite Outgrowth module.

Clearer imaging of whole organisms

- Large field of view enables imaging of an entire well of a 384-well plate with a single image
- Organisms remain in focus for the duration of the experiment
- Perform sophisticated analysis of images and create time-lapse videos

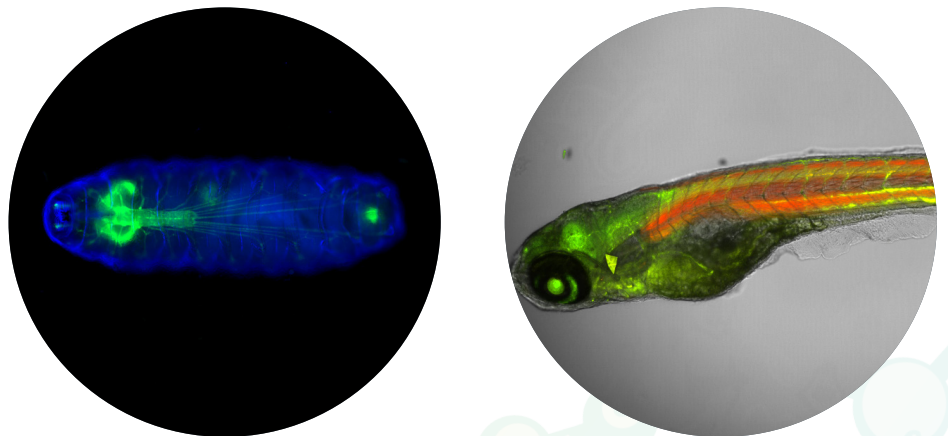

A complete solution for screening your most complex biological questions


Delivers seamless workflows from image acquisition to data analysis

Enjoy the benefits of a streamlined high-content screening (HCS) workflow in a fully integrated environment with our complete imaging solution for your most complex biological questions.


Acquire

MetaXpress High-Content Image Acquisition and Analysis Software powers our ImageXpress Micro Confocal system, giving you precise control over image acquisition and analysis, all within a unified interface.

- Acquisition wizard for the entire workflow guides protocol setup for ease of use in multi dimensional imaging
- Laser-based and software configurable image-based auto-focusing system ensures robust focus across a range of sample types
- Acquisition of live cell images enables monitoring of cell growth, death, differentiation, and migration; viral or bacterial invasion, cancer metastasis, chemotaxis, drug toxicity, or translocation

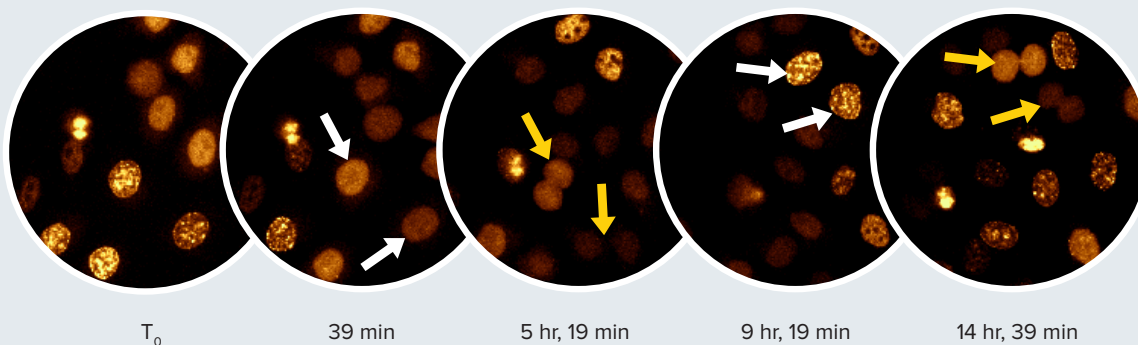

Live cell acquisition and analysis. HeLa cells expressing Cell Cycle Chromobody undergoing normal cell division while being imaged in confocal mode. In G1, the cells have a homogeneous fluorescence signal. During S phase, signal accumulates in the nucleus with formation of foci. In G2, the morphology returns to homogenous and the cell divides. White arrows indicates cell before cell division and yellow arrows indicates daughter cells after division event.

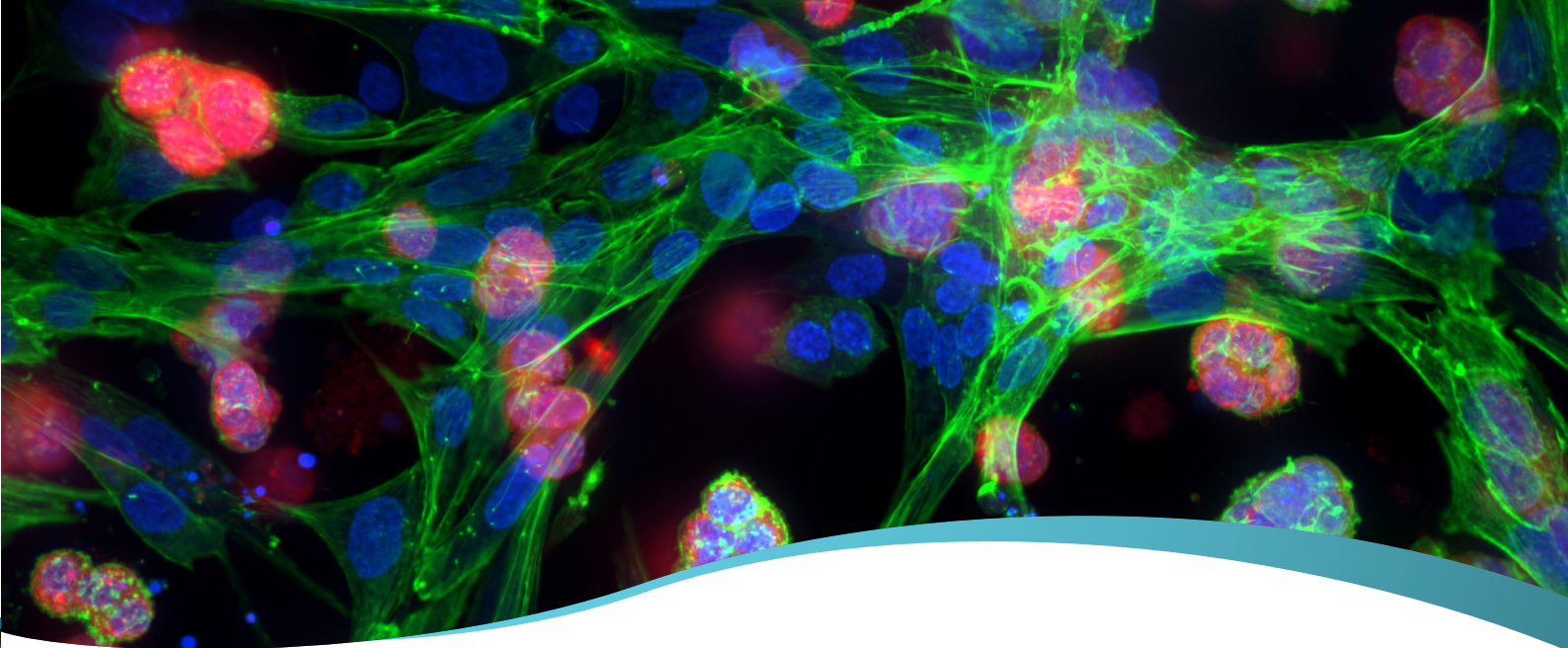

Analyze

Avoid delays in image analysis and data processing using our MetaXpress software with application modules that allow you to quickly and easily analyze your data.


- Plug-and-play application modules can be adapted to hundreds of image-based analysis workflows
- Custom module editor empowers you to further tailor your image analysis routines for a perfect fit
- Adaptive Background Correction adjusts image segmentation to the local intensity ranges and features within and between cells for better quantitation
- 2D projection algorithms include best focus, maximum and minimum, and sum projection for easy interpretation of 3D image data
- 3D volumetric analysis evaluates volume, XYZ location, distance to neighboring objects, diameter, depth, various intensity measurements, texture, or number of objects
- Save as cell-by-cell and/or image-by-image data

AcuityXpress Informatics Software, data visualization, mining, and hit selection are ready to use upon system installation.

Store

Regardless of the acquisition system used, images taken can be stored in the secure MDCStore Data Management Solution.

- Accessible for sophisticated analysis by the MetaXpress software
- Data migration portal for integration with third-party imaging systems or analysis tools to external host databases or third-party applications

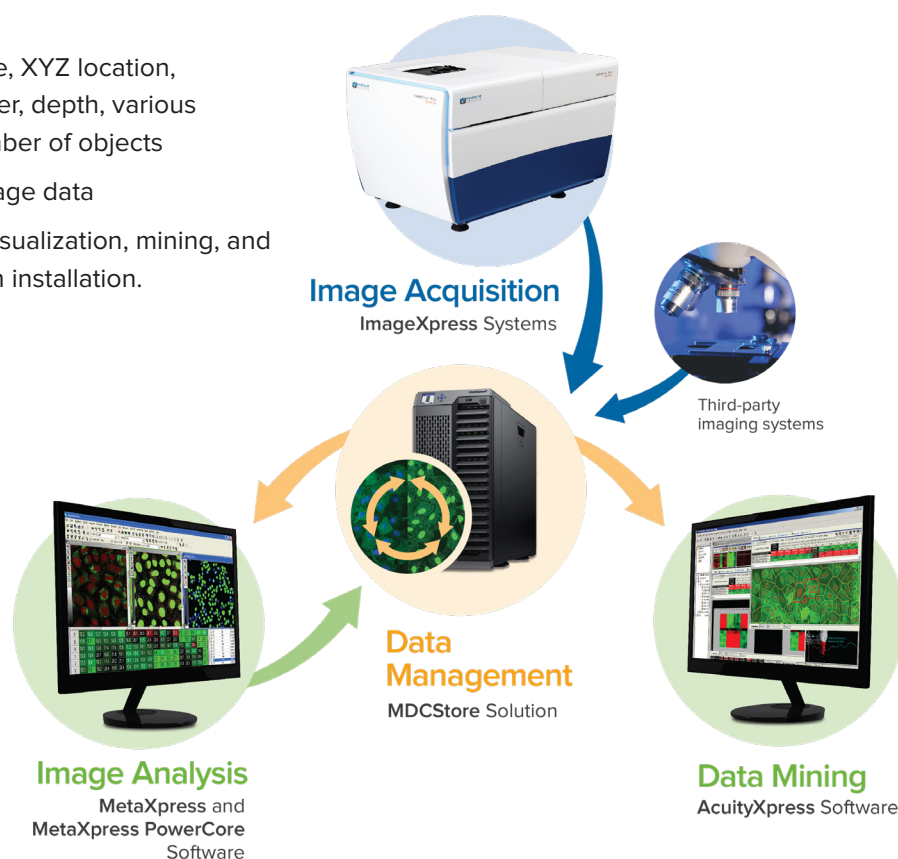

Specifications

System

- High-speed laser autofocus with integrated image autofocus option
- Linear encoded voice coil driven X, Y, and Z stages with better than 25 nm resolution
- 4-position automated objective changer*
- 5-position software selectable dichroic filter wheel*
- 8-position software selectable emission filter wheel*
- Sample compatibility: slides and one to 1536-well microplates, round or flat bottom, low to high profile

*User changeable

AgileOptix optical path

- AgileOptix technology enables the ImageXpress Micro Confocal system to deliver the sensitivity and throughput needed for demanding applications by combining a powerful solid-state light engine, high-quantum efficiency 16-bit, >4 megapixel scientific CMOS sensor, and selectable confocal geometries
- Large field of view (1.96 mm² at 10X) imaging maximizes collection of publication quality images and statistically relevant data
- >3 log dynamic range is available in both widefield and confocal modes
- Confocal can be purchased in one of the following three configurations:
 - Single-disk configuration with 60 µm confocal pinhole and widefield modes
 - High-throughput dual disk configuration with 60 µm confocal pinhole, unique and exclusive 50 µm slit confocal and widefield modes
 - High-resolution dual disk configuration with 60 µm and 42 µm confocal pinholes and widefield modes
- Laser illumination or custom light sources available

Option	Feature
Water Immersion Objectives	• 20X, 40X, and 60X (up to 1.2 NA) • Increase signal up to 4 times for brighter intensity at lower exposure times • Increase in penetration depths dependent on sample • Improve z-resolution and decrease optical aberrations • Auto water replenishment enables screening or imaging across a plate
Environmental Control	• Multi-day, live cell time-lapse imaging • Provides appropriate atmospheric conditions (e.g. 5% or 10% CO₂) • Mimics physiological environment (30–40 °C ± 0.5 °C) • Controls humidity and minimizes evaporation (0.5 µL/well/hour for 96- or 384-well formats)
Phase Contrast	• High contrast imaging where unstained cells are easily viewed or separated from background (4X–60X) • Ideal for non-fluorescent histochemically stained samples • Nikon 100W Pillar Diascopic Illuminator with TE-C ELWD Condenser • 0.3 NA with 65 mm WD and PhL, Ph1, and Ph2 selectable phase rings • Fluorophore-independent morphology visualization with fluorescent imaging overlay
On-board Fluidics	• Single-channel pipettor • Dispense volumes from 3 µL to 200 µL (±1 µL; ±5%) • Compatible with 96- or 384-well format FLIPR System pipette tips • Holds two plates for compound addition or media exchange • Optional plate heating • Environmental control

Note: all options, filters, and objectives are available at point of sale or as after market upgrades. Configurations shown herein do not encompass all configurations available. Contact your sales and support team today to identify the system configuration most suitable for your applications.

Austria VWR International GmbH Graumannsgasse 7 1150 Vienna Tel.: +43 1 97 002 0 Email: info.at@vwr.com	Belgium VWR International bvba Researchpark Haasrode 2020 Geldenaaksebaan 464 3001 Leuven Tel.: +32 (0) 16 385 011 Email: vwr.be@vwr.com	Denmark VWR International A/S Tobaksvæjen 21 2860 Søborg Tel.: +45 43 86 87 88 Email: info.dk@vwr.com	France VWR International S.A.S. Le Périgares – Bâtiment B 201, rue Carnot 94126 Fontenay-sous-Bois cedex Tel.: 0 825 02 30 30* (national) Tel.: +33 (0) 1 45 14 85 00 (international) Email: info.fr@vwr.com * 0.18 € TTC/min	Germany VWR International GmbH Hilpertstraße 20a D - 64295 Darmstadt Tel.: 0800 702 00 07* (national) Tel.: +49 (0) 6151 3972 0 (international) Email: info.de@vwr.com *Freecall
Ireland / Northern Ireland VWR International Ltd / VWR International (Northern Ireland) Ltd Orion Business Campus Northwest Business Park Ballycoolin Dublin 15 Tel.: +353 (0) 1 88 22 222 Email: sales.ie@vwr.com	The Netherlands VWR International B.V. Postbus 8198 1005 AD Amsterdam Tel.: +31 (0) 20 4808 400 Email: info.nl@vwr.com	Sweden VWR International AB Fagerstagatan 18a 163 94 Stockholm Tel.: +46 (0) 8 621 34 00 Email: kundservice.se@vwr.com	Switzerland VWR International GmbH Lerzenstrasse 16/18 8953 Dietikon Tel.: +41 (0) 44 745 13 13 Email: info.ch@vwr.com	UK VWR International Ltd Customer Service Centre Hunter Boulevard - Magna Park Lutterworth Leicestershire LE17 4XN Tel.: +44 (0) 800 22 33 44 Email: uksales@vwr.com