

OPTIONAL ACCESSORIES

MA9004	pH 4.01 buffer solution, 230 mL bottle
MA9007	pH 7.01 buffer solution, 230 mL bottle
MA9010	pH 10.01 buffer solution, 230 mL bottle
MA9015	Electrode storage solution, 230 mL bottle
MA9016	General cleaning solution, 230 mL bottle
M10000B	Rinse solution, 20 mL sachet (25 pcs.)
SE220	pH electrode with BNC connector and 1 m cable

SPECIFICATIONS

Range	0.00 to 14.00 pH
Resolution	0.01 pH
Accuracy (@ 25 °C)	±0.02 pH
Temperature	Manual setting
Compensation	0 to 50 °C
Calibration	Manual, 2-point with use of Offset and Slope trimmers
pH electrode	SE220 (included)
Environment	0 to 50 °C, 95 % RH max.
Battery Type	1 x 9 V alkaline (included)
Battery Life	approx. 70 hours of use
Dimensions	143 x 80 x 32 mm
Weight	220 g (meter with battery)

CERTIFICATION

Milwaukee Instruments conform to the CE European Directives.

Disposal of Electrical & Electronic Equipment. Do not treat this product as household waste. Hand it over to the appropriate collection point for the recycling of electrical and electronic equipment.

Disposal of waste batteries. This product contains batteries. Do not dispose of them with other household waste. Hand them over to the appropriate collection point for recycling.

Please note: proper product and battery disposal prevents potential negative consequences for human health and the environment. For detailed information, contact your local household waste disposal service or go to www.milwaukeeinstruments.com (USA & CAN) or www.milwaukeeinstruments.eu.



RECOMMENDATION

Before using this product, make sure it is entirely suitable for your specific application and for the environment in which it is used. Any modification introduced by the user to the supplied equipment may compromise the meter's performance. For your and the meter's safety do not use or store the meter in hazardous environment. To avoid damage or burn, do not perform any measurement in microwave ovens.

WARRANTY


This instrument is warranted against defects in materials and manufacturing for a period of 2 years from the date of purchase. Probe is warranted for 6 months. This warranty is limited to repair or free of charge replacement if the instrument cannot be repaired. Damage due to accidents, misuse, tampering or lack of prescribed maintenance is not covered by warranty. If service is required, contact your local Milwaukee Instruments Technical Service. If the repair is not covered by the warranty, you will be notified of the charges incurred. When shipping any meter, make sure it is properly packaged for complete protection.

Milwaukee Instruments reserves the right to make improvements in design, construction and appearance of its products without advance notice.



USER MANUAL

MW101 PRO
Portable pH Meter

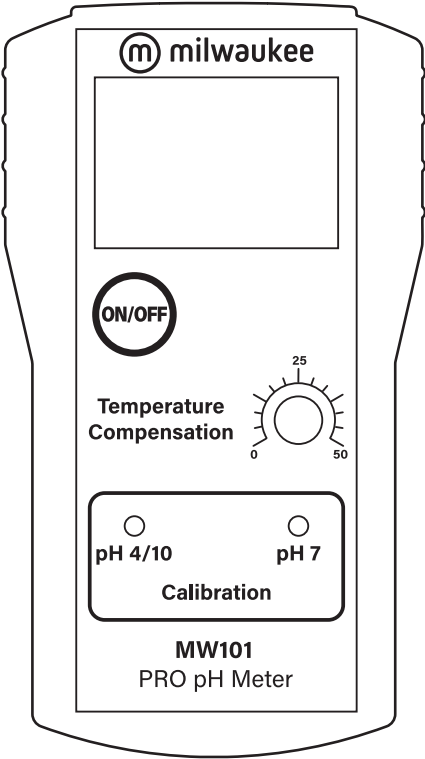

milwaukeeinstruments.com (USA & CAN)
milwaukeeinstruments.eu




PRELIMINARY EXAMINATION

Each meter is supplied complete with:


- SE220 pH electrode
- pH 7.01 Buffer solution (20 mL)
- Instruction manual
- 9 V battery (1 pc.)

OPERATION


- The meter is supplied complete with a 9 V battery. Slide off the battery compartment cover on the back of the meter. Install the battery into the battery clip connector while observing polarity.
- Always remove the electrode protective cap before taking any measurement. If the electrode has been left dry, soak the tip (≈ 2.5 cm) in **MA9015** Storage Solution for a few hours (or overnight) to reactivate.
- Connect the pH electrode to the BNC socket on the top of the meter.
- Turn the instrument on by pressing the ON/OFF key.
- Make sure that the meter has been calibrated before taking any measurements (see Calibration Procedure).

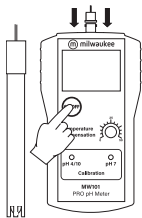

- Set the temperature knob to the value of testing solution (measured with help of an accurate thermometer).

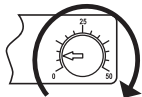

- Immerse pH electrode tip (2.5 cm) into the sample and stir gently.

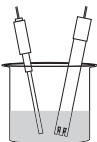

- After completing measurements, switch the meter off and store the electrode with a few drops of storage solution (**MA9015**) in the protective cap.



CALIBRATION PROCEDURE

A) Preparation:


Two calibration buffers are required

1. **pH 7.01 (MA9007)**
2. **pH 4.01 (MA9004)** if you are measuring in acid range (pH 0 - pH 7) or **pH 10.01 (MA9010)** if you are measuring in alkaline range (pH 7 - pH 14). Use two beakers for each

pH buffer. One beaker for rinsing the electrode, the other for calibration. Use a thermometer with 1 °C accuracy to measure the temperature of calibration solution.


B) Procedure:

- Remove the protective cap from the electrode. Rinse the tip of the electrode with some pH 7.01 solution, then immerse the pH electrode into a pH 7.01 buffer solution.
- Take the temperature of the buffer solution with a thermometer and set the temperature knob to the measured temperature (e.g. 15 °C).
- Adjust the OFFSET trimmer (pH 7) on the front panel, with a small screwdriver until the LCD shows the pH value at temperature of the buffer (see the pH versus temperature chart).

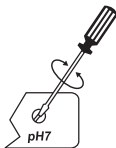

E.g. in this case, if the temperature is 15 °C, the meter display should be adjusted to read "pH 7.04".

- Now rinse the pH electrode in the first pH 4.01 beaker, then immerse it into the second pH 4.01 beaker or follow the same procedure if using pH 10.01 buffer.
- Adjust the SLOPE trimmer (pH 4/10) on the front panel, with a small screwdriver, until the LCD shows the pH value of the buffer at the temperature of measurement (see the pH versus temperature chart).

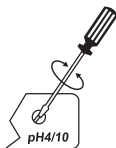

E.g. in this case, if the temperature is 15 °C, the meter display should be adjusted to read "pH 4.00" (or pH 10.01 would be adjusted to 10.12 pH). Calibration is now complete.

pH VERSUS TEMPERATURE CHART:

TEMP		pH VALUES		
°C	°F	MA9004	MA9007	MA9010
0	32	4.01	7.13	10.32
5	41	4.00	7.10	10.24
10	50	4.00	7.07	10.18
15	59	4.00	7.04	10.12
20	68	4.00	7.03	10.06
25	77	4.01	7.01	10.01
30	86	4.02	7.00	9.96
35	95	4.03	6.99	9.92
40	104	4.04	6.98	9.85
45	113	4.05	6.98	9.85
50	122	4.06	6.98	9.82
55	131	4.07	6.98	9.79
60	140	4.09	6.98	9.77
65	149	4.11	6.99	9.76
70	158	4.12	6.99	9.75

BATTERY REPLACEMENT

When the battery becomes weak the meter will display "⚡".

When the low battery indicator appears, only a few hours of battery life remain. A low battery will result in unreliable measurements. Prompt battery replacement is required.

Battery replacement must only take place in a non-hazardous area using an alkaline 9 V battery.

Turn the meter off, slide the battery compartment cover located at the rear of the meter off and replace the 9 V battery with a new one. Make sure the battery contacts are fully engaged in the connector, seat the battery in its compartment and replace the cover.

