


# Trusted Results at Your Fingertips

## Worry-free Weighing

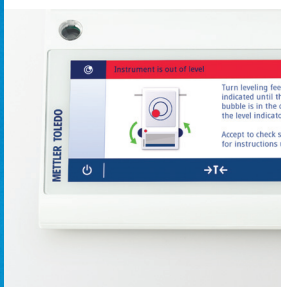

### Built-in function simplifies balance leveling

The built-in LevelControl function issues a warning when the balance is not level and provides onscreen guidance to help you level the balance correctly within seconds.



### MinWeigh function ensures process tolerances

The weight value remains red until the net sample is above the pre-programmed minimum value. A user-defined minimum weight can also be programmed into the balance to provide an additional safety factor.




### Passcodes protect balance settings

Create a passcode to restrict modification of balance settings to authorized personnel only – functions that affect metrological performance are protected against accidental alteration.




### Easy-cleaning features save time and effort

Remove, clean and replace all of the QuickLock draft shield glass panels in just a few simple steps, without tools and without moving the balance. The housing has smooth surfaces and rounded edges for easy cleaning.

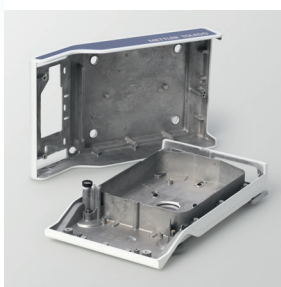

### Long balance lifetime

The full die-cast aluminum housing not only protects the weighing cell from environmental influences and impacts, it is also resistant to harsh chemicals, including acetone.

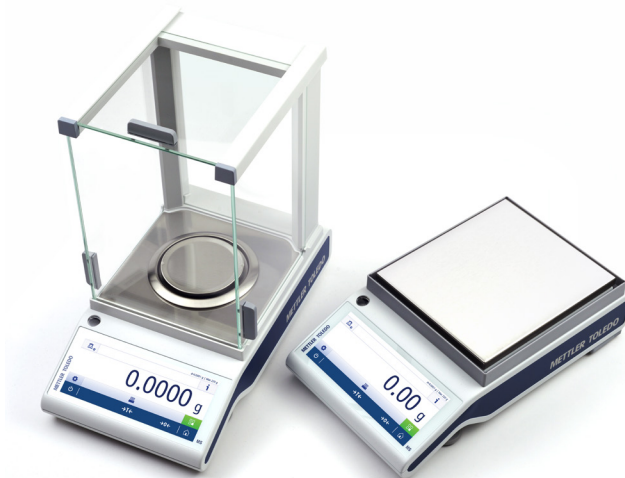

## MS-TS Analytical & Precision Balances

### Sophistication Made Simple


The built-in security features on MS-TS balances ensure results are always valid. Operators receive a warning if any problems are encountered or tolerances are not adhered to. Corrective actions can be taken immediately, avoiding costly reworking later.

Our renowned MonoBloc weighing cell, with proFACT automatic internal adjustment, delivers consistently reliable results. Built-in overload protection ensures a long balance lifetime. Features like LevelControl, MinWeigh, and passcode protection maintain security and compliance.

The 7" extra-large color TFT touchscreen display is operable through cotton, silicon and rubber gloves. An intuitive user interface and 18mm high digits bring comfort to your daily tasks.

## MS-TS Analytical Balances

	MS104TS	MS204TS	MS304TS
<b>Item Number</b>	30133522	30133523	30133524
<b>VWR Number</b>	<b>10753-566</b>	<b>10753-568</b>	<b>10753-570</b>
<b>Technical Specifications</b>			
Maximum capacity	120 g	220 g	320 g
Readability	0.1 mg	0.1 mg	0.1 mg
Tare range (from ... to)	0 ... 120 g	0 ... 220 g	0 ... 320 g
Repeatability (nominal) (sd)	0.1 mg (100 g)	0.1 mg (100 g)	0.1 mg (100 g)
Linearity deviation (test load)	0.2 mg (20 g)	0.2 mg (50 g)	0.2 mg (50 g)
Eccentricity deviation (test load) <sup>1</sup>	0.4 mg (50 g)	0.4 mg (100 g)	0.4 mg (100 g)
Sensitivity offset (test weight)	0.8 mg (100 g)	1 mg (200 g)	1.5 mg (300 g)
Sensitivity temperature drift <sup>2</sup>	1.5 ppm/°C	1.5 ppm/°C	1.5 ppm/°C
Sensitivity stability <sup>3</sup>	2.5 ppm/a	2.5 ppm/a	2.5 ppm/a
<b>Typical values</b>			
Repeatability (at low load)	0.08 mg	0.08 mg	0.08 mg
Linearity deviation	0.06 mg	0.06 mg	0.06 mg
Eccentricity deviation (test load) <sup>1</sup>	0.1 mg (50 g)	0.12 mg (100 g)	0.12 mg (100 g)
Sensitivity offset (test weight)	0.2 mg (100 g)	0.24 mg (200 g)	0.6 mg (300 g)
USP minimum sample weight (5% load, k=2, U=0.10%)	160 mg	160 mg	160 mg
Minimum sample weight (5% load, k=2, U=1%)	16 mg	16 mg	16 mg
Settling time	2 s	2 s	3 s

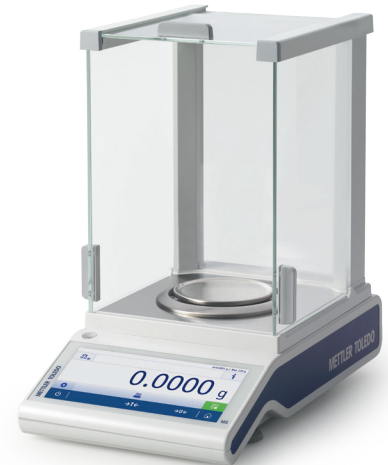

- 1) according to OIML R76;  
 2) in the temperature range 10 to 30°C,  
 3) Stability of sensitivity with proFACT self-adjustment switched on;  
 s: seconds; a: year (annum);  
 sd: standard deviation

## MS-TS Precision Balances, 1 mg

	MS303TS	MS403TS	MS603TS	MS1003TS
<b>Item Number</b>	30133525	30133526	30133527	30133528
<b>VWR Number</b>	<b>10753-558</b>	<b>10753-560</b>	<b>10753-562</b>	<b>10753-564</b>
<b>Technical Specifications</b>				
Maximum capacity	320 g	420 g	620 g	1020 g
Readability	1 mg	1 mg	1 mg	1 mg
Repeatability	1 mg	1 mg	1 mg	1 mg
Linearity deviation	2 mg	2 mg	2 mg	2 mg
<b>Typical values</b>				
Repeatability	0.7 mg	0.7 mg	0.7 mg	0.7 mg
Linearity deviation	0.6 mg	0.6 mg	0.6 mg	0.6 mg
Sensitivity offset (test weight)	3 mg (300 g)	3.5 mg (400 g)	3.5 mg (600 g)	4 mg (1000 g)
USP minimum sample weight (5% load, k=2, U=0.10%)	1.4 g	1.4 g	1.4 g	1.4 g
Minimum sample weight (5% load, k=2, U=1%)	0.14 g	0.14 g	0.14 g	0.14 g
Settling time	1.5 s	1.5 s	1.5 s	1.5 s
<b>Dimensions</b>				
Weighing pan size, WxD (mm)	127 x 127	127 x 127	127 x 127	127 x 127



## Embedded applications (All Models)

Weighing, statistics for all applications, check weighing, totaling, piece counting, formulation, percent weighing, factor weighing, dynamic weighing, density, dosing and PC Direct (for easy data transfer).

## MS-TS Precision Balances, 0.01 – 0.1 g

	MS1602TS	MS3002TS	MS4002TS	MS4002TS DR	MS6002TS	MS6002TS DR	MS12002TS	MS8001TS
<b>Item Number</b>	30133529	30133530	30133531	30133532	30133533	30133534	30133535	30133536
<b>VWR Number</b>	<b>10753-698</b>	<b>10753-700</b>	<b>10753-702</b>	<b>10753-704</b>	<b>10753-706</b>	<b>10753-708</b>	<b>10753-710</b>	<b>10753-712</b>
<b>Technical Specifications</b>								
Maximum capacity	1620 g	3200 g	4200 g	4200 g	6200 g	6200 g	12200 g	8200 g
Maximum capacity, fine range	–	–	–	820 g	–	1220 g	–	–
Readability	0.01 g	0.01 g	0.01 g	0.1 g	0.01 g	0.1 g	0.01 g	0.1 g
Readability, fine range	–	–	–	0.01 g	–	0.01 g	–	–
Repeatability	0.01 g	0.01 g	0.01 g	0.06 g	0.01 g	0.06 g	0.01 g	0.1 g
Repeatability, fine range	–	–	–	0.01 g	–	0.01 g	–	–
Linearity deviation	0.02 g	0.02 g	0.02 g	0.08 g	0.02 g	0.08 g	0.025 g	0.2 g
<b>Typical values</b>								
Repeatability	0.007 g	0.007 g	0.007 g	0.05 g	0.007 g	0.05 g	0.007 g	0.07 g
Repeatability, fine range	–	–	–	0.007 g	–	0.007 g	–	–
Linearity deviation	0.006 g	0.006 g	0.006 g	0.06 g	0.006 g	0.06 g	0.008 g	0.06 g
Sensitivity offset (test weight)	0.018 g (1600 g)	0.018 g (3000 g)	0.024 g (4000 g)	0.024 g (4000 g)	0.036 g (6000 g)	0.036 g (6000 g)	0.048 g (12000 g)	0.24 g (8000 g)
USP minimum sample weight (5% load, k=2, U=0.10%)	14 g	14 g	14 g	14 g	14 g	14 g	14 g	120 g
Minimum sample weight (5% load, k=2, U=1%)	1.4 g	1.4 g	1.4 g	1.4 g	1.4 g	1.4 g	1.4 g	12 g
Settling time	1.5 s	1.5 s	1.5 s	1.5 s	1.5 s	1.5 s	1.5 s	1 s
<b>Dimensions</b>								
Weighing pan size, WxD (mm)	170 × 200	170 × 200	170 × 200	170 × 200	170 × 200	170 × 200	170 × 200	190 × 226

## Features

<b>Accurate Results</b>	<ul style="list-style-type: none"> <li>MonoBloc weighing cell</li> <li>Strong overload protection</li> <li>proFACT internal adjustment</li> <li>MinWeigh warning function</li> <li>Full metal housing</li> </ul>
<b>Efficient Operation</b>	<ul style="list-style-type: none"> <li>7" Extra-large color TFT touchscreen</li> <li>18mm high digits</li> <li>LevelControl function</li> <li>Statistical data analysis</li> <li>Easy cleaning</li> <li>LevelLock</li> </ul>
<b>Quality Assurance</b>	<ul style="list-style-type: none"> <li>ISO-Log</li> <li>Sample ID input</li> <li>Passcode protection</li> </ul>
<b>Seamless Process</b>	<ul style="list-style-type: none"> <li>3 interfaces: USB device, USB host, RS232</li> <li>Bluetooth option</li> <li>PC Direct application</li> </ul>



## Accessories

### P-5x thermal printers



Fast, high quality printouts on paper, self-adhesive labels and continuous self-adhesive paper (including barcodes).

### Density kit



Dedicated kit for determining the density of solid and liquid substances with the built-in application.

### Bluetooth adaptors



Wirelessly send data between the balance and a PC, tablet or printer. No additional software needed.

### U-electrode




The freestanding U-Electrode removes electrostatic charges on samples and containers.


# Balance Service Qualification

Take away the guess work when you need to ensure quality and regulatory compliance. METTLER TOLEDO offers options for fast and professional balance qualification: EQ-Pac, IPac and GWP® Verification.

Benefits	Features
Save Time	No time wasted – a quick initial qualification allows you to start working with the balance – immediately.
Increase Productivity	Users receive product familiarization to ensure productive and safe working processes.
Ensure Compliance	EQ-Pacs and IPacs limit the risk of noncompliance during FDA audits and other quality management system assessments.

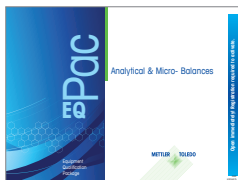

EQ-Pac covers the full range of Equipment Qualification (EQ) and is mainly used in the Pharmaceutical and Biotechnology industry. EQ-Pac is perfect for customers who prefer the entire qualification process is documented in a format designed and provided by the manufacturer. EQ-Pac includes Installation Qualification (IQ), Operational Qualification (OQ) and Performance Qualification (PQ).

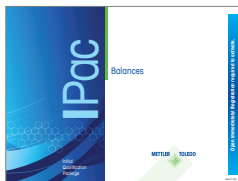

IPac ensures qualified IQ/OQ processes in all industries and is ideally suited for integration into an existing quality management system. These Pacs are designed for customers who require a professional equipment installation by the manufacturer.

Item	Catalog No.	Description
B39990001	10199-756	Basic StarterPac Installation
B39990002	10199-758	IPac Standard Qualification (IQ/OQ) – All Balance Models
B39990003	10199-760	EQPac Comprehensive Qualification (IQ/OQ/PQ) – Regulated Balances (USP Certificate)
B39990004	10199-762	EQPac Comprehensive Qualification (IQ/OQ/PQ) – Non Regulated Balances

[www.vwr.com](http://www.vwr.com)

For more information



1.800.932.5000  
vwr.com

Prices and product details are current when published; subject to change without notice. | Certain products may be limited by federal, state, provincial, or local regulations. | VWR makes no claims or warranties concerning sustainable/green products. Any claims concerning sustainable/green products are the sole claims of the manufacturer and not those of VWR International, LLC. All prices are in US dollars unless otherwise noted. Offers valid in US and Canada, void where prohibited by law or company policy, while supplies last. | VWR, the VWR logo and variations on the foregoing are registered (®) or unregistered trademarks and service marks, of VWR International, LLC and its related companies. All other marks referenced are registered by their respective owner(s). | Visit [vwr.com](http://vwr.com) to view our privacy policy, trademark owners and additional disclaimers. ©2015 VWR International, LLC. All rights reserved.