



VWR Collection Model 89511-416 Incubator cooled 0.65CF

50142486 Revision 00 December 2013



© 2013 VWR International, LLC. All rights reserved.

This Instruction Manual is copyright protected. Rights resulting thereof, particularly reprint, photomechanical or digital post-processing or reproduction, even in part, are only allowed with the written consent of VWR International, LLC.

This regulation does not apply to reproductions for in-plant use.

The contents of this document are subject to change without notice. Concerning translations into foreign languages, the English version of these operating instructions is binding.

Trademarks

All trademarks mentioned in the operating instructions are the exclusive property of the respective manufacturers.

VWR International, LLC
Radnor Corporate Center, Building One, Suite 200
100 Matsonford Road
Radnor, Pennsylvania 19087
United States

VWR International, LLC provides this document to its customers with a product purchase to use in the product operation. This document is copyright protected. This document is copyright protected and any reproduction of the whole or any part of this document is strictly prohibited, except with the written authorization of VWR International, LLC.

The content of this document are subject to change without notice.

All technical information in this document is for reference purposes only. System configurations and specifications in this document supersede all previous information received by the purchaser.

This document is not part of any sales contract between VWR International, LLC and a purchaser. This document shall in no way govern or modify any Terms and Conditions of Sale, which Terms and Conditions of Sale shall govern all conflicting information between the two documents.

Table of contents

Chapter 1	Safety Information.....	1-1
	Basic operating precautions	1-1
	Safe usage instructions	1-2
	Guarantee	1-2
	Explanation of safety information and symbols.....	1-2
	Intended use of the device.....	1-4
Chapter 2	Delivery of the device.....	2-1
	Packaging	2-1
	Acceptance inspection.....	2-1
	Scope of supply	2-1
Chapter 3	Installation.....	3-1
	Environmental conditions	3-1
	Room ventilation	3-2
	Installation.....	3-2
	Transport	3-2
Chapter 4	Description of unit	4-1
	Device overview.....	4-1
	Mains connection	4-2
	Useable space components.....	4-2
	Condensation tray.....	4-3
Chapter 5	Installation procedures.....	5-1
	Insert the perforated plates	5-1
	Preparing the useable space.....	5-1
	Mains connection	5-1
Chapter 6	Operation.....	6-1
	Preparing the device	6-1
	Starting Operation	6-1
Chapter 7	Handling and control.....	7-1
	Powering up and down.....	7-1
	Setting the temperature.....	7-2
	Interior light.....	7-2
Chapter 8	Shut-down.....	8-1
	Shutting the device down	8-1

Chapter 9	Cleaning and disinfection	9-1
	Empty the condensation tray	9-1
	Cleaning exterior surfaces.....	9-1
	Wipe / spray disinfection.....	9-1
Chapter 10	Maintenance	10-1
	Inspections and checks	10-1
Chapter 11	Disposal	11-1
	Overview of materials used.....	11-2
Chapter 12	Technical Data	12-1
Chapter 13	Accessories	13-1

Safety Information

Basic operating precautions

These operating instructions describe the VWR 18L personal incubator.

The VWR 18L personal incubator has been manufactured to the latest state of the art and been tested thoroughly for flawless functioning prior to shipping. However, the device may present potential hazards, particularly if it is operated by inadequately trained personnel or if it is not used in accordance with the intended purpose. Therefore, the following must be observed for the sake of accident prevention:

- The VWR 18L personal incubator is only to be used by trained and authorized personnel.
- The present operating instructions, applicable safety data sheets, plant hygiene guidelines and the corresponding technical rules issued by the operator shall be used to create written procedures targeted at personnel working with the subject matter device, detailing:
 - the disinfection measures to be employed for the device and the accessories used with it,
 - not to be used in a damp atmosphere,
 - the protective measures which are to be taken when processing certain agents,
 - action to be taken in case of accidents.
- Maintenance work on the device is only to be carried out by trained and authorized specialists. The device is not to be repaired. It will be repaired during the warranty period and then has to be disposed of by the customer; he having to purchase a new one should it still be required.
- The contents of these operating instructions are subject to change at any time without further notice.
- Concerning translations into foreign languages, the English version of these operating instructions is binding.
- Keep these operating instructions close to the device so that safety instructions and important information are always accessible.
- Should you encounter problems that are not detailed adequately in these operating instructions, please contact VWR International, LLC immediately for your own safety.

Safe usage instructions

- Observe the load of the stiffening plates (refer to chapter technical data).
- Do not charge the floor of the interior as there is a risk of the inserted item overheating.
- Distribute the items evenly and do not place them too close to the interior walls, so as to ensure a good temperature distribution.

Guarantee

VWR International, LLC warrants the operational safety and functions of the VWRs 18L personal incubator only under the condition that:

- the device is operated and serviced exclusively in accordance with its intended purpose and as described in these operating instructions,
- the device is not modified,
- only original spare parts and accessories that have been approved by VWR International, LLC are used.

The warranty is valid from the date of delivery of the device to the operator.

Explanation of safety information and symbols

Safety information and symbols used in the operating instructions

	 Indicates a hazardous situation which, if not avoided, will result in death or serious injuries.
	 Indicates a hazardous situation which, if not avoided, could result in death or serious injuries.
	 Indicates a situation which, if not avoided, could result in damage to property.

		NOTE Is used for useful hints and information regarding the application.
--	--	---

Additional symbols for safety information

		Wear safety gloves!
		Wear safety goggles!
		Harmful liquids!
		Electrical shock!
		Fire hazard!
		Contamination hazard!

Symbols on the device / Applicable standards



Observe operating instructions



VDE test mark

This device has been tested by the VDE on the basis of DIN EN/ISO/IEC 61010-1 and 61010-2-010 and bears the VDE GS symbol for tested safety. It is subjected to production monitoring by the VDE.



CE-conformity mark: confirms the conformity with the EU directives having validity for this product, this being confirmed in the EC Declaration of Conformity for this product.



Mark of conformity USA/Canada
This device has been tested by the Canadian Standards Association (CSA) on the basis of CAN/CSA-C22.2 No. 61010-1 and CAN/CSA-C22.2 No. 61010-2-010, UL 61010-1 and IEC 61010-2-010 and bears the cCSAus symbol for tested safety. It is subject to production monitoring by the CSA.

Intended use of the device

Correct use

The VWR 18L personal incubator is a laboratory device in the standard version and serves the preparation and cultivation of microbiological cultures and micro-organisms in safety levels L1 and L2. The device enables the special physiological ambient conditions which these cultures require to be simulated due to it exactly controlling the useable space temperature. The device is only to be used in closed spaces.

Incorrect use

To avoid the risk of explosion do not load the device with tissue, material, or liquids that:

- are easily flammable or explosive,
- release vapor or dust that forms combustible or explosive mixtures when exposed to air,
- release poisons,
- create a humid atmosphere,
- release dust,
- have an increased biological risk potential (safety level L3, L4).

The devices are not to be stacked.

	NOTE A use of the device which is not intended can be hazardous.
--	---

Delivery of the device

Packaging

The VWR 18L personal incubator are delivered in a rugged packaging box. All packaging materials can be separated and are reusable:

Packaging materials

Packaging carton: Recycled paper

Foam elements: Styrofoam (CFC-free)

Packaging film: Polyethylene

Acceptance inspection

After the device has been delivered, check the delivery immediately for:

- completeness,
- possible damage.

If components are missing or damage is found on the device or the packaging, in particular damage caused by humidity and/or water, please notify the carrier as well as VWR International, LLC Technical Support immediately.

Scope of supply

Quantity of components supplied	Pieces
Perforated plate	2
Power cord, version US/JP	1
Power cord, version EMEA/APAC	1
Mains adapter 12 V DC	1
Instruction Manual	1

Delivery of the device
Scope of supply

Installation

Environmental conditions

The VWR 18L personal incubator device must only be operated in a location that meets all of the specific ambient condition requirements listed below:

Location requirements

- Dry place of erection within a closed space which is free from draughts.
- Dust loading should be in keeping with the degree of contamination, which is Degree 2 according to IEC EN 61010-1. Using the incubator in an atmosphere with electrically conductive dust is prohibited.
- The minimal distance to adjacent surfaces must be observed on all sides (see Section 3.3).
- The operating room must be equipped with appropriate ventilation.
- Solid, level, fire-proof surface.
- Vibration-proof substructure (floor stand, lab table) capable of bearing the dead weight of the device and its accessories (particularly if two devices are stacked).
- The electrical circuitry of the device has been designed for an operating height of up to 2000 m above sea level.
- The permissible ambient temperature is between +18 °C and 32 °C (+64.4 °F and + 89.6 °F) for an incubation temperature of between 17 °C and 40 °C (62.6 °F and 104 °F).
- Relative humidity up to 80% (maximum; preferably 60-70%), non condensing.
- Dew is to be avoided, e.g. after changes of location or transportation. Should condensation exist, wait until the moisture has evaporated completely before connecting the device to a power source and powering up.
- Avoid direct exposure to sunlight.
- Devices that produce excessive amounts of heat must not be placed near the device.
- Storage temperature: When the device is placed in intermediate storage, make sure that the ambient temperature is between 20 °C and 60 °C (68 °F and 140 °F) and the maximum relative humidity does not exceed 90%, non-condensing.

Room ventilation

Heat dissipating from the device during continuous operation may cause a change in the room climate.

Therefore, the device must only be installed in rooms with sufficient ventilation.

Do not install the incubator in room recesses without ventilation.

When several devices are to be placed in the same room, additional ventilation may have to be provided as necessary.

The room ventilation is to be of a technical form which is designed in accordance with the thermal input.

Installation

When installing the device, make sure that the installation and supply connections remain freely accessible. The minimum wall clearance is 5 cm / 2 in on the sides and back. No other devices are to be placed on the device.

		NOTE The cooling fan on the back of the device is not to be covered. It is not to suck-in any loose items (e.g. paper).
--	--	---

Transport

The device is not to be lifted by the door or attachments when being transported.

Description of unit

Device overview

Front view



Figure 4-1: Front view VWR 18L personal incubator

- [1] Temperature display
- [2] Control panel
- [3] Door/door handle
- [4] Perforated plate

Rear view



Figure 4-2: Rear view VWR 18L personal incubator

- [1] Fan
- [2] Mains adapter

Mains connection

The mains connection of the device is done by connecting the enclosed mains adapter to the 100 V or 240V a.c. supply; mains fluctuation +/-10%, 50/60 Hz

Useable space components

The surface of the useable space of the VWRs 18L personal incubator is reduced to a minimum and therefore supports the protection against contamination and its easy and efficient elimination. All of the parts of the inner container of the VWR 18L personal incubator are made of plastic (the racks are made of stainless steel) and have a smooth and easy to clean surface.

A heat exchanger element [Fig. 4-3; 1] in the rear wall serves to heat or cool the useable space as required.

The useable space lighting can be switched on using the button  on the control panel.

Perforated plate

The metal plates have an anti-tipping device. The perforated plate is described in detail in the chapter entitled Commissioning.

Condensation tray

A removable tray [Fig. 4-3; 2] is provided underneath the device which collects the condensation which collects whilst the useable space is in use.

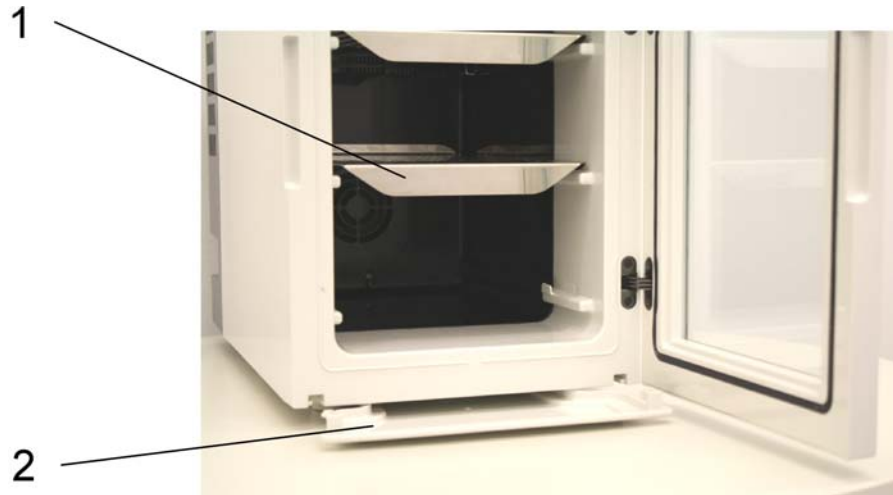


Figure 4-3: Cooling or heating elements and the condensation tray in the VWRs 18L personal incubator

Description of unit
Condensation tray

Installation procedures

Remove the packing foil.

Insert the perforated plates

Remove the transport securing devices from the perforated plates.

Insert the perforated plates in the grooved profiles on the side walls of the useable space.

Preparing the useable space

The VWR 18L personal incubator is not delivered in a sterile condition.

The following useable space components should be checked for cleanliness and disinfected prior to use:

- perforated plates,
- useable space surfaces,
- door seal,
- useable space door.

Mains connection

	 DANGER Electrical shock Contact with live electrical components may cause a lethal electric shock. Before connecting the device to the power supply, check the power cord and the plug for damage. Do not use damaged cables for connecting the device to the power supply!
---	--

Connect the device to the on site fused mains supply with the following values:

- Safety fuse (time-lag) T 16 A or circuit breaker B 16.
- Connection via a fault current circuit breaker (triggering current ≤ 30 mA) is recommended.

Connection to the power supply source

1. Before connecting the device to the power source, check to see if the power supply voltage corresponds with the specifications on the nameplate on the rear of the device. If the voltage (V) and current (A) ratings given are not as required, do not connect the device to the power source!
2. Connect the IEC connector to the socket on the rear of the device.
3. Insert the earthed plug [3] of the mains cable in a correctly earthed and fused socket.

	NOTE Power sockets The socket which the device is connected to is to be fully accessible when operating the device so as to ensure that the plug can be removed should danger occur!
--	--

4. Make sure the power cord is not subjected to tensile or compressive force.

	NOTE Condensation The device should be in the room for a sufficiently long duration before it is switched on, in order to ensure that it can reach the ambient temperature. This prevents condensation forming on live components.
--	--

Operation

Preparing the device

The device must not be released for operation before all major start-up activities have been completed (see chapter [Installation procedures](#)).

Device Check

Prior to starting operation, the following device components must be checked for their correct function:

- The door seal in the front frame must not be damaged.
- The perforated plates have to have been correctly inserted.

Starting Operation

1. Switch the device on with button .
2. The display shows the current useable space temperature.
3. Set the target values for the temperature using buttons  and  on the control panel.
4. Wait until the device has reached the set temperature. The temperature control regulates the temperature at the set target temperature value.
5. Load the useable space.

NOTE

Proper loading

To ensure sufficient air circulation and uniform heating of the samples, do not use more than 70% of the maximum surface area of the useable space. Bulky objects in the useable space that dissipate heat may impair heat distribution.

Handling and control

This chapter explains the operating elements of the VWRs 18L personal incubator in addition to providing operating instructions.

Figure 7-1: Control panel VWR 18L personal incubator

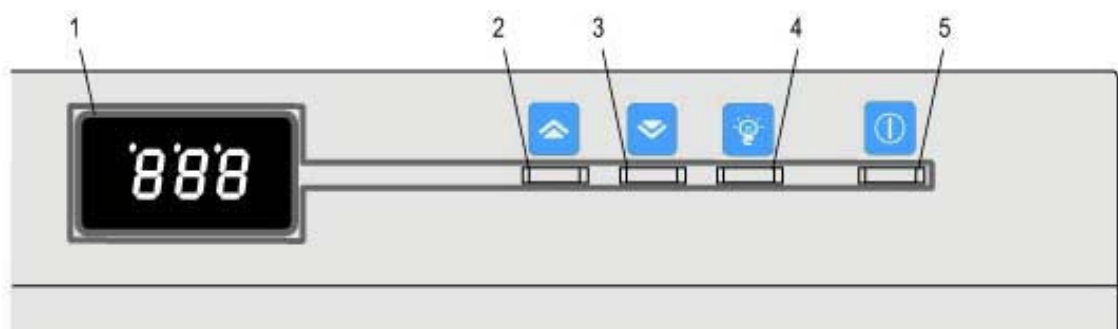


Table 7-1: Operating elements of the VWRs 18L personal incubator

No.	Designation	Function
1		Display with permanent indication of the current useable space temperature in °C in two figures (a switching over to other temperature units is not possible).
2		Increases the desired temperature by 1 °C (33.8 °F).
3		Decreases the desired temperature by 1 °C (33.8 °F).
4		Switches the interior light on or off.
5		Switches the VWR 18L personal incubator on or off.

Powering up and down

Switch the VWR 18L personal incubator on by pressing button  (Fig. 7-1; 5).

The display (Fig. 7-1; 5) is activated and shows the current useable space temperature.

To switch off the VWR 18L personal incubator, press button  (Fig. 7-1; 5) again. The device switches itself off and the display goes out.

Setting the temperature

You can preselect temperatures between 15 °C and 40 °C (+59 °F and +104 °F). The thermal output and cooling capacity can be influenced by:

- the ambient temperature,
- the number of samples which are to be heated up or cooled down,
- the frequency with which the door is opened.

Press button  (Fig. 7-1; 3), in order to reduce the required temperature by 1 °C (33.8 °F) each time or press button  (Fig. 7-1; 2), in order to increase the temperature by 1 °C (33.8 °F) each time.

The display (Fig. 7-1; 1) shows the set required temperature when you press button.

Interior light

If you wish to turn the interior light on or off, press the  (Fig. 7-1; 4) button.

Shut-down

Shutting the device down

	 WARNING	Contamination hazard
If the useable space surfaces are contaminated, germs may spread to the environment of the incubator.		

1. Remove all items and all auxiliary agents from the useable space.
2. Turn the device off using the control panel.
3. Unplug the power connector and secure it against accidental reconnection.
4. Until the incubator is shut down, the useable space must be continuously ventilated. Leave the door open and secure it against accidental closure.

Shut-down
Shutting the device down

Cleaning and disinfection

Empty the condensation tray

Empty the condensation tray from underneath the device. To do so, pull the condensation tray completely out of the device (see [Fig. 4-3](#)), when the door is open.

	 DANGER Electrical shock Contact with live electrical components may cause a lethal electric shock. Disconnect the device from the mains before carrying out cleaning and disinfection work! • Switch the device off with button  . • Remove the plug from the socket and secure it against reconnection. • Check whether the device is voltage-free.
---	---

Cleaning exterior surfaces

Remove dirt residues and depositions thoroughly using a solution of lukewarm water and commercial detergent.

Wipe the surfaces clean using a clean cloth and clear water.

Then, wipe the surfaces dry using a clean cloth.

Vacuum the fan guard on the back with a household vacuum cleaner.

Wipe / spray disinfection

The manual wipe and spray disinfection is a three-stage process:

- predisinfection,
- cleaning,
- final disinfection.

Space is to be made available for this disinfection work which is suitable for a treatment of the removed device components with cleaning agents and disinfectants.

	Incompatible cleaners CAUTION Some device components are made of plastics. Solvents may attack plastics. Strong acids or alkaline solutions may cause embrittlement of plastics. CAUTION Moisture-sensitive components Do not spray cleaning agent onto the display. When wiping over, ensure that no moisture is able to penetrate the display or run behind the door seal. Wipe the display dry with a rag made of 100% microfiber.
 	Alcoholic disinfectants! CAUTION Disinfectants having an alcohol content of more than 10% may form, in combination with air, easily combustible and explosive gas mixtures. When using such disinfectants, avoid open flames or exposure to excessive heat during the entire disinfection process! • Use such disinfectants only in adequately ventilated rooms. • After the disinfectant has been allowed to react, wipe the cleaned incubator components thoroughly dry. • Observe safety regulations to avoid fire and/or explosion hazard caused by alcohol-containing disinfectants. Chloride-containing disinfectants! CAUTION Chloride-containing disinfectants can corrode stainless steel. • Use only disinfectants that do not affect stainless steel!

 	 WARNING Health hazard The surfaces of the useable space may be contaminated. Contact with contaminated cleaning liquids may cause infections. Disinfectants may contain harmful substances. When cleaning and disinfecting, always observe the safety instructions and hygiene guidelines! • Wear safety gloves. • Wear safety goggles. • Wear mouth and respiratory system protection gear to protect your mucous membranes. • Observe the safety instructions of the disinfectant s manufacturer and the hygiene supervisor.
--	--

Predisinfection

1. Remove all samples from the useable space and store them in a safe place.
2. Spray or wipe over the surfaces of the useable space and or the perforated plates with a disinfectant.
3. Allow time for the disinfectant to act as specified by the manufacturer.

	NOTE Moisture-sensitive components Do not spray disinfectants onto the rear of the interior.
--	--

Cleaning

1. Remove the perforated plates from the useable space.
2. Wipe over the surfaces of the useable space and the removed perforated plates with tepid water containing household washing-up liquid. Also remove all stubborn stains with warm water and washing-up liquid without residues.
3. Rinse the cleaned surfaces 3-5 times with autoclaved water so as to fully remove cleaning agent residues.
4. Then rub the surfaces and cleaned perforated plates dry with a soft sterile cloth.

Final disinfection

1. Again spray or wipe over the surface of the useable space and the removed perforated plates with a disinfectant.
2. Allow time for the disinfectant to act as specified by the manufacturer.
3. Replace the perforated plates in the useable space.

Maintenance

Inspections and checks

The cooling circuit is maintenance-free and does not contain CFC.

To ensure the operational performance and safety of the device, their functions and the components listed below must be checked at regular intervals.

Daily check

Check the water level in the condensation tray underneath the device on a daily basis. In order to empty it, pull the condensation tray completely out of the device (see [Fig. 4-3](#)) when the door is open.

Annual inspection:

- Perform functional check of the control panel and of the device's temperature controller.
- Perform electrical safety check in accordance with the relevant national regulations.

Disposal

	 WARNING	Contamination hazard
The device may have been used for treating and processing infectious substances, which may have caused contamination of the device and its components. Prior to disposal, it is therefore mandatory that all incubator components be properly decontaminated! • Clean the incubator components thoroughly, then disinfect or decontaminate them (depending on application). • Attach a declaration of decontamination with details on decontamination activities performed to the items that are to be disposed of.		

Overview of materials used

Component	Material
Thermal insulation components	PU-foam
Printed circuit boards	Covered electrical components with diverse plastics, on epoxy resin-bonded printed circuit boards.
Plastic components, general	see material labeling
Exterior housing	Plastic
Device rear panel	Plastic
Outer door	Plastic
Heater	Peltier element
Inside chamber	Plastic
Perforated plates	Stainless steel
Seal	Cellular rubber
Fan impeller	Plastic
Cables	Plastic-sheathed stranded copper wire
Packaging	Corrugated board, polyethylene film, and styrofoam

Technical Data

Designation	Value
Volume	18 liter
Dimensions (W x H x D)	10.24 x 16.34 x 18.50 in
Interior dimensions	7.87 x 12.20 x 12.20 in
Weight	approx. 15.87 lb
Minimal wall distance	1.97 in on the sides and back
Number of stiffening plates	Standard 2, maximum 3
Current consumption	0.45 Amps
Protection/Installation class	II
Degree of protection	IP20
Max. mechanical load	4.41 lb per plate 11.02 lb per device
Connection voltage	100-240V AC +/-10% 50/60 Hz, Connection with the enclosed mains adapter
Power consumption	
- Maximum	45 W
- Energy consumption under no-load conditions at an ambient temperature of 23 °C (73.4 °F) and a useable space temperature of 37 °C (98.6 °F)	14 W
- Energy consumption under no-load conditions at an ambient temperature of 23 °C (73.4 °F) and a useable space temperature of 17 °C (62.6 °F)	20 W
Temperature deviations	
- spatial with an ambient temperature of 23 °C (73.4 °F)	+/- 1 °C at 37 °C (+/- 33.8 °F at 98.6 °F)
- temporal with an ambient temperature of 23 °C (73.4 °F)	+/- 0.2 °C at 37 °C (+/- 32.36 °F at 98.6 °F)
Preheat times (device empty, to 98% of the working temperature)	37 °C (98.6 °F) - 15 minutes

Designation	Value
Recovery times (device empty, door 30 s open, to 98 % of the working temperature)	37 °C (98.6 °F) - 5 minutes
On site single fuses*	- safety fuse (time-lag) T 16 A or circuit breaker B 16 - connection via a fault current circuit breaker (triggering current $\leq$ 30 mA) is recommended
Controllable temperature range in the useable space	+17...40 °C (+62.6...104 °F) with an ambient temperature of between 18 and 32 °C (64.4 and 89.6 °F)

* The national electrical engineering rules and corresponding technical connection conditions are to be adhered to when connecting to the mains supply.

Accessories

Material No.	Description
10031-744	Stainless steel perforated shelf for small cooled incubator, 18L

VWR International, LLC
Radnor Corporate Center, Building One, Suite 200
100 Matsonford Road
Radnor, Pennsylvania 19087
United States

(610) 386-1700

www.vwrsp.com