


VWR® SMART5 Networking- Quick Guide

■ Set-up Instructions

Note: It is recommended that the operating system be updated as part of this upgrade.

1. Insert the Dongle into a free USB slot (rear is recommended) and attach a network cable to the other end.

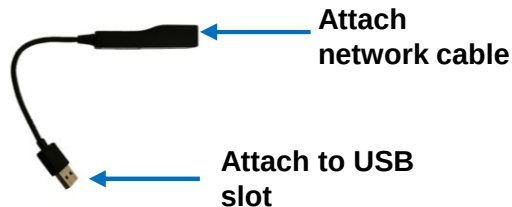


2. By default, the dongle will use DHCP to configure the network settings, but if you require a static IP, this can be set in the Settings screen.

Select




3. A pop-up window will now appear. Select 'Static' and press 'OK'.

Select

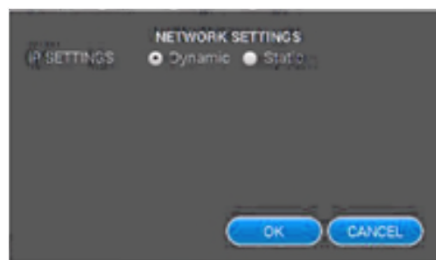


4. Once 'Static' has been selected the values can be added for the IP address, Subnet Mask and Gateway. To enter a value simply press on the red button and using the on screen keyboard type the value in. Once the values have been added press 'OK'.

Enter

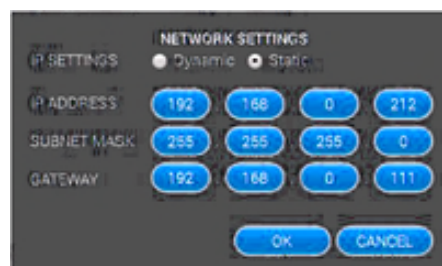


5. Now connect the Network cable.

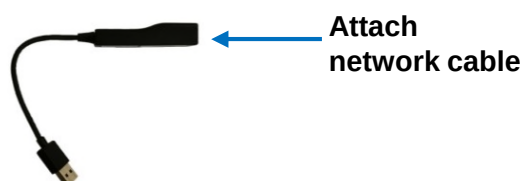

6. To save to a network location Select 'Network' in the Settings screen then press 'Location'.

Select

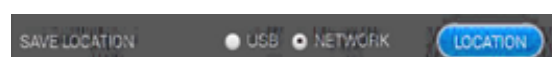

Enter



NETWORK LOCATION

SERVER *

SHARE *

USERNAME *

PASSWORD *

DOMAIN (none)

OK CANCEL

7. Fill out valid details and select 'OK'. This will attempt to connect to the server using the values you have provided.

Note: This will connect to Windows Servers. If you do not have a domain controller, leave DOMAIN as blank (will show (none)).

■ Updating Instructions

1. If downloading the update, unzip the contents of the zip file to a USB device that is formatted with FAT32.
2. Insert the USB device into the USB port by the side of the screen and turn the unit on.
3. When the system has booted, the user will be asked if they wish to update the software, click OK to update.
4. If the operating system can be upgraded and requires upgrading, a second dialog will ask if they wish to upgrade, click OK. (**N.B.** this upgrade takes ~5 minutes to complete).
5. Once the upgrade is complete VWR® SMART5 software will start.

CONTACT VWR®:
For further details , please consult the full Instruction Manual
Visit www.vwr.com to find your local VWR® contact