


170.0 mm



128.0 mm


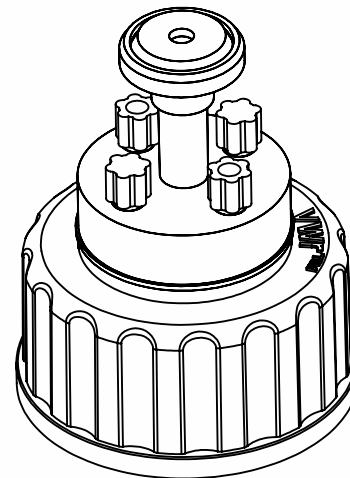


VWR®
Safety Cap

STEP 1

Port Configuration

Insert the solvent tubing through the safety cap fitting and then push the tubing through the port opening of Safety Cap. If available, continue to add as many solvent tubes through the available Safety cap port openings as needed.

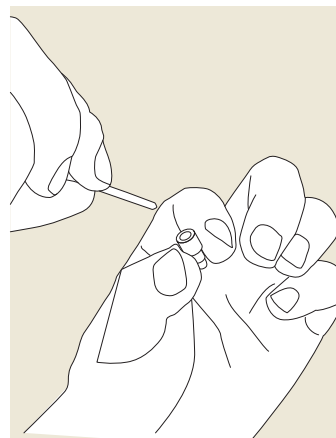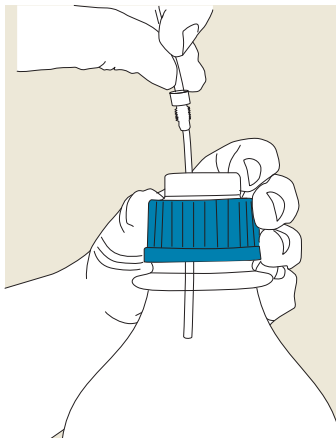

Install the Safety Cap

Hold the Safety Cap and insert the solvent tubing into solvent container. Rotate the Safety Cap clockwise onto the solvent bottle

Secure Tubing

Extend tubing so the solvent tubing reached the bottom of bottle. Rotate the fitting cap as tight as possible.

STEP 2


Safety Filter Valve Installation

Insert the Safety filter Valve along with blue sticker into the port opening on top of the Safety Cap. Turn the Screw by hand clockwise as far as it goes.

Mention the date on the given blank area (white).

This will serve as visual reminder of when to replace safety filter valve.

We recommend changing it every six months.

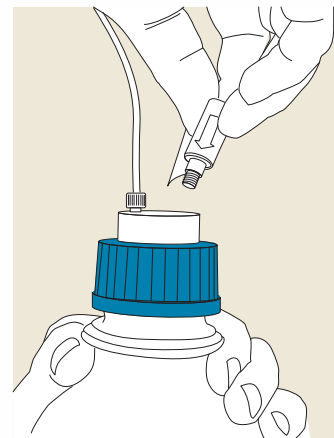

Ordering information

MEL Number	Pack Size	Pack Unit	Description
VWRA590-0057	4	EA	Safety Cap with filter valve 1 port & 1 filter
VWRA590-0058	4	EA	Safety Cap with filter valve 2 port & 1 filter
VWRA590-0059	4	EA	Safety Cap with filter valve 3 port & 1 filter
VWRA590-0060	4	EA	Safety Cap with filter valve 4 port & 1 filter
VWRA590-0061	4	EA	Safety Cap with filter valve 5 port & 1 filter
VWRA590-0062	4	EA	Filter valve with 1 filter & 4+4 blocks



Avantor Performance Materials India Limited
17th Floor, Building No. 5, Tower C
DLF Cyber City Phase - III,
Gurgaon - 122002, Haryana