
Bead Mill

Homogeniser

Quick Start Guide

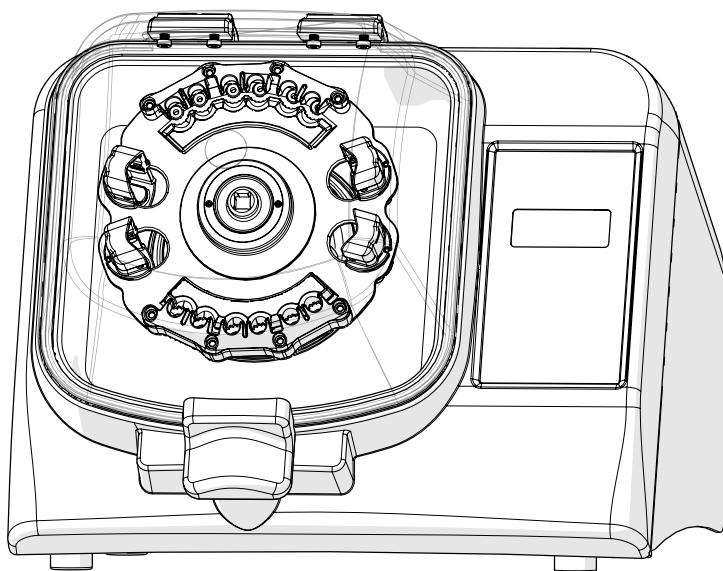
Catalog No	Description
432-0367	Bead Mill 230V, 300W


Warning / Safety Information

- ⚠ **DO NOT** attempt to modify any part of this product.
- ⚠ **DO NOT** allow the machine to be submerged in any liquid.
- ⚠ **DO NOT** use in any setting other than an indoor laboratory.
- ⚠ **DO NOT** plug power cord into an incorrect outlet or subject it to an incorrect voltage.
- ⚠ Use this product only for its intended purpose.
- ⚠ **DO NOT** use attachments not recommended by the manufacturer.
- ⚠ **DO NOT** operate the product if it is damaged in any way.
- ⚠ Keep this product away from heated surfaces.
- ⚠ **DO NOT** modify the plug or cord that is provided. If the plug will not fit the outlet, have the proper outlet installed by a qualified electrician.
- ⚠ **DO NOT** operate the product with the safety ground disconnected.
- ⚠ **WARNING:** Reduce the risk of unintentional starting; make sure the speed switch is in the OFF position before plugging in the motor.
- ⚠ **WARNING:** Damaged or worn power cords should be repaired or replaced immediately by a qualified electrician.
- ⚠ **WARNING:** Improper connection of the equipment can result in a risk of electric shock.

Unpacking

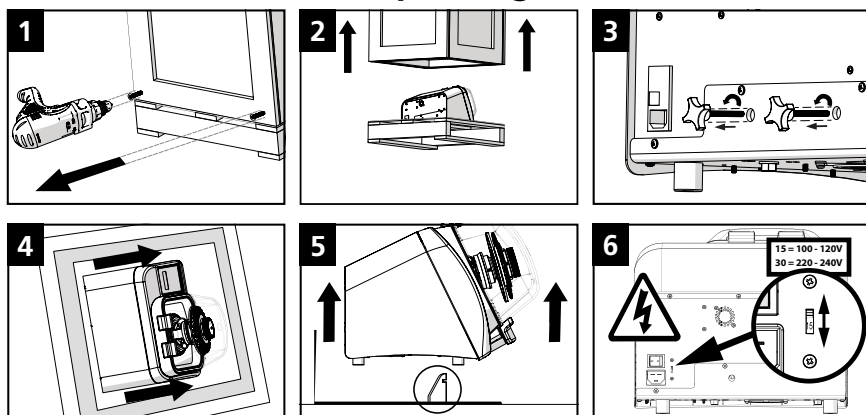
Voltage

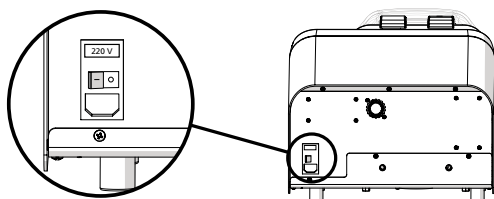Remove Warning stickers from the back of the unit and check that the voltage matches the voltage delivered by your electric power supply.

- ⚠ **CAUTION:** DO NOT connect the unit to the main power supply until correct voltage is selected. Plug the unit into power supply using the supplied power cord.
- ⚠ **CAUTION:** If the voltage supplied to the unit is higher than what is indicated on the voltage switch, the electronics can be irreversibly damaged, voiding the warranty.

General Operation

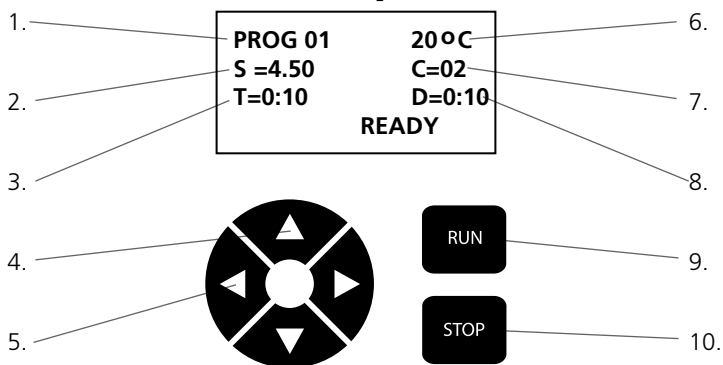- | | |
|----------------------------------|--|
| 1. Program Number | 6. Processing Chamber Temperature |
| 2. Speed | 7. Number of Cycles |
| 3. Cycle Time | 8. Dwell or Time interval between cycles |
| 4. Navigation Buttons | 9. Run Button |
| 5. Parameters Adjustment Buttons | 10. Stop Button |

Loading Sample Tubes

-
- A minimum of four tubes are required and must be spaced evenly*.
-
- Screw down to secure tubes.
- Move the locking switch to the "LOCK" position.
- Close the lid until the locking handle clicks into place.

⚠ CAUTION: DO NOT run without moving the switch to the "Lock" position.

⚠ CAUTION: The unit will not run without the lid fully closed.

For more details, please consult the full manual.

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VWR International AB
VWR International AG
VWR International Ltd