

Precision balance with Graphic display

User Manual

European Catalogue Numbers:

611-2266	ANAL. BALANCE LAG124 120G/0.1MG DISP	611-2304	ELECTR. PRECIS.BALANCE 1500G/0.01G DIP
611-2267	ANAL..BALANCE LAG124i 120G/0.1MG DISP	611-2305	ELECTR. PRECIS.BALANCE 2100G/0.01G DIS
611-2268	ANAL. BALANCE LAG164i 160G/0.1MG DISP	611-2306	ELECTR. PRECIS.BALANCE 3100G/0.01G DIS
611-2269	ANAL.BALANCE LAG 214i 220G/0.1MG DISP	611-2307	ELECTR. PRECIS.BALANCE 2100G/0.01G DIS
611-2270	ANAL. BALANCE LAG254i 250G/0.1MG DISP	611-2310	ELECTR. PRECIS.BALANCE 6500G/0.1G DIS
611-2271	ANAL. BALANCE LAG314i 310G/0.1MG DISP	611-2311	ELECTR. PRECIS.BALANCE 12000G/0.1G DIS
611-2280	ELECTR. PRECIS.BALANCE 520G/0.001G DIS	611-2371	ELECTR. PRECIS.BALANCE 4200G/0.01G DIS
611-2281	ELECTR. PRECIS.BALANCE 720G/0.001G DIS	611-2372	ELECTR. PRECIS.BALANCE 16000G/0.1G DIS
611-2282	ELECTR. PRECIS.BALANCE 520G/0.001G DIS	611-2373	ELECTR. PRECIS.BALANCE 20000G/0.1G DIS
611-2283	ELECTR. PRECIS.BALANCE 720G/0.001G DIS	611-3265	ANAL. BALANCE LAG124i 120G/0.1MG DISP
611-2284	ELECTR. PRECIS.BALANCE 520G/0.001G DIS	611-3266	ANAL. BALANCE LAG164i 160G/0.1MG DISP
611-2285	ELECTR. PRECIS.BALANCE 720G/0.001G DIS	611-3267	ANAL.BALANCE LAG 214i 220G/0.1MG DISP
611-2286	ELECTR.PRECIS.BALANCE 520G/0.001G DIS	611-3268	ANAL. BALANCE LAG254i 250G/0.1MG DISP
611-2287	ELECTR. PRECIS.BALANCE 720G/0.001G DIS	611-3273	ELECTR. PRECIS.BALANCE 520G/0.001G DIS
611-2290	ELECTR. PRECIS.BALANCE 2200G/0.01G DIS	611-3274	ELECTR. PRECIS.BALANCE 720G/0.001G DIS
611-2291	ELECTR. PRECIS.BALANCE 4200G/0.01G DIS	611-3275	ELECTR. PRECIS.BALANCE 520G/0.001G DIS
611-2295	ELECTR. PRECIS.BALANCE 210G/0.001G DIS	611-3276	ELECTR. PRECIS.BALANCE 720G/0.001G DIS
611-2297	ELECTR. PRECIS.BALANCE 210G/0.001G DIS	611-3278	ELECTR. PRECIS.BALANCE 4200G/0.01G DIS
611-2303	ELECTR. PRECIS.BALANCE 1000G/0.01G DIS	611-3297	ANAL.BALANCE LAG314i 310G/0.1MG CERT

Legal Address of Manufacturer**Europa**

VWR International bvba
Researchpark Haasrode 2020
Geldenaaksebaan 464
B-3001 Leuven
+ 32 16 385011
<http://be.vwr.com>

Country of origin

Italy

Intended use

The precision balances are used for precision scale in laboratories , industry , universities and wherever required weight measurement accurate and repeatable .

INDEX

1	INSTALLATION INSTRUCTIONS	3
2	STORAGE CONDITIONS	5
3	WEIGHING PLATE ASSEMBLY	6
4	KEYBOARD AND DISPLAY	7
5	OVERVIEW	8
5.1	FRONT PART OF THE BALANCE	8
5.2	REAR AND BOTTOM PART OF THE BALANCE	9
5.3	REAR AND BOTTOM PART OF THE BALANCE, MODELS VERSION PLANE SQUARE	10
6	WEIGHING	11
6.1	STAND BY	11
6.2	SIMPLE WEIGHING	11
7	CALIBRATION	12
7.1	BALANCES WITH EXTERNAL CALIBRATION	12
7.1.1	<i>External calibration</i>	12
7.2	BALANCES WITH INTERNAL CALIBRATION	14
7.2.1	<i>Internal calibration</i>	14
7.3	VERIFICATION	15
8	TARE FUNCTION	16
8.1	MANUAL TARE FUNCTION	17
8.2	OPTIONAL EXTERNAL ALPHANUMERIC KEYBOARD (CODE T201)	17
9	BALANCE PARAMETERS SETUP MENU	19
9.1	WEIGHT UNITS	20
9.2	SERIAL OUTPUT SETUP	21
9.3	TRANSMISSION SPEED SELECTION	22
9.4	AUTOZERO FUNCTION	23
9.5	FILTER SELECTION	24
9.6	STABILITY FUNCTION	25
9.7	CONTRAST ADJUSTMENT	26
9.8	BACKLIGHT REGULATION	27
9.9	TIMER-OFF FUNCTION	28
9.10	DATE AND TIME REGULATION	29
9.11	LANGUAGE SELECTION	30
9.12	CALIBRATION MODE SETTING	31
9.12.1	<i>Automatic Calibration (AUT-CAL)</i>	31
9.12.2	<i>Internal calibration (I-CAL) [locked for verified balances]</i>	32
9.12.3	<i>External calibration (E-CAL) [locked for verified balances]</i>	32
9.12.4	<i>Technical calibration (TEC-CAL) [locked for verified balances]</i>	33
9.13	CALIBRATION DATA	34
10	BALANCE PROGRAMS MENU	35
10.1	PIECE COUNTING FUNCTION	36
10.1.1	<i>Manual insertion of the average unit weight</i>	37
10.1.2	<i>Automatic updating of the average unit weight</i>	39
10.2	PROGRAM FOR THE DETERMINATION OF THE DENSITY OF A SOLID OR A LIQUID	40
10.2.1	<i>Solid density determination</i>	40
10.2.2	<i>Liquid density determination</i>	42
10.3	FORMULATION FUNCTION	44

10.3.1	Manual formulation.....	44
10.3.2	Formula saving.....	45
10.3.3	Formula recall.....	47
10.4	MAX-MIN THRESHOLDS FUNCTION.....	49
10.4.1	With both the limits set.....	50
10.4.2	With only the lower limit set.....	50
10.4.3	With only the upper limit set.....	50
10.5	PERCENTAGE WEIGHING FUNCTION.....	51
10.5.1	Automatic mode with reference weight.....	51
10.5.2	Mode with manual insertion of the reference weight.....	52
10.6	ANIMAL WEIGHING FUNCTION.....	53
10.7	MAXIMUM LOAD FUNCTION.....	54
10.8	GLP FUNCTION (GOOD LABORATORY PRACTICES).....	55
11	RS232 INTERFACE FEATURES.....	57
11.1	GENERAL FEATURES.....	57
11.2	MAP OF CONNECTOR.....	57
11.3	CONNECTION OF THE BALANCE TO THE OPTIONAL ALPHANUMERIC KEYBOARD.....	57
11.4	CONNECTION OF THE BALANCE TO COMPUTER.....	58
11.4.1	Continuous Transmission mode.....	59
11.4.2	On demand transmission mode.....	59
11.4.3	On demand transmission with G.L.P.....	61
11.5	CONNECTION OF BALANCE WITH SERIAL PRINTER.....	62
11.5.1	PRINT FORMATS.....	62
11.5.2	Generic Printer or TLP 50 printer with G.L.P.....	64
12	ERROR CODES.....	65
13	MAINTENANCE AND CARE.....	66
14	WARRANTY.....	66
15	QUICK GUIDE TO BALANCE PARAMETERS SETUP.....	67
16	BALANCE TECHNICAL FEATURES.....	69
17	OPTIONAL ACCESSORIES.....	70
18	TECHNICAL SERVICE.....	70
19	WARRANTY.....	71
20	COMPLIANCE WITH LOCAL LAWS AND REGULATIONS.....	71
21	EQUIPMENT DISPOSAL.....	71



1 Installation instructions



WARNING:

Please carefully read these installation and use instructions before you begin working with your new balance. The product's safety and accuracy can no longer be guaranteed if the instrument is used in a way that is different from that recommended in this manual.



- **Remove** the balance and its calibration weight (if supplied) from the carton and check for any visible damage to the instrument.
- **Do not install** the balance in environments in which there are drafts, strong temperature changes, or vibrations.
- **Do not use** the balance in potentially explosive environments.
- **The humidity of the environment** where the balance is used must be between 45% and 75%.

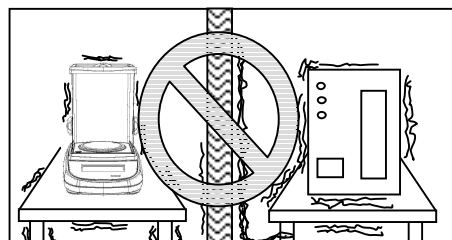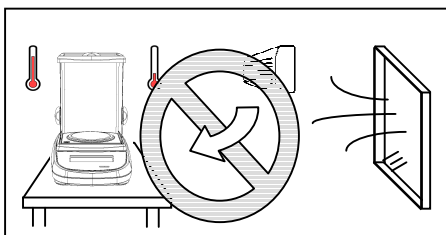

ATTENTION

In models with internal calibration remove the screw in the left lower part of the balance (see par5)

- **Place** the underplate and the plate on the balance (see page 6).
- **Level** the balance by adjusting the feet at its front.

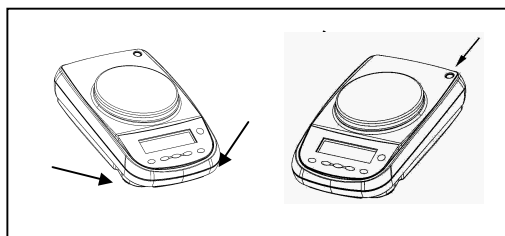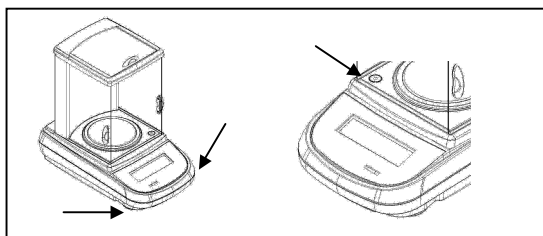

- **Connect** the power supply to connector 2 at the back of the balance (see page 7).
- **Connect** the power supply to a nearby outlet, **which must be easily accessible**; after a few seconds the balance will turn on automatically.

- **Balance Warm-up:**

- Wait 8 hours from switch one for Analytical balances (0,0001g)
- Wait 30 minutes for other models

Then calibrate the balance using the appropriate mass (if supplied), following instructions (par.7).


It is strongly suggested not to disconnect the balance from the electrical mains and to use the ON/OFF button to stand-by the balance when you finished using it

- **Calibrate** the balance every time that it is moved to another place.
- **Periodically check** the balance's calibration.
- **Do not** drop objects of an excessive weight on the balance's plate so as to avoid damaging it.
- **The assistance service** must be carried out by specialized personnel and the spare parts used must be original.
For this reason you must contact the retailer where the purchase was made.


2 Storage conditions

- **Storage temperature** +5 °C...+40°C
- **Storage humidity** 45% - 75%.
- **Keep the balance's packaging** in case it needs to be sent for servicing; disconnect all cables and any accessories to prevent unnecessary damage.
- **Do not expose** the balance to extreme temperatures and humidity unnecessarily, and avoid violent shocks.

3 Weighing plate assembly

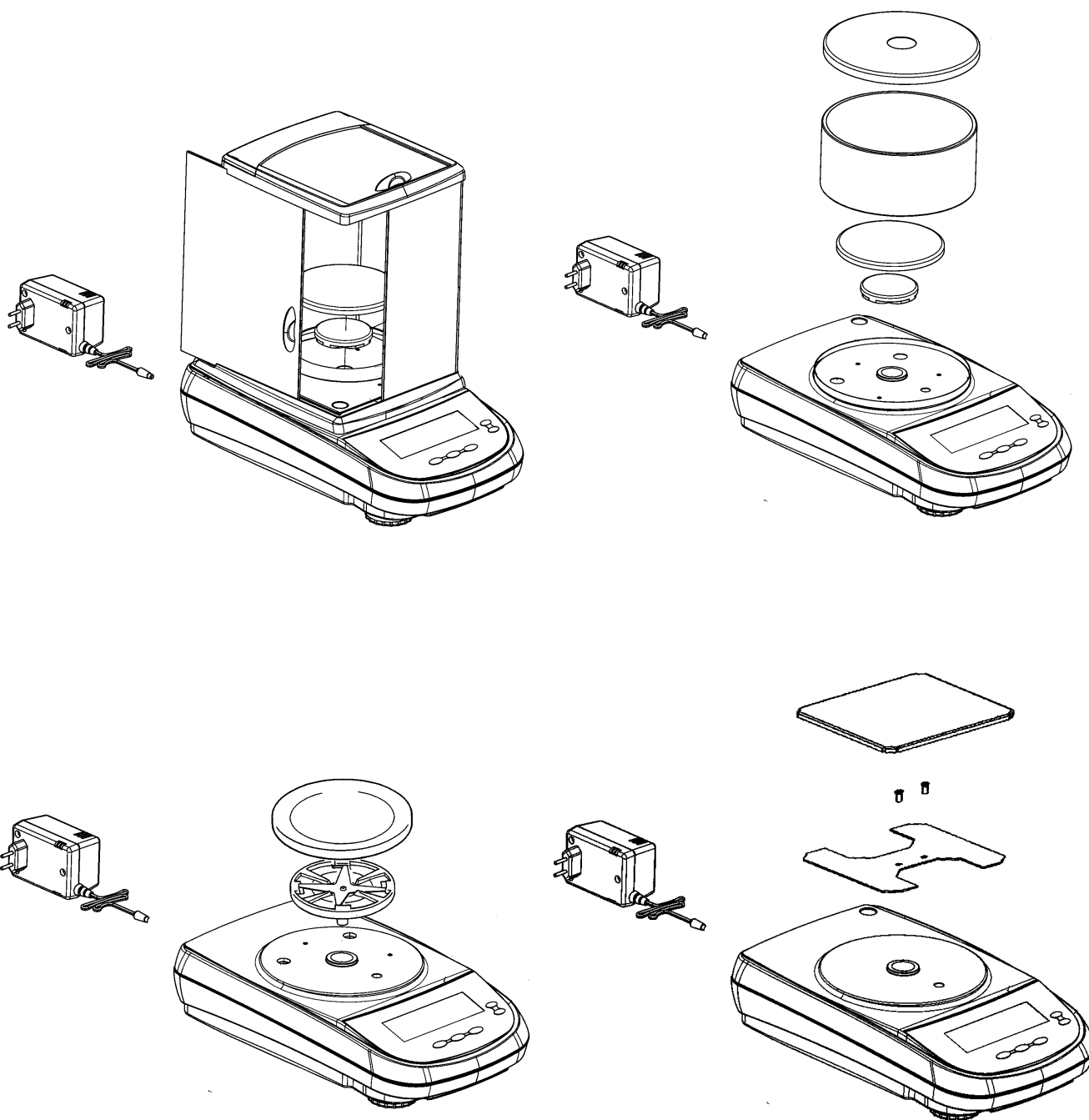

4 Keyboard and display

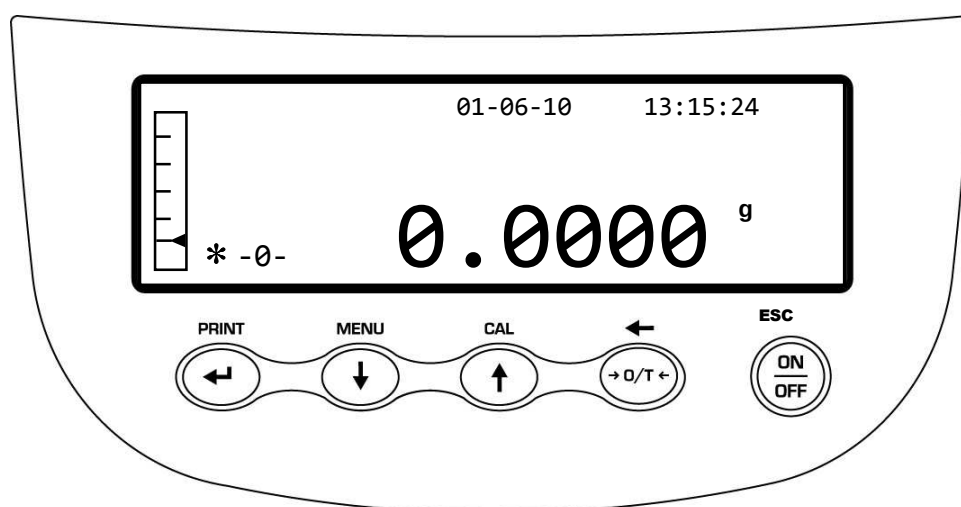

	Standby (OFF) or operating (ON) button or ESC		Stability indicator
	TARE and zeroing button		Zero indicator
	Selection CONFIRM or SEND data to the printer button		Percent weight
	Button for accessing the MENU to set the balance's parameters		Piece counting
	Balance CALIBRATION button		Battery indicator
			Data insertion mode
			Upper threshold
			Lower threshold
			Density measurement
			Unit of measurement

5 Overview

5.1 Front part of the balance

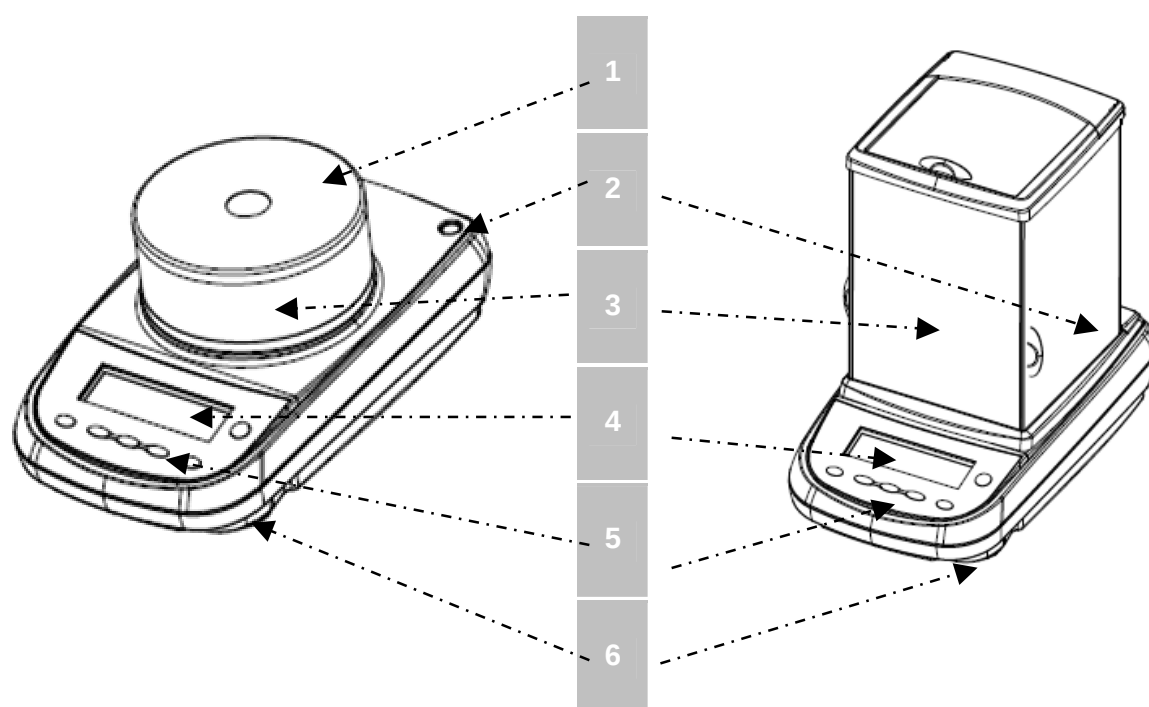


Fig.1

Balance with circular and three-panes windshield

Name and function	
1	Cap of windshield
2	Levelling bubble
3	Windshield
4	LCD display
5	Function buttons
6	Levelling feet

5.2 Rear and bottom part of the balance

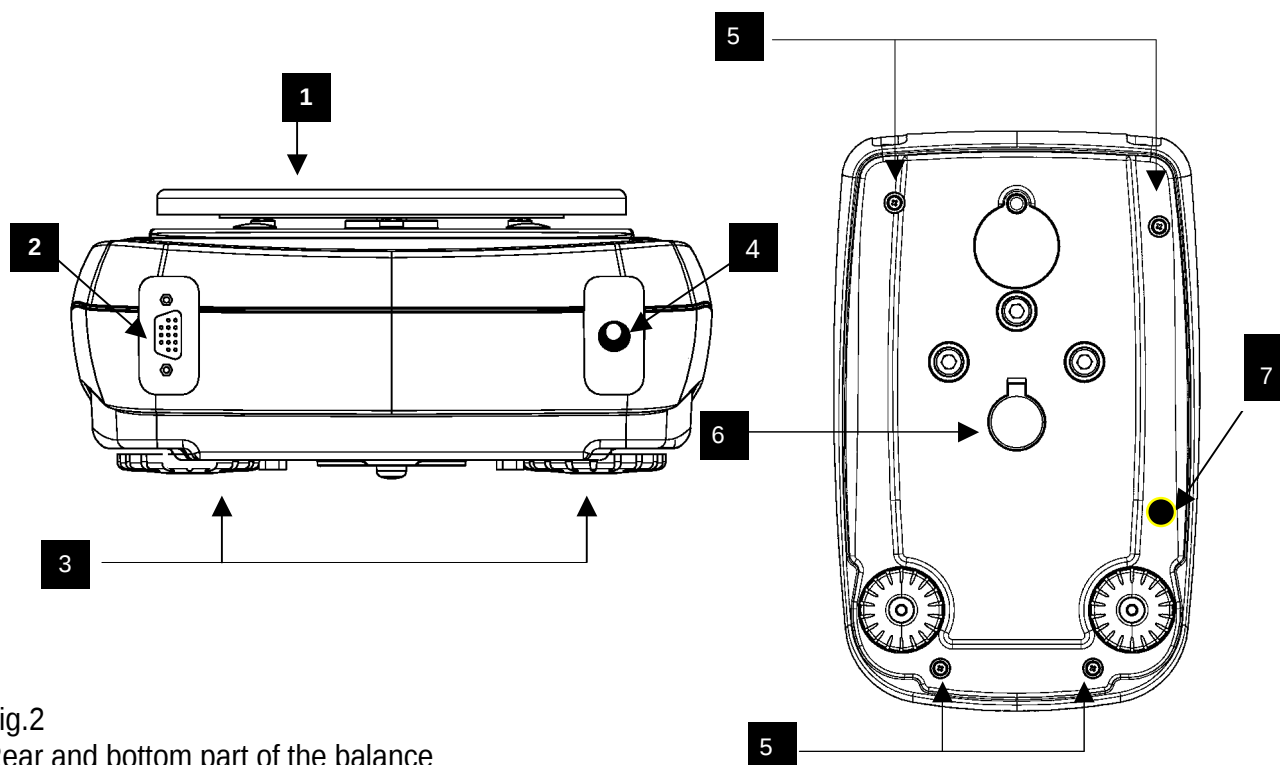


Fig.2
Rear and bottom part of the balance

Name and function (Fig.2)	
1	Weighing pan
2	9 poles (pin) female connector for RS232 interface for printer/PC/external keyboard
3	Levelling feet
4	connector for power supply
5	Screws for balance closing. In models with 4 feet it is necessary first to remove the two rear feet (one fixed and one moving) to access to the closing screws
6	Hook for under-balance weighing. Hook not available in internal calibration models and for verified balances. Note: please remove the cap to access the hook
7	Protection screw for autocalibration motor (<u>only for models with internal calibration</u>). ATTENTION: remove this screw by hand as soon as the balance is removed from its package

5.3 Rear and bottom part of the balance, models Version Plane Square

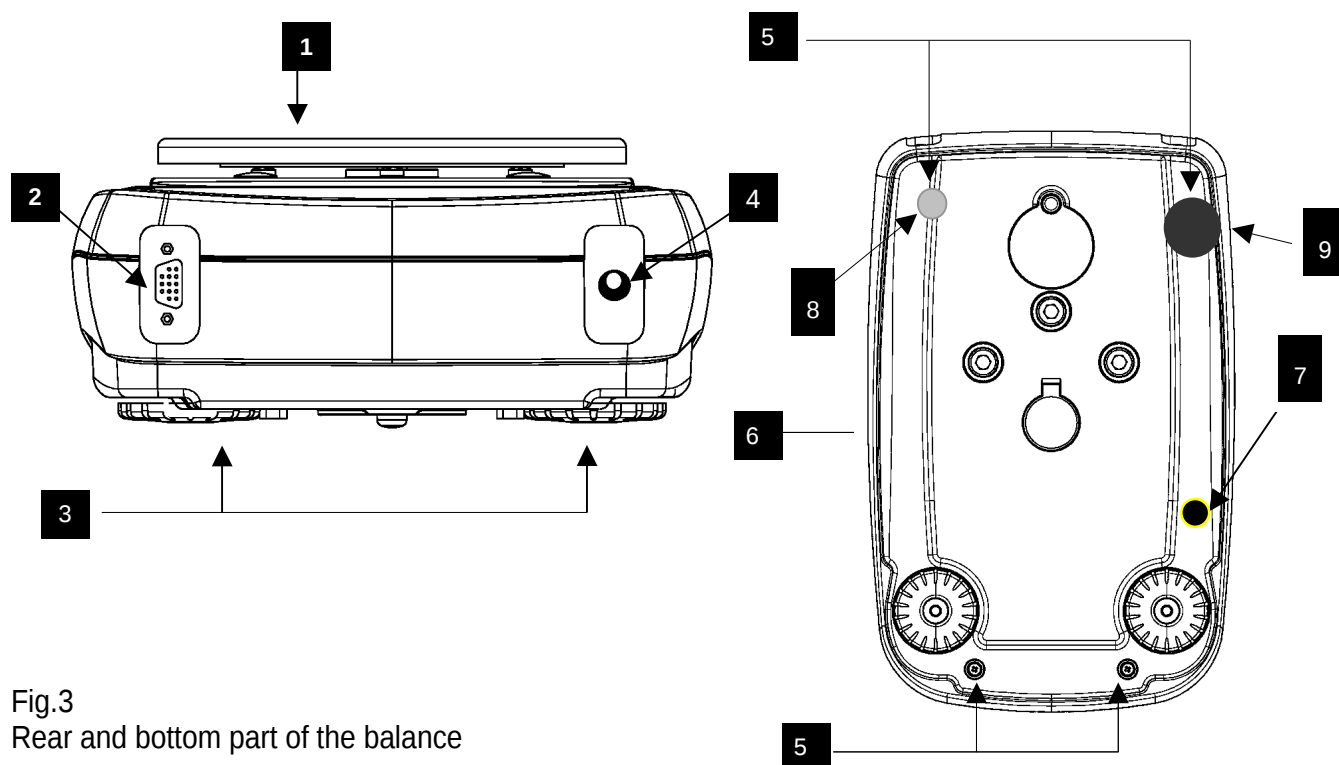


Fig.3
Rear and bottom part of the balance

Name and function (Fig.3)	
1	Weighing pan
2	9 poles (pin) female connector for RS232 interface for printer/PC/external keyboard
3	Levelling feet
4	connector for power supply
5	Screws for balance closing. In models with 4 feet it is necessary first to remove the two rear feet (one fixed and one moving) to access to the closing screws
6	Hook for under-balance weighing. Hook not available in internal calibration models and for verified balances. Note: please remove the cap to access the hook
7	Protection screw for autocalibration motor (<u>only for models with internal calibration</u>). ATTENTION: remove this screw by hand as soon as the balance is removed from its package
8	Adjustable rear
9	Fixed rear feet

6 Weighing



For verified Balances models at switch on it will appear “WARM UP” on display and a 30 minutes countdown will start

After having connected the balance to the power outlet, a self-diagnosis of the electronic circuits is automatically carried out that ends with the standby indication.



6.1 Stand By


From the “**STAND BY**” state:

- Press the **ON/OFF** button to bring the balance to working conditions.
- Press the **ON/OFF** button again to bring it back to the “**STAND BY**” state.



6.2 Simple weighing

Place the sample to be weighed on the plate and read the value of the weight on the display as soon as the asterisk stability symbol ✱ appears.

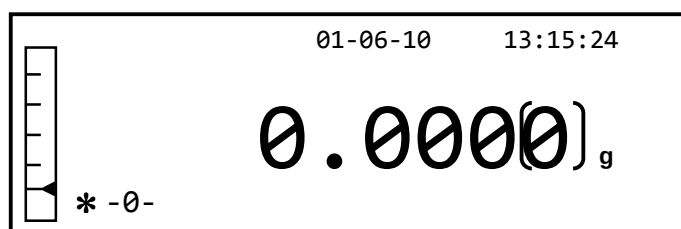

The digit that represents the verification division “e” is rounded with a frame

7 Calibration

The electronic balance carries out mass measurements using gravity (g). Differences in geographical regions and altitudes vary the gravitational acceleration (g).

The balance must therefore be adapted to environmental conditions to obtain accurate measurements. This regulation is carried out through the calibration function.

7.1 Balances with external calibration


7.1.1 External calibration




For verified Balances models the External calibration function is locked

The calibration is carried out through the CAL button.


1. Press the **CAL** button with the plate unloaded; the word CAL will be displayed.

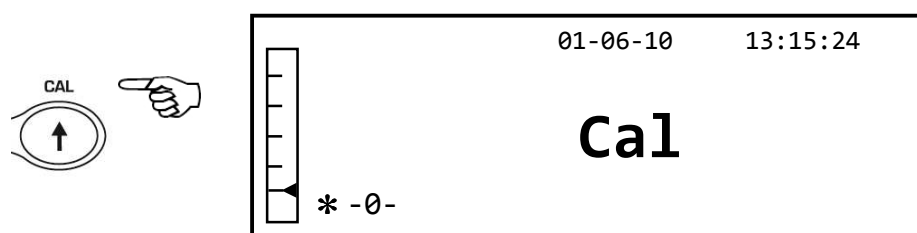

2. When the value of the calibration weight begins to flash, load the weight indicated by the display on the plate.

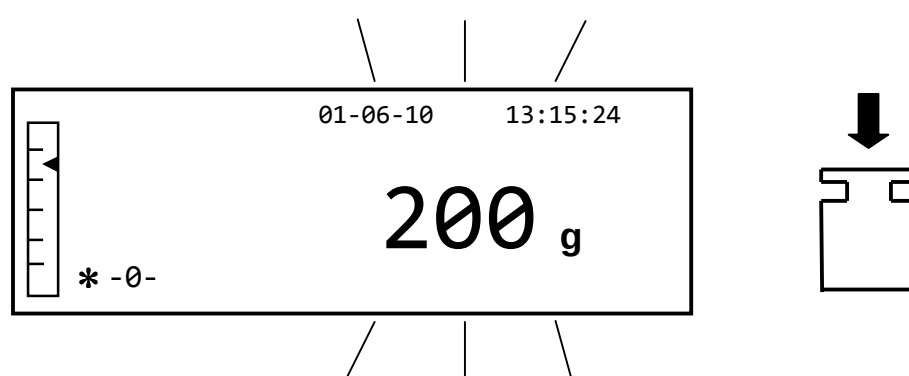


3. The display will stop flashing, indicating the value of the calibration weight with the stability indicator on.
Once the calibration has been carried out, the calibrated weight will be displayed with the indication of the current unit of measurement.
4. Remove the calibration weight.
The balance is ready for weighing operations.

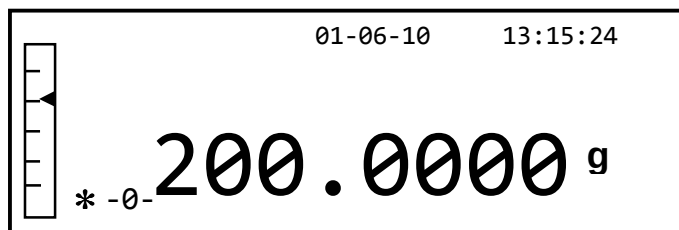

NOTE: an error message will be displayed if there is some interference during the calibration process. To interrupt the calibration process, press the ON/OFF button while the calibration weight indication flashes.

It is also possible to calibrate the balance with a calibration weight greater than the pre-set calibration weight:

1. Press and hold the **CAL** button with the plate empty until the beeping stops, and then release the button. The word “-CAL-” will be shown on the display, followed by the word “**LOAD**”, flashing.


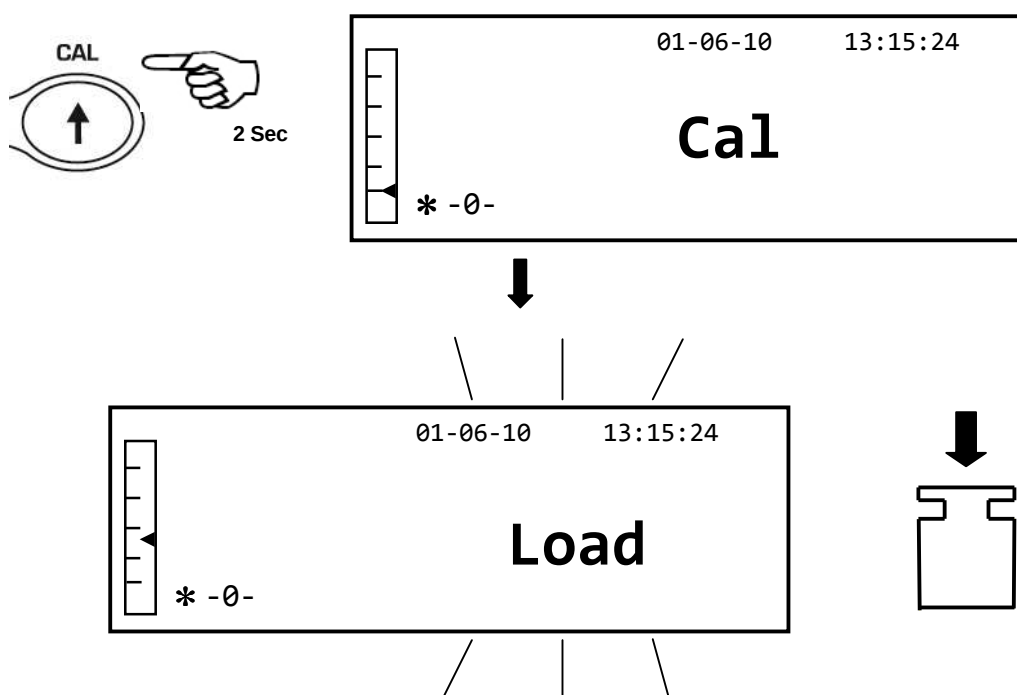


2. Load a weight that is equal to or greater than the pre-set calibration weight on the plate; the balance will recognize a weight that is equal to or greater than the calibration weight as valid provided that it is a whole weight with respect to the most significant figure of the calibration weight.

E.G.: if the calibration weight is 200 g, it will be possible to calibrate the balance with values that go from 200g, 300g, 400g up to the upper capacity limit of the balance.

The word “**LOAD**” on the display will stop flashing; once the calibration has been carried out, the value of the weight used will be displayed.

3. Remove the calibration weight; the balance is ready for the weighing operations.

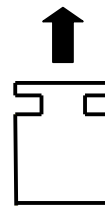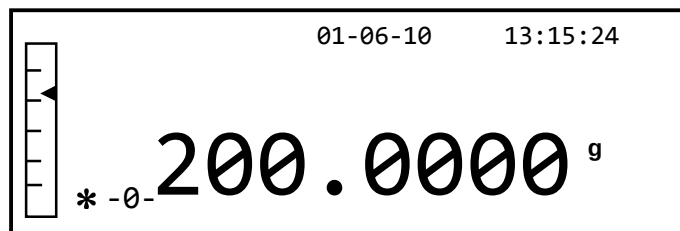

NOTE: an error message will be displayed if there is some interference during the calibration process. To interrupt the calibration process, press the ON/OFF button while the calibration weight indication flashes.

7.2 Balances with internal calibration



7.2.1 Internal calibration




This is Factory setting for verified balances models

In these models the calibration is carried out through an internal automatic system:

1. Press the **CAL** button with the plate empty.
The display will show the message "**CAL**" and the balance's calibration will be carried out automatically.

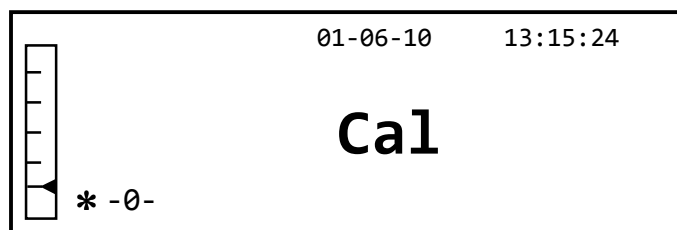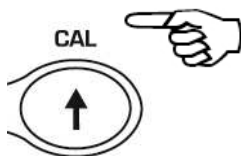

2. At the end of the calibration, the balance will return to normal weighing conditions.

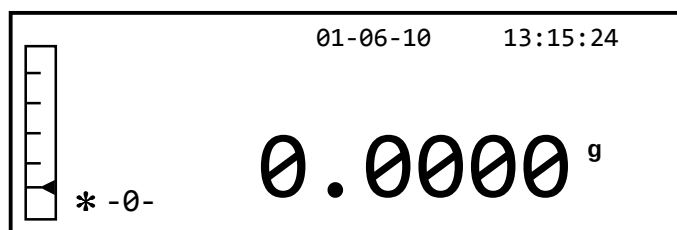

If the calibration is not completed due to vibrations or drafts, the message "**CAL BUT**" will be displayed. Press the CAL button again, and if the problem persists, select external calibration and contact the supplier.

To modify the calibration mode in these models with internal calibration, see section 9.12.1

7.3 Verification

General introduction:

According to EU directive 90/384/EEC balances must be verified if they are used for legally controlled area as follows:

- a) For commercial transactions if the price of goods is determined by weighing.
- b) For the production of medicines in pharmacies as well as for analyses in the medical and pharmaceutical laboratory.
- c) For official purposes.
- d) For manufacturing final packages.

In cases of doubt, please contact your local trade in standard.


Verification instructions

An EU type approval exists for balances described in their technical data as verifiable. If a balance is used where obligation to verify exists as described above, it must be verified and re-verified in regular intervals.

Re-verification of a balance is carried out according to the respective national regulations. The legal regulation of the country where the balance is used must be observed!

After verification the balance is sealed at the indicated positions.

Verification of the balance is invalid without the "seal".

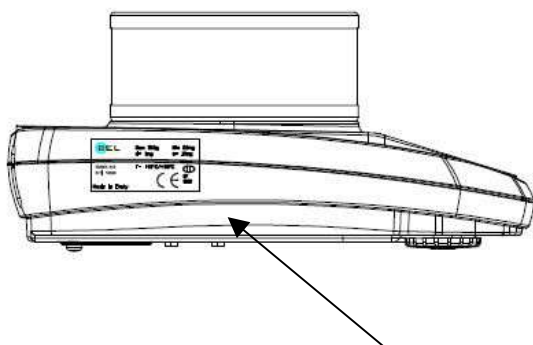

Position of the seal mark for Verified balances models



Balances with obligation to verify must be taken out of operation if:

- **The weighing result of the balance is outside the error limit.**
Therefore, in regular intervals load balance with known test weight (ca. 1/3 of the max. load) and compare with displayed value.
- **The reverification deadline has been exceeded.**





8 Tare function


1. The relative weight will be shown on the display.

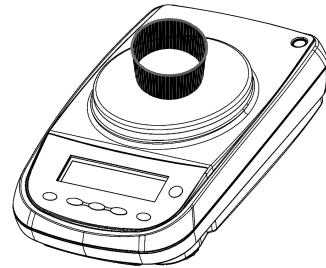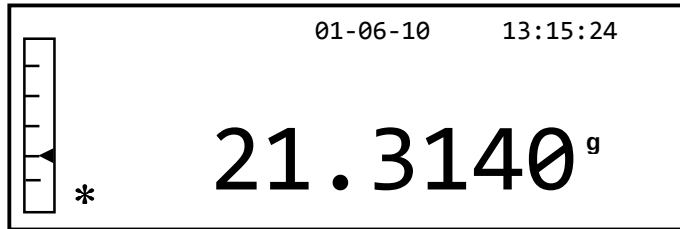



2. Press the **O/T** button. The word “Tare” will be displayed.

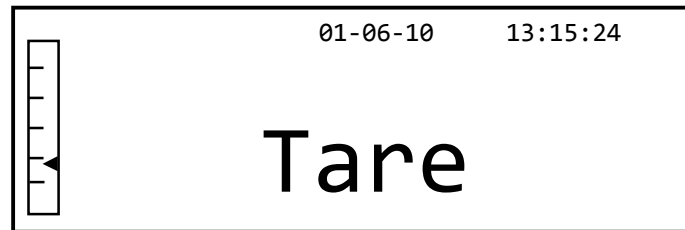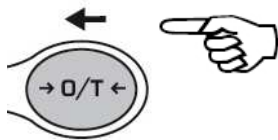

3. Once stability has been achieved, the zero value “**0.000**” will be displayed. In case stability is not achieved due to drafts, vibrations, or other types of disturbance, the dashes will continue to be displayed.

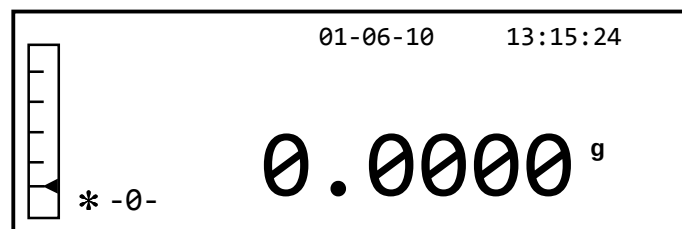

4. Put the objects to be weighed in the container. Read the value of the net weight on the display.

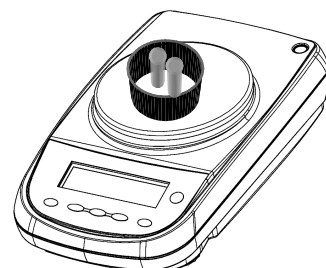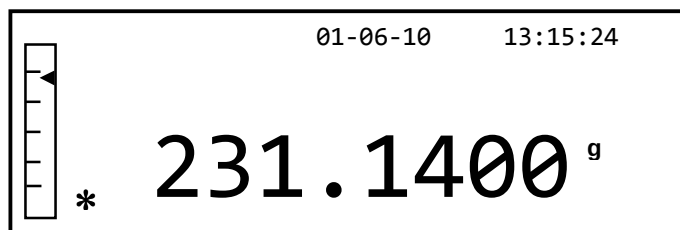


8.1 Manual tare function




Function not available in verified balances

This function allows a tare value to be entered manually.

1. Press and hold the **O/T** button with the plate empty until the beeping stops, and then release the button.
2. The following word will be shown on the display:

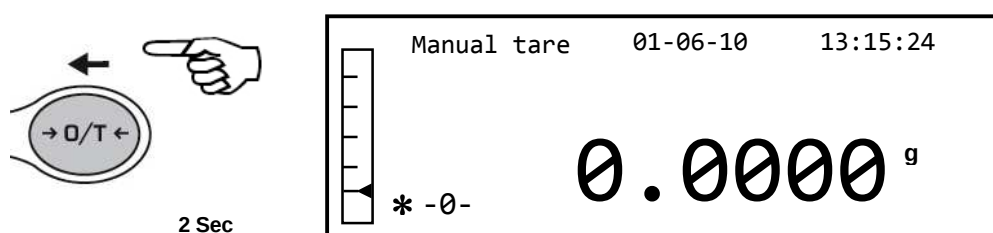

3. Now insert the desired tare value using the **CAL** and **MENU** buttons to increase and decrease the number, while pressing the **O/T** button to pass to the next number. During the entering phase, holding down the **O/T** button allows you to delete the inserted value.

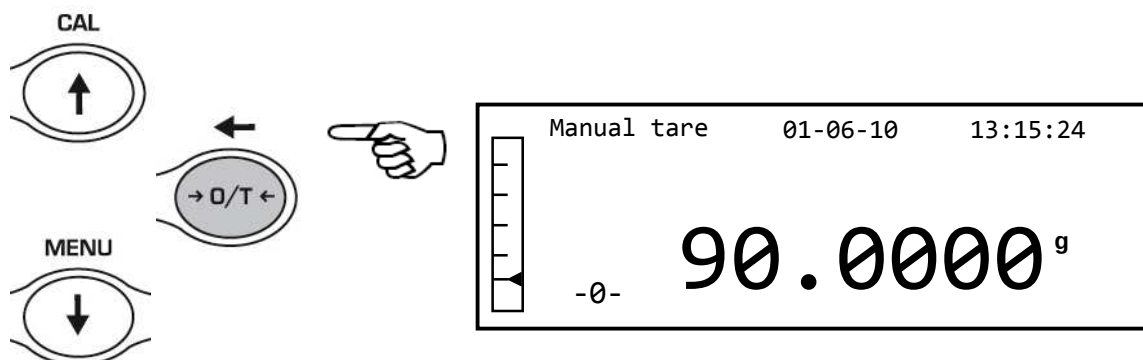


4. After having entered the desired value, press the **PRINT** button to confirm it. The value will remain in memory until the **TARE** button is pressed, or the instrument is disconnected from the power supply.

8.2 Optional external alphanumeric keyboard (Code T201)




Function locked in verified balances

If you have an optional alphanumeric keyboard, the tare can also be carried out by pressing the **TARE/DELETE** button of the keyboard itself, in the same way as described above.

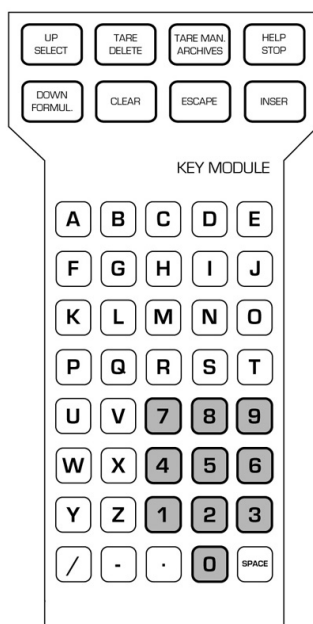


It is also possible to manually enter a known tare value from the keypad.

1. Press the **TARE MAN** button on the alphanumeric keypad. An arrow and the preceding manual tare value, if entered, will be shown on the display.

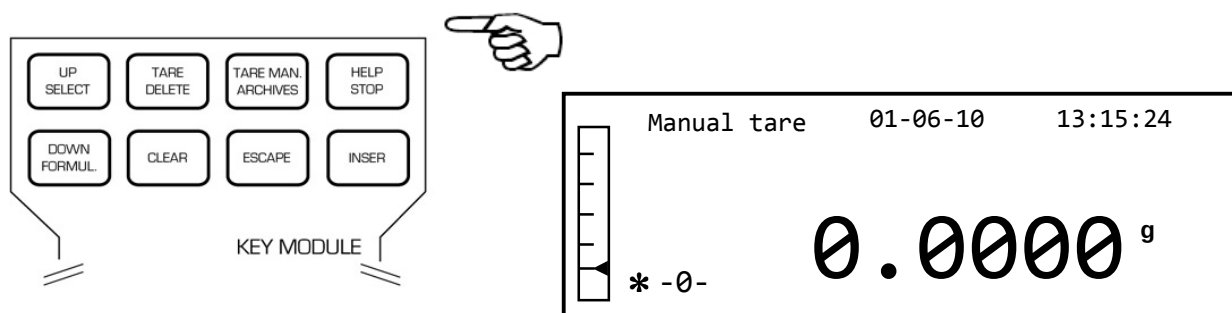

2. Press the **CLEAR** button to zero the preceding value and enter the new value, using the numerical keys on the lower part of the alphanumeric keypad.


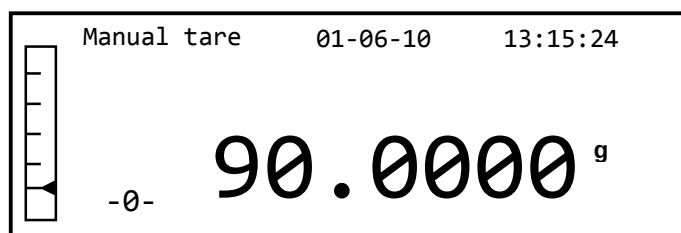


3. Press **INSER** to confirm.
4. Press the **ESCAPE** key to exit from the tare condition.


To zero the value of the manual tare inserted by keyboard, carry out a normal tare operation with the **O/T** button on the balance, or with the **TARE/DELETE** button on the optional alphanumeric keyboard.

9 Balance parameters setup menu

1. Press and hold the **MENU** button with the plate empty until the beeping stops, and then release the button.
2. The following writing will be shown on the display:


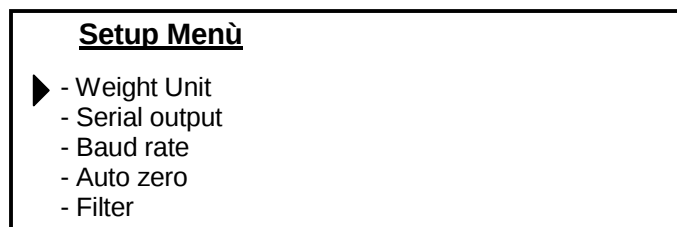
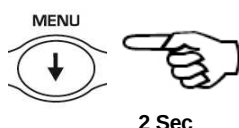


3. Now use the CAL and MENU buttons to navigate forward or backward in the parameters menu.

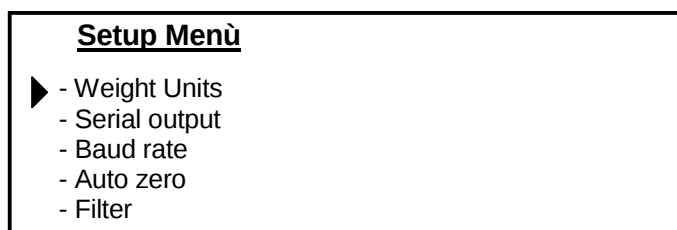

...


...

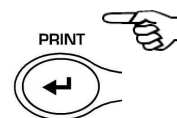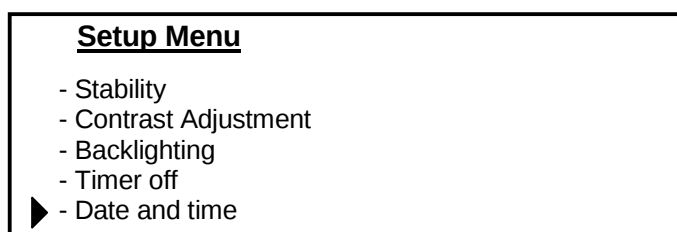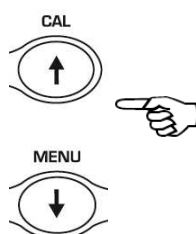

...

...

...

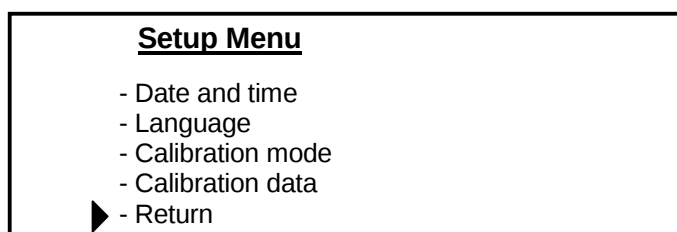


4. Position the cursor on the desired parameter and press the PRINT button to confirm the selection.
5. Press the ON/OFF button to exit from the menu or select the return function and press the PRINT button.


9.1 Weight Units.

The weight units with which the balance will weigh can be selected.

1. Set the weight units parameter as described in chapter 10.
The currently set weight units will be shown on the display:

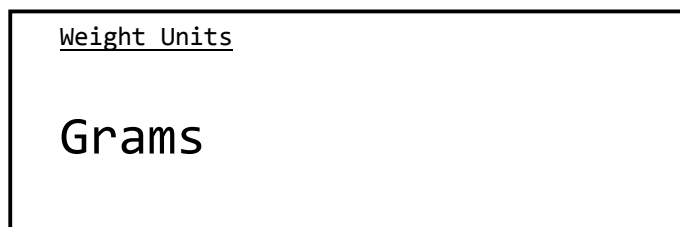

2. Now by pressing the **MENU** or **CAL** button it will be possible to scroll through the weight units forward or backward.
3. Then press the **PRINT** button to confirm the desired weight units.

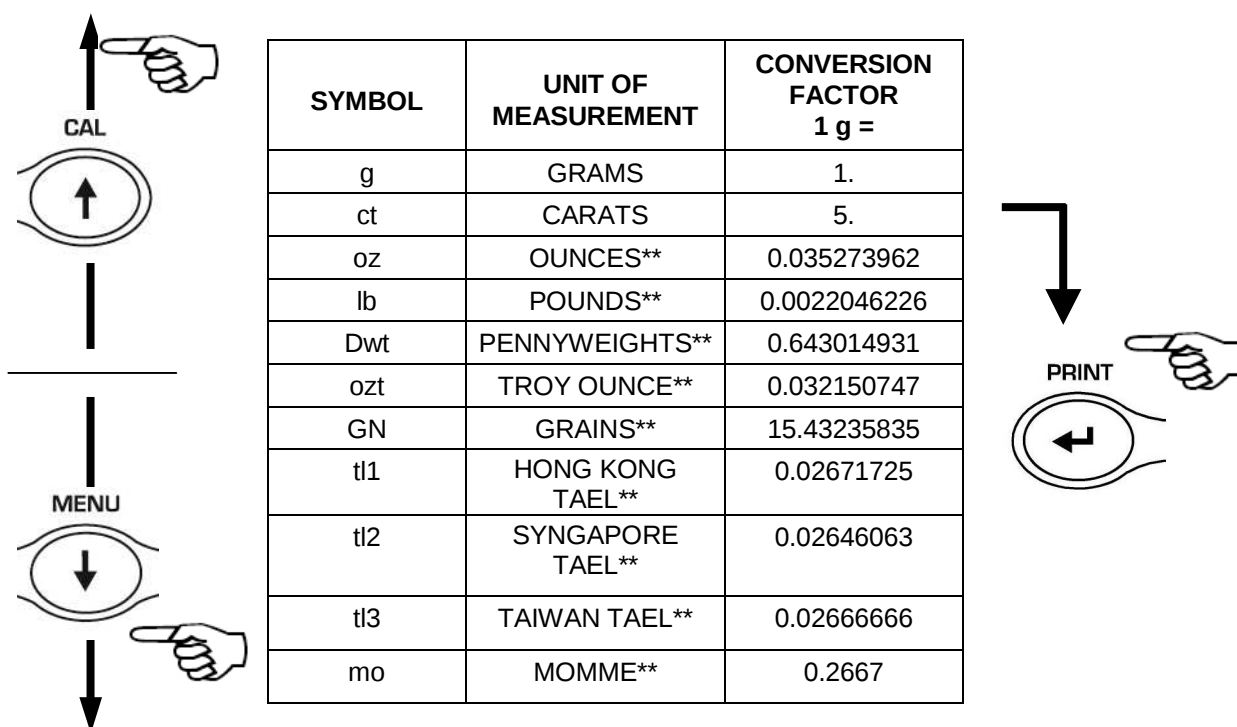

****not available in verified balances**


4. After having selected the desired weight unit, the screen relative to the menu of balance parameters will be displayed again. It will now be possible to select another parameter or return to weighing mode by pressing the **ON/OFF** button.

By returning to the weighing mode in the right part of the display, the symbol of the selected weight unit will be displayed.

9.2 Serial output setup

Different data transmission devices and modes can be selected.

1. Select the serial output parameter as described in paragraph 10.
The currently set transmission mode will be shown on the display:

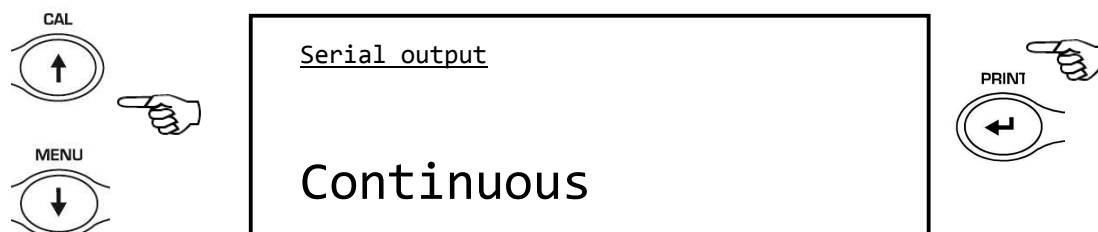

2. Now by pressing the **MENU** or **CAL** button it will be possible to scroll through the serial output **MENU** forward or backward.
3. Then press the **PRINT** button to confirm the desired transmission mode.

The different transmission modes are illustrated below:

TRANSMISSION MODE	FEATURES
Continuous	Transmits the weight data in a continuous way
On demand	Transmits the weight data only when the PRINT button is pressed
Generic printer	The weight data is printed only when the Busy command is active
Tlp50 printer	The weight data is printed only if the Tlp50 model printer is connected
Upon request - Glp	Transmits the weight data and the Glp information only when the PRINT button is pressed
Generic printer - Glp	The weight data and the Glp information are printed only when the Busy command is active
Tlp – Glp printer	The weight data and the Glp information are printed only if the Tlp50 model printer is connected

NOTE: transmission speed selection (paragraph 9.3)

4. After having selected the desired transmission mode, the screen relative to the balance parameters menu will be displayed again. It will now be possible to select another parameter or return to weighing mode by pressing the **ON/OFF** button.

9.3 Transmission speed selection

Different data transmission speeds can be selected.

1. Select the baud rate parameter as described in paragraph 10.
The currently set transmission speed will be shown on the display:




2. Select the serial data transmission speed (1200-2400-4800-9600 baud). By pressing the **MENU** or **CAL** button you can scroll the different transmission speeds forward or backward; then confirm the choice with the **PRINT** button.
3. After having selected the desired transmission speed, the screen relative to the balance parameter menu will be displayed again. It will now be possible to select another parameter or return to weighing mode by pressing the **ON/OFF** button.

9.4 Autozero function

Autozero is a correction from a possible deviation from zero.

Different autozero levels can be selected.

1. Select the autozero parameter as described in chapter 10.
The currently set autozero parameter will be shown on the display:

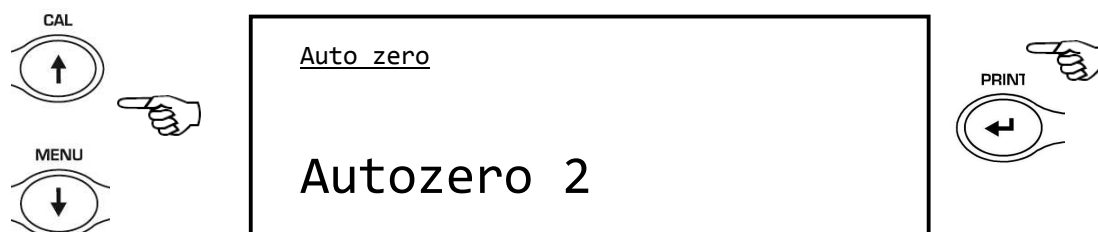

2. Select the desired autozero level. By pressing the **MENU** or **CAL** button it will be possible to scroll through the various levels forward or backward; then confirm your choice with the **PRINT** button.

AUTOZERO MENU	AUTOZERO LEVEL
Autozero off	Autozero off
Autozero 1	Light autozero
Autozero 2*	Average autozero
Autozero 3*	Heavy autozero
Autozero 3E*	Heavy full-scale autozero


* locked in verified balances

3. After having selected the desired autozero, the screen relative to the balance parameters menu will be displayed again. It will now be possible to select another parameter or return to weighing mode by pressing the **ON/OFF** button.

9.5 Filter selection

The balance can be adapted to different environmental conditions thanks to the selection of three different filters:

1. Select the filter parameter as described in paragraph 10.
The currently set filter type will be shown on the display:

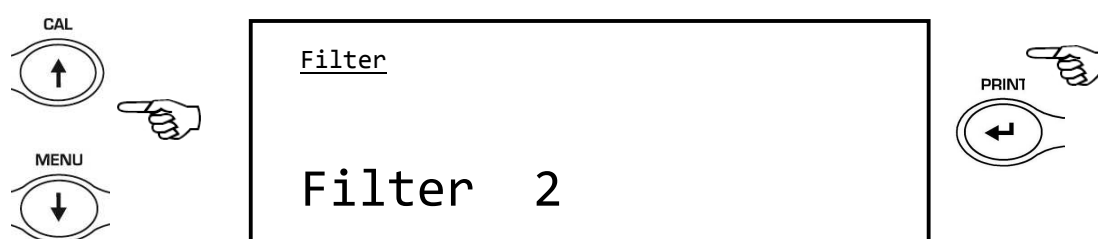

2. Select the desired filter level. Pressing the **MENU** or **CAL** button it will be possible to scroll through the various levels forward or backward; then confirm your choice with the **PRINT** button.


FILTER MENU	FILTER LEVEL
Filter 1	Use this filter level in stable environmental conditions and for use of the instrument in filling or dosing mode
Filter 2	Use this filter level when the environmental conditions are not stable
Filter 3	Use this filter level when the environmental conditions are particularly unstable

3. After having select the desired filter level, the screen relative to the balance parameters menu will be displayed again. It will now be possible to select another parameter or return to weighing mode by pressing the **ON/OFF** button.

9.6 Stability function

The stability symbol will appear on the display when the weight is stable within a defined interval.

1. Select the stability parameter as described in paragraph 10.
The currently set type of stability will be shown on the display:

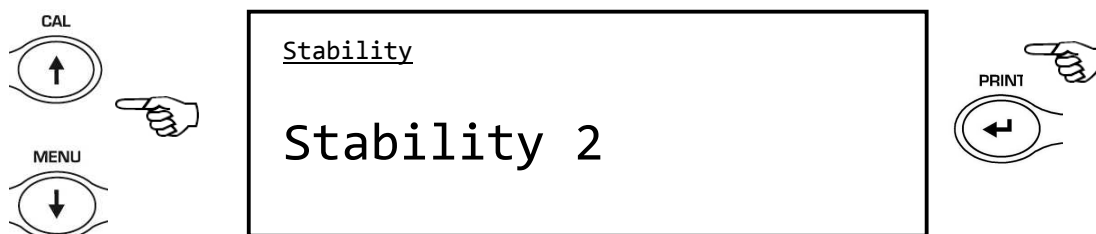

2. Select the desired stability level. By pressing the **MENU** or **CAL** button you can scroll through the various levels forward or backward; then confirm your choice with the **PRINT** button.

STABILITY MENU	LEVEL OF STABILITY
Stability 1	Use this level of stability when the environmental conditions are stable
Stability 2	Use this level of stability when the environmental conditions are less stable
Stability 3	Use this level of stability when the environmental conditions are unstable

3. After having selected the desired level of stability, the screen relative to the balance parameter menu will be displayed again. It will now be possible to select another parameter or return to weighing mode by pressing the **ON/OFF** button.


9.7 Contrast Adjustment

The balance's display is equipped with an LCD display; the contrast can be regulated in order to make the indication as visible as possible from different angles.

1. Select the contrast adjustment parameter as described in paragraph 10.
The currently set contrast value will be shown on the display:



2. Select the desired contrast value. By pressing the **MENU** or **CAL** button it will be possible to increase or decrease the value; then confirm the choice with the **PRINT** button.

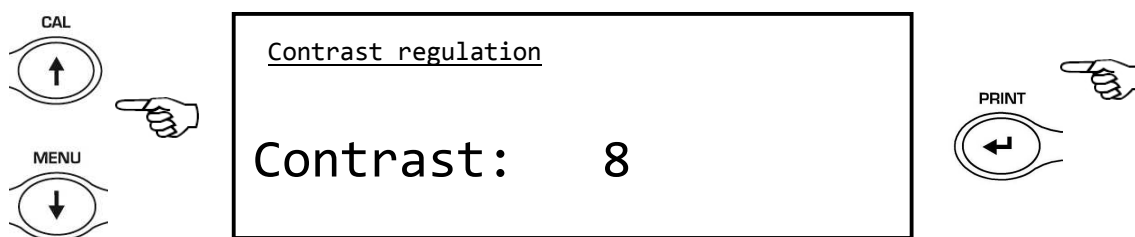


3. After having selected the desired contrast level, the screen relative to the balance parameters menu will be displayed again. It will now be possible to select another parameter or return to weighing mode by pressing the **ON/OFF** button.

9.8 Backlight regulation

The balance's display is equipped with a backlight to make the indication visible even in low light conditions.

1. Select the backlight parameter as described in paragraph 10.
The currently set mode will be shown on the display:

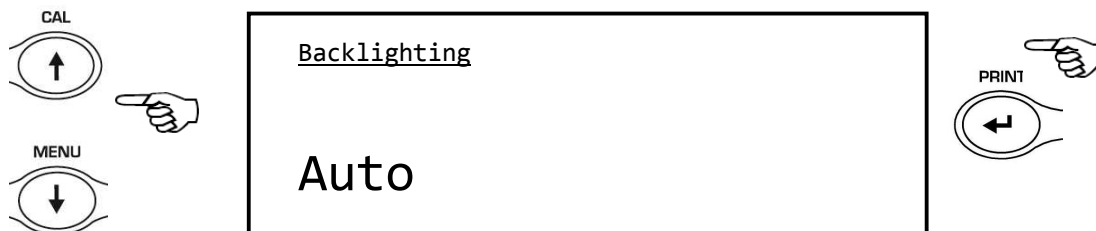

2. Select the desired modality. By pressing the **MENU** or **CAL** button it will be possible to scroll through the various levels forward or backward; then confirm the choice with the **PRINT** button.


BACKLIGHT MENU	BACKLIGHT MODE
Auto	Backlight automatically active during the weighing phases
On	Backlight always on
Off	Backlight always off

3. After having selected the desired mode, the screen relative to the balance parameters menu will be displayed again. It will now be possible to select another parameter or return to weighing mode by pressing the **ON/OFF** button.

9.9 Timer-off function

This function allows you to activate the automatic turn-off of the balance after a preset time of inactivity.

1. Select the Timer off parameter as described in paragraph 10.
The currently set mode will be shown on the display:

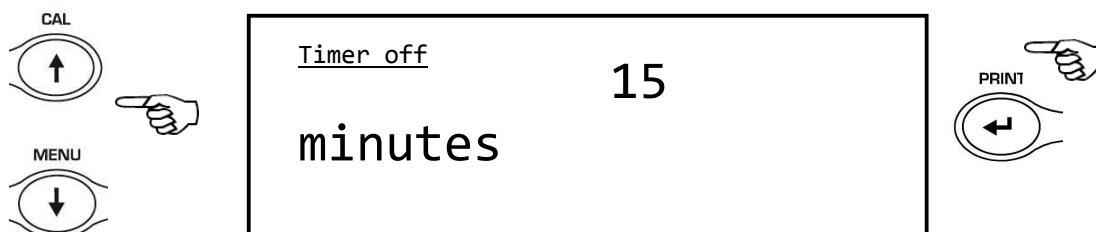

2. Select the desired auto-off modality. By pressing the **MENU** or **CAL** button it will be possible to scroll the various levels forward or backward; then confirm the choice with the **PRINT** button.

TIMER-OFF MENU	AUTO OFF MODE
Disabled	Timer-off disabled
2 minutes	Timer-off after 2 minutes of inactivity
5 minutes	Timer-off after 5 minutes of inactivity
15 minutes	Timer-off after 15 minutes of inactivity


NOTE: In balances equipped with an internal battery, it is recommended not to disable the Timer-off mode in order to prolong battery life.

3. After having selected the desired mode, the screen relative to the balance parameters menu will be displayed again. It will now be possible to select another parameter or return to weighing mode by pressing the **ON/OFF** button.


9.10 Date and time regulation

This function allows you to regulate the date and time, and to modify the date display format.


1. Select the date and time parameter as described in paragraph 10.
The current set date and time will be shown on the display:

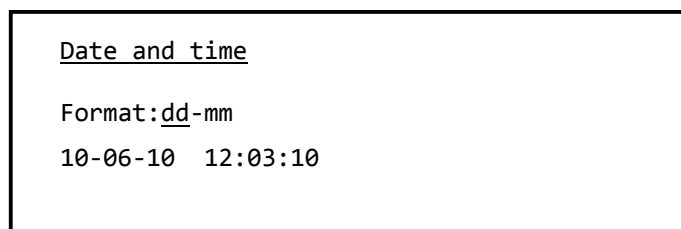

2. Select the desired format of the date. Pressing the **MENU** or **CAL** button it will be possible to modify the format dd-mm or mm-dd; then confirm the choice with the **PRINT** button.

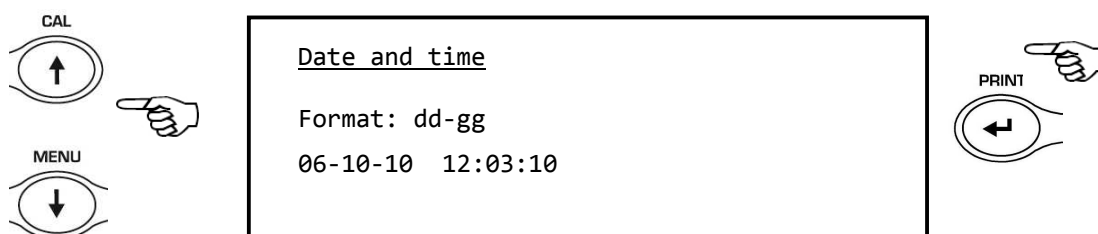

3. Set the desired date and time by using the **MENU** and **CAL** buttons to increase and decrease the number and the **PRINT** button to pass to the next date.

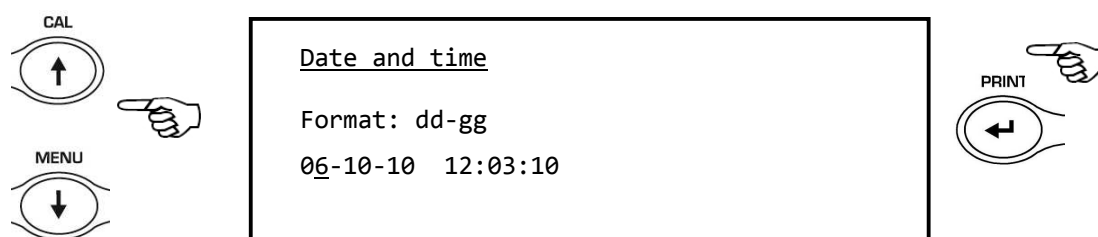

4. After having regulated the date and time, press and hold the **PRINT** button until the beeping stops and then release the button to save the settings.
5. The screen relative to the balance parameters menu will then be displayed. It will now be possible to select another parameter or return to weighing mode by pressing the **ON/OFF** button.

9.11 Language selection

This function allows you to set the desired usage language.

1. Select the language parameter as described in paragraph 10.
The currently set language will be shown on the display:



2. Select the desired language. By pressing the **MENU** or **CAL** button it will be possible to scroll the various levels forward or backward; then confirm your choice with the **PRINT** button.

LANGUAGE MENU	LANGUAGE
Italian	Italian language
English	English language
Português	Portuguese language
Deutsch	German language
Français	French language
Español	Spanish language


- After having selected the desired language, the screen relative to the balance parameters menu will be displayed again. It will now be possible to select another parameter or return to weighing mode by pressing the **ON/OFF** button.

9.12 Calibration mode setting

This function allows you to set the calibration mode.

* This function is only available for the models equipped with internal calibration (not available for all models). To check if the instrument is equipped with internal calibration mode, refer to the table of technical features, chapter 18.

1. Select the calibration mode parameter as described in paragraph 10.
The currently set calibration mode will be shown on the display:

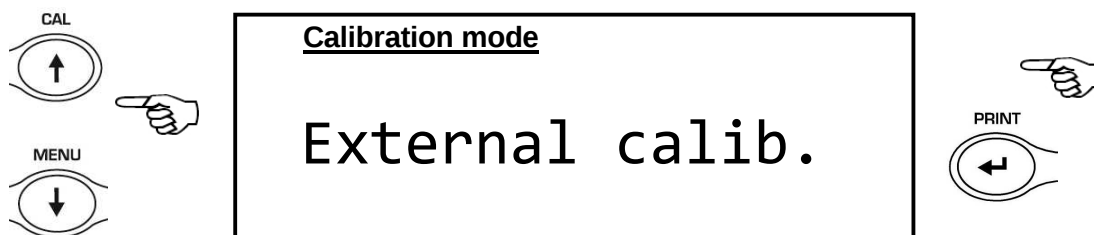

2. Select the desired mode. By pressing the **MENU** or **CAL** buttons it will be possible to scroll the different calibration modes forward or backward:
 - External calibration*
 - Internal calibration*
 - Automatic calibration
 - Technical calibration****Locked for verified balances models**
3. Press the **PRINT** button to confirm “**AUT-CAL**”, “**I-CAL**”, “**E-CAL**”.
To confirm “**TEC-CAL**”, keep the **PRINT** button pressed until the beeping stops.
4. The screen relative to the balance parameters menu will then be displayed. It will now be possible to select another parameter or return to weighing mode by pressing the **ON/OFF** button.


9.12.1 Automatic Calibration (AUT-CAL)



This is Factory setting for verified balances models

The balance self-calibrates when the temperature variation exceeds the factory preset value and at factory preset time intervals, through the internal reference mass, and only if the balance pan is empty.

When the balance needs to perform the Automatic calibration, the display will show the following message:

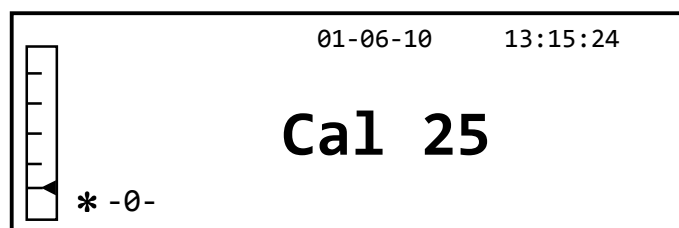

A 25-seconds countdown will start during which you can decide if:

- Stop the automatic-calibration procedure by pressing the “**ON/OFF**” button that will be delayed of 5 minutes


or

- Let the countdown finish so that the automatic calibration starts


Note: during the countdown do NOT load anything on the plate!

In this mode, it is also possible to carry out the calibration with the internal reference mass by pressing the **CAL** button at any moment, first ensuring that no weight is loaded on the plate.

1. Press the **CAL** button with the plate empty. The display will show the message “**CAL**” and the balance’s calibration will be carried out automatically.

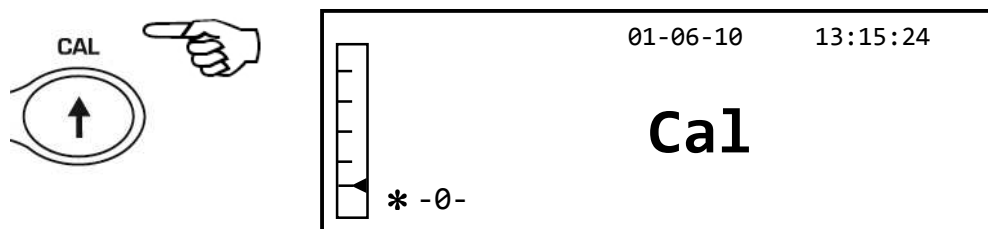

2. At the end of the calibration, the balance will return to normal weighing conditions.

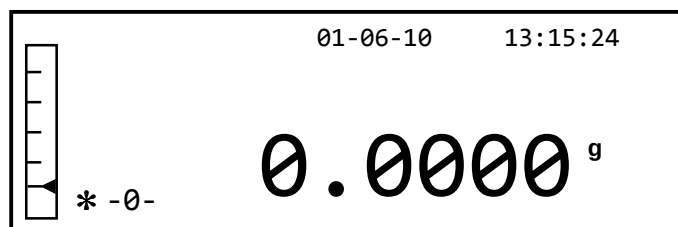

If the calibration is not completed due to vibrations or drafts, the message “**CAL but**” will be displayed. Press the **CAL** button again, and if the problem persists, select external calibration and contact the supplier.

9.12.2 Internal calibration (I-CAL) [locked for verified balances]

The balance calibrates itself through the internal reference mass **ONLY** upon the request of the user by pressing the **CAL** button.

Before carrying out the internal calibration, ensure that no weight is loaded on the plate.

9.12.3 External calibration (E-CAL) [locked for verified balances]


The balance will be calibrated by using the external reference mass.

(Follow the procedures described in paragraph 8.1.2)


9.12.4 Technical calibration (TEC-CAL) [locked for verified balances]

This function allows the internal reference mass to be calibrated whenever assistance-control-maintenance interventions make this necessary.


1. After having selected the **TEC-CAL** calibration mode, press the **CAL** button with the plate empty. The word “CAL” will be displayed.

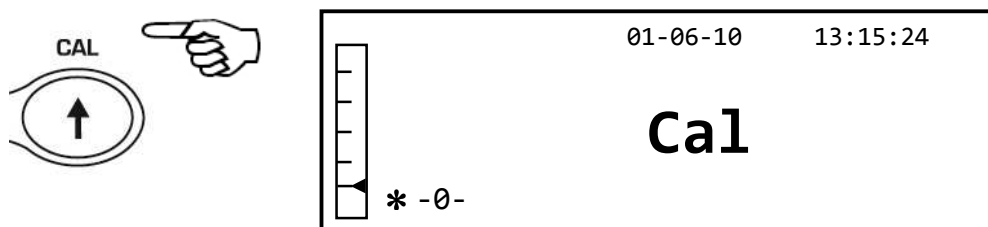

2. When the value of the calibration weight begins to flash, load the calibration weight on the plate.

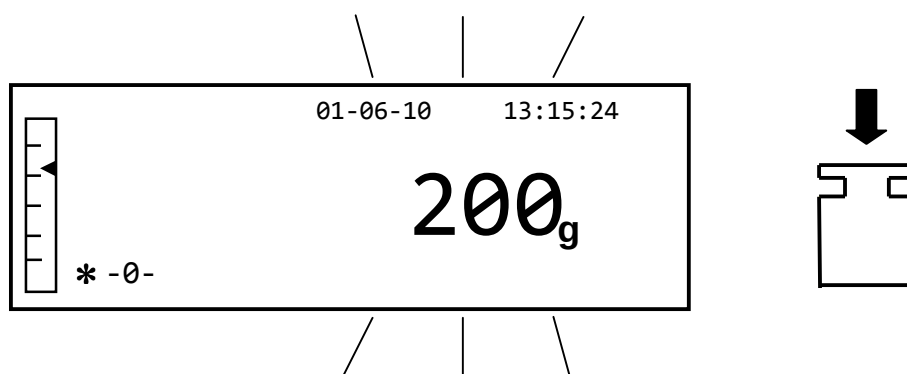

3. Wait for the calibrated weight to be displayed and the stability symbol to turn on, and then remove the weight from the plate.
4. When “0.000” is shown on the display, press the **PRINT** button in a prolonged manner until the beeping stops. The acquisition and automatic storage of the internal weight will now begin. During the acquisition cycle, the display will show the following flashing writing:

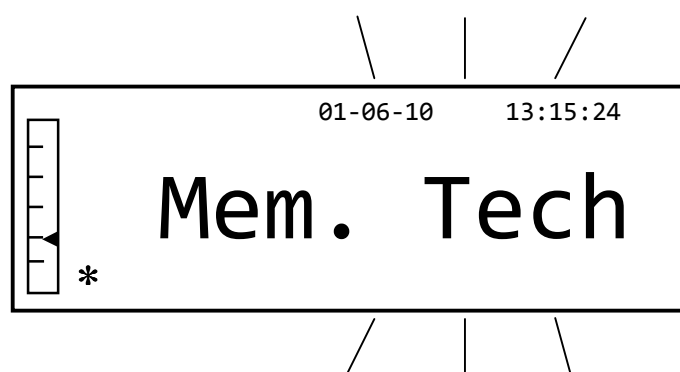

5. Once the internal calibration has been stored, the balance will return to the normal weighing condition.
6. Now re-enter the calibration menu as described in paragraph 9.1.12 and set the desired internal, automatic, or external calibration modality.



WARNING: this procedure must be carried out only using class E2 reference masses.

9.13 Calibration data

This function allows you to display the data relative to the last calibration carried out.

- Date
- Calibration mode
- Correction



1. Select the calibration data parameter as described in paragraph 10.
The data relative to the last calibration carried out will be shown on the display:

Calibration data	
10-06-10	
External calib.	2000.0g
Corr.:	1.2g




2. Press the **PRINT** button to print the calibration data.
3. Press the **ON/OFF** button to exit from the screen and return to the balance parameters menu. It will now be possible to select another parameter or return to weighing mode by pressing the **ON/OFF** button.

10 Balance programs menu

1. Press the **MENU** button with the plate empty.
2. The following writing will be shown on the display:

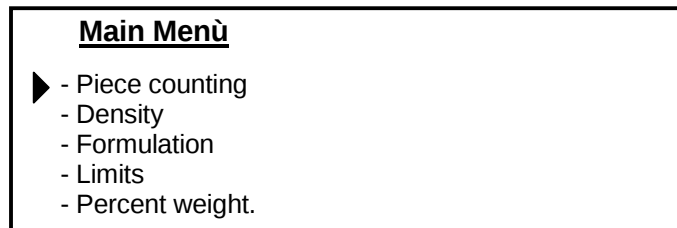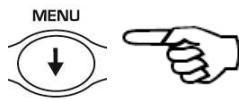

3. Now use the CAL and MENU buttons to navigate forward or backward in the menu of parameters.

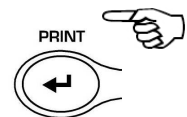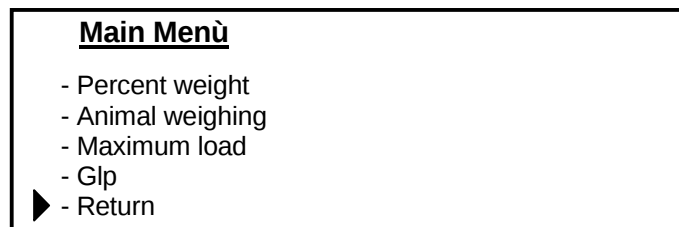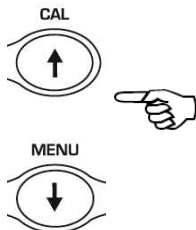


4. Position the cursor on the desired parameter and press the PRINT button to confirm the selection.
5. Press the ON/OFF button to exit from the menu or select the return function and press the PRINT button.


10.1 Piece counting function

The piece counting program allows you to carry out a total count of the pieces after having carried out a sampling of pieces or having inserted the average unit weight of the piece.


1. Select the piece counting program as described in paragraph 10.
The following screen will be shown on the display:

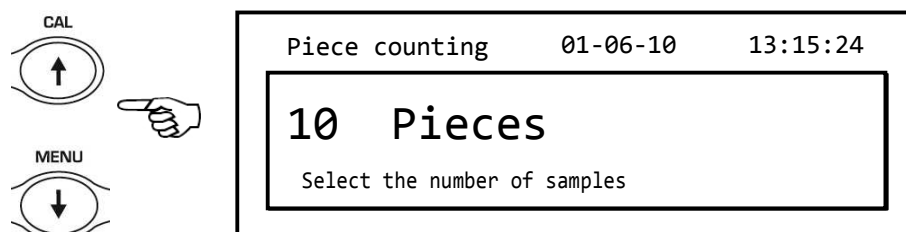


1. Select the number of pieces to put on the plate as a sample, pressing in sequence the **MENU** button to increase and the **CAL** button to decrease.
2. Load an empty container, if used, then press the **PRINT** button to confirm.
The choice of the number of pieces (10, 25, 50, 100, manual, see chapter 10.3) is a function of the weight of an individual piece. Load the empty container, if used.

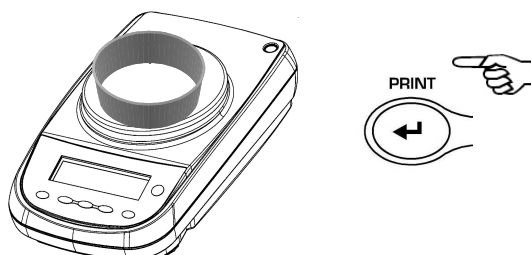

3. Load the number of pieces indicated on the display on the plate and press the **PRINT** button.

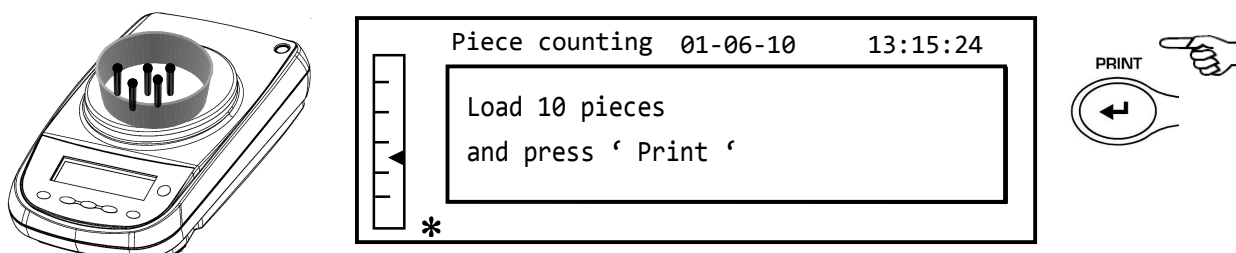




If there are enough samples (e.g. 10 as in the figure), the number of pieces loaded will appear on the display. It will now be possible to proceed with the counting of the pieces.

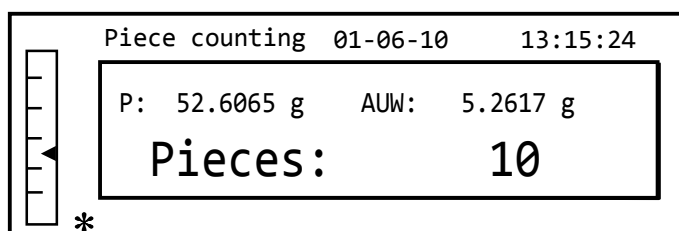

If the pieces to be counted have a weight that is too little with respect to the balance's resolution, an error message will be displayed.


In this case it will be necessary to turn to a balance with greater resolution.

If the weight of the samples is acceptable but not sufficient, the following message will be displayed: Add enough pieces so as to approximately double the quantity loaded on the plate, then press the **PRINT** button.

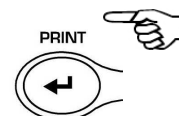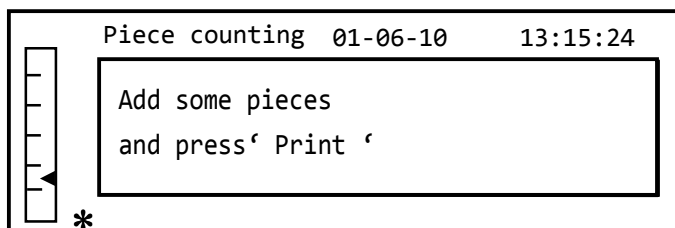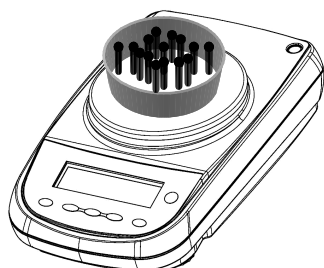

If the number of pieces is still insufficient, the message indicated above will be displayed again. Double the quantity of pieces loaded again.

Once a sufficient number of pieces has been reached, their number will be displayed and it will be possible to proceed with the counting, loading the pieces to be counted on the plate.

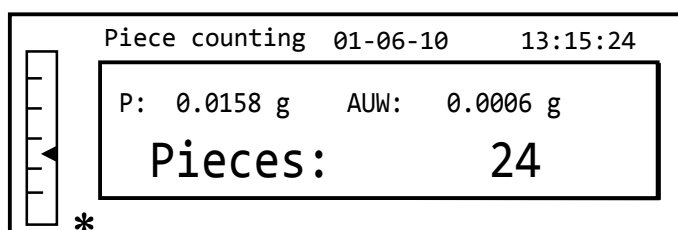


4. To exit from piece counting mode, press the **ON/OFF** button and the balance will return to the normal weighing conditions.


10.1.1 Manual insertion of the average unit weight

This function allows the user to enter, when known, the average unit weight of the piece, thus avoiding the sampling of the pieces.


1. Select the piece counting program as described in paragraph 10.
The following screen will be shown on the display:

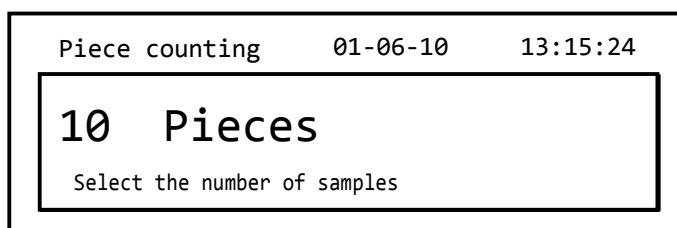


- Press the **MENU** button until the following message appears on the display:

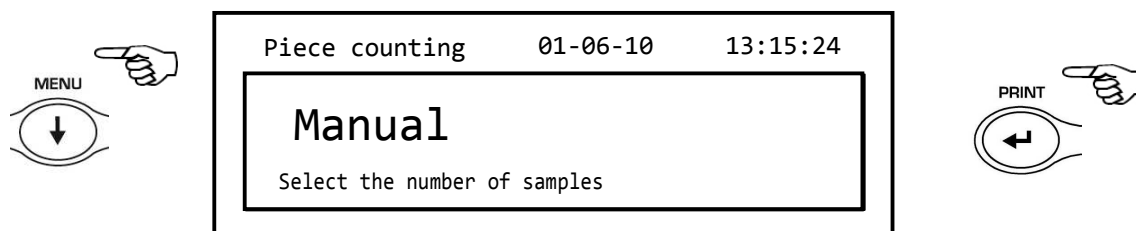


- Then press the **PRINT** button to confirm.

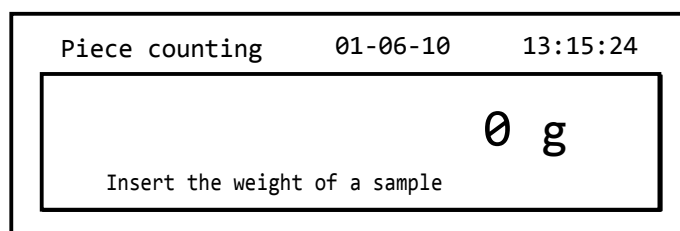

- Insert the piece's unit weight in grams using the **CAL** and **MENU** buttons to increase and decrease the value, while pressing the **O/T** button to pass to the next value. To insert a decimal point, hold down the **CAL** button for a prolonged time. During the entering phase, the prolonged pressing of the **O/T** button allows you to delete the inserted value.

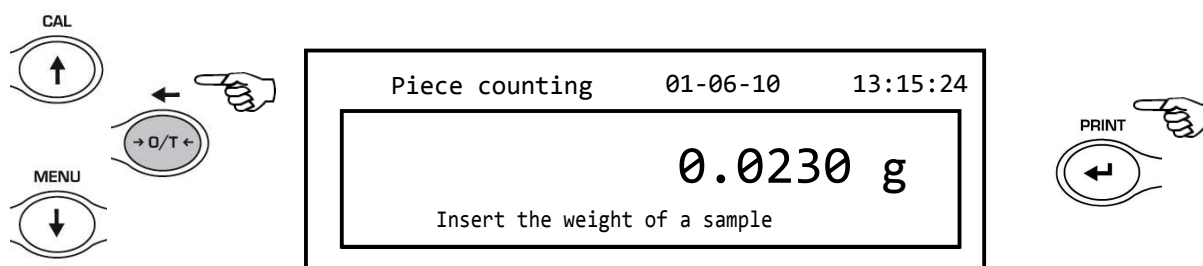


- Press the **PRINT** button to confirm.
 If the piece entered is less than 100 times the resolution of the balance, an error message will be displayed.
 To exit without entering the weight, press the **ON/OFF** button.
- If the weight is sufficient, "0" will be shown on the display; it is now possible to proceed with the counting, loading the pieces on the plate.

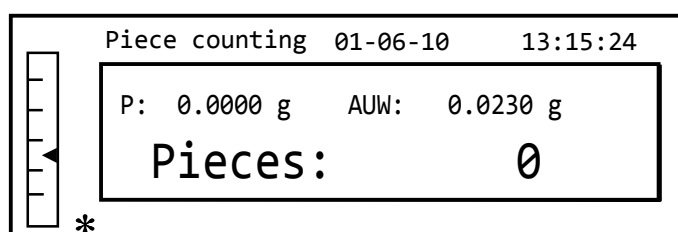

- To exit from the piece counting function, press the **ON/OFF** button.

It is also possible to use the optional alphanumeric keyboard to insert the average unit weight of the sample. In this case, carry out the same procedure described above to enter manual insertion mode.

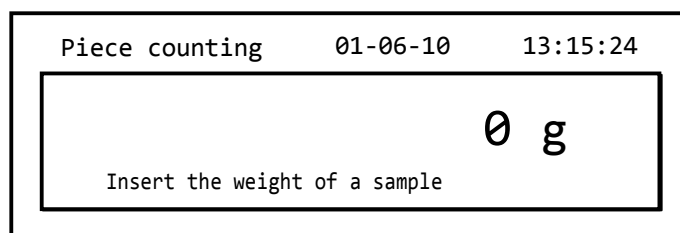

1. Insert the unit weight in grams of the sample by using the numeric keys from 0 to 9 and the decimal point.
In case of error, press the **CLEAR** button and restart.
2. Press the **INSER** button to confirm.
3. If the piece entered is less than 100 times the resolution of the balance, an error message will be displayed.
To exit without entering the weight, press **ESCAPE** (on the alphanumeric keyboard) or **ON/OFF**.
4. If the weight is sufficient, "0" will be shown on the display; it is now possible to proceed with the counting, loading the pieces on the plate.
5. To exit from the piece counting function, press the **ON/OFF** button.

10.1.2 Automatic updating of the average unit weight

After having carried out the sampling, the average unit weight can be updated in the following way.

1. Instead of loading all of the pieces to be counted, load a number of pieces approximately double that of those loaded on the plate and wait for the beep.
2. This procedure can now be repeated up to a maximum of 255 pieces or you can proceed with the normal counting of the pieces.
This mechanism allows for a more accurate estimate of the average unit weight and a better precision in the counting of the pieces.

NOTE: the automatic updating mechanism is not active if the sampling has been carried out through insertion of the average unit weight.


10.2 Program for the determination of the density of a solid or a liquid

The density calculation program allows the determination of the density of a solid or liquid through the use of the lower weighing hook or the hydrostatic kit (optional Code T027).


10.2.1 Solid density determination

1. Select the density program as described in paragraph 10.


The following screen will be shown on the display:

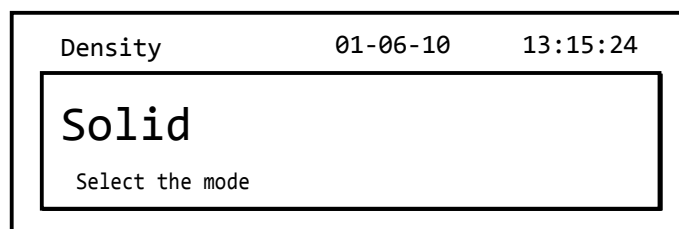

2. Then press the **PRINT** button to confirm the selection.
3. The density value of the liquid to be used will be displayed. The default value is equal to 1.0000 (distilled water at 20°C).

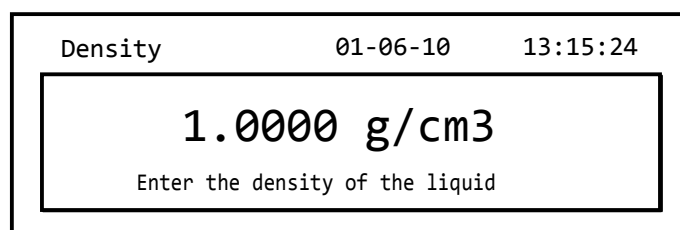


4. It is possible to insert a different value using the **CAL** and **MENU** buttons to increase and decrease the value, while pressing the **O/T** button to pass to the next value. During the entering phase, prolonged pressure on the **O/T** button allows you to cancel the inserted value.

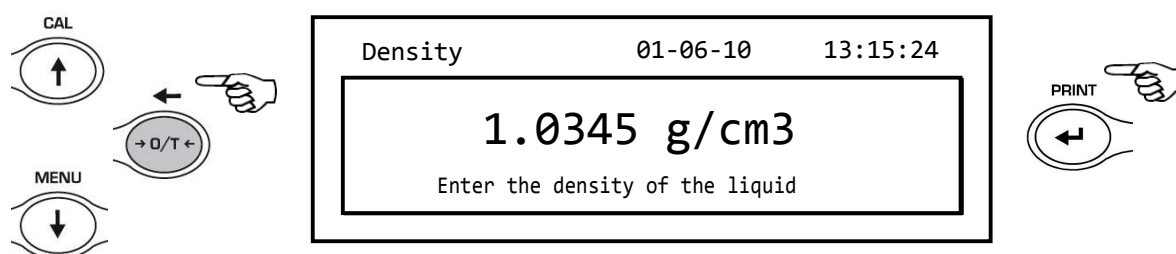


Note: It is also possible to set the value using the optional alphanumeric keypad.


5. Once the desired value has been set, press the **PRINT** button.
6. It will now ask you to weigh the solid in the air.

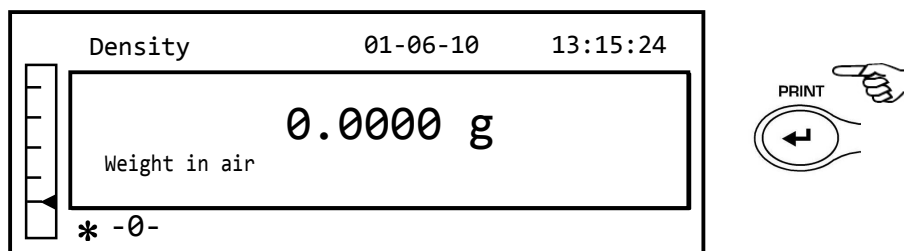


7. Carry out a tare if necessary and load the solid. Wait for the stability symbol to appear and press the **PRINT** button to acquire the value. The word 'wait...' will appear during the acquisition of the value.
8. The weight of the solid in the liquid will then be requested. Carry out the tare of the drum in the liquid. Put the solid in the drum, immerse the solid, and wait for the stability indicator to appear. Then press the **PRINT** button. The word 'wait...' will be displayed during the acquisition of the value.

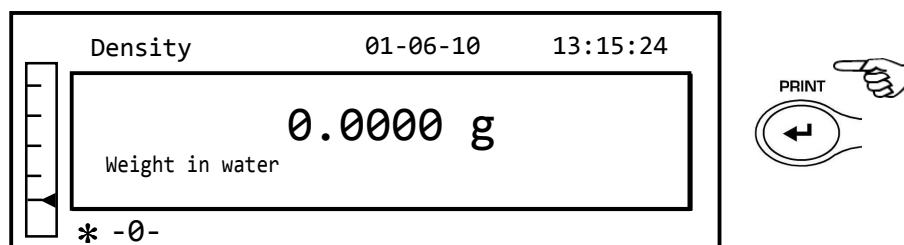

9. The result of the density calculation of the solid will now be displayed. If the balance is equipped with a printer, it will be possible to print the density value by pressing the **PRINT** button.

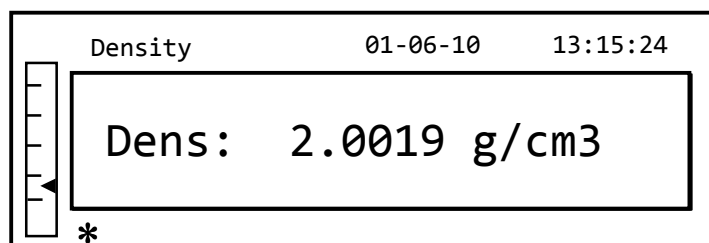


10. The following string will be shown on the display in case of error:

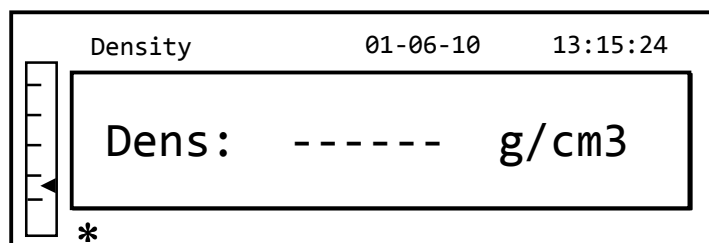



11. Now press the **ON/OFF** button to exit from the density function, or the **MENU** button to carry out the density measurement for another solid.

10.2.2 Liquid density determination


1. Select the density program as described in paragraph 10.
The following screen will be shown on the display:

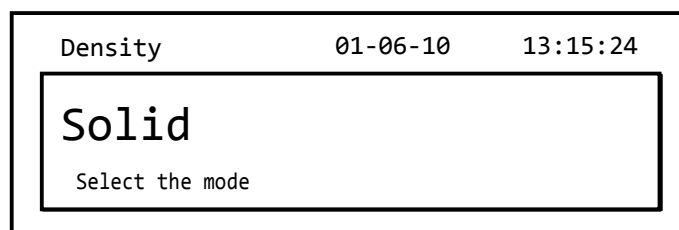




2. Press the **MENU** button select liquid mode. Then press the **PRINT** button to confirm.

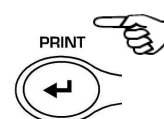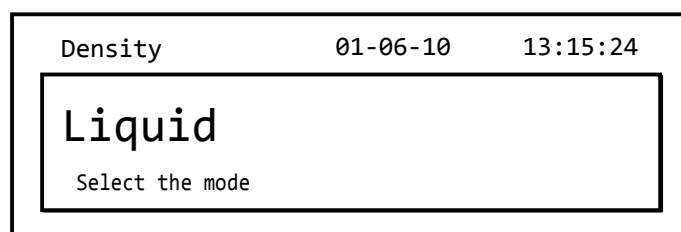

3. The default value of the solid's density will be displayed. The default value is equal to 3.0000 g/cm^3 .


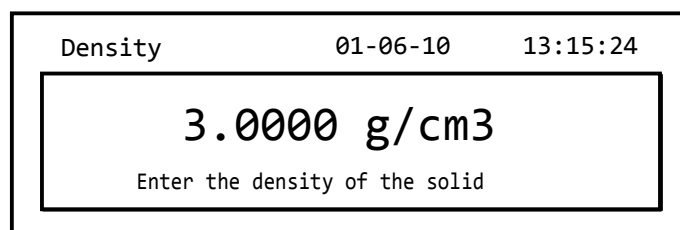


4. A different value can be entered by using the **CAL** and **MENU** buttons to increase and decrease the value, while pressing the **O/T** button to pass to the next value. During the entering phase, prolonged pressure on the **O/T** button allows you to delete the inserted value.

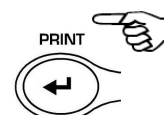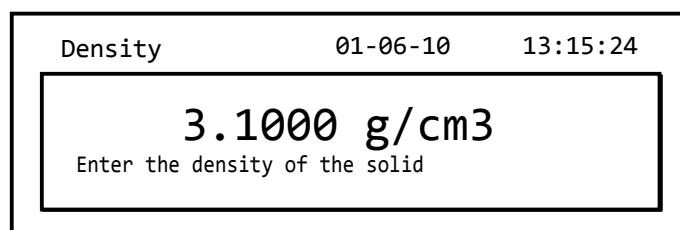
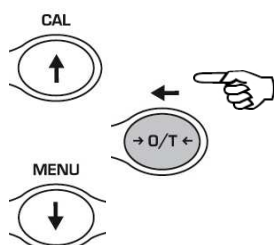


Note: It is also possible to set the value using the optional alphanumeric keypad.


5. Once the desired value has been set, press the **PRINT** button.
6. It will now ask you to weigh the dipstick in the air.

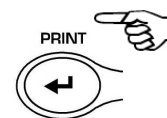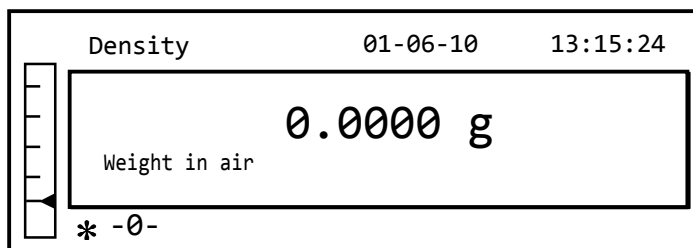


7. Carry out a tare if necessary and load the dipstick. Wait for the stability symbol to appear and press the **PRINT** button to acquire the value. The word 'wait...' will appear during the acquisition of the value.
8. The weighing of the dipstick immersed in the liquid will then be requested. Then immerse the solid in the liquid, wait for the stability indicator to appear, and then press the **PRINT** button. The word 'wait...' will be displayed during the acquisition of the value.

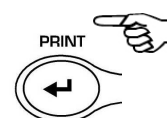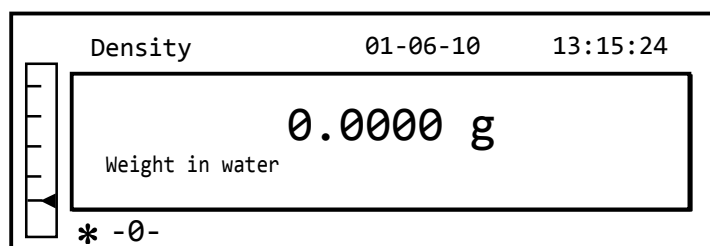

9. The result of the density calculation of the liquid will now be displayed. If the balance is equipped with a printer, it will be possible to print the density value by pressing the **PRINT** button.

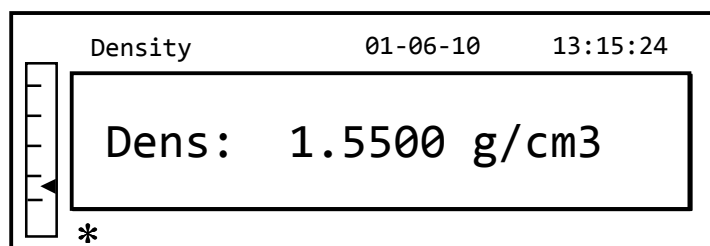

10. The following string will be shown on the display in case of error:

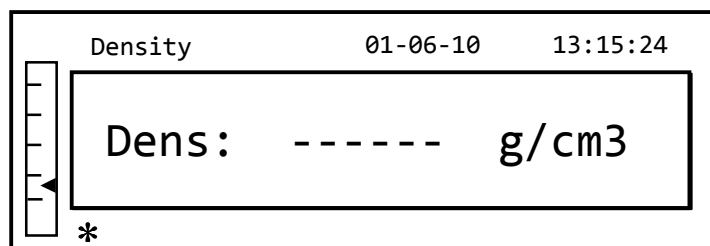


11. Now press the **ON/OFF** button to exit from the density function, or the **MENU** button to carry out the density measurement for another liquid.


10.3 Formulation function

10.3.1 Manual formulation


1. Select the formulation program as described in paragraph 10.
The following screen will be shown on the display:

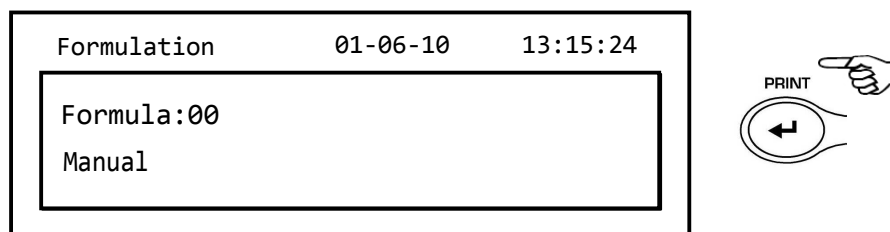


2. Then press the PRINT button to confirm the selection.

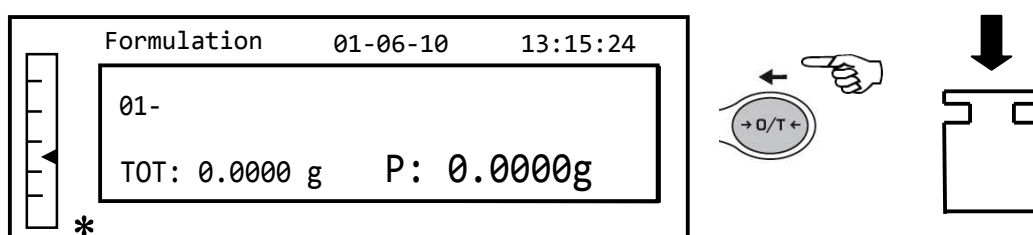


3. Carry out a tare operation if necessary, and load the first ingredient.

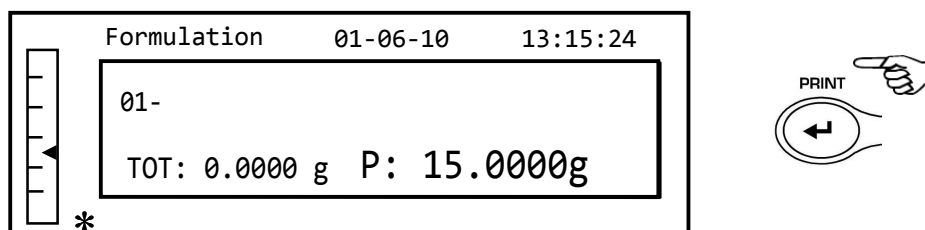

4. Then press the **PRINT** button to confirm.

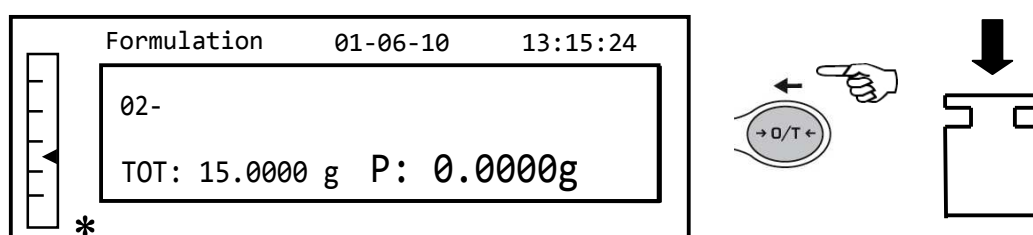


5. Carry out a tare operation if necessary, and load the second ingredient.

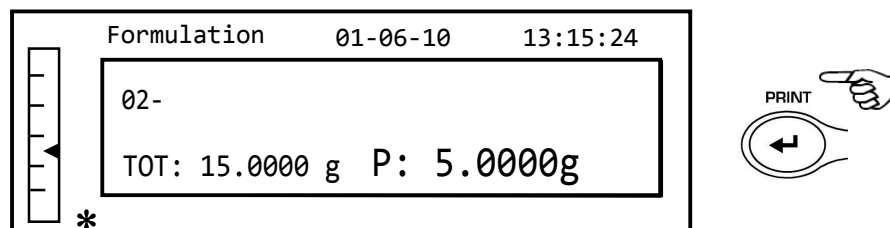

6. Then press the **PRINT** button to confirm.
7. Repeat the operation for a maximum number of 99 ingredients.

Note: During the acquisition of the ingredient, the display of Err10 indicates a negative weight value. Check not to have made a mistake with the ingredient loading and zeroing procedure.


- To end, print the value of the individual components and the total value, and press and keep pressed the **PRINT** button until the beeping stops. The display will show the following screen:

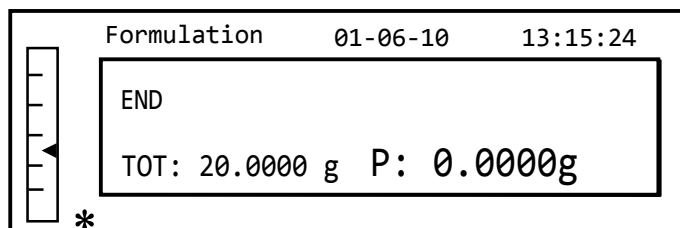


- To exit from the screen and carry out a new formulation, press the **ON/OFF** button once. To exit from the program and return to the weighing screen, press the **ON/OFF** button two consecutive times.

10.3.2 Formula saving

- Select the formulation program as described in paragraph 10. The following screen will be shown on the display:

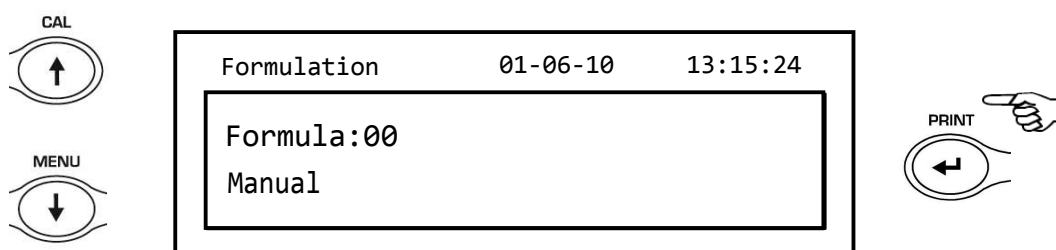


- Select the number of the formula to save or modify using the **CAL** and **MENU** buttons to increase and decrease the value, after the word 'Formula'. Then press and keep pressed the **PRINT** button until the beeping stops to confirm the selection and enter the 'setup formula' menu.

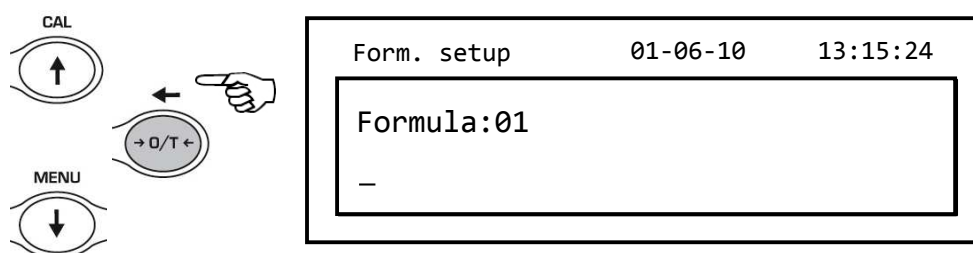

- Enter the name of the formula (it can be a series of numbers or letters, max 20 characters) using the **MENU** or **CAL** buttons to scroll all of the available characters, and the **TARE** button to pass to the next character. To select the uppercase or lowercase character, press and hold the **MENU** button until the beeping stops.

Note: It is also possible to set the value using the optional alphanumeric keypad.

Form. setup	01-06-10	13:15:24
Formula:01 Torta		

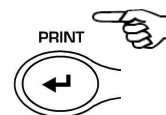

- Press the **PRINT** button to confirm.

CAL
↑

→ 0/T ←

MENU
↓

Form. setup	01-06-10	13:15:24
Cake 01 _		

PRINT

- Enter the name of the first component (it can be a series of numbers or letters, max 11 characters) using the **MENU** or **CAL** buttons to scroll through the available characters.
- Then press the **PRINT** button to confirm and save the value.

CAL
↑

→ 0/T ←

MENU
↓


Form. setup	01-06-10	13:15:24
Cake 01- Flour 10.0000g T- =- 1.0 % T+ = + 2.0 %		

PRINT


- Now enter the quantity of the component using the **MENU** or **CAL** buttons to increase or decrease the value while pressing the **0/T** button to pass to the next value and the **PRINT** button to pass to the next parameter
- Now enter the negative tolerance and press **PRINT** button to pass next parameter
- Now enter the positive tolerance.
- Then press the **PRINT** button to confirm and save the value.
- Repeat the operation described from point 5 to point 10 to enter all of the desired components up to a maximum of 20.
- After having entered all of the desired components press the **ON/OFF** button to exit from the formula saving procedure.

10.3.3 Formula recall


1. Select the formulation program as described in paragraph 10.
The following screen will be shown on the display:

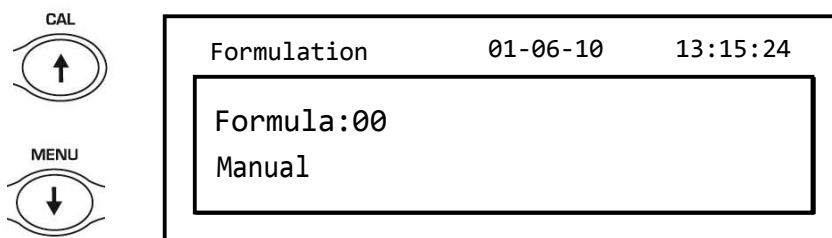

2. Choose the name of the formula (previously saved) using the **CAL** and **MENU** keys to scroll through the various formulas inserted.

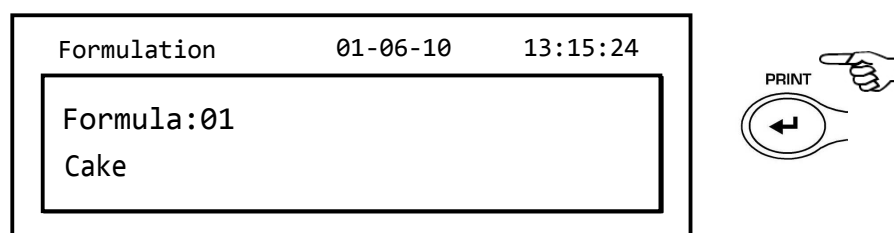


3. Then press the **PRINT** button to confirm the selection.

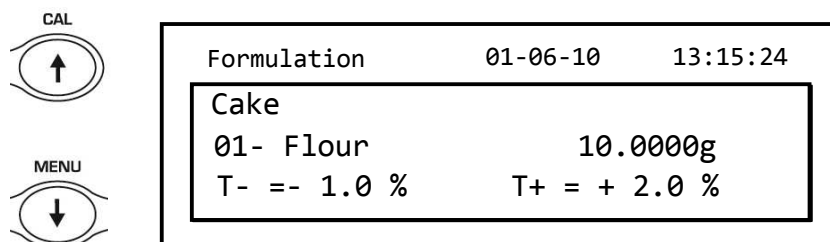

4. It will now be possible to display the various components and the relative quantities of the selected formula using the **MENU** and **CAL** buttons.


5. Press the **PRINT** button again to insert the constant multiplier.

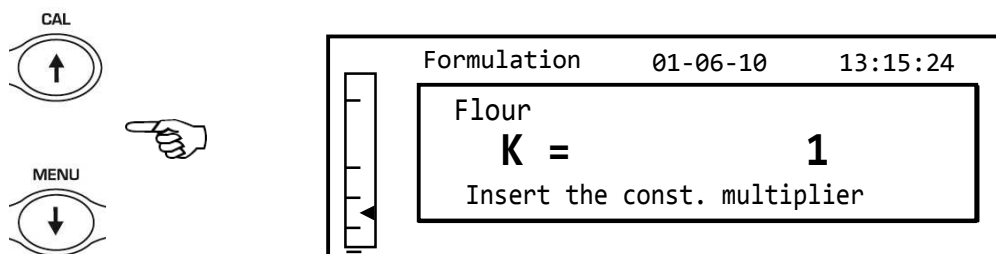

6. Insert now the multiplicative constant K to determine the desired amount of product. Use the **MENU** or **CAL** buttons to increase or decrease the value.


Example: if the entered formula is for 100g of product, inserting $K = 2$ the values of all components will be recalculated to obtain a total amount of product equal to 200g.

7. Press the **PRINT** button again to begin weighing the various components. If necessary, carry out the tare operation before measuring out the quantity of component indicated at the top right of the display

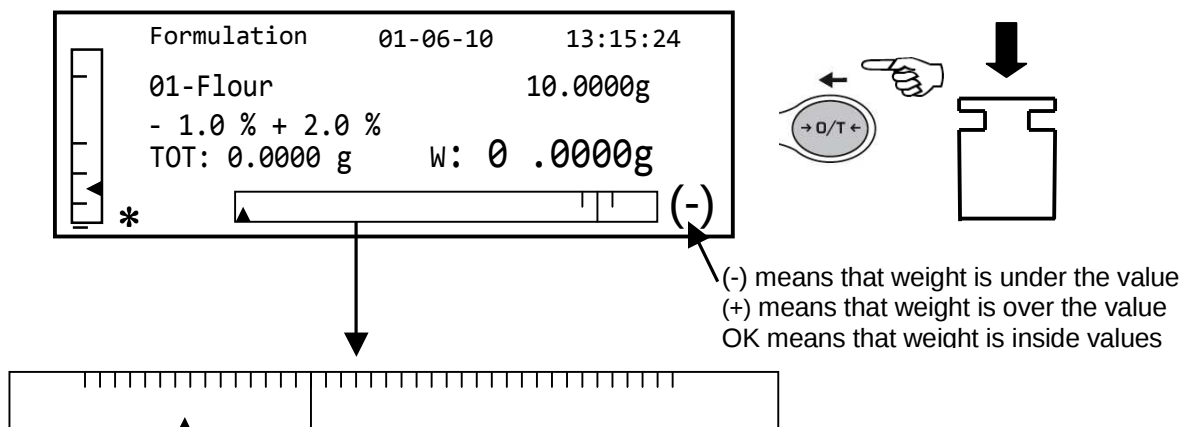


To facilitate the dosing operation, when the value of the component is approaching the threshold of the acceptable value, the dosing bar will automatically zoom.

8. Then press the **PRINT** button to pass to the next component.

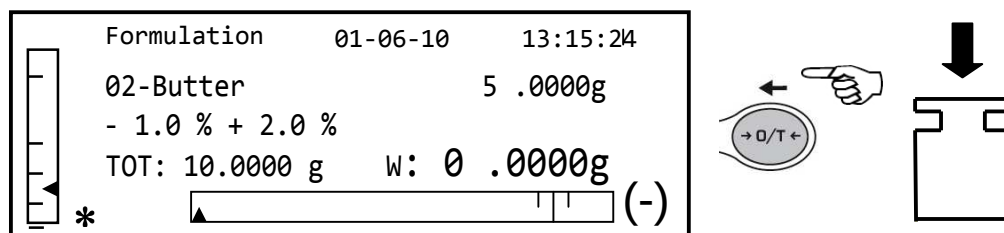

9. If necessary, carry out the tare operation before measuring out the quantity of component indicated at the top right of the display.
10. Then press the **PRINT** button to pass to the next component.
11. Repeat the procedure until the last component, after which the weights of the single components measured and the total weight will be printed if the instrument is equipped with a printer. The display will show the following screen:

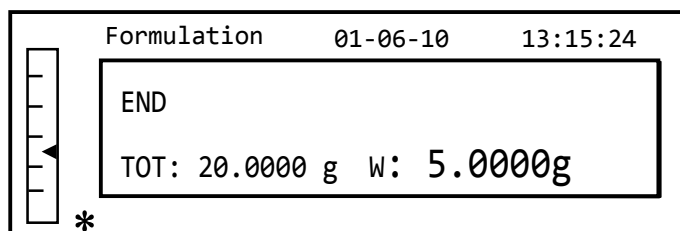


12. To exit from the screen and carry out a new formulation, press the **ON/OFF** button once.
To exit from the program and return to the weighing screen, press the **ON/OFF** button two consecutive times.

To interrupt and exit from the formulation function at any time, press the **ON/OFF** button.


10.4 Max-Min thresholds function.

The thresholds function allows you to determine if the weight loaded on the plate is above or below two thresholds pre-set by the user.

1. Select the thresholds function as described in paragraph 10.
The following screen will be shown on the display:

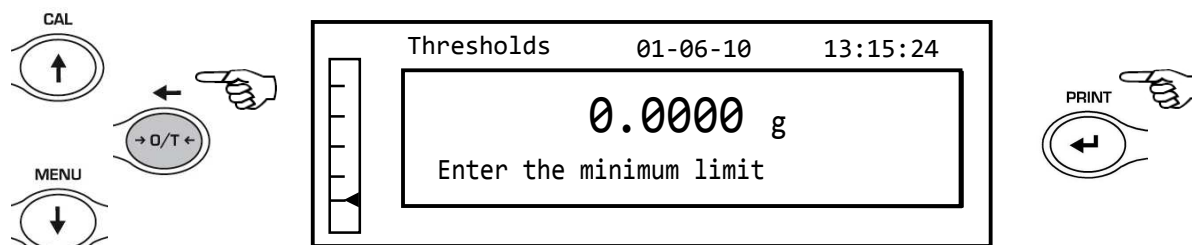


2. Enter the MINIMUM limit value by using the **CAL** and **MENU** buttons to increase and decrease the value, while pressing the **O/T** button to pass to the next number. During the entering phase, prolonged pressure on the **O/T** button allows you to delete the entered value.
3. Then press the **PRINT** button to confirm. The entered value will remain in memory until the balance is turned off.
4. The following screen will then be displayed.

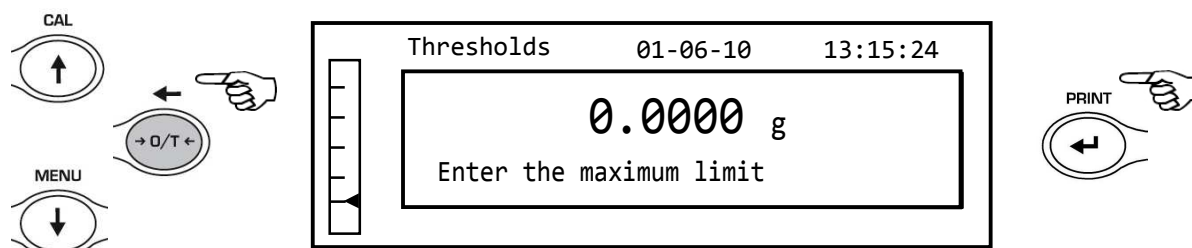


5. Now insert the MAXIMUM limit using the same procedure described for the insertion of the MINIMUM limit.
6. Then press the **PRINT** button to confirm. The entered value will remain in memory until the balance is turned off.

Note: It is also possible to set the value by using the optional alphanumeric keypad.

7. The following screen will then be displayed.

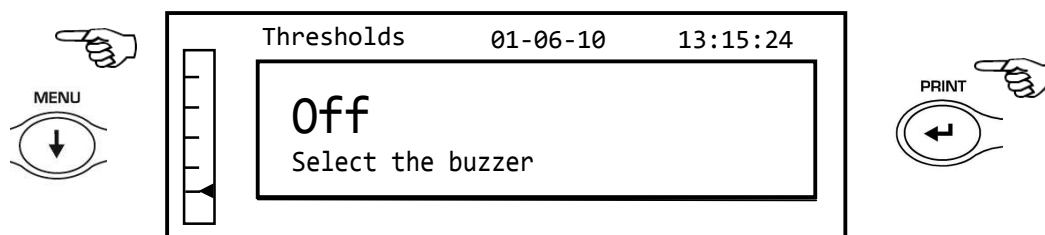

8. Through the **MENU** key, select the activation or not of the acoustic signal when the weight is within the two set limits. Then confirm the selection by pressing the **ENTER** button.
9. If the thresholds have been inserted correctly, the balance will return to weighing mode with an indication of the threshold status (+ MAX threshold, - MIN threshold, **OK** within the two limits sets).

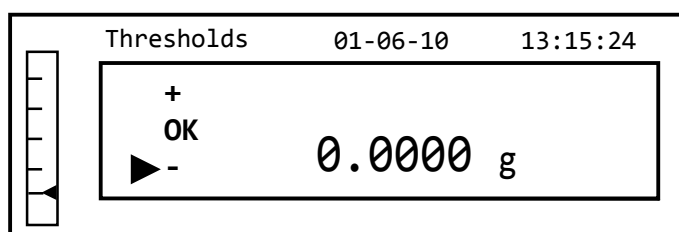

NOTE: If the values have not been set correctly, the word **ERROR 07** will be displayed.

The thresholds function has three operating modes.

10.4.1 With both the limits set

This mode allows to identify an acceptance range, inserting a lower limit and an upper limit, in which the value of weight is considered ok, identified by the "OK" symbol that is visualized on the display together with the acoustic signal, if activated. When the weight is under the value of the lower limit set, the symbol "L" is visualized on display, while if the value is over the upper limit set, the symbol "H" is visualized on display.

10.4.2 With only the lower limit set

When only the lower limit is set and the upper limit is left to zero, the weight is considered ok each time the value of weight is over the lower limit set, identified by the "OK" symbol that is visualized on the display together with the acoustic signal, if activated. When the weight is under the value of the lower limit set, the symbol "L" is visualized on display.

10.4.3 With only the upper limit set

When only the upper limit is set and the lower limit is left to zero, the weight is considered ok each time the value of weight is under the upper limit set, identified by the "OK" symbol that is visualized on the display together with the acoustic signal, if activated. When the weight is over the value of the upper limit set, the symbol "H" is visualized on display.

10.5 Percentage weighing function


This function allows you to read the weight as a percentage of a reference weight. The reference weight is assumed as the 100% value (factory setting).

There are two modes for the acquisition of the reference weight – an automatic one (with reference weight), and a manual one (with the manual entry of the value of the reference weight).

10.5.1 Automatic mode with reference weight


1. Select the percentage weighing function as described in paragraph 10.

The following screen will be shown on the display:

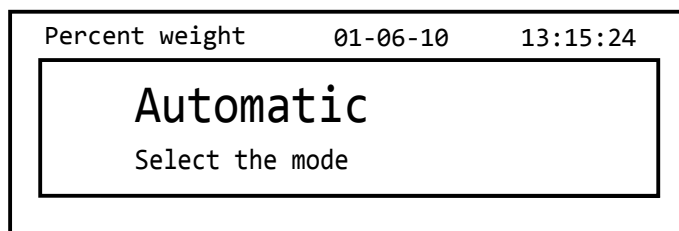


2. Confirm automatic mode by pressing the **PRINT** button.

3. The tare will be carried out and you will be asked to load the reference weight on the plate.

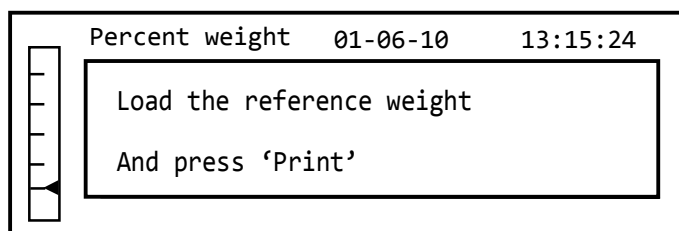

4. Load the reference weight on the plate and then press the **PRINT** button; the word "Wait" will be shown. Once the weight is acquired, a screen with an indication of the weight loaded, reference weight, and percentage weight will be shown.

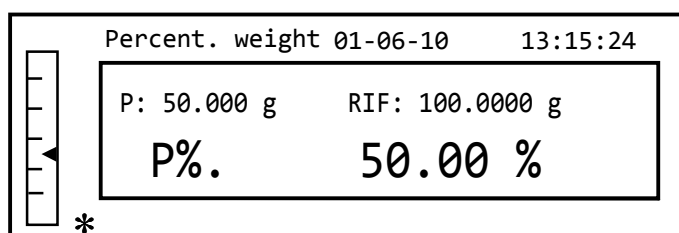


5. Now remove the reference weight, load the sample and read the percentage weight.

6. Press the **ON/OFF** button to exit from the percentage weighing function.


NOTE: If the reference weight entered is less than 10 displayed digits, the word **ERROR 07** will be shown.

10.5.2 Mode with manual insertion of the reference weight.


1. Select the percentage weight function as described in paragraph 10.
The following screen will be shown on the display:

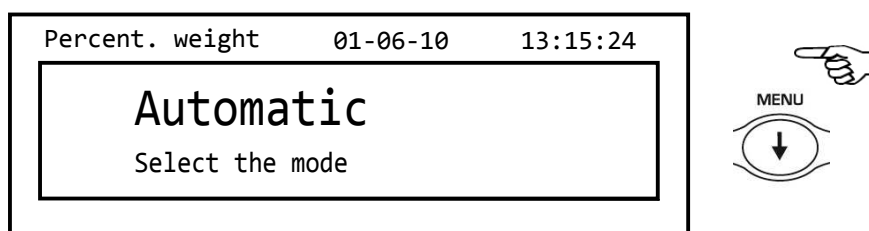


2. Press the **MENU** button to select manual mode

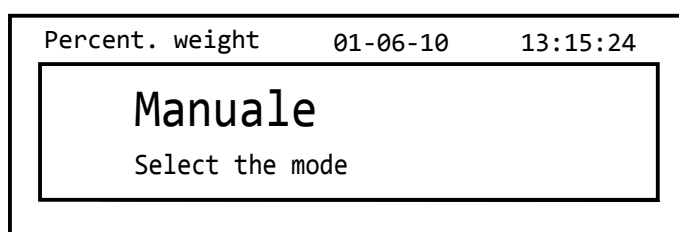

3. Confirm manual mode by pressing the **PRINT** button.

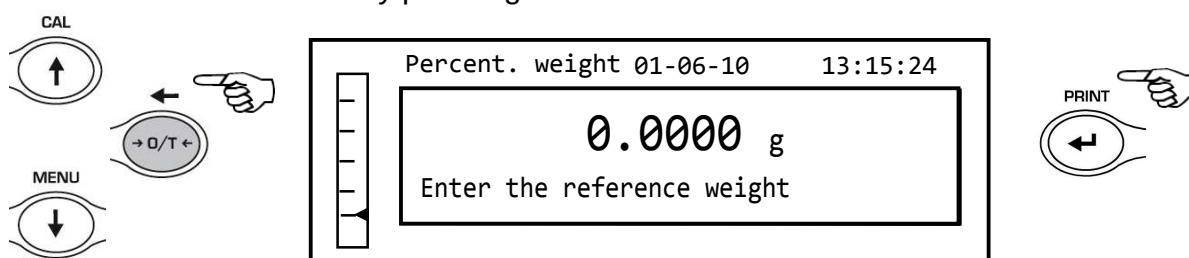

4. You can now enter the reference weight value, using the **CAL** and **MENU** keys to increase and decrease the value, while pressing the **O/T** button to pass to the next value. During the entry phase, holding down the **O/T** button allows you to delete the value entered. The value entered will remain in the memory until the balance is turned off. It is also possible to enter the value using the optional alphanumeric keypad.
5. After having inserted the desired reference weight value, press the **ENTER** key.

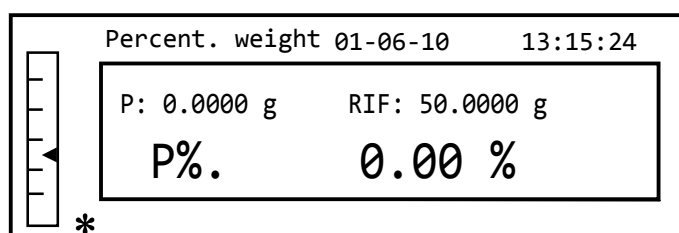


6. Now load the sample and read the percentage value.
7. Press the **ON/OFF** button exit from the percentage weighing function.

NOTE: If the reference weight entered is less than 10 displayed digits, the word **ERROR 07** will be shown.


10.6 Animal weighing function

Thus function allows you to acquire an averaged weight of moving objects or animals for a settable period of time.


1. Select the animal weighing function as described in paragraph 10.
The following screen will be shown on the display:

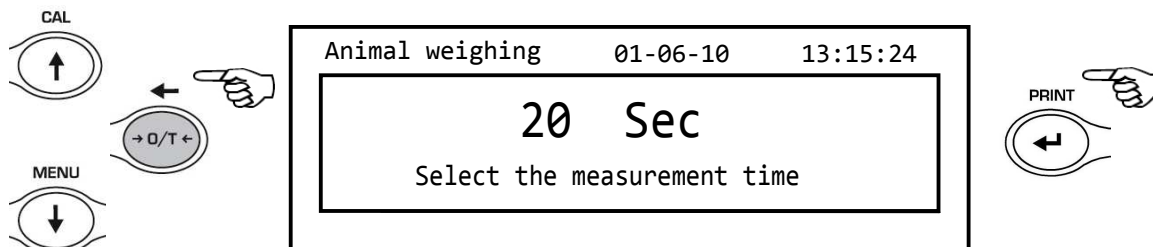


2. Set the desired time from 5 to 90 seconds using the **MENU** key to decrease and **CAL** to increase. Then confirm by pressing the **PRINT** button.

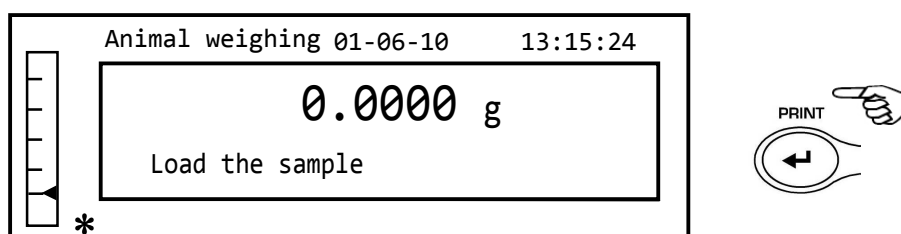

3. Load the sample to be weighed on the plate and press the **PRINT** button; the value of the current weight and the set sampling countdown time will be displayed.

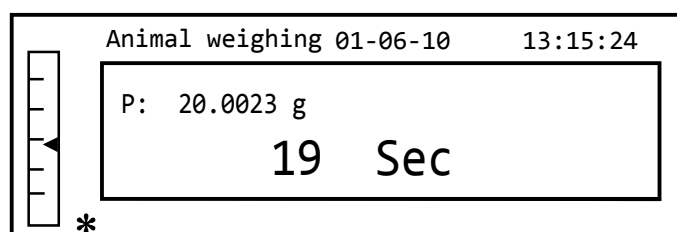

4. Once acquired, the weight will be shown on the display with an indication of the average weight detected.

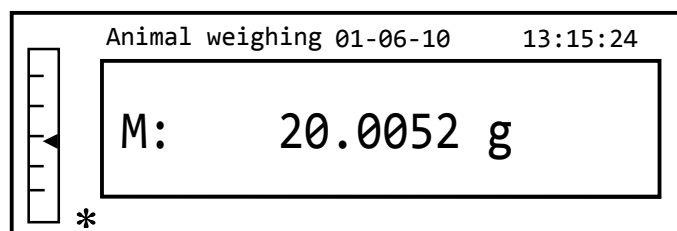


5. Press the **ON/OFF** button once to carry out another measurement, or twice to exit from the function.


10.7 Maximum load function

The “maximum load” function allows you to measure the maximum breakage load of a solid.

1. Select the maximum load function as described in paragraph 10.
A tare will automatically be carried out and the following screen will be shown on the display with an indication of the maximum load function at the top left:

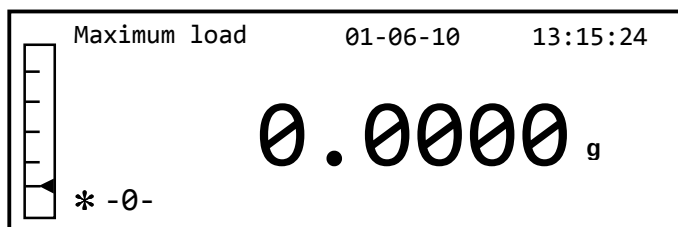

2. The breakage weight can now be detected.

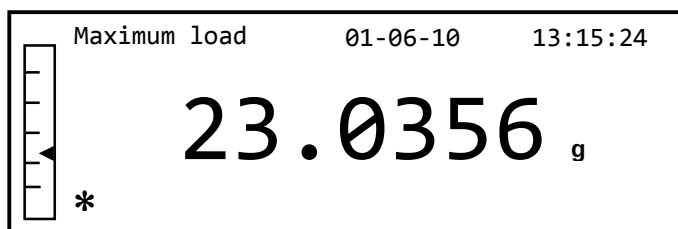


3. Press the **TARE** button to carry out another measurement.
4. Press the **ON/OFF** button to exit from the maximum load function.

10.8 GLP function (Good Laboratory Practices)


The “GLP” function allows you to save the identifying parameters of the instrument and operator to be able to print them along with the value of the test results.

1. Select the GLP function as described in paragraph 10.
The following screen will be shown:

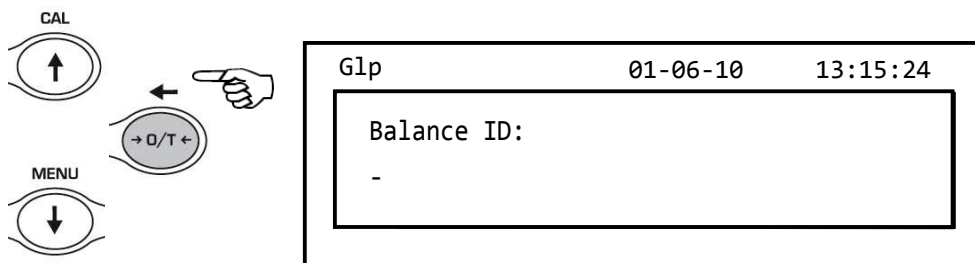


2. Enter the balance ID (it can be a series of numbers or letters, max 18 characters) using the **MENU** and **CAL** buttons to scroll through all of the available characters. To select uppercase or lowercase characters, press and hold the **MENU** button until the beeping stops.

Note: It is also possible to set the value by using the optional alphanumeric keypad.

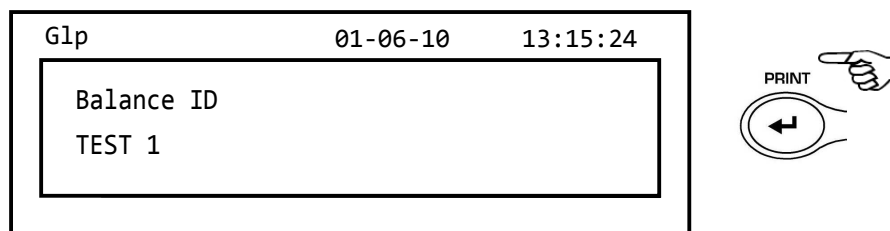


3. Enter the user ID (it can be a series of numbers or letters, max 18 characters) using the **MENU** and **CAL** buttons to scroll through all of the available characters.

Note: It is also possible to set the value by using the optional alphanumeric keypad.

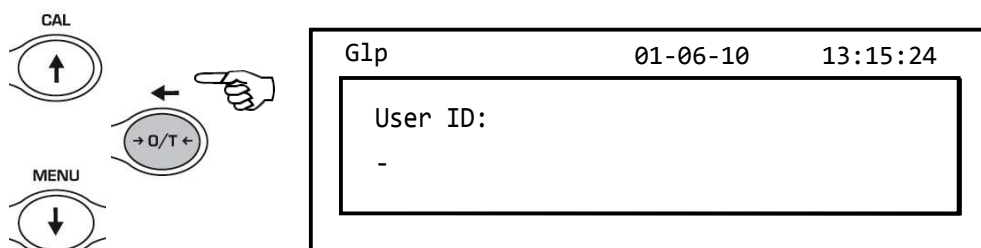


4. Confirm by pressing the **PRINT** button.

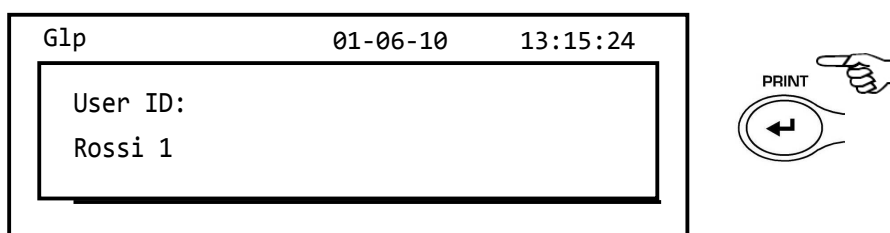


5. Enter the project identifier (it can be a series of numbers or letters, max 18 characters) using the **MENU** and **CAL** buttons to scroll through all of the available characters.

Note: It is also possible to set the value by using the optional alphanumeric keypad.

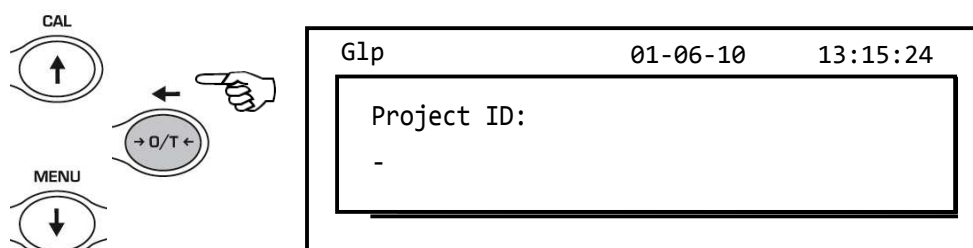

6. Then confirm all of the data entered by pressing and holding the **PRINT** button until the beeping stops.

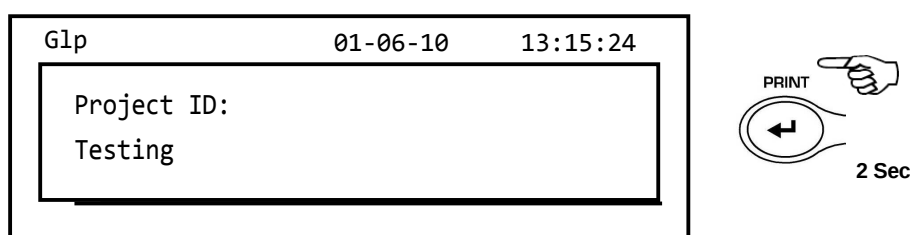


7. The balance will automatically return to the weighing screen.

11 RS232 Interface features

11.1 General features

The balance transmits the value visualized on the display following serial RS232C standard, allowing to print the value of weight to a PC monitor or to a serial printer. In the case of connection to a PC, it will be possible to select the transmission in continuous mode or transmission at user command through pressing of the PRINT button (as described at par.9.2). The balance is also capable of receiving commands, always through the standard RS232C, that allow performing all the functions available through the keyboard of PC itself. The speed of transmission and reception can be selected, as described previously (see par.9.3), to 1200, 2400, 4800, e 9600 baud. The character format is of 8 bit preceded by one bit of start and followed by a bit of stop. Parity is not considered.

11.2 Map of connector

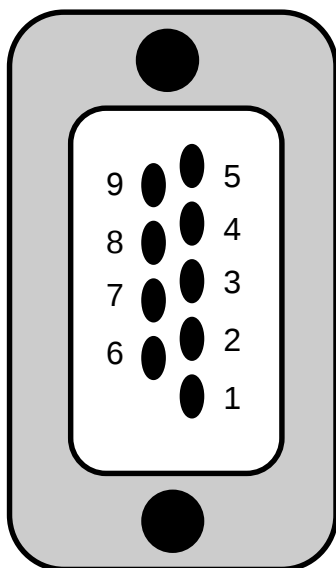

CONNECTOR MAP CONNECTION FOR KEYPAD, PC AND PRINTER


- pin 1 = Power +5v for keyboard
- pin 2 = signal Tx
- pin 3 = signal Rx
- pin 4 = busy signal
- pin 5 = Gnd
- pin 4-6 = connected eachother for transmission to PC

11.3 Connection of the balance to the optional alphanumeric keyboard

To connect the optional external alphanumeric keyboard to the balance, it must be used the same connector that is used to link the balance to PC. In this case the connection to PC or to printer must be effected through the connector placed on the optional alphanumeric keyboard.

11.4 Connection of the Balance to computer

To receive/transmit data, link the connector (n.2 of fig.1 at par.**Errore. L'origine riferimento non è stata trovata.**) of the balance to the serial port of your Personal Computer as shown below:

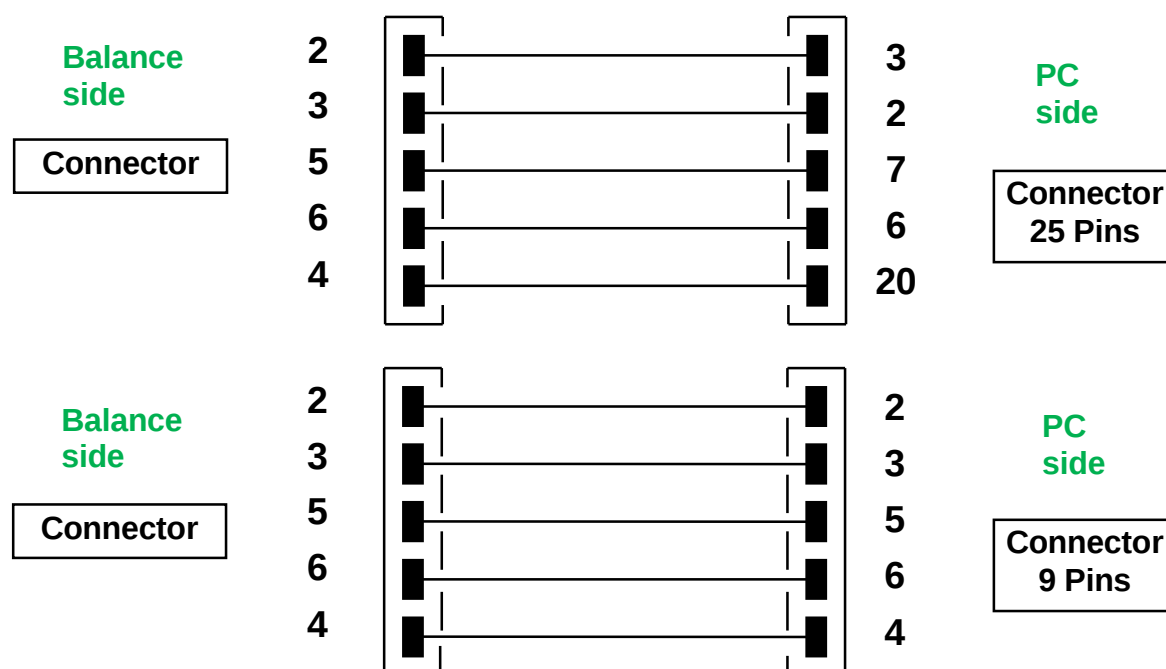

There are three ways of transmission in which the Balance and the computer can be interfaced:

- Continuous transmission of weight data (continuous mode must be set from the menu as explained in the paragraph 9.2).
- On demand transmission of weight data (on demand mode must be set from the menu as explained in the paragraph 9.2).
- On demand transmission with GLP of weight data (on demand mode with GLP must be set from the menu as explained in the paragraph 9.2).

In all the modes it is possible to execute all the balance's functions directly from the computer's keyboard, transmitting to balance the ASCII codes as shown in the table below.

CODE	1 st FUNCTION (SINGLE PRESS)
"T" = H54	TARE
"C" = H43	CALIBRATION
"E" = H45	ENTER
"M" = H4D	MENU
"O" = H4F	ON/OFF

CODICE	2 nd FUNCTION (PROLONGED PRESS)
"t" = H74	TARE
"c" = H63	CALIBRATION
"e" = H65	ENTER
"m" = H6D	MENU
"o" = H6F	ON/OFF

11.4.1 Continuous Transmission mode

String transmitted is composed by the following 14 characters:

- First character: weight sign (blank or -)
- Second to ninth character: weight or other data
- Tenth to twelfth character: weight unit symbol
- thirteenth character: stability indicator
- fourteenth character: carriage return
- fifteenth character: line feed

Eventual non-significative zero are put as spaces.

In the following table the various transmission formats are shown:

Weight mode (valid for both continuous and on demand transmission)

1°	2°	3°	4°	5°	6°	7°	8°	9°	10°	11°	12°	13°	14°	15°
Sign	weight								measure unit			Stability	CR	LF

11.4.2 On demand transmission mode

When in on demand mode, the transmitted data to computer do not include only informations of the weight value but also date/time and other informations that depend on the function you're currently using are transmitted to computer

Below are shown the data transmitted in each situation possible:

WEIGHT:

03-04-11 10:13:44

Weight: 0.00 g

PIECE COUNTING:

03-04-11 10:49:28

Pcs.: 10
Weight: 100.02 g
MPW: 10.00 g

DENSITY:

03-04-11 10:51:15

d: 1.4504 g/cm3

FORMULATION:

03-04-11 10:54:57

Manual

1. 31.05 g
2. 100.02 g
3. 26.89 g

T = 157.96 g

NOTE: To transmit the print of total of weights, press and keep pressed the PRINT button

THRESHOLDS:

Value under threshold

Value inside thresholds

Value over threshold

03-04-11 11:02:19

03-04-11 11:01:50

03-04-11 11:01:50

Lim.1 : 10.00 g
Lim.2 : 100.00 g
Weight: -0.01 g
TEST: KO! ---

Lim.1 : 10.00 g
Lim.2 : 100.00 g
Weight: 31.08 g
TEST: OK!

Lim.1 : 10.00 g
Lim.2 : 100.00 g
Weight: 131.10 g
TEST: KO! +++

PERCENTUAL WEIGHT

03-04-11 11:58:39

Perc. 100.0 %
Weight: 18.69 g
Refer.: 18.69 g

ANIMAL WEIGHING:

03-04-11 12:01:06

Time = 20 Sec
M: 56.53 g


MAXIMUM LOAD:

03-04-11 12:01:57

Max.: 2.76 g

11.4.3 On demand transmission with G.L.P.


In the on demand transmission with G.L.P., the data transmitted to computer are the same as described as in the on demand transmission without G.L.P. mode but with the add of G.L.P. parameters before of each transmission, as described below:

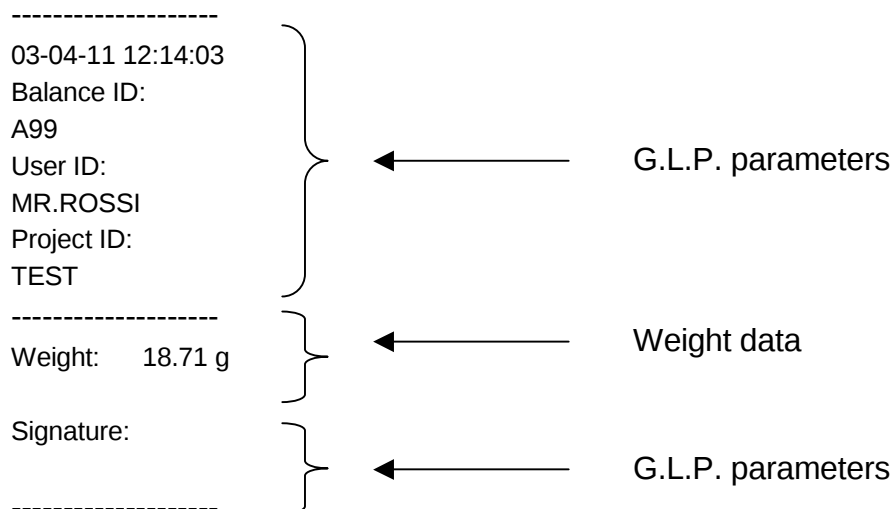

11.5 Connection of balance with serial printer

It is possible to connect the balance to a printing peripheral.

To print the weight, connect the connector of the balance to the serial printer as shown in the scheme below:

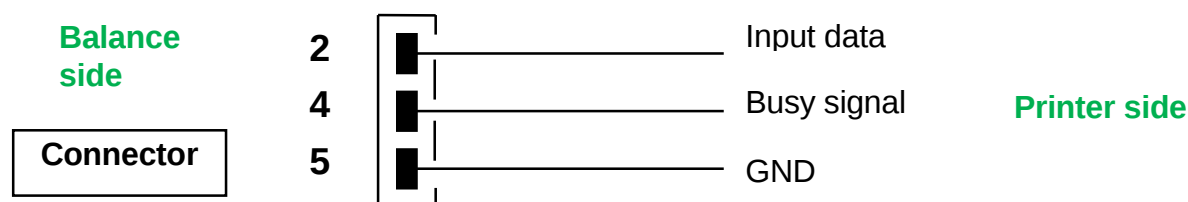

Here you can find description of the several modes of printing that can be selected:

- Print of weight data with a generic serial printer (from the menu, set the generic printing mode as described in the paragraph 9.2 and manage the busy signal).
- Print of weight data together with GLP indications with generic serial printer (from the menu, set the generic printing-GLP mode as described in the paragraph 9.2 and manage the busy signal)
- Print of weight data with printer model TLP50 (from the menu, set the printer TLP mode as described in the paragraph 9.2).
- Print of weight data together with GLP indications with printer model TLP50 (from the menu, set the printer TLP - GLP mode as described in the paragraph 9.2).

Note: In all different printing modes just described, if the weight is not stable during transmission of data to printer, an acoustic signal is emitted and ERR05 is displayed and weight is not printed.

11.5.1 PRINT FORMATS

Here are described the different types of print, depending on the print mode and on the function selected:

Generic printing or TLP 50 printer

Weighing mode:

03-04-11 10:13:44

Weight: 0.00 g

Piece counting:

03-04-11 10:49:28

Pcs.: 10
Weight: 100.02 g
MPW: 10.00 g

Density:

03-04-11 10:51:15

d: 1.4504 g/cm³

Formulation:

03-04-11 10:54:57

Manual

1. 31.05 g
2. 100.02 g
3. 26.89 g

T = 157.96 g

NOTE: To transmit the print of total of weights, press and keep pressed the PRINT button

Thresholds:

Value under threshold

Value inside thresholds

Value over threshold

03-04-11 11:02:19

03-04-11 11:01:50

03-04-11 11:01:50

Lim.1 : 10.00 g

Lim.2 : 100.00 g

Weight: -0.01 g

TEST: KO! ---

Lim.1 : 10.00 g

Lim.2 : 100.00 g

Weight: 31.08 g

TEST: OK!

Lim.1 : 10.00 g

Lim.2 : 100.00 g

Weight: 131.10 g

TEST: KO! +++

Percentual weight:

03-04-11 11:58:39

Perc. 100.0 %

Weight: 18.69 g

Refer.: 18.69 g

Pesata animali:

03-04-11 12:01:06

Time = 20 Sec

M: 56.53 g


Maximum load:

03-04-11 12:01:57

Max.: 2.76 g

11.5.2 Generic Printer or TLP 50 printer with G.L.P.

In the print mode with G.L.P. the printed data are the same as shown in the print mode without G.L.P. but with the add of G.L.P. paramters as shown below:

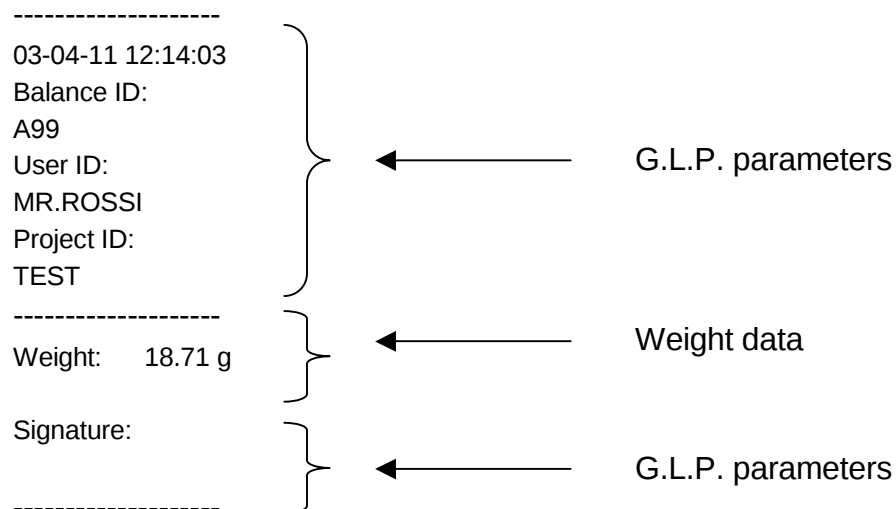

12 Error codes

ERROR MESSAGE ON DISPLAY	MEANING	POSSIBLE SOLUTIONS
ERR01	Weight not stable after operation of tare	Protect the balance from air flows or from vibrations of the working table
ERR02	impossible to start the calibration due to instability of the balance	Protect the balance from air flows or from vibrations of the working table.
ERR03	calibration weight not correct or balance unstable	Calibrate with correct weight or protect the balance from environment disturbs. In models with internal calibration remove the screw in the left lower part of the balance (see par5)
ERR04	weight of samples for the piececounting function not adequate or unstable	Select a bigger number of samples or protect the balance from vibrations.
ERR05	impossible to print because of weight unstable	Protect the balance from environment disturbs
ERR06	Weight cannot get stable in density mode	Protect the balance from environment disturbs
ERR07	Weight cannot get stable in percentual weighing mode	Protect the balance from environment disturbs
ERR08	Anomaly on autocalibration motor	Contact service staff
ERR09	Weight cannot get stable in formulation mode	Protect the balance from environment disturbs
ERR10	Weight of component out of tolerance in formulation mode	Reduce quantity
ERR F	Flash memory damaged	Ask assistance to authorized service staff
“UNLOAD”	weight loaded on the pan or pan not positioned properly	Remove the weight from the pan or position properly the pan and underpan.
“CAL But”:	the balance requires to be re-calibrated	Unload weights, if any, on the pan, and press the CAL button
	Over-range condition	Unload the weights loaded on the pan
	Under-range condition	Place properly pan and underpan

13 Maintenance and care

Regular maintenance of yours balance guarantee accurate measurements.

- **Cleaning**

Before cleaning the balance unplug the power supply of the balance from the voltage supply of your room. Do not use aggressive cleaning product (as solvents or similar), use a humid towel with soft detergent, Avoid liquids to go inside the instrument during the cleaning. Wipe the balance with a soft towel. Parts of samples or powder can be removed using a brush or vacuum cleaner.

- **Safety checks**

Safety of the instrument is no more guaranteed when:

- balance power supply is clearly damaged
- balance power supply is not working anymore
- balance power supply is stored for long time in hard environment conditions.

In these instances refer to the assistance centre where specialized technician will make reparations to bring back the instrument in the safety conditions eventually.


14 Warranty

- Duration of warranty is of 24 months from the date of purchase shown in the invoice concerning the product or by delivery note.
- Warranty covers all parts resulting defective at the origin. It does not cover mechanical or electronic parts damaged by wrong installation, tampering or incorrect use.
- Warranty does not cover damages caused by impacts, balance drops or drop of objects on weighing pan.
- Shipment to and from service centre is at customer charge.

15 Quick guide to balance parameters setup

- To enter the balance parameters setup menu, press and keep pressed the **MENU** button until the acoustic alarm gets mute.
- Use then the **MENU** button to go to next parameter, use the **CAL** button to go to previous and the **PRINT** button to confirm the choice.
- To escape from menu, press and keep pressed the **MENU** button until the acoustic alarm gets mute.

Menu settaggi Parametri bilancia	Measuring units			Grams
	Serial output			Carats
	Baud rate	1200 Baud rate	Tlp Printer	Ounces
		2400 Baud rate	Glp on demand	Pounds
		4800 Baud rate	Generic Printer-	Penny Weights
		9600 Baud rate	Glp	Troy Ounces
	Auto zero	Autozero Off	Generic Printer	Grains
	Filter	Filter 1	Glp on demand	HongKong Tael
		Filter 2	Continuous	Singapore Tael
		Filter 3	Tlp-Glp Print	Taiwan Tael
	Stability	Stability 1		Momme
		Stability 2		
		Stability 3		
	Contrast regulation	Contrast 0-15		
	Backlight	Disabled		
		2 Minutes		
		5 Minutes		
		15 Minutes		

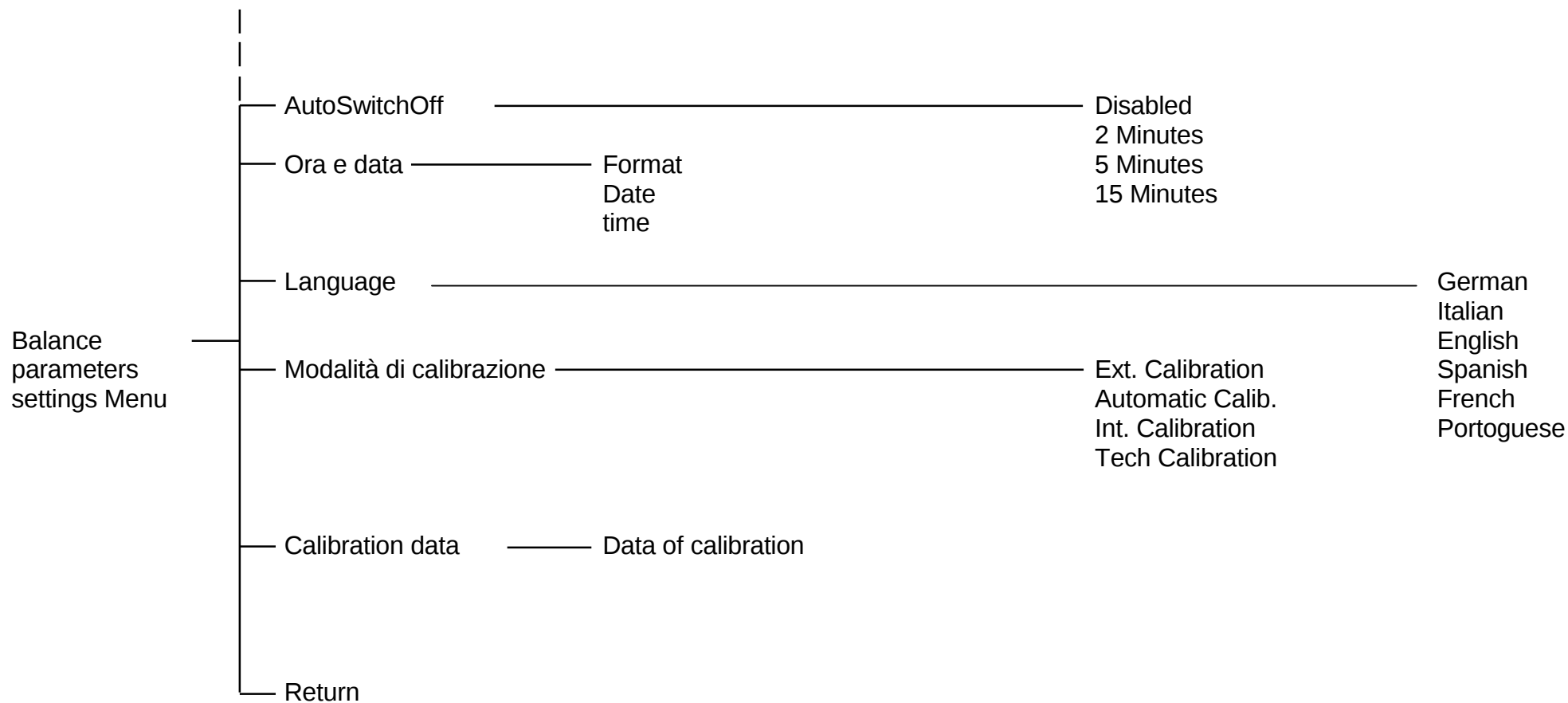

16 Balance Technical features

For models with power supply 9 V:

All the models listed are only for internal use. Maximum altitude using limit: 4000m. Pollution level: 2. Over voltage category: II

Power supply provided:	INPUT: 230V ~ 50Hz o 115V ~ 60Hz, OUTPUT: 9V DC 1000mA, Max power absorbed 9VA
Enviroment conditions adaption:	Filters selection
Autozero:	Selectable from Menu
Serial output:	RS232C
Operating temperature:	+5°C - +35°C

For models with power supply 24 V:

All the models listed are only for internal use. Maximum altitude using limit: 4000m. Pollution level: 2. Over voltage category: II

Power supply provided:	INPUT: 230V ~ 50Hz o 115V ~ 60Hz, OUTPUT: 24V DC 550mA, Max power absorbed 13.2VA-
Enviroment conditions adaption:	Filters selection
Autozero:	Selectable from Menu
Serial output:	RS232C
Operating temperature:	+5°C - +35°C

17 Optional Accessories

ACCESSORIES	
Code	DESCRIPTION
630-1485	TLP-50 Serial Printer, with date/time (connection cable included)
630-1486	TLP-50 Serial Printer, with date/time (connection cable included) Interface UK
611-2794	DPP-250 Serial Printer (can work also with batteries, connection cable included)
611-2793	BLUETOOTH printer DPP-250-BT (works also w/ battery) for Tablet system
630-1488	Roll of paper for TLP-50 printer
630-1487	Roll of 54x24mm labels (500) for TLP-50 printer
611-2795	Roll of paper for DPP250 printer
630-1488	Alphanumeric external keyboard mod. Keymodule
611-2786	RS232/USB CONVERTER
611-2787	Serial cable 9pins M/M for serial output printer/PC
611-2553	Density KIT for solids and liquids for balances 0.001g and 0.0001g
611-2552	Density KIT for solids and liquids for balances 0.01g
611-8684	Portable rechargeable battery pack

18 Technical service

Web Resources

Visit the VWR's website at www.vwr.com for:

- Complete technical service contact information
- Access to VWR's Online Catalogue, and information about accessories and related products
- Additional product information and special offers

Contact us For information or technical assistance contact your local VWR representative or visit. www.vwr.com.

19 Warranty

VWR International warrants that this product will be free from defects in material and workmanship for a period of two (2) years from date of delivery. If a defect is present, VWR will, at its option and cost, repair, replace, or refund the purchase price of this product to the customer, provided it is returned during the warranty period. This warranty does not apply if the product has been damaged by accident, abuse, misuse, or misapplication, or from ordinary wear and tear. If the required maintenance and inspection services are not performed according to the manuals and any local regulations, such warranty turns invalid, except to the extent, the defect of the product is not due to such non-performance.

Items being returned must be insured by the customer against possible damage or loss. This warranty shall be limited to the aforementioned remedies. IT IS EXPRESSLY AGREED THAT THIS WARRANTY WILL BE IN LIEU OF ALL WARRANTIES OF FITNESS AND IN LIEU OF THE WARRANTY OF MERCHANTABILITY.

20 Compliance with local laws and regulations

The customer is responsible for applying for and obtaining the necessary regulatory approvals or other authorizations necessary to run or use the Product in its local environment. VWR will not be held liable for any related omission or for not obtaining the required approval or authorization, unless any refusal is due to a defect of the product.

21 Equipment Disposal



This equipment is marked with the crossed out wheeled bin symbol to indicate that this equipment must not be disposed of with unsorted waste.

Instead it's your responsibility to correctly dispose of your equipment at lifecycle -end by handing it over to an authorized facility for separate collection and recycling. It's also your responsibility to decontaminate the equipment in case of biological, chemical and/or radiological contamination, so as to protect from health hazards the persons involved in the disposal and recycling of the equipment.

For more information about where you can drop off your waste of equipment, please contact your local dealer from whom you originally purchased this equipment.

By doing so, you will help to conserve natural and environmental resources and you will ensure that your equipment is recycled in a manner that protects human health.

Thank you

Austria

VWR International GmbH
Graumannsgasse 7
1150 Wien
Tel.: 01 97 002 0
Fax: 01 97 002 600
E-mail: info@at.vwr.com

Belgium

VWR International bvba
Researchpark Haasrode 2020
Geldenaaksebaan 464
3001 Leuven
Tel.: 016 385 011
Fax: 016 385 385
E-mail: customerservice@be.vwr.com

Czech Republic

VITRUM VWR s. r. o.
a VWR International Company
Pražská 442
CZ - 281 67 Stříbrná Skalice
Tel.: +420 321 570 321
Fax: +420 321 570 320
E-mail: info@cz.vwr.com

Denmark

VWR - Bie & Berntsen
Transformervej 8
2730 Herlev
Tel.: 43 86 87 88
Fax: 43 86 87 90
E-mail: info@dk.vwr.com

Finland

VWR International Oy
Valimotie 9
00380 Helsinki
Tel.: 09 80 45 51
Fax: 09 80 45 52 00
E-mail: info@fi.vwr.com

France

VWR International S.A.S.
Le Périgares – Bâtiment B
201, rue Camot
94126 Fontenay-sous-Bois cedex
Tel.: 0 825 02 30 30 (0,15 € TTC/min)
Fax: 0 825 02 30 35 (0,15 € TTC/min)
E-mail: info@fr.vwr.com

Germany

VWR International GmbH
Hilpertstraße 20a
D - 64295 Darmstadt
Freecall: 0800 702 00 07
Fax: 0180 570 22 22*
Email: info@de.vwr.com
*0,14 €/Min. aus d. dt. Festnetz

Hungary

VWR International Kft.
Simon László u. 4.
4034 Debrecen
Tel.: (52) 521-130
Fax: (52) 470-069
E-mail: info@hu.vwr.com

Ireland / Northern Ireland

VWR International Ltd /
VWR International (Northern Ireland) Ltd
Orion Business Campus
Northwest Business Park
Ballycoolin
Dublin 15
Tel.: 01 88 22 222
Fax: 01 88 22 333
E-mail: sales@ie.vwr.com

Italy

VWR International PBI S.r.l.
Via San Giusto 85
20153 Milano (MI)
Tel.: 02-3320311/02-487791
Fax: 800 152999/02-40090010
E-mail: info@it.vwr.com
info@internationalpbi.it

The Netherlands

VWR International B.V.
Postbus 8198
1005 AD Amsterdam
Tel.: 020 4808 400
Fax: 020 4808 480
E-mail: info@nl.vwr.com

Norway

VWR International AS
Haavard Martinsens vei 30
0978 Oslo
Tel.: 02290
Fax: 815 00 940
E-mail: info@no.vwr.com

Poland

VWR International Sp. z o.o.
Limbowa 5
80-175 Gdansk
Tel.: 058 32 38 200 do 204
Fax: 058 32 38 205
E-mail: labart@pl.vwr.com

Portugal

VWR International -
Material de Laboratório, Lda
Edifício Neopark
Av. Tomás Ribeiro, 43- 3 D
2790-221 Camaxide
Tel.: 21 3600 770
Fax: 21 3600 798/9
E-mail: info@pt.vwr.com

Spain

VWR International Eurolab S.L.
C/ Tecnología 5-17
A-7 Llinars Park
08450 - Llinars del Vallès
Barcelona
Tel.: 902 222 897
Fax: 902 430 657
E-mail: info@es.vwr.com

Sweden

VWR International AB
Fagerstagatan 18a
163 94 Stockholm
Tel.: 08 621 34 00
Fax: 08 621 34 66
E-mail: kundservice@se.vwr.com

Switzerland

VWR International GmbH
Lerzenstrasse 16/18
8953 Dietikon
Tel.: 044 745 13 13
Fax: 044 745 13 10
E-mail: info@ch.vwr.com

Turkey

Pro Lab Laboratuvar Teknolojileri Ltd.Şti.
a VWR International Company
Orta Mah. Cemal Gürsel Caddesi
Ördekcioglu İşmerkezi No.32/1
34896 Pendik - İstanbul
Tel.: +90216 598 2900
Fax: +90216 598 2907
Email: info@pro-lab.com.tr

UK

VWR International Ltd
Customer Service Centre
Hunter Boulevard - Magna Park
Lutterworth
Leicestershire
LE17 4XN
Tel.: 0800 22 33 44
Fax: 01455 55 85 86
E-mail: uksales@uk.vwr.com

Australia

VWR International, Pty Ltd.
Unit 1/31 Archimedes Place
Murarie, Queensland 4172
Tel.: 1300 727 696
Fax: 1300 135 123

China

VWR (Shanghai) Co., Ltd
2nd Floor, Building 4,
Lane 998, Halei Rd,
Zhangjiang Hi-tech Park
Shanghai, 201203
Tel.: +86-21-5898 6888
Fax: +86-21-5855 8801
E-mail: info_china@vwr.com

India

VWR Lab Products Private Limited
135/12, Brigade Towers, 2nd Floor
Front wing, Brigade Road,
Bengaluru, India – 560 025
Tel.: +91-80-41117125/26 (Bengaluru)
Tel.: +91-2522-647911/922 (Mumbai)
Fax: +91-80-41117120
E-mail: vwr_india@vwr.com

New Zealand

Global Science - A VWR Company
241 Bush Road
Albany 0632, Auckland
Tel.: 0800 734 100
Fax: 0800 999 002
E-mail: sales@globalscience.co.nz

Singapore

VWR Singapore Pte Ltd
18 Gul Drive
Singapore 629468
Tel.: +65 6505 0760
Fax: +65 6264 3780
E-mail: sales@sg.vwr.com

GO TO VWR.COM FOR THE LATEST
NEWS, SPECIAL OFFERS AND
DETAILS OF YOUR LOCAL VWR
DISTRIBUTOR
