VWR

Microscope  **VisiScope®** series 200

INSTRUCTION MANUAL

Model	European Catalogue Number
BL254	630-1928
TL254	630-1929
BL254PI	630-2684
TL254PI	630-2688
BL254T1	630-2673

Version: 5
Issued: 29, 02, 2016

Legal Address of Manufacturer

Europe

VWR International bvba
Researchpark Haasrode 2020
Geldenaaksebaan 464
B-3001 Leuven
+ 32 16 385011
<http://be.vwr.com>

Country of origin: **ITALY**

Table of Contents

Warning

Safety information

Package contents

Unpacking

Intended use

Symbols and conventions

Product specifications

Overview

Operation

Adjust the observation head

Place the specimen on the stage

Lamp settings

Adjust interpupillary distance (binocular models)

Focus and diopter adjustment

Condenser adjustment

Set the numerical aperture

Photo and video capture

Troubleshooting

Repair and maintenance

User replaceable accessories and spare parts

Technical service

Warranty

Compliance with local laws and regulations

Disposal

Warning

This microscope is a scientific precision instrument designed to last for many years with a minimum of maintenance. It is built to high optical and mechanical standards and to withstand daily use.

We remind you that this manual contains important information on safety and maintenance, and that it must therefore be made accessible to the instrument users.

We decline any responsibility deriving from incorrect instrument use that does not comply with this manual.

Safety Information

Avoiding Electrical Shock

Before plugging in the power supply, make sure that the supplying voltage of your region matches with the operation voltage of the equipment and that the lamp switch is in off position.

Users should observe all safety regulations of the region. The equipment has acquired the CE safety label. However, users have full responsibility to use this equipment safely.

Please follow the guidelines below, and read this manual in its entirety to ensure safe operation of the unit.

Package Contents

DESCRIPTION	QUANTITY
Microscope stand with 4 position nosepiece, stage, condenser	1
Optical head (BL254: binocular; TL254: trinocular)	1
Objective 4x	1
Objective 10x	1
Objective 40x	1
Objective 100x	1
Eyepiece WF10x/20mm	2
Immersion oil	1
Power supply 6Vdc	1

Unpacking

The microscope is housed in a moulded Styrofoam container. Remove the tape from the edge of the container and lift the top half of the container. Take care to ensure that the optical items (objectives and eyepieces) do not fall out and get damaged. Using both hands (one around the neck and one around the base), lift the microscope from the container and put it on a stable desk.

Place the observation head onto the top of the neck and tighten the lock-screw. Insert the eyepieces into the eye tubes.

Connect the provided power supply to the power supply input jack on the rear of the microscope.

Intended use

For research and teaching and teaching use only. Not intended for any animal or human therapeutic or diagnostic use.

Symbols and conventions

The following chart is an illustrated glossary of the symbols that are used in this manual.


CAUTION

This symbol indicates a potential risk and alerts you to proceed with caution

Product Specifications

Model	Head	Eyepiece	Nosepiece	Objectives	Stage	Focusing	Illuminator	Condenser	Power supply
BL254	Binocular, 360° rotating, 30° inclined.	Wide Field 10x/20mm	Quadruple, reversed	DIN E-PLAN Achromatic 4x/0.10, 10x/0.25, 40x/0.65, 100x/1.25	Double layer, 139x133 mm with X-Y 75x33mm moving range	Coaxial coarse and fine focusing	P-LED system with adjustable intensity	Sliding-in, 1.25 N.A. Abbe type with adjustable height and iris diaphragm	100-240 Vac 50-60 Hz external power adapter 6Vdc output
TL254	Trinocular, 360° rotating, 30° inclined	Wide Field 10x/20mm	Quadruple, reversed	DIN E-PLAN Achromatic 4x/0.10, 10x/0.25, 40x/0.65, 100x/1.25	Double layer, 150x139 mm with X-Y 75x33mm moving range	Coaxial coarse and fine focusing	P-LED system with adjustable intensity	Sliding-in 1.25 N.A. Abbe type with adjustable height and iris diaphragm	100-240 Vac 50-60 Hz external power adapter 6Vdc output
BL254T1	Binocular, 360° rotating, 30° inclined.	Wide Field 10x/20mm	Quadruple, reversed	DIN E-PLAN Achromatic 4x/0.10, 10x/0.25, 40x/0.65, 100x/1.25	Double layer, 150x139 mm with X-Y 75x33mm moving range	Coaxial coarse and fine focusing	P-LED system with adjustable intensity	Sliding-in, 1.25 N.A. Abbe type with adjustable height and iris diaphragm	100-240 Vac 50-60 Hz external power adapter 6Vdc output
TL254T1	Trinocular, 360° rotating, 30° inclined	Wide Field 10x/20mm	Quadruple, reversed	DIN E-PLAN Achromatic 4x/0.10, 10x/0.25, 40x/0.65, 100x/1.25	Double layer, 150x139 mm with X-Y 75x33mm moving range	Coaxial coarse and fine focusing	P-LED system with adjustable intensity	Sliding-in, 1.25 N.A. Abbe type with adjustable height and iris diaphragm	100-240 Vac 50-60 Hz external power adapter 6Vdc output
BL254PI	Binocular, 360° rotating, 30° inclined.	Wide Field 10x/20mm	Quadruple, reversed	Infinity corrected E-PLAN 4x/0.10, 10x/0.25, 40x/0.65, 100x/1.25	Double layer, 150x139 mm with X-Y moving range 75x33mm. Belt-drive in X direction	Coaxial coarse and fine focusing	P-LED system with adjustable intensity	Sliding-in 1.25 N.A. Abbe type with adjustable height and iris diaphragm	100-240 Vac 50-60 Hz external power adapter 6Vdc output

Overview

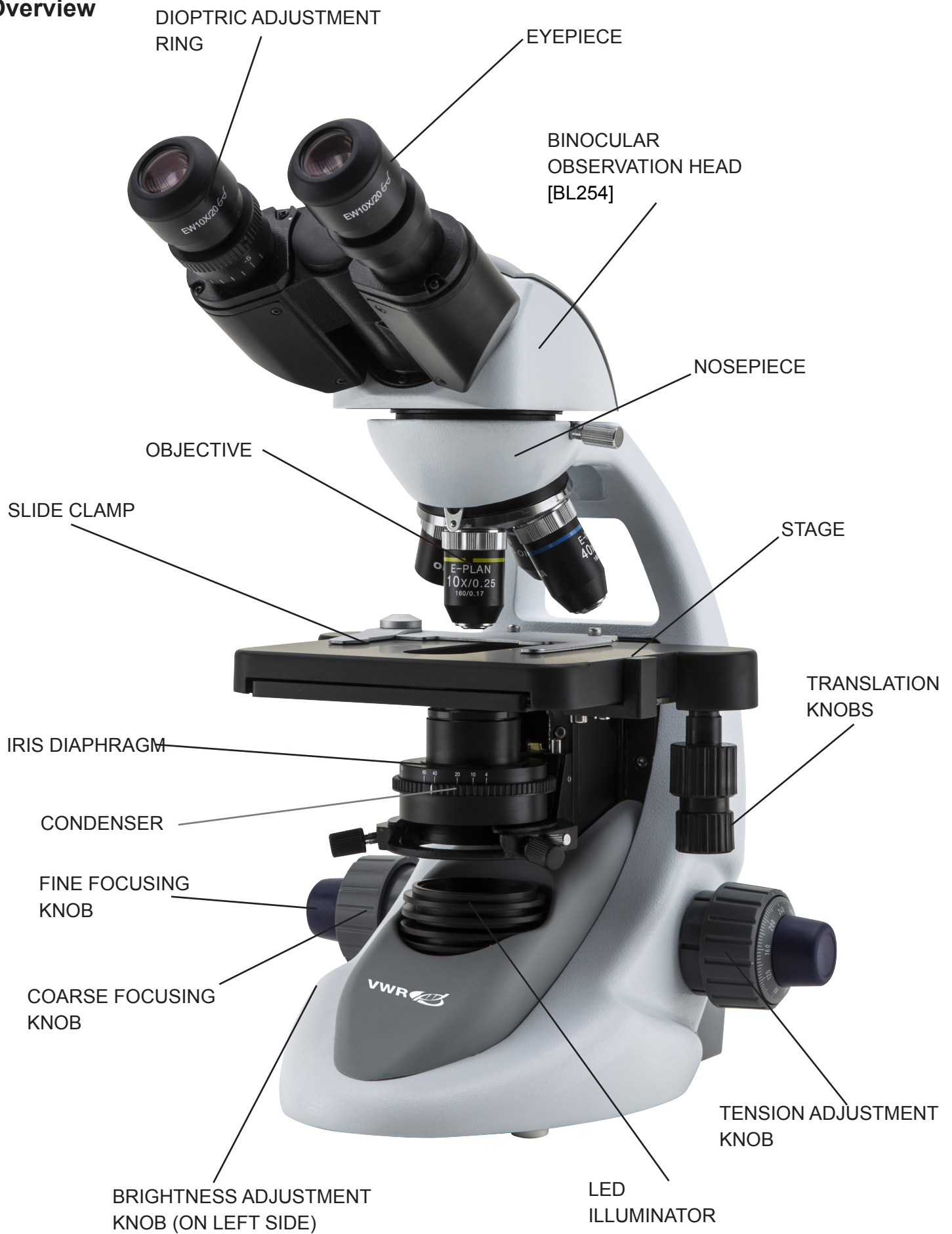TRINOCULAR
OBSERVATION
HEAD [TL254]

PHOTO PORT


Rear:

TENSION
ADJUSTMENT KNOB

ON/OFF SWITCH

POWER SUPPLY
INPUT JACK

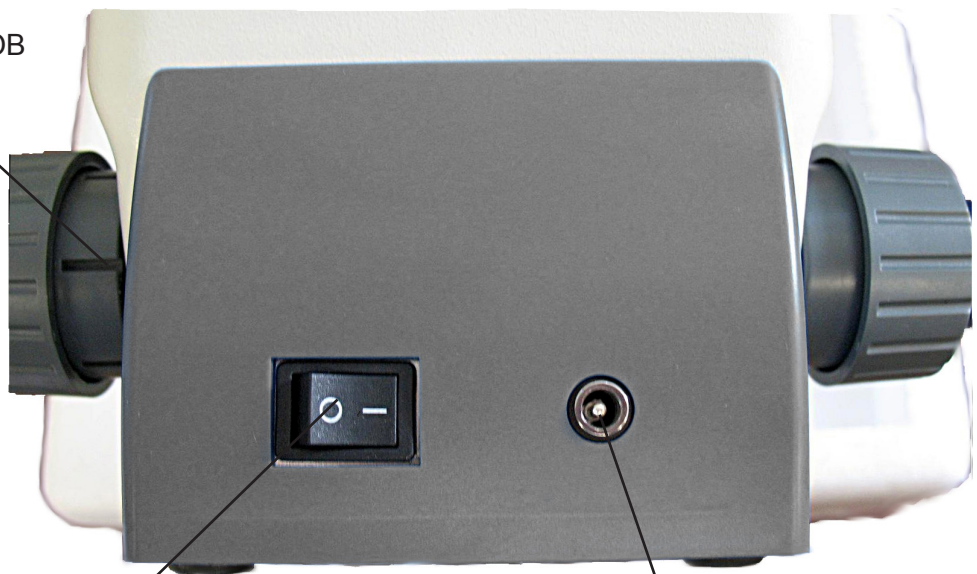Operation

Adjust the observation head

Loosen the lock-screw, turn the observation head to a comfortable position for observation, and then lock the lock-screw.

Place the specimen on the stage

Lock the specimen slide on the mechanical stage using the slide holder. Ensure that the specimen is centred over the stage opening by adjusting the coaxial stage controls.

Lamp settings

The microscope is fitted with a white LED illuminator. Before turning on the light, read the section about electrical safety precautions. Insert the plug of the power supply into the jack on the rear and turn on the ON/OFF switch. Turn the brightness adjustment knob to a brightness suitable for observation.

Adjust interpupillary distance (binocular models)

Hold the right and left side of the observation head with both hands and adjust the interpupillary distance by turning the two parts until one circle of light can be seen, when looking down the microscope with both eyes.

Focus and diopter adjustment

Rotate the coarse focusing knob to bring the slide into focus using a low magnification objective. Adjust the fine focusing knob to get the image sharp and clear while observing with your stronger eye, then turn the dioptic adjustment ring on the opposite eyepiece until a sharp and clear image can be observed also with the other eye. Turn the tension adjust knob to achieve a suitable focus tension without the specimen dropping out of focus. When the image is in focus, turn the nosepiece to choose the objective you need.

Condenser adjustment

Raise or lower the condenser through the knob to obtain a clear and uniform illumination of the sample. To center the condenser: completely close the iris diaphragm. Using the condenser centering screws, move the diaphragm in the center of the field of view. Then gradually expand the diaphragm until it is tangent to the edges of the field of view. If necessary, you can perform an additional adjustment. The condenser is centered when the edges of iris diaphragm are tangent to the field of view.

Set the numerical aperture

Adjust the aperture of the iris diaphragm to set the numerical aperture of the illuminator, thus controlling image contrast and resolution.

Photo and video capture

To use a camera, for photo or video capturing, connect it to the trinocular tube using the proper adaptor. Please refer to the camera manuals for further details.

Troubleshooting

Review the information in the table below to troubleshoot operating problems.

PROBLEM	CAUSE	SOLUTION
LIGHT DOESN'T TURN ON	Power supply not connected	Check that the power supply jack is well inserted on the rear of the microscope
	Connected	Rotate the brightness adjustment control and check if there is an increase in the light output
IMAGE CANNOT BE SEEN OR IS DARK.	The iris diaphragm aperture is not completely opened	Completely open the iris diaphragm aperture
	Brightness level is low	Rotate the brightness adjustment potentiometer
	Objective is not aligned with the optical axis	Rotate the nosepiece until an objective is well inserted in the optical path (it "clicks")
IMAGE IS UNCLEAR, BLURRED OR HAS INSUFFICIENT CONTRAST.	Objectives or filters are dirty	Wipe them clean
	The iris diaphragm aperture is not opened correctly	Open the iris diaphragm aperture completely
	Condenser at wrong height	Rotate the condenser knob until you see a uniform illumination

Repair and maintenance

Microscopy environment

This microscope is recommended for use in a clean, dry and shock free environment with a temperature of 0-40°C and a maximum relative humidity of 85 % (non condensing). Use a dehumidifier if needed.

After using the microscope

The microscope should always be kept vertical when moving it so that no moving parts, such as the eyepieces, fall out.

Never mishandle or impose unnecessary force on the microscope.

Never attempt to service the microscope yourself.

After use, turn down the illumination intensity control and turn the light off. Cover the microscope with included the dust cover, and keep it in a dry and clean place.

Electrical safety precautions

Before plugging in the power supply, make sure that the supplying voltage of your region matches with the operation voltage of the equipment and that the lamp switch is in off-position.

Users should observe all safety regulations of the region. The equipment has acquired the CE safety label. However, users do have full responsibility to use this equipment safely.

Cleaning the optics

If the optical parts require cleaning first use compressed air.

If that is not sufficient use a soft lint-free cloth with water and a mild detergent.

And as a final option use the piece of cloth moistened with a 3:7 mixture of ethanol and ether.

Note: ethanol and ether are highly flammable liquids. Do not use them near a heat source, near sparks or near electric equipment. Use these chemicals in a well ventilated room.

Remember to never wipe the surface of any optical items with your hands. Fingerprints can damage the optics. Do not disassemble objectives or eyepieces in an attempt to clean them.

For the best results, use the VWR cleaning kit (see catalogue number below).

If you need to send the microscope to manufacturer for maintenance, please use the packaging if possible .

User replaceable accessories and spare parts

DESCRIPTION	QUANTITY	CAT. NO.
Eyepiece WF10x/20mm	2	VWRI630-2000
Eyepiece WF15x	2	VWRI630-2001
Eyepiece WF20x	2	VWRI630-2002
Micrometer eyepiece WF10x/20mm	1	VWRI630-2003
Slide micrometric 26x76mm 1mm	1	VWRI630-1650
Objective 4x/0,10 E-PLAN	1	VWRI630-2004
Objective 10x/0,25 E-PLAN	1	VWRI630-2005
Objective 20x/0,40 E-PLAN	1	VWRI630-2006
Objective 40x/0,65 E-PLAN	1	VWRI630-2007
Objective 60x/0,80 E-PLAN	1	VWRI630-2008
Objective 100x/1,25 E-PLAN (Oil)	1	VWRI630-2009
Objective E PLAN IOS 4X	1	VWRI630-1639
Objective E PLAN IOS 10X	1	VWRI630-1640
Objective E PLAN IOS 20X	1	VWRI630-1641
Objective E PLAN IOS 40X	1	VWRI630-1642
Objective E PLAN IOS 100X	1	VWRI630-1643
Objective 60x/0,80 IOS E-PLAN	1	VWRI630-2049
Polarizing set (filters only)	1	VWRI630-2010
Rotating table for polarizing set	1	VWRI630-2011
Battery Solar Pack	1	VWRI630-1637
Phototube adapter for APS-C sensor	1	VWRI630-1645
Camera adapter CCD 0,5x	1	VWRI630-1644
Cleaning kit	1	VWRI630-2081

Technical service

Web Resources

Visit the VWR's website at **www.vwr.com** for:

- Complete technical service contact information
- Access to VWR's Online Catalogue, and information about accessories and related products
- Additional product information and special offers

Contact us For information or technical assistance contact your local VWR representative or visit. **www.vwr.com**.

Warranty


VWR International warrants that this product will be free from defects in material and workmanship for a period of five (5) years from date of delivery. If a defect is present, VWR will, at its option and cost, repair, replace, or refund the purchase price of this product to the customer, provided it is returned during the warranty period. This warranty does not apply if the product has been damaged by accident, abuse, misuse, or misapplication, or from ordinary wear and tear. If the required maintenance and inspection services are not performed according to the manuals and any local regulations, such warranty turns invalid, except to the extent, the defect of the product is not due to such non-performance.

Items being returned must be insured by the customer against possible damage or loss. This warranty shall be limited to the aforementioned remedies. IT IS EXPRESSLY AGREED THAT THIS WARRANTY WILL BE IN LIEU OF ALL WARRANTIES OF FITNESS AND IN LIEU OF THE WARRANTY OF MERCHANTABILITY.

Compliance with local laws and regulations

The customer is responsible for applying for and obtaining the necessary regulatory approvals or other authorizations necessary to run or use the Product in its local environment. VWR will not be held liable for any related omission or for not obtaining the required approval or authorization, unless any refusal is due to a defect of the product.

Disposal

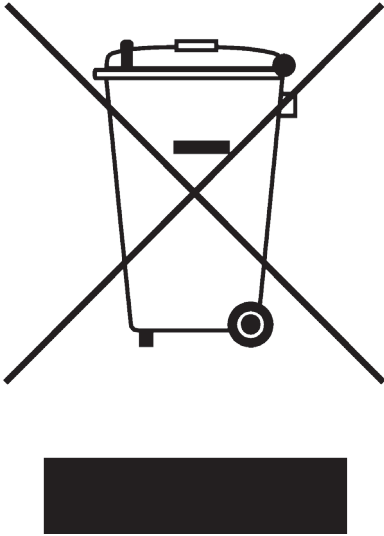This equipment is marked with the crossed out wheeled bin symbol to indicate that this equipment must not be disposed of with unsorted waste.

Instead it is your responsibility to correctly dispose of your equipment the end of its life cycle by handling it over to an authorized facility for separate collection and recycling. It is also your responsibility to decontaminate the equipment in case of biological, chemical and/or radiological contamination, so as to protect from health hazards the persons involved in the disposal and recycling of the equipment.

For more information about where you can drop off your waste equipment, please contact your local dealer from whom you originally purchased this equipment.

By doing so, you will help to conserve natural and environmental resources and you will ensure that your equipment is recycled in a manner that protects human health.

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Ireland / Northern Ireland

VWR International Ltd /
VWR International (Northern Ireland) Ltd
Orion Business Campus
Northwest Business Park
Ballycoolin
Dublin 15
Tel.: 01 88 22 222
Fax: 01 88 22 333
E-mail: sales@ie.vwr.com

Italy

VWR International PBI S.r.l.
Via San Giusto 85
20153 Milano (MI)
Tel.: 02-3320311/02-487791
Fax: 800 152999/02-40090010
E-mail: info@it.vwr.com

The Netherlands

VWR International B.V.
Postbus 8198
1005 AD Amsterdam
Tel.: 020 4808 400
Fax: 020 4808 480
E-mail: info@nl.vwr.com

New Zealand

Global Science - A VWR Company
241 Bush Road
Albany 0632, Auckland
Tel.: 0800 734 100
Fax: 0800 999 002
E-mail: sales@globalscience.co.nz

Norway

VWR International AS
Haavard Martinsens vei 30
0978 Oslo
Tel.: 02290
Fax: 815 00 940
E-mail: info@no.vwr.com

Poland

VWR International Sp. z o.o.
Limbowa 5
80-175 Gdansk
Tel.: 058 32 38 200 do 204
Fax: 058 32 38 205
E-mail: info@pl.vwr.com

Portugal

VWR International -
Material de Laboratório, Lda
Edifício Neopark
Av. Tomás Ribeiro, 43- 3 D
2790-221 Carnaxide
Tel.: 21 3600 770
Fax: 21 3600 798/9
E-mail: info@pt.vwr.com

Singapore

VWR Singapore Pte Ltd
18 Gul Drive
Singapore 629468
Tel.: +65 6505 0760
Fax: +65 6264 3780
E-mail: sales@sg.vwr.com

Spain

VWR International Eurolab S.L.
C/ Tecnología 5-17
A-7 Llinars Park
08450 - Llinars del Vallès
Barcelona
Tel.: 902 222 897
Fax: 902 430 657
E-mail: info@es.vwr.com

Sweden

VWR International AB
Fagerstagatan 18a
163 94 Stockholm
Tel.: 08 621 34 00
Fax: 08 621 34 66
E-mail: kundservice@se.vwr.com

Switzerland

VWR International GmbH
Lerzenstrasse 16/18
8953 Dietikon
Tel.: 044 745 13 13
Fax: 044 745 13 10
E-mail: info@ch.vwr.com

Turkey

Pro Lab Laboratuvar Teknolojileri Ltd.Şti.
a VWR International Company
Orta Mah. Cemal Gürsel Caddesi
Ördekcioglu İşmerkezi No.32/1
34896 Pendik - İstanbul
Tel.: +90216 598 2900
Fax: +90216 598 2907
Email: info@pro-lab.com.tr

UK

VWR International Ltd
Customer Service Centre
Hunter Boulevard - Magna Park
Lutterworth
Leicestershire
LE17 4XN
Tel.: 0800 22 33 44
Fax: 01455 55 85 86
E-mail: uksales@uk.vwr.com

**GO TO VWR.COM FOR THE
LATEST NEWS, SPECIAL OFFERS
AND DETAILS OF YOUR LOCAL
VWR DISTRIBUTOR**

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