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
Homogeniser

User Manual

en	User Manual
fr	Manuel d'utilisation
it	Manuale d'uso
es	Manual del usuario
de	Benutzerhandbuch
po	Manual do usuário

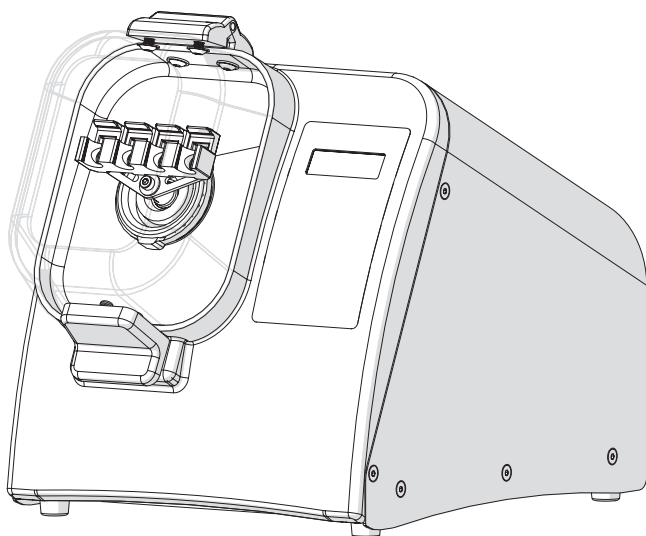Catalog No	Description
432-0366	Mini Bead Mill 230V, 300W

Legal Address of Manufacturer

United States Europe

VWR International, LLC
100 Matsonford Rd
Radnor, PA 19087
800-932-5000
<http://www.vwr.com>

VWR International bv
Researchpark Haasrode 2020
Geldenaaksebaan 464
3001 Leuven
+ 32 16 385011
<http://be.vwr.com>

Country of Origin

Made in USA

Intended Use

The VWR Mini Bead Mill is intended for milling and homogenising of a number of substances. Not intended for medical or diagnostics purposes.

Technical Service

Web Resources

Visit the VWR website at www.vwr.com for:

- Complete technical service contact information
- Access to the VWR Online Catalogue, and information about accessories and related products
- Additional product information and special offers

Contact Us

For information or technical assistance, contact your local VWR representative or visit www.vwr.com.

Warranty

VWR warrants that this product will be free from defects in material and workmanship for a period of one year from date of delivery. If a defect is present, VWR will, at its option and cost, repair, replace, or refund the purchase price of this product to the customer, provided it is returned during the warranty period. This warranty does not apply if the product has been damaged by accident, abuse, misuse, or misapplication, or from ordinary wear and tear. If the required maintenance and inspection services are not performed according to the manuals and any local regulations, such warranty turns invalid, except to the extent, the defect of the product is not due to such non-performance.

Items being returned must be insured by the customer against possible damage or loss. This warranty shall be limited to the aforementioned remedies. IT IS EXPRESSLY AGREED THAT THIS WARRANTY WILL BE IN LIEU OF ALL WARRANTIES OF FITNESS AND IN LIEU OF THE WARRANTY OF MERCHANTABILITY.

Compliance with Local Laws and Regulations

The customer is responsible for applying for and obtaining the necessary regulatory approvals or other authorizations necessary to run or use the Product in its local environment. VWR will not be held liable for any related omission or for not obtaining the required approval or authorization, unless any refusal is due to a defect of the product.

Equipment Disposal

This equipment is marked with the crossed out wheeled bin symbol to indicate that this equipment must not be disposed of with unsorted waste. Instead it is your responsibility to correctly dispose of your equipment at life cycle - end by handing it over to an authorized facility for separate collection and recycling. It is also your responsibility to decontaminate the equipment in case of biological, chemical and/or radiological contamination, so as to protect from health hazards the persons involved in the disposal and recycling of the equipment. For more information about where you can drop off your waste of equipment, please contact your local dealer from whom you originally purchased this equipment. By doing so, you will help to conserve natural and environmental resources and you will ensure that your equipment is recycled in a manner that protects human health.


Thank you.

Table of Contents

- 4** Unpacking
- 5** Overview
- 6** General Operation
- 7** Loading Sample Tubes
- 8** Changing the Tube Holder
Changing the Neck Gasket
- 9** Troubleshooting
- 10** Transporting the Mini Bead Mill


Unpacking

1.

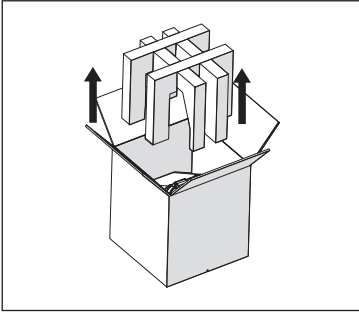
Lift foam from box.

2.

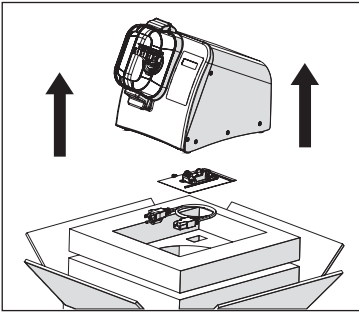
Lift unit up and off of the box bottom.
DO NOT LIFT THE UNIT BY HOLDING THE LID.
Remove Tool Kit and Power Cord

3.

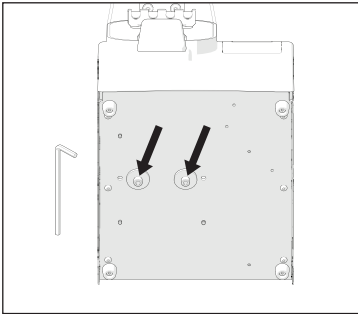Remove the 2 screws from the bottom of the unit.



4.

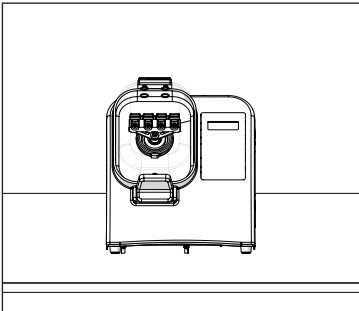Place the unit on a stable surface.

5. **Connect power cord**

Overview

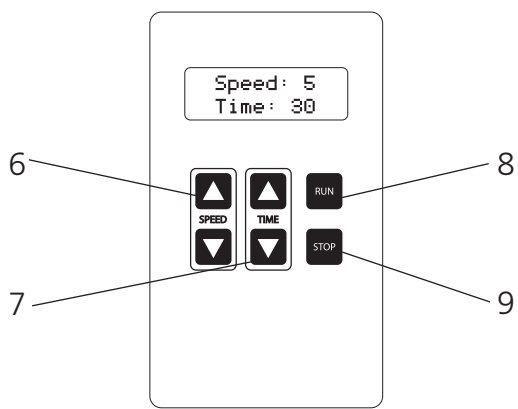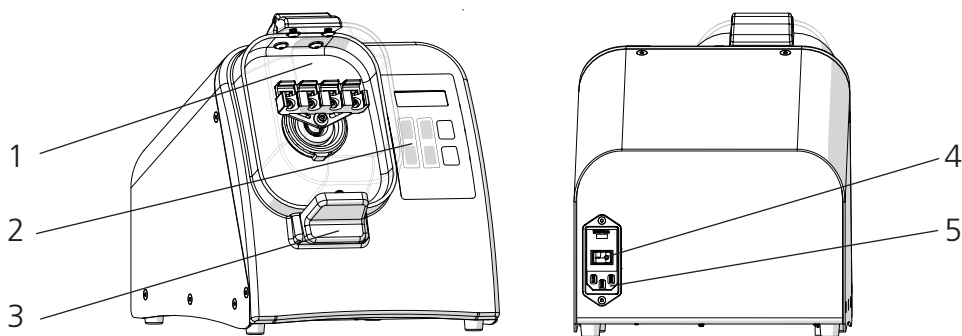
- 1. Lid
- 2. Keypad
- 3. Lid Latch
- 4. ON/Off Switch
- 5. Power Cord Port
- 6. Speed Control
- 7. Timer Control
- 8. Run Button
- 9. Stop Button

The Mini Bead Mill consists of the following:

Description	Quantity
Mini Bead Mill	1
Power Cord	1
Spare Fuse	1
7 mL Tube Holder Assembly	1
2 mL Tube Holder Assembly	1
Allen Key	1
User Manual	1

General Operation

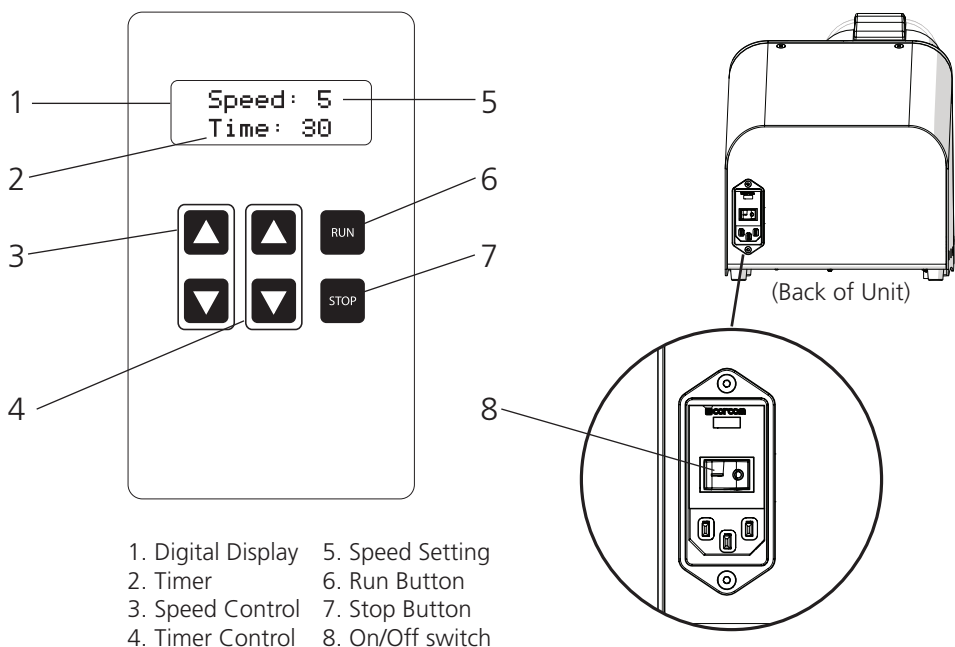1. Power

- Plug in the power cord and move the On/Off switch on back of the unit to the On position.

2. Set Speed

- Adjust speed using up and down buttons.

3. Set Timer

- Adjust run time using up and down buttons from 1 sec. to 5 min.

4. Start

- Press Run button.



5. Pause/Stop

- Press Stop Button to Pause.


- Press Stop again to reset timer and cancel run

Speed Setting	Processing Power in G's	Processing Speed in m/s
1	30	1.5
2	60	2.5
3	90	3.7
4	120	4.4
5	150	5.0


Loading Sample Tubes

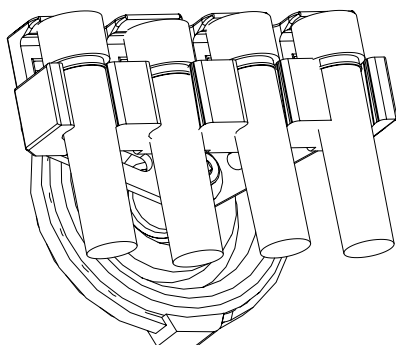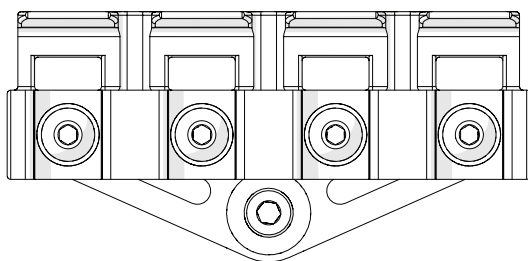
1.5 mL & 2 mL Tubes

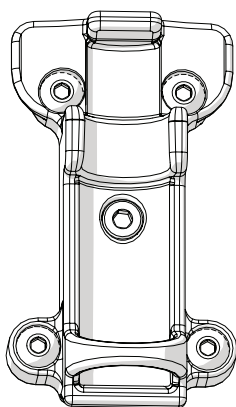
7 mL Tubes

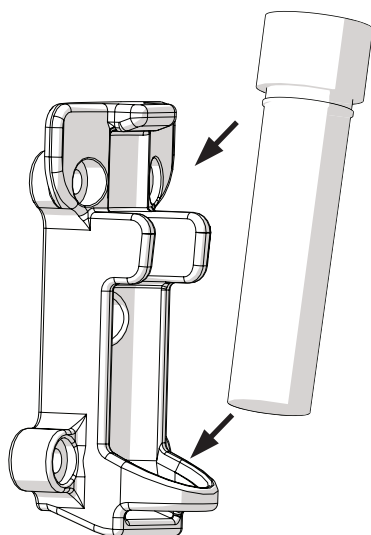
Changing the Tube Holder

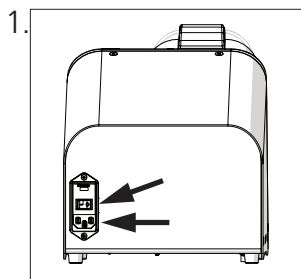Caution: Turn unit off and unplug power cord.

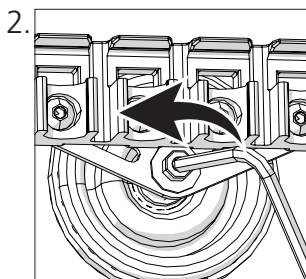
Use the supplied 3 mm Allen wrench to loosen tube carriage holder screw

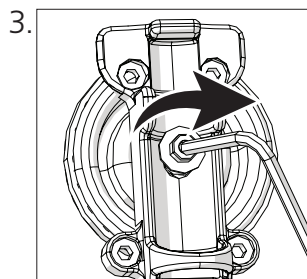
Use the supplied 3 mm Allen wrench to tighten tube carriage holder screw


Changing the Neck Gasket

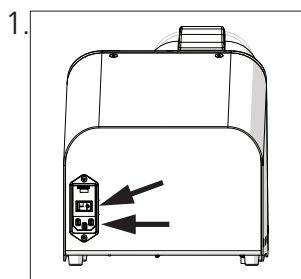
Caution: Turn unit off and unplug the power cord.

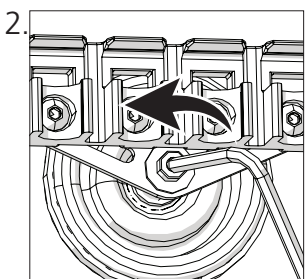
Use the supplied 3 mm Allen wrench to loosen the central screw and remove restriction tie and old neck gasket.

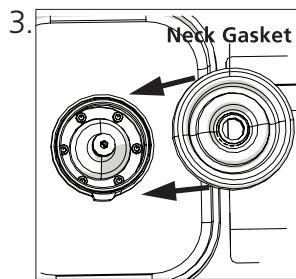Slide the new neck gasket onto the pivot shaft and housing

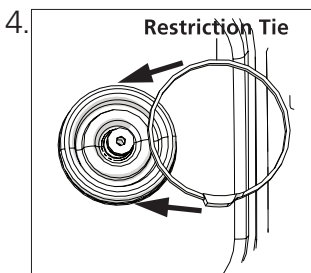Secure the new neck gasket using the restriction tie with the clamp positioned toward the bottom of the instrument. Tighten central screw.

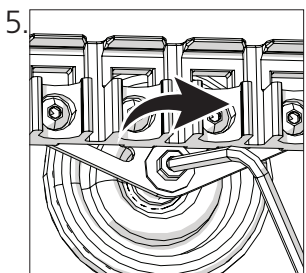Use the supplied 3 mm Allen wrench to tighten the central screw and secure the tube carriage.


Troubleshooting


Error Codes

Error Code	Definition	Cause	Possible Solution
Error 1	Display screen shows incomplete or unusual characters	CPU needs to be reset	Turn off the unit. Wait 30 seconds. Turn the unit back on.
Error 2	No display on the screen	No Power	Check main voltage. Check that the voltage of the unit (100V-240V) matches what is being delivered by the main power supply. Make sure the unit is plugged in properly. Check and replace fuses if blown.
		Faulty display screen	Turn off the unit. Contact manufacturer for assistance.
Error 3	Unit powers on but motor does not turn	Faulty motor	Contact manufacturer for assistance.
Error 4	Close Door	Lid is not seated properly.	Check Lid latch for obstruction Press lid and ensure the handle is closed
		Detection system is faulty.	

Transporting the Mini Bead Mill

**The Mini Bead Mill must be shipped in its original packaging.
Before transporting please follow these steps.**

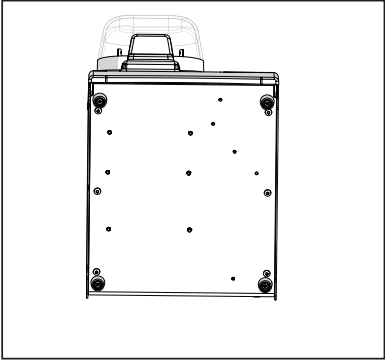1. 

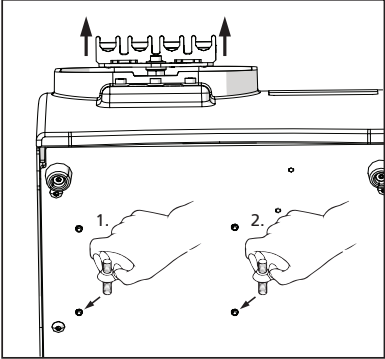
Remove all tubes, and cover the unit in its plastic bag.
2. 

The tube carriage will slide back. Lift the tube carriage to align the motor mount holes to the holes in the base plate.

Insert screws with Allen key until you feel them engage (It may take some adjusting to get them aligned with the threads in the unit).

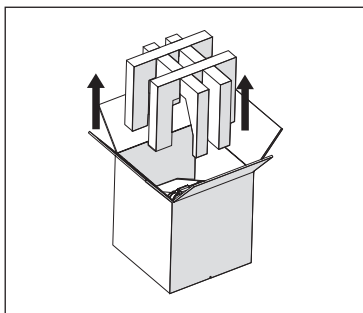Tighten screws until completely seated.

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VWR International AB
VWR International AG
VWR International Ltd