

Instruction manual

VWR® PREMIUM CONTROL PC SOFTWARE

EU cat. no: 390 - 1323



Legal address of Manufacturer

Europe

VWR International bv
Researchpark Haasrode 2020
Geldenaaksebaan 464 B-3001 Leuven
+3216385011
be.vwr.com

Country of Origin

POLAND

Table of Contents

1. Main description	4
2. Hardware requirements.....	5
3. Program installation.....	5
3.1 Ports used.....	7
4. First start-up	7
4.1 Application boot	7
4.2 Logging in	7
5. Operating the application	8
5.1 Main Menu	8
5.2 List of units.....	10
5.3  Status	10
5.4  Programs.....	11
5.4.1. Preview of the program list	12
5.4.2. Creating / edition of program	13
5.4.3. Creating an offline program	14
5.5  Graph.....	15
5.6  Import of PLKX file	15
5.7  Register data	16
5.8  Statistics.....	17
5.9  Manage the users	18
5.9.1. User creation / edition	18
5.9.2. Assigning units	19
5.10  Manage units	20
5.10.1. Adding a unit	20
5.10.2. Connecting to the unit	21
5.11 Settings and information.....	21
5.11.1. Settings	22
5.12 Sterilisation (applies only to CO ₂ incubators)	23

1. Main description

Premium Control Software was designed to meet customers' needs to facilitate operation and control the VWR PREMIUM models. The application is the integral part of PREMIUM series units. It can be downloaded from unit's internal memory to a USB stick.

You may purchase the Premium Control in SET (a program on a USB memory stick, a special key- the licence for Premium Control, and a LAN cable) as additional equipment for:

Description	Cat. No.
IL 56 Prime	390-1346
IL 112 Prime	390-1349
IL 180 Prime	390-1350
DL 56 Prime	466-0600
DL 112 Prime	466-0601
DL 180 Prime	466-0602
VL 56 Prime	466-0603
VL 112 Prime	466-0604
VL 180 Prime	466-0605
IL 68R	390-1341
IL 150R	390-1342
IL 250R	390-1343



Two-way (to the unit and from the unit) **communication is possible only for PREMIUM models**. You can send to and download from the unit. **For the other models you can only download data from the unit.**

Application functionalities:

- simultaneously connect several PREMIUM units
- control units remotely
- overview of running program status
- save real-time running program data to the .plx¹⁾ file
- alarm information
- option to create programs and upload them remotely
- modify existing programs
- create programs offline
- start / stop programs
- set a delayed start for a program
- open registry data file / events downloaded from the unit
- **generate reports** /graphs from registry or events data file
- overview of current data statistics
- generate reports from current statistics
- download statistics of the launched program
- download registered data / events
- produce charts
- user management panel
- change time zone
- change temperature correction
- set alarms

¹⁾ plx. file is a file format saved by the Premium Control application or saved to a USB memory stick (pendrive) when downloading data/events on from the unit. These files may contain registration data (date/time, temperature) or event data (date/time, event name, type of event, code). It is not possible to open these files in another application.

2. Hardware requirements

Below there is the specification of the computer on which the Premium Control program can be installed:

- Windows 8, 10 operating system,
- min 300 MB free space on the hard disc,
- minimum screen resolution: HDTV 720p 1280×720 32-bit,
- LAN network card,
- NET Framework 4.6.1 installed (you can download it on www.microsoft.com),
- Microsoft SQL Server Compact 3.5 SP2 ENU installed,
- Microsoft SQL Server Compact 4.0 installed.

3. Program installation

The installation program can be downloaded to a USB memory stick from the unit in the PREMIUM version. To do so:


- insert a USB memory stick into the slot,
- in main menu go to **info** and press the icon 
- .



A folder 'Download' will be saved onto the USB flash drive in which you can find the following folders:


- Premium Control – includes an installation program and the instruction manual of Premium Control application
- Manual – contains a manual of the unit

From the Download\Premium Control folder run the file setup.exe.

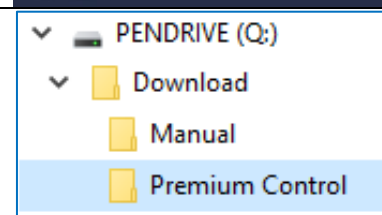


During the installation of Premium Control two additional applications: SQL Server CE 4.0 and SQL Server CE 3.5 will be automatically installed. They are necessary for a correct operation of the program.


Choose the target location of the program using **Browse** button or confirm the default one by pressing **Next**. If the target folder does not exist, the program will ask if you want to create it. Press the button: **Yes**.

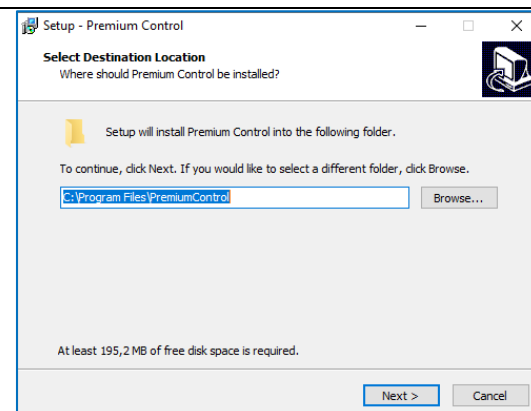


In this window confirm the installation by pressing **Next**.

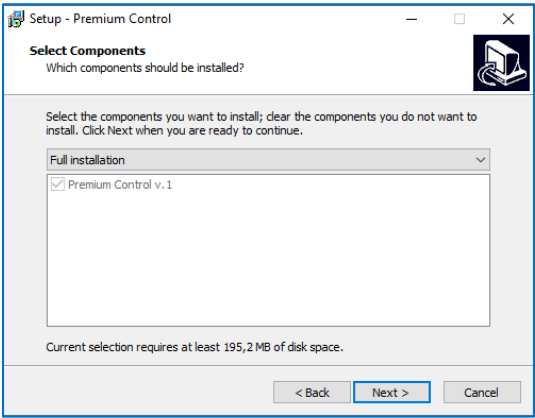


Choose Menu Start folder, in which to create shortcuts to the Premium Control application. You can skip creating shortcuts in the Start Menu by selecting the "Don't create a Start Menu folder" option at the bottom of the installation window. After finding the folder or accepting the default folder, click **Next**.

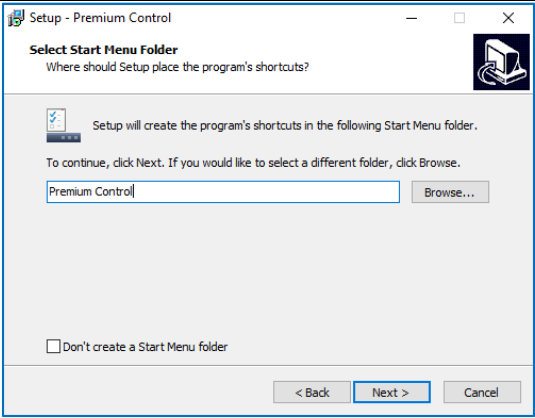


You can create an additional icon on your desktop by selecting (recommended) "Create a desktop shortcut". Click Next to continue.

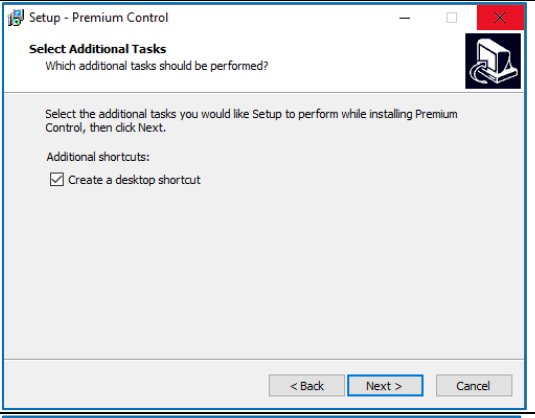


Confirm the start of installation by pressing **Install**.

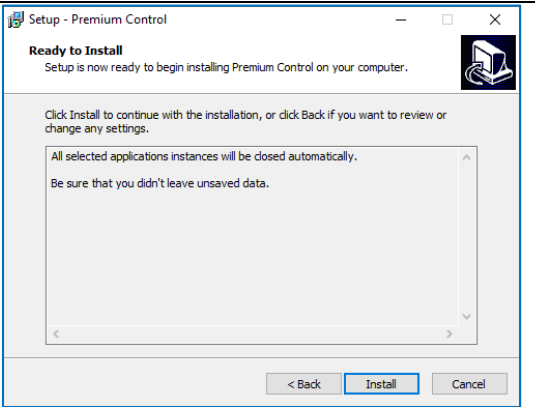

The installation progress is shown graphically.

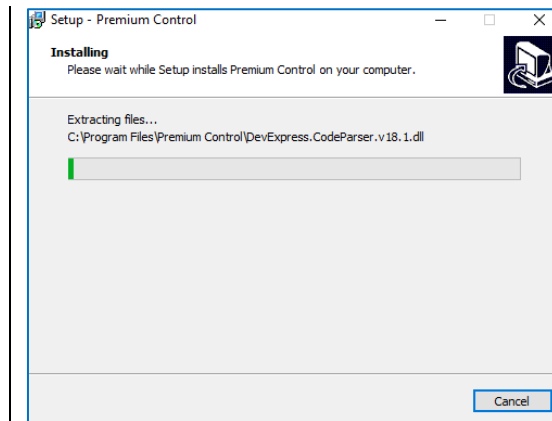

If the installation process run correctly, a summary window will be displayed. Press **Finish**.

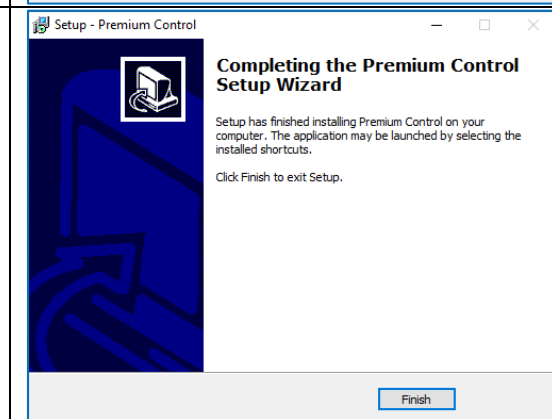

3.1 Ports used

The port used to communicate the Premium Control application with VWR units is : **56506**

4. First start-up

4.1 Application boot

- If **Create a desktop shortcut** was chosen during the installation process, the following shortcut will display on the desktop:




click on it to run the installed program.

- If the option "Do not create a folder in the Start Menu" was not selected, the program can be started by clicking Start→Programs→Premium Control.
- The program can also be started by clicking on the icon shown above and located in the installation folder. Default: "C: \ Program Files \ Premium Control.

4.2 Logging in

Working with the application is possible only by the logged-in user. To do this, click on the "Login" button.

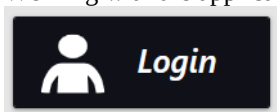

Then, enter a name of a user and a password in the login window (Figure 1). During the first start-up you should enter:

login: admin

password: admin

It is recommended to change the password as soon as possible. It should contain uppercase and lowercase, symbols and numbers.
Information about changing the password is on [page 18](#).



Enter your user data

Username:

Password

Login

Figure 1 – Logging in

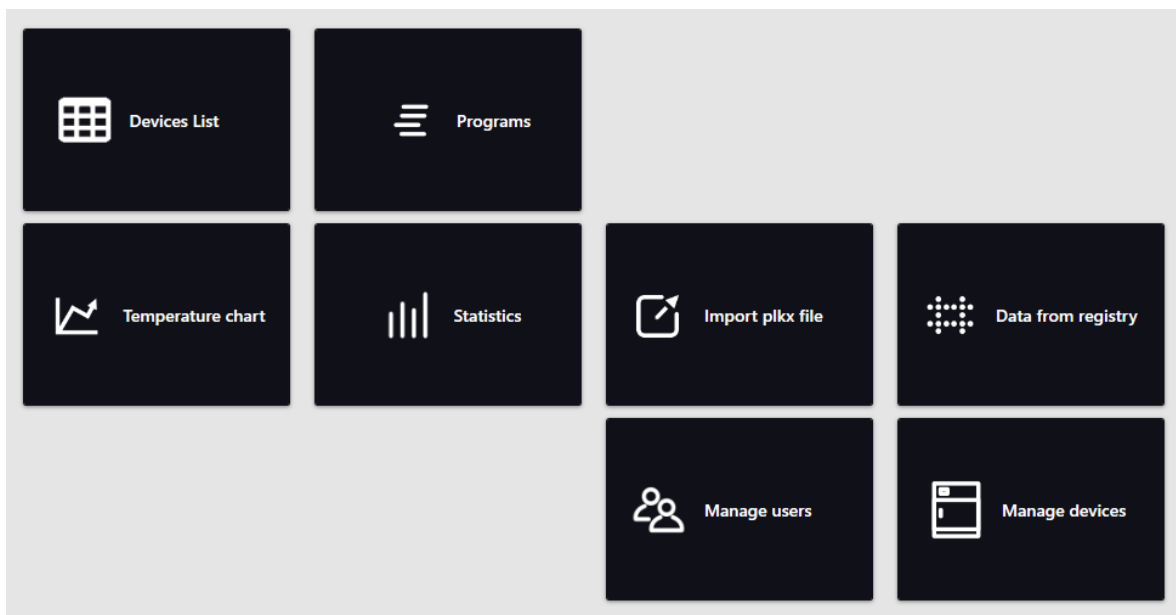


Figure 2 – Main menu

5. Operating the application

5.1 Main Menu

After the user has been successfully logged in, the main application menu will be displayed ([Figure 2](#)). The available bars depend on the rights of the logged in user.

All links to the program's function panels are available in the main menu. The same links are available in the sidebar ([Figure 3](#)).

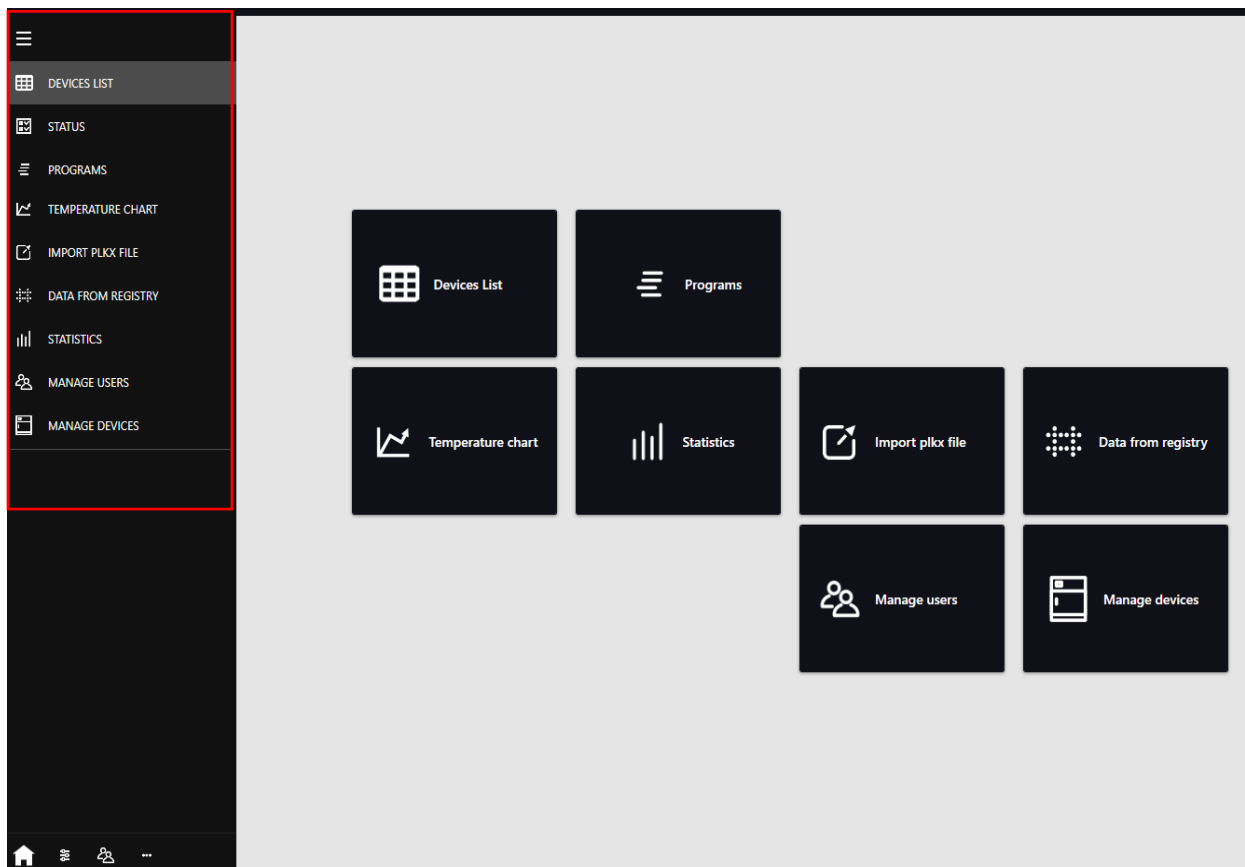


Figure 3 – Extended sidebar

The sidebar (Figure 3) enables quick moving between panels. It is always available in every panel. The extended list of panels is given below:

	Extended sidebar		Statistics
	List of units		Manage users
	Status		Manage devices
	Programs		Main menu
	Temperature chart		Settings
	Import of .plkx file		Logged in user
	Data from registry		

The panels are described in details further in the manual.

5.2 List of units


This panel displays units available to the logged-in user.

Save Layout		Restore Layout					
*	Serial number	Device name	Temp[1]	Program	Door	Alarm	Save to file
	IW1TF 197777	IW1TF 197777	25,7°C				
	IW1TF 196666	IW1TF 196666	21,4°C				

Figure 4 – List of units


Available functions:

- **Save the diagram** – saves the current column layout
- **Restore layout** – restores default column layout
- **Save to file** – clicking the icon starts recording on a regular basis to a file
- **Concise header menu** - available by right-clicking on the table header. Enables to quickly move to a chosen panel.
- **Concise device menu** - available after right-clicking on the unit's row. It allows you to move quickly to the selected panel.

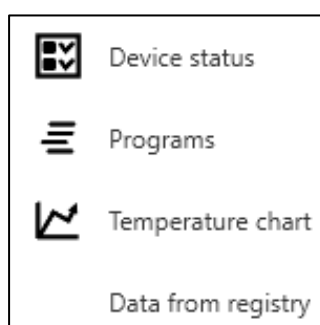

5.3 Status

The "status" panel is used to monitor the currently running program and the status of the equipment. It consists of independent windows that can be moved and closed. To save the layout of windows, save it by clicking on **Save diagram**. The Restore Layout button restores the default layout of the windows. For other units (not Premium) a key is required.

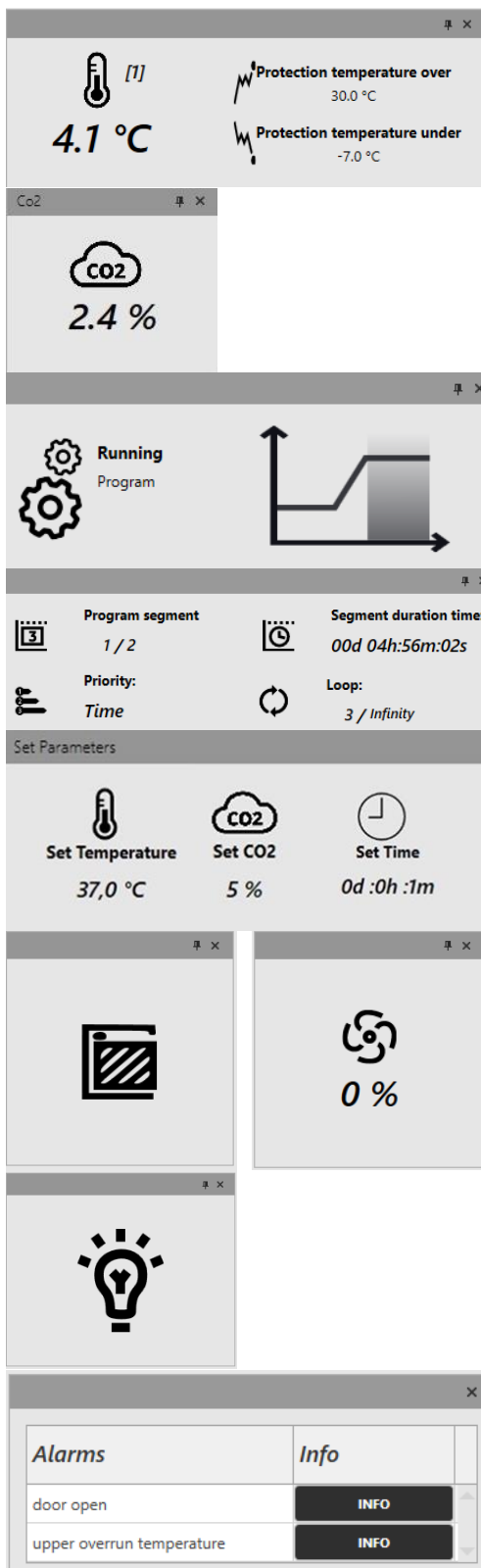

Current temperature,
Over and under temperature protection

Current CO₂ concentration

Program status: works or is turned off
Program stage: heating, cooling or maintaining the temperature

The program segment currently being executed, the selected priority, the duration of the segment, the number of the loop being executed

Set temperature, CO₂ concentration, time

Door status: open or closed
Fan performance

Indicates that the internal lightning is on

List of active alarms.
button „**INFORMATION**” – displays a preview of all instances of the selected alarm

5.4 Programs

In this panel (Figure 5) you can browse programs, add a new one or edit an existing program. The top menu is available:



D. START – planning the program launch according to set date and time.

START – starts the chosen program.

STOP – stops the program.

CREATE PROGRAM – launches the new program creator.

EDIT PROGRAMU – edits the chosen program.

PROGRAM OFFLINE – enables to create a program without connecting to the unit (more information [page 14](#)).

DELETE PROGRAM – deletes the selected program.

Available only for Premium units. For other equipment a key is required and only a program preview is available.

5.4.1. Preview of the program list

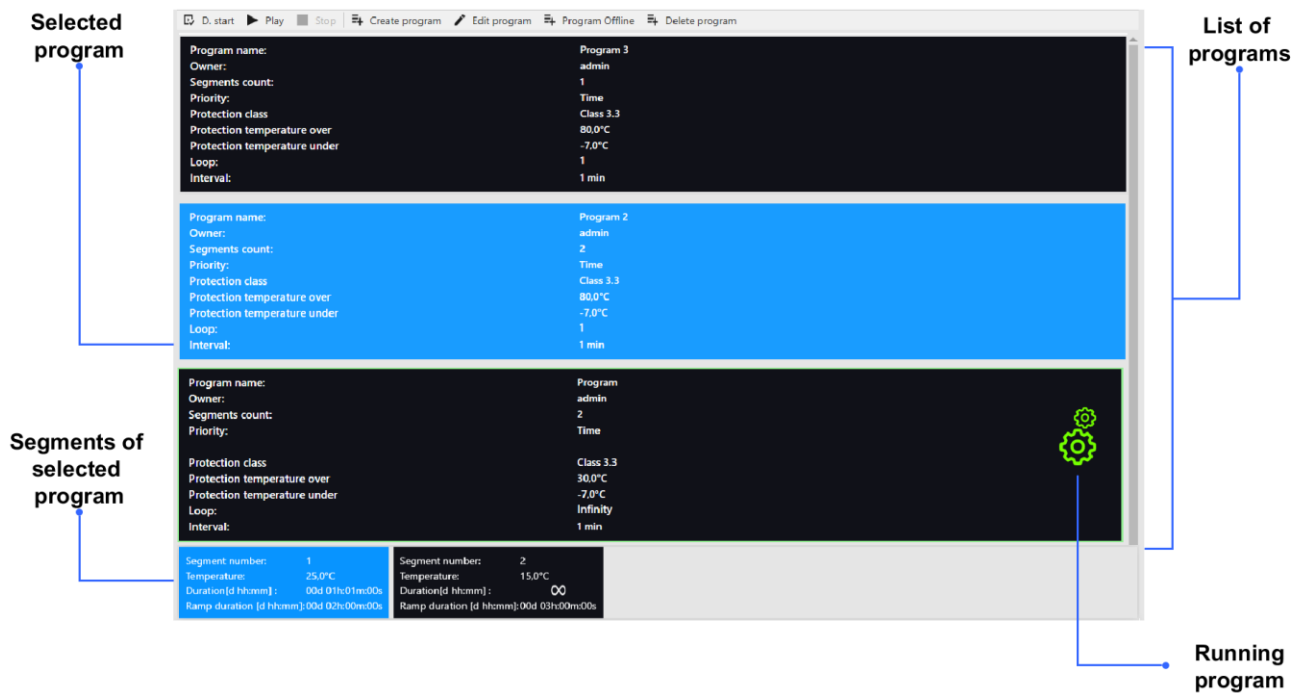


Figure 5 – Program list

5.4.2. Creating / edition of program



A detailed description of the set parameters in the program can be found in the manual of PREMIUM equipment.

Available only for PREMIUM units.

Click „add a program” or „program edition”, and program creator will be displayed.

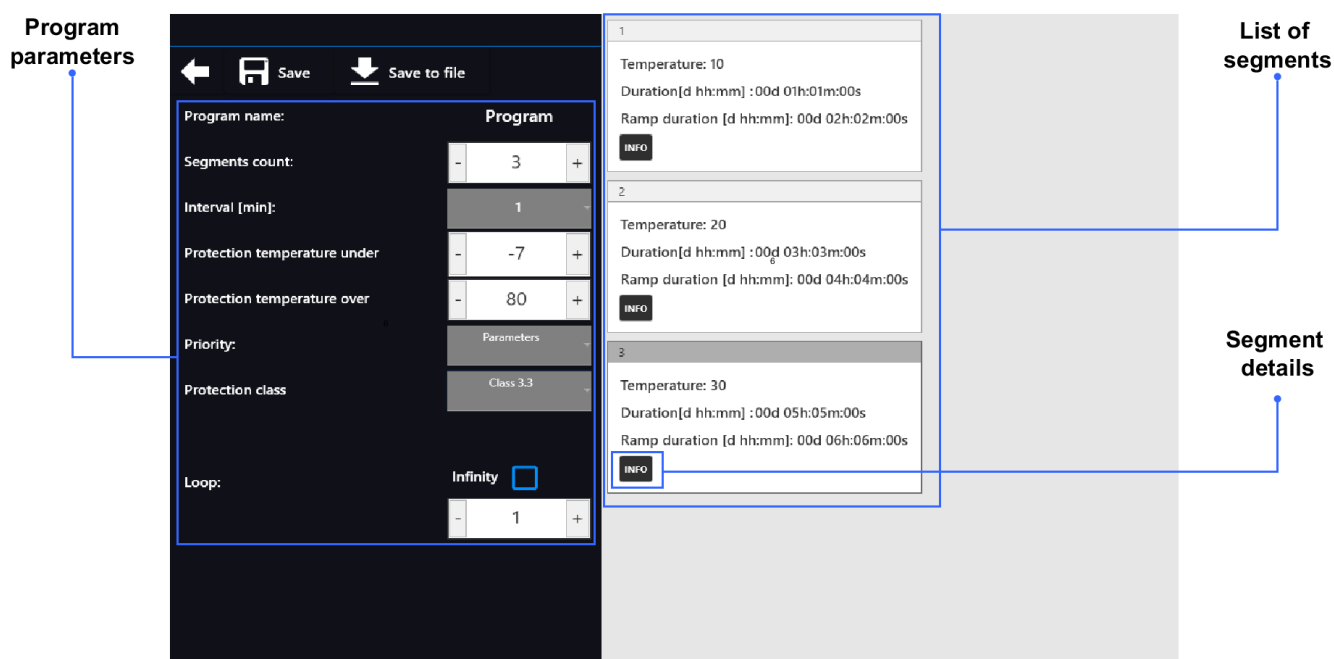


Figure 6 - Adding / edition of program

In the left part of the window, in the program parameters, you can set:

- The name of a program
- Number of segments
- **Interval [min]** – frequency of data saving to the data register
- Under temperature protection
- Over temperature protection
- **Priority** – time or parameters priority, more information
- **Protection mode** – the security class of the program
- **Loop** – the number of program repetitions, more information.

Once all the above program parameters have been set, go to the List of segments in the right part of the window where you can edit individual segments.

By clicking **INFO** you can display / edit the detailed information about the segment (Figure 7):

- **temperature** – set temperature which the unit is to achieve in the segment (cannot be lower than the protection temperature (under temperature) +2°C and higher than the temperature of the over-temperature -2°C),
- **duration of the segment** – time of maintaining the set temperature ([d hh:mm]) in days, hours and minutes; in the last segment you can choose continuous work,
- **ramp duration** – time of reaching the set temperature ([d hh:mm]) marked in days, hours and minutes.

The following parameters are available on selected models or as options available on request:

- **fan** – fan efficiency in %,
- **ramp fan** – efficiency of the fan while maintaining the set temperature,
- **CO₂** – set CO₂ concentration which the unit is to achieve in the segment,
- **ventilation pipe** – the level of ventilation pipe opening,
- **ramp pipe** – the degree of ventilation pipe opening during maintaining the set temperature.



- saves the introduced changes in the program, available after passing through all segments.

Figure 7 – Segment edition

5.4.3. Creating an offline program

Available only for PREMIUM units.

The offline program gives opportunity to create a program without connecting the unit to the network or a computer. You need a file with a configuration of the equipment which the program is created for. It can be created by following steps below:

- 1) Insert the USB stick to the selected unit.
- 2) Log in.
- 3) Go to the „Info” panel.
- 4) Confirm the window asking if you want to save the unit ID to the data carrier, by pressing .
- 5) Remove the USB stick and connect it to the computer with the Premium Control application installed.
- 6) Go to the „Programs” panel and click „Offline Program”.
- 7) In the file opening window select the USB stick and browse the file saved on/downloaded from the equipment. The file has a type, model and serial number in its name.
- 8) Create a program in the creator.
- 9) Click „Save to file” – program will be saved in the location specified in „Settings” in the section „Storage of unit records”.
- 10) You can create several programs.
- 11) Transfer the saved file to the USB stick.
- 12) Remove the USB stick and insert to the equipment.
- 13) In the equipment enter „Main Menu” and the “Programs” panel.
- 14) Click the icon and then .
- 15) Select the saved program from the list.
- 16) The program will appear on the list.



The saved program file is only visible in the equipment based on which the configuration was created.

5.5 Graph

For other units (not Premium) a key is required.

In this panel you have a preview of the temperature and CO₂ concentration that prevailed in the unit at a given time in the form of a graph. In this panel, a graph is drawn from the data currently downloaded from the unit. Drawing the graph starts when connected to the unit and continues even if the graph panel is not displayed. Presented series of the data depend on the type and equipment of the unit with which we are connected.

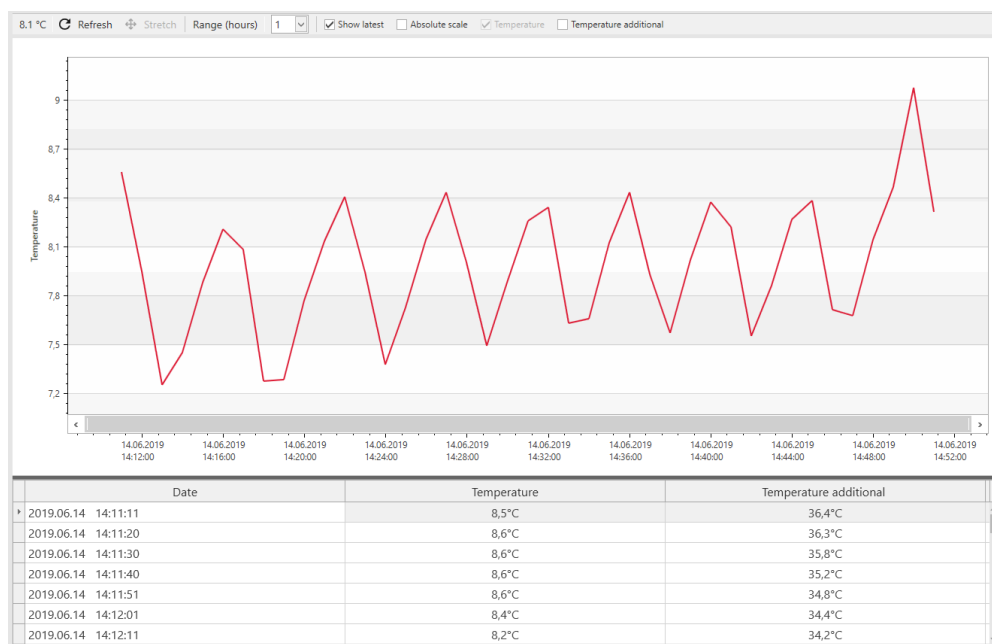


Figure 8 – Temperature graph

In the upper left corner the current temperature from the main sensor is displayed.

This panel can be refreshed on a regular basis using the "**refresh**" option on the top taskbar.

The data for the chart is saved every defined period of time (in this case, 5 seconds), however, you can turn on or off the chart refreshing with the "**show latest**" button. If the button is blue, it means that the chart is updated on a regular basis. However, if it is white then the graph will be in place.

The "**stretch**" option displays the entire graph. Available if the option to update data on the chart is off.

Selecting the "**absolute scale**" scales the y axis so that it starts from 0 °C.

The "**Temperature**", "**Additional sensor**", "**CO₂**" options display or hide the data series on the chart.


5.6 Import of PLKX file

For other units (not Premium) a key is required.

Panel (**Błąd! Nie można odnaleźć źródła odwołania.**) enables the opening of *.plx file which contains data registered in the unit:

- data: date and time, temperature of the main sensor, optionally the temperature of the additional sensor, CO₂ concentration
- events: date and time, event description, event priority, service code.

The following main menu is available:

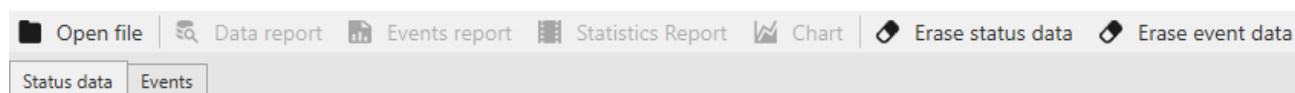

Open file – opens a file with .plx extension.

Data report – opens a window with data report.

Events report – opens a window with events report.

Statistics report – opens a window with a report summarizing each segment of the program. It contains information about registration time, the number of samples, the registered minimum and maximum value and the average temperature value.

Chart – opens a window with a chart – available only for data.

Erase status data – closes a file with data.

Erase event data – closes a file with events.

Status data and Events – buttons that enable to switch between open files with data and events.

Saving the PLKX file on the unit:

- 1) connect a USB stick directly to the unit.
- 2) enter the main menu and then to „**data register**” panel where you can download the registered data or to „**event register**” panel where you can download register data on the unit.
- 3) Press the icon  choose *.plx type and confirm.
- 4) Connect the USB stick to the computer and in Premium Control application click „**open file**” in the top left corner.

Date	Temperature	Program phase	CO2
2021-07-19 05:21:47	24,5°C	set temp.	4,5%
2021-07-19 05:22:47	24,5°C	set temp.	4,5%
2021-07-19 05:23:47	24,5°C	set temp.	4,5%
2021-07-19 05:24:47	24,5°C	set temp.	4,5%
2021-07-19 05:25:47	24,5°C	set temp.	4,5%
2021-07-19 05:26:47	24,5°C	set temp.	4,5%
2021-07-19 05:27:47	24,5°C	set temp.	4,5%
2021-07-19 05:28:47	24,5°C	set temp.	4,5%
2021-07-19 05:29:47	24,5°C	set temp.	4,5%
2021-07-19 05:30:47	24,5°C	set temp.	4,5%
2021-07-19 05:31:47	24,5°C	set temp.	4,5%
2021-07-19 05:32:47	24,5°C	set temp.	4,5%
2021-07-19 05:33:47	24,5°C	set temp.	4,5%
2021-07-19 05:34:47	24,5°C	set temp.	4,5%
2021-07-19 05:35:47	24,5°C	set temp.	4,5%
2021-07-19 05:36:47	24,5°C	set temp.	4,5%
2021-07-19 05:37:47	24,5°C	set temp.	4,5%
2021-07-19 05:38:47	24,5°C	set temp.	4,5%
2021-07-19 05:39:47	24,5°C	set temp.	4,5%
2021-07-19 05:40:47	24,5°C	set temp.	4,5%


Figure 9 – Import of PLKX file

5.7 Register data

For other units (not Premium) a key is required.

Panel (Figure 9) enables to download registered data and events from the unit.

The following menu is available:

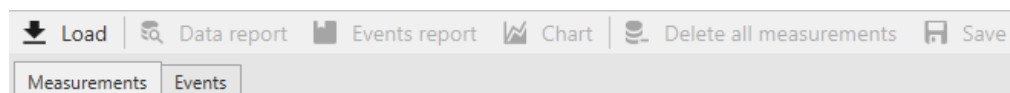

Load – downloads registered data and events from the unit.

Data report – opens a window with data report.

Events report – opens a window with events report.

Chart – opens a window with a graph – available only for data.

Delete all measurements – deletes all the measurements on the unit.

Delete data register – deletes all the events on the unit.

Save – saves the downloaded data in the location specified in the settings.

To download data/events in the right panel with the list of units, select the unit and click Connect. Then click **"LOAD"**, the recorded data and events will be downloaded.

Clicking one of the DELETE options erases the data on the unit irretrievably. Therefore, it is recommended to download and save data and events before this.

Measurements		Events		
	Date °	Temperature °	Program phase °	Additional °
↑				
▶	2019-06-14 14:59:24	7,5°C	set temp.	29,9°C
	2019-06-14 14:58:24	8,3°C	set temp.	30,0°C
	2019-06-14 14:57:24	8,4°C	set temp.	29,7°C
	2019-06-14 14:56:24	8,0°C	set temp.	27,8°C
	2019-06-14 14:55:24	7,6°C	set temp.	29,2°C
	2019-06-14 14:54:24	8,2°C	set temp.	29,5°C
	2019-06-14 14:53:24	8,4°C	set temp.	28,8°C
	2019-06-14 14:52:24	8,0°C	set temp.	27,0°C
	2019-06-14 14:51:24	7,6°C	set temp.	28,6°C
	2019-06-14 14:50:24	8,1°C	set temp.	30,2°C
	2019-06-14 14:49:24	9,1°C	set temp.	30,2°C
	2019-06-14 14:48:24	8,5°C	set temp.	29,8°C
	2019-06-14 14:47:24	8,1°C	set temp.	28,2°C
	2019-06-14 14:46:24	7,6°C	set temp.	28,4°C
	2019-06-14 14:45:24	7,8°C	set temp.	30,8°C
	2019-06-14 14:44:24	8,5°C	set temp.	30,0°C
	2019-06-14 14:43:24	8,3°C	set temp.	29,1°C
	2019-06-14 14:42:24	7,8°C	set temp.	27,7°C
	2019-06-14 14:41:24	7,5°C	set temp.	29,8°C
	2019-06-14 14:40:24	8,3°C	set temp.	30,1°C

Figure 9 – Downloaded data

5.8 Statistics

For other units (not Premium) a key is required.

In this panel (Figure 10) statistics are downloaded from the unit and displayed, on the currently running program or program that has ended. Statistics are calculated separately for each segment and separately for each cycle. Data logging starts for calculation after 30 seconds from reaching the set temperature in the segment. The next data is registered every 1 minute. Statistical data is retrieved from the unit.


Data report			
Loop number: 1			
Segment number	Average temperature	Minimal temperature	Maximal temperature
1	6.10	3.60	21.49
2	8.00	4.97	8.52
Loop number: 2			
Segment number	Average temperature	Minimal temperature	Maximal temperature
1	4.42	3.18	9.38
2	7.99	4.04	8.52
Loop number: 3			
Segment number	Average temperature	Minimal temperature	Maximal temperature
1	4.41	3.28	16.22
2	7.80	4.53	9.05

Figure 10 – Statistics

5.9 Manage the users

Panel (Figure 11) is used to manage users of Premium Control application and allows to connect them with users available in units

The top menu is available:

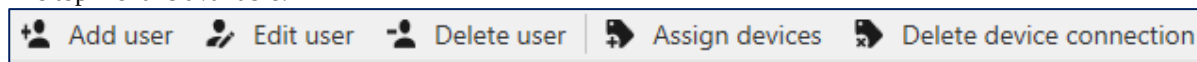

Add user – adds a new Premium Control user.


Edit user – allows to change data.

Delete user – deletes Premium Control user.

Assign devices – allows you to assign the Premium Control application to the user of the connected unit.

Delete device connection – allows you to delete the user's connection to the unit.

After entering the panel, users have a panel divided into two windows and a taskbar at the top. The upper window contains information about the created users and the bottom window contains information about the units assigned to them.

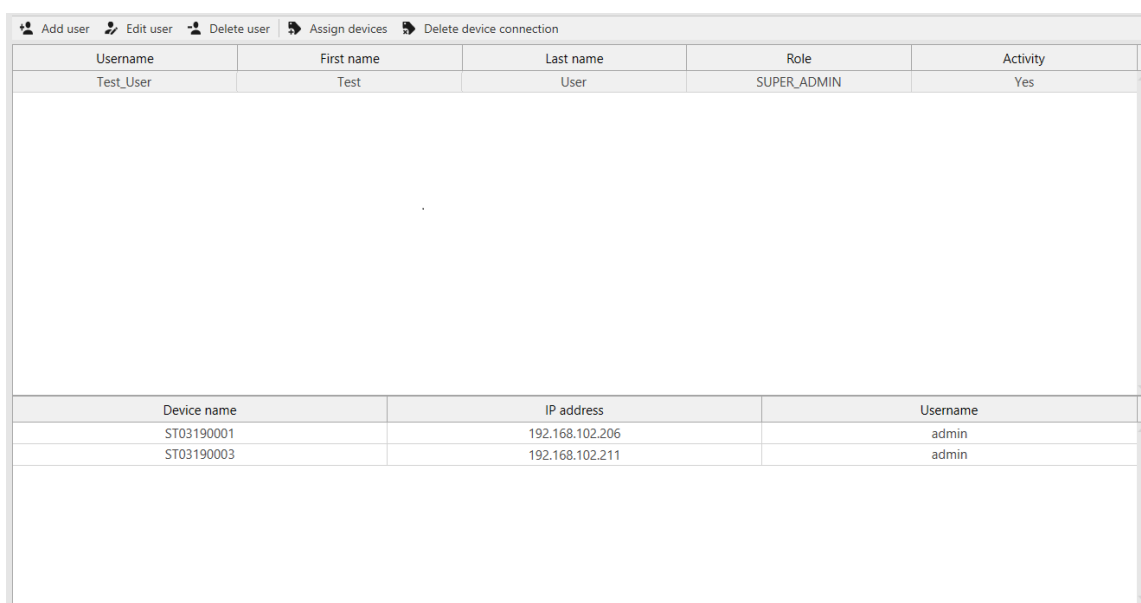

Username	First name	Last name	Role	Activity
Test_User	Test	User	SUPER_ADMIN	Yes

Device name	IP address	Username
ST03190001	192.168.102.206	admin
ST03190003	192.168.102.211	admin

Figure 11 – Manage the users

5.9.1. User creation / edition

Click Add a user or Edit a user, a creator will display (Figure 12) where in the following windows you can enter the name and user data, select the role and enter the password. All fields must be filled out.

Figure 12 – User edition

5.9.2. Assigning units

In order to enable a user connection to a given unit, it must be associated with the user on the unit. To do this, go to the **"manage users"** panel and then click **"Assign devices"**. A wizard will be launched where in the subsequent windows select the unit and the user with whom you want to link it.

Name	Programs limit	Role	Assigned to
admin	10	SUPER_ADMIN	admin

Figure 13 – Connecting to a unit

5.10 Manage units

Panel (Figure 14) enables managing the units available in the application (maximum 10 units).

The following top menu is available:



Refresh – refreshes the list of available units,

Add device – adds a new unit to the Premium Control application (the unit must be running and connected to the network),

Edit device – enables to change the address or description of the unit,

Remove device – removes the unit from the list.

Device name – the name given in the unit (Main menu -> Panel),

Description – the description of the unit provided in the Premium Control application during adding or edition of the unit,

Application connection status – „Connected” means that the Premium Control application is currently connected to the unit,

IP address – IP address of the unit,

Network connection status – „Connected” means that the unit is running and connected to the network.

Device name	Serial number	Description	Application connection status	IP address	Network connection status
ST03190001	ST03190001		Disconnected	192.168.102.206	Connected
IW1TF 197777	IW1TF 197777		Disconnected	192.168.102.138	Connected

Figure 14 – Managing units

5.10.1. Adding a unit

- 1) Click „add a unit” on the toolbar.
- 2) Enter IP address and description (optionally) in the new window (Figure 15).
- 3) Click „confirm” to add the unit to the list.
- 4) After adding, a wizard will be displayed that facilitates linking the logged user with the user on the unit (Figure 16).

 Add device

Device ID: 0

IP address
192.168.102.138

Description
New_device

Confirm

Cancel

Figure 15 – Adding a unit

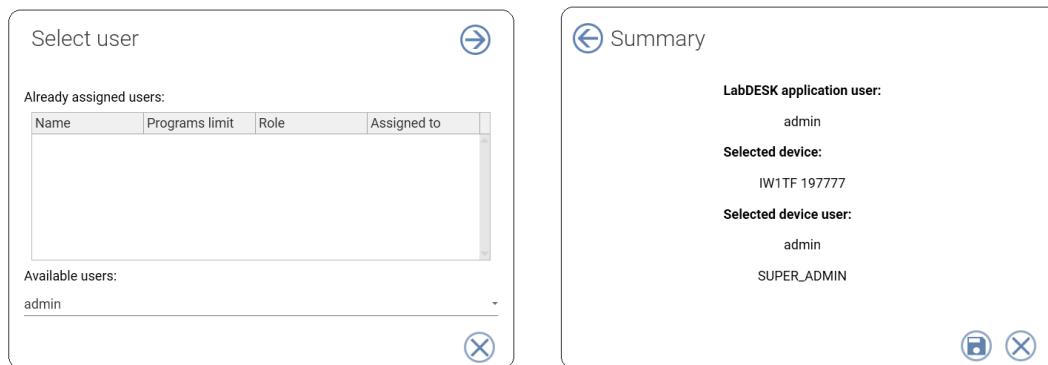


Figure 16 – Assigning a user

5.10.2. Connecting to the unit

When you're done assigning the unit to the user, go to the connection with the previously assigned unit. To do this, click on "**Devices**" on the right side of the menu (you can do this from each panel). During the first call you will be asked to enter the login and password from the unit.

This tab can be permanently attached to the program using this icon

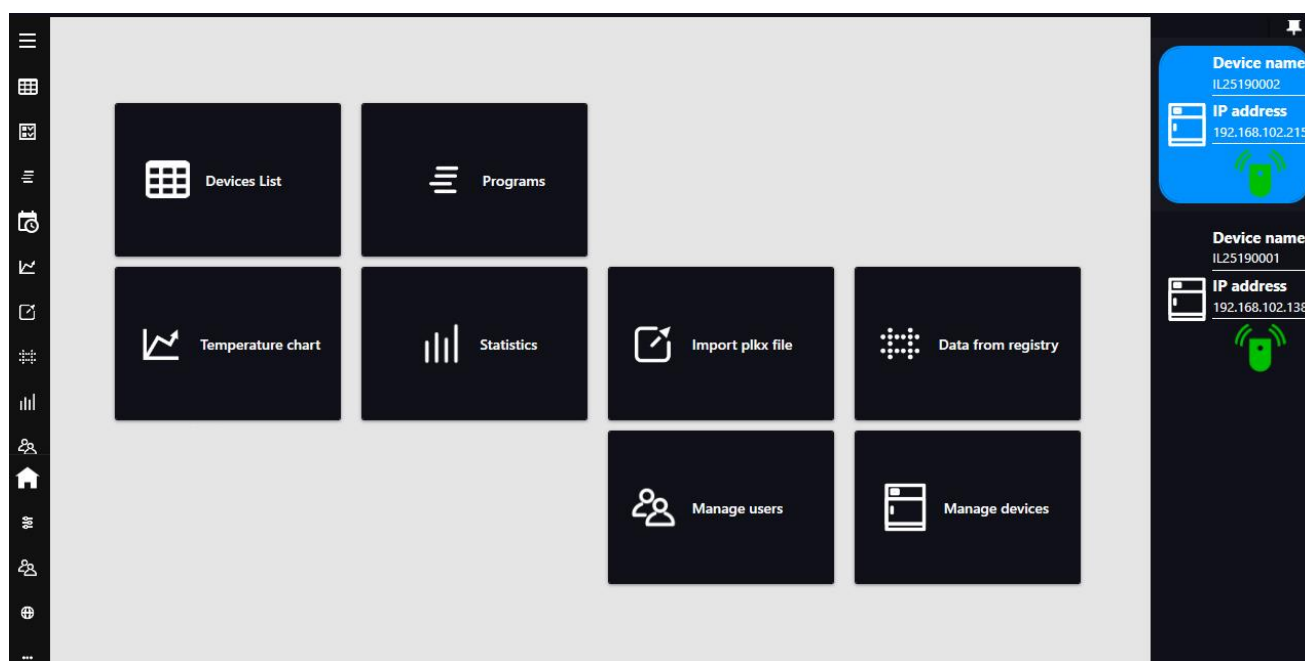


Figure 17 – Connecting to the unit

To connect with a selected unit click this icon . When the program connects with the unit, the icon will change into green .

5.11 Settings and information

The settings and information bar can be found in the lower left corner in the main menu.



displays main menu.



displays application settings panel.



displays information about the logged user (name, surname, assigned role).

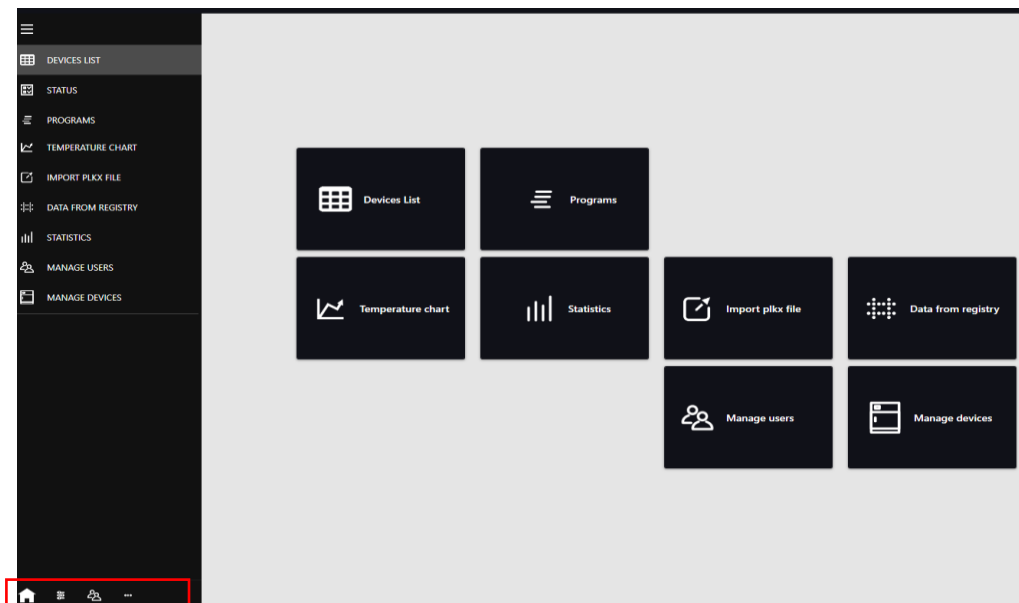


Figure 18 - Menu

5.11.1. Settings

Allows you to change Premium Control settings. Here you can set the following things:

- Program language,
- Enable or disable deletion of data after disconnection with the unit,
- Start window – the default window displayed after running the application,
- Database –database file, default location after installation: C:\PremiumControlDataBase,
- The frequency of unit temperature refresh,
- Device record storing – a folder in which data registered on a regular basis or downloaded from the unit will be saved.

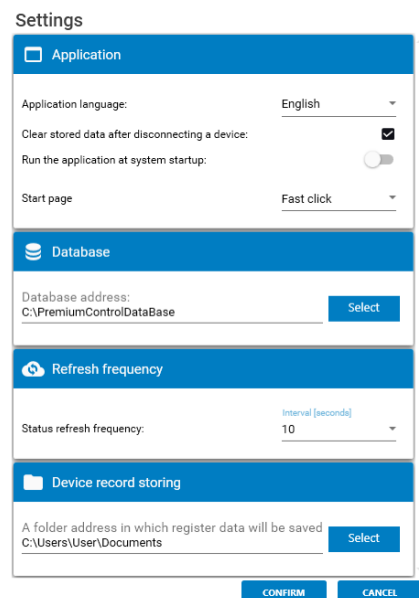


Figure 19 – Settings

5.12 Sterilisation (applies only to CO₂ incubators)

The CO₂ incubator is equipped with an automatically controlled hot-air sterilization program. To start the program, browse the STERILIZATION program from the programs list and press the PLAY button in the top menu.

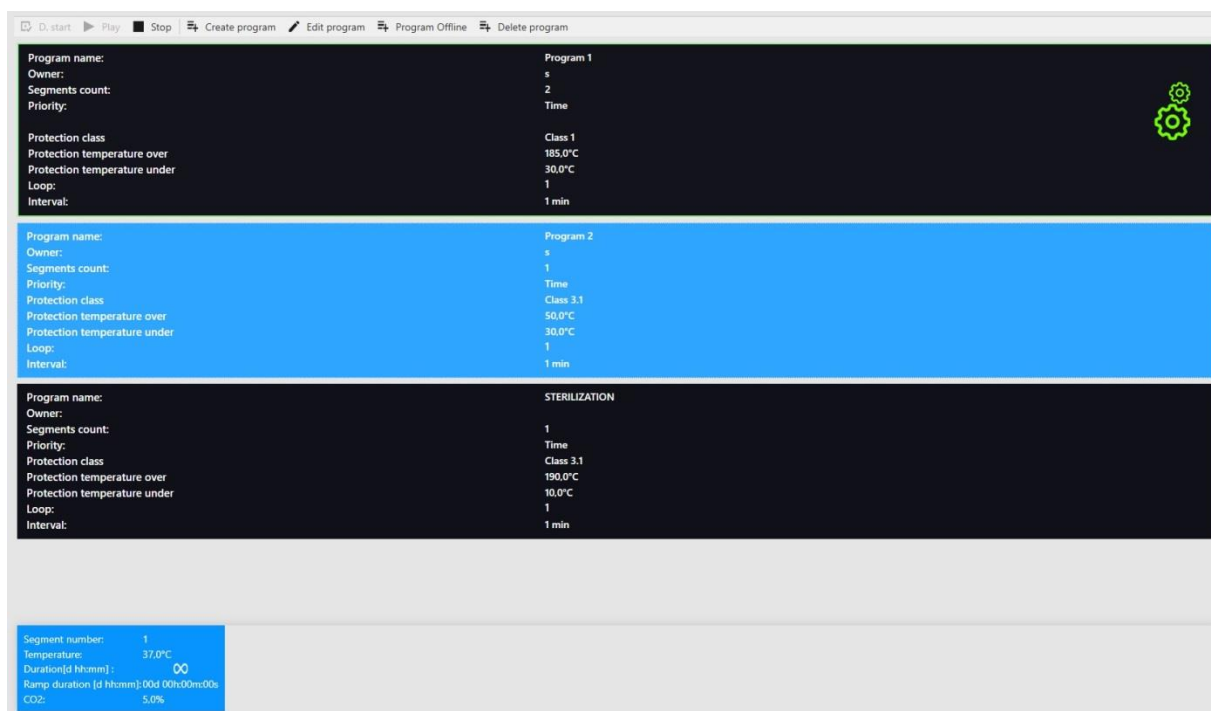


Figure 20 – Program list

Before running the sterilization program a confirmation prompt is displayed ([Figure 21](#)). Make sure that the water pan was emptied and all samples were removed from the chamber. Confirm the message. The program starts immediately or when the program delay time is over.



Figure 21 – Confirmation prompt

The effective sterilization process (heating up and holding) is automatically finished after holding 180°C inside the chamber for 2 hours. The information of successfully performed process pops up on the controller display ([Figure 22](#)).



Figure 12 – Successful sterilization

The total cycle of sterilization finishes when the temperature inside the chamber cools down to 37°C. Then, the information of completed sterilization occurs (Figure 23). Confirm the message.

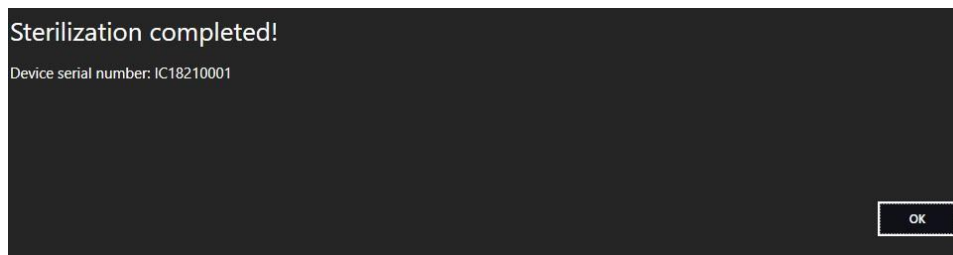


Figure 23 – Sterilization completed

If during the heating up and holding phase:

- the program is stopped,
- the external door is opened,
- the unit is turned off with the main power switch or in case of power failure,
- the unit's error occurs,

the sterilization is prematurely terminated. The information of incompletely performed sterilization pops up on the display (Figure 24) and the event is registered in the Event log ("Sterilization failed").

If the sterilization cycle was terminated prematurely, the germs inside the chamber might not be killed totally. It is recommended to repeat the sterilization.

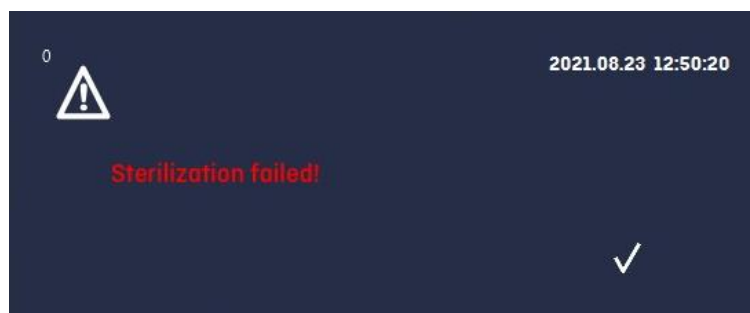


Figure 24 - Sterilization failed

If the door is open when the sterilization program is starting, the program won't start. The information of cancelled sterilization pops up on the controller display (Figure 25). Close the door and run the sterilization program again.

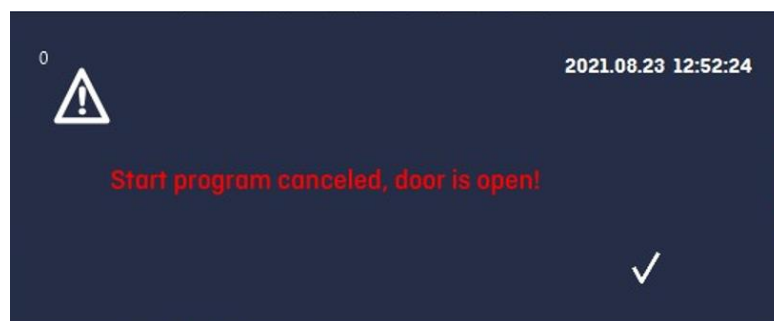


Figure 25 – Sterilization start cancelled

Your Distributor

Austria

VWR International GmbH
Graumannsgasse 7
1150 Vienna
Tel.: +43 01 97 002 0
Email: info.at@vwr.com

Belgium

VWR International bv
Researchpark Haasrode 2020
Geldenaaksebaan 464
3001 Leuven
Tel.: +32 016 385 011
Email: vwr.be@vwr.com

China

VWR (Shanghai) Co., Ltd
Bld.No.1, No.3728 Jinke Rd,
Pudong New District
Shanghai, 201203- China
Tel.: 400 821 8006
Email: info_china@vwr.com

Czech Republic

VWR International s. r. o.
Veetee Business Park
Pražská 442
CZ - 281 67 Stříbrná Skalice
Tel.: +420 321 570 321
Email: info.cz@vwr.com

Denmark

VWR International A/S
Tobaksvejen 21
2860 Søborg
Tel.: +45 43 86 87 88
Email: info.dk@vwr.com

Finland

VWR International Oy
Valimotie 17-19
00380 Helsinki
Tel.: +358 09 80 45 51
Email: info.fi@vwr.com

France

VWR International S.A.S.
Immeuble Estréo
1-3 rue d'Aurion
93114 Rosny-sous-Bois cedex
Tel.: 0 825 02 30 30* (national)
Tel.: +33 (0) 1 45 14 85 00 (international)
Email: info.fr@vwr.com
* 0,18 € TTC/ min + prix appel

Germany

VWR International GmbH
Hilpertstraße 20a
D - 64295 Darmstadt
Freecall: 0800 702 00 07
Tel.: +49 (0) 6151 3972 0 (international)
Email: info.de@vwr.com

Hungary

VWR International Kft.
Simon László u. 4.
4034 Debrecen
Tel.: +36 (52) 521-130
Email: info.hu@vwr.com

India

Avantor Performance Materials India
Limited
17th Floor, Building No. 5, Tower C
DLF Cyber City Phase – III
Gurgaon - 122002, Haryana
Tel.: +91-1244-65-6700
help@avantorinc.com

Ireland

VWR International Ltd
Orion Business Campus
Northwest Business Park
Ballycoolin
Dublin 15
Tel.: +353 01 88 22 222
Email sales.ie@vwr.com

Italy

VWR International S.r.l.
Via San Giusto 85
20153 Milano (MI)
Tel.: +39 02-3320311
Email: info.it@vwr.com

Korea

VWR International ~
17, Daehak 4-ro, Yeongtong-gu,
Suwon-si, Gyeonggi-do
Tel.: +82 31 645 7256
saleskorea@avantorsciences.com

Middle East & Africa

VWR International FZ-LLC
DSP Laboratory Complex
125, Floor 01
Dubai, United Arab Emirates
Tel.: +971 4 5573271
info.mea@vwr.com

The Netherlands

VWR International B.V.
Postbus 8198
1005 AD Amsterdam
Tel.: +31 020 4808 400
Email: info.nl@vwr.com

Norway

VWR International AS
Brynsalleen 4
0667 Oslo
Tel.: +47 22 90 00 00
Email: info.no@vwr.com

Poland

VWR International Sp. z o.o.
Limbowa 5
80-175 Gdansk
Tel.: +48 058 32 38 200
Email: info.pl@vwr.com

Portugal

VWR International –
Material de Laboratório, Lda
Centro Empresarial de Alfragide
Rua da Indústria, nº 6
2610-088 Amadora
Tel.: +351 21 3600 770
Email: info.pt@vwr.com

Singapore

VWR Singapore Pte Ltd
18 Gul Drive
Singapore 629468
Tel: +65 6505 0760
Email: sales.sg@vwr.com

Spain

VWR International Eurolab S.L.U.
C/ Tecnología 5-17
A-7 Llinars Park
08450 - Llinars del Vallès
Barcelona
Tel.: +34 902 222 897
Email: info.es@vwr.com

Sweden

VWR International AB
Fagerstagatan 18a
163 94 Stockholm
Tel.: +46 08 621 34 00
Email: info.se@vwr.com

Switzerland

VWR International AG
Lerzenstrasse 16/18
8953 Dietikon
Tel.: +41 044 745 13 13
Email: info.ch@vwr.com

UK

VWR International Ltd
Customer Service Centre
Hunter Boulevard
Magna Park
Lutterworth
Leicestershire
LE17 4XN
Tel.: +44 (0) 800 22 33 44
Email: uksales@vwr.com