

# Product Information

## VivoBrite™ Rapid Antibody Labeling Kit for Small Animal *In Vivo* Imaging, CF™ 680 SE

Catalog Number: 92160

### Components and Storage:

Material	Amount	Storage Condition*	Stability*
CF680, SE (#92160A)	3 X 0.1 umole	4 °C and moisture-free	If stored as directed, stable for at least 3 months from time of reception
DMSO, anhydrous (#99953)	150 uL		
Sodium bicarbonate solution, 1 M, pH 8.3 (#99954)	1 mL		
1X PBS, pH 7.4 (#99955)	50 mL		
Ultrafiltration vial, 10K MWCO (#99956)	3 X 1		
Reaction vial (#99957)	3 X 1		
Storage vial, sterile (#99958)	3 X 1		
Syringe, 1 mL (#99959)	3 X 1		
Syringe filter, 4 mm, 0.2 um (#99960)	3 X 1		
*If the reactive dye (component #92160A) is stored separately at -20°C, the shelf-life of the kit may be extended to at least 6 months.			

### Spectral Property

$\lambda_{abs}/\lambda_{em} = 681/698$  nm (protein conjugate in pH 7.4 buffer; see Figure 1 for spectra);  
 $\epsilon = \sim 210,000$ ;  
 $A_{280}/A_{max}$ , or CF = 0.09 (correction factor for estimating degree of protein labeling)

### Product Application

VivoBrite Rapid Antibody Labeling kits are designed for preparing fluorescently labeled antibodies for *in vivo* imaging in small animals. Each kit contains one of our superior near-IR CF dyes and everything else you need for carrying out the labeling reaction and purifying the labeled product. (For more details on near-IR CF dyes, see near-IR CF dye flyer at the Biotium website.)

The reactive dye has a succinimidyl ester group, which reacts with an amine group of the protein (i.e., lysine side-chain amine) to form a stable amide linkage. The dye is supplied in three separate vials, each containing sufficient dye for labeling 1-2 mg of IgG antibody. Following the labeling reaction, unconjugated dye is conveniently and rapidly removed by ultrafiltration using the provided membrane filtration vials.

### Protocol for Labeling IgG antibodies

The protocol below is for labeling 1 mg of an IgG antibody. The procedure may be scaled up or down for any amount of protein as long as the ratios of the reagents are maintained.

Note: Warm all reagents to room temperature before use.

#### 1. Prepare the Antibody for Labeling

If the antibody is already in solution at  $\geq 1$  mg/mL in PBS or a similar buffer free of any amine-containing chemicals or preservatives, such as Tris, ammonium or amino acids, it is ready for the labeling reaction (See next step). However, if any amine-containing chemical is present, it must be removed as it will interfere with the labeling reaction. This can be accomplished by dialysis against 1X PBS (#99955) using a micro-dialysis device (not provided), or by ultrafiltration using a membrane filtration vial provided in the kit (see procedure below). Presence of sodium azide does not affect the labeling.

To remove amine contaminants by ultrafiltration:

- Load the antibody solution (0.3 mL max each time- do not overfill or you may experience some loss of antibody) in the upper chamber of a filtration vial (#99956) and centrifuge at 14,000 rpm until nearly all the liquid is in the lower collection tube;
  - Empty the collection tube and add any additional antibody solution to the upper chamber. Repeat steps a & b until all the antibody solution has been centrifuged.
  - Dilute the concentrated antibody in the upper chamber to  $\sim 0.3$  mL with 1X PBS (#99955) and centrifuge to complete the second round of ultrafiltration;
  - Repeat Step c.
  - Transfer the ultrafiltered antibody to the provided reaction vial (#99957) with a total volume of  $\sim 0.9$  mL 1X PBS.
- Important: Save the used filtration vial for dye removal below.

Caution! Avoid touching the membrane of the filtration vial during liquid transfer using a pipet. Any damage to the membrane may result in serious product loss.

Note: The filtration vial contains a membrane that is permeable only to molecules with a molecular weight less than ~10,000 Daltons. Thus, small molecules freely pass through the membrane into the collection tube while the antibody, which has a molecular weight of ~150,000 Daltons, stays above the membrane.

## 2. Carry Out the Labeling Reaction

a) Add 0.1 mL 1M pH 8.3 sodium bicarbonate (#99954) to the antibody solution. If necessary, add 1X PBS so that the antibody concentration is ~1 mg/mL.

Note: Labeling efficiency varies with antibody concentration. In general, the higher the antibody concentration, the higher the labeling efficiency of the dye. This protocol is optimized for 1 mg/mL antibody concentration, at which the labeling efficiency of the dye can be expected to be ~30%. Thus, the degree of labeling (DOL) can be predicted using  $(n_d/n_p) \times 30\%$ , where  $n_d/n_p$  is the molar ratio of dye to antibody used in the labeling reaction.

b) Allow a vial of CF™680 SE (0.1 umole) to warm up to room temperature, and then add 25 uL anhydrous DMSO (#99953). Vortex to dissolve the dye and then centrifuge down the dye solution to the bottom of the vial.

c) Add ~12 uL of the dye stock from Step 2b to the antibody solution prepared in Step 2a. Protect the antibody/dye solution from light by wrapping the vial in aluminum foil and incubate the reaction for 1 hr at room temperature on a slow rocker.

Note: i) At this dye-to-protein molar ratio (i.e., ~7.2), the degree of labeling (DOL) is expected to be 2~2.2 dye/protein. Animal studies showed that antibodies labeled with near-IR CF dyes at a DOL of ~2 produced excellent results. However, because our near-IR CF dyes are engineered to be biocompatible and less prone to fluorescence quenching, antibodies labeled with a higher DOL (i.e., 2-4) may produce even better in vivo results, although this has yet to be confirmed. If a higher DOL is desired, increase the dye-to-protein ratio. The optimal DOL will need to be determined empirically.

ii) Any left-over dye stock solution may be kept for later use. If stored at -20 °C, the dye solution should be stable for at least one month.

## 3. Purify the Labeled Antibody

a) Transfer ~0.3 mL of the reaction solution from Step 2c to the upper chamber of a filtration vial (#99956) and centrifuge.

Caution! Avoid touching the membrane of the filtration vial during liquid transfer using a pipet. Any damage to the membrane may result in serious product loss.

b) Centrifuge the vial at 14,000 rpm until nearly all of the liquid is in the collection tube below (5-10 min). A small amount of 1X PBS may be used to rinse the reaction vial and complete the solution transfer. Repeat until all of the antibody/dye solution has been centrifuged.

c) Empty the collection tube, which contains the unconjugated free dye. Add ~0.3 mL 1X PBS (#99955) to the upper chamber of the vial, which now contains the labeled antibody. Flick or shake the tube gently to fully dissolve the labeled protein. Centrifuge again to ultrafilter the antibody solution.

Note: The filtration vial contains a separation membrane that is permeable only to molecules with a molecular weight less than ~10,000 Daltons. Thus, small molecules, such as phosphate and the unconjugated free dye, freely pass through the membrane into the collection tube while the antibody, which has a molecular weight of ~150,000 Daltons, stays above the membrane. The process is opposite to gel filtration, where small molecules are trapped in the gel matrix while proteins pass through the gel.

d) Repeat Step 3c two more times. By the third ultrafiltration, the color of the solution in the collection tube should be almost clear, indicating complete removal of the free dye from the labeled antibody. If some free dye is still detected in the third wash, one or two more washes may be necessary to completely remove the unconjugated dye.

e) Add an appropriate amount of 1X PBS to reach the desired antibody concentration.

f) Filter the labeled antibody solution from Step 3e into a sterile storage vial (#99958) using the provided syringe (#99959) and syringe filter (#99960). Remove the cap from the syringe tip and firmly attach the filter to the end of the syringe with a twist. Remove the plunger and add the labeled antibody carefully to the syringe. This may be best done with a pipette. Replace the plunger, and filter the antibody through in a slow, dropwise fashion.

Caution: Do not plunge the antibody through the filter too quickly or with too much force as the filter may disconnect from the syringe.

Note: Typical yield for preparing CF dye-labeled antibody using this protocol is 80-90%.

## 4. Determine the Degree of Labeling (DOL)

a) The final concentration of the antibody conjugate can be calculated from the formula:

$$[\text{conjugate}] \text{ (mg/mL)} = \{[A_{280} - (A_{\text{max}} \times \text{CF})]/1.4\} \times \text{dilution factor}$$

In the above formula:

[conjugate] is the concentration (in mg/mL) of the antibody conjugate solution prepared in Step 3f;

"dilution factor" is the fold of dilution used for spectral measurement (See Note below);

$A_{280}$  and  $A_{\text{max}}$  are the absorbance readings of the conjugate at 280 nm and the absorption maximum (~681 nm for CF™680), respectively; CF is the absorbance correction factor (0.09 for CF™680);

and the value 1.4 is the extinction coefficient of whole (H+L) IgG. Proteins other than whole IgG may have very different extinction coefficients, and thus, the formula will have to be adjusted for accurate determination of DOL.

Note: the antibody solution prepared in Step 3f is typically too concentrated for accurate absorbance measurement and should be diluted to approximately ~0.1 mg/mL. For example, if the labeled antibody is in ~0.5 mL PBS, which would roughly give an antibody concentration of ~2 mg/mL, you will need to perform a 1:20 dilution (i.e., dilution factor = 20) for spectral measurement.

b) Calculate the DOL according to the formula:

$$\text{DOL} = (A_{\text{max}} \times \text{Mwt} \times \text{dilution factor}) / (\epsilon \times [\text{conjugate}])$$

In the above formula:

$A_{\text{max}}$ , "dilution factor" and [conjugate] are as defined in Step 4a;

Mwt is the molecular weight of IgG: ~150,000;

and  $\epsilon$  is the molar extinction coefficient of CF™680: ~210,000.

## 5. Storage of Labeled Antibody

a) Store the labeled antibody at 4 °C and protect from light.

Note: Long term storage of antibodies less than 0.1 mg/mL may result in binding to plastic vials and degradation of product.

## Normalized Absorption and Emission Spectra of CF680 Conjugate

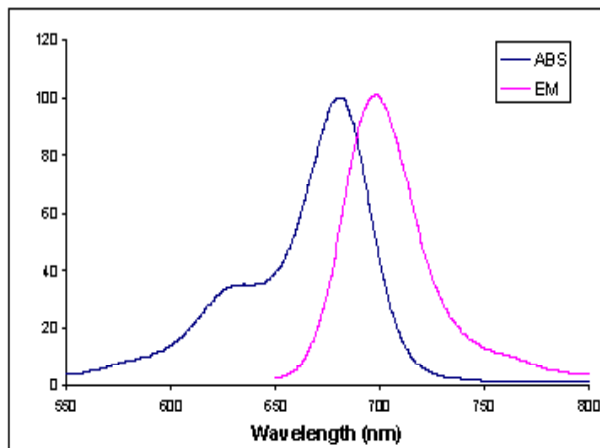


Figure 1. Absorption and emission spectra of CF680 conjugated goat anti-mouse IgG in pH 7.4 PBS buffer.

## Related Products

The following is a complete list of near-IR CF dye products:

Table 1. Near-IR CF™ Dyes and Related Products

Cat. #	Product Name	Abs/Em (nm)	Unit Size
92137	CF660, succinimidyl ester	667/683*	1 umole
92139	CF680, succinimidyl ester	681/698*	1 umole
92142	CF750, succinimidyl ester	755/775*	1 umole
92150	CF770, succinimidyl ester	770/797*	1 umole
90082	DMSO, anhydrous	N/A	10 mL
92160	VivoBrite rapid antibody labeling kit for small animal in vivo imaging, CF680 SE	681/698	3 labelings
92161	VivoBrite rapid antibody labeling kit for small animal in vivo imaging, CF750 SE	755/775	3 labelings
92162	VivoBrite rapid antibody labeling kit for small animal in vivo imaging, CF770 SE	770/797	3 labelings
20065	CF680 goat anti-mouse IgG(H+L), highly cross-absorbed, 2 mg/mL	681/698	0.25 mL
20067	CF680 goat anti-rabbit IgG(H+L), highly cross-absorbed, 2 mg/mL	681/698	0.25 mL
20070	CF750 goat anti-mouse IgG(H+L), 2 mg/mL	755/775	0.5 mL
20073	CF750 goat anti-rabbit IgG(H+L), 2 mg/mL	755/775	0.5 mL
20077	CF770 goat anti-mouse IgG(H+L), highly cross-absorbed, 2 mg/mL	770/797	0.25 mL
20078	CF770 goat anti-rabbit IgG(H+L), highly cross-absorbed, 2 mg/mL	770/797	0.25 mL

\*Wavelengths are for antibody conjugates in pH 7.4 PBS buffer.

Biotium products are high-quality reagents and materials intended for research purposes only. Some products are potentially hazardous chemicals - please read the Material Safety Data Sheet for additional information regarding handling potentially hazardous chemicals. Several of Biotium products and product applications are covered by U.S. and international patents and pending patents. Our products are not available for resale or other commercial uses without a specific agreement from Biotium, Inc. We welcome inquiries about licensing the use of our dyes, trademarks or technologies. Please submit inquiries by e-mail to [btinfo@biotium.com](mailto:btinfo@biotium.com). All names containing the designation ™ are trademarks of Biotium, Inc.

The manufacture, use, sale, offer for sale, or import of this product may be subject to one or more patents or pending applications owned or licensed by Biotium, Inc. The purchase price of this product includes a limited, non-transferable immunity and corresponding patent claims outside of the United States, to use only the purchased amount of the products and components of the product solely for the buyer's own internal research (whether the buyer is an academic or for-profit entity), in a manner consistent with the accompanying product literature. No right to use, sell or otherwise transfer (a) this product (b) its components or (c) materials made using this product or its components to a third party, or otherwise use this product or its components or materials made using this product or its components for any other purpose or a Commercial Purpose, is hereby granted expressly, by implication or by estoppel. Commercial Purpose means any activity by a party for consideration and include, but is not limited to: (1) use of the product or its components in manufacturing; (2) use of the product or its components to provide a service including without limitation, reporting information or data resulting from use of the product or its components for a fee or other commercial consideration; (3) use of the product or its components for therapeutic, diagnostic or prophylactic purposes; or (4) resale of the product or its components, whether or not such product or its components are resold for use in research. For products that are subject to multiple limited use label licenses, the most restrictive terms apply. This product is for research use only. Diagnostic uses require a separate license from Biotium, Inc. For information on purchasing a license to this product including for purposes other than research, contact Biotium, Inc., 3159 Corporate Place, Hayward, CA 94545, Tel: (510) 265-1027. Fax: (510) 265-1352.

\*CF™ dye technology is covered by pending US and international patents. ™ Alexa® is a registered trademark of Invitrogen, and Cy is a trademark of GE Healthcare.