




**Thermo Scientific Orion**  
ROSS pH Electrodes

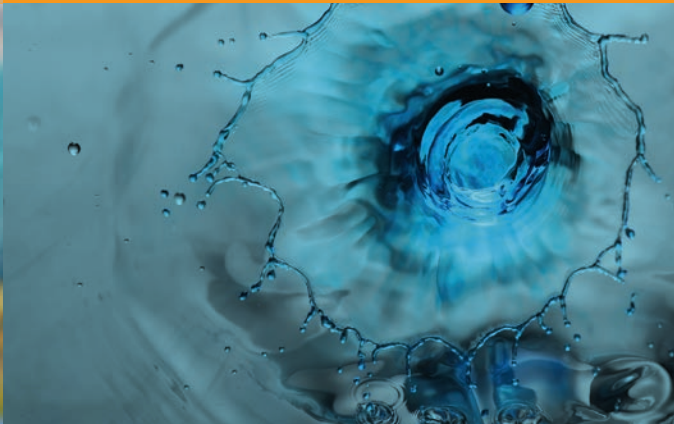

## Have confidence

in your pH measurements

Proven electrodes for accurate and reproducible pH measurements

**Thermo**  
SCIENTIFIC

# Thermo Scientific Orion ROSS pH Sensors

## Not all pH Electrodes are Created Equal

It is critical that the electrochemical measurements you perform day-to-day are accurate and reproducible. You rely on your pH electrodes to measure your samples quickly and precisely, making them an essential part of your lab. Whether you're testing samples in pharmaceuticals, food or general lab work, one incorrect reading can have devastating effects. Why take chances with a low quality electrode?



### Consumer Health Products

People's safety is of the utmost importance, so you need to be confident in the validity of your pH measurements.

### Pharmaceuticals

You need highly reproducible pH measurements to comply with strict GLP and GMP operating procedures.

### Food/Beverage Quality Tests

The critical control point (CCP) pH measurement needs to keep up with your manufacturing process.

### Field Water Testing

Efficient field testing requires rapid and reproducible pH measurements at any temperature to keep up with the workload and keep your facility EPA compliant.


### Quality Control Purposes

When budgets are constrained and throughput is critical, you want a reliable and long-lasting electrode that does not require frequent calibration and maintenance.

### Contract Testing Laboratory Work

When you need to run multiple tests but have a limited amount of sample to work with, you need to retain the integrity of your sample.

Chemists have been relying on Thermo Scientific™ Orion™ ROSS™ pH electrodes to provide critical measurements for over 35 years. Other electrodes come and go, but the ROSS reference system continues to be the foundation of the most stable and responsive electrodes available.

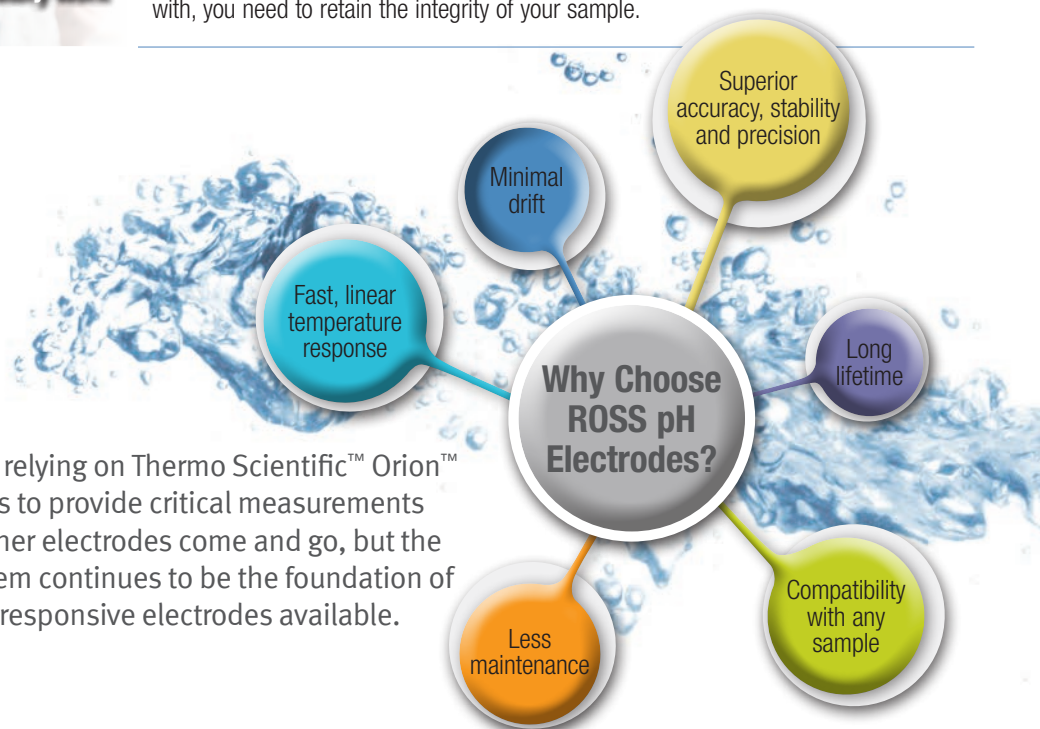



## Stable reference system and unique design features

### Orion ROSS pH Electrodes Deliver Best-in-Class Performance

Choose from a wide variety of Thermo Scientific™ Orion™ ROSS™ glass or epoxy body electrodes, such as: refillable, gel-filled, Sure-Flow™, Triode™, semi-micro, micro, flat bulb, spear tipped, and Ultra™ for extended warranty (18 months for gel / 24 months for refillable). ROSS electrode connectors are also available to fit most titrators.

#### ROSS Reference System

- Superior measurement stability
- Fast response
- High accuracy and precision in samples with varying temperatures
- No long-term drift

#### Double Junction Design

- Reference system is isolated from the sample by an intermediate fill solution to make the electrode compatible with samples such as those that contain TRIS, proteins and sulfides
- Electrode can be filled with a variety of solutions which are compatible with various sample types

#### Proprietary pH Sensitive Glass Bulb

- Optimized for sensitive response and wide dynamic range at all sample temperatures



#### Unique Coil Reference Design

- Path from junction to reference wire is 3 to 4 times longer than in conventional sensors minimizing dispersion of outer fill solution to the reference wire which increases reference stability
- Reference part of the electrode is maximally protected from the sample and lasts longer than in electrodes with conventional design

#### Inert Platinum Wire

- Does not introduce metal ions into sample
- Does not deplete over time

#### The Difference is in the Unique ROSS Reference System

$[I_3^-]/[I^-]$ iodide/triiodide ion pair chemistry with platinum wire	<b>Produces no metal ions to contaminate samples</b>
Attains thermal and chemical equilibrium quickly	<b>Fast to stabilize</b>
Ions are soluble over a broad temperature range	<b>Stable at all temperatures</b>
Chemical equilibrium is completely reversible	<b>Does not deplete reference material</b>





ROSS Electrode Specifications	
<b>Slope</b>	92-102 % of theoretical Nernst slope
<b>Isopotential Point</b>	pH 7
<b>Accuracy of Measuring a pH 6.86 Buffer after Standardization at 25 °C</b>	±0.03 pH for buffer at any temperature between 0-100 °C using automatic temperature compensation
<b>Speed of Response in 6.86 Buffer Going from 25 to 75 °C</b>	Stable to 0.01 pH within 30 seconds
<b>Speed of Response between 6.86 and 4.01 Buffers at 25 °C</b>	Stable to 0.005 pH within 30 seconds
<b>Considerations</b>	Silver free, TRIS, protein and sulfide compatible
<b>Temperature Ranges</b>	0 °C up to 100 °C 0 °C to 80 °C for gel filled triodes models
<b>Junction Types</b>	Sure-Flow (clog-free), ceramic, glass capillary, glass fiber, sleeve
<b>Body Materials</b>	Glass or epoxy
<b>Body Styles</b>	Standard, semi-micro, micro, rugged bulb, spear tip, flat surface
<b>Dimensions</b>	Total length: 120 mm, 155 mm (micro), 165 mm (semi-micro) Diameter: 12 mm Tip dimensions (D x L): 3 mm x 40 mm (micro); 6 mm x 95 mm (glass semi-micro), 8 mm x 95 mm (epoxy semi-micro)
<b>Connectors</b>	BNC, BNC waterproof, EDIN waterproof, MiniDin, screw cap, RCA
<b>Cable Lengths</b>	1 m, 3 m (select models)



For help in choosing a ROSS electrode, refer to our interactive Thermo Scientific Electrode Selection Guide.  
[www.thermoscientific.com/orionelectrodes](http://www.thermoscientific.com/orionelectrodes)



Select a Star-performance Thermo Scientific™ Orion™ pH meter to display your ROSS electrode measurements.  
[www.thermoscientific.com/orionmeters](http://www.thermoscientific.com/orionmeters)

[thermoscientific.com/water](http://thermoscientific.com/water)

© 2014 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific and its subsidiaries.

Water Analysis Instruments

**North America**  
 Toll Free: 1-800-225-1480  
 Tel: 1-978-232-6000  
[info.water@thermo.com](mailto:info.water@thermo.com)

**Netherlands**  
 Tel: (31) 020-4936270  
[info.water.uk@thermo.com](mailto:info.water.uk@thermo.com)

**India**  
 Tel: (91) 22-4157-8800  
[wai.asia@thermofisher.com](mailto:wai.asia@thermofisher.com)

**Japan**  
 Tel: (81) 045-453-9175  
[wai.asia@thermofisher.com](mailto:wai.asia@thermofisher.com)

**China**  
 Tel: (86) 21-68654588  
[wai.asia@thermofisher.com](mailto:wai.asia@thermofisher.com)

**Singapore**  
 Tel: (65) 6778-6876  
[wai.asia@thermofisher.com](mailto:wai.asia@thermofisher.com)

**Australia**  
 Tel: (613) 9757-4300  
 in Australia (1300) 735-295  
[InfoWaterAU@thermofisher.com](mailto:InfoWaterAU@thermofisher.com)

**Thermo**  
 SCIENTIFIC

A Thermo Fisher Scientific Brand