


MIOX Antibody (N-term)

Peptide Affinity Purified Rabbit Polyclonal Antibody (Pab)

Catalog #	Applications:	Reactivity:	Accessions:
AP13229a	WB, IHC, E	H	Q9UGB7 , NP_060054.4
Concentration:	Size:	Isotype:	Clone Name:
0.25 mg/ml	0.1 mg	Rabbit Ig	RB33052

Application Data:

Calculated MW: 33010 Da

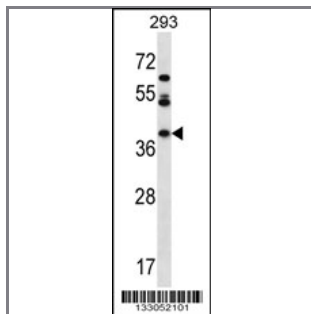

MIOX Antibody (N-term) (Cat. #AP13229a) western blot analysis in 293 cell line lysates (35ug/lane). This demonstrates the MIOX antibody detected the MIOX protein (arrow).



MIOX Antibody (N-term) (Cat. #AP13229a) immunohistochemistry analysis in formalin fixed and paraffin embedded human kidney tissue followed by peroxidase conjugation of the secondary antibody and DAB staining. This data demonstrates the use of MIOX Antibody (N-term) for immunohistochemistry. Clinical relevance has not been evaluated.

Gene ID:	Gene Symbol:
55586	MIOX

Other Names:

Inositol oxygenase; MIOX; ALDRL6; KSP32; RSOR

Target/Specificity:

This MIOX antibody is generated from rabbits immunized with a KLH conjugated synthetic peptide between 1-30 amino acids from the N-terminal region of human MIOX.

Application Notes:

The suggested dilution is:

ELISA 1:1,000
Western blotting 1:100~500
Immunohistochemistry 1:10~50

Format:

Purified polyclonal antibody supplied in PBS with 0.09% (W/V) sodium azide. This antibody is purified through a protein A column, followed by peptide affinity purification.

Storage:

MIOX Antibody (N-term) can be refrigerated at 2-8°C for up to 6 months. For long term storage, place the at -20°C.

Precautions:

MIOX Antibody (N-term) is for research use only and not for use in diagnostic or therapeutic procedures.
