ADP600 Series Polarimeter

QUICK START GUIDE
37-212'02 REV A, OCT-2023
Bellingham + Stanley, a Xylem brand, Xylem Analytics Germany GmbH, AM Achalaich 11, 82362 Weilheim, Germany
Email: sales.bs.uk@xylem.com

Thank you for purchasing this instrument. In order to ensure that this product provides many years of service please follow the guidance in this document. To register your instrument with the manufacturer and print a warranty certificate visit: www.bellinghamandstanley.com

Unpacking

WARNING: Risk of injury. This product is heavy. Single person lift could cause injury. Use assistance when moving or lifting.

Carefully remove all of the packing material. It is recommended that the box and other packing materials are retained so that, should the need arise, the instrument can be safely returned to the supplier for service. Check that all parts listed below are present and that no transit damage has occurred. If any are damaged or missing, contact the supplier immediately.


Contents list

Quantity	Item	Code
1	ADP600 Series Polarimeter	See below
1	Mains lead set	55-109
1	Power supply	55-106
1	Quick Start Guide	37-212
2	Card type ID tag	37-530
1	Lamp alignment target	37-510
1	Spare lamp	37-511
1	Touchscreen protector	19-204

ADP600 polarimeter code numbers:

Model	Complete polarimeter	Polarimeter module only	Wavelengths
ADP610	37-61	37-M61	589nm
ADP620	37-62	37-M62	546 & 589nm
ADP622	37-63	37-M63	365 & 589nm
ADP640	37-64	37-M64	405, 436, 546 & 589nm
ADP650	37-65	37-M65	365, 405, 436, 546 & 589nm
ADP660	37-66	37-M66	325, 365, 405, 436, 546 & 589nm
ADP670	37-67	37-M67	325, 365, 405, 436, 546, 633 & 589nm
ADS620P	37-72	37-M72	880 & 589nm

Place the adaptor on the instrument rails so that the flat surface is in contact with the temperature controlled plate.

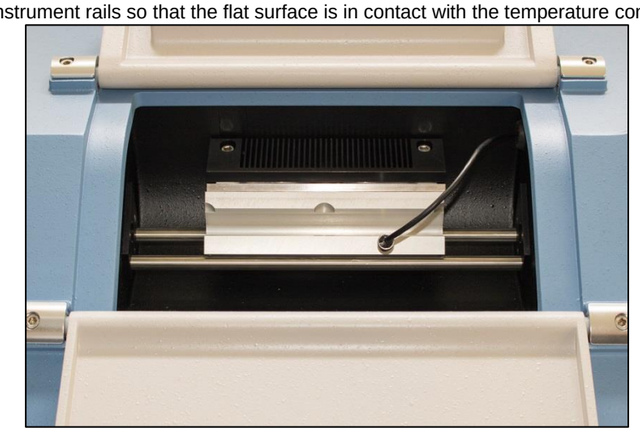Positioning the system

- Dry and indoors.
- Capable of holding the instrument's weight.
- Away from draughty or hot equipment like fans or heaters.
- Out of direct sunlight or strong ambient light.
- Away from potential sources of interference, such as RFI generating equipment.
- Within reach of a power point.
- Not using a power circuit that also has large motors or noise generating equipment connected to it.

Do not block the air vents.

Mains connection

The power supply adapter is supplied with a moulded mains cord and plug to suit one of several socket types. For UK leads, replace the fuse only with the type indicated on the plug.

- Voltage 110 to 230 V ~ ±10%
- Frequency 50 to 60 Hz
- Maximum current 2 A

RISK OF ELECTRIC SHOCK:

- For indoor use only.
- Must be kept dry.
- Disconnect the equipment from the mains supply before unplugging the mains lead from the instrument.
- Do not cover, designed to operate with free air convection.
- No cleaning required.
- Do not jet wash

Switching on and off

Plug in the power supply and switch on the mains supply, after startup the instrument will show a loading screen indicating the software version. Once loaded, the instrument will be in standby mode. To switch on, swipe a finger or stylus across the display. After a few moments the instrument will show the Measurement Display (or the Setup Wizard if switching on for the first time).


To switch to standby mode press and hold the Menu Button for several seconds and select Yes on the confirmation screen. It is then safe to turn off the mains supply and disconnect the power lead.

Manoeuvring through menus

The instrument features a capacitive, high-resolution touchscreen display. On-screen buttons can be easily used by simply touching the display. Use the supplied stylus and screen protector to keep the touchscreen clean and avoid damage from improper use. The touchscreen should only be used with a finger or stylus designed for use with a capacitive touchscreen.

Alternatively, a USB keyboard, available from Bellingham + Stanley, can be used to select items and navigate menus.

Many buttons feature a quick select option, indicated by the two arrows on the right of the button. Simply touch either arrow to cycle through the possible settings.

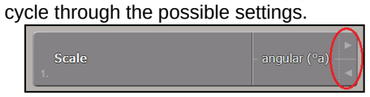When in a menu, the Quit and Home buttons become available. Press the Quit Button to return to the previous menu. Any changes made will be saved. The Home Button will exit to the Measurement Display, again any changes made will be saved.

Setup Wizard

When switching on for the first time, after the instrument has completed its startup procedures, the Setup Wizard will appear. First select the desired language. A list of different options will then appear:

Copy user guides and PC software to USB: Makes copies of the full User Guides and PC software compatible with the instrument to a removable USB storage device. Insert a USB removable storage device before selecting. Full User Guides provide detailed instructions on instrument operation.
Set clock: Set the instrument to local time and set the date format.
Set reading mode: Make changes to the measurement settings.
Set record settings: Specify how reading results will be printed or saved.
Set user settings: Alter the security settings to a level suitable for the area of operation.

Select each in turn and alter the settings as desired, then press Quit to return to the Setup Wizard Menu. Each menu will be marked by a tick after the settings are altered. Once all required settings are changed press Quit on the Setup Wizard Menu to show the Measurement Display.

This symbol is an internationally agreed indicator that the product bearing it should not be disposed of as general waste or garbage which might end up in landfill sites, but should instead be sent for special processing and/or recycling in those countries where appropriate legislation and facilities are in place.

Full User Guides can be downloaded from the instrument's Help Menu or Setup Wizard

Caution or warning.


Electrical hazard.

ADP600 Series Polarimeter

QUICK START GUIDE
37-212'02 REV A, OCT-2023
Bellingham + Stanley, a Xylem brand, Xylem Analytics Germany GmbH, AM Achalaich 11, 82362 Weilheim, Germany
Email: sales.bs.uk@xylem.com

Auto Zero


When the instrument is turned on a Zero calibration will be requested. This is to ensure that readings are as accurate as possible. Whilst the temperature stabilises the Reading Result section of the screen will display a warning message.

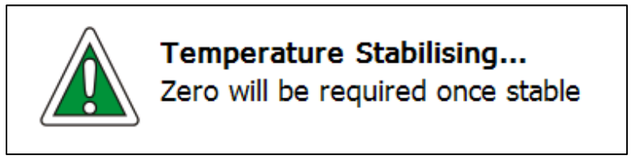The Zero can be performed with an empty chamber or adaptor, or a tube filled with the solvent to be used to dissolve the sample.

To Zero the instrument press the Zero button on the Measurement Display or access Zero via the Calibration Menu in the Main Menu. Readings for the Zero calibration will be taken with the settings in the Basic measurement settings of Mode or the currently selected Method. If a stability condition is set, this will also be used. The required number of scans will then be taken before the Zero point is set.

Measurement Display

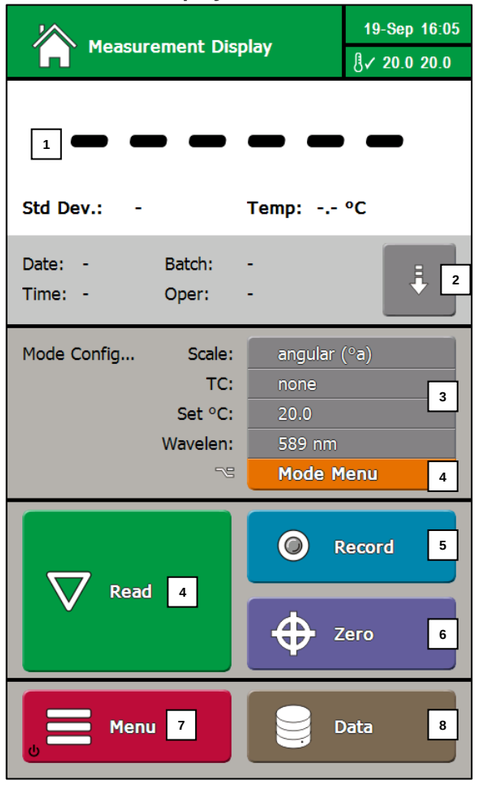From the Measurement Display readings can be taken, viewed and saved. Quick access to all the frequently used features of the instrument is also available.

1. Reading result: After taking a reading, the final result will be displayed here. Warning and alert messages will also appear here. For example, after the temperature set-point is changed, a Zero calibration will be requested.
2. Display Scan results: Press the button to reveal the individual result from each scan, the reading result is an average of these individual scans. Press again to view the reading configuration.
3. Reading Configuration Panel: Displays the measurement settings for the currently displayed reading or, if a result is not currently displayed, the settings to be used for the next reading. It also allows for quick changes to the Basic reading settings.
4. Mode button: Gives access to edit all Basic and Advanced reading settings.
5. Read button: Initiates a reading.
6. Record button: Pushing the button will record the result to the Reading Log or printer depending on the print settings.
7. Zero button: Initiates a Zero calibration routine. A Zero can be performed with any sample. Follow on screen instructions.

Holding down the above buttons for several seconds will open an associated Menu.

8. Menu button: Opens the instrument's Main Menu from where all instrument settings can be changed. Hold for several seconds to switch the instrument to standby.
9. Data button: Opens the Data Menu, this includes viewers for saved results and the System Log as well as the ability to download the results and any saved PDFs.

Reading

Lift open the chamber covers, place the sample tube in the adaptor and place the temperature probe in the tube neck. Ensure that the temperature probe is suitably immersed. With the sample tube and probe in place close the chamber lids, allow some time for thermal stabilisation and press the read button. The instrument will display a progress bar whilst data is collected and the measurement calculated.

The reading will be displayed using the scale and temperature compensation set in the Mode menu or current Method, shown on the Reading Configuration Panel.

If more than one scan was requested the displayed value is the average of multiple measurements.

Menu Structure

- **1. Calibration**
 - 1. Zero
 - 2. Span (Top Calibration)
 - 3. Undo
 - 4. Report
 - 5. Warning Frequency
- **2. Data**
 - 1. Saved Results
 - 2. System Log
 - 3. Copy Files to USB
 - 4. Move Files to USB
 - 5. Export Logs to USB
- **3. Setup**
 - 1. Setup Wizard…
 - 2. Methods
 - N. Add new…
 - S. Settings…
 - 1. Use Methods
 - 2. Hide preset
 - 3. Selection Type
 - 3. Users
 - N. Add new…
 - 1. administrator
 - 2. calibrator
 - 3. operator
 - S. Settings…
 - 1. Protection Level: 1. basic 5. cfr; 2. low; 3. medium; 4. high
 - 2. User Types: N. Add new… 1. administrator 2. calibrator 3. operator
 - 4. Reading
 - 1. Scales
 - N. Add new…: 1. From Library 2. From Data
 - 2. Temp. Comps.
 - N. Add new…: 1. From Data
 - 5.Recording
 - 1. Mode: 1. none 5. statistics 2. single 3. save 4. single + save
 - 2. Printer: 1. none 5. network 2. serial 3. usb 4. pdf
 - 3. LIMS Output
 - S. Select Fields
 - 1. none
 - 2. csv
 - 3. xml
 - 1. serial 2. network 3. usb file 4. ftp file
 - 4. Record After Read: 1. no 2. yes
 - 5. Traceability
 - 1. Location
 - 2. Batch Codes: 1. off 2. number 3. text 4. dateinc
 - 3. Auto Increment: 1. off 2. on
 - 4. Log Users: 1. off 2. number 3. text 4. login
 - 5. Remember Text: 1. none 2. batch 3. user 4. batch and user
 - 6. Ports
 - 1. Serial Port: 1. Baud Rate 2. Word Length
 - 2. Network: 1. DHCP 2. IP Address 3. Subnet Mask 4. Default Gateway
 - 6. System
 - 1. Language
 - 2. Remote Display: 1. off 2. on
 - 3. Display Brightness
 - 4. Screensaver
 - 5. Time / Date
 - 1. Date Format: 1. dd/mm/yy 2. mm/dd/yy
 - 2. Summer Time: 1. no 2. yes
 - 3. Time Server
 - 4. Set Clock…
 - 4. UTC Offset
 - 6. Temp Control: 1. Tick Criteria
 - 7. Clone From Backup
- **4. Maintenance**
 - 1. Maintenance Report
 - 2. Lamp Details
 - 3. Backup: 1. Backup Now 2. Schedule 3. Destination 4. Last Backup
 - 4. Clean Screen
 - 5. Update Software
- **5. Help**
 - 1. Information
 - 2. Contact
 - 3. User Guide
 - 4. PC software

This symbol is an internationally agreed indicator that the product bearing it should not be disposed of as general waste or garbage which might end up in landfill sites, but should instead be sent for special processing and/or recycling in those countries where appropriate legislation and facilities are in place.

Full User Guides can be downloaded from the instrument's Help Menu or Setup Wizard

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