


TEX261 Antibody (C-term)

Peptide Affinity Purified Rabbit Polyclonal Antibody (Pab)


Catalog #	Applications:	Reactivity:	Accessions:
AP13251b	WB, IHC, E	H	NP_653183.2 , Q6UWH6
Concentration:	Size:	Isotype:	Clone Name:
0.25 mg/ml	0.1 mg	Rabbit Ig	RB33130

Application Data:

Calculated MW: 22538 Da

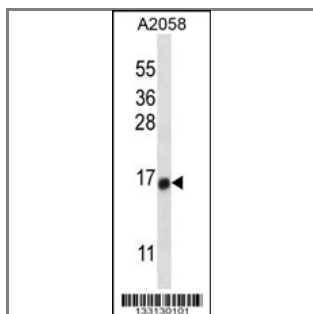

TEX261 Antibody (C-term) (Cat. #AP13251b) western blot analysis in A2058 cell line lysates (35ug/lane). This demonstrates the TEX261 antibody detected the TEX261 protein (arrow).

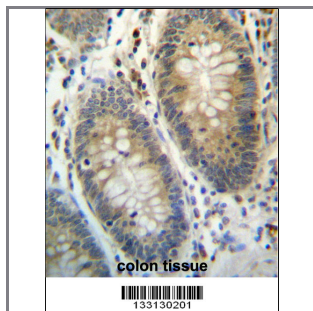

TEX261 Antibody (C-term) (Cat. #AP13251b) immunohistochemistry analysis in formalin fixed and paraffin embedded human colon tissue followed by peroxidase conjugation of the secondary antibody and DAB staining. This data demonstrates the use of TEX261 Antibody (C-term) for immunohistochemistry. Clinical relevance has not been evaluated.

Gene ID:	Gene Symbol:
113419	TEX261

Other Names:

Protein TEX261; TEX261;

Target/Specificity:

This TEX261 antibody is generated from rabbits immunized with a KLH conjugated synthetic peptide between 147-176 amino acids from the C-terminal region of human TEX261.

Application Notes:

The suggested dilution is:

ELISA 1:1,000
Western blotting 1:100~500
Immunohistochemistry 1:10~50

Format:

Purified polyclonal antibody supplied in PBS with 0.09% (W/V) sodium azide. This antibody is purified through a protein A column, followed by peptide affinity purification.

Storage:

TEX261 Antibody (C-term) can be refrigerated at 2-8°C for up to 6 months. For long term storage, place the at -20°C.

Precautions:

TEX261 Antibody (C-term) is for research use only and not for use in diagnostic or therapeutic procedures.
