

 DELTALAB GMDN code 14175	TECHNICAL DATA SHEET	429910en Revision 2	Data: 22/06/07 Page: 1/2

Codes: 429910, 429920, 429946
Screw cap conical tube. 15 ml (17x122 mm)

Catalogue 2007: page 27


Sanitary Product. IVDD 98 / 79 / CE

SPECIFIC USE AND DESCRIPTION

Code	429910	429920	429946
Case quantity	500	500	500
Case weight	3,8 kg	3,9 kg	3,9 kg
Case volume	0,029 m ³	0,029 m ³	0,038 m ³
Cases per pallet	60	60	54
Material	Transparent polypropylene	Transparent polypropylene	Transparent polypropylene
External mouth diameter	17 mm	17 mm	17 mm
External diameter with cap	20,9 mm	20,9 mm	20,9 mm
Total length with cap	122 mm	122 mm	122 mm
Maximum volume	15,4 ml	15,4 ml	15,4 ml
Maximum recommended volume	14 ml	14 ml	14 ml
Graduations	0,5 to 14 ml, graduations in 0,5 ml increments	0,5 to 14 ml, graduations in 0,5 ml increments	0,5 to 14 ml, graduations in 0,5 ml increments
Supplied in	Bag of 500 pieces	5 bags of 100 pieces each	Individual bag
Sterile (radiation)	NO	YES	YES
Autoclavable up to 121°C (not screwed cap)	YES	NO	NO
Withstand centrifugation:	7.000 G	5.000 G	5.000 G

- Conical bottom tube suitable for the centrifugation tests in immunology and microbiology labs.
- Screw on and splitted tube. Features external moulded graduations that guarantee non-existence of particles inside. Includes a matte writing surface of 55x10mm on side of tube.
- Caps ensure a positive hermetic seal from outside as well as from inside thanks to their inner ridge that perfectly fit the rim of containers.
- Blue ribbed high density polyethylene cap.

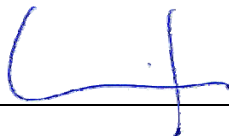

Elaborated	Verified	Authorised
		

 DELTALAB	TECHNICAL DATA SHEET	429910en	Date: 22/06/07
		Revision 2	Page: 2/2

- **Recommendation:** During autoclaving it is necessary to place the cap on the tube without screwing it. This allows the vapour to penetrate inside the tube so that avoids therefore deformation of the tube.

Code **429920** : Sterile, supplied in 100 pieces bags made of polyethylene (25 microns). On the bag are printed: IVD, CE, lot n°, expiry date, sterile, name of manufacturer.

Code **429946** : Sterile, supplied in individual bags made of polyethylene. On the bag are printed: IVD, CE, lot n°, expiry date, sterile, name of manufacturer.



Elaborated	Verified	Authorised
