INSTALLATION GUIDE

RephiDuo P Pack S2

RephiLe Catalogue Number:

RR800CP01

Corresponding Millipore Part Number:

PR0GT0S02

Suitable Millipore Systems：

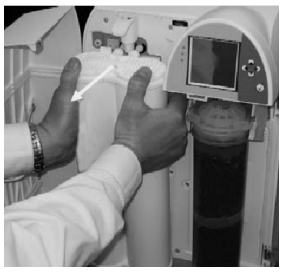
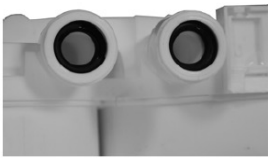
Milli-Q Integral

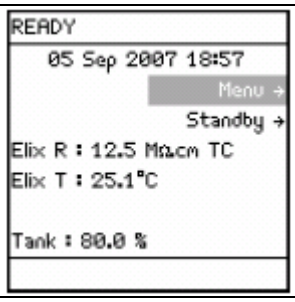
Elix Advantage

This product is comparable to the corresponding Millipore cartridge. Please refer to Millipore's User Manual for a detailed instruction on installation.

A brief description of the installation process is described below:

Step	Action	Result
1	Set the system to **STANDBY** Mode	
2	● Open the system left door. ● Lift up the cartridge Locking Handle.	

3	Remove the used cartridge. The system indicates that the cartridge is removed.	
4	• Remove the covers on the 2 ports of the new cartridge. • Look inside the ports to make sure the rubber O-rings are firmly in place at the bottom of the ports. • Wet the O-rings with pure water.	
5	• Push the top of the cartridge into the ports on the system first till it stops. • Then push on the bottom of the cartridge into its position.	
6	• Push the cartridge Locking Handle down. • Close the left door. The system indicates that a new cartridge is installed.	
7	Press the  button. System goes into flush mode	

8	System enters READY mode automatically after the flush step.	

The system is ready for use.

RephiSolo, RephiDuo and RephiQuatro are trademarks of RephiLe Bioscience, Ltd.

© 2010. RephiLe Bioscience, Ltd. All rights reserved.

All other trademarks are trademarks of their respective manufacturers.