INSTALLATION GUIDE

RephiDuo U Pack TL

RephiLe Catalogue Number:

RR500Q101

Comparable Millipore Part Number:

QGARDTL01

Suitable Millipore Systems:


RiOs 30 - 200

Elix 20 – 100, Elix Clinical, Elix Gulfstream


This product is a comparable to the corresponding Millipore cartridge. Please refer to Millipore's User Manual for a detailed instruction on installation.

A brief description of the installation process is described below:


1. Place the system into **Standby** by pressing the STANDBY button for about 2 seconds

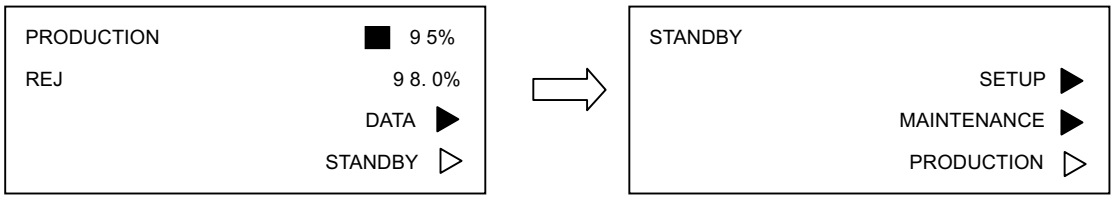


2. Press the MAINTENANCE button

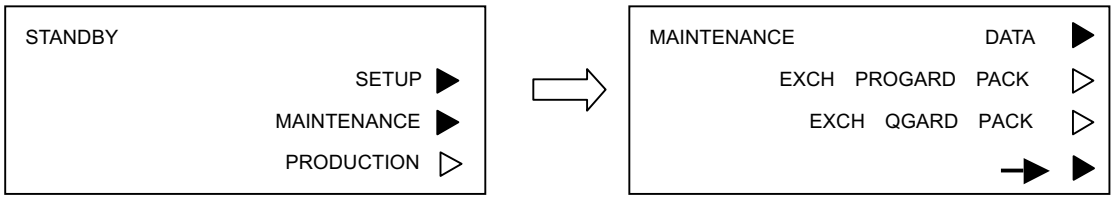

3. Press the EXCH QGARD PACK button for 2 seconds.

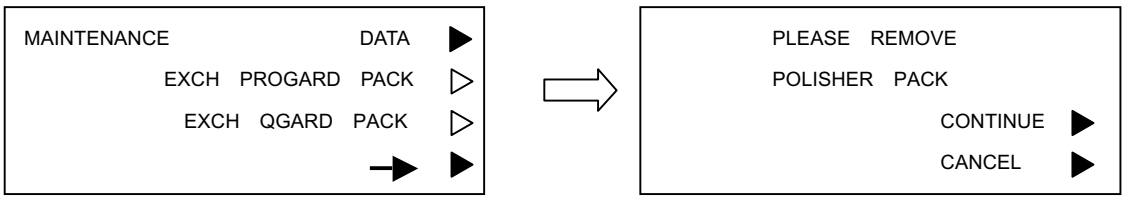

4. Remove the exhausted cartridge by lifting up the lock clip on top of the cartridge then remove the cartridge.
5. Press CONTINUE
6. Remove the new cartridge from the package. Verify product name and part number. Then remove the plastic covers on the two ports, check O-rings inside ports sit firm and flat against the bottom.
7. Wet the O-rings with water, install it onto the system.
8. Push the cartridge firmly in, then lock it into its position by pushing down the lock clip.
9. Press CONTINUE to go to PRODUCTION Mode

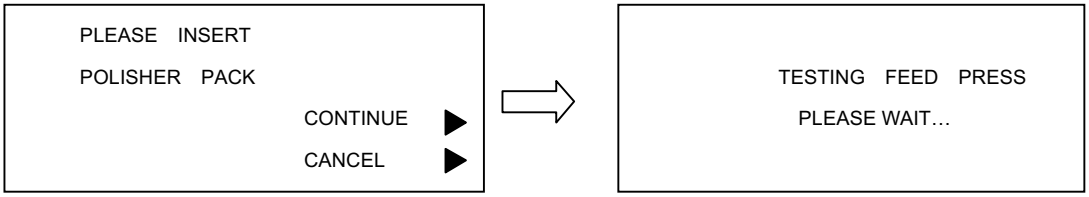

The LCD will display TESTING FEED PRESS for 10 seconds before entering PRODUCTION mode.

System is ready for use

RephiSolo, RephiDuo and RephiQuatro are trademarks of RephiLe Bioscience, Ltd.

© 2011. RephiLe Bioscience, Ltd. All rights reserved.

All other trademarks are trademarks of their respective manufacturers.