

NEW PPE REGULATION



In April 2018, a new PPE Regulation (EU) 2016/425 became applicable replacing the existing PPE Directive 89/686/EEC.

In addition, requirements for chemical protective gloves have been revised and a new EN ISO 374-1:2016 has been published which will imply several changes in the testing and labelling of our VWR Collection gloves.

In future, the Declaration of Conformity for PPE-labelled gloves will be accessible on the vwr.com website and packaging will include a direct link (QR Code) to access the Declaration of Conformity for each product.

1

THE NEW EN ISO 374-1:2016 STANDARD FOR CHEMICAL PROTECTIVE GLOVES


The new standard considers three types of chemical protective gloves:



TYPE A GLOVES - need to achieve at least a permeation level of 2 (longer than 30 minutes) against six of the chemicals listed in EN ISO 374-1. As a minimum, the six tested chemicals will be identified by their code letters under the flask pictogram.



TYPE B GLOVES - need to achieve at least a permeation level of 2 (longer than 30 minutes) against at least three of the chemicals listed in EN ISO 374-1. The tested chemicals will be identified by their code letter under the flask pictogram.

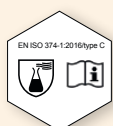

TYPE C GLOVES - need to achieve at least a permeation level 1 (longer than 10 minutes) against one of the chemicals listed in EN ISO 374-1.

2

The list of challenge chemicals in EN ISO 374-1 includes 6 additional chemicals

Code letter	Chemical	CAS no.	Class
A	Methanol	67-56-1	Primary alcohol
B	Acetone	67-64-1	Ketone
C	Acetonitrile	75-05-8	Nitrile compound
D	Dichloromethane	75-09-2	Chlorinated paraffin
E	Carbon disulphide	75-15-0	Organic compound containing sulphur
F	Toluene	108-88-3	Aromatic hydrocarbon
G	Diethylamine	109-89-7	Amine
H	Tetrahydrofuran	109-99-9	Heterocyclic and ether compound
I	Ethyl acetate	141-78-6	Ester
J	n-Heptane	142-82-5	Saturated hydrocarbon
K	Sodium hydroxide (40%)	1310-73-2	Inorganic base
L	Sulphuric acid (96%)	7664-93-9	Inorganic mineral acid
M	Nitric acid (65%)	7697-37-2	Inorganic mineral acid
N	Acetic acid (99%)	64-19-7	Organic acid
O	Ammonium hydroxide (25%)	1336-21-6	Organic base
P	Hydrogen peroxide (30%)	7722-84-1	Peroxide
S	Hydrofluoric acid (40%)	7664-39-3	Inorganic mineral acid
T	Formaldehyde (37%)	50-00-0	Aldehyde

3

As a newly added part to the standard series of chemical protective gloves, EN 374-4 specifies the determination of resistance to degradation by chemicals.

Coming into contact with a chemical can have an impact on the material properties of a protective glove, not only during contact with the chemical, but also subsequently. Degradation testing, therefore, includes the testing of the change of material's properties after exposure to chemicals.



Although EN ISO 374-1 contains no performance requirements for degradation, the EN 374-4 result will need to be included in the user instructions, and is indicated as a percentage degradation.

4

EN ISO 374-5:2016 is a newly introduced standard and specifies the requirements and test methods for protective gloves intended to protect the user against microorganisms.




It refers to resistance testing to **penetration of virus, fungi and bacteria**, and relates to the testing standard ISO 16604:2004.

VWR will continue to test and certify its gloves to the more stringent ASTM standard for microorganisms, which is FDA-approved and also accepted by the PPE Notified Body as part of the EU type examination procedure.

As such, VWR will continue to label in accordance with the internationally accepted ASTM F1671/F1671M-13 testing standard.

VWR, together with its trusted Notified Bodies, are jointly working on a transition plan in order to issue new EU type examination certificates and respective Declarations of Conformity in accordance with PPE Regulation (EU) 2016/425. In the course of this process, all Category III marked protective gloves will be tested and certified according to the newest EN ISO 374-1:2016 standard.

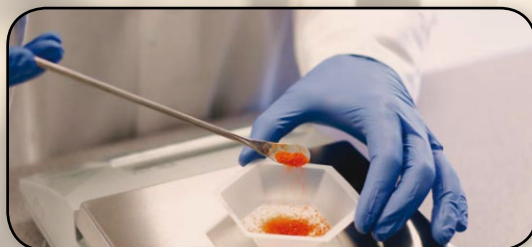

GO TO VWR.COM FOR THE LATEST NEWS, SPECIAL OFFERS AND DETAILS FROM YOUR LOCAL VWR SUPPORT TEAM

Austria

VWR International GmbH
Tel.: +43 1 97 002 0
Email: info.at@vwr.com

Belgium

VWR International bvba
Tel.: +32 (0) 16 385 011
Email: vwr.be@vwr.com

Czech Republic

VWR International s. r. o.
Tel.: +420 321 570 321
Email: info.cz@vwr.com

Denmark

VWR International A/S
Tel.: +45 43 86 87 88
Email: info.dk@vwr.com

Finland

VWR International Oy
Tel.: +358 (0) 9 80 45 51
Email: info.fi@vwr.com

France

VWR International S.A.S.
Tel.: 0 825 02 30 30* (national)
Tel.: +33 (0) 1 45 14 85 00 (international)
Email: info.fr@vwr.com
* 0,18 € TTC/min

Germany

VWR International GmbH
Tel.: 0800 702 00 07* (national)
Tel.: +49 (0) 6151 3972 0 (international)
Email: info.de@vwr.com
*Freecall

Hungary

VWR International Kft.
Tel.: +36 52 521130
Email: info.hu@vwr.com

Ireland / Northern Ireland

VWR International Ltd / VWR International (Northern Ireland) Ltd
Tel.: +353 (0) 1 88 22 222
Email: sales.ie@vwr.com

Italy

VWR International S.r.l.
Tel.: +39 02 3320311
Email: info.it@vwr.com

The Netherlands

VWR International B.V.
Tel.: +31 (0) 20 4808 400
Email: info.nl@vwr.com

Norway

VWR International AS
Tel.: +47 22 90 00 00
Email: info.no@vwr.com

Poland

VWR International Sp. z o.o.
Tel.: +48 58 32 38 200
Email: info.pl@vwr.com

Portugal

VWR International - Material de Laboratório, Lda
Tel.: +351 21 3600 770
Email: info.pt@vwr.com

Spain

VWR International Eurolab S.L.
Tel.: +34 902 222 897
Email: info.es@vwr.com

Sweden

VWR International AB
Tel.: +46 (0) 8 621 34 00
Email: kundservice.se@vwr.com

Switzerland

VWR International GmbH
Tel.: +41 (0) 44 745 13 13
Email: info.ch@vwr.com

UK

VWR International Ltd
Tel.: +44 (0) 800 22 33 44
Email: uksales@vwr.com