Operating Instructions

Original Operating Instructions

Cubis®

MCA Models
Electronic Mass Comparator



1000132460

SARTORIUS

Contents

1 About These Instructions 7
1.1 Validity 7
1.2 Related Documents 7
1.3 Target Groups 7
1.4 Symbols Used 8
1.4.1 Warnings in Operation Descriptions 8
1.4.2 Other Symbols Used 8

2 Safety Instructions 9
2.1 Intended Use 9
2.1.1 Modifications to the Device 9
2.1.2 Repairs and Maintenance on the Device 10
2.2 Qualifications of Personnel 10
2.3 Functionality of the Device Parts 10
2.4 Safety Information on the Device 10
2.5 Electrical Equipment 11
2.5.1 Power Supply Unit and Power Supply Cable 11
2.5.2 Connection Point for Power Supply Unit and Power Supply Cable 11
2.6 Conduct in an Emergency 11
2.7 Accessories 11
2.8 Risk of Injury During Transporting 12
2.9 Glass Breakage 12
2.10 Trip hazard from connection cable 12

3 Device Description 13
3.1 Device Overview 13
3.2 Draft Shield 14
3.3 Components in the weighing chamber 15
3.3.1 Devices with scale interval 1 µg 15
3.3.2 Devices with scale interval 0.01 mg and 0.1 mg 16
3.4 Connections 17
3.4.1 Back of the Device 17
3.4.2 Rear wall of weighing compartment 18
3.4.3 Control Unit 18
3.5 Transport Lock (Only for Models MCA1005MC | MCA605MC | MCA1004MC) 19
3.6 Accessories 19
3.7 Below-balance Weighing 19
3.8 Protective Caps and Covers for Connections 20

4 Operating Design 21
4.1 Operating Elements in the Main Menu 21
4.2 Operating Elements in Task Management 22
4.3 Operating Elements in the Weighing Display 23
4.4 User Guide 24
4.4.1 Messages 24
4.4.2 Help Function 24
4.4.3 Advanced User Guide 25
4.5 Status Center 26
4.6 Keypad 27

4.7 Status Display of the Buttons 27
4.8 Buttons in the Operating Display 28
4.8.1 Buttons for Navigation or Organization in Displays 28
4.8.2 Buttons for Editing or Managing Entries 29
4.8.3 Buttons for Weighing, Printing, and Export Functions 31
4.9 Displays in the Operating Display 32
4.10 Applications and Tasks 33
4.11 Licensable Device Functions 33
4.12 Navigating the Menus 34
4.13 Menu Structure 36
4.13.1 Main Menu 36
4.13.2 "Settings" Menu 37

5 Installation 43
5.1 Scope of Delivery 43
5.2 Selecting an Installation Site 44
5.3 Unpacking 44
5.4 Attaching or Removing the Control Unit 44
5.5 Placing the Device on its Side and Installing 45
5.6 Inserting the Connection Cable for the Control Unit 45
5.7 Inserting the Connection Cable for Ethernet 46
5.8 Positioning the Weighing Pan and Associated Components 46
5.8.1 Devices with scale interval 1 µg 46
5.8.2 Devices with scale interval 0.01 mg and 0.1 mg 47
5.9 Installing the Draft Shield 47
5.10 Removing the Transport Lock (Only for Models MCA1005MC | MCA605MC | MCA1004MC) 48
5.11 Installing the Control Unit 49
5.12 Acclimatizing 49

6 Commissioning 50
6.1 Connecting the Connection Cable for the Control Unit 50
6.2 Attaching the Connection Cable for Ethernet 50
6.3 Connecting the Power Supply Unit to the Device 50
6.4 Connecting the Power Supply 51
6.5 Connecting Accessories 51
6.6 Attaching Protective Caps and Covers 51

7 System Settings 52
7.1 Switching the Device On or Off 52
7.2 User Login or Logout 52
7.3 Performing System Settings 53
7.4 Using the Help Function 53
7.5 Activating Applications (QAPPs) 54
7.5.1 Activating QAPP Package With All Applications 54
7.5.2 Activating Individual Applications From QAPP Package 55
7.6 Switching Off the isoCAL Function 55
7.7 Configuring Motorized Opening and Closing of the Draft Shield 55
7.8 Managing User Data 56
7.8.1 Assigning a Password 56
7.8.2 Changing a User Profile or User Role 56
7.8.3 Creating a User Profile 56
7.8.4 Creating a User Role 57

7.9 Setting up Device for Network Printer 57
7.9.1 Selecting Configuration for Network Printing 57
7.9.2 Setting up Device for Network Printer via Independent Wi-Fi Network . 58
7.9.3 Setting up Device for Network Printer via Company Network 58
7.9.4 Configuring print profiles 59
7.10 Downloading Additional Information 60

8 Operation 61
8.1 Observing Warm-up Time 61
8.2 Manually Opening or Closing a Draft Shield 61
8.3 Motorized Opening or Closing of a Draft Shield 61
8.3.1 Opening or Closing on the Control Unit 61
8.3.2 Opening or Closing via Proximity Sensors 62
8.3.3 Opening or Closing a Motorized Draft Shield at the Door Handle 62
8.4 Leveling the Device 62
8.5 Calibration, Adjustment or Linearization 63
8.5.1 Adjusting With the isoCAL Function 63
8.5.2 Adjusting the Device Internally 64
8.6 Creating a Task 65
8.7 Adding Weighing and Print Profiles to a Task 65
8.8 Weighing 65
8.9 Weighing and Printing with ID Marking 66
8.9.1 Saving Values for Printouts 66
8.9.2 Marking Saved Values as Invalid 67
8.9.3 Marking Saved Values as Valid 67
8.9.4 Printing Saved Values 68
8.9.5 Ending the Task 68
8.10 Viewing the Data Storage Device 68
8.11 Running Applications (Examples) 69
8.11.1 Executing the “Unit change” function 69
8.11.2 Running the “Statistics” Application (from “Essentials” QAPP Package) 70

9 Cleaning and Maintenance 71
9.1 Preparing the Device for Cleaning 71
9.2 Cleaning the Device 71
9.3 Maintenance Schedule 72
9.4 Recommissioning 72
9.5 Performing a Software Update 72
9.6 Performing a QAPP Center Update 73

10 Malfunctions 75
10.1 Malfunctions on the Control Unit or During Weighing Processes 75

11 Decommissioning 76
11.1 Decommissioning the Device 76
11.2 Dismantling Device Parts 76
11.2.1 Removing the Weighing Pan and Associated Components 76
11.2.2 Dismounting the Draft Shield 77
11.2.3 Device with scale interval 1 µg: Dismounting the Adapter Module 78
11.2.4 Dismounting the Connection Cable 78

12 Transport . . . 79
12.1 Producing the Transport Lock (Only for Models MCA1005MC | MCA605MC | MCA1004MC) . . . 79
12.2 Transporting the Device . . . 80

13 Storage and Shipping . . . 81
13.1 Storing . . . 81
13.2 Returning the Device and Parts . . . 81

14 Disposal . . . 82
14.1 Disposing of the Device and Parts . . . 82

15 Technical Data . . . 83
15.1 Dimensions and Weights . . . 83
15.2 Ambient Conditions During Storage and Transport . . . 83
15.3 Installation Conditions . . . 84
15.3.1 Installation Site . . . 84
15.3.2 Ambient Conditions at the Installation Site . . . 85
15.4 Acclimatization Before Power Supply . . . 85
15.5 Electrical Data . . . 85
15.5.1 Power Supply . . . 85
15.5.2 Safety of Electrical Equipment . . . 86
15.5.3 Electromagnetic Compatibility . . . 86
15.5.4 Interfaces . . . 86
15.6 Warm-up Time to Reach Operating Temperature . . . 87
15.7 Model-specific and Licensable Device Functions . . . 87
15.8 Recommended Calibration Weight . . . 87
15.8.1 Models MCA206MC | MCA106MC | MCA66MC | MCA36MC . . . 87
15.8.2 Models MCA1005MC | MCA605MC | MCA1004MC . . . 87
15.9 Conditions for isoCAL Function . . . 88
15.10 Data Storage Device . . . 88
15.11 Integrated Clock . . . 88
15.12 Backup Battery . . . 88
15.13 Materials . . . 89
15.14 Cleaning Agents and Cleaning Procedures . . . 90
15.14.1 Approved Cleaning Agents . . . 90
15.14.2 Approved Cleaning Procedures . . . 90
15.15 Metrological Data . . . 91
15.15.1 Models MCA206MC | MCA106MC | MCA66MC | MCA36MC . . . 91
15.15.2 Models MCA1005MC | MCA605MC | MCA1004MC . . . 92

16 Accessories and Spare Parts . . . 93
16.1 Accessories . . . 93
16.2 Spare Parts . . . 94

17 Sartorius Service . . . 94

18 Conformity Documents . . . 95

19 Trademark Information . . . 95

20 Index . . . 98

1 About These Instructions

1.1 Validity

These instructions are part of the device; they must be read in full and kept in a safe place. These instructions apply to the device in the following versions:

Device	Model
Cubis® Mass Comparator	
Models with scale interval 1 µg:	MCA206MC | MC106MC | MCA66MC | MCA36MC
Models with scale interval 0.01 mg and 0.1 mg:	MCA1005MC | MCA605MC | MCA1004MC

1.2 Related Documents

▶ In addition to these instructions, observe the following documents:
- Only for models with scale interval 1 µg: Instructions for the motorized inner draft shield YDS125A
- Only for models with scale interval 0.01 mg and 0.1 mg: Instructions for the glass insert to reduce height YDSHR
- Instructions for the accessories used, e.g. YCM20MC climate module or weighing table
- Optional: Additional information about cleaning the device (Best Cleaning Practices for Cubis II Ultra-High Resolution Balances)

1.3 Target Groups

These instructions are addressed to the following target groups. The target groups must possess the specified knowledge.

Target Group	Knowledge and Qualifications
Operator	The operator is familiar with the device and the associated work processes. The operator understands the hazards which may arise when working with the device and knows how to prevent them.*

* If a person in the target group operates the software interface of the device, they are also the “user”.

1.4 Symbols Used

1.4.1 Warnings in Operation Descriptions

⚠ WARNING
Denotes a hazard that may result in death or serious injury if it is **not** avoided.

⚠ CAUTION
Denotes a hazard that may result in moderate or minor injury if it is **not** avoided.

NOTICE
Denotes a hazard that may result in property damage if it is **not** avoided.

1.4.2 Other Symbols Used

▶	Required action: Describes activities that must be carried out. The activities in the sequence must be carried out in succession.
▷	Result: Describes the result of the activities carried out.
[ ]	Refers to operating and display elements. Indicates status, warning, and error messages.

Figures in These Instructions

Depending on the device configuration, the figures depicting the device and operating display may differ slightly from the supplied device. The variants shown in these instructions are examples.

2 Safety Instructions

2.1 Intended Use

The device is a high-resolution mass comparator, which can be used in testing institutes and calibration services. The device is intended for accurate mass determination of test weights.

The device can be operated as follows:

- In individual operation
- Connected to a PC
- Integrated into a network

The device is intended solely for use in accordance with these instructions. Any other use is considered **improper** and can impair the protection supported by the device, e.g. protection against mechanical dangers.

Foreseeable Misuse

The following applications are **not** permitted: Operation under any atmosphere other than the normal atmosphere.

Operating Conditions for the Device

Do **not** use the device in potentially explosive environments. Only use the device indoors.

Do **not** use the device in a vacuum or negative pressure environment (permissible air pressure range see Chapter “15.3 Installation Conditions”, page 84).

Do **not** use physical measures to alter the as-delivered state of the device, and only connect approved accessories (see Chapter “16 Accessories and Spare Parts”, page 93).

Adhere to the maintenance and calibration intervals for the device (see Chapter “9.3 Maintenance Schedule”, page 72).

The device may only be used with the equipment and under the operating conditions described in the Technical Data section of these instructions.

2.1.1 Modifications to the Device

If the device is modified: Persons may be put at risk. Device-specific documents and product approvals may lose their validity.

For queries regarding modifications to the device, contact Sartorius.

2.1.2 Repairs and Maintenance on the Device

Device repairs and maintenance may only be carried out by persons with specialized knowledge of the device. If the device is **not** repaired or maintained by a specialist: Persons may be put at risk. Device-specific documents and product approvals may lose their validity.

Sartorius recommends that any repair work, even that carried out after the end of the warranty period, is carried out by Sartorius Service or after consulting with Sartorius Service.

Only the maintenance tasks described in these instructions should be carried out. For maintenance tasks that need to be carried out by Sartorius Service, contact Sartorius Service.

2.2 Qualifications of Personnel

Persons who do not possess adequate knowledge about how to use the device safely may injure themselves and other persons.

If a particular qualification is required for an activity: The target group is specified. If **no** qualification is specified: The activity can be performed by the "Operator" target group.

2.3 Functionality of the Device Parts

Non-functioning device parts, e.g., as a result of damage or wear, can cause malfunctions. There is a risk of injury to persons.

- If device parts are **not** functioning: Do **not** use the device.
- Comply with the maintenance intervals (intervals and maintenance tasks see Chapter "9.3 Maintenance Schedule", page 72).

2.4 Safety Information on the Device

Symbols, e.g., warnings and safety stickers, are safety information for handling the device. Missing or illegible safety information may result in this information **not** being observed. There is a risk of injury to persons.

- Do **not** conceal, remove, or modify the symbols.
- Replace the symbols if they become illegible.

2.5 Electrical Equipment

2.5.1 Power Supply Unit and Power Supply Cable

The use of an **unauthorized** power supply unit or power supply cable may cause life-threatening injuries as a result of electric shocks, for example.

- ▶ Only use the original power supply unit and power supply cable provided with the device.
- ▶ If the power supply unit or power supply cable needs to be replaced: Contact Sartorius Service. Do **not** repair or modify the power supply unit or power supply cable.

2.5.2 Connection Point for Power Supply Unit and Power Supply Cable

If the connection point for the power supply unit and power supply cable is not suitable: There is a risk of people becoming seriously injured, e.g., due to electric shocks.

- ▶ Protect the power supply unit and the power supply cable from liquids.
- ▶ Do **not** use a damaged power supply unit or power supply cable.

2.6 Conduct in an Emergency

If an emergency occurs, e.g., due to malfunctions of the device or dangerous situations: Persons may be injured. The device must be immediately taken out of operation:

- ▶ Disconnect the device from the power supply.
- ▶ Prevent the device from recommissioning.

2.7 Accessories

The use of unsuitable accessories can affect the functionality and safety of the device and has the following consequences:

- — Risk of injury to persons
- — Damage, malfunctions, or failure of the device

- ▶ Only use accessories that have been approved by Sartorius and are safe for use.

2.8 Risk of Injury During Transporting

If the device is **not** transported correctly: The device could fall and cause an injury, e.g. a foot injury. If the device is **not** placed down correctly, e.g., on a laboratory table: Fingers may be crushed.

- ▶ Disconnect the device from all connections at the installation site.
- ▶ Use both hands when transporting the device and setting it down. To do this, reach sideways under the device with both hands in the rear part of the device.
- ▶ Do **not** carry the device on the draft shield or the control unit.

2.9 Glass Breakage

Glass components can break if they fall or are handled incorrectly. Glass fragments can cause cuts.

- ▶ Do **not** use sharp or hard objects on the operating display.
- ▶ Do **not** allow items to fall onto the operating display.
- ▶ In the event of damage to the operating display or draft shield, do **not** use the device. Contact Sartorius Service.

2.10 Trip hazard from connection cable

If the device's connection cables, e.g. the power supply cable, are laid haphazardly: People may trip over the connection cables and injure themselves.

- ▶ Lay all connection cables in such a way that trip hazards are avoided.

3 Device Description

3.1 Device Overview

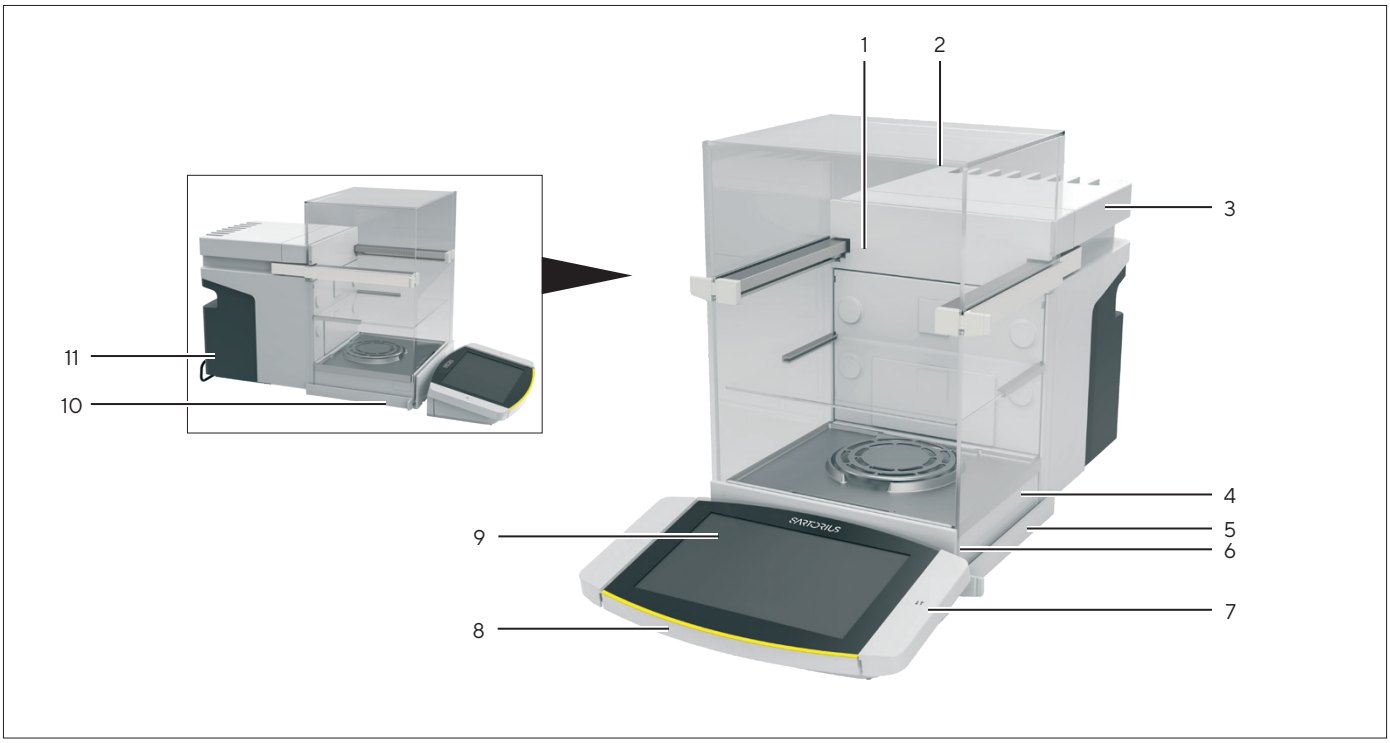


Fig. 1: Cubis® Mass Comparator (example)

Pos.	Name	Description
1	Weighing chamber	
2	Draft shield	
3	Cooling element	
4	Weighing module	Contained in the housing.
5	Leveling foot	Motorically adjustable. Located on the bottom of the device.
6	Proximity sensor	The proximity sensor is intended for contact-free opening or closing of the following components: – Doors of the motorized draft shield – Device with scale interval 1 µg: Motorized inner draft shield
7	Palm key	Enables opening of the doors of a motorized draft shield.
8	Control Unit	Removable.
9	Operating display	
10	Transport lock	Only for models MCA1005MC | MCA605MC | MCA1004MC: For securing the weighing system during transport (see Chapter "3.5 Transport Lock (Only for Models MCA1005MC | MCA605MC | MCA1004MC)", page 19).
11	Manufacturer's ID label	

3.2 Draft Shield

The draft shield can be opened in a motorized manner or manually.

When opening the draft shield, the top and side doors slide backwards. The top and side doors can be opened together or independently.

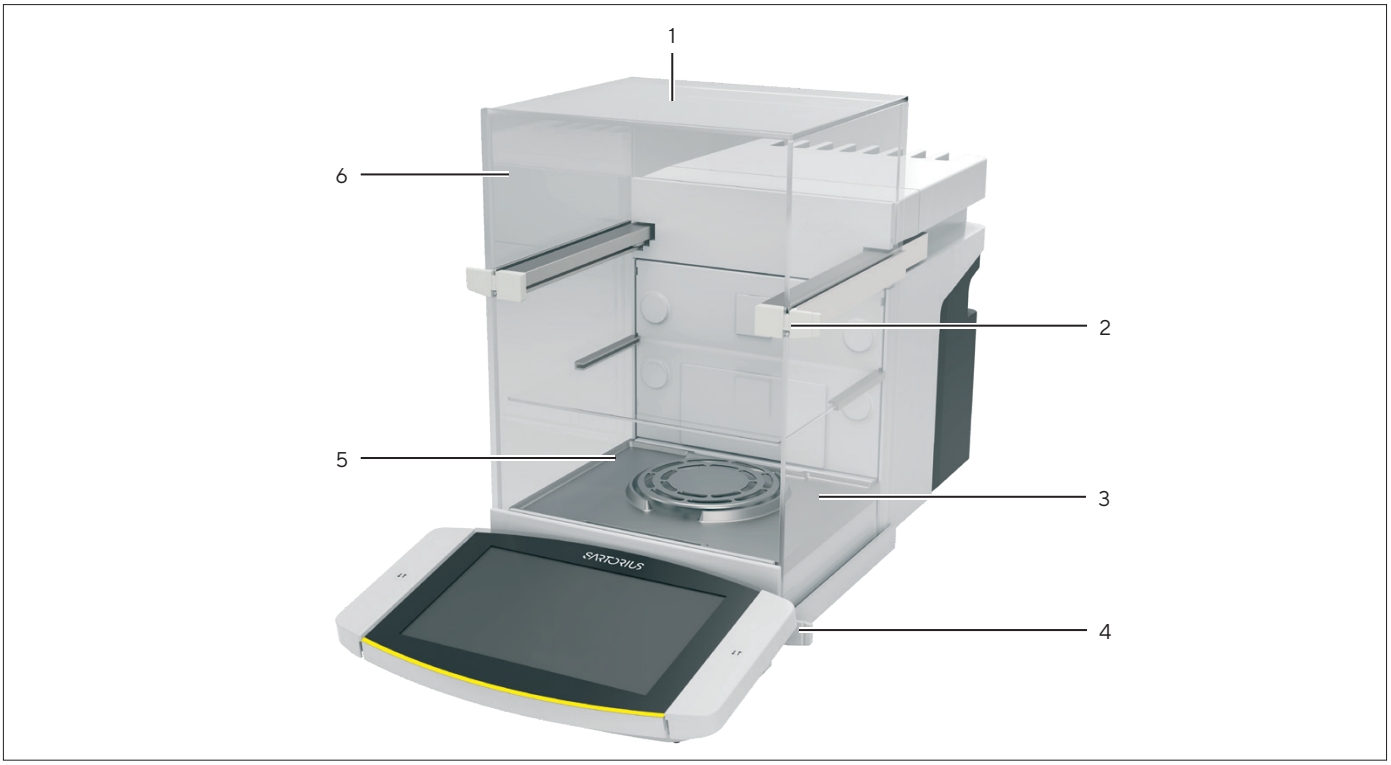


Fig. 2: Draft shield, motorized, all doors closed (example)

Pos.	Name	Description
1	Rear panel	Is installed on delivery.
2	Top door handle	Opens the top door.
3	Side door	
4	Bottom door handle	Opens the side door.
5	Front panel	
6	Top door	Consists of one top panel and two side panels.

3.3 Components in the weighing chamber

3.3.1 Devices with scale interval 1 µg

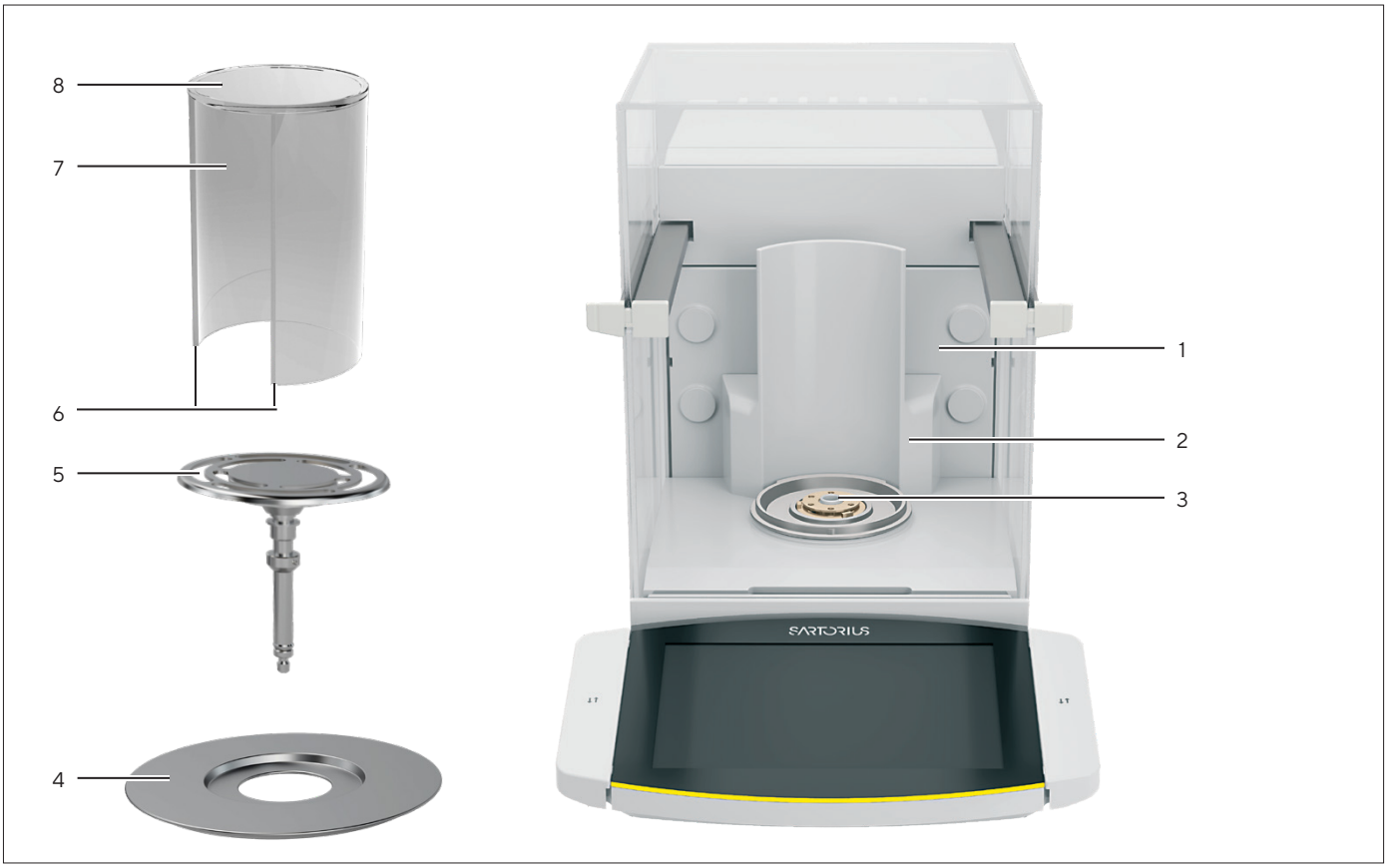


Fig. 3: Components in the weighing chamber (example)

Pos.	Name	Description
1	Rear wall of weighing compartment	
2	Adapter module	Can be taken out, e.g. for cleaning.
3	Pan retainer	For inserting the weighing pan.
4	Shield plate	
5	Weighing pan	
6	Glass cylinder opening	
7	Glass cylinder	
8	Glass lid	

3.3.2 Devices with scale interval 0.01 mg and 0.1 mg

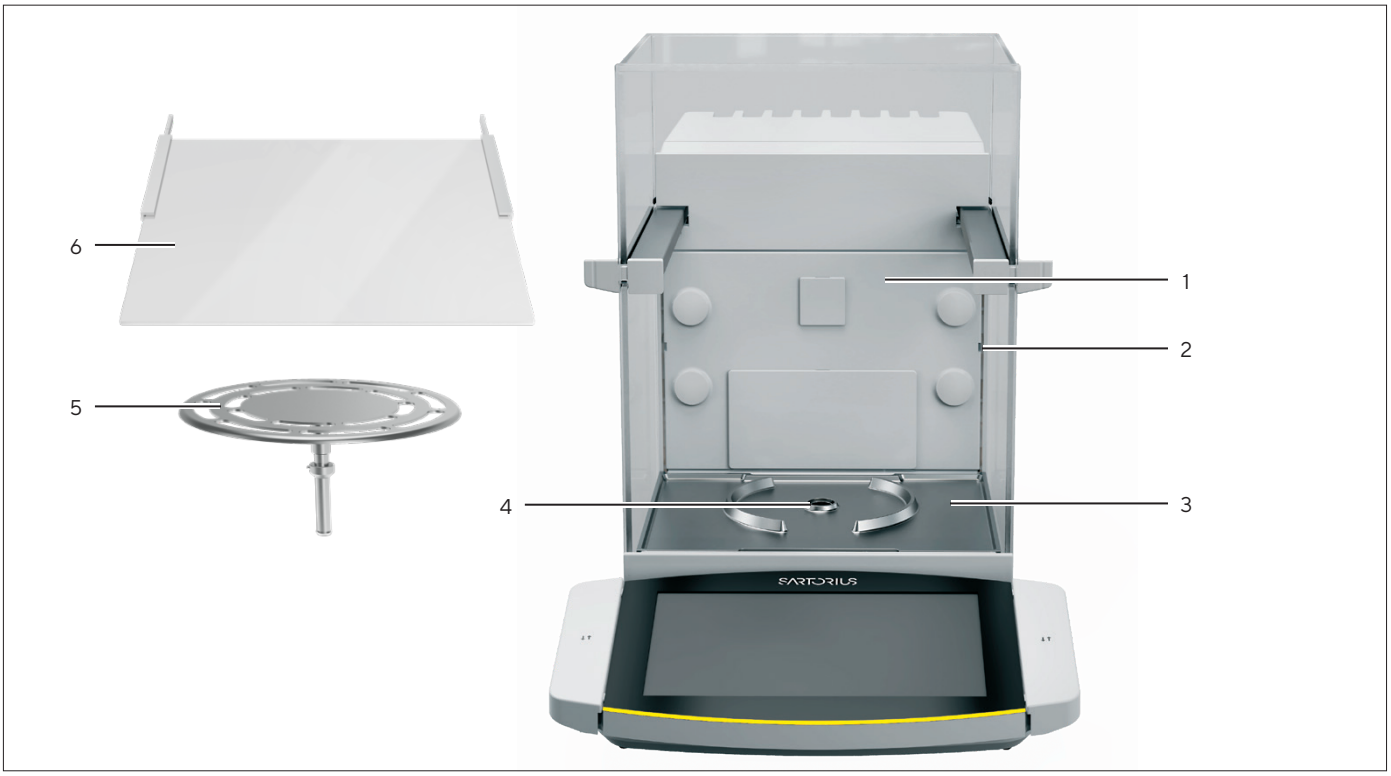


Fig. 4: Components in the weighing chamber (example)

Pos.	Name	Description
1	Rear wall of weighing compartment	
2	Guide for glass insert	
3	Base plate of weighing compartment	Can be taken out, e.g. for cleaning.
4	Pan retainer	Holds the weighing pan. Prevents the weighing pan from turning.
5	Weighing pan	
6	Glass insert	For reducing the height of the weighing chamber. This reduces drafts.

3.4 Connections

3.4.1 Back of the Device

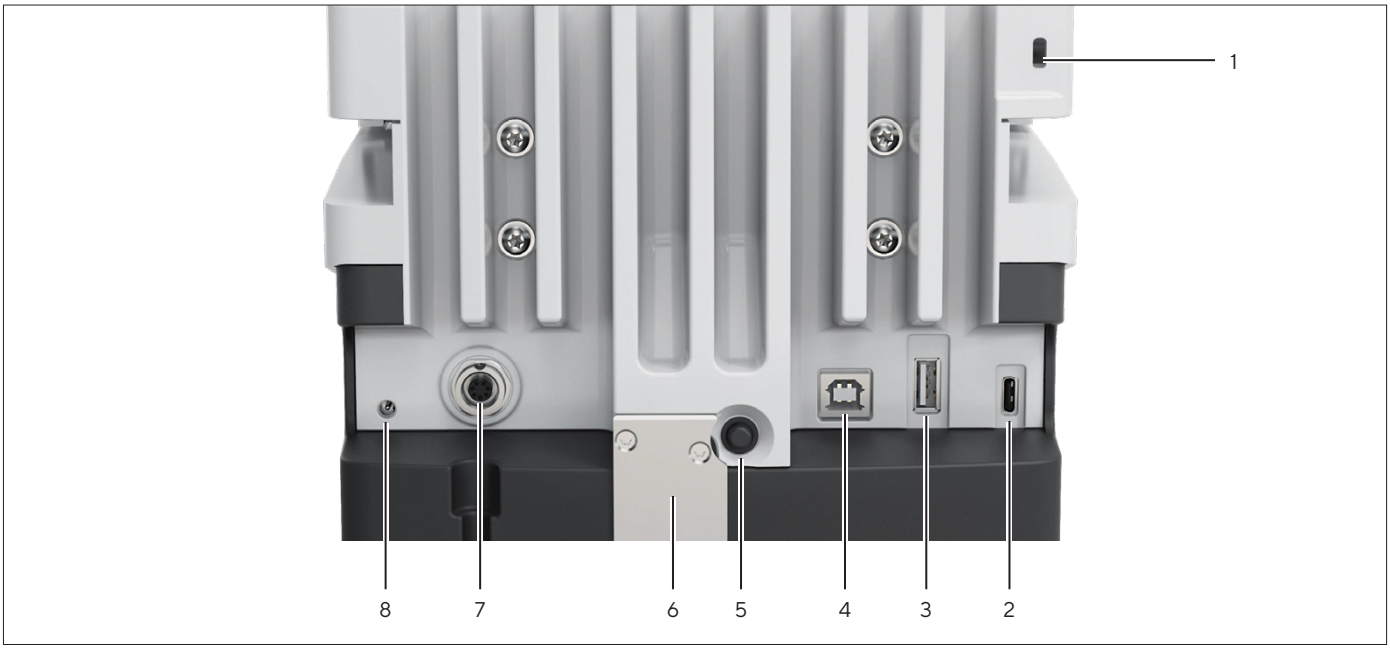


Fig. 5: Connections on the back of the device

Pos.	Name	Description
1	Slot	For attaching a “Kensington” anti-theft device.
2	USB-C connection	For connecting accessories.
3	USB-A connection	For connecting accessories.
4	USB-B connection	For connection to a PC.
5	On key	When using the customized switch-on settings: Switches the device on.
6	Access switch	Protects the device from changes to the device settings.
7	Control unit connection	Receives the connection cable for the control unit.
8	Power supply	

3.4.2 Rear wall of weighing compartment



Fig. 6: Connections on the rear wall of the weighing compartment corresponding bezels removed

Pos.	Name	
1	Connection for an internal climate module	Available as an accessory.
2	Connection for adapter module	For connecting a motorized inner draft shield.

3.4.3 Control Unit

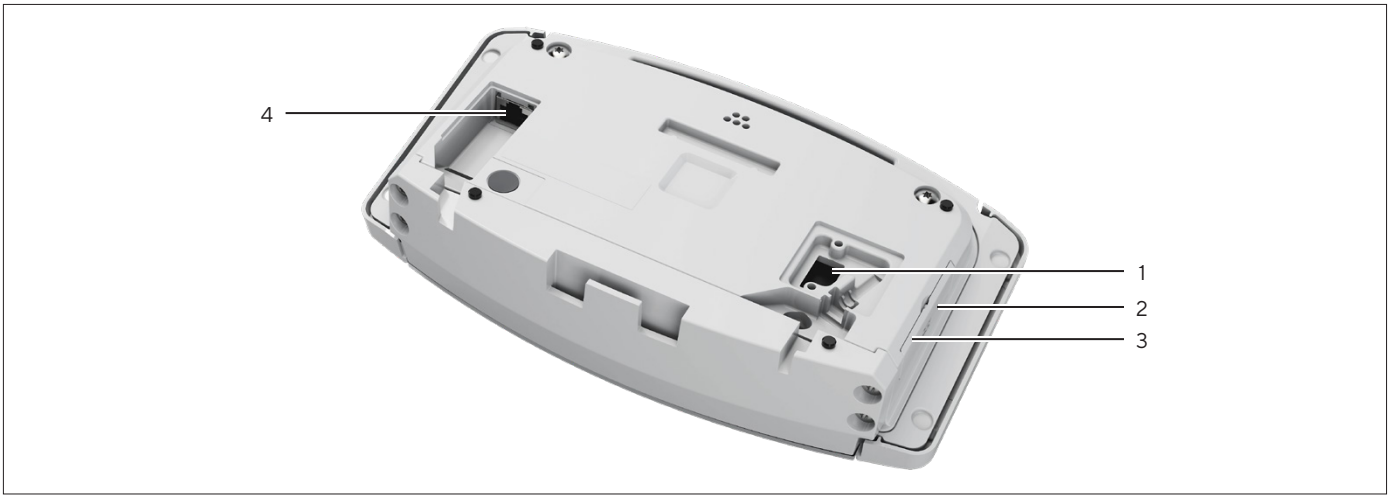


Fig. 7: Connections on the control unit, some covers removed

Pos.	Name	Description
1	Weighing module connection	Receives the connection cable for the control unit.
2	USB-C connection	For connecting accessories (**not** shown).
3	USB-A connection	For connecting accessories (**not** shown).
4	Ethernet connection	For connection to a network.

3.5 Transport Lock (Only for Models MCA1005MC | MCA605MC | MCA1004MC)

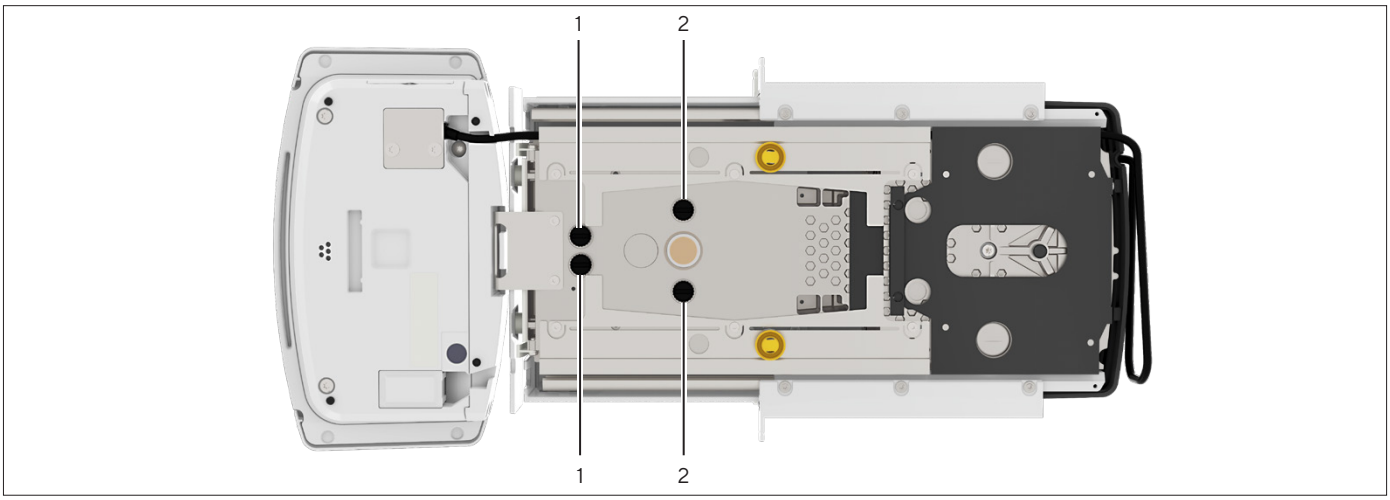


Fig. 8: Transport lock on the underside of the device

Pos.	Name	
1	Short sealing screw	For sealing the weighing chamber during operation.
2	Long safety screw	For securing the weighing system during transport.

3.6 Accessories

Accessories are available for the device. These enable the device to be customized to the specific conditions during weighing procedures, e.g., climate module or weighing table.

3.7 Below-balance Weighing

The device is suitable for below-balance weighing. Samples can be suspended for weighing using below-balance weighing, e.g. samples, which do **not** fit on the weighing pan. Below-balance weighing is possible under the following conditions:

- The device must be set up on a weighing table with recess.
- For below-balance weighing, a below-balance weighing hook must be inserted in the device base (see Chapter “16 Accessories and Spare Parts”, page 93).

3.8 Protective Caps and Covers for Connections

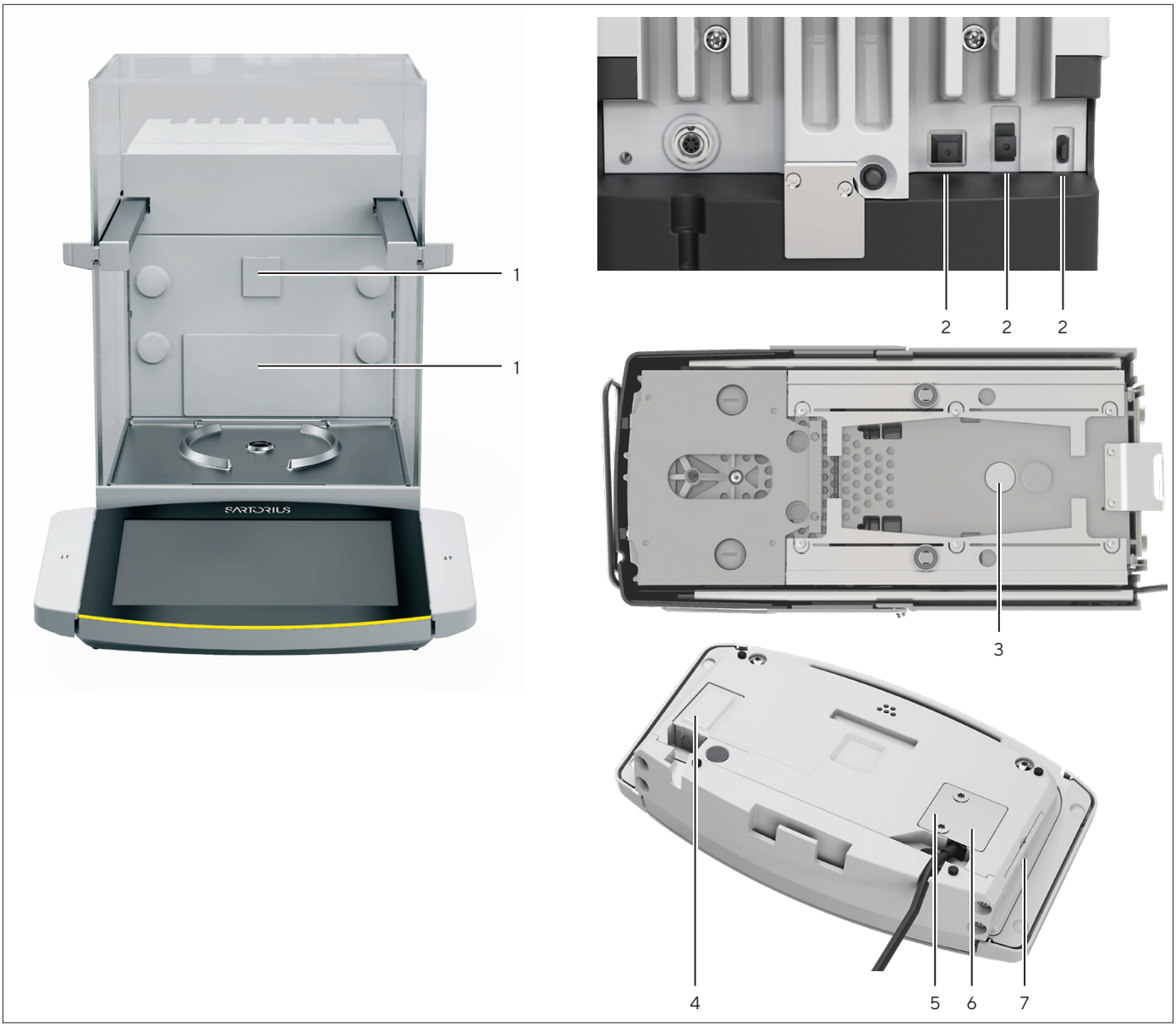


Fig. 9: Protective caps and covers on the rear wall of the weighing compartment, the back of the device, and the control unit

Pos.	Name	Description
1	Cover for rear wall of weighing compartment	Installed in the rear wall of the weighing compartment.
2	Protective cap for USB connection	Removable.
3	Cover for below-balance weighing	Screwed on.
4	Cover for Ethernet connection	Removable.
5	Cover for weighing module connection	Screwed on.
6	Seal	Screwed on **(not** depicted).
7	Cover for USB connection	Removable.

4 Operating Design

4.1 Operating Elements in the Main Menu

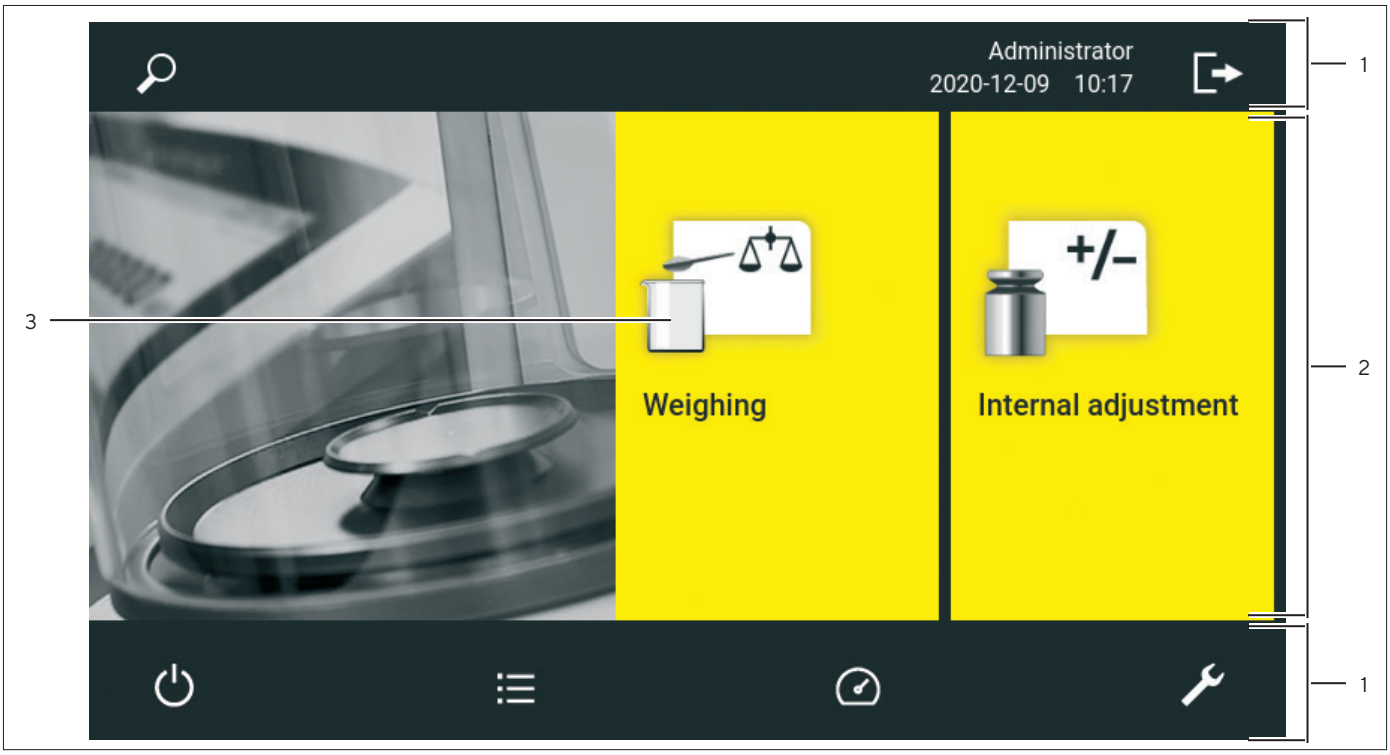


Fig.10: Operating elements in the main menu (example)

Pos.	Name	Description
1	Function bar	– Displays available submenus and operating functions for the current display and current user. – Enables navigation and searching in menus and lists. – In the “Settings” Menu: Displays the name of the menu.
2	Available tasks	Displays all tasks available for the active user.
3	Task	Starts the described task.

4.2 Operating Elements in Task Management

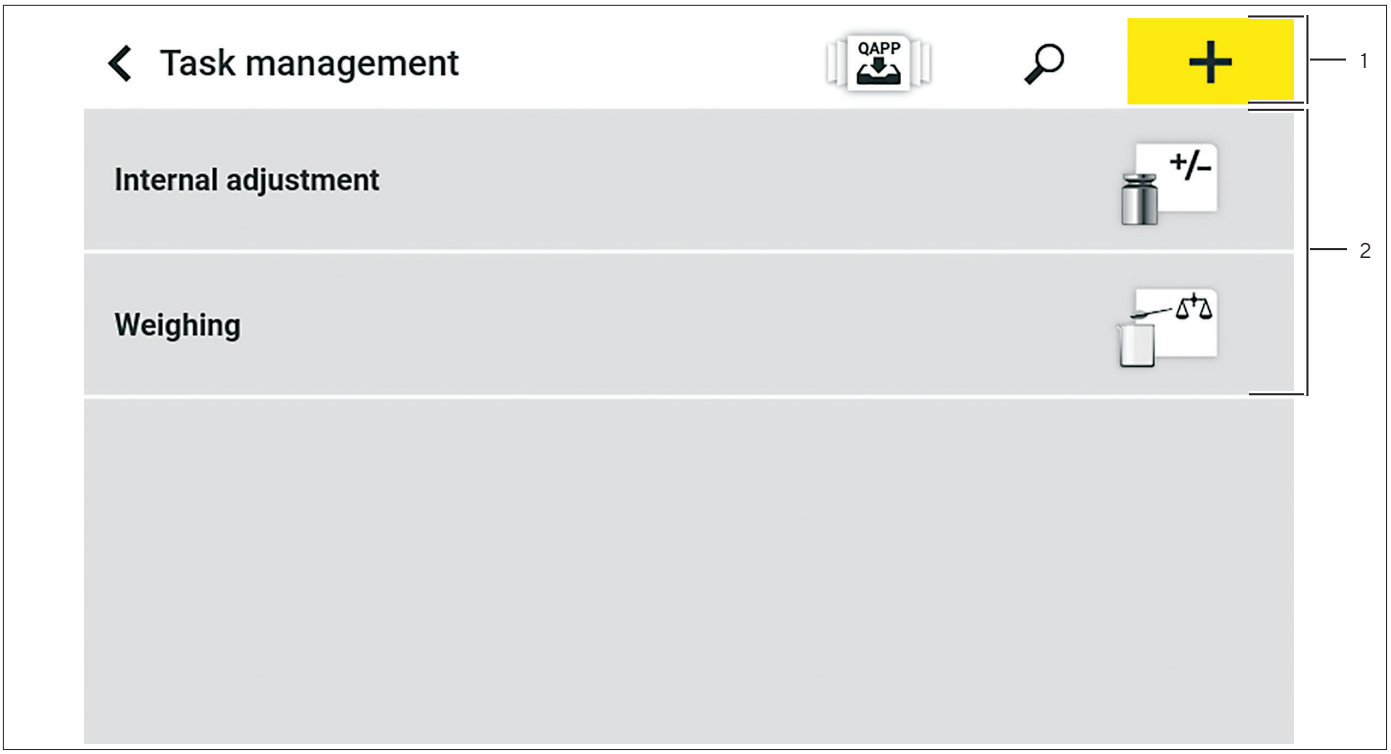


Fig. 11: Operating elements in Task Management (example)

Pos.	Name	Description
1	Function bar	– Enables navigation and searching in menus and lists. – Enables the addition of tasks. – Opens the QAPP Center. – Displays the name of the menu.
2	Available tasks	– Displays all available tasks. – Opens a summary of the properties for the displayed task.

4.3 Operating Elements in the Weighing Display

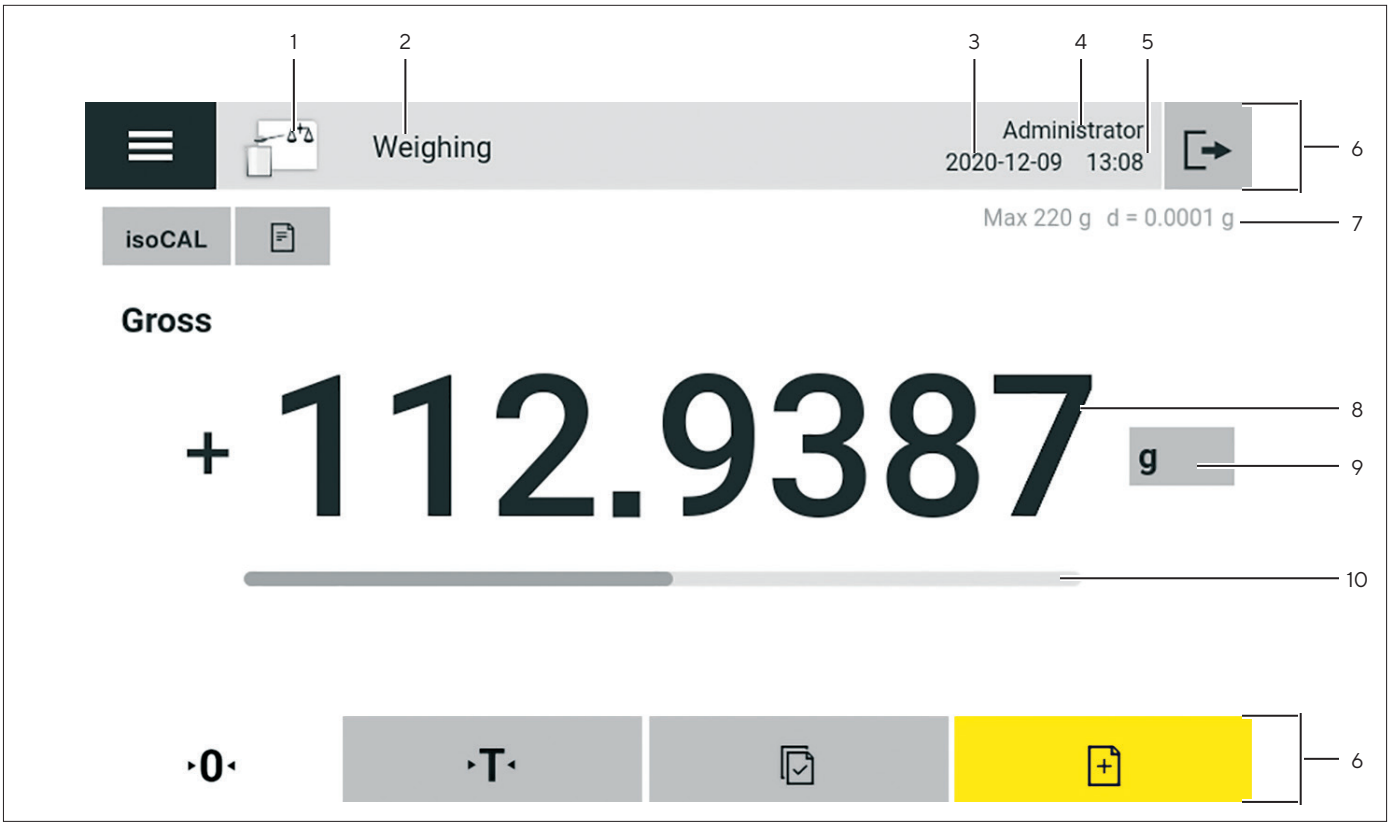


Fig. 12: Weighing display (example)

Pos.	Name	Description
1	Application symbol	Displays the symbol for the active application.
2	Task name	Displays the name of the active task.
3	Date display	Displays the current date.
4	User name	Displays the name of the active user profile.
5	Time display	Displays the current time.
6	Function bar	– Enables navigating the menus. – Displays available operating functions for the current display.
7	Metrological Data	Only for models MCA206MC | MCA1005MC | MCA1004MC: Also displays the variable weighing range with associated scale interval.
8	Weight value display	In the selected unit and resolution
9	Weighing unit	– Displays the selected unit, e.g. grams, [g]. – Enables the unit and resolution to be selected.
10	Capacity display	Displays the measured value as a percentage of weighing capacity utilization.

4.4 User Guide

4.4.1 Messages

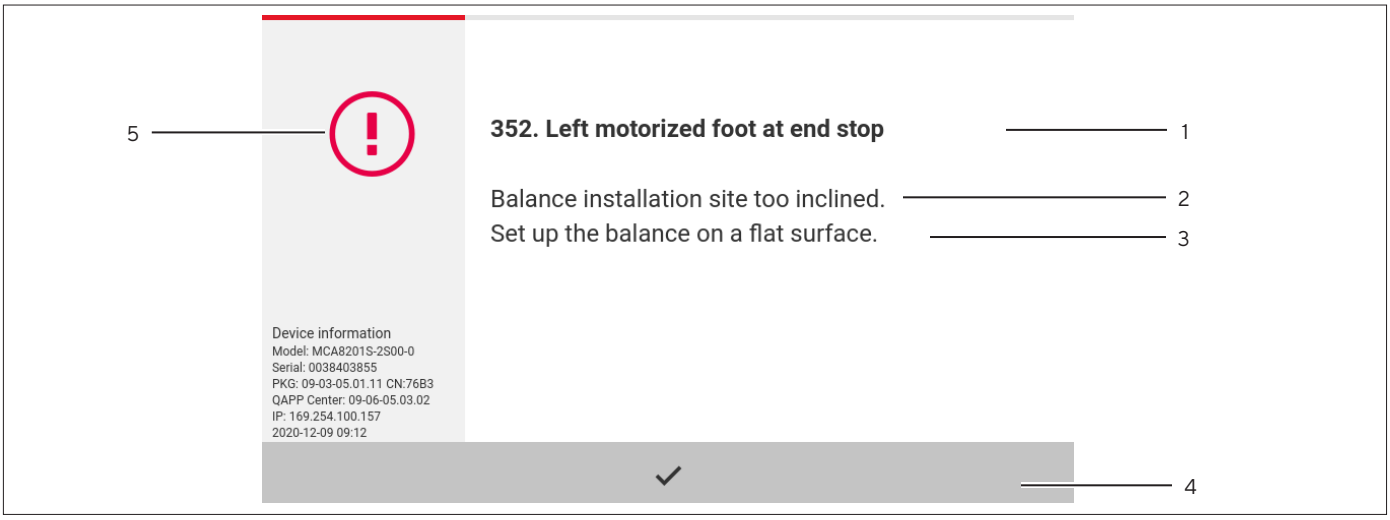


Fig. 13: Error message (example)

Pos.	Name	Description
1	Title	Specifies the type of message, with or without error number.
2	Description	Specifies the cause.
3	Solution	Specifies the measures necessary to eliminate the cause of the message.
4	Confirm	Confirms and closes the message.
5	Message type	Indicates that the message is a status message, warning message, or an error message.

4.4.2 Help Function

For some menus, additional descriptions of menu or parameter contents can be accessed, e.g. printing concept or user management. The help texts have a full text search and are equipped with links. The help texts can be exported.



Fig. 14: Help function in the “Settings” menu (example)

Pos.	Name	Description
1	[Help] button	
2	[Export] button	

4.4.3 Advanced User Guide

Some applications have an advanced user guide. The advanced user guide guides the user through the active task and comprises multiple switchable displays depending on the configured application:

- First display: Shows the step currently being carried out as a graphic display.
- Second display: Shows the current parameters for the task.
- Third display, only for applications with statistics function: Represents statistical information as a bar chart.

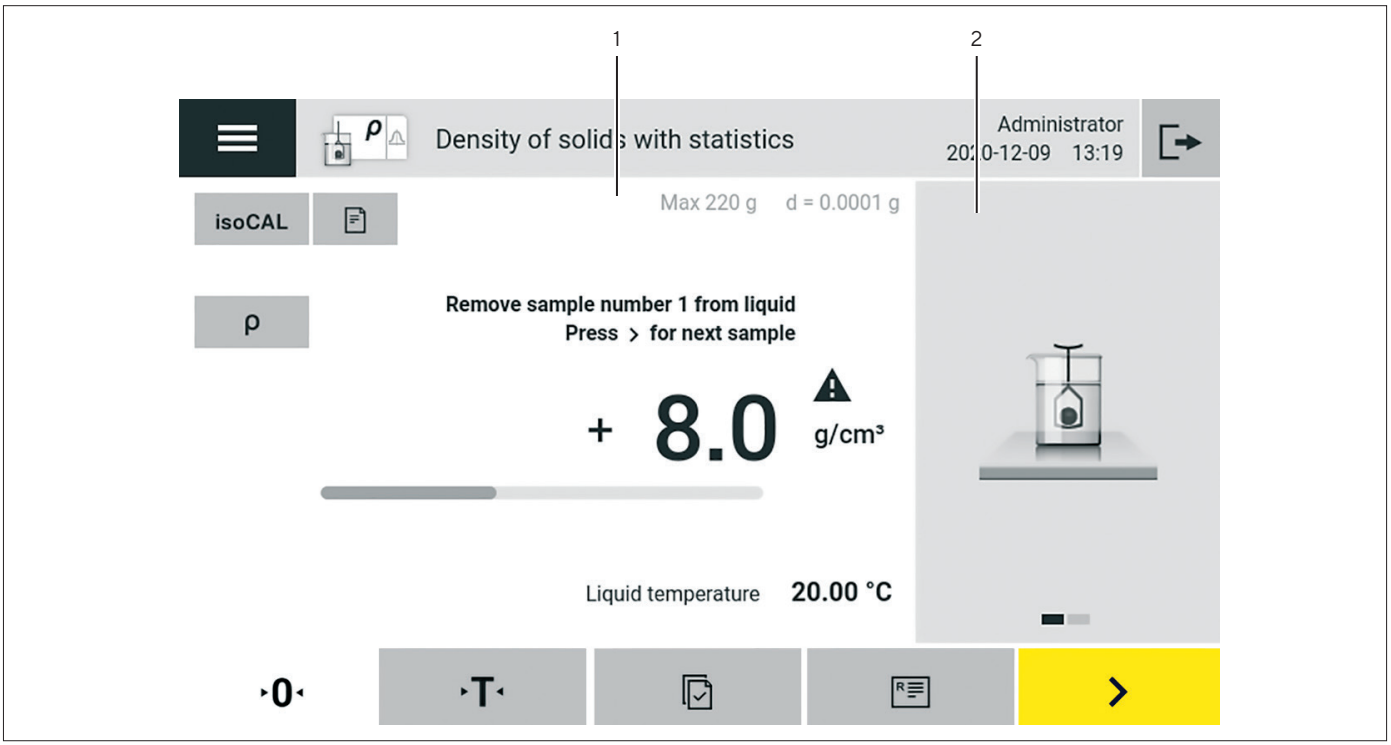


Fig. 15: Advanced user guide (example)

Pos.	Name	Description
1	Weighing display with operator guidance	
2	Advanced User Guide	

4.5 Status Center

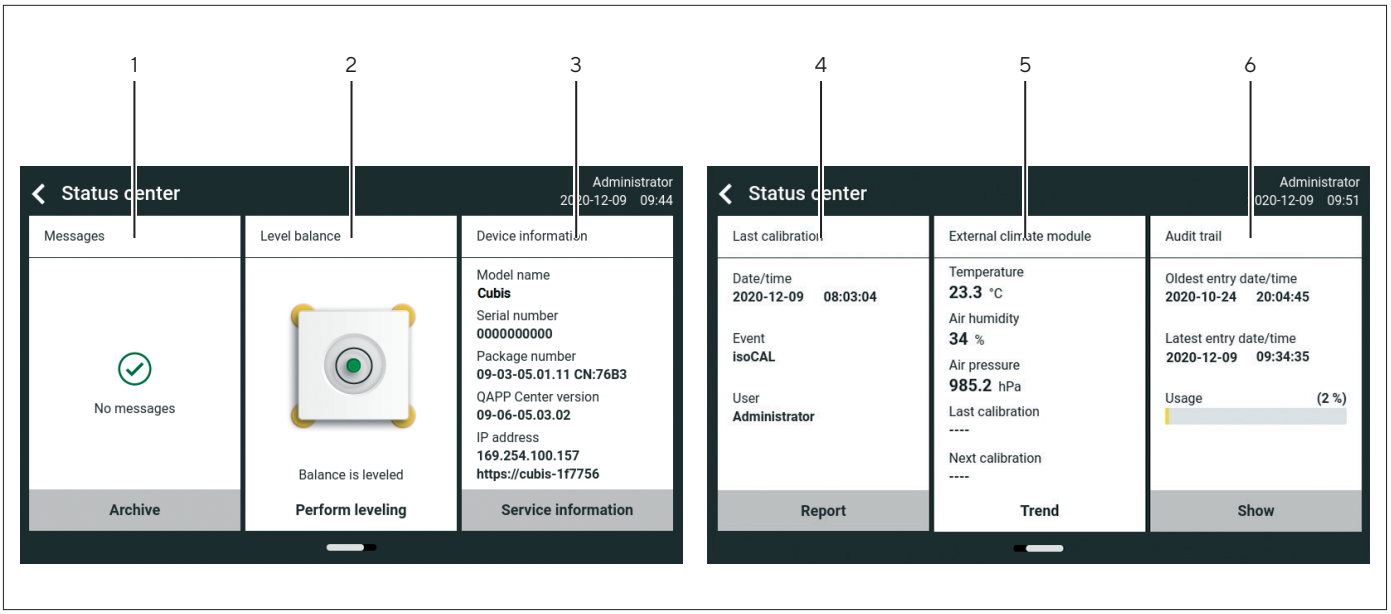


Fig. 16: Status Center (example)

Pos.	Name	Description
1	Messages	Displays information, warning, and error messages.
2	Leveling status	Displays the status of the level.
3	Status for the device	Displays general device information.
4	Calibration and adjustment report	Displays the data for the last adjustment and calibration.
5	Climate data status	Displays the data from the optional climate module.
6	Audit trail	Displays the oldest and newest audit trail entry and the memory usage in percent.
7	Cleaning status	Displays the cleaning status of the device (when using the Cleaning QAPP, **not** shown).

4.6 Keypad

The keypad is intended for entering values in input fields and is adapted to the selected language. If an input field is activated: The alphanumeric keypad or numerical keypad is displayed.

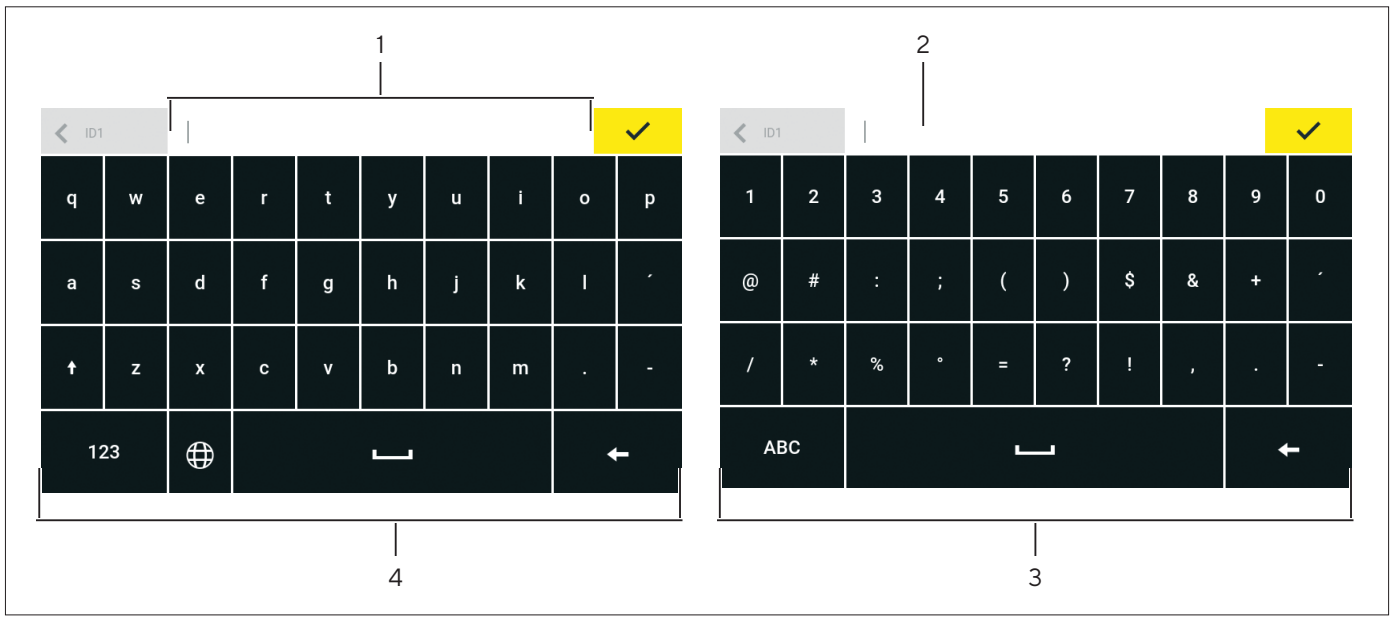


Fig. 17: Alphanumeric keypad and numerical keypad (example)

Pos.	Name	Description
1	Input field	
2	Input assistance	Indicates which values can be entered in the input field, e.g. only numbers.
3	Numeric keypad	
4	Alphanumeric keypad	

4.7 Status Display of the Buttons

Symbol	Name	Description
	Predominant button	Indicates that the function must be executed. The button is color-highlighted.
	Secondary button	Indicates that the function can be executed. The button is highlighted gray.
	Inactive button	Indicates that the function currently **cannot** be executed. The button is grayed out.

4.8 Buttons in the Operating Display

4.8.1 Buttons for Navigation or Organization in Displays

Symbol	Name	Description
	[Menu] button	Quits the active task and opens the main menu.
	[Back] button	– Returns to the previous display. – In the main menu: Accesses the last-performed task.
	[Search] button	Displays options for browsing tasks and list elements.
	[Filter] button	Displays options for filtering tasks and list elements.
	[Sort] button	Displays options for sorting tasks and list elements.
	[Standby] button	Switches the operating display to standby mode.
	[Task Management] button	Opens the task management.
	[Status center] button	Opens the status center.
	[Setup] button	Opens the “Settings” menu.
	[User Logout] button	Logs the currently active user out and accesses the login display.
	[User login] button	Accesses the login display.
	[Next] button	Accesses the next display, e.g., the next processing step for a wizard.
	[Active search] button	Displays an active search for tasks or list elements.
	[Active Filter] button	Displays an active filter for tasks or list elements.
	[Sort Ascending] button	Sorts tasks or list elements in ascending order.
	[Sort descending] button	Sorts tasks or list elements in descending order.
	[10 positions forward] button	Jumps forwards 10 positions in tasks or list elements.

Symbol	Name	Description
10>	[10 Positions Right] button	Jumps 10 positions to the right in tasks or list elements.
Service Information	[Service information] button	Opens the “Settings” / “Device Settings” / “Device Information” / “Service” menu.
Archive	[Status archive] button	Opens an overview of all status messages, warning messages, and error messages.
Perform leveling	[Perform leveling] button	Opens the Leveling Wizard.
	[Warning] button	Opens the list of current warning messages.
	[Error] button	Opens the list of current error messages.
	[Info] button	– Opens a display with additional information about the current menu. – In the status center: Opens the list of current status messages.
	[QAPP Center] button	Opens the QAPP Center.
	[Details] button	Shows additional information about an element.
	[Display Image] button	Displays available images or videos.
	[Close Menu] button	Closes the menu.
	[**Invalid** weight value] button	Opens a display in which all the pending messages are listed.

4.8.2 Buttons for Editing or Managing Entries

Symbol	Name	Description
	[New] button	Starts the wizard for creating new entries, e.g., a new task, or a new weighing profile or user.
	[Activate QAPP] button	Opens the display for activating an application.
	[OK] button	Saves a selection or entry.
	[More] button	Shows the keypad for entering a user-defined value.
	[Select all] button	Selects all elements of a list.

Symbol	Name	Description
	[Deselect all] button	Deselects the selection for all elements of a list.
	[Cancel] button	— Cancels the current process without saving the changed settings or values. — In the display for editing the print memory: Marks the selected value as invalid.
	[Edit] button	Accesses the Adjustment Wizard for editing the displayed element, e.g. for setting the parameters for a task.
	[Delete] button	Deletes the selected element.
	[User Activated] button	Opens an input for deactivating the user.
	[User Deactivated] button	Opens an input for activating the user.
LICENSE	[License] button	Opens the input field for the license key for activating a QAPP package or application.
	[Keypad] button	Shows the keypad.
	[Select character set] button	Switches over the keypad character set for the language version.
▸ PT1	[Show Preset Tare 1 Function] button	Shows the preset tare 1 function.
◂ PT1	[Hide Preset Tare 1 Function] button	Hides the preset tare 1 function.
	[Edit password] button	Opens the input field to change a user's password.
	[Deactivate entry] button	Deactivates an entry.
	[Activate entry] button	Activates an entry.

4.8.3 Buttons for Weighing, Printing, and Export Functions

Symbol	Name	Description
	[Leveling] button	Opens the Leveling Wizard.
isoCAL	[isoCAL] button	Starts the isoCAL function.
Prc	[Result] button	Toggles between the result and weight value display for the current application, e.g. weighing in percent.
	[Zero] button	Starts zeroing.
	[Tare] button	Starts taring.
	[Tare 1] button	Stores the current weight value in the tare 1 memory.
Tare 1 0.00 g	[Tare 1 Entry] button	Opens an input field for manually entering the tare 1 value.
	[Delete Tare 1] button	Deletes the tare 1 memory.
	[Start] button	Starts the selected application.
	[End] button	Quits the active application and opens the Print Memory display.
	[Confirm] button	Confirms the current display and initializes the next step.
	[Save] button	Saves the weight value and sends it to the print memory.
	[Print memory] button	Opens the Print Memory display.
	[Print] button	If the Print Memory is displayed: Exports the saved print orders via the integrated data interfaces.
	[Connect] button	Checks the connection to a connected printer and prints a test page.
	[Connector] button	Checks the connection to the respective connector.
g	[Unit change] button	– Toggles between the weight units and the resolution. – Accesses the “Unit change” function menu.
	[Restart] button	If an application is active: Deletes the saved values and re-starts the application.
	[Report] button	If an application is active, e.g. density determination: Displays a report about the progress of the application.
	[Result report] button	If an application is active, e.g. density determination: Displays a report about the result of the application.

Symbol	Name	Description
	[Export] button	– Creates files depending on the context. – Releases the exported files, e.g., PDF or HTmL files, via a connector.
	[Help] button	Displays descriptions of the menu and parameter contents, e.g., weighing characteristics and weighing functions, printing concept, user management

4.9 Displays in the Operating Display

Symbol	Name	Description
	[Leveling] display	Indicates that the device is **not** leveled.
	[Print memory] display	Indicates that elements are located in the print memory.
g	[Unit symbol] display	Indicates the set weight unit, e.g. [g] for “grams”.
Gross	[Result] display	Indicates whether the weight value being displayed is a gross value or the result of an application.
	[**Invalid** weight value] display	– Together with the unit indicates: – The display shows the calculated weight value of an application, e.g. for the “Totalizing” application. – The displayed weight value is less than the minimum initial weight. – For conformity-assessed devices without unit display: – Negative gross value – The device has **not** yet reached the necessary operating temperature to accurately determine the weight. – The device is **not** leveled. – The device requires the isoCAL function to be executed.
	[Sign] display	Indicates whether the value being displayed is positive or negative.
	[Current page] display	– Indicates that the current display is a multi-page display. – Indicates which page of a multi-page display is being displayed currently.
	[Copy] button	Indicates that data are being copied.
	[Selection] display	Indicates that an entry is selected in a list.
	[Target value] display	Shows the target value in the capacity display.

4.10 Applications and Tasks

QAPP applications (applications) are grouped together in QAPP packages. The device is supplied with some freely accessible applications in the “Essentials” QAPP package. These applications can be used to carry out the most important functions, such as weighing and calibration.

Other applications may be activated in the QAPP Center for a fee. The applications contain additional or specific functions, such as user management or specific weighing functions.


An application must be configured as a task before it can be used. Specific settings need to be applied for this, with the help of the assistant function. Tasks are visible to all users who have the required role to perform the task.

4.11 Licensable Device Functions

Some device functions are licensable. In order for these device functions to be usable, they must be activated in the QAPP Center for a fee. The device functions may be activated upon delivery or they can be activated later on in the QAPP Center (for licensable device functions, see Chapter 15.7, page 87).

4.12 Navigating the Menus

Procedure

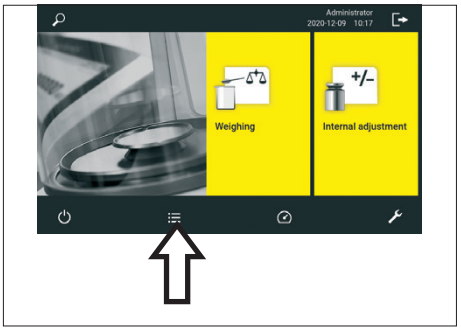


▶ To open a menu: Press the button for the preferred menu, e.g., the [Task management] button.


▷ The menu is opened and the name of the open menu is opened in the function bar.




▶ To return to the main menu from other displays: Press the [Menu] button or press the [Back] button (multiple times) until the main menu is opened.

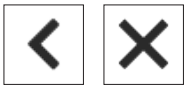


▶ To exit a display: Press the [Back] or [Cancel] button.

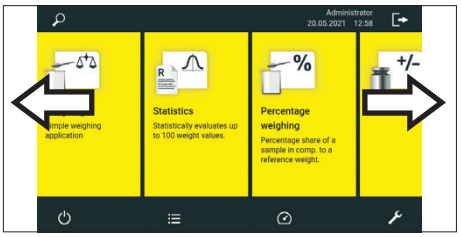

▶ To scroll through the available tasks: Swipe the task bar to the left or right.

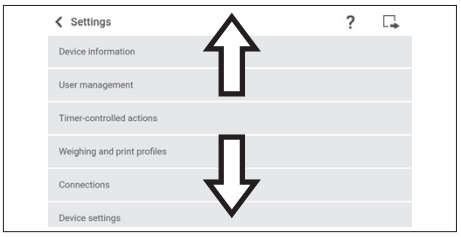


▶ TO SCROLL THROUGH THE LISTS IN A MENU: SWIPE THE LIST DOWNWARDS OR UPWARDS.

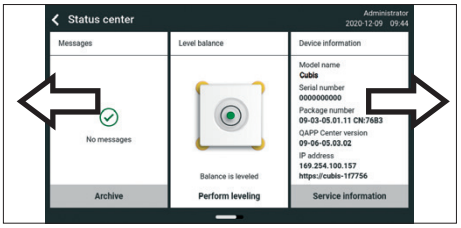


▶ To display the next or previous page in a display containing several pages, e.g. in the Status center: Swipe the display to the left or right.

- ▶ If a value needs to be selected from a list:
  - ▶ Scroll to the preferred value in the display. To do this, swipe the display upwards or downwards.
  - ▶ Press the preferred value.
  - ▶ To confirm the selection: Press the [OK] button.
- ▷ The selected value is saved and the list closes.

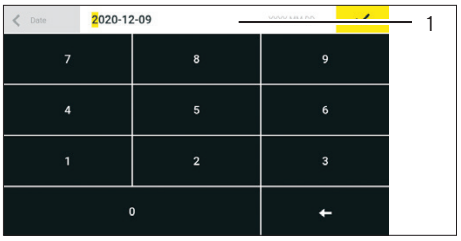

- ▶ If elements from a display need to be filtered or a display needs to be browsed:
  - ▶ Press the [Search] or [Filter] button.
  - ▷ The keypad is displayed.
  - ▶ Type the search value or value to be filtered into the input field (1) using the keypad.
  - ▶ Press the [OK] button.
- ▶ To close the input field for searching and filtering without starting a search or filter operation: Do **not** enter a value in the input field or delete the entered value.
  - ▶ Press the [OK] button.

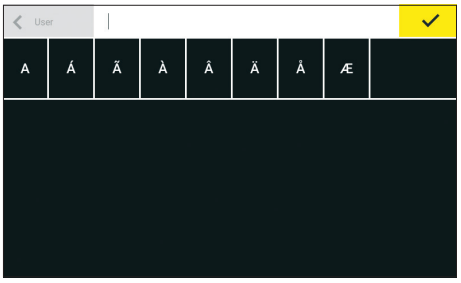

- ▶ If language-specific characters need to be entered using the keypad:
  - ▶ Press and hold a letter on the keypad.
  - ▷ If language-specific characters are available for the letter being pressed: A display opens containing all the language-specific characters available for the letter being pressed.
  - ▶ To select a language-specific character and return to the keypad display: Press the preferred language-specific character.

4.13 Menu Structure

4.13.1 Main Menu

▶ Navigating the menus (see Chapter 4.12, page 34).

Level 1	Level 2	Level 3	Description
Task management			Displays all available tasks. Opens a summary of the properties for the displayed task.
	QAPP Center	QAPP package	– Lists the QAPPs contained in a QAPP package. – Displays the licensing status of the QAPPs. – Leads to the input screen for the licensing of the QAPP package.
Status center			Display information about the device status, e.g. level status, amount of storage or messages.
Settings	Device information		
	User management		
	Access management*		
	Timer-controlled actions		
	Weighing and Print Profiles		
	Connections		
	Device settings		
	Device maintenance		

* If a related QAPP is activated: The menu is opened.

4.13.2 “Settings” Menu

▶ Navigating the menus (see Chapter 4.12, page 34).

Level 1	Level 2	Level 3	Description
Device Information	General device information		Displays general device information, e.g. model name, serial number, QAPP Center version.
	Service		Displays general information on service-related topics, e.g. telephone number, last maintenance, or qualification data.
	View data storage device		Display, sort, browse, or export the contents of the data storage device.
	View audit trail		If an audit trail license has been purchased: Display, sort, browse, or export the contents of the audit trail memory.
	Software licenses		Display the list of all open source software modules used.
User management			Display or edit the 4 standard user profiles, e.g., password, language. Create or edit additional user profiles, user roles, and rules*.
Access management*	Role management		Display or edit the 4 standard user roles. Create or edit additional user roles.
	Rules		Define the rules for failed attempts at password entry.
	Local password rules		Define the permissible characters, length, and validity time for the password.
	LDAP server configuration		Configure the data for device communication with an LDAP server.
Timer-controlled actions			Display the created timer-controlled actions, e.g., “Start task”.
Weighing and print profiles	Weighing		Display and manage the saved weighing profiles.
	YDP30		Display and manage the saved print profiles, e.g. print profiles for the printer YDP30 or print profiles for PDF or CSV.
	PDF		
	CSV		
	PC direct		
	SBI direct		
	Web service report		
	Data file		

* If a related QAPP is activated: The menu is opened.

Level 1	Level 2	Level 3	Description
Connections	Connectors	USB	Display and manage the saved connectors, e.g., a USB stick.
		YDP30-NET	
		FTP	If this extension has been configured in the QAPP Center: The connection to the FTP or FTPS servers facilitates file transfer using the FTP or FTPS protocol.
		FTPS	
		Network printer	Display and manage the saved network printers.
		SMB	If the "Windows File Server" extension has been configured in the QAPP Center: The connection to a Windows file server facilitates file transfer using a Server Message Block (SMB protocol).
		HTTP	If the "Windows File Server" extension has been configured in the QAPP Center: Enables data connection to a server with service end points (web services, ReST).
		HTTPS	
		HTTP (dynamic)	
		HTTPS (dynamic)	
	Network	General settings	Define the device host name.
		Ethernet	Display and manage settings for the device's Ethernet or Wi-Fi connection.
		Wi-Fi	
		Diagnostic functions	Performs a network test to check the connection.
	Website / web services	Website access	Determines the settings for the display of the website for the device. Depending on the settings for the website access, a user log-in may be required. The user has the same access rights for the website as on the device.
		Web service password	Define the password for the web services.
	Certificates	Trusted certification authorities	Display and export the trusted certification authorities (CA).
		Device certificates	Display and export the device certificates.
		Trust unknown certificate authorities	Specifies the security for unknown certification authorities.
	Interfaces	Serial communication over network	Display and edit the profile for the Ethernet connection.
		COM-RS232 interface	Display and edit the profile for the RS232 connection.
		USB-B interface	Display and edit the profile for the USB-B connection.

* If a related QAPP is activated: The menu is opened.

Level 1	Level 2	Level 3	Description
Connections	SBI protocol	Format	Configure the settings for the data output and data output format.
		Output	Define whether the output takes place with or without stability.
		Automatic output	Activate or deactivate the output rate for automatic data output.
		Display SBI error messages	Activates or deactivates the error messages for the SBI protocol.
		Formatting output values	Defines the line length for the data output.
	Climate modules		View the settings for the connected climate module, e.g. the internal climate module.
	Connected devices	Motion sensor	If a motion sensor is connected to the device: Define the number and functions of the gestures.
		External USB keys	If an external USB peripheral device is connected to the device, e.g. a keyboard or foot switch: Define the functions of the USB keys. In this way, key functions for the device can be triggered via the external USB peripheral device.

* If a related QAPP is activated: The menu is opened.

Level 1	Level 2	Level 3	Description
Device settings	Date and time	NTP configuration	Activate or deactivate time synchronization via NTP.
		Set date and time	Enter the date and time. Determine the time zone.
	Safe Weighing	isoCAL execution mode	Set the isoCAL function.
		isoCAL limit rule	Define rules for triggering the isoCAL function.
		User-defined warning limit	Specify the warning limit for the maximum percentage deviation compared to the reference weight. When the warning limit is reached, a message is displayed.
		Safety level	Set the alert level for the message indicating that the isoCAL function must be executed.
	Preload	Show current preload	Displays the currently set preload.
		Set preload	Sets a preload for technical equipment that needs to be constantly on the device. In this case the maximum capacity is reduced by the set preload.
		Delete preload	Deletes the set preload.
	Calibration/adjustment report	Automatic printout for isoCAL	Activates the data output for the isoCAL function.
		Automatic printout for cal./adj. Task	Activates the data output for the calibration and adjustment function.
		Print profile	Defines the output device on which a printout must be carried out.
	Minimum weight*		Functions for monitoring the permissible operating range.
	Measurement uncertainty*		Functions for dynamic display of the measurement uncertainty for the weight value.
	Device ID		Determine the device identifier.
	Electronic signature*	Sign report	Activate or deactivate the electronic signature.

* If a related QAPP is activated: The menu is opened.

Level 1	Level 2	Level 3	Description
Device settings	Draft shield*	Left/right palm key	Set the function of the left and right palm key.
		Automatic mode	If enabled: Closes the draft shield, carries out the function, and opens the draft shield again.
		Sensor mode*	Activates or deactivates the proximity sensors.
		Left proximity sensor*	Set the sensitivity of the left proximity sensor.
		Right proximity sensor*	Set the sensitivity of the right proximity sensor.
		Opening angle for left button*	Set the opening angle of the motorized inner draft shield for the left palm key.
		Opening angle for right button*	Set the opening angle of the motorized inner draft shield for the right palm key.
	Start-up behavior	Initial zeroing/taring	Activate or deactivate automatic zeroing and taring when starting the device.
		Automatic login	Activate or deactivate automatic log-in of the last user when starting the device.
		Service user identification	Activates a user identification for the service.
		Automatic start of task	Select or deactivate automatic start of task when staring the device. All tasks that are configured on the device can be selected.
	Device cleaning*	Cleaning support	Functions to support device cleaning

* When the following conditions have been met: The menu or parts of the menu are open:
— The motorized draft shield is activated.
— The device is fitted with a proximity sensor.
— A motorized inner draft shield is mounted.

Level 1	Level 2	Level 3	Description
Device settings	Display settings	Display brightness	Define the brightness of the operating display.
		Power button	Select the “Energy saving” or “Switch-off” function.
		Color scheme*	Define the settings for the operating display color scheme.
	Sound (Loudspeaker)	Touch sounds	Activate or deactivate the acoustic signal when operating a button.
		Message sounds	Activate or deactivate the acoustic signal for messages.
		Sound for end-of-action	Activate or deactivate the acoustic signal for the end of an action.
Device maintenance	Update firmware	List of connectors and then list of firmware	— Starts a firmware update. — **Not** for conformity-assessed devices.
	Install QAPP Center	List of connectors and then list of QAPP Centers	Updates the QAPP Center using the selected connector.
	Install customized QAPP	List of connectors and then list of customized QAPPs	Updates a customized QAPP using the selected connector.
	Backup device settings	List of connectors	Saves the device data on the selected connector.
	Restore device settings	List of connectors	Loads saved device data from the selected connector.
	Export device settings	List of device settings	Save the preferred settings on a selected connector.
	Import device settings	List of connectors for import	Load selected settings, e.g. from a server or connected USB stick.
	Languages	Import language	Loads the language file from the selected connector.
		Delete imported languages	
	Reset to factory settings	Restore factory settings	
	Create diagnostic file		

* If a related QAPP is activated: The menu is opened.

5 Installation

5.1 Scope of Delivery

Item	Quantity
Device	1
Base plate of weighing compartment	1
Models MCA1005MC | MCA605MC | MCA1004MC:	
90 mm weighing pan	1
Glass insert to reduce height YDSHR	1
Models MCA206MC | MC106MC | MCA66MC | MCA36MC:	
50 mm weighing pan	1
Shield plate for 50 mm weighing pan	1
Adapter module for inner draft shield YDS125A	1
Glass cylinder	1
Glass lid	1
Power supply unit	1
Country-specific power supply cable with test seal	1
USB connection cable	1
Control unit for Cubis® MCA	1
Connection cable for control unit, mounted on delivery	1
Draft shield	
Front panel	1
Top door	1
Left door	1
Right door	1
Rear panel, installed on delivery	1
Palm key for motorized draft shield, mounted on the control unit on delivery	2
In-use dust cover for display and control unit	1
Models MCA1005MC | MCA605MC | MCA1004MC:	
Transport lock	1
Operating Instructions	1
Specific certificates for the device	1
‘Cleaning case’ with cleaning utensils, e.g. brush for cleaning	1

5.2 Selecting an Installation Site

Procedure

- ▶ Ensure that the installation conditions have been met (see Chapter “15.3 Installation Conditions”, page 84).
- ▶ **NOTICE** Risk of damage to the power supply unit from argon! Observe the instructions for using argon (see Chapter “15.3 Installation Conditions”, page 84).


5.3 Unpacking

The device is packed in multilayer foam packaging. The device components, e.g. the weighing pan, are inserted into the second and third layers of the foam packaging.

The individual layers of the foam packaging must be removed one after the other.

Procedure

- ▶ Open the package.
- ▶ Undo the tension belt.
- ▶ Remove the top layer of foam packaging.
- ▶ Remove the second and third layers of foam packaging.
- ▶ Grip the device by the grooves (1) on the left and right-hand side and lift the device out of the bottom layer of foam packaging.
- ▶ Place the device down on a stable, completely flat surface.
- ▶ Sartorius recommends keeping the original packaging to return the device appropriately, e.g., for repairs.

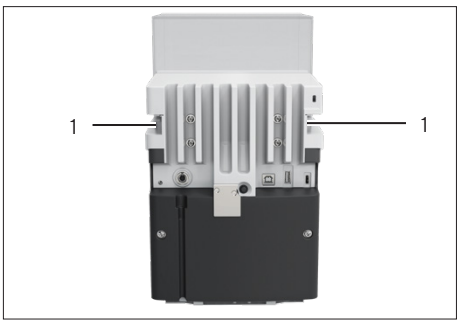

5.4 Attaching or Removing the Control Unit

Procedure

- ▶ In order to attach the control unit to the device: Hook the control unit into the control unit holder (1). The control unit must be positioned entirely on the base.
- ▶ In order to remove the control unit from the device: Pull the control unit up off of the control unit holder.



5.5 Placing the Device on its Side and Installing

For some installation work, the device has to be placed on its side, for example when inserting connection cables.

Material: 1 soft support base, for placing down the device

Requirements


- There are **no** components installed in the pan retainer.
- All outer draft shield doors are removed: Front panel, upper panel, side panels.
- Only for models MCA1005MC | MCA605MC | MCA1004MC: All safety screws of the transport lock have been mounted. The transport lock is in place (see Chapter 12.1, page 79).

⚠ CAUTION

Risk of injury when lifting or transporting!

- ▶ Use both hands when transporting the device and setting it down. To do this, reach sideways under the device with both hands in the rear part of the device.


Procedure

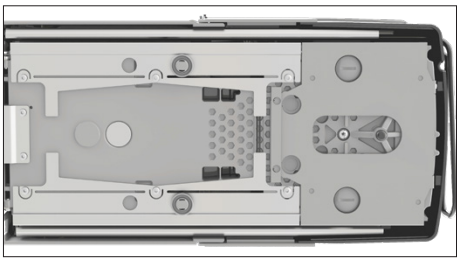

- ▶ Remove the control unit from the device.
- ▶ If the device is to be placed on its side:
  - ▶ Reach sideways under the device with both hands in the rear part of the device.
  - ▶ Turn the device on its side and place it on the soft support base.
- ▶ If the device is to be returned to the upright position:
  - ▶ Reach sideways under the device with both hands in the rear part of the device.
  - ▶ Place the device back onto its base.

5.6 Inserting the Connection Cable for the Control Unit

Procedure

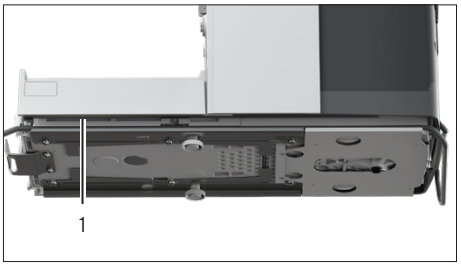


- ▶ Place the device on its side (see Chapter 5.5, page 45).
- ▶ Fully press the control unit connection cable into the cable channel (1) along the side of the device. The control unit connection cable must **not** come into contact with the side panels as this can cause device malfunctions.
- ▶ Place the device back onto its base.

5.7 Inserting the Connection Cable for Ethernet

An Ethernet connection cable can be attached to the device. The Ethernet connection cable must be inserted into the cable channel on the side of the device.

Material: 1 Ethernet connection cable

Procedure

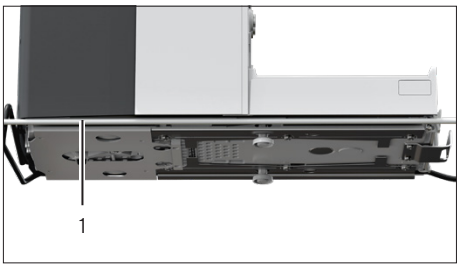

- ▶ Place the device on its side (see Chapter 5.5, page 45).
- ▶ Fully press the Ethernet connection cable into the cable channel (1) along the side of the device. The Ethernet connection cable must **not** come into contact with the side panels as this can cause device malfunctions.
- ▶ Place the device back onto its base.

5.8 Positioning the Weighing Pan and Associated Components

5.8.1 Devices with scale interval 1 µg

Requirements

– The front panel of the draft shield is removed.

Procedure




- ▶ Remove the bottom cover for the rear wall of the weighing compartment.
- ▶ Mount motorized inner draft shield in the weighing chamber, see enclosed Installation Instructions for the motorized inner draft shield YD-S125A.

5.8.2 Devices with scale interval 0.01 mg and 0.1 mg


Requirements

– The front panel of the draft shield is removed.

Procedure

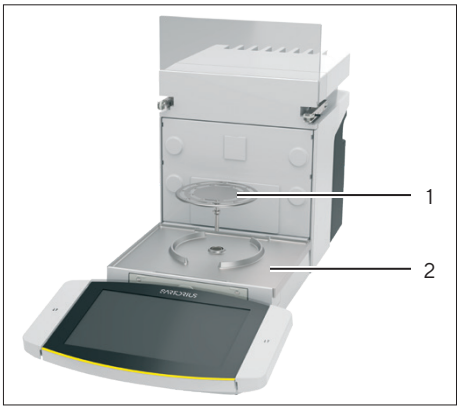


▶ Insert the base plate of the weighing compartment (2) into the device.
▶ Insert the weighing pan (1) into the pan retainer.
▶ If the weighing pan needs to be turned, e.g. with sample holders:
  ▶ Lift the weighing pan and turn it to the left or right.
  ▷ The weighing pan snaps into the next position of the pan retainer. This avoids accidental turning of the weighing pan.

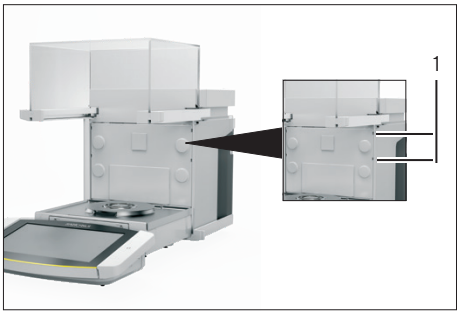


▶ Install the glass insert to reduce height in the weighing chamber, see enclosed Installation Instructions YDSHR.

5.9 Installing the Draft Shield


Procedure

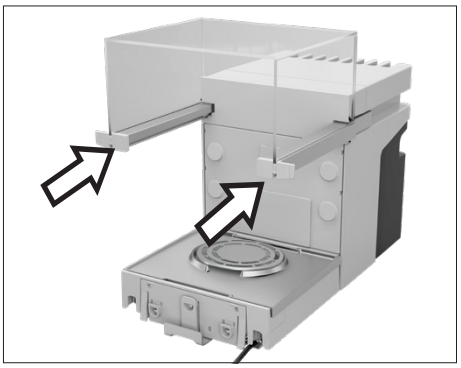

▶ Remove the control unit from the device.
▶ Install the upper door. To do this, insert the two door panel mounts into the two guide rods on the side of the device and push backwards.

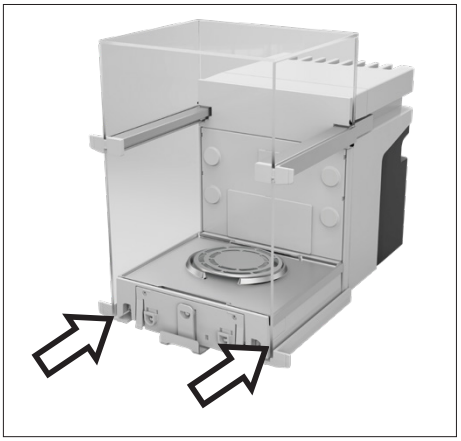

▶ Install the right and left doors. To do this, insert the door panel mounts, one after the other, into the guide rods at the bottom on the side of the device, and push backwards.

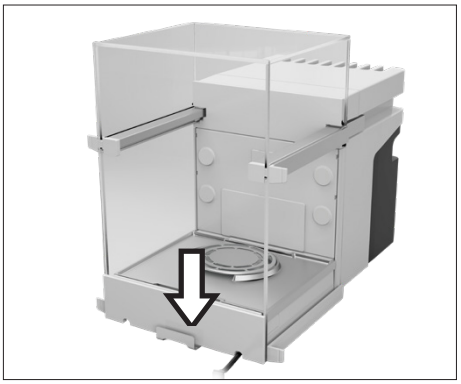


▶ Insert the panel mount for the front panel into the recesses in the front of the device.

5.10 Removing the Transport Lock (Only for Models MCA1005MC | MCA605MC | MCA1004MC)

Removing the transport lock will release the weighing system and mean that the device can be used.


Procedure

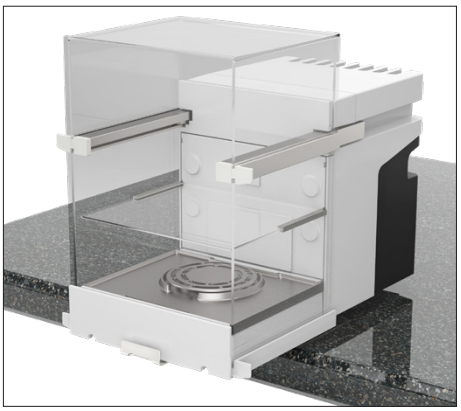

▶ Remove the control unit from the device.
▶ Push the front of the device over the edge of the table so that the safety and sealing screws on the underside of the device are accessible.



- ▶ Unscrew the short sealing screws (2).
- ▶ Unscrew the longer safety screws (1).

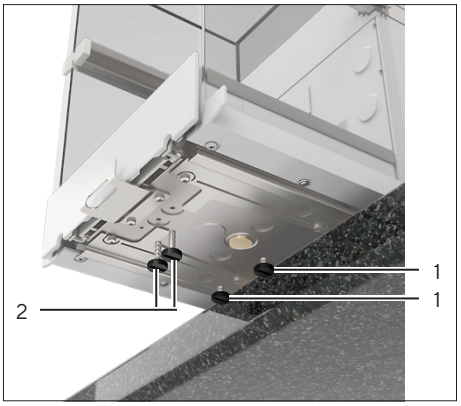

- ▶ Change over the safety and sealing screws for storage:
  - ▶ Screw the short sealing screws (1) into the rear position.
  - ▶ Screw the longer safety screws (2) into the front position.
- ▶ Carefully push the device back so that it is fully on the weighing table.

5.11 Installing the Control Unit

The control unit can be installed in front of or next to the device.

Procedure

- ▶ Remove the control unit from the device.
- ▶ Place the control unit in the preferred position (for the dimensions for positioning the control unit, see Chapter “15.1 Dimensions and Weights”, page 83). The control unit must be positioned entirely on the base.

5.12 Acclimatizing

When a cold device is brought into a warm environment: The temperature difference can lead to condensation from humidity in the device (moisture formation). Moisture in the device can lead to malfunctions.


Acclimatize the weights to the temperature of the load cell. This is the only way to avoid measurement errors caused by air buoyancy and fluctuations resulting from convection currents across the surface of the sample.

Procedure

- ▶ Allow the device to acclimatize at the installation site (for the acclimatization duration see Chapter “15.4 Acclimatization Before Power Supply”, page 85). The device must be disconnected from the power supply during that time.

6 Commissioning

6.1 Connecting the Connection Cable for the Control Unit

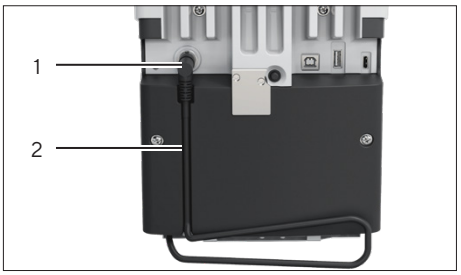

Procedure

- ▶ Insert the connection cable for the control unit into the recess in the rear side of the device (2).
- ▶ Plug the connection cable for the control unit into the “Control unit connection” (1) and tighten by hand.


6.2 Attaching the Connection Cable for Ethernet



Procedure


- ▶ Remove the control unit from the device.
- ▶ Remove the cover for the Ethernet connection from the control unit.
- ▶ Connect the Ethernet connection cable to the Ethernet connection (1).
- ▶ Plug the other end of the Ethernet connection cable into the Ethernet connection at the installation site.

6.3 Connecting the Power Supply Unit to the Device

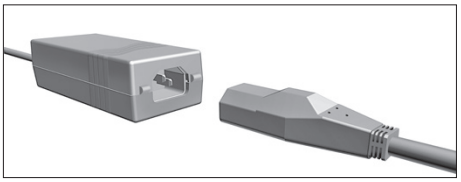

Procedure

- ▶ Plug the power supply cable into the power supply unit connection.

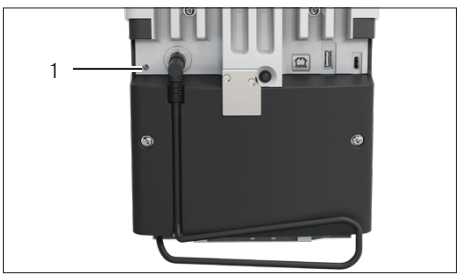

- ▶ Connect the connection cable from the power supply unit to the “Power supply” connection (1) on the back of the device.

6.4 Connecting the Power Supply

Procedure

- ▶ Check whether the country-specific power plug matches the power supplies at the installation site.
  - ▶ If required: Contact Sartorius Service.
- ▶ Check whether the voltage specifications on the manufacturer's ID label match those of the power supply at the installation site.
  - ▶ If the input voltage is too high or too low: Do **not** connect the device to the power supply.
  - ▶ Contact Sartorius Service.
- ▶ Connect the device to the power supply at the installation site. To do this, connect the power plug from the power supply cable into the power socket.
- ▷ The device is switched on and performs the initial functions for device startup.


6.5 Connecting Accessories

Accessories can be connected to the device. For certain accessory parts, the covers must be removed in the wall of the weighing compartment, e.g. for connecting a sample holder module.

Requirements

The accessories must be suitable for the device (see instructions for the accessories).

Procedure

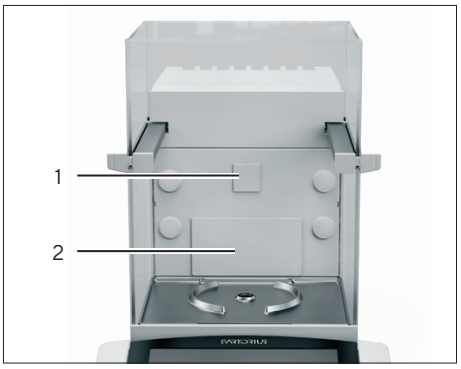

- ▶ If a cover for the rear wall of weighing compartment needs to be removed from the device: Remove the top (1) and bottom (2) covers for the wall of weighing compartment from the device.
- ▶ Connect the accessories to the appropriate connections on the device (for connections for accessory parts, see instructions for accessories).

6.6 Attaching Protective Caps and Covers

If connections of the device are **not** being used during operation: We recommend using the protective caps and covers provided to seal the connections on the rear wall of the weighing compartment and on the back of the device.

Procedure

- ▶ Check whether all unused connections have been sealed.
  - ▶ If required: Seal the unused connections on the device using the corresponding covers or protective caps.

7 System Settings

7.1 Switching the Device On or Off

If the device is being connected to the power supply for the first time, or after restoring the factory settings: The device is switched on and the Setup Wizard opens. All steps in the Setup Wizard must be completed.


If the device is switched on after the Setup Wizard is complete and **no** passwords have been assigned for the user profiles: The user profile for the last user to be logged in is loaded.

Procedure

- ▶ If the Setup Wizard is shown: Follow the instructions in the Setup Wizard in the operating display.
- ▷ If the Setup Wizard is **not** shown: The login display is shown or the user profile of the last user to be logged in is loaded.
- ▶ To switch the device off: Disconnect the device from the power supply.

7.2 User Login or Logout

Procedure

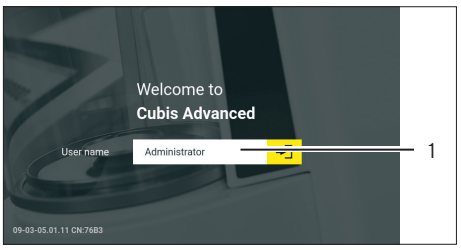

- ▶ If the name of the preferred user profile is displayed in the input field (1) of the login display: Press the [Login] button.
- ▶ If the name of the preferred user profile is **not** displayed in the input field of the login display:
  - ▶ Press the input field in the login display.
  - ▷ The user selection is opened.
  - ▶ Press the name of the preferred user profile.
  - ▶ If a password has been assigned for the preferred user: Enter the password in the input field and press the [OK] button.
- ▷ The user profile opens and the main menu or the last task being performed by the active user is displayed.
- ▶ To log out the active user profile from the device: Press the [Logout] button.

7.3 Performing System Settings

Default settings can be adjusted for the device and the applications in order to align with the ambient conditions and individual operating requirements.

The following settings are necessary to operate the device together with connected components:

- Set up the communication of the connected devices
- Set up additional components

The following settings are recommended to set up the device:

- Set the behavior of the isoCAL function
- Set the behavior of the motorized draft shield (only for devices with a motorized draft shield)
- If the associated QAPP is activated and the LDAP server is configured: Enter a password.

Procedure

▶ Open the main menu.
▶ Press the [Setup] button.
▶ To adjust settings: Open the preferred submenu.
▶ Select the preferred configuration value.
▶ Exit the menu.
▷ With some settings, the [Booting device] display appears in the operating display and the device restarts.

7.4 Using the Help Function

If help texts are available in a menu: The [Help] button is displayed.

Procedure



▶ Press the [Help] button.
▶ The help texts are displayed.
▶ To navigate through the help text: Swipe the text downwards or upwards.



▶ If the full text search is to be used:
  ▶ Press the [Search] button and enter a search term.



  ▷ The [Active search] button is displayed. The search hits are highlighted in yellow.



  ▶ To switch between the search hits: Press the [Next] and [Back] buttons.



- ▶ To call up a structure in a help text: Press the [Info] button.
- ▷ The structure for the help text is displayed. The structure headings are links.
- ▶ Press on the preferred heading.




- ▶ To export the help text: Press the [Export] button.
- ▷ The available connectors are displayed, e.g. the connector [USB stick].
- ▶ Select the connector and follow the instructions in the operating display.

7.5 Activating Applications (QAPPs)

All applications from the QAPP package “Essentials” are activated for the device at the factory. Additional applications may be activated in the QAPP Center. These applications and QAPP packages can be tested for 30 days free of charge, and after that require a license.

Procedure

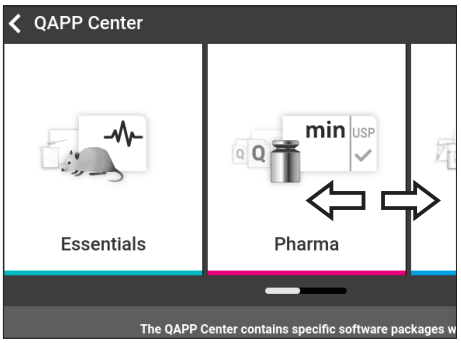

- ▶ Open the task.
- ▶ Press the [QAPP Center] button.
- ▷ An overview of the available QAPP packages is displayed.
- ▶ Select the preferred QAPP application package e.g. “pharmaceuticals” or “hardware”.
- ▷ A list of the applications contained in the QAPP package is displayed.

7.5.1 Activating QAPP Package With All Applications

Procedure

- ▶ Press the [License] button.
- ▷ The input field for the license key is displayed.
- ▶ If an additional cost is associated with the QAPP package: Enter the license key in the input field and press the [OK] button.
- ▶ If no additional cost is associated with the QAPP package: Press the [OK] button.

7.5.2 Activating Individual Applications From QAPP Package

Procedure

- ▶ Press the preferred application.
- ▷ A display opens containing details about the selected application.
- ▶ Press the [License] button.
- ▷ The input field for the license key is displayed.
- ▶ If an additional cost is associated with the application: Enter the license key in the input field and press the [OK] button.
- ▶ If no additional cost is associated with the application: Press the [OK] button.

7.6 Switching Off the isoCAL Function

Procedure

- ▶ In the “Settings / Device settings / Safe Weighing / isoCAL execution mode” menu, for the “isoCAL function” parameter, select the “Off” configuration value.

7.7 Configuring Motorized Opening and Closing of the Draft Shield

Pressing the palm key on the control unit of the device enables motorized opening and closing of the draft shield doors. The draft shield has a learning capability which enables the following opening parameters to be saved:

- It is possible to control all doors together or individually.
- The opening width of the doors can be adjusted.

Requirements

The application for use of the motorized draft shield is activated.

Procedure

- ▶ Close all draft shield doors.
- ▶ To determine how far a door is opened by pressing the palm key: Manually push the door open to the preferred position.
- ▶ If several doors are to be controlled by the palm key simultaneously: Manually push the preferred doors open to the preferred position.
- ▶ Press the preferred palm key.
- ▷ All open doors are closed.
- ▷ The settings for motorized opening and closing of the draft shield are saved.

7.8 Managing User Data

7.8.1 Assigning a Password

Procedure

- ▶ Log into the device using the user profile for which a password needs to be assigned.
- ▶ Open the “Settings / User management” menu.
- ▶ Press the [Edit password] button.
- ▷ The user password input field is displayed.
- ▶ Enter the preferred password in the input field and confirm with the [OK] button.

7.8.2 Changing a User Profile or User Role

For certain functions, the “User management” extension must be licensed in the QAPP Center.

Procedure

- ▶ To change a user profile: Open the “Settings / User management” menu.
- ▶ To change a user role: Open the “Settings / Access management / Role management” menu.
- ▷ An overview of the saved user profiles or user roles is displayed.
- ▶ Press the preferred user profile or user role.
- ▷ An overview of the settings is displayed.



- ▶ Press the [Edit] button.
- ▷ The “User settings” or “Role parameters” display is opened.
- ▶ Specify the preferred settings in the “User management” menu and the “Access management” menu.
- ▶ Confirm the entries. To do this, press the [OK] button.

7.8.3 Creating a User Profile

Requirements

The “User management” extension is licensed in the QAPP Center.

Procedure

- ▶ Open the “Settings / User management” menu.
- ▷ An overview of the saved user profiles is displayed.



- ▶ Press the [New] button.
- ▷ The “User settings” display is opened.
- ▶ Specify the preferred settings in the “User management” menu and the “Access management” menu.
- ▶ Confirm the entries. To do this, press the [OK] button.

7.8.4 Creating a User Role

Requirements

The “User management” extension is licensed in the QAPP Center.

Procedure

▶ Open the “Settings / Access management / Role management” menu.
▷ An overview of the saved user role is displayed.
▶ If an existing entry is to be used as a template:
  ▶ Press the user role.
  ▶ Press the [Select all] button.
  ▷ The settings are applied and it is opened in a new entry.





▶ If a new entry needs to be created: Press the [New] button.
▷ The “Role management” display is opened.
▶ Specify the preferred settings in the “User management” menu and the “Access management” menu.
▶ Confirm the entries. To do this, press the [OK] button.

7.9 Setting up Device for Network Printer

7.9.1 Selecting Configuration for Network Printing

Configuration Options	Description
Network Printing via Independent Wi-Fi Network	– Device and network printer communicate via an independent Wi-Fi network – Standard Wi-Fi components can be used – Special security measures or network settings are **not** required – Recommended configuration for easy printing, without the use of additional network functions

Setting up Network Printing via Company Network	– Device and network printer communicate via the company network
	– All network functions of the device are available in the company network
	– The company network must meet the following requirements:
	  – DHCP server present
	  – Standard authentication methods allowed, **no** special security measures are required to log in via Wi-Fi

Procedure

▶ Check which configuration is suitable for network printing.

7.9.2 Setting up Device for Network Printer via Independent Wi-Fi Network

Requirements

- The individual user profile has the role “Administrator” or “Service”.
- The network printer is connected to a Wi-Fi router.

Procedure

- ▶ Test the suitability of the Wi-Fi router and Wi-Fi stick (for suitability, see Chapter “16 Accessories and Spare Parts”, page 93).
- ▶ Insert the Wi-Fi stick into the device in a USB connection.
- ▶ Open the “Settings / Connections / Network / Wi-Fi” menu.
- ▶ Press the [Edit] button.
- ▶ Enter the service set identifier for the router under “Wi-Fi SSID”. The service set identifier is specified on the back of the Wi-Fi router.
- ▶ Enter the router password under “Wi-Fi password”. The password is specified on the back of the Wi-Fi router.
- ▶ Press the [Back] button.
- ▷ In the “Wi-Fi” overview display, the status “Ready” is displayed within 10 seconds.

Creating a Network Printer

Procedure

- ▶ Open the “Settings / Connections / Connectors / YDP30-NET” menu.
- ▶ Press the [New] button.
- ▶ Enter a name for the network printer under “Connector name”, e.g., printer YDP30-NET.
- ▶ Enter the network printer IP address under “IP or Host”. The IP address is specified in the network printer “Network Settings” printout.
- ▶ Confirm the entries. To do this, press the [OK] button.

7.9.3 Setting up Device for Network Printer via Company Network

Requirements

- The individual user profile has the role “Administrator” or “Service”.
- The company network has a DHCP server.
- The network printer is connected to the company network.

Procedure

- ▶ Open the “Settings / Connections / Network / Wi-Fi” menu.
- ▶ Press the [Edit] button.
- ▶ Enter the service set identifier for the company network under “Wi-Fi SSID”. The service set identifier is available from the company network administrator.

- ▶ Enter the company network password under “Wi-Fi password”. The password is available from the company network administrator.
- ▶ Press the [Back] button.
- ▷ In the “Wi-Fi” overview display, the status “Ready” is displayed within 10 seconds.

Creating a Network Printer

Procedure

- ▶ Open the “Settings / Connections / Connectors / YDP30-NET” menu.
- ▶ Press the [New] button.
- ▶ Enter a name for the network printer under “Connector name”, e.g., printer YDP30-NET.
- ▶ Enter the network printer “HOST NAME” under “IP or Host”. The HOST NAME is specified in the network printer “Network Settings” printout.
- ▶ If connection with the host name is delayed: Enter the network printer IP address under “IP or Host”. The IP address is specified in the network printer “Network Settings” printout.
- ▶ Confirm the entries. To do this, press the [OK] button.
- ▶ If the IP address changes depending on the DHCP network settings: Enter the new IP address of the network printer again.

7.9.4 Configuring print profiles

Procedure

- ▶ Open the “Settings / Weighing and print profiles / YDP30” menu.
- ▶ Edit a pre-configured print profile or create a new print profile by pressing the [New] button, e.g., network printer YDP30-NET.
- ▶ Follow the wizard’s instructions in the operating display for more settings.
- ▶ Confirm the entries. To do this, press the [OK] button.
- ▶ Add the created print profile to a task (see Chapter “8.7 Adding Weighing and Print Profiles to a Task”, page 65).

7.10 Downloading Additional Information

As part of the MCA firmware package on the Sartorius website, additional information is available for the device, e.g. a description of interface protocols, or a set of installation instructions for a website certificate. The information is available as a PDF file, partially in English.

Procedure

- ▶ Download the "Cubis® MCA Firmware" from the Sartorius website (www.sartorius.com).
- ▶ Find the additional information you require, e.g. the description of interface protocols.

8 Operation

8.1 Observing Warm-up Time

After connecting to the power supply, the warm-up time must be observed. This enables the device to reach its required operating temperature and ensures accurate values during weighing processes.

Procedure

- ▶ Ensure that the warm-up time is observed (see Chapter “15.6 Warm-up Time to Reach Operating Temperature”, page 87).

8.2 Manually Opening or Closing a Draft Shield

All doors can be fully or partially opened.

Procedure

- ▶ To open the draft shield manually, e.g. the right-hand door: Take the corresponding door handle and push backwards.
- ▶ To close the draft shield manually, e.g. the right-hand door: Take the corresponding door handle and push it fully forwards.

8.3 Motorized Opening or Closing of a Draft Shield

8.3.1 Opening or Closing on the Control Unit

Procedure

- ▶ Press the palm key. This enables motorized opening or closing of the motorized draft shield according to the stored setting.

8.3.2 Opening or Closing via Proximity Sensors

The proximity sensor works in “crossover mode”:

- Left proximity sensor: Opens and closes the right-hand and the upper door
- Right proximity sensor: Opens and closes the left-hand and the upper door

The sensitivity of the proximity sensors can be adjusted (see Chapter “4.13.2 “Settings” Menu”, page 37). Wearing safety gloves may inhibit the use of the proximity sensor.

If a motorized inner draft shield is mounted: The motorized inner draft shield also opens.

Procedure

▶ Hold your hand in front of the left or right proximity sensor. This enables the door to be fully opened or closed.

8.3.3 Opening or Closing a Motorized Draft Shield at the Door Handle

Procedure



▶ To open or close a door: Press on the door handle (1) of a door. The door is automatically fully opened or closed.

8.4 Leveling the Device

Leveling compensates any inclines at the device’s installation site. If leveling is necessary: The [Leveling] button is displayed in the weighing display and a message is displayed in the status center.

Procedure

▶ If the weighing display is opened: Press the [Leveling] button.
▶ When the status center is displayed: Press the [Perform leveling] button.
▷ The Leveling Wizard opens.
▶ Follow the wizard’s instructions.

8.5 Calibration, Adjustment or Linearization

Function	Description
Calibration	The device checks how much the displayed value deviated from the specified setpoint.
Adjustment	The device corrects the deviation to the setpoint.
Linearization	The device corrects the deviation from the ideal characteristic curve and the setpoint (only for models MCA206MC | MC106MC | MCA66MC | MCA36MC | MCA605MC).

Only for devices MCA206MC | MCA1005MC | MCA1004MC: The electrical weighing range is substituted for the full specified weighing range in a fully automatic, mechanical manner. Absolute weight values can only be calibrated, adjusted, and linearized within the electrical weighing range.

The device needs to be calibrated and adjusted regularly. Various methods can be selected for this:

- Adjusting With the isoCAL Function
- Internal Calibration and Adjustment
- External Adjustment
- Internal linearization

Only internal adjustment is described below.

Procedure

▶ Where one of following conditions occurs, calibrate and adjust the device using the preferred method:
  - Daily, every time the device is switched on
  - After every leveling
  - After changing the ambient conditions (temperature, humidity, or air pressure)
  - After setting the device up at a new installation site

8.5.1 Adjusting With the isoCAL Function

The device can be automatically internally calibrated and adjusted using the isoCAL function.

Requirements

- The isoCAL function is configured in the “Safe Weighing” menu, e.g. “On, automatic execution”.
- The conditions for triggering and executing the isoCAL function are met (see Chapter “15.9 Conditions for isoCAL Function”, page 88).

Procedure

- ▶ If the isoCAL function is set to automatic start and the isoCAL function is triggered:
  - ▷ The [isoCAL] button flashes in the operating display.
  - ▶ Wait until the isoCAL function is executed.
  - ▷ In the operating display, a time display counts down from 15 seconds to 0.
  - ▷ If **no** load change or **no** operation takes place on the device before the expiration of the time display: The isoCAL function starts.
- ▶ If the isoCAL function is set to manual start and the isoCAL function is triggered:
  - ▷ The [isoCAL] button flashes in the operating display.
  - ▶ Press the [isoCAL] button.
  - ▷ The isoCAL function starts.
- ▷ If the isoCAL function is complete: The device confirms the completion of the calibration or adjustment process with an acoustic signal, and the calibration report is displayed.
- ▶ To output the calibration report via a connector: Press the [Print Memory] button.
- ▶ To close the calibration report and return to the previous display: Press the [OK] button.

8.5.2 Adjusting the Device Internally

Requirements

The weighing pan is unloaded.

Procedure

- ▶ Open the main menu.
- ▶ Press the "Internal adjustment" task.
- ▷ The internal adjustment function is executed.
- ▷ If automatic leveling is configured: The device levels itself automatically.
- ▷ When the adjustment function is complete: The device confirms the completion of the adjustment process with an acoustic signal, and the adjustment report is displayed.
- ▶ To output the adjustment report via a connector: Press the [Print Memory] button.
- ▶ To close the adjustment report and return to the previous display: Press the [OK] button.

8.6 Creating a Task

The QAPP applications must be configured in a task so that they can be executed, e.g. defining weighing and print profiles.

Procedure

- ▶ Open the task.
- ▶ Press the [New] button.
- ▷ A list of the activated applications is displayed.
- ▶ To select an application: Press the preferred application.
- ▷ The wizard for creating a new task starts.
- ▶ Follow the wizard's instructions in the operating display.

8.7 Adding Weighing and Print Profiles to a Task

Weighing or print profiles must be added to a task in order to be used. The weighing and print profiles must be configured in the Settings menu.

Procedure

- ▶ Open the task.
- ▶ Start the wizard to create or edit a task and follow the wizard's instructions in the operating display.

8.8 Weighing

Handling the mass comparator requires steady and continuous work. Observe the recommendations from OIML R111 during this process.

Perform a few trial weighing procedures before you begin the actual weighing to get a feel for the handling method.

When weighing chemicals, suitable vessels must be used for the material to be weighed. This makes it possible to prevent damage to the device or its accessories.

Requirements

The device is leveled and adjusted.

Procedure

- ▶ Avoid measurement errors due to moisture or warmth from the hand: Always wear gloves or use forceps or similar utensils.
- ▶ Start a task with weighing function.
- ▶ Press the [Zero] button.
- ▶ Carefully center the weights on the weighing pan.
- ▶ If magnetized/magnetizable samples need to be identified:
  - ▶ To keep such effects to a minimum, increase the distance between the sample and the weighing pan using a non-magnetic material.
  - ▶ Shield from external magnetic fields to some extent using soft magnetic plates.
- ▶ If below-balance weighing is being carried out: Suspend the sample on the hook for below-balance weighing, e.g. with a wire.
- ▶ If using a vessel for the sample:
  - ▶ Place the vessel for the material to be weighed onto the weighing pan.
  - ▶ Press the [Tare] button. This compensates for the weight of the vessel.
  - ▶ Place the sample in the vessel.
- ▶ Once the weight value is displayed in black and the weighing unit is displayed: Read off the measured value.

8.9 Weighing and Printing with ID Marking

8.9.1 Saving Values for Printouts

Every lot and every sample can be assigned an ID number. The ID numbers are saved in the print memory and exported during the printing process.

Requirements

The prompt for the batch IDs and sample IDs is activated for the “Weighing” task.

Procedure

- ▶ Open the main menu.
- ▶ Press the “Weighing” task.
- ▷ The input field for the lot ID is displayed.
- ▶ If a barcode scanner is connected to the device: Scan or type the lot ID into the input field.
- ▶ If no barcode scanner is connected to the device: Type the lot ID into the input field.
- ▶ Press the [OK] button.
- ▷ The weighing display is opened.

Entering a Sample ID

Procedure

- ▶ Press the [Zero] button.
- ▶ Place the sample on the weighing pan.
- ▶ Press the [Save] button.
- ▷ The input field for the sample ID is displayed.
- ▶ Type or scan the sample ID into the input field.
- ▶ Press the [OK] button.
- ▷ The [Print memory] button is displayed in the operating display.
- ▷ The weight value and entered IDs are saved.
- ▶ If additional values are to be saved:
  - ▶ Remove the sample being weighed.
  - ▶ Place the next sample on the weighing pan and press the [Save] button.
  - ▶ Type or scan the sample ID into the input field.
  - ▶ Press the [OK] button.

8.9.2 Marking Saved Values as Invalid

Procedure

- ▶ Press the [Print Memory] button.
- ▷ The print memory opens and a list of all saved values is displayed.
- ▶ Press the preferred items.
- ▷ An overview of all data to be printed is displayed for the selected items.
- ▶ Press the [Yes] button.
- ▶ If a reason for the invalidity of the value needs to be displayed on the printout: Enter a reason in the input field and press the [OK] button.
- ▶ If no reason for the invalidity of the value needs to be displayed on the printout: Press the [No] button.

8.9.3 Marking Saved Values as Valid

Procedure

- ▶ Press the [Print Memory] button.
- ▷ The print memory opens and a list of all saved values is displayed.
- ▶ Press the items which are preferred to be marked as invalid.
- ▷ An overview of all data to be printed is displayed for the selected items.
- ▶ Press the [Yes] button.

8.9.4 Printing Saved Values

Procedure

▶ Press the [Print Memory] button.
▷ The print memory opens and a list of all saved values is displayed.
▶ To start the print process using the print profile embedded in the current task: Press the [Print] button.
▷ A file, e.g. a PDF file, is produced and sent to the printer defined in the print profile. For two print profiles with different connectors: The print processes will be started in succession.

8.9.5 Ending the Task

Procedure

▶ Press the [End] or [Menu] button.
▷ If additional values are to be saved in the print memory:
  ▷ A dialog for prematurely ending the task is opened.
  ▶ To return to the weighing display and print the saved values: Press the [Yes] button and print the saved values.
  ▶ To end the task and delete the saved values from the print memory: Press the [No] button.
▷ The task ends and the main menu is opened.

8.10 Viewing the Data Storage Device

The weight values are stored in the data storage device along with the date, time, and process number. The content of the data storage device can be searched and sorted using a filter. The data storage device has a ring buffer. This means that the oldest data records can be overwritten by new ones. The data storage device is designed for a specific number of data records (for the number of data records see Chapter “15.10 Data Storage Device”, page 88).

The operator is responsible for regularly exporting the entries for archiving in the data storage device. The available storage capacity is displayed to the operator.

Procedure

- ▶ Open the “Settings” / “Device settings” / “Device information” menu.
- ▶ Press the “View data storage device” menu item.
- ▷ A list of all the values saved in the data storage device is displayed.
- ▶ If only the values for a preferred date need to be displayed:
  - ▶ Press the [Filter] button.
  - ▶ Type the preferred date into the input field.
  - ▶ To filter the list with the date entered: Press the [OK] button.
- ▶ If the list needs to be searched for a particular ID:
  - ▶ Press the [Search] button.
  - ▶ Type the preferred ID into the input field.
  - ▶ To start the search with the entered ID: Press the [OK] button.

8.11 Running Applications (Examples)

8.11.1 Executing the “Unit change” function

The “Unit change” function enables switching between the different units and resolutions defined in the weighing profile of the active task. The units and resolutions can be set at the beginning of the weighing process.

Procedure

- ▶ Start the preferred task.
- ▶ Press the [Unit change] button.
- ▷ All units defined in the weighing profile for the active task are displayed in a list.
- ▷ All resolutions for the weight value defined in the weighing profile for the active task are displayed in a list.
- ▶ Press the preferred unit.
- ▶ To set the resolution for the selected unit: Press the preferred resolution.
- ▶ To confirm the selection and return to the weighing display: Press the [OK] button.
- ▷ The current weight value is displayed in the selected unit and resolution.

8.11.2 Running the “Statistics” Application (from “Essentials” QAPP Package)

The “Statistics” application saves up to 1000 weight values and evaluates these statistically.

The following values are saved and exported by the statistics application:

- Number of components
- Mean value
- Standard deviation
- Variation coefficient
- Sum of all values
- Lowest value (minimum)
- Highest value (maximum)
- Spread: Difference between maximum and minimum

The “Statistics” application can be combined with the following functions:

- Unit change, only available in the weighing display, before saving the initial weight value
- ID marking
- Automatic taring

Procedure

- ▶ Open the main menu.
- ▶ Start the task for the “Statistics” application.
- ▶ Press the [Zero] button.
- ▶ Place a sample on the weighing pan.
- ▶ To start the recording of the statistics: Press the [Confirm] button.
- ▷ The current weight value is saved and the device is automatically tared.
- ▶ If the next value is to be saved:
  - ▶ Remove the sample on the weighing pan.
  - ▶ Place a new sample on the weighing pan and press the [Confirm] button.
- ▶ To display an overview of the saved data: Press the [Report] button.
- ▶ To print the current statistics: Press the [Print] button.
- ▶ To end the current statistics, and to delete the saved values: Press the [End] button.

9 Cleaning and Maintenance

9.1 Preparing the Device for Cleaning

Procedure

- ▶ If an accessory is connected to the device: Disconnect the accessory from the device (see instructions for the accessory).
- ▶ If the upper or lower cover for the rear wall of the weighing compartment is removed: Insert the cover for the weighing compartment rear wall into the rear wall of the weighing compartment (see Chapter 6.5, page 51).

9.2 Cleaning the Device

Sartorius recommends cleaning the device at regular intervals, e.g. weekly. There must be **no** foreign substances, e.g., particles, fibers or liquids, in the area of the weighing pan.

To clean the device, you can use the cleaning utensils supplied or a moistened cloth. Cleaning of the device is described in the QAPPs for device cleaning, e.g., “Advanced Cleaning”.

Procedure

- ▶ ⚠ **WARNING** Risk of injury due to electrical current! Protect the power supply unit and the power supply cable from liquids.
- ▶ Only use suitable cleaning agents and cleaning procedures and observe the product information for the cleaning agent used (cleaning agents see Chapter “15.14 Cleaning Agents and Cleaning Procedures”, page 90).
- ▶ Open the QAPP for cleaning the device and follow the instructions in the operating display.
- ▶ If the dirt is **not** coming away easily:
  - ▶ Lightly spray the surfaces of the device with the cleaning agent.
  - ▶ Leave the cleaning agent to work (for recommended contact time, see Chapter “15.14 Cleaning Agents and Cleaning Procedures”, page 90).
  - ▶ Dry the surfaces with a **non**-abrasive cloth.

9.3 Maintenance Schedule

Interval	Component	Activity	Chapter, Page
12 to 24 months, depending on operating conditions	Device	Contact Sartorius Service.	

9.4 Recommissioning

Procedure

- ▶ Reinstall the removed components back onto the device (see Chapter 5, page 43):
  - — Draft shield
  - — Weighing Pan and Associated Components
  - — Accessories
- ▶ Connect the preferred accessories (see Chapter 6.5, page 51).
- ▶ Connect the device to the power supply (see Chapter 6.4, page 51).

9.5 Performing a Software Update

A software update can be installed (software package) from a USB mass storage device using the device's USB-A or USB-C connection. An update can also be carried out from a server via other connectors of the device. The following describes the installation of a USB mass storage device, taken from the Sartorius website.

A software update can extend or change the functionality of the device. Sartorius recommends carrying out software updates:

- — Prior to starting the software updates, save the device data to a USB mass storage device.
- — If a QAPP Center update is to be carried out in addition: Perform the software update for the device first.

Two files are required for the software update: Firmware file with file ending ".upd", and checksum file with file ending ".upd.md5".

Execution of the software update and troubleshooting during this process are described in the help text "Device maintenance" (for accessing help texts, see Chapter "7.4 Using the Help Function", page 53).

Requirements

The device is connected to the power supply.

Procedure

- ▶ Download the software package from the Sartorius website onto the USB mass storage device. To do this, download the "Cubis® MCA Firmware" file.
- ▶ If this relates to a zip file: Unpack the software package on the USB mass storage device. The files must be saved in the main directory (at root level). The files must **not** be copied into a folder.
- ▶ Insert the USB mass storage device with the software package into the device's USB-A connection.
- ▶ Press the "Update firmware" menu entry in the "Settings / Device maintenance" menu.
- ▶ Press the "USB stick" as the connector and select the preferred software version.
- ▷ The software update takes approx. 5 minutes.
- ▷ Once the software update is complete: The software version number is updated in the login display.

9.6 Performing a QAPP Center Update

The QAPP Center package can be installed from a USB mass storage device using the device's USB-A or USB-C connection. An update can also be carried out from a server via other connectors of the device. The following describes the installation of a USB mass storage device, taken from the Sartorius website.

Sartorius recommends the following when carrying out QAPP Center updates:

- Prior to starting the QAPP Center update, save the device data to a USB mass storage device.
- If a software update is also to be carried out: Perform the software update for the device first.

Two files are required for the QAPP Center update: QAPP Center with file ending ".appcenter" and checksum file with file ending ".qappcenter.md5.md5".

Execution of the QAPP Center update and troubleshooting during this process are described in the help text "Device maintenance" (for accessing help texts, see Chapter "7.4 Using the Help Function", page 53).

Requirements

- The device is switched on.
- The QAPP Center package is saved on a USB mass storage device or on a server via a connector.

Procedure

- ▶ Download the QAPP Center package from the Sartorius website onto the USB mass storage device. To do this, download the "Cubis® MCA Firmware" file.
- ▶ If this relates to a zip file: Unpack the QAPP Center package on the USB mass storage device. The files must be saved in the main directory (at root level). The files must **not** be copied into a folder.
- ▶ Insert the USB mass storage device with the QAPP Center package into one of the device's USB-A connections.
- ▶ Press the "Install QAPP Center" menu entry in the "Settings" / "Device maintenance" menu.
- ▶ Press "USB stick" as the connector.
- ▶ Press the preferred package.
- ▶ Once the QAPP Center update is complete: Confirm successful installation with the [OK] button.
- ▷ Existing tasks remain unchanged after the QAPP Center update. In the existing tasks, the original QAPP versions are used.
- ▶ In order to use the new QAPP version: Create a new task with the new QAPP version. Existing tasks are **not** automatically adapted via the QAPP Center update.

10 Malfunctions

10.1 Malfunctions on the Control Unit or During Weighing Processes

Malfunction	Cause	Solution	Chapter, Page
The operating display is blank.	The device is disconnected.	Check whether the power supply unit is connected to the device and the power supply at the installation site.	6.3, 50; 6.4, 51
The operating display is blank or an error is displayed.	The control unit is **not** connected.	Check if the connection cable for the control unit is connected to the device.	6.1, 50
The displayed weight value changes constantly.	The device installation site is unstable.	Adjust the parameters in the “Ambient conditions” submenu.	
		Change the installation site.	
	A foreign object is positioned between the weighing pan and the housing.	Remove the foreign object.	
The weight readout displayed by the device is obviously wrong.	The device has **not** been adjusted.	Adjust the device.	
	The device was **not** tared before weighing.	Tare the device.	
The side or upper panel of the motorized draft shield does not open or close.	The operating range of the side or upper panel of the motorized draft shield is not functional or is blocked.	Hold down the palm key. If the problem continues to occur: Contact Sartorius Service.	

11 Decommissioning

11.1 Decommissioning the Device

Procedure


- ▶ Disconnect the device from the power supply.
- ▶ Disconnect the device from all connections at the installation site.
- ▶ If an accessory is connected to the device: Disconnect the accessory from the device (see instructions for the accessory).
- ▶ If the upper or lower cover for the rear wall of the weighing compartment is removed: Insert the cover for the weighing compartment rear wall into the rear wall of the weighing compartment (see Chapter 6.6, page 51).
- ▶ Clean the device.

11.2 Dismantling Device Parts

11.2.1 Removing the Weighing Pan and Associated Components


Procedure

- ▶ ⚠ **CAUTION** Risk of cuts due to broken glass! Handle the product with care.
- ▶ Open the draft shield, e.g. the right-hand door.
- ▶ Remove the associated components from the device:
  - — Device with scale interval 1 µg: Glass cylinder and lid
  - — Weighing pan (1)
  - — Base plate of weighing compartment (2)

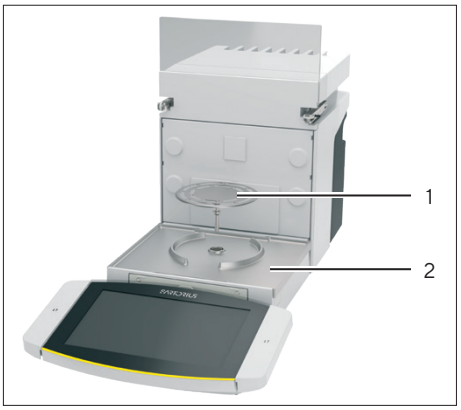


11.2.2 Dismounting the Draft Shield


Procedure

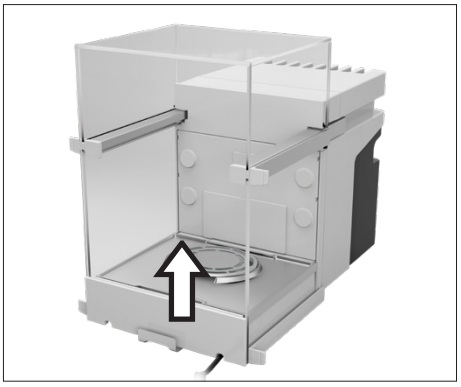

- ▶ Remove the control unit from the device.
- ▶ Pull off the front panel of the device upwards.

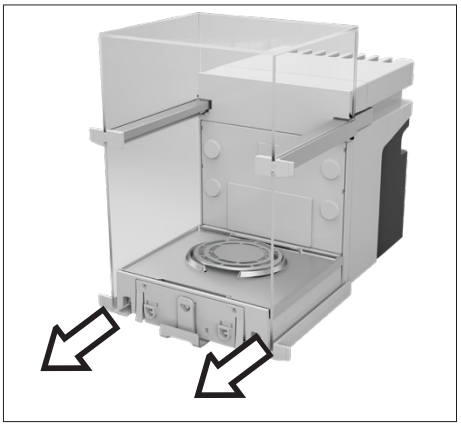

- ▶ Remove the right and left doors. To do this, pull each door forwards to pull it out of the guide rods on the device.

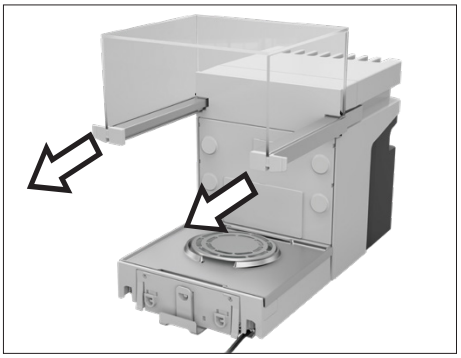


- ▶ Pull the upper door forwards on both sides and remove it from the guide rods on the device.

11.2.3 Device with scale interval 1 µg: Dismounting the Adapter Module

Requirements

- The device is **not** connected to the power supply.
- The device draft shield is fully opened.
- The front panel of the draft shield is removed.


Procedure

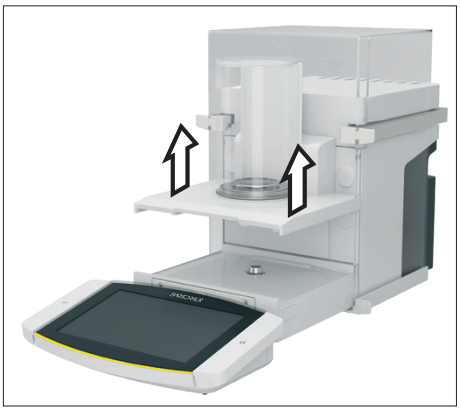

▶ Lift the adapter module out of the device.

11.2.4 Dismounting the Connection Cable

Procedure

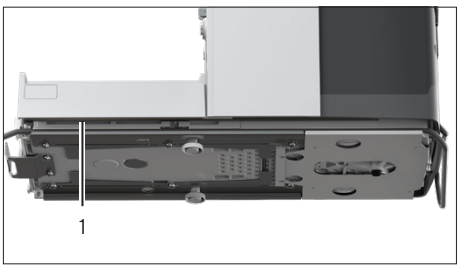

▶ Place the device on its side (see Chapter 5.5, page 45).
▶ Remove the connection cable for the control unit (1) out of the cable channel along the side of the device.
▶ Remove the Ethernet connection cable out of the cable channel along the side of the device.
▶ Place the device back onto its base.

12 Transport

12.1 Producing the Transport Lock (Only for Models MCA1005MC | MCA605MC | MCA1004MC)


Procedure



- ▶ **NOTICE** The weighing system may be damaged if the transport lock is not mounted!
- ▶ Remove the control unit from the device.
- ▶ Push the front of the device over the edge of the table so that the safety and sealing screws are accessible.



- ▶ Unscrew the 2 safety and 2 sealing screws (1).

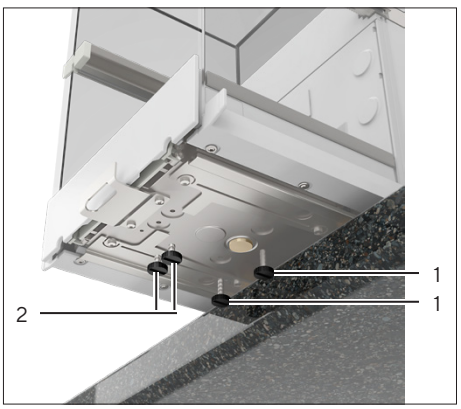


- ▶ If the weighing system needs to be secured with the safety screws:
  - ▶ Screw the short sealing screws (2) into the front position.
  - ▶ Screw the longer safety screws (1) into the rear position.

12.2 Transporting the Device

Requirements

- The device has been decommissioned.
- The control unit is secured to the device.
- Only for models MCA1005MC | MCA605MC | MCA1004MC: The transport lock is in place (see Chapter 12.1, page 79).

Procedure

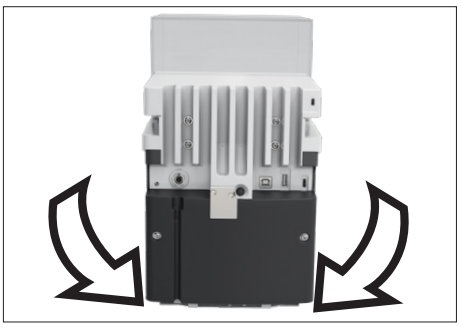

- ▶ ⚠ **CAUTION** Risk of injury when lifting or transporting!
  - ▶ Disconnect the device from all connections at the installation site.
  - ▶ Use both hands when transporting the device and setting it down. To do this, reach sideways under the device with both hands in the rear part of the device.
  - ▶ Do **not** carry the device on the draft shield or the control unit.
- ▶ For transporting over longer distances, use a trolley with soft mats. The control unit must be positioned entirely on the base.

13 Storage and Shipping

13.1 Storing

Requirements

The device has been decommissioned.

Procedure

- ▶ Store the device according to the ambient conditions (see Chapter “15.2 Ambient Conditions During Storage and Transport”, page 83).

13.2 Returning the Device and Parts

Defective devices or parts can be returned to Sartorius. Returned devices must be clean and packed in their original packaging.

Transport damage as well as measures for subsequent cleaning and disinfection of the device or parts by Sartorius are charged to the sender.

Devices contaminated with hazardous materials, e.g., harmful biological or chemical substances, will **not** be accepted for repair or disposal.

Procedure

- ▶ Decommission the device.
- ▶ Contact Sartorius Service for instructions on how to return devices or parts (see www.sartorius.com).
- ▶ Pack the device and its parts in their original packaging for return.

14 Disposal

14.1 Disposing of the Device and Parts

The device and the device accessories must be disposed of properly by disposal facilities.

A lithium cell battery, type CR2032, is installed inside the device. Batteries must be disposed of properly by disposal facilities.

Procedure

- ▶ Dispose of the device in accordance with local government regulations. Inform the disposal facility that there is a lithium cell battery, type CR2032, installed inside the device.
- ▶ Dispose of the packaging in accordance with local government regulations.

15 Technical Data

15.1 Dimensions and Weights

	Unit	Value
Dimensions		
Dimensions (L x W x H)	mm	501 x 240 x 301
Dimensions (L x W x H), with control unit removed	mm	376.5 x 240 x 301
With control unit removed: Max. distance between device and control unit	cm	25
Width of door opening, with fully opened side door	mm	145
Width of door opening, with fully opened upper door	mm	145
Weighing pan dimensions (diameter)		
Models MCA206MC | MC106MC | MCA66MC | MCA36MC	mm	50
Models MCA1005MC | MCA605MC | MCA1004MC	mm	90
Models MCA206MC | MC106MC | MCA66MC | MCA36MC): Glass cylinder dimensions		
Internal diameter	mm	80
Inner height	mm	115
Weight, with draft shield mounted, approx.	kg	16

15.2 Ambient Conditions During Storage and Transport

	Unit	Value
Temperature		
During storage and transport	°C	-20 – +60
Storage, dry		

15.3 Installation Conditions

15.3.1 Installation Site

	Unit	Value
Height above sea level, maximum	m	3000
Conventional laboratory rooms		
For indoor use only		
Permissible air pressure range	mbar	600 – 1200
No potentially explosive atmospheres		
Pollution level according to IEC 61010-1		2
The installation site meets the requirements relating to the ambient conditions.		
Suitable for protection class		
IP protection: Protected against dust and water according to DIN EN 60529-1		IP30
Access to operation-relevant parts is guaranteed.		
Space Requirements		
Suitable for the dimensions of the device and the associated components.		
Setup surface		
Suitable for the weight of the device and the associated components.		
Stable, fully flat, even, low vibrations		
Additional properties		
Not directly against a wall		
No heat from heating systems or direct sunlight		
No direct drafts from open windows, AC systems, or doors		
No vibrations		
No “heavy traffic” areas (personnel)		
No electromagnetic fields		
No dry air		
For air conditioned rooms, bear in mind that the humidity and temperature control of the air conditioning system works with time constants that are as long as possible and small amplitudes. Avoid erratic changes.		
Comply with the recommendations for ambient conditions from OIML R111.		
When setting up the device in an isolator with protective argon atmosphere		
Power supply unit is positioned outside the isolator and the argon atmosphere		

15.3.2 Ambient Conditions at the Installation Site

	Unit	Value
Temperature		
In operation	°C	+5 – +40
In operation, with isoCAL function	°C	+10 – +30
In operation, without isoCAL function	°C	+17 – +27
Relative humidity		
For temperatures up to 31°C	%	80
Then linear decrease from 80 % at 31°C to 50% at 40°C		

15.4 Acclimatization Before Power Supply

	Unit	Value
Time period between unpacking and connecting to a power supply	h	2

15.5 Electrical Data

15.5.1 Power Supply

	Unit	Value
Power supply only permitted via the power supply cable and power supply unit supplied by Sartorius		
Sartorius power supply unit, type 1000099844		
Primary (power supply unit)		
AC voltage	V	100 – 240 ± 10%
Frequency	Hz	50 – 60
Current consumption, maximum	A	0.8
Secondary (device)		
Max. DC voltage at 4.3 A output current	V	15 ±15%
Power consumption, maximum	W	64.5
Overvoltage category according to IEC 61010-1		II
Power supply cable		
Power supply cable according to IEC 60320-1 C13 | C14, with IEC plug, 3-pin, and with country-specific power plug		

15.5.2 Safety of Electrical Equipment

Safety provisions according to IEC 61010-1 Safety requirements for electrical equipment for measurement, control, and laboratory use – Part 1: General requirements

15.5.3 Electromagnetic Compatibility

Electromagnetic compatibility according to DIN EN 61326-1: Electrical equipment for measurement, control and laboratory use – EMC requirements – Part 1: General requirements (IEC 61326-1:2013)
Interference immunity: Suitable for use in industrial areas (Table 2 of the standard)
Transient emissions: Class B: Suitable for use in residential areas and areas that are directly connected to a low voltage network that (also) supplies residential buildings
NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures: — Reorient or relocate the receiving antenna — Increase the separation between the equipment and receiver — Connect the equipment into an outlet on a circuit different from that to which the receiver — Consult the dealer or an experienced radio/TV technician for help
Shielded Cables: Connections between the devices and peripherals must be made using shielded cables in order to maintain compliance with FCC radio frequency emission limits.
Modifications: Any modifications made to this devices that are not approved by Sartorius may void the authority granted to the user by the FCC to operate this equipment.

15.5.4 Interfaces

USB-A connection
Communication: USB host (master)
USB-B connection
Communication: USB device (slave)
Type of interface: Virtual serial interface (virtual COM-port, VCP) and “PC direct” communication
USB-C connection
Communication: Downstream-facing port (DFP), USB host (Master)
Communication: RS232 connection with accessory YCC-USB-C-D09M

15.6 Warm-up Time to Reach Operating Temperature

	Unit	Value
Time period between connecting to the power supply and carrying out weighing tasks	h	1

15.7 Model-specific and Licensable Device Functions

	Models with scale interval 1 µg	Models with scale interval 0.01 mg and 0.1 mg
Model-specific device functions		
Proximity sensor	Available	Available
Use of motorized draft shield	Available	Available
Motorized inner draft shield	Available	Accessories
Use of QAPP-controlled functions for specific tasks on the device, e.g., access management, electronic signature, device cleaning	Licensable	Licensable

15.8 Recommended Calibration Weight

15.8.1 Models MCA206MC | MCA106MC | MCA66MC | MCA36MC

		MCA206MC	MCA106MC	MCA66MC	MCA36MC
	Unit	Value	Value	Value	Value
Resolution	mg	0.005	0.001	0.001	0.001
Supported load	g	211	111	61	32
Recommended accuracy class		E2	E2	E2	E2
External test weight	g	100	100	50	20

15.8.2 Models MCA1005MC | MCA605MC | MCA1004MC

		MCA1005MC	MCA605MC	MCA1004MC
	Unit	Value	Value	Value
Resolution	mg	0.001	0.001	0.01
Supported load	g	1110	610	1110
Recommended accuracy class		E2	E2	E2
External test weight	g	500	500	500

15.9 Conditions for isoCAL Function

	Unit	Value
Possible Conditions for Triggering the isoCAL Function		
In the event of a temperature change	K	1.5
After a time interval	h	12
After successful leveling		
Required conditions for executing the isoCAL function		
Operating display is in weighing mode (**not** in menu)		
Alphanumeric inputs are **not** active.		
Min. time period without entry on the device	min	2
Min. time period with unchanged loading of the weighing pan	min	2
Max. loading of weighing pan, high-capacity	%	2

15.10 Data Storage Device

	Value
Maximum number of data records	150,000

15.11 Integrated Clock

	Unit	Value
Maximum deviation per month (RTC)	s	30

15.12 Backup Battery

	Unit	Value
Lithium battery, type CR2032		
Service life at room temperature, minimum	Years	10

15.13 Materials

Housing
Stainless steel 1.4401 | 1.4404, Aluminum
Plastic PBT | PA
Float glass Optiwhite
Control unit:
Aluminum, painted
Plastic PBT | PP
Float glass
Weighing pan: Titanium

15.14 Cleaning Agents and Cleaning Procedures

15.14.1 Approved Cleaning Agents

	Cleaning agents and concentration					
Device components	Ethanol, 70%	Isopropanol, 70%	Citric acid, 10%	Diluted hydrogen peroxide, 3.5%	Sodium hydroxide, 32%	Ecolab™ Klercide Sporicidal Chlorine
Draft shield						
Top door	x	x	xx	xx	-	xx
Left and right doors	x	x	xx	xx	-	x
Front panel	x	x	x	x	-	x
Components in the weighing compartment						
Weighing pan	x	x	x	x	xx	x
Shield plate	x	x	x	x	xx	x
Base plate of weighing compartment (removable)	x	x	x	x	xx	x
Rear wall of weighing compartment	xx	x	x	x	x	x
Weighing compartment support (holds the base plate of the weighing compartment)	x	x	x	x	x	x
Control unit, including operating display, without palm key	x	x	x	x	x	x
Palm key	x	x	x	x	-	x
Back of the device						
Plastic surfaces	x	xx	x	x	x	x
Cooling element	x	xx	x	x	x	x

x Suitable
xx Suitable, may affect visual appearance of the device, does **not** impact mechanical stability
- **Not** suitable

15.14.2 Approved Cleaning Procedures

Wiping the device surfaces with a lightly moistened cloth		
Spraying the device surfaces with cleaning fluid, contact time	Min	5–10

15.15 Metrological Data

15.15.1 Models MCA206MC | MCA106MC | MCA66MC | MCA36MC

		MCA206MC	MCA106MC	MCA66MC	MCA36MC
	Unit	**Value**	**Value**	**Value**	**Value**
Scale interval (d)	mg	0.001	0.001	0.001	0.001
Maximum capacity (Max)	g	211	111	61	32
Repeatability according to the “ABA” process					
Number of “ABA” cycles		5	5	5	5
Standard deviation of the differences	mg	0.004	0.002	0.002	0.0015
Standard deviation of the differences, typical value	mg	0.0025	0.0015	0.001	0.001
Repeatability according to the “ABA” process with defined test weight					
Test weight	g	100	10	2	2
Standard deviation of the differences	mg	0.003	0.0015	0.0007	0.0007
Electronic weighing/taring range	g	111	111	61	32
Substitution weights	g	100	–	–	–
Linearity deviation					
Tolerance	mg	0.01	0.008	0.008	0.006
Deviation when load is off-center, positions according to OIML R76					
Test weight	g	100	50	20	10
Tolerance	mg	0.05	0.03	0.02	0.020
Typical value	mg	0.04	0.02	0.01	0.008
Sensitivity drift of +10°C to +30°C	ppm/K	1	1	1	1
Typical stabilization time	s	3.5	3.5	3.5	3.5
Cycle time (ABA)	s	100	90	90	90

15.15.2 Models MCA1005MC | MCA605MC | MCA1004MC

		MCA1005MC	MCA605MC	MCA1004MC
	Unit	**Value**	**Value**	**Value**
Scale interval (d)	mg	0.01	0.01	0.1
Maximum capacity (Max)	g	1110	610	1110
Repeatability according to the “ABA” process				
Number of “ABA” cycles		5	5	5
Standard deviation of the differences	mg	0.02	0.02	0.07
Standard deviation of the differences, typical value	mg	0.015	0.01	0.05
Repeatability according to the “ABA” process with defined test weight				
Test weight		100	50	100
Standard deviation of the differences	mg	0.015	0.01	0.05
Electronic weighing/taring range	g	610	610	610
Substitution weights	g	500	–	500
Linearity deviation				
Tolerance	mg	0.1	0.1	0.15
Typical value	mg	0.06	0.06	0.1
Deviation when load is off-center, positions according to OIML R76				
Test weight	g	500	200	500
Tolerance	mg	0.35	0.25	0.4
Typical value	mg	0.25	0.2	0.3
Sensitivity drift of +10°C to +30°C	ppm/K	1	1	1
Typical stabilization time	s	3.5	3.5	3.5
Cycle time (ABA)	s	90	90	90

16 Accessories and Spare Parts

16.1 Accessories

This table contains a selection of the accessories that can be ordered. For information on other products, contact Sartorius.

Item	Quantity	Order number
Thermal transfer	thermal printer for GLP	GMP printouts on continuous paper and labels	1	YDP30	YDP30-NET
Wireless nano USB adapter for company network or independent Wi-Fi network, e.g., operation with a Sartorius network printer YDP30-NET (only for Europe)	1	YWLAN01MS
Wireless nano router, e.g., for Sartorius network printer YDP30-NET for operation in an independent Wi-Fi network (only for Europe)	1	YWLAN02MS
Connection cable for control unit, length 3 m	1	YCC01-MCD3-3
Connection cable with RS232 adapter, USB-C to RS232, 9-pin	1	YCC-USB-C-D09M
Ethernet extension, 1 m	1	YCC-RJ-CAT7
Motion sensor with USB connection cable	1	YHS02USB
For models with scale interval 0.01 mg and 0.1 mg: Inner draft shield		
Motorized	1	YDS125A
Manual	1	YDS125U
Climate module (connection to rear wall of weighing compartment)		
For determining temperature, air humidity and air pressure	1	YCM20MC
For determining temperature, air humidity and air pressure, with DAkkS calibration certificate	1	YCM20MC-DAKKS
Climate module tower (USB connection), for determining temperature, air humidity and air pressure	1	YCM20MC tower
Weighing pan		
For models with scale interval 0.01 mg and 0.1 mg: 50 mm weighing pan, slotted, with shield plate for 50 mm weighing pan	1	YWP09-3
Dust cover	1	YDCC2MCA
Under-pan hook for below-balance weighing	1	On request
Density determination set for solids and liquids	1	YDK03MS
Sartorius Wedge, software for data communication between the PC and device	1	YSW02
PC software “ScalesNet-M” for highly accurate mass determination	1	YSN03C
Additional control unit (second display)	1	YSD01

Item	Quantity	Order number
Weighing table		
Weighing table made from synthetic stone, with vibration dampening	1	YWT03
Weighing table made from wood with synthetic stone	1	YWT09
Wall console	1	YWT04
Cleaning kit	1	YCK01MC
Second draft shield	1	YDS24C
Calibration weights		
Design: Knob weight, 20 g, accuracy class E2 with DAkkS calibration certificate	1	YCW422-02
Design: Knob weight, 50 g, accuracy class E2 with DAkkS calibration certificate	1	YCW452-02
Design: Knob weight, 100 g, accuracy class E2 with DAkkS calibration certificate	1	YCW512-02
Design: Knob weight, 500 g, accuracy class E2 with DAkkS calibration certificate		YCW552-02

16.2 Spare Parts

This table contains a selection of the spare parts that can be ordered. For information on other products, contact Sartorius.

Item	Quantity	Order number
Left door for draft shield (external draft shield)	1	YCCDSL
Right door for draft shield (external draft shield)	1	YCCDSR
Upper door for draft shield (external draft shield)	1	YCCDSU
Front panel for draft shield (external draft shield)	1	YCCDSF

17 Sartorius Service

Sartorius Service is available for queries regarding the device. Visit the Sartorius website (www.sartorius.com) for information about the service addresses, services provided, or to contact a local representative.

For inquiries about the system or when contacting Sartorius Service in the event of a malfunction, ensure that you have the device information, e.g., serial number, hardware, firmware, configuration, close at hand so you can pass it on to Sartorius Service. To do this, refer to the information on the manufacturer's ID label and in the "Device Information" menu.

18 Conformity Documents

The attached documents declare the conformity of the device with the designated directives or standards.

19 Trademark Information

Ecolab™ Klercide is a registered trademark of Ecolab Europe GmbH.

Original

EG-/EU-Konformitätserklärung
EC / EU Declaration of Conformity

Hersteller
Manufacturer

Sartorius Lab Instruments GmbH & Co. KG
37070 Goettingen, Germany

erklärt in alleiniger Verantwortung, dass das Betriebsmittel
declares under sole responsibility that the equipment

Geräteart
Device type

Elektronischer Massekomparator | Netzgerät
Electronic Mass Comparator | Power Supply

Modelle
Models

MCAvMC, | ZAG65US15
v = 36, 66, 106, 206, 605, 1004, 1005

in der von uns in Verkehr gebrachten Ausführung allen einschlägigen Bestimmungen der folgenden Europäischen Richtlinien entspricht und die anwendbaren Anforderungen folgender harmonisierter Europäischer Normen einschließlich deren zum Zeitpunkt der Erklärung geltenden Änderungen erfüllt:

in the form as delivered fulfils all the relevant provisions of the following European Directives and meets the applicable requirements of the harmonized European Standards including any amendments valid at the time this declaration was signed listed below:

	EMV | EMC	RoHS	Maschinen | Machines
Richtlinie *Directive*	2014/30/EU	2011/65/EU	2006/42/EG 2006/42/EC
Norm(en) *Standard(s)*	EN 61326-1:2013 EN IEC 61326-1:2021	EN IEC 63000:2018	EN ISO 12100:2010 EN 61010-1:2010 +A1:2019 +AC2014-04

Die Person, die bevollmächtigt ist, die technischen Unterlagen zusammenzustellen:
The person authorised to compile the technical file: Sartorius Lab Instruments GmbH & Co. KG
Electronics & Product Compliance
37070 Goettingen, Germany

Sartorius Lab Instruments GmbH & Co. KG
Goettingen, **2024-05-21**

Mario Hespe
Managing Director

Halil Yildirim
Senior Compliance Manager (SLI)

*: angewandte, jedoch für Maschinen nicht harmonisierte Norm /
applied standard, which however is not harmonized for machines

Doc: 3331587-00 SLI24CE003-00.de,en 1 / 1 PMF: 2369813 OP-113_fo1_2020.03.03

Original

UK Declaration of Conformity

Manufacturer Sartorius Lab Instruments GmbH & Co. KG
37070 Goettingen, Germany

declares under sole responsibility that the equipment

Device type **Electronic Mass Comparator | Power Supply**

Models **MCAvMC, | ZAG65US15**
v = 36, 66, 106, 206, 605, 1004, 1005

in the form as delivered fulfils all the relevant provisions of the following British Regulations and meets the applicable requirements of the British Designated Standards including any amendments valid at the time this declaration was signed listed below:

The Electromagnetic Compatibility Regulations 2016
UK Statutory Instruments 2016 No. 1091

BS EN 61326-1:2013; BS EN IEC 61326:2021

The Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012; UK Statutory Instruments 2012 No. 3032

BS EN IEC 63000:2018

The Supply of Machinery (Safety) Regulations 2008
UK Statutory Instruments 2008 No. 1597

BS EN ISO 12100:2010, BS EN 61010-1:2010+A1:2019*)

The person authorised to compile the technical file:

Sartorius UK Ltd.
Longmead Business Centre, Blenheim Road
KT19 9 QQ Epsom, Surrey, UK

Sartorius Lab Instruments GmbH & Co. KG
Goettingen, 2023-05-23

Mario Hespe
Managing Director

Halil Yildirim
Product Compliance Officer (SLI)

*: applied standard, which however is not harmonized for machines

20 Index

A

Accessories....11, 19, 93
  Connecting....51
Acclimatizing....49, 85, 86
Adjustment....63
  Calibration weight....87
  Internal....64
  With isoCAL Function....63
Ambient conditions....83, 85
Application....33
  Examples....69
Assigning a Password....56

B

Backup Battery....88

C

Cable. Siehe Connection cable
Calibration....63
Calibration weight....87
Cleaning
  Cleaning the Device....71
  Preparation....71
Cleaning Agents and Cleaning Procedures....90
Clock....88
Commissioning....50
Conduct in an Emergency....11
Configuring Print Profiles....59
Conformity....95
Connection cable
  For control unit....45, 50
  For Ethernet....46, 50
  Risk of tripping....12
Connections....17
Control unit
  Attaching or removing....44
  Connection cable for control unit....45, 50
  Setting up....49
Covers....20

D

Data storage device....88
Decommissioning....76
Device....13
  Components in the weighing chamber....15, 16
  Connections....17
  Draft shield....14
  Switching on or off....52
Dimensions and Weights....83
Disposal....82
Draft shield
  Dismantling....77
  Mounting....47
  Overview....14
Draft shield, manual
  Opening and Closing....61
Draft shield, motorized
  Opening or Closing....61
  Saving Functions....55

E

Electrical Equipment
  Connection Point for Power Supply Unit and Power Supply Cable....11
  Interfaces....86
  Power consumption....86
  Power supply....85
  Power Supply Unit and Power Supply Cable....11
  Safety and Electromagnetic Compatibility....86
Electromagnetic Compatibility....86
Emergency....11

F

Foreseeable Misuse....9
Functionality of the Device Parts....10

G

Glass Breakage....12

I

Installation....43
Installation Conditions....84
Intended Use....9
  Modifications to the Device....9
  Operating Conditions for the Device....9
  Repairs and Maintenance on the Device....10
Interfaces....86
isoCAL Function
  Adjustment....63
  Conditions....88
  Switch-off....55

L

Leveling....62
Licensable device functions....33
Linearization....63

M

Maintenance....10
Maintenance Schedule....72
Malfunctions....75
Materials....89
Metrological Data....91
Misuse....9
Models....7
Modifications to the Device....9

O

Opening the Help Function....53
Operating Conditions for the Device....9
Operation....61

P

Performing a QAPP Center Update 73
Performing a Software Update 72
Performing System Settings 53
Placing the device on its side and installing 45
Power consumption 86
Power supply 51, 85
Power supply unit and power supply cable
  Mounting 50
  Power supply 51
Power Supply Unit and Power Supply Cable 11
  Connection Point 11
Print
  Marking Values as Invalid 67
  Marking Values as Valid 67
  Printing Values 68
  Saving Values for Printout 66
Protective Caps 20

Q

Qualifications of Personnel 10

R

Recommissioning 72
Related Documents 7
Repairs and Maintenance on the Device 10

S

Safety Information on the Device ... 10
Safety of Electrical Equipment 86
Sartorius Service 94
Scope of Delivery 43
Selecting an Installation Site 44, 83, 84, 85
Setting up a Network Printer 57
Spare Parts 94
Storage and Shipping 81
Switching off 52
Switching on 52
Symbols on the Device 10
Symbols Used 8

T

Target Groups
  Description 7
  Qualifications of Personnel 10
Task 33
  Creating or Adding 65
  Ending 68
Technical data 83
Transport 12, 79
  Transporting 80
Transport lock 19
  Manufacturing 79
  Removing 48

U

Unpacking 44
Update
  QAPP Center 73
  Software 72
User Login or Logout 52

V

Viewing the Data Storage Device .. 68

W

Wait times 85, 86
Warm-up Time 61, 85, 86
Weighing 65
Weighing pan 15, 16, 83, 93
  Inserting 46
  Removing 76

Sartorius Lab Instruments GmbH & Co. KG
Otto-Brenner-Strasse 20
37079 Goettingen, Germany

Phone: +49 551 308 0
www.sartorius.com

The information and figures contained in these instructions correspond to the version date specified below.
Sartorius reserves the right to make changes to the technology, features, specifications and design of the equipment without notice.
Masculine or feminine forms are used to facilitate legibility in these instructions and always simultaneously denote all genders.

Copyright notice:
These instructions, including all components, are protected by copyright.
Any use beyond the limits of the copyright law is not permitted without our approval.
This applies in particular to reprinting, translation and editing irrespective of the type of media used.

Last updated:

09 | 2024

© 2024

Sartorius Lab Instruments GmbH & Co. KG
Otto-Brenner-Str. 20
37079 Goettingen, Germany

SP | Publication No.: WMC6032-e240903