


SLC22A4 Antibody (C-term)

Peptide Affinity Purified Rabbit Polyclonal Antibody (Pab)

Catalog #	Applications:	Reactivity:	Accessions:
AP5986b	WB, E	H	Q9H015 , NP_003050.2
Concentration:	Size:	Isotype:	Clone Name:
0.25 mg/ml	0.1 mg	Rabbit Ig	RB27501

Application Data:

Calculated MW: 62155 Da

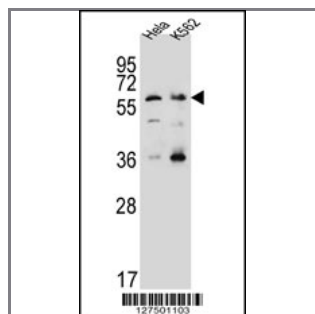

SLC22A4 Antibody (C-term) (Cat. #AP5986b) western blot analysis in HeLa, K562 cell line lysates (35ug/lane). This demonstrates the SLC22A4 antibody detected the SLC22A4 protein (arrow).

Gene ID:	Gene Symbol:
6583	SLC22A4

Other Names:

SLC22A4;ETT; OCTN1; UT2H

Target/Specificity:

This SLC22A4 antibody is generated from rabbits immunized with a KLH conjugated synthetic peptide between 521-550 amino acids from the C-terminal region of human SLC22A4.

Application Notes:

The suggested dilution is:

ELISA 1:1,000

Western blotting 1:100~500

Format:

Purified polyclonal antibody supplied in PBS with 0.09% (W/V) sodium azide. This antibody is purified through a protein A column, followed by peptide affinity purification.

Storage:

Maintain refrigerated at 2-8°C for up to 6 months. For long term storage store at -20°C in small aliquots to prevent freeze-thaw cycles.

Precautions:

This product is for research use only. Not for use in diagnostic or therapeutic procedures.
