INSTALLATION GUIDE

RephiQuatro U Pack 1（RO or distilled feed）

RephiLe Catalogue Number:
RR100Q101
Corresponding Millipore Part Numbers:
CPMQ004R1


Suitable Millipore Systems:
Milli-Q Plus
Milli-Q PF Plus
Milli-Q UF Plus
Milli-Q UV Plus
Milli-Q RG
Alpha-Q

This product is comparable to the corresponding Millipore cartridge. Please refer to Millipore's User Manual for a detailed instruction on installation.

A brief description of the installation process is described below:

1. Press "OPERATE/STANDBY" to let the system into Standby mode.
2. Turn the RECIRCULATION/PRODUCTION valve to PRODUCTION and allow any excess water to drain from the system. Then return the valve to the RECIRCULATION position; this will depressurize the purification cartridge.
3. Unlock the cartridge security bar by pressing down on the bar and swinging it to the left. Remove the exhausted cartridge from the system.
4. Open the new cartridge package; remove the two blue protective covers from the new cartridge, then wet the O-rings on the ports with pure water.
5. Snap-in the new cartridge by pushing it into its position; close the security bar to lock the cartridge to its position.
6. Degas the cartridge: Remove the final filter. Connect the degas connector to the water outlet and the sink. If cannot locate the degas connector, put a large container underneath the water dispense outlet.
7. Turn the RECIRCULATION / PRODUCTION valve to PRODUCTION. Press "OPERATE/STANDBY" to get the system into Operate mode and to start the pump. Allow the system to purge water and air until any trapped air is removed. Then allow the system to run to drain for 5 minutes.
8. Turn the RECIRCULATION/PRODUCTION valve to RECIRCULATION. Press "OPERATE/STANDBY" to put the system into Standby mode.
9. Remove the degas connector or beaker, allow the system to stand undisturbed overnight in Standby. This allows the new cartridge to hydrate.
10. Next morning, open the point-of use valve. Press "OPERATE/STANDBY" to get the system into Operate mode. Allow the system to run to drain for 5 minutes.
11. Close the point-of-use valve. Press "OPERATE/STANDBY" to put the system into Standby mode.
12. Proceed to the installation of a new final filter where required.

The system is ready for use.

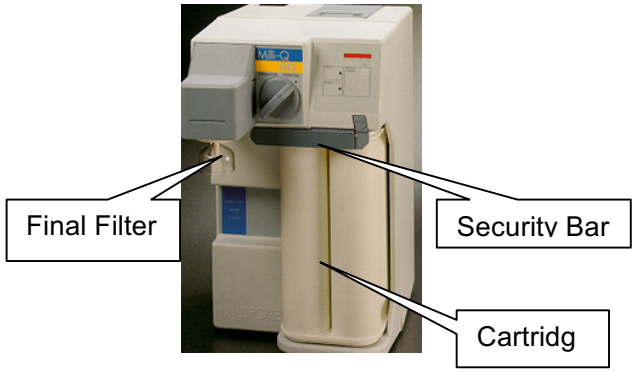

RephiSolo, RephiDuo and RephiQuatro are trademarks of RephiLe Bioscience, Ltd.
© 2010. RephiLe Bioscience, Ltd. All rights reserved.
All other trademarks are trademarks of their respective manufacturers.