INSTALLATION GUIDE

RephiDuo P Pack L1

RephiLe Catalogue Number:

RR500CPP1

Corresponding Millipore Part Number:

PRPK0L0S1

Suitable Millipore Systems:

RiOs 30 - 200
Elix 20 – 100. Elix Clinical


This product is comparable to the corresponding Millipore cartridge. Please refer to Millipore's User Manual for a detailed instruction on installation.

A brief description of the installation process is described below:

Attention: Do not remove the Cartridge without following the system instructions shown below.

Removing the used cartridge

1. Place the Water System into STANDBY Mode.

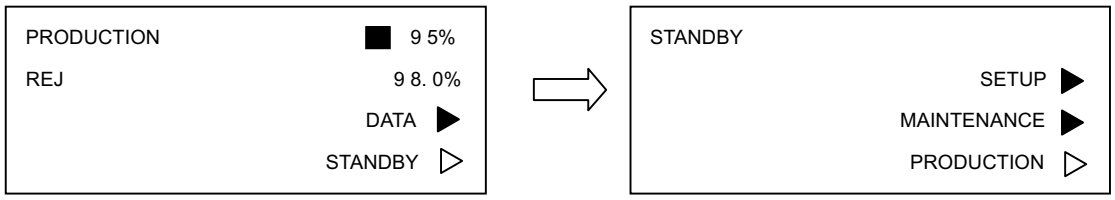



If an inlet solenoid valve is not installed before the cartridge, then close the manual feed water valve before the cartridge.

2. Make sure the 2 pieces of drain tubing are secured in the drain or a large container in placed directly under them to collect the water in the cartridge. The container needs to be large enough to hold the water in the cartridge.

Caution: Look at the pressure gauges at the top of the cartridge holder. If they are not 0 psi, perform the next step with caution.


3. SLOWLY open the Drain valve marked “DRAIN” located on the bottom of the cartridge holder. This will release any residual pressure as well as empty the cartridge.

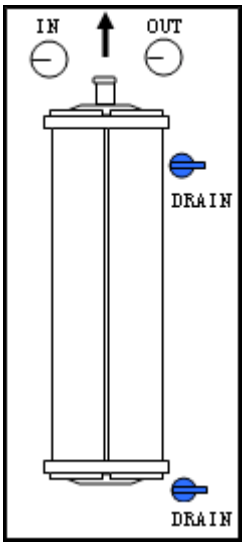


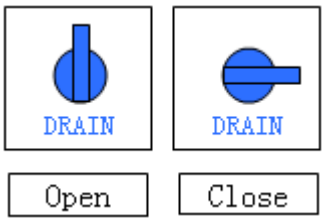

4. Open the valve marked “DRAIN” on the top of the cartridge holder. Allow all the water to drain from the cartridge.

5. Unlock the cartridge by lifting up the **Locking Handle** that holds the cartridge in place. Gently remove the used cartridge from the Water System. Close both of the valves marked “DRAIN”.

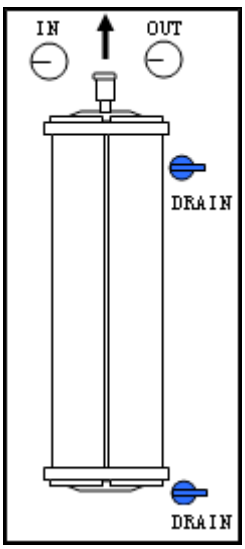

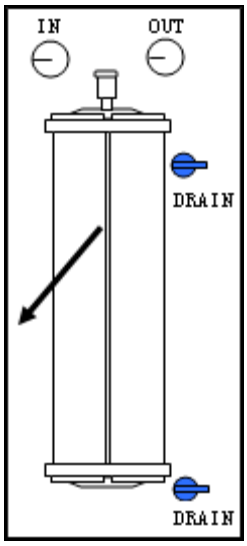

Installing the new cartridge


6. Remove the plastic covers on the 4 ports of the new cartridge.

7. Look inside the ports to locate the black rubber O-rings. Make sure these are pushed in so they are firmly in-place. Wet the O-rings with water.

8. Install the new cartridge with the **This End Up** arrow pointing up. **Note: Cartridge is directional. It MUST be installed with the "This End Up" arrow pointing up.**

9. Make sure all top and bottom ports are securely pushed in to their positions.

10. Lock cartridge into its place by pushing down the Locking Handle to its original position.

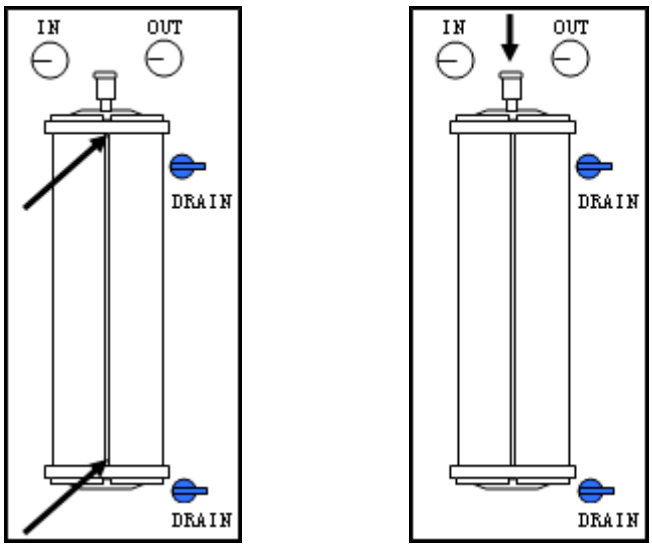

11. Make sure the two Drain valves are in closed position (in horizontal position).

No special rinsing is needed.

The system is ready for use.

RephiSolo, RephiDuo and RephiQuatro are trademarks of RephiLe Bioscience, Ltd.

© 2010. RephiLe Bioscience, Ltd. All rights reserved.

All other trademarks are trademarks of their respective manufacturers.