

HemoCue® Hb 201+ System



The Standard for Uncompromising Accuracy

Providing lab accuracy and ease of use, the HemoCue® Hb 201+ system has become a standard in Hb point-of-care testing. Healthcare providers around the world rely on the immediate results so they can make the right decisions when they need them most — right at the point of care.

With dedicated support and service as well as unmatched training and education based on 30 years of experience, you can count on HemoCue for the right solutions for all your needs.

Accuracy Starts With Us



Have Confidence in Your Answers at the Point of Care

- ▶ Precise factory calibration against the ICSH reference method
- ▶ Patented microcuvette technology with excellent lot-to-lot reproducibility
- ▶ Built-in self-test
- ▶ Blood-based liquid controls available
- ▶ Industry leading 3 year warranty


Get Easy Access to Lab-Quality Accuracy

- ▶ Capillary, venous or arterial whole blood sample
- ▶ Brief training with virtually no maintenance
- ▶ Link result with patient ID for medical record integration
- ▶ Printer interface



HemoCue® Hb 201+ System


Principle	Modified azidemethemoglobin reaction; dual wavelengths (570 nm and 880 nm) for compensation of turbidity
Calibration	Factory calibrated against the ICSH reference method; needs no further calibration
Sample Material	Capillary, venous or arterial whole blood
Measurement Range	0–25.6 g/dL (0–256 g/L, 0–15.9 mmol/L)
Results	Within 60 seconds
Sample Volume	~10 µL
Dimensions	160×43×85 mm (6.30×1.69×3.35 inches)
Weight	350 g (0.77 pounds) with batteries installed
Storage Temp.	Analyzer: 0–50 °C (32–122 °F) Microcuvettes: unopened 15–30 °C (59–86 °F); three-month open vial stability
Operating Temp.	15–30 °C (59–86 °F)
Power	AC Adapter or 4 AA batteries
Interface	Printer and HemoCue® Basic Connect™ including optional barcode scanner
Quality Control	Built-in self-test ; optional liquid controls
CPT Code	85018

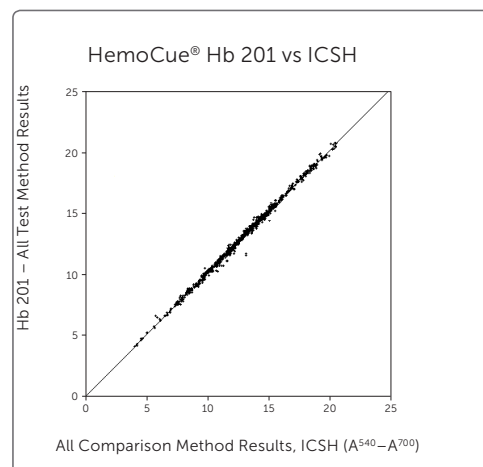


The CPT codes provided are based on AMA guidelines and are for informational purposes only. CPT coding is the sole responsibility of the billing party. Please direct any questions regarding coding to the payer being billed.

Product Item Number	VWR Catalog No.	Supplier Part No.
Hb 201+ Analyzer	10842-886	121721
Hb 201+ Microcuvettes (200/box)	10842-854	111716
Hb 201+ Individually Packaged Microcuvettes (100/box)	10842-836	111715
Basic Connect Software	10842-868	131074
HemoCue Hb 201+ Cuvette Holder	10842-832	131717
Adapter (12 volt power supply) (US Standard-all HemoCue devices except DM)	10842-874	130658
HemoCue Cleaner (5/box)	10842-876	139123

Three Simple Steps

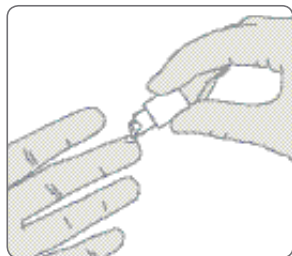

1 Fill microcuvette.

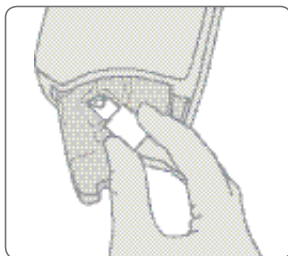

2 Place microcuvette into analyzer.



3 View results.


1.800.932.5000
vwr.com

Prices and product details are current when published; subject to change without notice. [Certain products may be limited by federal, state, provincial, or local regulations.] VWR makes no claims or warranties concerning sustainable/green products. Any claims concerning sustainable/green products are the sole claims of the manufacturer and not those of VWR International, LLC. All prices are in US dollars unless otherwise noted. Offers valid in US, void where prohibited by law or company policy, while supplies last. [VWR, the VWR logo and variations on the foregoing are registered (®) or unregistered trademarks and service marks, of VWR International, LLC and its related companies. All other marks referenced are registered by their respective owner(s).] Visit vwr.com to view our privacy policy, trademark owners and additional disclaimers. ©2016 VWR International, LLC. All rights reserved.

1116 Lit. No. 081858W

HemoCue has been a leader in point-of-care medical diagnostics for over 30 years. We specialize in giving healthcare providers lab-quality accuracy with results comparable to that of a clinical lab.

