INSTALLATION GUIDE

RephiQuatro P Pack

RephiLe Catalogue Number:
RR100CP01
Corresponding Millipore Part Numbers:
CPR0P0402


Suitable Millipore Systems:
Milli-RO 5/8/10 Plus
Milli-RO 30/60/90
Milli-RX 20/45/75
AFS 200/300

This product is comparable to the corresponding Millipore cartridge. Please refer to Millipore's User Manual for a detailed instruction on installation.

A brief description of the installation process is described below:

1. Press "OPERATE/STANDBY" for 2 seconds to let the system into **Standby** mode.
2. Push down and open the cartridge security bar.
3. Remove the exhausted cartridge from the system.
4. Open the new cartridge package, remove the two blue protective covers from the new cartridge, then wet the O-rings with pure water.
5. Snap-in the new cartridge by pushing it into its position, close the security bar and lock the cartridge to its position.
6. Press "OPERATE/STANDBY" for 2 seconds to get the system into Operate mode.

The system is ready for use.

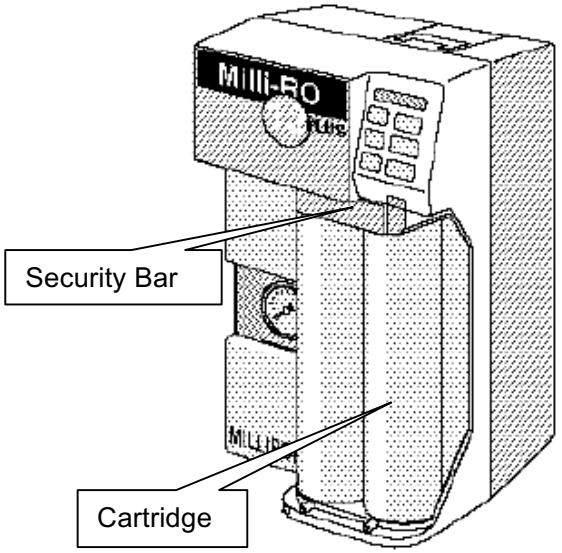

RephiSolo, RephiDuo and RephiQuatro are trademarks of RephiLe Bioscience, Ltd.
© 2010. RephiLe Bioscience, Ltd. All rights reserved.
All other trademarks are trademarks of their respective manufacturers.