RFM300-M Series Refractometer

QUICK START GUIDE
19-407'04 REV A, OCT-2023
Bellingham + Stanley, a Xylem brand, Xylem Analytics Germany GmbH, AM Achalaich 11, 82362 Weilheim, Germany
Email: sales.bs.uk@xylem.com

Thank you for purchasing this instrument. In order to ensure that this product provides many years of service please follow the guidance in this document. To register your instrument with the manufacturer and print a warranty certificate visit:
www.bellinghamandstanley.com

Unpacking
Carefully remove all of the packing material. It is recommended that the box and other packing materials are retained so that, should the need arise, the instrument can be safely returned to the supplier for service. Check that all parts listed below are present and that no transit damage has occurred. If any are damaged or missing, contact the supplier immediately.

Contents list

Quantity	Item	Code
1	RFM300-M Series Refractometer	See below
1	Mains lead set	55-109
1	Power supply	55-105
1	Quick Start Guides	19-407
2	Card type ID tag	37-530
2	Spare air vent filter	22-482

RFM300-M Refractometer code numbers:

Model	Complete refractometer	Refractometer module only
RFM330-M	19-35	19-335
RFM340-M	19-45	19-345

Positioning the system
- Dry and indoors.
- Capable of holding the instrument's weight.
- Away from draughty or hot equipment like fans or heaters.
- Out of direct sunlight or strong ambient light.
- Away from potential sources of interference, such as RFI generating equipment.
- Within reach of a power point.
- Not using a power circuit that also has large motors or noise generating equipment connected to it.

Do not block the air vents.

Mains connection
The power supply adapter is supplied with a moulded mains cord and plug to suit one of several socket types. For UK leads, replace the fuse only with the type indicated on the plug.

- Voltage 110 to 230 V ~ ±10%
- Maximum current 2 A
- Frequency 50 to 60 Hz

RISK OF ELECTRIC SHOCK:

- For indoor use only.
- Must be kept dry.
- Disconnect the equipment from the mains supply before unplugging the mains lead from the instrument.
- Do not cover, designed to operate with free air convection.
- No cleaning required.
- Do not jet wash

Switching on and off
Plug in the power supply and switch on the mains supply, after startup the instrument will show a loading screen indicating the software version. Once loaded, the instrument will be in standby mode. To switch on, press the green Power softkey.

To switch to standby mode press and hold the green Power softkey for several seconds and select Yes on the confirmation screen. It is then safe to turn off the mains supply and disconnect the power lead.

Manoeuvring through menus
The softkey functions change depending upon what inputs are required at that time. If a softkey has no text then that button is disabled.

Menu items can be selected by either:

- Pressing the Up/Down softkeys until it is highlighted and then pressing Select.
- Pressing the item number key.

Alternatively, a USB keyboard, available from Bellingham + Stanley, can be used to select items and navigate menus.

If the menu item contains a list of options to choose from then either:

- Select the menu item and then choose the option from the displayed list.
- Press the Forward/Back arrow keys to step through the options.

Menu items that can be modified by using the arrow keys are identified by arrows being displayed next to the menu items value.


Return to the Measurement Display by pressing Quit to step back through the menus.

Setup Wizard
When the instrument is switched on for the first time, a Setup Wizard will step through a number of questions so that the system can be configured for a particular application. First select the desired language by pressing the associated number key. Press the Next softkey to access each setting screen:

Copy user guides and PC software to USB: Makes copies of the full User Guides and PC software compatible with the instrument to a removable USB storage device. Insert a USB removable storage device before selecting. Full User Guides provide detailed instructions on instrument operation.
Set the standard measurement conditions: Make changes to the Measurement Settings.
Set clock: Set the instrument to local time.
Set the printer options: Specify how reading results will be printed.

When the desired conditions are set in each screen press the Quit softkey to proceed.

Auto Zero
When the instrument is turned on a Zero calibration will be requested by the instrument. This is to ensure that readings are as accurate as possible. Whilst the temperature stabilises the Reading Result section of the screen will display a warning message.

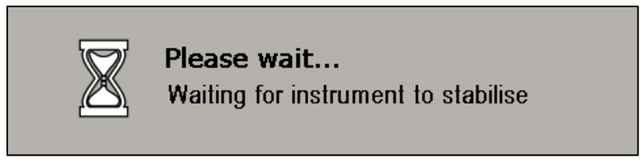Follow the onscreen instructions – ensure that the prism is thoroughly clean and that the presser is then lowered.

The instrument stores the details of the last sample used for Zero calibration – typically distilled water. If a different sample is to be used press the Change Sample button and enter the data as requested, the sample value must be entered in the scale used for calibration.

Once the sample details are correct apply the Zero sample to the instrument prism. The complete prism surface must be covered, 1ml is generally sufficient to achieve this. Lower the presser and press ok. Continuous readings will then be taken until the sample stabilises The Zero point will then be set.

This symbol is an internationally agreed indicator that the product bearing it should not be disposed of as general waste or garbage which might end up in landfill sites, but should instead be sent for special processing and/or recycling in those countries where appropriate legislation and facilities are in place.

Full User Guides can be downloaded from the instrument's Help Menu or Setup Wizard

Caution or warning.

Electrical hazard.


RFM300-M Series Refractometer

Bellingham + Stanley, a Xylem brand, Xylem Analytics Germany GmbH, AM Achalaich 11, 82362 Weilheim, Germany
Email: sales.bs.uk@xylem.com

QUICK START GUIDE
19-407'04 REV A, OCT-2023

Measurement Display

From the Measurement Display readings can be taken, viewed and saved.

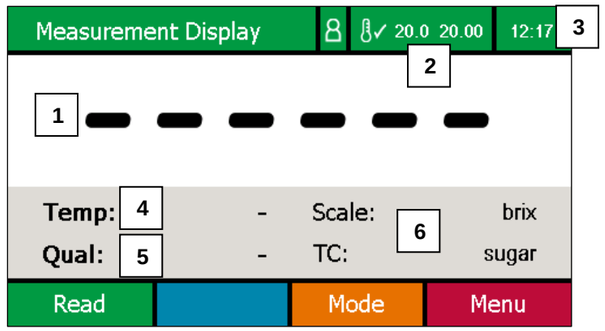1. Reading result: After taking a reading, the final result will be displayed here. Warning and alert messages will also appear here. For example, after the temperature set-point is changed, a Zero calibration will be requested.
2. Temperature control: the stability status, prism set temperature and current prism temperature.
3. Time: The current time.
4. Temperature: The temperature that the prism was at when the reading was recorded.
5. Quality: The Quality figure is derived from the optical pattern that falls upon the instrument detection system. A high Quality value indicates that the optical pattern is well-defined, making the reading signal easy to resolve. A low Quality value, caused for instance by an opaque sample or poor sample application, means a less well defined optical pattern so the reading will be less reliable. A Quality value of 100 is set when zeroing the instrument, which can then be used as a reference with which to compare other measured samples.
6. Reading Configuration: Displays the measurement settings for the currently displayed reading or, if a result is not currently displayed, the settings to be used for the next reading.

Softkeys


Read button: Initiates a reading. The complete prism surface must first be covered, 1ml is generally sufficient to achieve this.
Record button: Pushing the button will record the result to the Reading Log or printer depending on the print settings.
Mode button: Gives access to edit all Basic and Advanced reading settings e.g. scale and temperature compensation.
Menu button: Opens the instrument's Main Menu from where Calibration can be performed, saved results downloaded and all instrument settings can be changed.

Taking a reading

- Raise the presser and thoroughly clean the prism surface.
- Place a small amount of the sample to be measured on the prism. The complete prism surface must be covered, 1ml is generally sufficient.
- Lower the presser.
- Readings will start or press the Read button.

The instrument will display a progress bar. The reading will then be displayed using the scales and temperature compensations specified in the Mode Menu or current Method as shown on the Reading Configuration Panel.

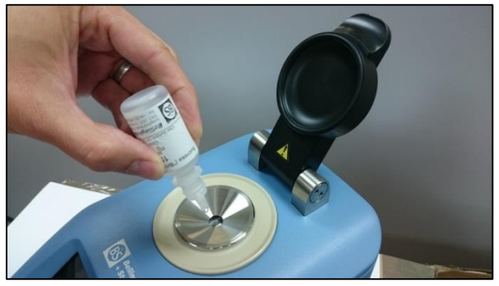Menu Structure

1. Calibration			
1. Zero			
2. Span (Top Calibration)			
3. Undo			
4. Report			
2. Data			
1. Saved Results			
2. System Log			
3. Copy PDFs to USB			
4. Move PDFs to USB			
5. Export Logs to USB			
3. Setup			
1. Methods	1. Add new…		
	2. Settings…	1. Use Methods	
		2. Hide preset	
		3. Selection Type	
2. Users	1. Add new…		
	2. Protection	1. low	
		2. medium	
		3. high	
		4. custom	
	3. administrator – 000		
	4. calibrator – 355		
	5. operator – 123		
3. Reading	1. Recurrence	1. single	
		2. continuous	
		3. auto single	
	2. Appl. Correction	1. none	
		2. offset	
		3. citric acid	
	3. Scales	1. Add new…	1. From Library
			2. From Data
	4. Temp. Comps.	1. Add new…	1. From Data
4.Recording	1. Mode	1. none 5. statistics	
		2. single 6. multi	
		3. save	
		4. single + save	
	2. Printer	1. none 5. network	
		2. serial	
		3. usb	
		4. pdf	
	3. Record After Read	1. no	
		2. yes	
	4. Traceability	1. Batch Codes	1. off
			2. number
			3. text
			4. dateinc
		2. Auto Increment	1. off
			2. on
		3. Log Users	1. off
			2. number
			3. text
			4. login
		4. Remember Text	1. none
			2. batch
			3. user
			4. batch and user
	5. Ports	1. Serial Port	1. Baud Rate
			2. Word Length
		2. Network	1. DHCP
			2. IP Address
			3. Subnet Mask
			4. Default Gateway
		3. CSV (Lims)	1. none
			2. serial
			3. network
			4. usb
5. System	1. Language		
	2. Time / Date	1. Date Format	1. dd/mm/yy
			2. mm/dd/yy
		2. Summer Time	1. no
			2. yes
		3. Set Clock…	
	3. Display	1. Brightness	
		2. Remote Display	
		3. Setup	
	4. Temperature Control	1. Mode	1. off
			2. on
		2. Disable Presser	1. no
			2. yes
	5. Clone From Backup		
4. Maintenance			
1. Maintenance Report			
2. Filter Change Reminder		1. Frequency	
		2. Last Changed	
		3. New Filters Fitted	
3. Backup to USB			
4. Update Software			
5. Help			
1. Setup Wizard			
2. Information			
3. Contact			
4. User Guide			
5. PC software			

This symbol is an internationally agreed indicator that the product bearing it should not be disposed of as general waste or garbage which might end up in landfill sites, but should instead be sent for special processing and/or recycling in those countries where appropriate legislation and facilities are in place.

Full User Guides can be downloaded from the instrument's Help Menu or Setup Wizard

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