


MGLL Antibody (N-term)

Peptide Affinity Purified Rabbit Polyclonal Antibody (Pab)

Catalog #	Applications:	Reactivity:	Accessions:
AP16876a	WB, E	H	Q99685 , NP_001003794.1 , NP_009214.1
Concentration:	Size:	Isotype:	Clone Name:
0.25 mg/ml	0.1 mg	Rabbit Ig	RB30365

Application Data:

Calculated MW: 33261 Da

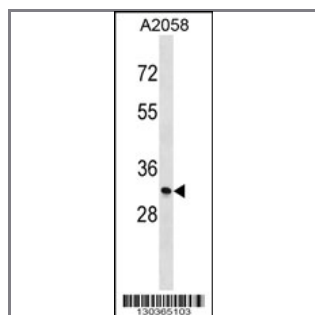

MGLL Antibody (N-term) (Cat. #AP16876a) western blot analysis in A2058 cell line lysates (35ug/lane). This demonstrates the MGLL antibody detected the MGLL protein (arrow).

Gene ID:	Gene Symbol:
11343	MGLL

Other Names:

Monoglyceride lipase;MGLL;

Target/Specificity:

This MGLL antibody is generated from rabbits immunized with a KLH conjugated synthetic peptide between 33-63 amino acids from the N-terminal region of human MGLL.

Application Notes:

The suggested dilution is:

ELISA 1:1,000

Western blotting 1:100~500

Format:

Purified polyclonal antibody supplied in PBS with 0.09% (W/V) sodium azide. This antibody is purified through a protein A column, followed by peptide affinity purification.

Storage:

MGLL Antibody (N-term) can be refrigerated at 2-8°C for up to 6 months. For long term storage, place the at -20°C.

Precautions:

MGLL Antibody (N-term) is for research use only and not for use in diagnostic or therapeutic procedures.
