VWR®

Data Acquisition Software

INSTRUCTION MANUAL

European Catalogue Number:

705 - 1086

Version: 1
Issued: 25. January 2018

Legal Address of Manufacturer

Europe

VWR International bvba

Researchpark Haasrode 2020

Geldenaaksebaan 464

B-3001 Leuven

+ 32 16 385011

http://be.vwr.com

Country of origin

Made in Germany

Table of contents

Warning 4
Safety Information 4
Delivery contents. 4
Intended use 4
Product Specifications 5
Software Installation 6
Data transmission 12
 The option window 12
Selecting the Saved Data Format 12
 Exel Selection 12
 Notepad Selection 13
 Availability of the IRIM 13
 Store location. 14
 User Interface Window with IRIM Connection 15
Transmission of Data 16
Technical service 16
Warranty. 17
Compliance with local laws and regulations 17

Warning

CAUTION

This manual is a reference guide for users of the Data Acquisition Software. The user should have prior knowledge of using Windows-based applications. The screenshots provided herein have been created using Windows XP and Excel 2007. If other operating systems or Office versions are being used, differences may occur.

Safety Information

Read this entire manual carefully before use!
The system may only be used by qualified persons!
Keep this manual for later reference!
Also observe important notices in order to avoid malfunctions and faults.

Delivery contents

Standard contents for VWR IRIM:

- ❑ 1 Instruction manual Infrared Interface Module
- ❑ 1 Instruction manual Data Acquisition Software
- ❑ 1 Instruction manual Device Firmware Update (DFU)
- ❑ 4 batteries (Type AA)
- ❑ 4 screws
- ❑ 1 Screwdriver
- ❑ 1 CD-ROM
- ❑ 1 USB-cable

Intended use

The Data Acquisition Software is designed to receive data from the instrument and store it either in an Excel-Spreadsheet or as a .txt-file. Depending on the instrument being used for transmission, an IRIM (InfraRed interface Module) or an RS232 computer interface is required. If the latter is not available, an RS232 to USB adapter can also be used.

Always adhere to the limits specified in chapter "Product Specifications" (page 5). To ensure correct and intended usage, the instructions have to be read and understood.

The instrument is not used as intended if

- it is used in application areas not stated in these instructions.
- it is used in conditions which deviate from the ones described in these instructions.

Product Specifications

Processor	Pentium® 4/M or similar
RAM	512 MB
Screen resolution	1280 x 1024 pixels
Operating system	Windows® XP, Windows Vista, Windows 7
Excel version *	Excel® 2003 or higher
Disk Space	90 MB

* If an Excel version older than 2003 is being used or Excel software has not been installed, no Excel option will be available when the Software is started (see chapter 'The option window' page 12).

Software installation

The Data Acquisition Software is available on CD or from the download area of the VWR homepage: www.vwr.com. If loading via the CD, the installation starts automatically when the CD is inserted into the CD-ROM drive of the computer.

Note: If an older version of the Data Acquisition Software is already installed, please remove the program before installing the new software. Please remove National Instrument Software and Photometer Data Acquisition.

On automatic start-up of the setup.exe, the location for program installation has to be confirmed. It is recommended to confirm the default file path by pressing the "Next" button.

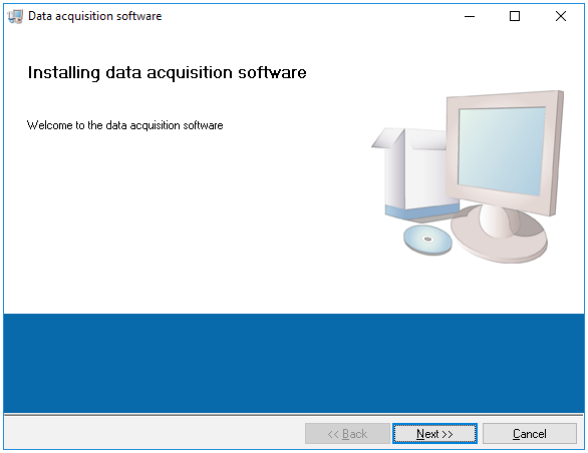
Image 1

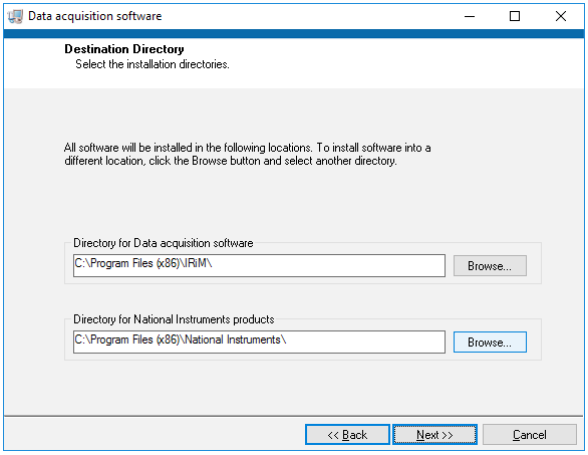
Image 2

On confirmation, two license agreements will appear which must be accepted and confirmed by pressing the "Next" Button.

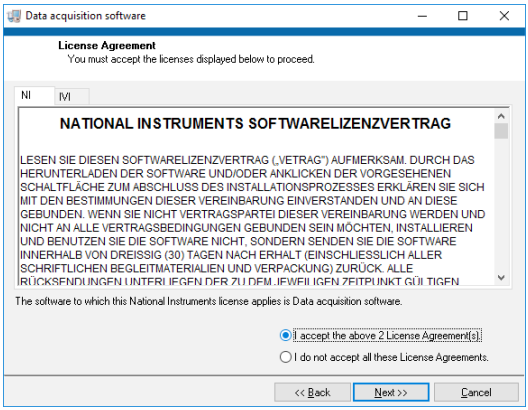
Image 3

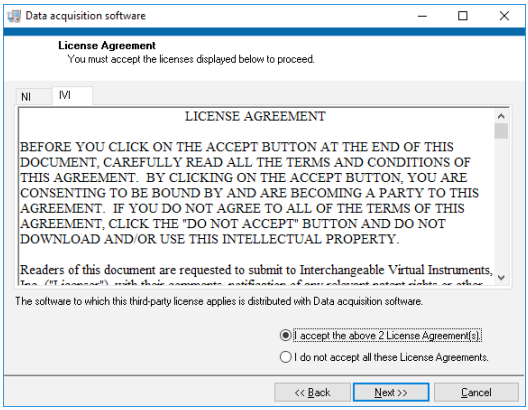
Image 4

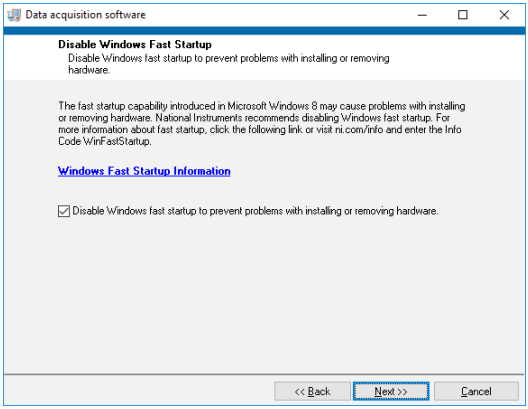
Image 5

Finally, a summary of the program to be installed is displayed. To start the installation press the "Next" button.

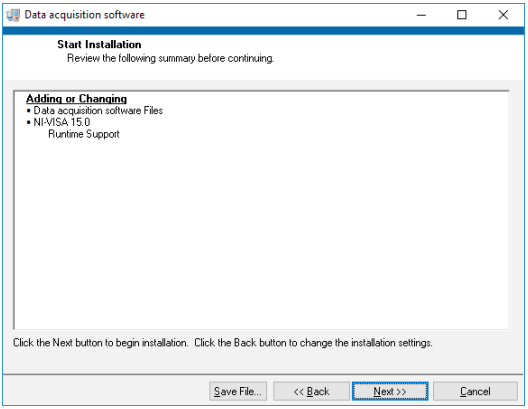
Image 6

The installation starts.

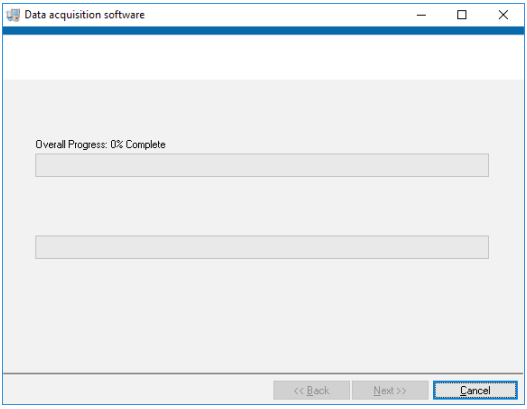
Image 7

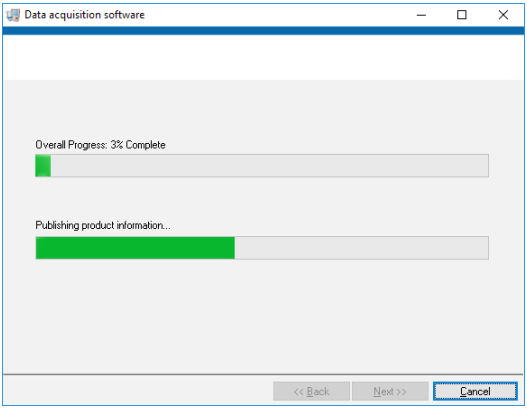
Image 8

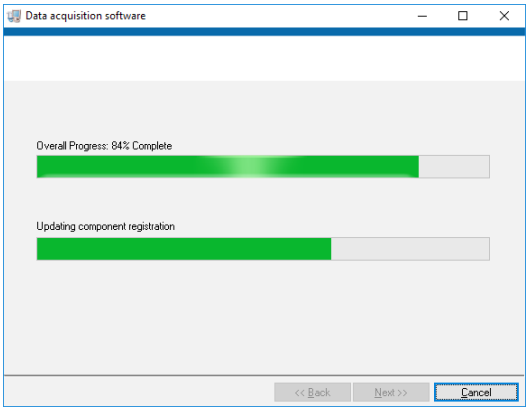
Image 9

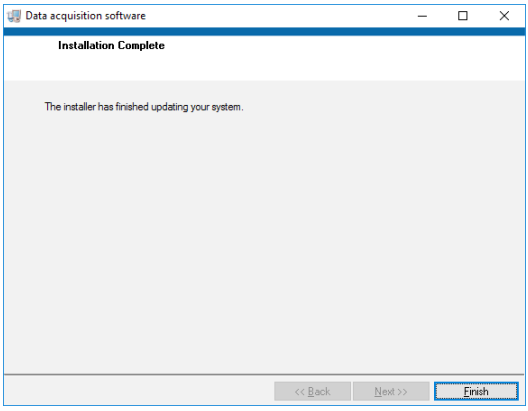
Image 10

Once the installation has been successfully carried out, pressing the "Finish" button ends the procedure.

After the installation, a Restart of the PC is necessary.

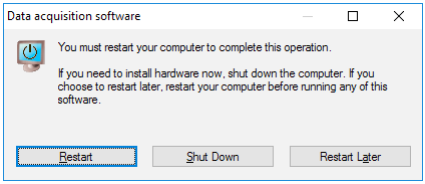
Image 11

On restarting the PC, a link will be created on the desktop which can, in future, be used to start the software.

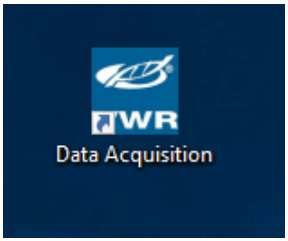
Image 12

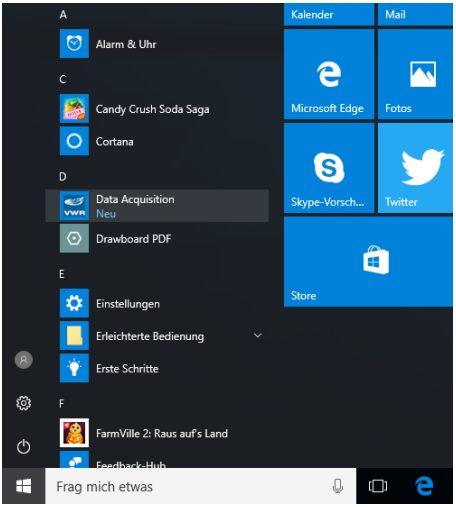
Image 13

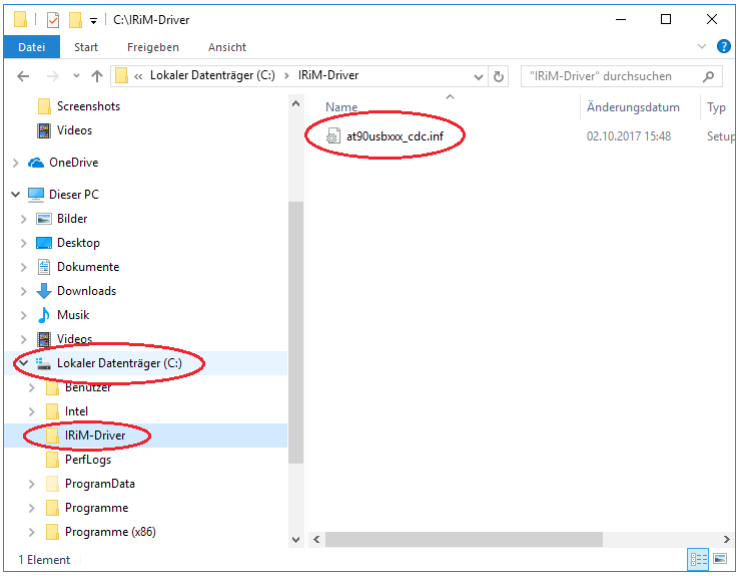
Image 14

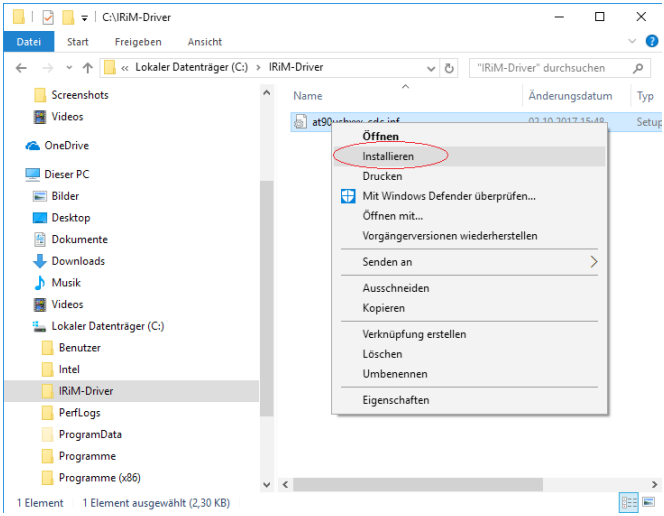
Image 15

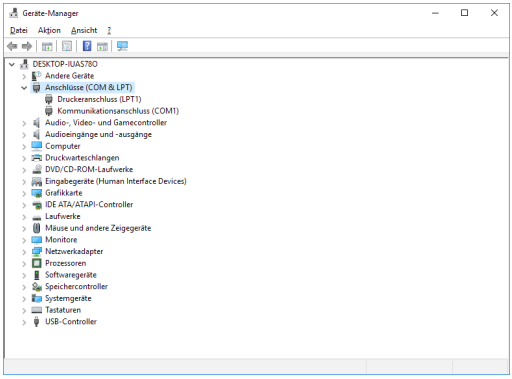
Image 16

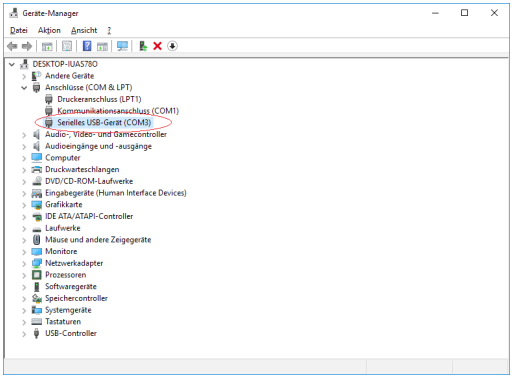
Image 17

Data transmission

The option window

Once the software has been started, an option window appears. Here you may enter your User Name, select your preferred language and request which type of report is required. You must also select which instrument you are using for the data transmission.

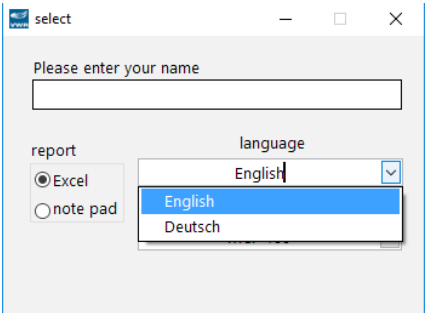
Image 18

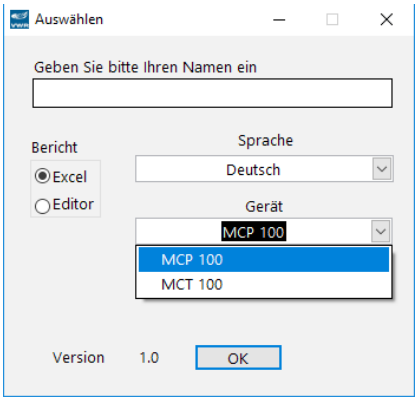
Image 19

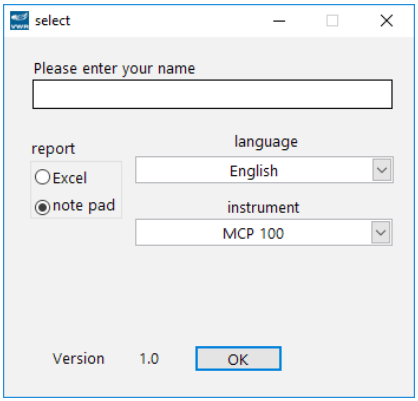
Image 20

Selecting the Saved Data Format

The downloaded data can be stored in two kinds of file format: Excel and Notepad as selected under chapter 'The option window'see page 5.

Excel Selection

If the Excel format option was selected in Section 2.1, an Excel Template automatically opens when the Save button is pressed. This template will be used to store the transmitted data.

Notepad Selection

If the Notepad format option was selected in chapter 'The option window'see page 5, the data will automatically be stored as a .txt file.

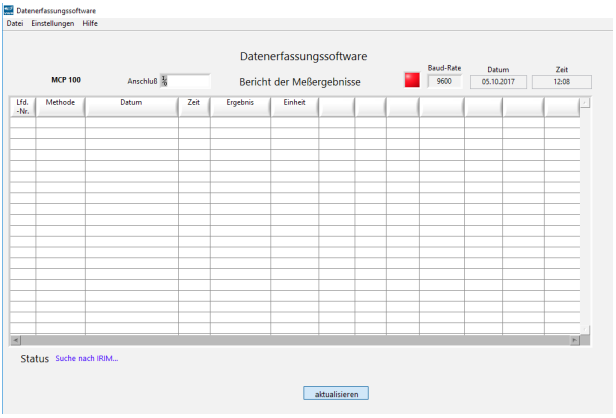
Image 21

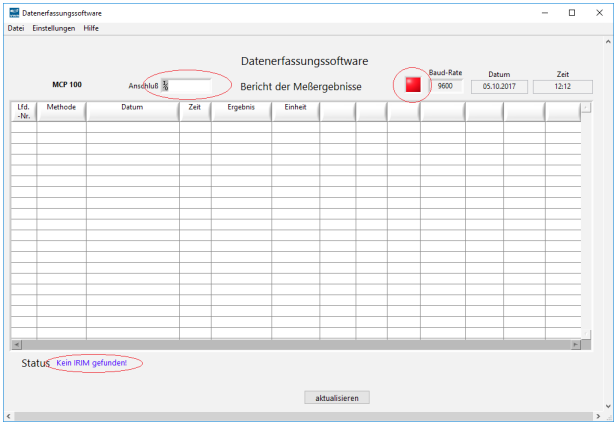
Image 22

Availability of the IRIM

After confirming the selection as shown in chapter 'The option window'see page 5, the software will automatically verify the IRIM's availability.

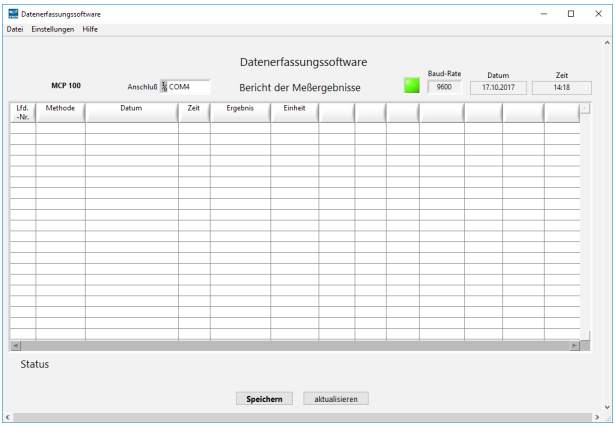
Image 23

If this verification fails, the following messages will appear in the status box.

In this case, the IRIM must be reconnected. To establish this has been successful, a new search must be carried out by clicking the Refresh button. Once the IRIM has been found to be correctly connected, the storage location window opens. Please confirm the storage location or select an alternative accordingly and press “save”.

Store location

Once the selection has been confirmed, a window appears where the file name and type may be entered. Please note, you may save data either as an Excel Spreadsheet or a Notepad *.txt file.

In Excel, if an existing file name is selected, a new worksheet will be created within this file. If you select to save to an existing Notepad file, the newly transmitted data will be attached to the already existing .txt file.

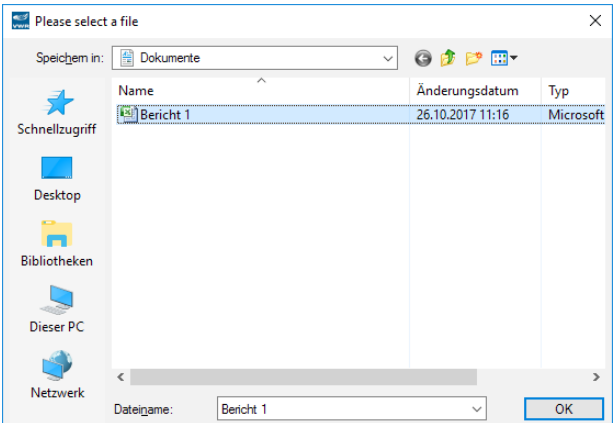
Image 24

User Interface Window with IRIM Connection

The green status LED indicates an established and correct connection of the IRIM. The COM interface to which the IRIM is attached is also displayed. If the Excel Spreadsheet Data Storage Option was selected, then the "save" button will also be shown. This will not appear if you chose to save the data in Notepad as the data will be saved automatically (Example: MCP 100).

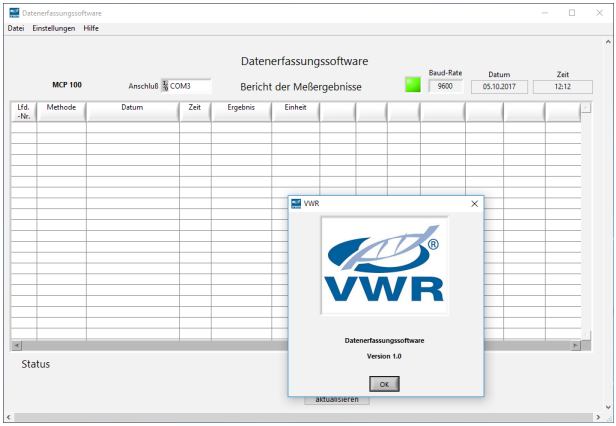
Image 25

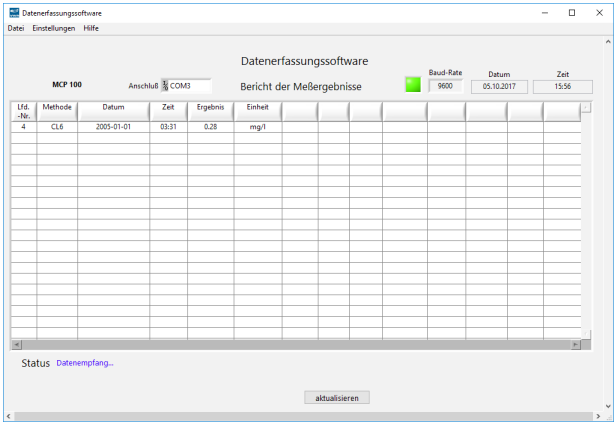
Image 26

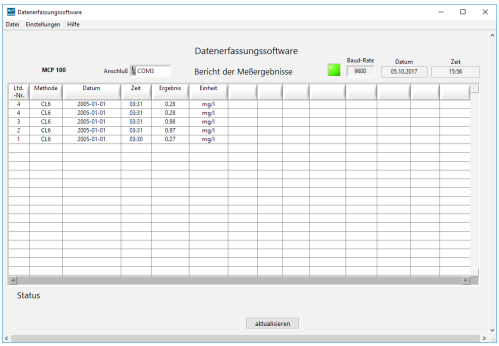
Image 27

Transmission of Data

The data can now be transmitted from the instrument to the PC. For details on how to do this, please refer to the instrument's Operation Manual. Once transmission has begun, the status row displays "Data receiving". The transmitted data will be displayed in the User Interface Window (Example: MCP 100).

On pressing the "save" button, the data will be transmitted to an Excel Spreadsheet where it can be stored.

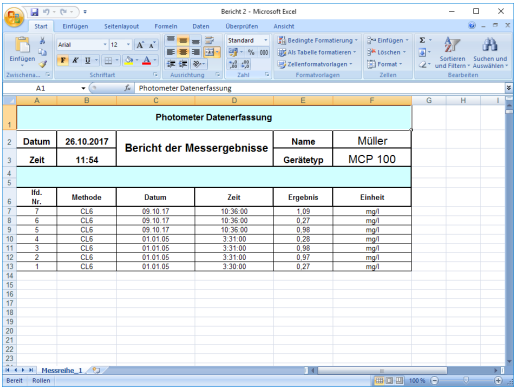
Image 28

In the case of Notepad, the data is automatically stored as a .txt-file.

Technical service

Web Resources

Visit the VWR website at www.vwr.com for:

- Complete technical service contact information
- Access to the VWR Online Catalogue, and information about accessories and related products
- Additional product information and special offers

Contact us For information or technical assistance contact your local VWR representative or visit. www.vwr.com.

Warranty

VWR warrants that this product will be free from defects in material and workmanship for a period of two (2) years from date of delivery. If a defect is present, VWR will, at its option and cost, repair, replace, or refund the purchase price of this product to the customer, provided it is returned during the warranty period. This warranty does not apply if the product has been damaged by accident, abuse, misuse, or misapplication, or from ordinary wear and tear. If the required maintenance and inspection services are not performed according to the manuals and any local regulations, such warranty turns invalid, except to the extent, the defect of the product is not due to such non-performance.

Items being returned must be insured by the customer against possible damage or loss. This warranty shall be limited to the aforementioned remedies. IT IS EXPRESSLY AGREED THAT THIS WARRANTY WILL BE IN LIEU OF ALL WARRANTIES OF FITNESS AND IN LIEU OF THE WARRANTY OF MERCHANTABILITY.

Compliance with local laws and regulations

The customer is responsible for applying for and obtaining the necessary regulatory approvals or other authorisations necessary to run or use the Product in its local environment. VWR will not be held liable for any related omission or for not obtaining the required approval or authorisation, unless any refusal is due to a defect of the product.

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Local VWR offices in Europe and Asia Pacific

Austria
VWR International GmbH
Graumanngasse 7
1150 Vienna
Tel.: +43 1 97 002 0
Fax: +43 1 97 002 600
Email: info.at@vwr.com

Belgium
VWR International bvba
Researchpark Haasrode 2020
Geldenaaksebaan 464
3001 Leuven
Tel.: 016 385 011
Fax: 016 385 385
Email: vwr.be@vwr.com

China
VWR International China Co., Ltd
Shanghai Branch
Room 256, No. 3058 Pusan Road
Pudong New District
Shanghai, 200123
Tel.: +86-21 589 868 88
Fax: +86-21 585 588 01
Email: info_china@vwr.com

Czech Republic
VWR International s. r. o.
Veetee Business Park
Pražská 442
CZ - 281 67 Stříbrná Skalice
Tel.: +420 321 570 321
Fax: +420 321 570 320
Email: info.cz@vwr.com

Denmark
VWR International A/S
Tobaksvejen 21
2860 Søborg
Tel.: 43 86 87 88
Fax: 43 86 87 90
Email: info.dk@vwr.com

Finland
VWR International Oy
Valimotie 9
00380 Helsinki
Tel.: +358 9 80 45 51
Fax: +358 9 80 45 52 00
Email: info.fi@vwr.com

France
VWR International S.A.S.
Le Périgares – Bâtiment B
201, rue Carnot
94126 Fontenay-sous-Bois cedex
Tel.: 0 825 02 30 30 (0,18 EUR TTC/min)
Fax: 0 825 02 30 35 (0,18 EUR TTC/min)
Email: info.fr@vwr.com

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India
VWR Lab Products Private Limited
No.139. BDA Industrial Suburb,
6th Main, Tumkur Road,
Peenya Post,
Bangalore, India – 560058
Tel.: +91-80-28078400
Fax +91-80-41117120
Email: vwr_india@vwr.com

Ireland / Northern Ireland
VWR International Ltd /
VWR International (Northern Ireland) Ltd
Orion Business Campus
Northwest Business Park
Ballycoolin
Dublin 15
Tel.: 01 88 22 222
Fax: 01 88 22 333
Email sales.ie@vwr.com

Italy
VWR International S.r.l.
Via San Giusto 85
20153 Milano (MI)
Tel.: 02-3320311/02-487791
Fax: 02-332031307/02-40090010
Email: info.it@vwr.com

The Netherlands
VWR International B.V.
Postbus 8198
1005 AD Amsterdam
Tel.: 020 4808 400
Fax: 020 4808 480
Email: info.nl@vwr.com

Norway
VWR International AS
Haavard Martinsens vei 30
0978 Oslo
Tel.: 22 90 00 00
Fax: 815 00 940
Email: info.no@vwr.com

Poland
VWR International Sp. z o.o.
Limbowa 5
80-175 Gdansk
Tel.: 058 32 38 210
Fax. 058 32 38 205
Email: info.pl@vwr.com

Portugal
VWR International - Material de Laboratório, Lda
Centro Empresarial de Alfragide
Rua da Indústria, nº 6
2610-088 Alfragide
Tel.: 21 3600 770
Fax: 21 3600 798/9
Email: info.pt@vwr.com

Singapore
VWR Singapore Pte Ltd
18 Gul Drive
Singapore 629468
Tel: +65 6505 0760
Fax: +65 6264 3780
Email: sales.sg@vwr.com

Spain
VWR International Eurolab S.L.
C/ Tecnología 5-17
A-7 Llinars Park
08450 - Llinars del Vallès
Barcelona
Tel.: 902 222 897
Fax: 902 430 657
Email: info.es@vwr.com

Sweden
VWR International AB
Fagerstagatan 18a
163 94 Stockholm
Tel.: 08 621 34 00
Fax: 08 621 34 66
Email: info.se@vwr.com

Switzerland
VWR International GmbH
Lerzenstrasse 16/18
8953 Dietikon
Tel.: 044 745 13 13
Fax: 044 745 13 10
Email: info.ch@vwr.com

UK
VWR International Ltd
Customer Service Centre
Hunter Boulevard
Magna Park
Lutterworth
Leicestershire
LE17 4XN
Tel.: 0800 22 33 44
Fax: 01455 55 85 86
Email: uksales@vwr.com