



Carboy - Narrow Mouth



Carboy - Wide Mouth



NA Cat#	Europe Cat#	Description	Capacity (Liter)	Brin Capacity (Liter)	Material	Autoclavable	Sterile (Gamma)	Cap Size (mm)	Neck Internal Dia (mm)	General wall Thickness (mm)	Product Size mm (±2MM)	Height With out Cap mm (±2MM)
10755-100	216-1708	Carboy, Round, Wide Mouth, 10 L - LDPE	10 L	11.5	LDPE	NO	-	Ø100	Ø85.5	2.5	Ø249	331
10755-102	216-1730	Carboy, Round, Wide Mouth, 10 L - PP	10 L	11.5	PP	Yes	-	Ø100	Ø85.5	2.5	Ø249	331
75795-690	216-1919	Carboy, Round, Wide Mouth, 15 L - LDPE	15 L	16.74	LDPE	NO	-	Ø100	Ø85.5	2.5	Ø285	372
10755-122	216-1733	Carboy, Round, Wide Mouth, 20 L - PP	20 L	22	PP	Yes	-	Ø100	Ø85.5	2.5	Ø285	472
10755-098	216-1736	Carboy, Round, Narrow Mouth, 10 L - HDPE	10 L	11.5	HDPE	NO	-	Ø83	Ø65	2.5	Ø249	376
10755-084	216-1704	Carboy, Round, Narrow Mouth, 10 L - LDPE	10 L	11.5	LDPE	NO	-	Ø83	Ø65	2.5	Ø249	376
10755-086	216-1703	Carboy, Round, Narrow Mouth, 10 L - PP Autoclavable	10 L	11.5	PP	Yes	-	Ø83	Ø65	2.5	Ø249	376
75795-688	216-1918	Carboy, Round, Narrow Mouth, 15 L - LDPE	15 L	17.45	LDPE	NO	-	Ø83	Ø65	2.5	Ø285	415
10755-110	216-1737	Carboy, Round, Narrow Mouth, 20 L, HDPE	20 L	22	HDPE	NO	-	Ø83	Ø65	2.5	Ø285	518
10755-112	216-1706	Carboy, Round, Narrow Mouth, 20 L - LDPE	20 L	22	LDPE	NO	-	Ø83	Ø65	2.5	Ø285	518
10755-114	216-1701	Carboy, Round, Narrow Mouth, 20 L - PP Autoclavable	20 L	22	PP	Yes	-	Ø83	Ø65	2.5	Ø285	518
75795-674	216-1908	Carboy, Round, Narrow Mouth, 25 L - LDPE	25 L	27.6	LDPE	NO	-	Ø83	Ø65.5	2.5	Ø285	583
10755-124	216-1738	Carboy, Round, Narrow Mouth, 50 L - LDPE	50 L	52	LDPE	NO	-	Ø83	Ø65	2.5	Ø376	665
10755-126	216-1702	Carboy, Round, Narrow Mouth, 50 L - PP Autoclavable	50 L	52	PP	Yes	-	Ø83	Ø65	2.5	Ø376	665
10755-121	216-1729	Carboy, Round, Narrow Mouth, 20 L, HDPE - Sterile	20 L	22	HDPE	NO	Yes	Ø83	Ø65	2.5	Ø285	518
75795-776	Not Applicable	Carboy, Round, Narrow Mouth, 50 L, HDPE - Sterile	50 L	52	HDPE	NO	Yes	Ø83	Ø65	2.5	Ø376	665
10755-092	216-1709	Carboy Round, Heavy Duty Round, 10L - PP Autoclavable	10 L	11.25	PP	Yes	-	Ø83	Ø65	3	Ø249	376
10755-104	216-1707	Carboy Round, Heavy Duty Round, 20L - PP Autoclavable	20 L	22	PP	Yes	-	Ø83	Ø65	3	Ø285	518
76210-660	216-2180	Carboy, Round, Narrow Mouth, 50 L, HDPE	50 L	52	HDPE	NO	-	Ø83	Ø65	2.5	Ø376	665

Spare : Refer NA Cat# 10820-450 / EU Cat# 216-1840 for Narrow Mouth Cap.

Refer NA Cat# 10820-436 / EU Cat# 216-1841 for Wide Mouth Cap.

Refer NA Cat# 10820-438 / EU Cat# 216-1842 for Narrow Mouth Filling and Vent Cap.

Material

Carboy : Refer above table for material type, 100 % Virgin, FDA 21 CFR Compliance, USP Class VI passed.

Closure : Polypropylene, 100 % Virgin, FDA 21 CFR Compliance, USP Class VI passed With Food grade colour.

BSE/TSE Statement: These product is manufactured without the use of Raw materials of animal origin.

Latex Free Statement: No materials containing latex or natural rubber are used in the manufacturing, handling and packaging processes for this product.

Caution: Applicable only for PP Products

1. Before Autoclaving separate the caps and Spigot / Stopcock threads.
2. Do not stack products, this may cause vacuum and may lead to physical damage / collapse of product.

Material Recycle Code



Description: Carboy without Spigot

Drawing:
Not to Scale



Release: 12/23/16

Revision: 3



NA Cat#	Europe Cat#	Description	Capacity (Liter)	Brim Capacity (Liter)	Material	Autoclavable	Sterile (Gamma)	Cap Size (mm)	Neck Internal Dia (mm)	General wall Thickness (mm)	Product Size mm (±2MM)	Height With out Cap mm (±2MM)
75795-676	216-1909	Carboy Round, Spigot, 5L - LDPE	5L	6L	LDPE	NO	-	Ø53	Ø38.8	2.0	Ø180	350
10755-088	216-1731	Carboy Round, Spigot, 10L - LDPE	10L	11	LDPE	NO	-	Ø83	Ø65	2.5	Ø249	376
10755-090	216-1710	Carboy Round, Spigot, 10L - PP	10L	11	PP	Yes	-	Ø83	Ø65	2.5	Ø249	376
10755-116	216-1717	Carboy Round, Spigot, 20L - LDPE	20L	21.5	LDPE	NO	-	Ø83	Ø65	2.5	Ø285	518
10755-118	216-1705	Carboy Round, Spigot, 20L - PP	20L	21.5	PP	Yes	-	Ø83	Ø65	2.5	Ø285	518
75795-698	216-1928	Carboy Round, Spigot, 25L - HDPE	25L	26.7	HDPE	NO	-	Ø83	Ø65	2.5	Ø285	583
75795-682	216-1899	Carboy Round, Spigot, 25L - LDPE	25L	27.6	LDPE	NO	-	Ø83	Ø65	2.5	Ø285	583
10755-128	216-1734	Carboy Round, Spigot, 50L - LDPE	50L	52	LDPE	NO	-	Ø83	Ø65	2.5	Ø376	665
10755-130	216-1718	Carboy Round, Spigot, 50L - PP	50L	52	PP	Yes	-	Ø83	Ø65	2.5	Ø376	665
10766-762	216-1745	Spare stopcock / Spigot for Carboy	-	-	-	Yes	-	-	-	-	-	-

Spare : Refer NA Cat# 10820-450 / EU Cat# 216-1840 for Narrow Mouth Cap.

Refer NA Cat# 10820-436 / EU Cat# 216-1841 for Wide Mouth Cap.

Refer NA Cat# 10820-438 / EU Cat# 216-1842 for Narrow Mouth Filling and Vent Cap.

Refer NA Cat# 10766-762 / EU Cat# 216-1745 for Stopcock / Spigot / Tap Set

Material

Carboy : Refer above table for material type, 100 % Virgin, FDA 21 CFR Compliance, USP Class VI passed.

Closure : Polypropylene, 100 % Virgin, FDA 21 CFR Compliance, USP Class VI passed
With Food grade colour.

BSE/TSE Statement: These product is manufactured without the use of Raw materials of animal origin.

Latex Free Statement: No materials containing latex or natural rubber are used in the manufacturing, handling and packaging processes for this product.

Spare Spigot Set



Spigot



Use internal Diameter 16mm Tube for without Connector application

Tabulation Connector



- Use internal Diameter 6mm – 8mm Tube for with Connector application.
- Connector Size towards Tap fitting 13.5mm

Stopper



Caution: Applicable only for PP Products

1. Before Autoclaving separate the caps and Spigot / Stopcock threads (If Applicable).
2. Do not stack products, this may cause vacuum and may lead to physical damage / collapse of product.

Material Recycle Code



Description: Carboy with Spigot

Drawing:
Not to Scale



Release: 9/26/2018

Revision: 4



NA Cat#	Europe Cat#	Description	Capacity (Liter)	Brim Capacity (Liter)	Material	Autoclavable	Neck Internal Diameter (mm)	General wall Thickness (mm)	Product Size mm (±2MM)	Height With out Cap mm (±2MM)
75795-692	216-1920	Carboy Round with Sanitary neck-76mm, PP, 10L	10 L	12.4L	PP	Yes	Ø65	2.5	Ø249	376
75795-678	216-1897	Carboy Round with Sanitary Neck-76mm, PP, 20L	20L	23.5	PP	Yes	Ø65	2.5	Ø285	518
75795-684	216-1910	Carboy Round with Sanitary neck-76mm, PP, 50L	50L	54	PP	Yes	Ø65	2.5	Ø376	665
76000-976	-	Carboy Round with Sanitary neck-76mm, PP, 50L - Sterile	50L	54	PP	Yes	Ø65	2.5	Ø376	665

Material

Carboy : Refer above table for material type, 100 % Virgin, FDA 21 CFR Compliance, USP Class VI passed.

BSE/TSE Statement: These product is manufactured without the use of Raw materials of animal origin.

Latex Free Statement: No materials containing latex or natural rubber are used in the manufacturing, handling and packaging processes for this product.

Caution: Applicable only for PP Products

1. Before Autoclaving separate the caps and Spigot / Stopcock threads / Connectors (If Applicable).
2. Do not stack products, this may cause vacuum and may lead to physical damage / collapse of product.

Material Recycle Code



Description: Carboy with Sanitary Neck

Drawing:
Not to Scale



Release: 12/23/16

Revision: 1



NA Cat#	Europe Cat#	Description	Capacity (Liter)	Brim Capacity (Liter)	Material	Autoclavable	Cap Size (mm)	Neck Internal Dia (mm)	General wall Thickness (mm)	Product Size mm (±2MM)	Height With out Cap mm (±2MM)
75795-694	216-1926	Carboy, Narrow mouth, with Cap & Sanitary Flange 41mm connection, PP, 10L	10 L	12.5L	PP	Yes	Ø83	Ø65	2.5	295 X Ø 249	376
75795-686	216-1916	Carboy, Narrow mouth, with Cap & Sanitary Flange 41mm connection, PP, 50L	50L	54	PP	Yes	Ø83	Ø65	2.5	423 X Ø376	665

Spare : Refer NA Cat# 10820-450 / EU Cat# 216-1840 for Narrow Mouth Cap.
Refer NA Cat# 10820-438 / EU Cat# 216-1842 for Narrow Mouth Filling and Vent Cap.

Material

Carboy : Refer above table for material type, 100 % Virgin, FDA 21 CFR Compliance, USP Class VI passed.

BSE/TSE Statement: These product is manufactured without the use of Raw materials of animal origin.

Latex Free Statement: No materials containing latex or natural rubber are used in the manufacturing, handling and packaging processes for this product.

Caution: Applicable only for PP Products

1. Before Autoclaving separate the caps and Spigot / Stopcock threads / Connectors (If Applicable).
2. Do not stack products, this may cause vacuum and may lead to physical damage / collapse of product.

Material Recycle Code



Description: Carboy with Sanitary Flange

Drawing:
Not to Scale



Release: 2/27/18

Revision: 2



NA Cat#	Europe Cat#	Description	Capacity (Liter)	Brim Capacity (Liter)	Material	Autoclavable	Sterile (Gamma)	Cap Size (mm)	Neck Internal Dia (mm)	General wall Thickness (mm)	Product Size mm (±2MM)	Height With out Cap mm (±2MM)
10755-094	216-1719	Carboy Rectangular Spigot, HDPE, 10L	10 L	11.5 L	HDPE	NO	-	Ø100	Ø87.3	2.5	284 X 176	389
10755-096	216-1735	Carboy Rectangular Spigot - PP, 10L	10 L	11.5 L	PP	Yes	-	Ø100	Ø86.3	2.5	284 X 176	389
10755-106	216-1700	Carboy Rectangular Spigot - HDPE, 20L	20 L	22 L	HDPE	NO	-	Ø100	Ø86.3	2.5	340 X 209	486
10755-108	216-1732	Carboy Rectangular Spigot , PP, 20L	20 L	22 L	PP	Yes	-	Ø100	Ø85.7	2.5	340 X 209	486
75795-680	216-1898	Carboy Rectangular, PP, 10L	10 L	12.4L	PP	Yes	-	Ø100	Ø86.3	2.5	284 X 176	389
75795-672	216-1907	Carboy Rectangular, PP, 20L	20L	22 L	PP	Yes	-	Ø100	Ø85.7	2.5	340 X 209	486
75795-696	216-1927	Carboy Rectangular, HDPE, 20L	20L	22 L	HDPE	NO	-	Ø100	Ø86.3	2.5	340 X 209	486

Spare : Refer NA Cat# 10766-762 / EU Cat# 216-1745 for Spare Stopcock / Spigot / Tap set
Refer NA Cat# 10820-436 / EU Cat# 216-1841 for Spare Cap.

Material

Carboy : Refer above table for material type, 100 % Virgin, FDA 21 CFR Compliance, USP Class VI passed.

Closure : Polypropylene, 100 % Virgin, FDA 21 CFR Compliance, USP Class VI passed
With Food grade colour.

BSE/TSE Statement: These product is manufactured without the use of Raw materials of animal origin.

Latex Free Statement: No materials containing latex or natural rubber are used in the manufacturing, handling and packaging processes for this product.

Caution: Applicable only for PP Products

- Before Autoclaving separate the caps and Spigot / Stopcock threads.
- Do not stack products, this may cause vacuum and may lead to physical damage / collapse of product.

Material Recycle Code



Spare Spigot Set

Rectangular Carboy With Spigot



Spigot



Use internal Diameter 16mm Tube for without Connector application

Tabulation Connector



- Use internal Diameter 6mm – 8mm Tube for with Connector application.
- Connector Size towards Tap fitting 13.5mm

Stopper



Rectangular Carboy With out Spigot



Description: Carboy Rectangular

Drawing:
Not to Scale



Release: 09/26/2018

Revision: 4



NA Cat#	Europe Cat#	Description	Capacity (Liter)	Brim Capacity (Liter)	Material	Autoclavable	Cap Size (mm)	Neck Internal Dia (mm)	General wall Thickness (mm)	Product Size mm (±2MM)	Height With out Cap mm (±2MM)
75836-684	-	Jerrican Rectangular, HDPE, 20L	10 L	12.5L	HDPE	NO	Ø53	Ø39.0	2.5	351 X 245	463

Material

Carboy : Refer above table for material type, 100 % Virgin, FDA 21 CFR Compliance, USP Class VI passed.

BSE/TSE Statement: These product is manufactured without the use of Raw materials of animal origin.

Latex Free Statement: No materials containing latex or natural rubber are used in the manufacturing, handling and packaging processes for this product.

Caution: Applicable only for PP Products

1. Before Autoclaving separate the caps and Spigot / Stopcock threads / Connectors (If Applicable).
2. Do not stack products, this may cause vacuum and may lead to physical damage / collapse of product.

Material Recycle Code



Description: Jerry Can

Drawing:
Not to Scale



Release: 2/27/18


Revision: 1



Carboy - Narrow Mouth



NA Cat#	Europe Cat#	Description	Capacity (Liter)	Brim Capacity (Liter)	Material	Autoclavable	Sterile (Gamma)	Cap Size (mm)	Neck Internal Dia (mm)	General wall Thickness (mm)	Product Size mm (±2MM)	Height With out Cap mm (±2MM)
10755-121	216-1729	Carboy, Round, Narrow Mouth, 20 L, HDPE - Sterile	20 L	22	HDPE	NO	Yes	Ø83	Ø65	2.5	Ø285	518
75795-776	Not Applicable	Carboy, Round, Narrow Mouth, 50 L, HDPE - Sterile	50 L	52	HDPE	NO	Yes	Ø83	Ø65	2.5	Ø376	665

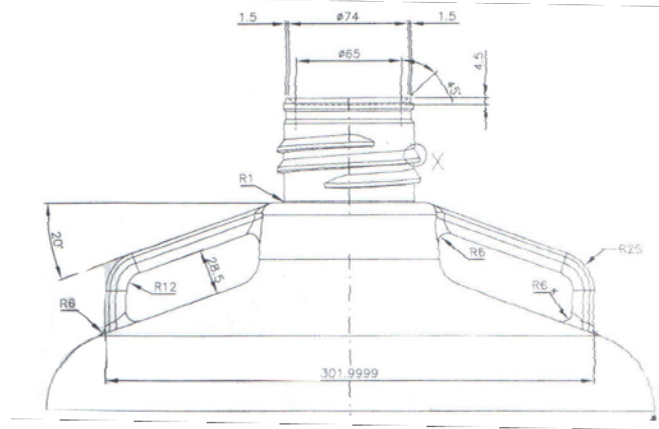

Spare : Refer NA Cat# 10820-450 / EU Cat# 216-1840 for Narrow Mouth Cap.
Refer NA Cat# 10820-438 / EU Cat# 216-1842 for Narrow Mouth Filling and Vent Cap.

Material

Carboy : Refer above table for material type, 100 % Virgin, FDA 21 CFR Compliance, USP Class VI passed.

Closure : Polypropylene, 100 % Virgin, FDA 21 CFR Compliance, USP Class VI passed With Food grade colour.

BSE/TSE Statement: These product is manufactured without the use of Raw materials of animal origin.

Latex Free Statement: No materials containing latex or natural rubber are used in the manufacturing, handling and packaging processes for this product.

Caution: Applicable only for PP Products

1. Before Autoclaving separate the caps and Spigot / Stopcock threads.
2. Do not stack products, this may cause vacuum and may lead to physical damage / collapse of product.

Material Recycle Code



Description: Carboy without Spigot - Sterile

Drawing:
Not to Scale



Release: 5/21/19

Revision: 1