Tempshield, Inc.

P.O. Box 199 Mt. Desert, ME 04660 USA

Notice of Use: **Cryo-Apron**®

Issue Date: 10/25/99 Rev.: C

Statement of Use: Cryo-Aprons® are multi-layer protective cryogenic aprons designed to provide a high level of protection to the torso and upper legs from extremely cold temperatures, both in cold atmospheres, from transmission of cold through direct contact with cold objects, and from splashes from cold liquids. The aprons are intended for use in a variety of applications and under a wide range of conditions.

- When handling cryogenic liquids, Waterproof Cryo-Gloves® should be used along with the apron.
- Select an appropriate apron length for your application.
- Periodically inspect the condition of your apron - replace aprons that show excessive wear.

Limits of Use: The maximum duration of exposure at ultra cold temperatures is dependent on many variables, including atmospheric conditions, the task being performed and the users physiology. The thermal flux transmitted through the apron is sufficiently low to allow the user adequate time to safely remove themselves from the hazard.


Care instructions: Store in a clean dry space.
Do not machine wash or dry.*




* Cryo-Aprons® were not submitted for CE wash testing.

Packaging: Cryo-Aprons® are individually packaged in poly bags.
For shipping: Individually packaged aprons shall be additionally packaged in a suitable packing carton constructed to meet all applicable freight requirements.

Cryo-Aprons® satisfy all the basic health and safety requirements of the PPE Directive. They are ergonomically designed, absent of risks and nuisance factors, contain no materials known to cause allergies or health hazards, and are comfortable to wear for extended periods.

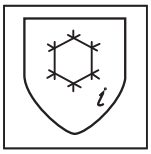

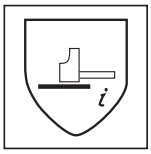

Cryo-Aprons® Available sizes and styles: EN 340

Flexibility Behavior ..... Pass
Contact Cold 2

Model	Size: Width	Length
CA36	61cm / 24”	91cm / 36”
CA42	61cm / 24”	107cm / 42”
CA48	61cm / 24”	122cm / 48”
CA54	61cm / 24”	137cm / 54”

Abrasion Resistance.....1
Tear Resistance.....2

Additional copies of this document can be obtained from Tempshield, Inc. P.O. Box 199 Mt. Desert, ME 04660 USA

CE EC type-examination was conducted by INSTITUT TEXTILE DE FRANCE
Avenue Guy de Collongue B.P. 60 69132 ECULLY CEDEX France Notified Body # 0072