
Milli-Q®
Lab Water Solutions

Milli-Q® SQ 2Series

USER MANUAL


MANUEL D'UTILISATION, MANUALE D'USO, BEDIENUNGSANLEITUNG,
MANUAL DE USUARIO, MANUAL DE INSTRUÇÕES, ユーザーマニュアル, 用户手册,
РУКОВОДСТВО ПОЛЬЗОВАТЕЛЯ

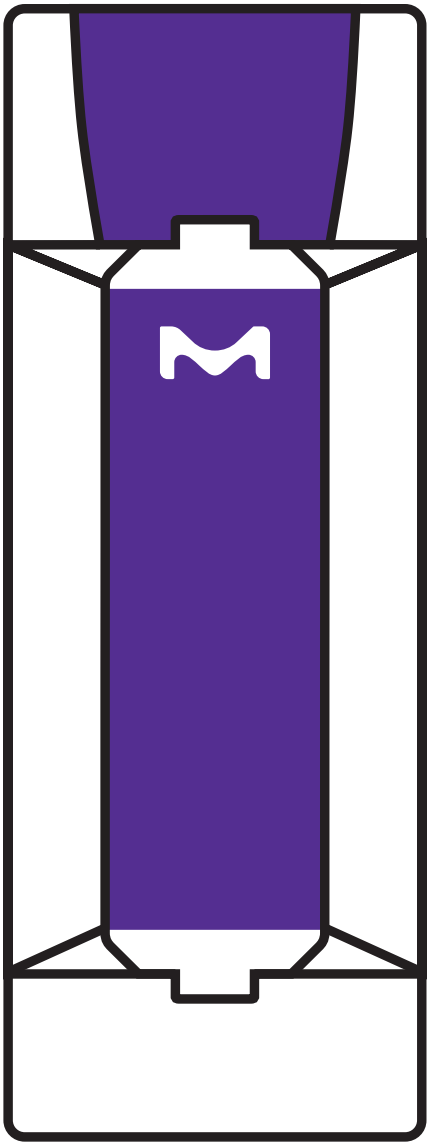
Prefilter
ZSR240PFT0

Milli-Q® SQ 2Series
Welcome to a modular universe

Bienvenue dans un univers modulaire / Benvenuto in un universo modulare / Herzlich Willkommen in einem modularen Universum / Bienvenido a un universo modular / Bem-vindo a um universo modular / モジュール式の世界にようこそ / 欢迎来到一个模块化世界 / Добро пожаловать в мир модульных решений

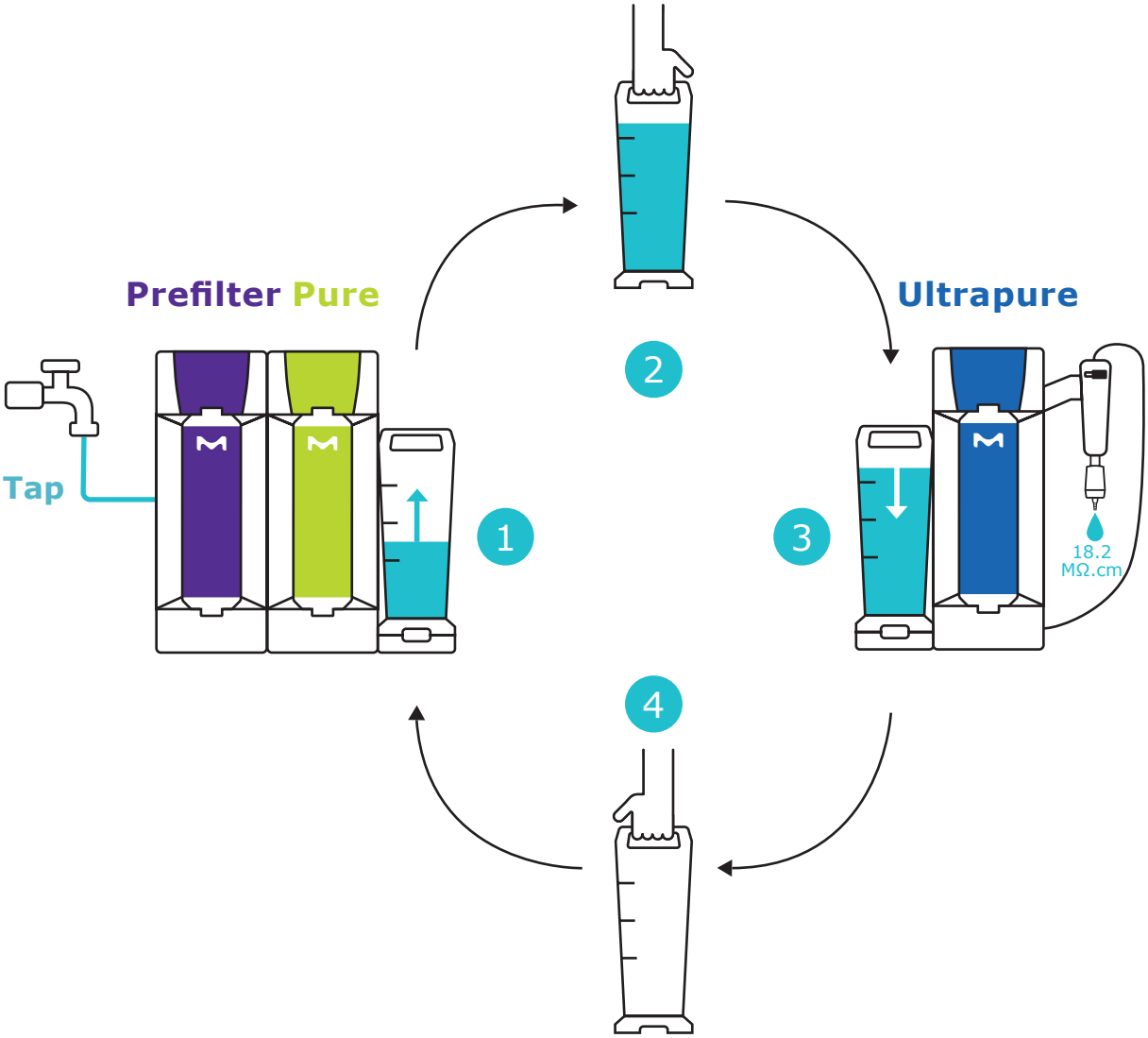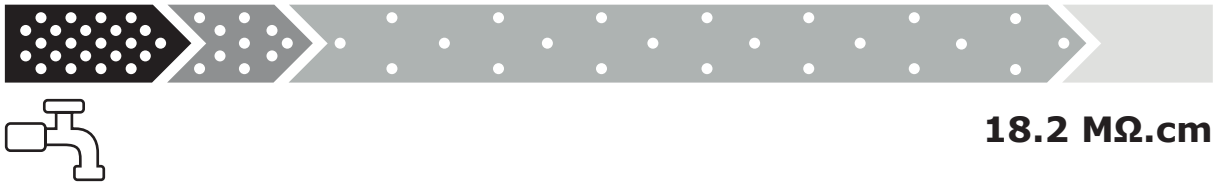
18.2 MΩ.cm

Adapt your setup to your evolving needs

Adaptez votre installation à l'évolution de vos besoins / Potrai adattare il tuo impianto all'evolversi delle tue necessità / Passen Sie das Setup an Ihre Anforderungen an / Adapte su configuración a sus necesidades en constante evolución / Adapte sua configuração às suas necessidades em evolução / 変化するニーズに設定を適合させましょう / 让您的设置适应您不断变化的需求 / Меняйте структуру в соответствии с растущими потребностями

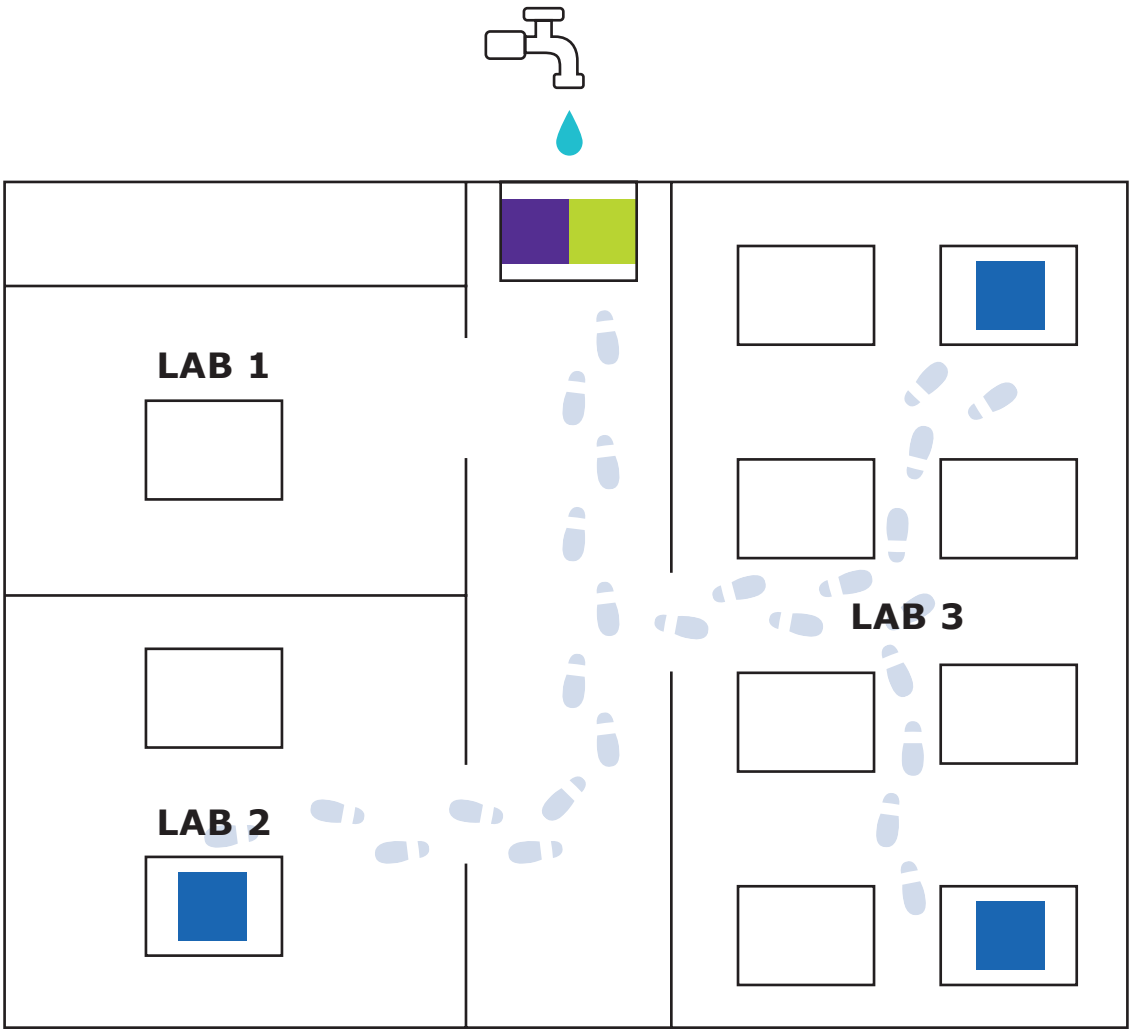
Prefilter

Pure

Ultrapure

Connection Pre-requisite & Environmental


Prérequis de connexion et environnementaux / Prerequisiti ambientali e per le connessioni / Voraussetzungen für Anschlüsse & Umfeld / Requisitos previos de conexión y entorno / Pré-requisitos de conexão e ambientais / 接続、環境要件 / 连接及环境要求 / Подсоединение и необходимые условия окружающей среды

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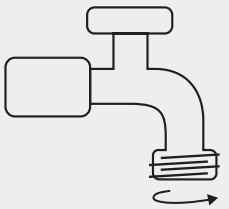
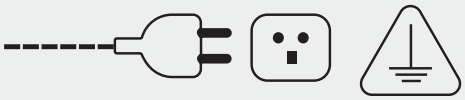
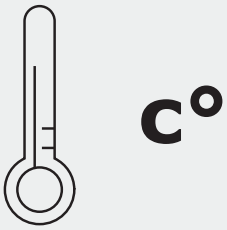
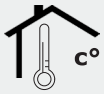
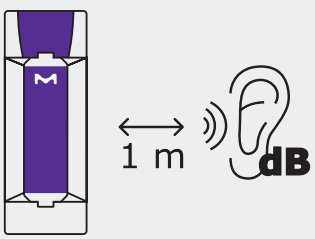
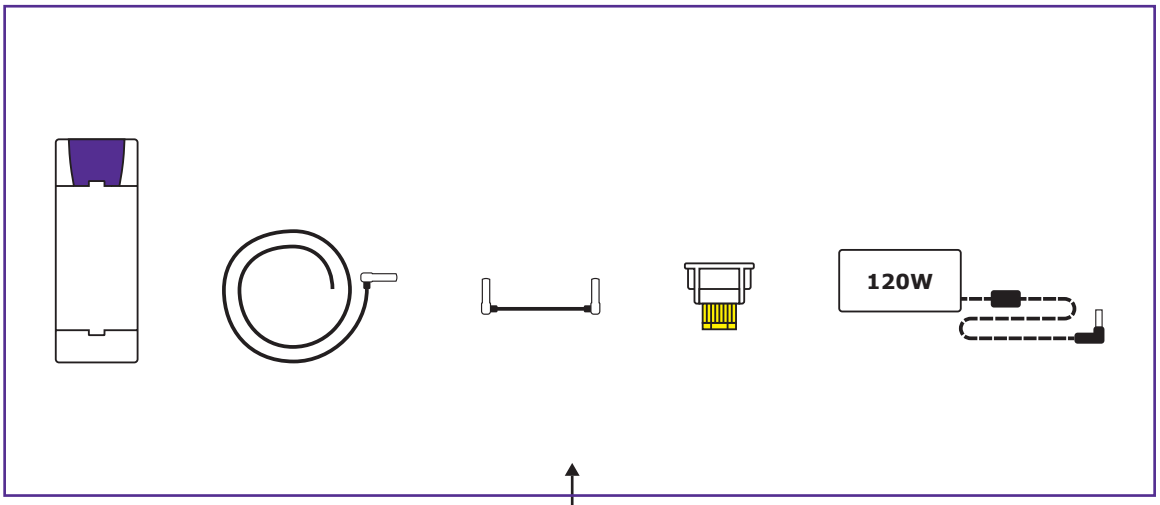
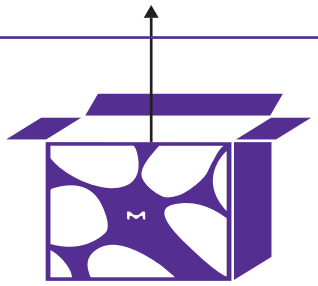
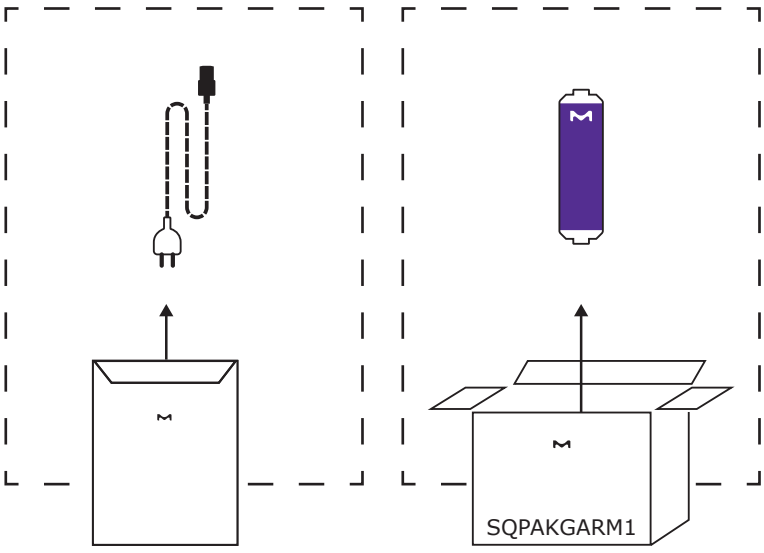
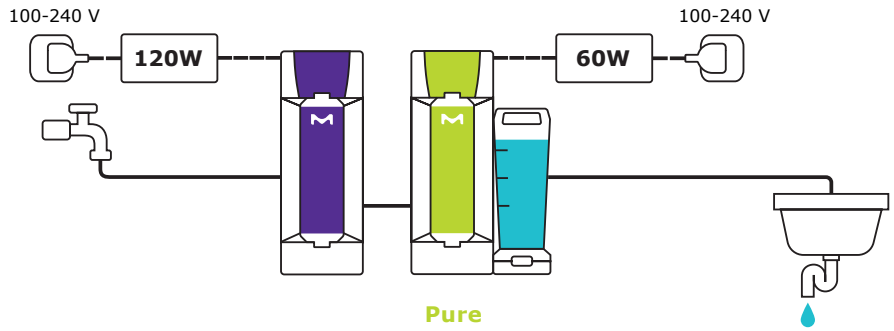
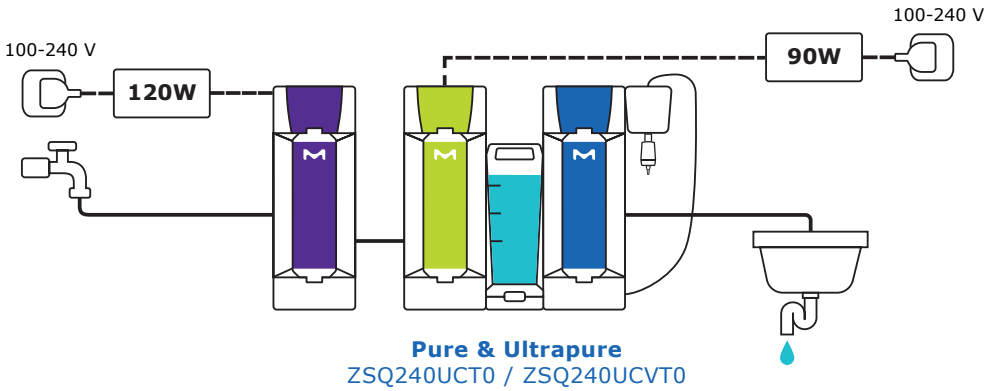
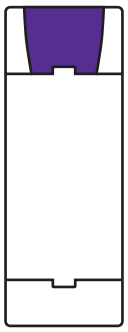
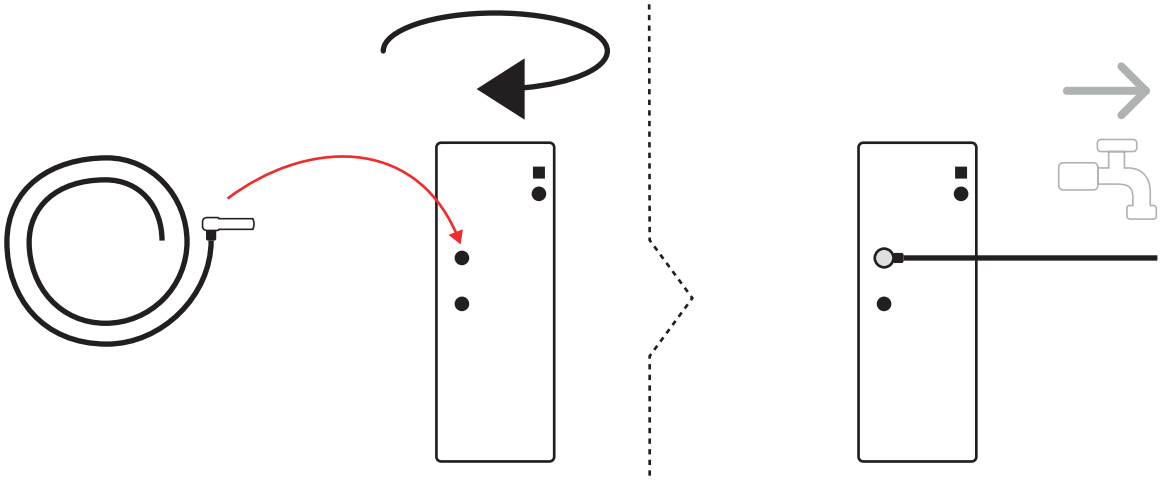
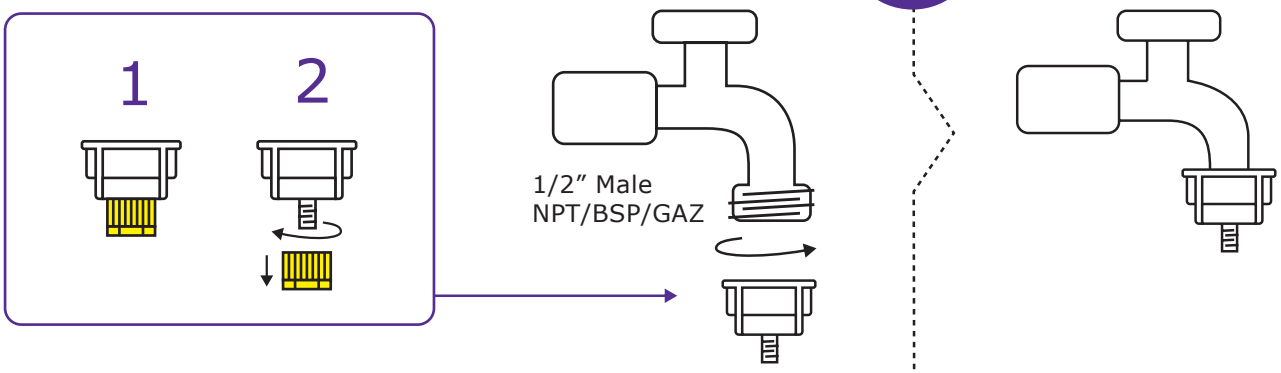
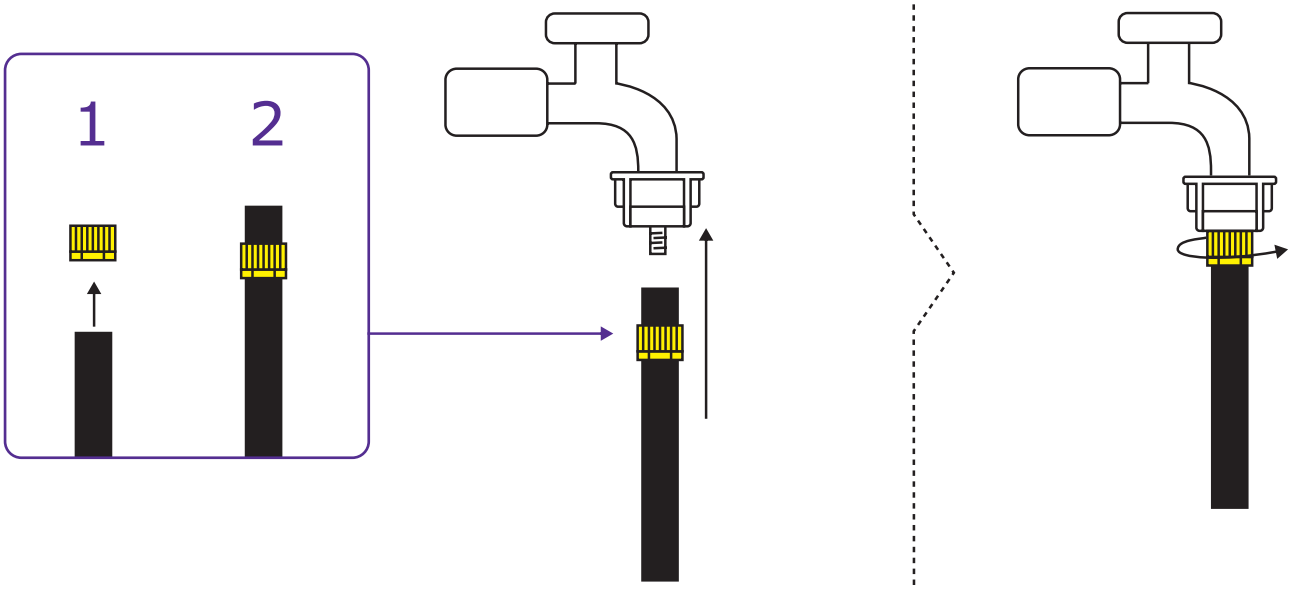
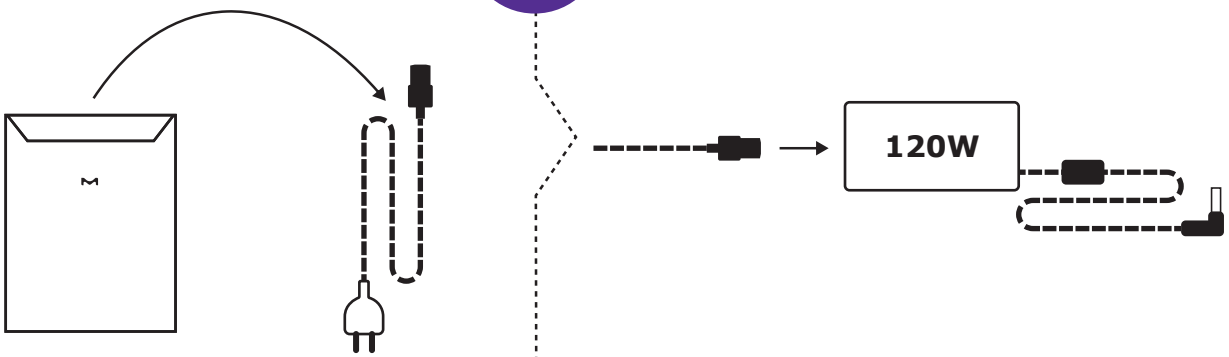
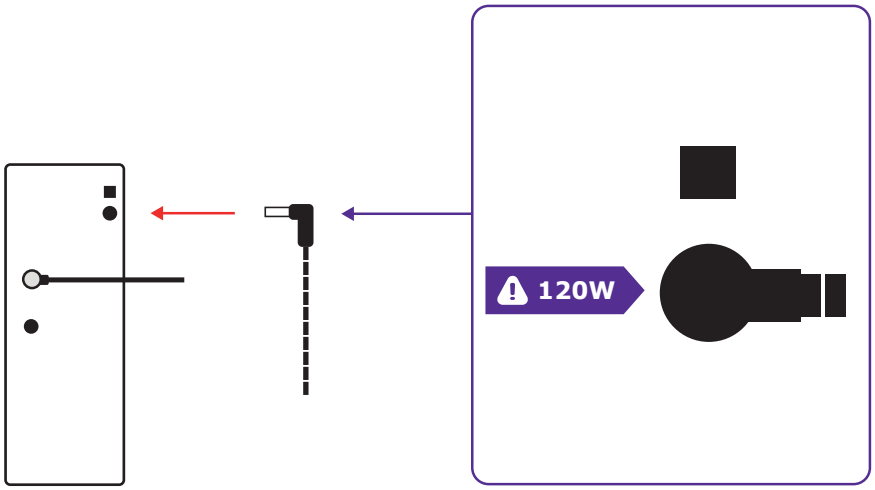
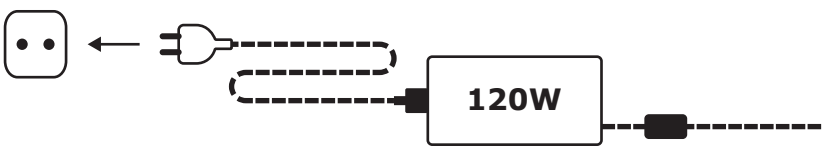
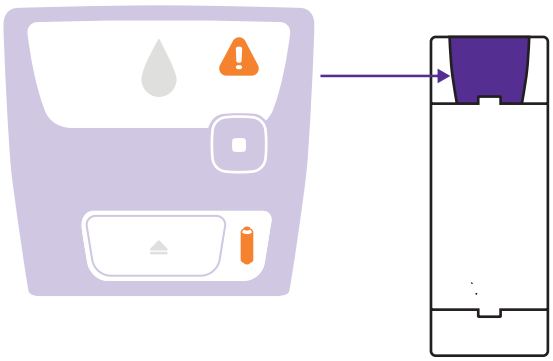
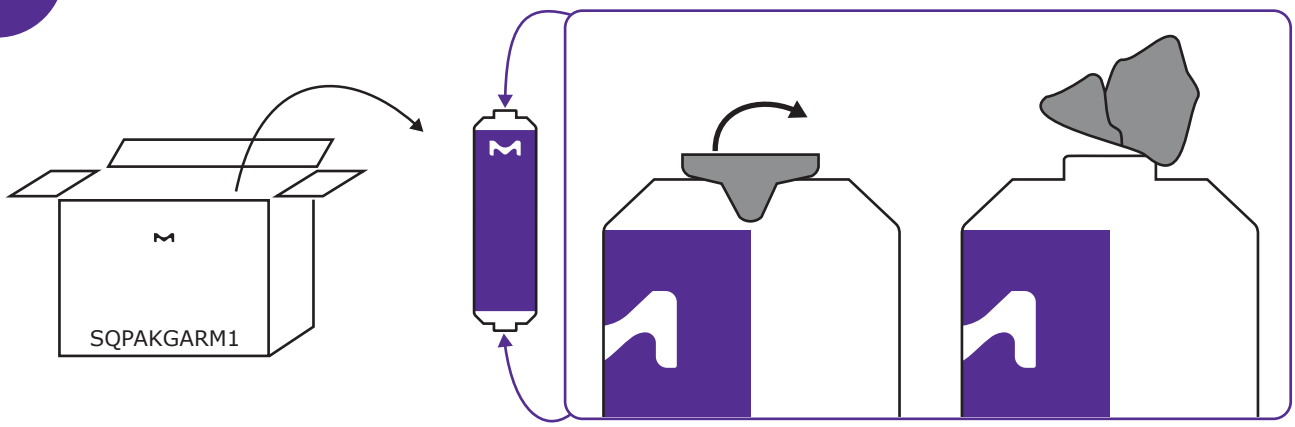
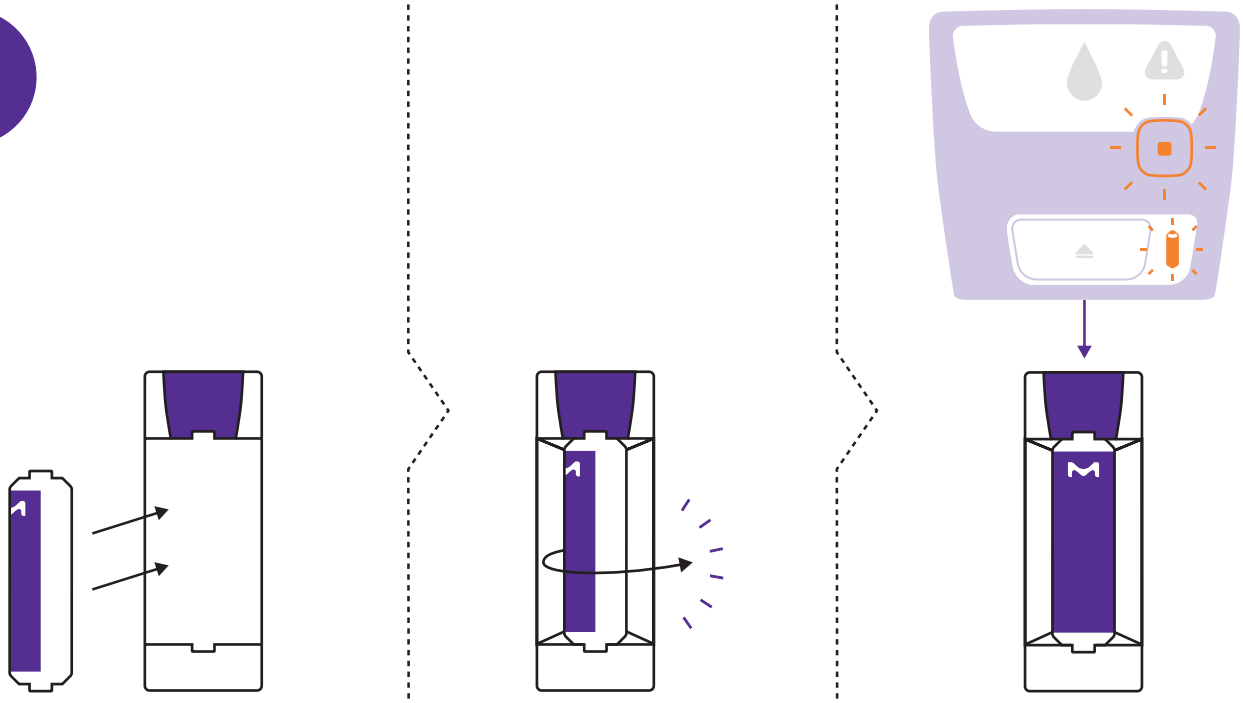
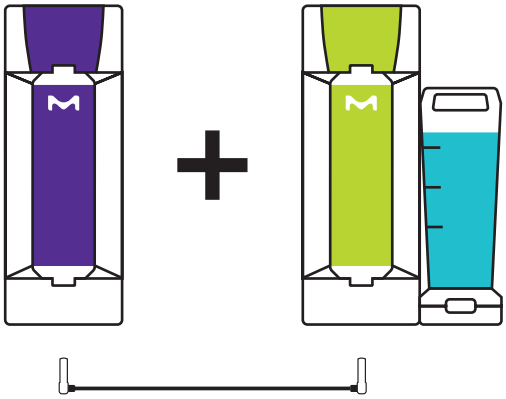
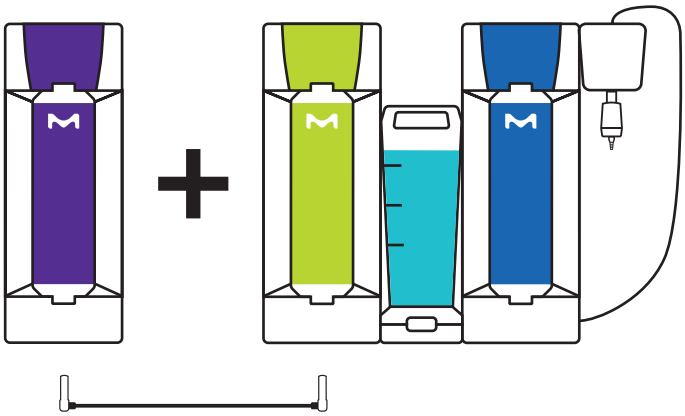
Accessory: Water Sensor Installation

Accessoire : Installation du détecteur d'eau / Accessorio: Installazione del sensore antiallagamento /
Zubehör: Installation des Wasserwächters / Accesorio: Instalación del sensor de agua /
Acessório: Instalação do sensor de água / アクセサリー:漏水センサーの設置 / 配件:漏水检测器安装 /
Дополнительный материал: Установка датчика воды


1

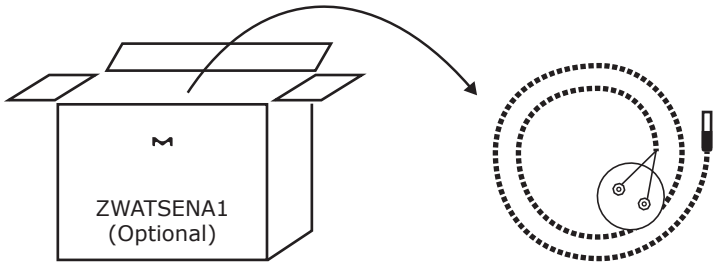2

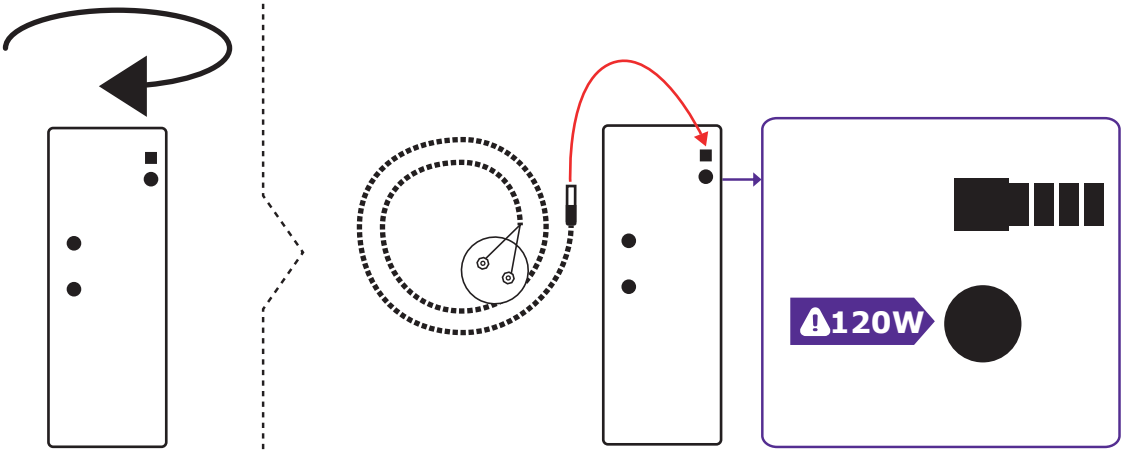3

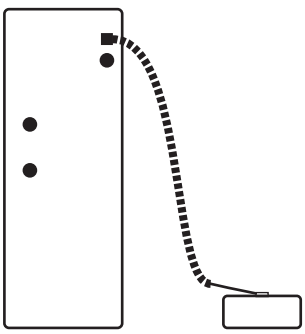Milli-Q®
Lab Water Solutions


System Presentation

PRÉSENTATION DU SYSTÈME / PRESENTAZIONE DEL SISTEMA / SYSTEMDARSTELLUNG /
PRESENTACIÓN DEL SISTEMA / APRESENTAÇÃO DO SISTEMA / システムの説明 / 系统介绍 / ПРЕЗЕНТАЦИЯ СИСТЕМЫ

SYSTEM PRESENTATION

System Overview


Aperçu du système / Panoramica del sistema / Systemüberblick / Descripción del sistema / Visão geral do sistema / システムの概要 / 系统概述 / Обзор системы


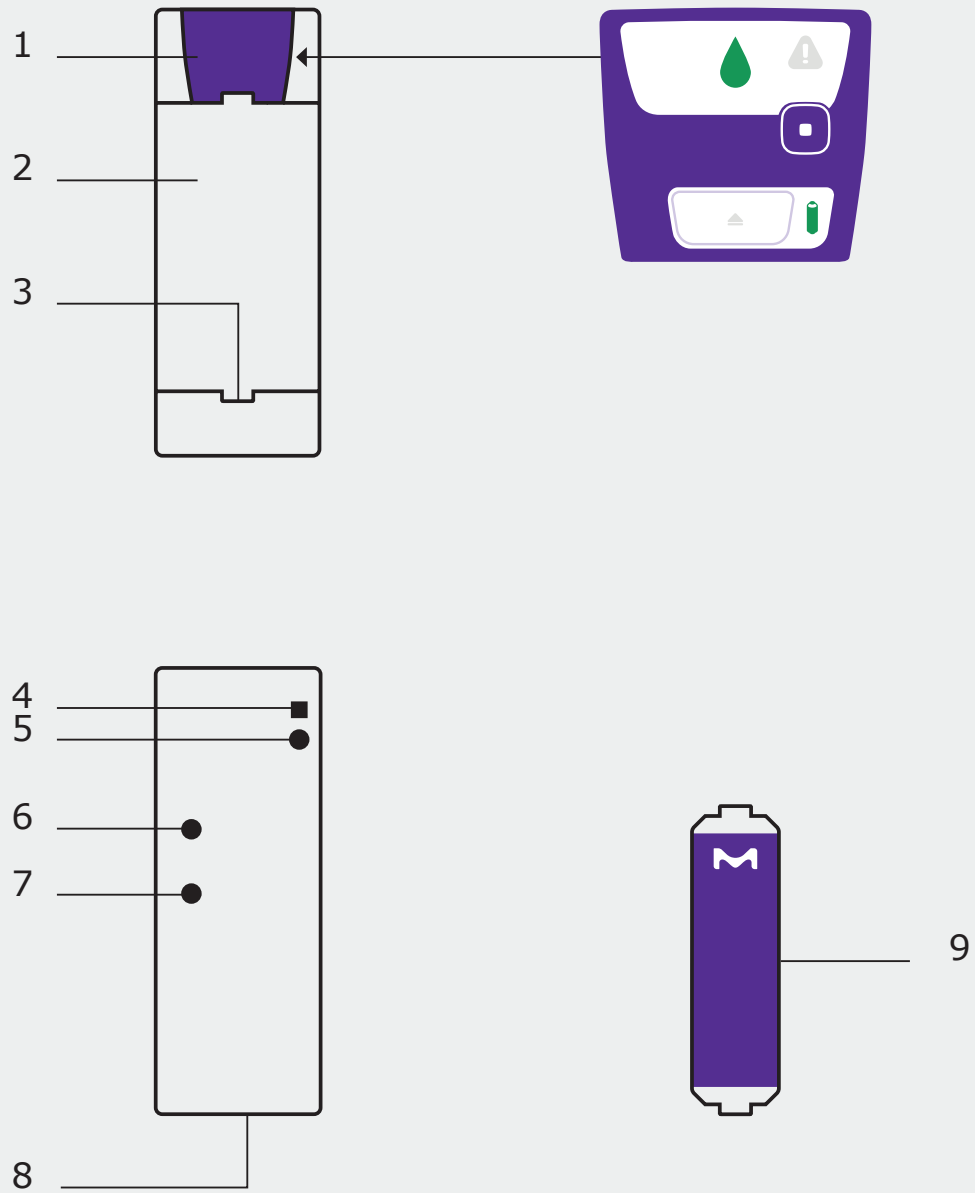
1. LED Interface	LED d'interface / LED dell'interfaccia / Schnittstellen-LED / LED del interfaz / LED de interface / インターフェイスのLED / 界面 LED 指示灯 / Светодиодные индикаторы интерфейса
2. System	Système / Sistema / System / Sistema / Sistema / システム / 系统 / Система
3. Cartridge holder	Support de cartouche / Portacartuccia / Filterelementhalter / Portacartuchos / Porta-cartucho / カートリッジホルダー / 滤柱座 / Держатель картриджа
4. Water sensor socket	Prise du détecteur d'eau / Presa per il sensore antiallagamento / Wasserwächterbuchse / Toma del sensor de agua / Soquete do sensor de água / 漏水センサーソケット / 漏水检测器插口 / Гнездо датчика воды
5. Power supply socket	Prise de l'alimentation électrique / Presa per l'alimentatore / Stromversorgungsbuchse / Toma de alimentación / Soquete da fonte de alimentação / 電源ソケット / 电源插座 / Гнездо источника питания
6. Feed water connection	Raccordement à l'eau d'alimentation / Connessione con l'alimentazione dell'acqua / Speisewasseranschluss / Conexión al agua de alimentación / Conexão de água de alimentação / 供給水接続口 / 进水连接 / Соединение для подачи воды
7. Prefilter to Pure connection	Raccordement Préfiltre-Module pur / Prefiltro - Connessione con il modulo per la produzione di acqua pura / Vorfilter – Anschluss des Reinwassermoduls / Prefiltro: conexión módulo Pure / Pré-filtro - conexão de módulo puro / プレフィルター - Pureモジュールの接続 / 预过滤器至纯化模块的连接 / Фильтр предварительной очистки — подсоединение модуля Pure
8. Firm label	Plaque de firme / Piastra fissa / Feste Platte / Placa de marca / Placa firme / 銘板 / 铭牌(Firm plate) / Фирменная табличка
9. Prefilter Cartridge	Préfiltre / Cartuccia prefiltro / Filterelement mit Vorfilter / Cartucho de prefiltro / Cartucho de pré-filtro / プレフィルターカートリッジ / 预过滤器滤柱 / Картридж с фильтром предварительной очистки

Dimensions & Weights

Dimensions et poids / Dimensioni e pesi / Abmessungen & Gewichte / Dimensiones y pesos / Dimensões e pesos / 寸法と重量 / 尺寸和重量 / Размеры и вес

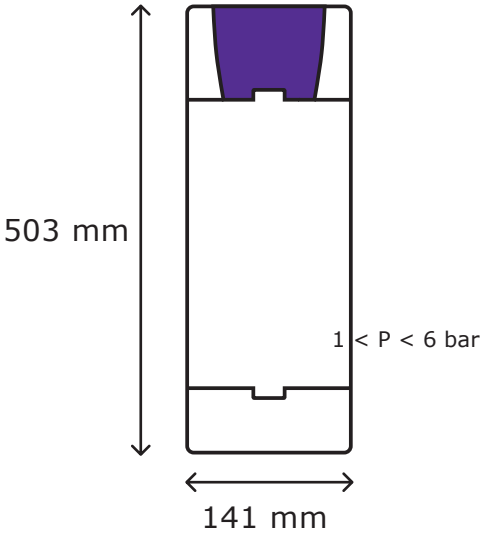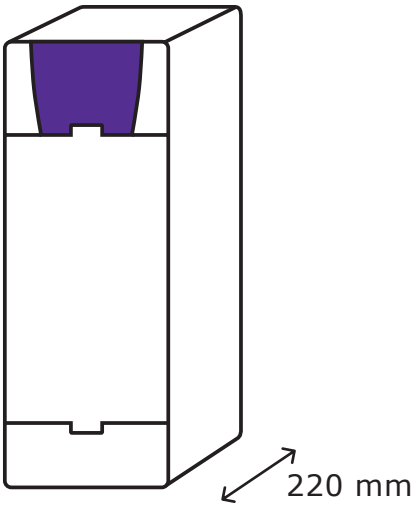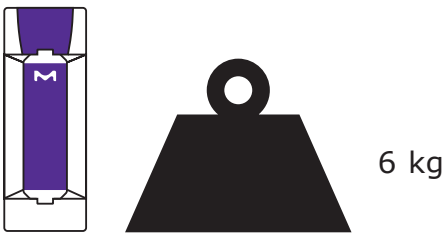
Connection Pre-requisite & Environmental

Prérequis de connexion et environnementaux / Prerequisiti ambientali e per le connessioni / Voraussetzungen für Anschlüsse & Umfeld / Requisitos previos de conexión y entorno / Pré-requisitos de conexão e ambientais / 接続、環境要件 / 连接及环境要求 / Подсоединение и необходимые условия окружающей среды

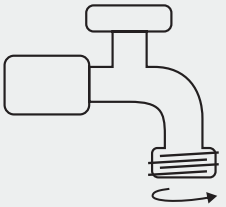
Q_{max} ≥ 120 L.h^{-1} @ 1 bar


1 < P < 6 bar

1/2” Male NPT/BSP/GAZ

50/60 Hz 10%	120 W	100 - 240 V 10%

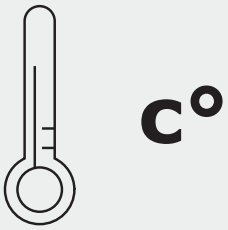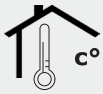4°C ≤ T ≤ 40°C

5°C ≤ T ≤ 35°C

Air and water / Air et eau / Aria e acqua / Luft und Wasser / Aire y agua / Ar e água / 空気と水 / 空气和水 / Воздух и вода


< 80% @ 31°C,
< 50 % @ 40°C

Relative humidity / Humidité relative / Umidità relativa / Relative Luftfeuchtigkeit / Humedad relativa / Umidade relativa / 相对湿度 / 相对湿度 / Относительная влажность

< 3000 m


Above sea level / Au-dessus du niveau de la mer / Sopra il livello del mare / Über dem Meeresspiegel / Sobre el nivel del mar / Acima do nível do mar / 海面以上 / 高于海平面 / Высота над уровнем моря
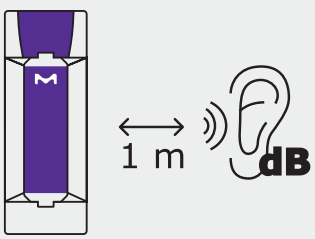
< 60 dB @ 1 m @ 20 kPa

Technical Specifications


Spécifications techniques / Specifiche tecniche / Technische Daten / Especificaciones técnicas /
Especificações técnicas / 技術仕様 / 技术规格 / Технические характеристики

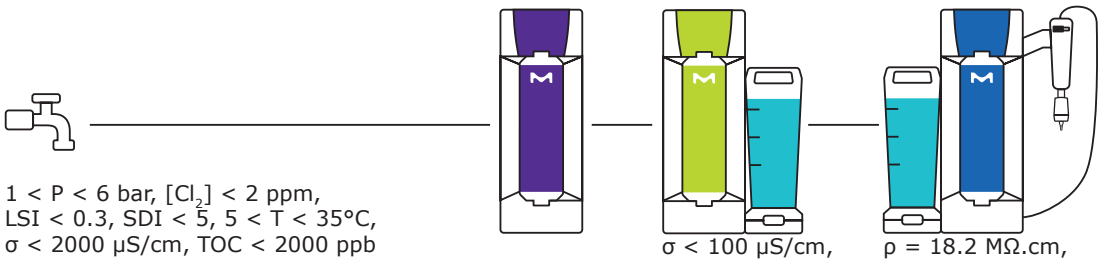1 < P < 6 bar, [Cl_2] < 2 ppm,
LSI < 0.3, SDI < 5, 5 < T < 35°C,
σ < 2000 µS/cm, TOC < 2000 ppb

σ < 100 µS/cm,
30 ≤ Q ≤ 40 L/h


ρ = 18.2 MΩ.cm,
TOC ≤ 5 ppb*

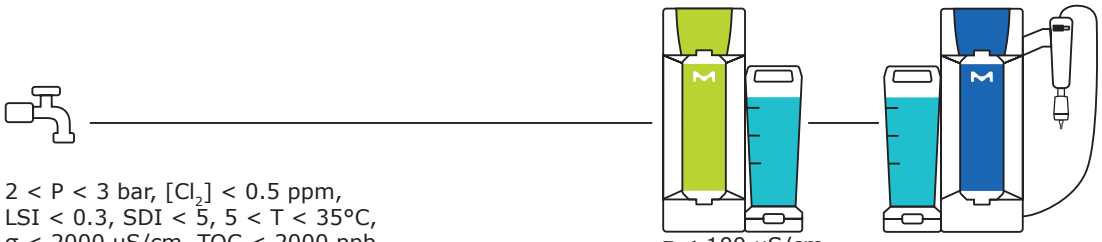2 < P < 3 bar, [Cl_2] < 0.5 ppm,
LSI < 0.3, SDI < 5, 5 < T < 35°C,
σ < 2000 µS/cm, TOC < 2000 ppb

σ < 100 µS/cm,
5 ≤ Q ≤ 40L/h

ρ = 18.2 MΩ.cm

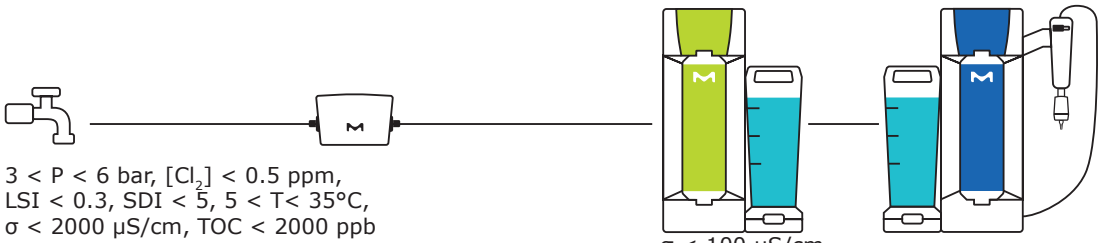3 < P < 6 bar, [Cl_2] < 0.5 ppm,
LSI < 0.3, SDI < 5, 5 < T< 35°C,
σ < 2000 µS/cm, TOC < 2000 ppb

σ < 100 µS/cm,
5 ≤ Q ≤ 40 L/h

ρ = 18.2 MΩ.cm

All parameters corrected @ 25°C /
Tous les paramètres corrigés à 25 °C /
Tutti i parametri corretti a 25 °C /
Alle Parameter korrigiert @ 25 °C /
Todos los parámetros corregidos a 25 °C /
Todos os parâmetros corrigidos a 25 °C /
全パラメーターは25 °Cに補正 /
所有参数均在 25 °C 下校正 /
Все параметры скорректированы при температуре 25 °C

SQPAK TOC QUANTA Cartridge TOC ≤ 5 ppb*


SQPAK™ Final Filter

Q ≤ 1.6 L/min

Particulates Particules / Particelle / Partikel / Partículas / Particulados / 粒子/ 颗粒 / Частицы	≤ 0.22 µm
Bacteria Bactéries / Batteri / Bakterien / Bacterias / Bactérias / 生菌数/ 细菌 / Бактерии	≤ 10 CFU/L

Organization Organisme / Ente / Organisation / Organización / Organização / 組織 / 组织 / Организация	Water quality / grade Qualité/Grade de l'eau / Qualità/grado dell'acqua / Wasserqualität / Calidad o grado del agua / Qualidade/grau da água / 水質グレード / 水质/级别 / Качество/класс воды
European Pharmacopeia	Purified water in bulk / Eau purifiée en vrac / Acqua purificata in bulk / Gereinigtes Wasser in Großmengen / Agua purificada a granel / Água purificada a granel / バルク精製水 / 散装纯水 / Очищенная вода нерасфасованная
U.S. Pharmacopeia	
Japanese Pharmacopoeia / 日本薬局方	Purified water / Eau purifiée / Acqua purificata / Gereinigtes Wasser / Agua purificada / Água purificada / 精製水 / 纯水 / Очищенная вода
Chinese Pharmacopoeia / 中国药典	
ASTM® D1193	Type I water / Eau de Type I / Acqua di Tipo I / Wasser vom Typ I / Agua de tipo I / Água tipo I / タイプ Ⅰ 水 / Ⅰ类水 / Вода типа I
ISO® 3696	Grade 1 water / Eau de classe 1 / Acqua di Grado 1 / Wasser der Qualität 1 / Agua de grado 1 / Água de grau 1 / グレード1水 / 1 级水 / Вода класса 1
Chinese National Standard GB/T 6682 / 中国国家标准 GB/T 6682	
Chinese National Standard GB/T 33087 / 中国国家标准 GB/T 33087	Ultrapure water / Eau ultra pure / Acqua ultrapura / Reinstwasser / Agua ultrapura / Água ultrapura / 超純水 / 超纯水 / Сверхчистая вода
JIS K 0557	A4 water / Eau A4 / Acqua A4 / A4-Wasser / Agua A4 / Água A4 / A4水 / A4 水 / Вода A4
Clinical and Laboratory Standards Institute (CLSI®)	Clinical Laboratory Reagent Water (CLRW) / Eau de qualité CLRW / Acqua a grado reagente per il laboratorio clinico CLRW / Reinwasser für das klinische Labor CLRW / Agua de calidad reactiva para laboratorios clínicos CLRW / Água reagente para laboratórios clínicos CLRW / Clinical Laboratory Reagent Water (CLRW) / 临床实验室试剂水(CLRW) / Вода клинического лабораторного назначения CLRW

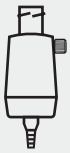SQPAK™ Bio Final Filter

Q ≤ 0.7 L/min

Bacteria Bactéries / Batteri / Bakterien / Bacterias / Bactérias / 生菌数 / 细菌 / Бактерии	< 10 CFU/L
Pyrogens (endotoxins) Pyrogènes (endotoxines) / Pirogeni (endotossine) / Pyrogene (Endotoxine) / Pirógenos (endotoxinas) / Pirogênios (endotoxinas) / パイロジェン（エンドトキシン） / 热原（内毒素） / Пирогены (эндотоксины)	< 0.001 EU/mL
RNases RNases / RNasi / RNasen / ARNasas / RNases / RNase / 核糖核酸酶（RNase） / РНКазы	< 1 pg/mL
DNases DNases / DNasi / DNasen / ADNasas / DNases / DNases / 脱氧核糖核酸酶（DNases） / ДНКазы	< 5 pg/mL
Proteases Protéases / Proteasi / Proteasen / Proteasas / Proteases / プロテアーゼ/ 蛋白酶 / Протеазы	< 0.15 µg/mL

Milli-Q®
Lab Water Solutions


User Manual


MANUEL D'UTILISATION / MANUALE D'USO / BEDIENUNGSANLEITUNG / MANUAL DE USUARIO /
MANUAL DO USUÁRIO / ユーザーマニュアル / 用户手册 / РУКОВОДСТВО ПОЛЬЗОВАТЕЛЯ


USER MANUAL

How To Use

Mode d’emploi / Istruzioni per l’uso / Anwendung / Cómo utilizar / Como usar /
ご利用方法 / 如何使用 / Способ применения

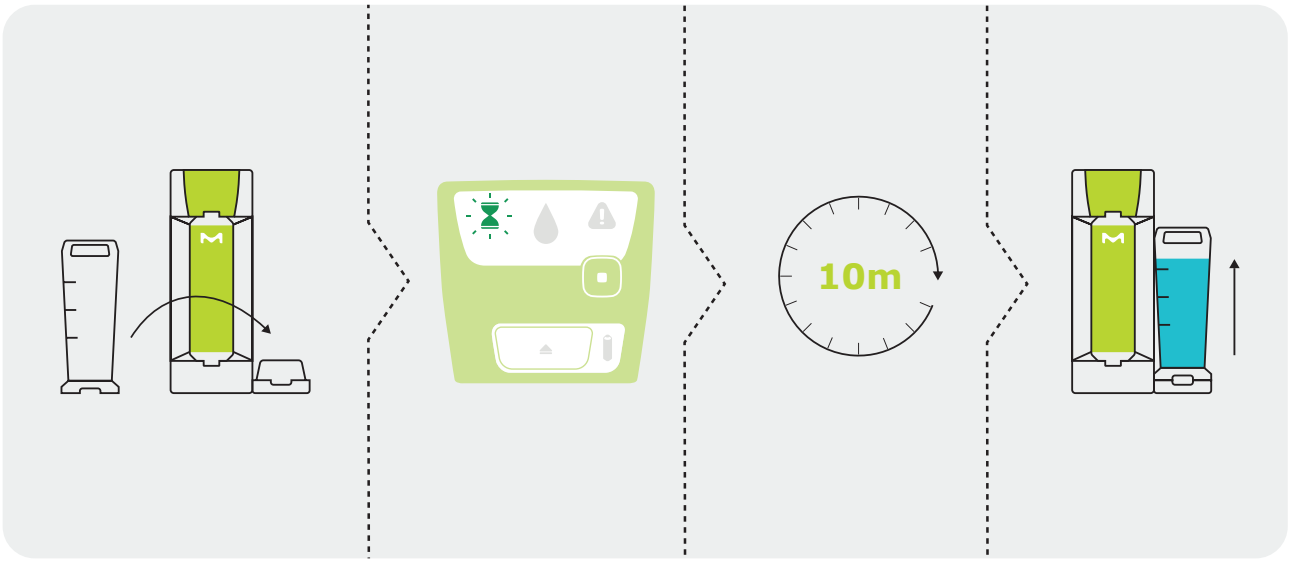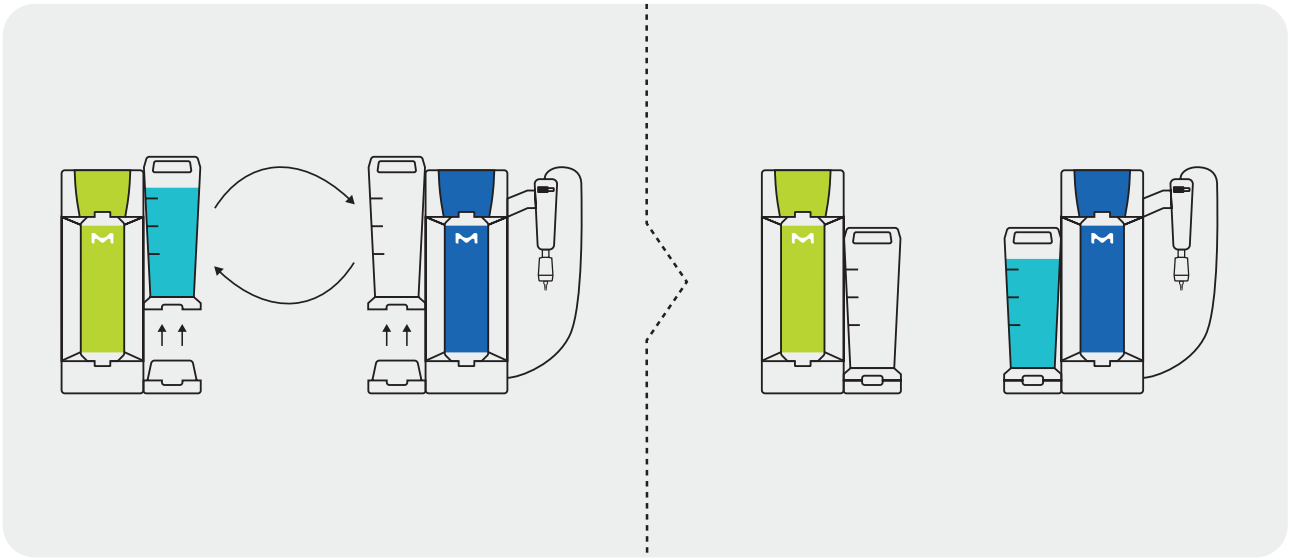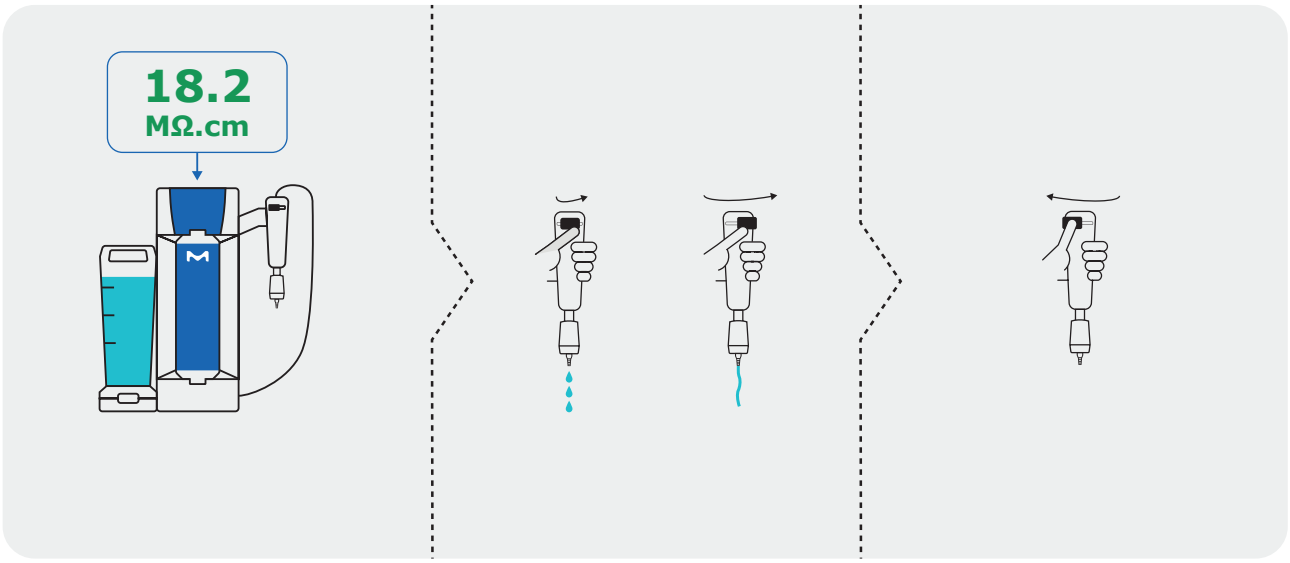
Signals & Messages


Signaux et messages / Segnali e messaggi / Signale & Mitteilungen / Señales y mensajes / Sinais e mensagens / シグナルとメッセージ / 信号和消息 / Сигналы и сообщения

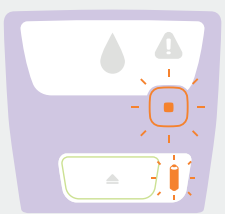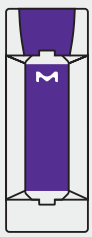
Troubleshooting


Dépannage / Risoluzione dei problemi / Fehlersuche / Resolución de problemas / Resolução de problemas / トラブルシューティング / 故障排除 / Устранение неполадок

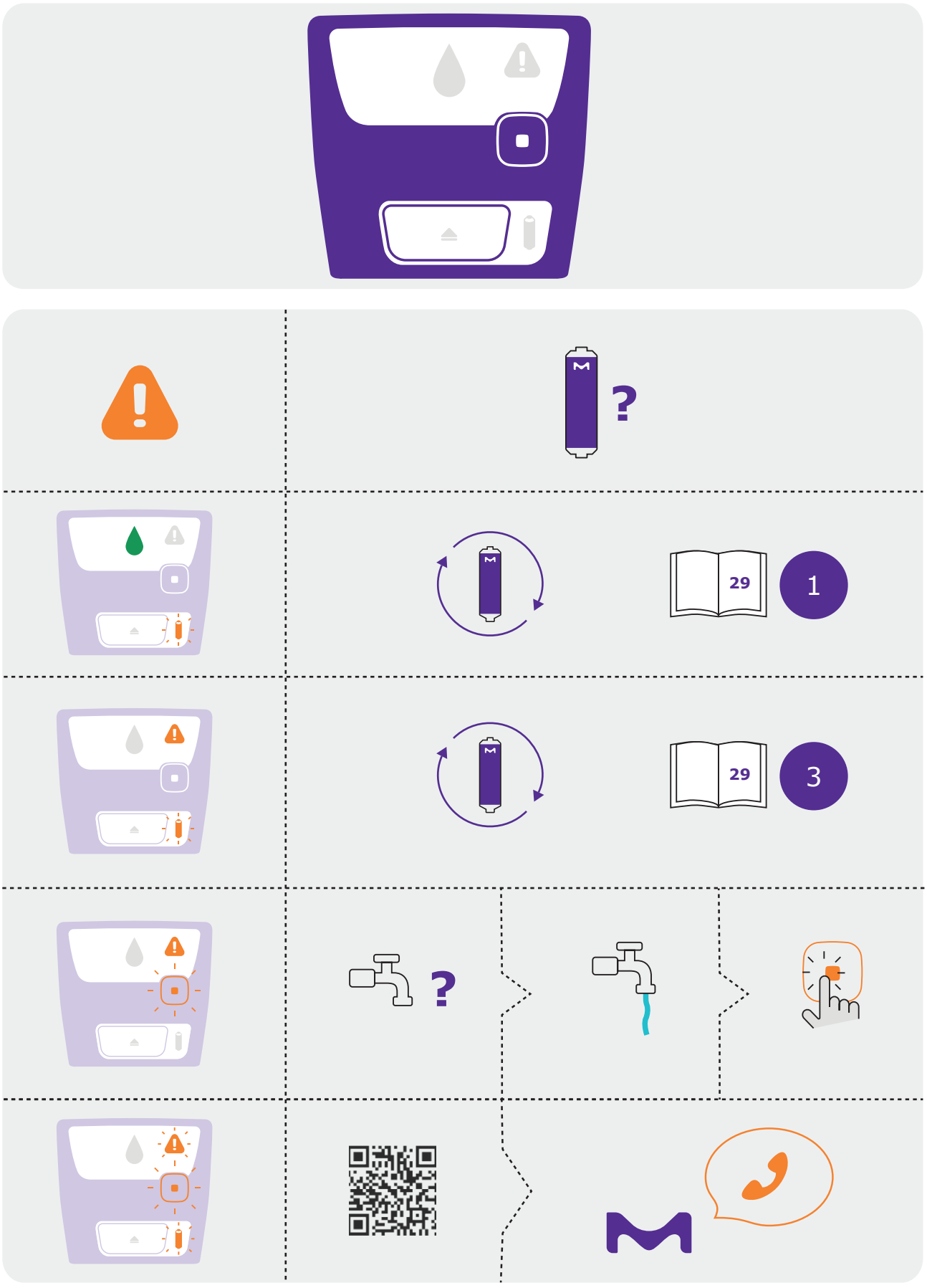
Cartridge Replacement

Remplacement des cartouches / Sostituzione delle cartucce / Ersetzen von Filterelementen / Sustitución del cartucho / Substituição do cartucho / カートリッジの交換 / 滤柱更换 / Замена картриджа

1

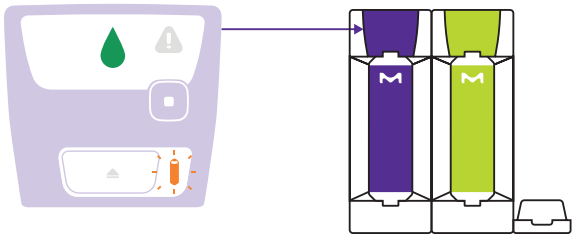
2


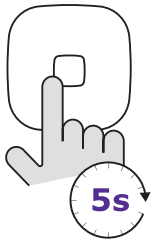
3

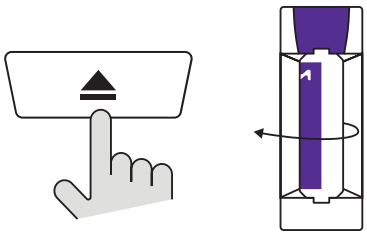
4

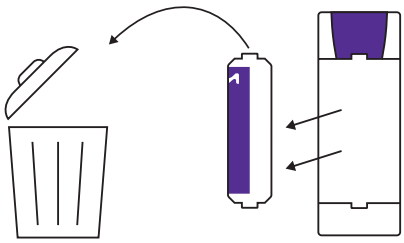
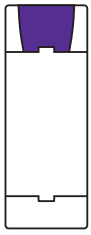
5

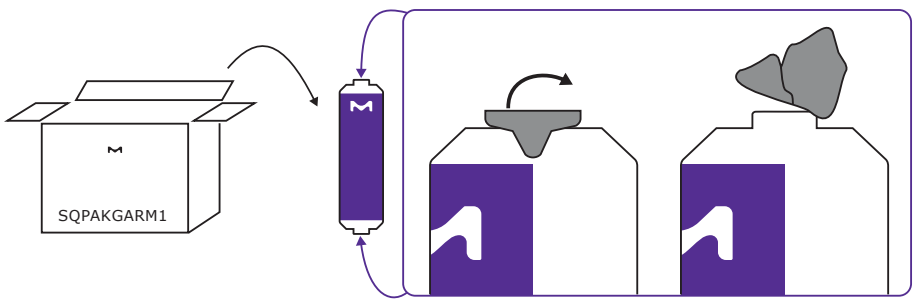
6

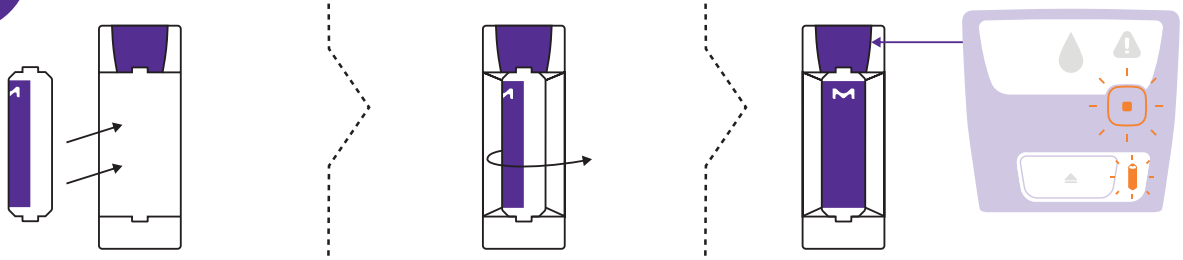7

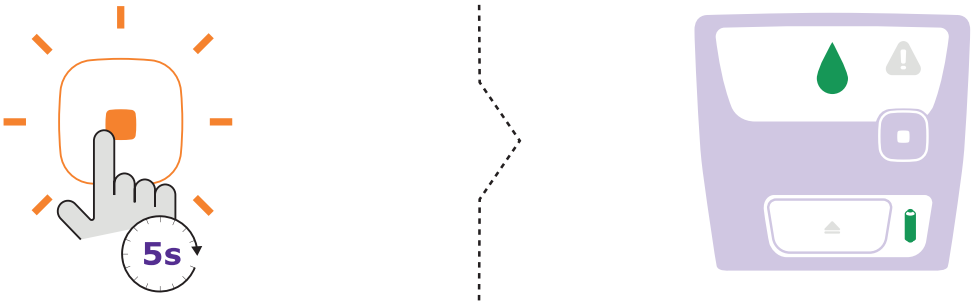8

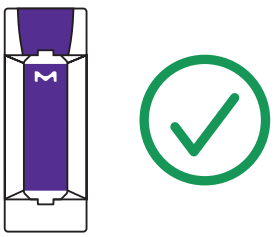Milli-Q®
Lab Water Solutions

Ordering Information

GUIDE D'ACHAT / INFORMAZIONI PER GLI ORDINI / BESTELLINFORMATIONEN / INFORMACIÓN PARA PEDIDOS / INFORMAÇÕES PARA PEDIDOS / ご注文情報 / 订货信息 / ИНФОРМАЦИЯ ДЛЯ ЗАКАЗА


ORDERING INFORMATION

Ordering Information


Guide d'achat / Informazioni per gli ordini / Bestellinformationen / Información para pedidos / Informações para pedidos / ご注文情報 / 订货信息 / Информация для заказа


ZSR240PFT0
Milli-Q® SP 240 Prefilter Module

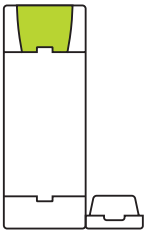ZSR240R0T0
Milli-Q® SR 240 Purification System

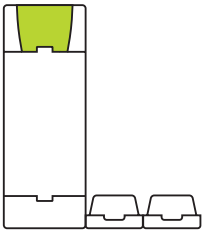
ZSR240RLT0
Milli-Q® SR 240L Purification System

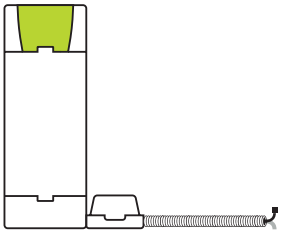
ZSR240RXT0
Milli-Q® SR 240XL Purification System


TANKSQ003
Milli-Q® SQ Switch 3.5L Tank

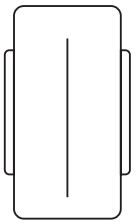
TANKSQ050
Milli-Q® SQ XL 50L Tank

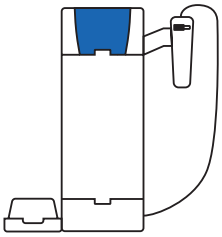
ZSQ200U0T0
Milli-Q® SQ 200 Purification System

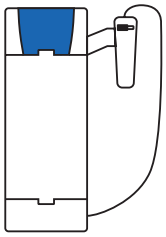ZSQ200UPT0
Milli-Q® SQ 200P Purification System

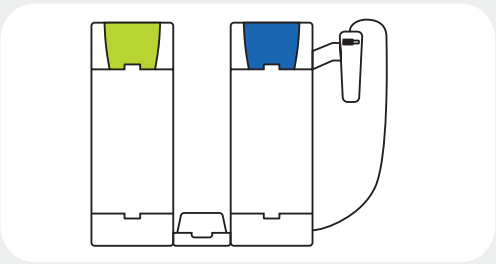ZSQ240UCT0
Milli-Q® SQ 240C Purification System

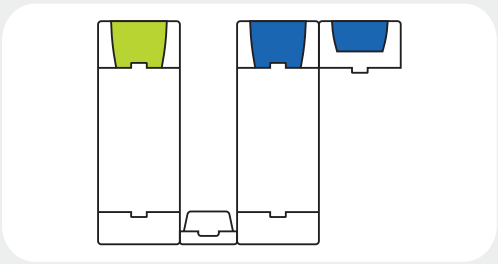ZSQ240UCVT0
Milli-Q® SQ 240CV Purification System

SQPAKGARM1
SQPAK™ Gard Pack

SQPAKPR0M1
SQPAK™ RO Cartridge


SQPAKQUAM1
SQPAK™ Quanta Cartridge


SQPAKT0CM1
SQPAK™ TOC Quanta Cartridge


SQPAKFINM1
SQPAK™ Final Filter

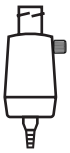
SQPAKBI0M1
SQPAK™ Bio Final Filter

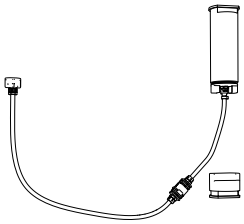
ZSQ240CTT0
Milli-Q® SQ Cleaning Tool

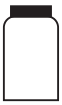
ZWCL01F50
ROProtect C

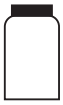
5874316024
EfferSan Effervescent Tablets (USA only) /
Pastilles effervescentes EfferSan (USA uniquement) /
Compresse effervescenti EfferSan (USA) /
EfferSan-Brausetabletten (nur USA) /
Comprimidos efervescentes EfferSan (EE.UU. solo) /
Tabletes efervescentes EfferSan (EUA apenas) /
EfferSan 発泡錠（米国のみ） /
EfferSan 泡腾片（仅在美国）) /
Пенящиеся таблетки EfferSan (только в США)

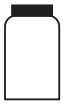5874316024C
EfferSan Effervescent Tablets (CA only) /
Pastilles effervescentes EfferSan (CAN uniquement) /
Compresse effervescenti EfferSan (Canada) /
EfferSan-Brausetabletten (nur CA) /
Comprimidos efervescentes EfferSan (CAN solo) /
Tabletes efervescentes EfferSan (Canadá apenas) /
EfferSan 発泡錠（カナダのみ） /
EfferSan 泡腾片（仅在加拿大）/
Пенящиеся таблетки EfferSan (только в Канаде)

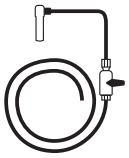ZSQ240GVT0
Milli-Q® SQ Gravity Valve

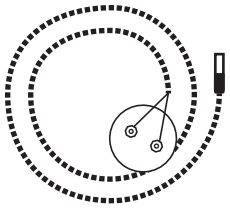ZWATSENA1
Water sensor

ZFMQPR3BT0
Milli-Q® Pressure Regulator

SQPAKACARE
SQPAK™ Care A Cleaner

Legal Notices

Mentions légales / Note legali / Rechtliche Hinweise / Avisos legales / Avisos Legais /
法的記載事項 / 法律声明 / Правовые уведомления

Legal information & warranty
It has always been Millipore SAS policy to continuously improve its products. The information in this document is subject to change without notice and should not be construed as a commitment by Millipore SAS. Millipore SAS assumes no responsibility for any errors that might appear in this document. This user manual is believed to be complete and accurate at the time of publication. In no event shall Millipore SAS be liable for incidental or consequential damages in connection with or arising from the use of this user manual.

Product warranty and limitation of liability
The applicable warranty and limitation of liability for the products listed in this publication may be found at www.SigmaAldrich.com within the "Conditions of Sale" applicable to your purchase transaction.

Copyright
© Millipore SAS 2023.

All rights reserved. This document or parts thereof may not be reproduced in any form without the written permission of Millipore SAS. The photographs illustrating the products are non-contractual.

Trademarks
New trademark names

The initial M, Millipore, Milli-Q, Q-POD, SQPAK, Millipak (…) and Biopak, are trademarks of Merck KGaA, Darmstadt, Germany.Millipore SAS is an affiliate of Merck KGaA, Darmstadt, Germany.All other trademarks are trademarks of their respective manufacturers.The Life Science Business of Merck KGaA, Darmstadt, Germany operates as MilliporeSigma in the US and Canada.

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All rights reserved.この文書はその全体または一部分でも、Millipore SASの書面による許可無くしてはいかなる形体でも複製を禁じます。製品を図示する写真は契約の一部を成しません。