# INSTRUCTION MANUAL

## VWR® Hotplate W10

**European Catalogue Number:**
720-2422

**Instructions pour l'utilisation (FR) 6-9 (27)**
**Gebrauchsanweisung (DE) 10-13**
**Istruzioni per l'uso (IT) 14-17**
**Instrucciones para el uso (ES) 18-21**
**Instruções para el uso (PT) 22-25**

Version: 1
Issued: 23122014

CE IVD

Please follow the guidelines below, and read this manual in its entirety to ensure safe operation of the unit.

The unit can have hot surfaces and hot water

Do not use the unit in an explosive environment

## Avoiding Electrical Shock

Make sure that the voltage rating on the unit corresponds to the mains power voltage

Make sure the mains power socket (wall socket) has a functioning electrical ground connection

Make sure the units power plug corresponds to the mains power socket

Always use the units grounded power cord between the unit and mains power socket

Do not use the unit if the power cord is damaged or if the unit is damage

Do not open the unit, all repairs should be done by an authorized technician

---

## Compliance with local laws and regulations

The customer is responsible for applying for and obtaining the necessary regulatory approvals or other authorisations necessary to run or use the Product in its local environment. VWR will not be held liable for any related omission or for not obtaining the required approval or authorisation, unless any refusal is due to a defect of the product.

## Equipment disposal

This equipment is marked with the crossed out wheeled bin symbol to indicate that this equipment must not be disposed of with unsorted waste.

Instead it's your responsibility to correctly dispose of your equipment at lifecycle -end by handling it over to an authorized facility for separate collection and recycling. It's also your responsibility to decontaminate the equipment in case of biological, chemical and/or radiological contamination, so as to protect from health hazards the persons involved in the disposal and recycling of the equipment.

For more information about where you can drop off your waste of equipment, please contact your local dealer from whom you originally purchased this equipment.

By doing so, you will help to conserve natural and environmental resources and you will ensure that your equipment is recycled in a manner that protects human health.

Thank you

## Package Contents

| Description | Quantity |
|---|---|
| Hotplate W10 | 1 |
| Power cord (EU-plug, UK-plug) | 2 |
| Instruction manual | 1 |

## Installation

Unpack the delivery. Check that the unit is not damaged and that all parts are included.

Place the unit on the work area. Connect the unit to an electrical grounded mains power socket, use the units power cord (which has a with ground connection).

Check that the unit is OFF; no lights on the power switch or display. (If the unit is ON, green light on the power switch and red light on the display, then press the down side of the power switch, the lights goes OFF.)


## Intended use and overview

The product is for use at a pathology/histology laboratory for in vitro diagnostics, for preparation of human tissue sections on glass slides for microscopy. The units casing is made of white powder coated aluminium. The heated surface is black anodized.

The unit is equipped with a power switch with green light, a temperature control unit with digital display, internal heating element and an overheat protection fuse.

## Product Specifications

| | |
|---|---|
| Mains power: | 230 VAC 50/60Hz |
| Power consumption: | 260 W |
| Temperature range: | +20 to +89°C |
| Thermostat: | Digital |
| Temp. accuracy: | ±1°C |
| External dimensions: | W:146 x D:486 mm x H:70 mm |
| Weight: | 3,1 kg |

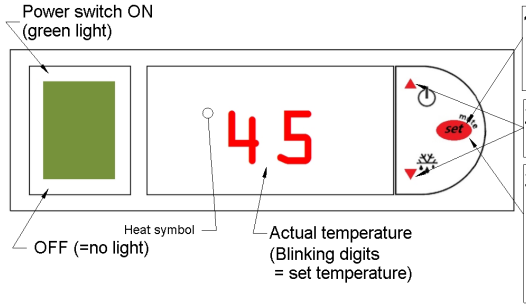**1** Press firmly with a finger on "set" for 1 second until the temperature digits starts to blink. Remove your finger.

**2** Press firmly on up or down arrow until the display shows the desired temp.

**3** Press firmly on "set" for ½ second. Remove your finger. Now the display shows the actual temperature. Wait some time for the water to reach the desired temperature.

## Getting Started

Turn the unit ON by pressing the upper part of the power switch; the power switch shows green light and the display shows actual temperature.

To change the desired temperature, see instructions **1-2-3** in the figure above.

If you press longer than 3 seconds and the display shows **P5** then reset the unit by pressing the power switch OFF and then ON. Start at 1 again, see above.

## Important Guidelines

First turn the unit OFF. Then disconnect the power cord first from the mains supply and then from the unit.

Please observe that the heated area can be hot.

Service and repair should only be made by an authorized technician.

# Troubleshooting

Review the information in the table below to troubleshoot operating problems.

| Problem | Cause | Solution |
|---|---|---|
| No heat | Too low set temperature | Increase the desired (set) temperature |
| | Heating element problem | Repair by an authorized technician |
| No power, no light | Power cord not firmly attached | Check the connection to mains power and to the unit |
| | The units fuse has blown | Replace the units fuse, located under the power cord socket |
| | The internal thermo fuse has blown | Repair by an authorized technician |
| Power switch has green light but no light on the display | Control unit not functioning | Repair by an authorized technician |

# Technical service

**Web Resources**

Visit the VWR website at **www.vwr.com** for:

- Complete technical service contact information
- Access to the VWR Online Catalogue, and information about accessories and related products
- Additional product information and special offers

**Contact us** For information or technical assistance contact your local VWR representative or visit. **www.vwr.com**.

# Warranty

**VWR International** warrants that this product will be free from defects in material and workmanship for a period of two (2) years from date of delivery. If a defect is present, VWR will, at its option and cost, repair, replace, or refund the purchase price of this product to the customer, provided it is returned during the warranty period. This warranty does not apply if the product has been damaged by accident, abuse, misuse, or misapplication, or from ordinary wear and tear. If the required maintenance and inspection services are not performed according to the manuals and any local regulations, such warranty turns invalid, except to the extent, the defect of the product is not due to such non-performance.

Items being returned must be insured by the customer against possible damage or loss. This warranty shall be limited to the aforementioned remedies. IT IS EXPRESSLY AGREED THAT THIS WARRANTY WILL BE IN LIEU OF ALL WARRANTIES OF FITNESS AND IN LIEU OF THE WARRANTY OF MERCHANTABILITY.

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# Declaration of Conformity

Undersigned representative for the following Manufacturer:

| | |
|---|---|
| Manufacturer: | KUNZ INSTRUMENTS AB |
| Address: | Tradgardsgatan 9, SE-14934 Nynashamn, Sweden |

Declare hereby that the following product: 720-2422 Hotplate W10

Comply with the following EC Directive and standards:

| Directive | Title |
|---|---|
| 98/79/EC | In Vitro Diagnostic Medical Devices (IVD), Annex III: other |
| 2011/65/EU | RoHS2 |

| Standard | Title |
|---|---|
| IEC/SS-EN 61010-1:2001 | Safety requirements for electrical equipment for measurement, control and laboratory use Part 1: General requirements |
| IEC/SS-EN 61010-2-101:2002 | Safety requirements for electrical equipment for measurement, control and laboratory use Part 2-101: Particular requirements for in vitro diagnostic (IVD) medical equipment |
| IEC/SS-EN 61326-1:1997 | Electrical equipment for measurement, control and laboratory use - EMC requirements |
| IEC/SS-EN 61326-2-6:2006 | Electrical equipment for measurement, control and laboratory use - EMC requirements Part 2-6: Particular requirements - In-vitro diagnostic (IVD) medical equipment |
| SS-EN ISO 14971:2012 | Medical devices - Application of risk management to medical devices |

Nynashamn December 5, 2014

Anders Olsson, Manufacturers representative

# Local VWR offices in Europe and Asia Pacific

**Australia**
VWR International Pty.LTD
Unit 1/31 Archimedes Place
Murarrie, Queensland 4172
Tel.: 1300 727 696
Fax: 1300 135 123

**Austria**
VWR International GmbH
Graumanngasse 7
1150 Vienna
Tel.: +43 1 97 002 0
Fax: +43 1 97 002 600
E-mail: info@at.vwr.com

**Belgium**
VWR International bvba
Researchpark Haasrode 2020
Geldenaaksebaan 464
3001 Leuven
Tel.: 016 385 011
Fax: 016 385 385
E-mail:
customerservice@be.vwr.com

**China**
VWR (Shanghai) Co., Ltd
Rm.219
No. 2100, Dongming Road
Pudong New District
Shanghai, 200126
Tel.: +86-21 589 868 88
Fax: +86-21 585 588 01
E-mail: info_china@vwr.com

**Czech Republic**
VWR International s. r. o.
Veetee Business Park
Pražská 442
CZ - 281 67 Stříbrná Skalice
Tel.: +420 321 570 321
Fax: +420 321 570 320
E-mail: info@vitrum.cz

**Denmark**
VWR - Bie & Berntsen
Transformervej 8
2730 Herlev
Tel.: 43 86 87 88
Fax: 43 86 87 90
E-mail: info@dk.vwr.com

**Finland**
VWR International Oy
Valimotie 9
00380 Helsinki
Tel.: +358 9 80 45 51
Fax: +358 9 80 45 52 00
E-mail: info@fi.vwr.com

**France**
VWR International S.A.S.
Le Périgares – Bâtiment B
201, rue Carnot
94126 Fontenay-sous-Bois cedex
Tel.: 0 825 02 30 30 (0,15 EUR TTC/min)
Fax: 0 825 02 30 35 (0,15 EUR TTC/min)
E-mail: info@fr.vwr.com

**Germany**
VWR International GmbH
Hilpertstrasse 20a
D - 64295 Darmstadt
Freecall: 0800 702 00 07
Fax: 0180 570 22 22*
Email: info@de.vwr.com
*0,14 €/Min. aus d. dt. Festnetz

**Hungary**
VWR International Kft.
Simon László u. 4.
4034 Debrecen
Tel.: (52) 521-130
Fax: (52) 470-069
E-mail: info@hu.vwr.com

**India**
VWR Lab Products Private Limited
135/12, Brigade Towers, 2nd Floor
Front wing, Brigade Road,
Bengaluru , India – 560 025
Tel: +91-2522-647911/922 (Mumbai)
Tel: +91-80-41117125/26 (Bangalore)
Fax +91-80-41117120
E-mail: vwr_india@vwr.com

**Ireland / Northern Ireland**
VWR International Ltd / VWR International
(Northern Ireland) Ltd
Orion Business Campus
Northwest Business Park
Ballycoolin
Dublin 15
Tel.: 01 88 22 222
Fax: 01 88 22 333
E-mail sales@ie.vwr.com

**Italy**
VWR International PBI S.r.l.
Via San Giusto 85
20153 Milano (MI)
Tel.: 02-3320311/02-487791
Fax: 02-332031307/02-40090010
E-mail: info@it.vwr.com

**The Netherlands**
VWR International B.V.
Postbus 8198
1005 AD Amsterdam
Tel.: 020 4808 400
Fax: 020 4808 480
E-mail: info@nl.vwr.com

**New Zealand**
VWR International LP
241 Bush Road
Albany
0632 New Zealand

**Norway**
VWR International AS
Haavard Martinsens vei 30
0978 Oslo
Tel.: 02290
Fax: 815 00 940
E-mail: info@no.vwr.com

**Poland**
VWR International Sp. z o.o.
Limbowa 5
80-175 Gdansk
Tel.: 058 32 38 210
Fax. 058 32 38 205
E-mail: labart@pl.vwr.com

**Portugal**
VWR International - Material de
Laboratório, Lda
Edifício Neopark
Av. Tomás Ribeiro, 43- 3 D
2790-221 Carnaxide
Tel.: 21 3600 770
Fax: 21 3600 798/9
E-mail: info@pt.vwr.com

**Singapore**
VWR Singapore Pte Ltd
18 Gul Drive
Singapore 629468
Tel: +65 6505 0760
Fax: +65 6264 3780
E-mail: sales@sg.vwr.com

**Spain**
VWR International Eurolab S.L.
C/ Tecnología 5-17
A-7 Llinars Park
08450 - Llinars del Vallès
Barcelona
Tel.: 902 222 897
Fax: 902 430 657
E-mail: info@es.vwr.com

**Sweden**
VWR International AB
Fagerstagatan 18a
163 94 Stockholm
Tel.: 08 621 34 00
Fax: 08 621 34 66
E-mail: info@se.vwr.com

**Switzerland**
VWR International GmbH
Lerzenstrasse 16/18
8953 Dietikon
Tel.: 044 745 13 13
Fax: 044 745 13 10
E-mail: info@ch.vwr.com

**Turkey**
Pro Lab Laboratuar Teknolojileri Ltd.Şti.
a VWR International Company
Orta Mah. Cemal Gürsel Caddesi
Ördekcioglu Işmerkezi No.32/1
34896 Pendik - Istanbul
Tel.: +90216 598 2900
Fax: +90216 598 2907
Email: info@pro-lab.com.tr

**UK**
VWR International Ltd
Customer Service Centre
Hunter Boulevard
Magna Park
Lutterworth
Leicestershire
LE17 4XN
Tel.: 0800 22 33 44
Fax: 01455 55 85 86
E-mail: uksales@uk.vwr.com