INSTALLATION GUIDE

RephiDuo P Pack TL1 with RO cleaning media

RephiLe Catalogue Number:

RR500CPC1

Comparable Millipore Part Number:

PR0GTLCS1


Suitable Millipore Systems:

RiOs 30 - 200
Elix 20 – 100. Elix Clinical


This product is a comparable to the corresponding Millipore cartridges. Please refer to Millipore's User Manual for a detailed instruction on installation.

A brief description of the installation process is described below:


1. Place the system into **Standby** by pressing the STANDBY button for about 2 seconds

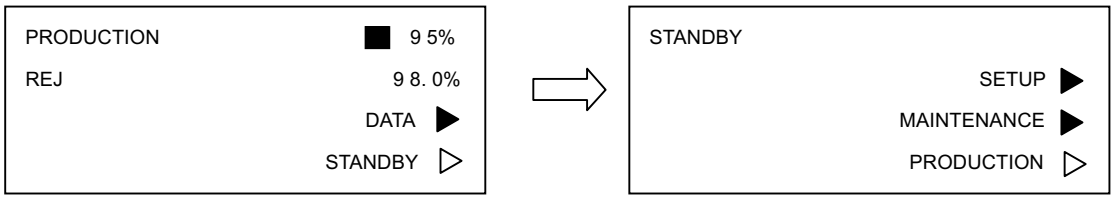


2. Press the MAINTENANCE button

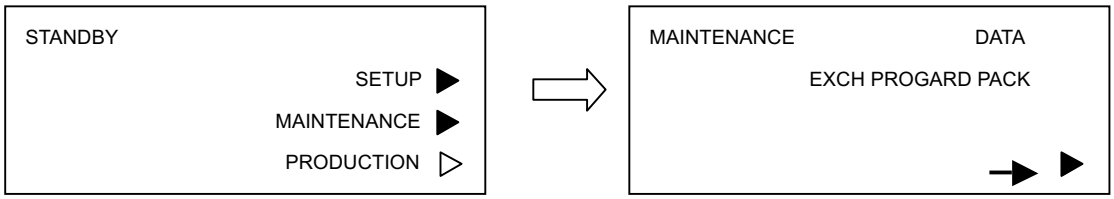


3. Press the EXCH PROGARD PACK button for 2 seconds to get the system into depressurize and purge mode. System will purge for 12 seconds during which time the pump will operate.

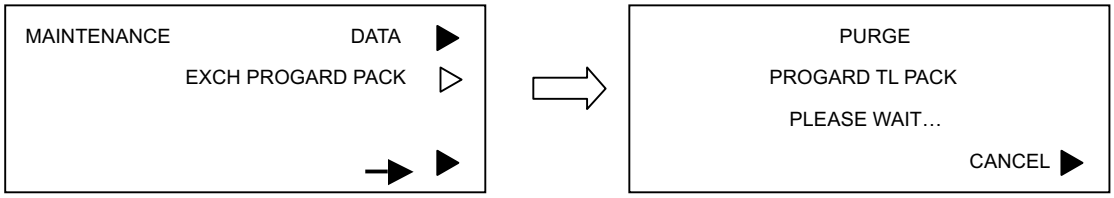


4. Wait till the purge completes. The system is now ready to have cartridge replaced.

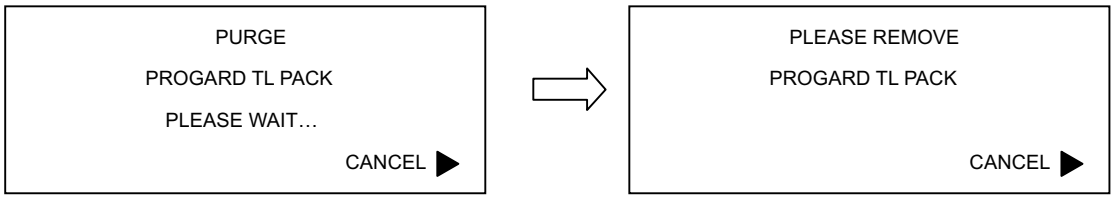


4.1 Lift up the cartridge Locking Handle as shown in the picture.

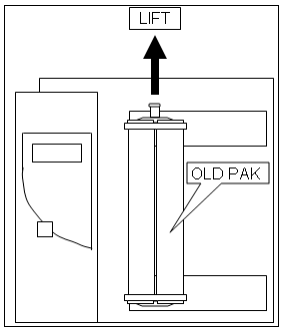

4.2 Gently pull the used cartridge out of the system. Make sure no old O-rings are left on the four ports of the water system.

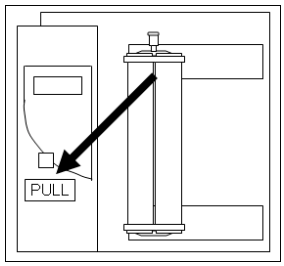


4.3 Take out the new cartridge; remove the four blue protective caps on the cartridge nozzles. Make sure each O-ring inside nozzles is firmly pushed in and seated flat.

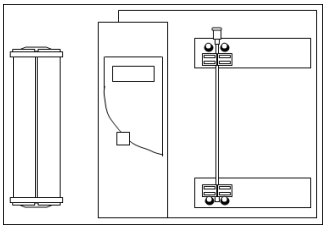

4.4 Install the cartridge onto the system with the **THIS END UP** arrow on the cartridge pointing up and all four nozzles are firmly pushed in.


Note: Cartridge is directional. It must be installed with the THIS SIDE UP arrow pointing up.

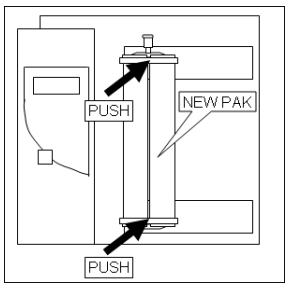

4.5 Push down the cartridge Locking Handle until it is securely locked.


If the pin does not lock all the way down, the cartridge probably has not been pushed all the way into ports.


Reinstall the cartridge if necessary till the lock pin is fully locked.

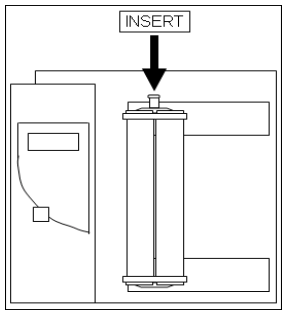

5 Press the CONTINUE a few times till the system goes into cleaning mode

Note: If you have a Storage Reservoir with a Vent Filter, the LCD will remind you to change the Vent Filter at this time. You can press the CONTINUE button one more time to skip it.

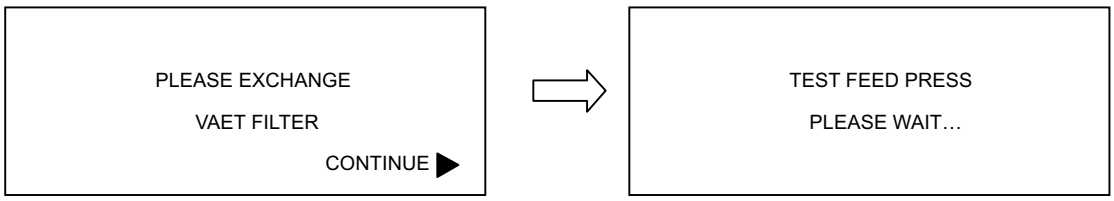

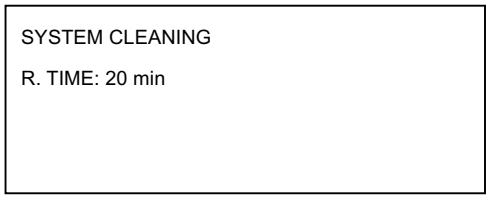

System goes into Cleaning for 20 minutes. Once the cleaning cycle completes, system will produce water automatically.

RephiSolo, RephiDuo and RephiQuatro are trademarks of RephiLe Bioscience, Ltd.

© 2010. RephiLe Bioscience, Ltd. All rights reserved.

All other trademarks are trademarks of their respective manufacturers.