


Test Chromatogram

PRODUCT DETAILS

Catalogue Number:	HIRPB-250A	Serial Number:	M2402263
Material:	Hichrom 5µm RPB	Dimensions:	250 x 4.6 mm
Batch:	B1223	Data File:	A40007.dat

TEST CONDITIONS

Mobile Phase:	MeOH:H ₂ O (85:15)	Storage Solvent:	Mobile Phase
Flow Rate:	1.0 mL/min	Wavelength:	254 nm
Pressure:	104 bar	Temperature:	22°C

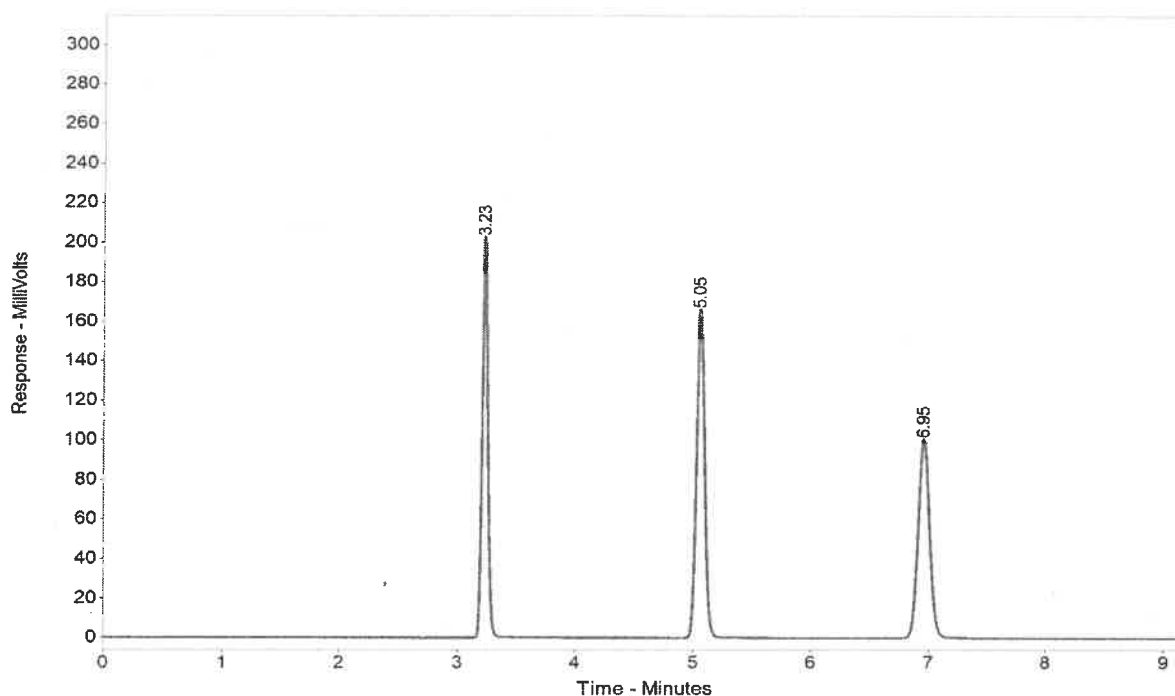

COLUMN PERFORMANCE

Peak	Name	t _r (minutes)	Efficiency N _{theo}	Asymmetry*	
				A _{s1}	A _{s2}
1	Dimethyl Phthalate	3.23	23300	92	1.24
2	Toluene	5.05	27200	95	1.17
3	Biphenyl	6.95	24200	95	1.12

*See overleaf for an explanation of asymmetry measurements.