

V. Moving the Rack

Caution: Do not carry the rack by the Knob (D).
Always use main column on the Base (A) (see Fig. 5).



Fig. 5



620 Lakeview Parkway • Vernon Hills, IL 60061
 Phone: 1-847-816-5070 • Fax: 1-847-816-5072
 Email: info@heathrowscientific.com
 Web: www.heathrowscientific.com

V. Moving the Rack

Caution: Do not carry the rack by the Knob (D).
Always use main column on the Base (A) (see Fig. 5).




Fig. 5



620 Lakeview Parkway • Vernon Hills, IL 60061
 Phone: 1-847-816-5070 • Fax: 1-847-816-5072
 Email: info@heathrowscientific.com
 Web: www.heathrowscientific.com

Flip & Grip™ Pipettor Rack Instructions

Part listing:

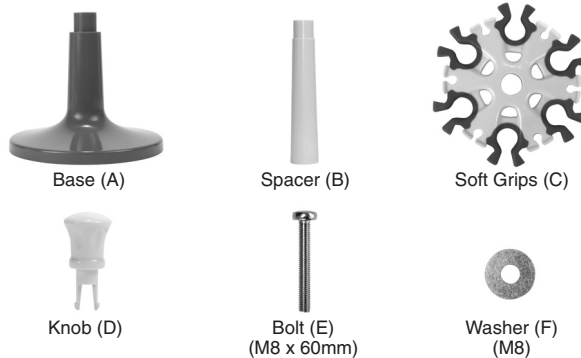

Tools needed:

- Phillips head #2 or #3 screw driver
- Drill (bit size) $\frac{5}{16}$ ", 8mm or O

I. Bench top use assembly


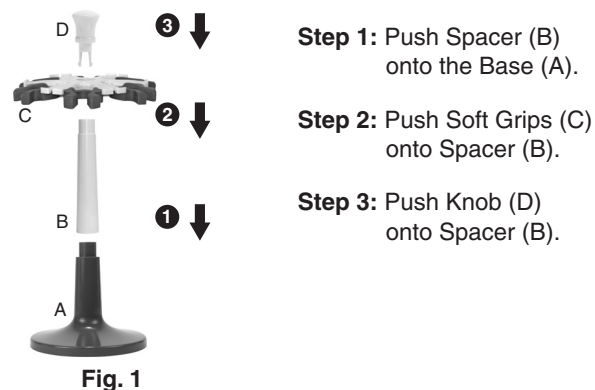


Fig. 1

Flip & Grip™ Pipettor Rack Instructions

Part listing:

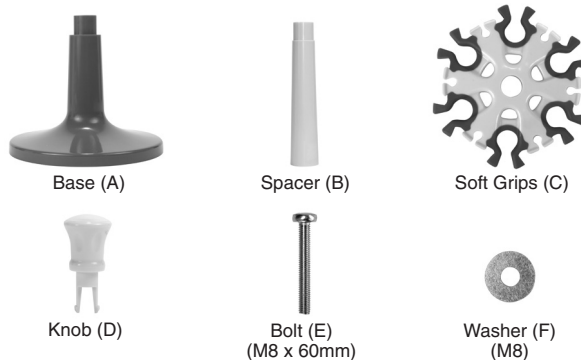

Tools needed:

- Phillips head #2 or #3 screw driver
- Drill (bit size) $\frac{5}{16}$ ", 8mm or O

I. Bench top use assembly

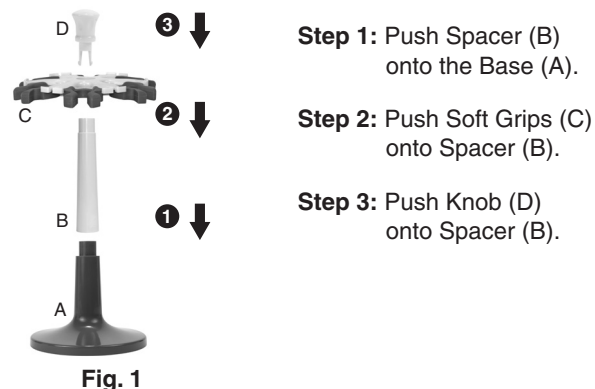


Fig. 1

II. Shelf mount assembly

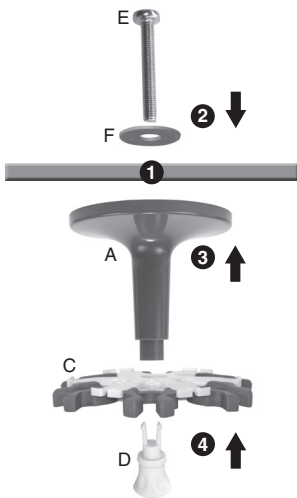


Fig. 2

Step 1: Drill $\frac{5}{16}$ " or 8mm hole in shelf.

Step 2: Insert the Bolt (E) and Washer (F) as shown in Fig. 2.

Step 3: Flip the Base (A) and screw the Bolt (E) into the Base (A).

Step 4: Assemble the Soft Grips (C) to the Base (A) with the Knob (D) as shown in Fig. 2.

III. Changing Rack Storage Mode

- Remove Knob (D) by pulling firmly.
- Refer to the steps in previous section I or II.

IV. Rack use

- Hold the Knob (D) to steady the rack.
- Load the pipettor into the Soft Grip (C) as shown in Fig. 3.
- Hang the pipettor on the rack as shown in Fig. 4.
- Hold Knob (D) firmly when removing the pipettor.

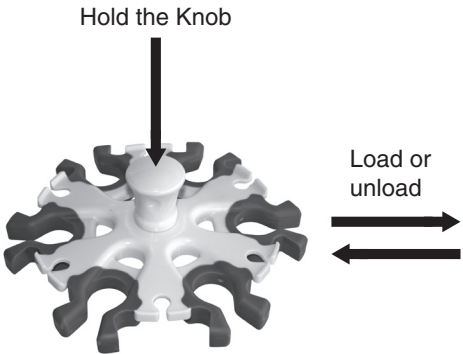


Fig. 3



Fig. 4

II. Shelf mount assembly

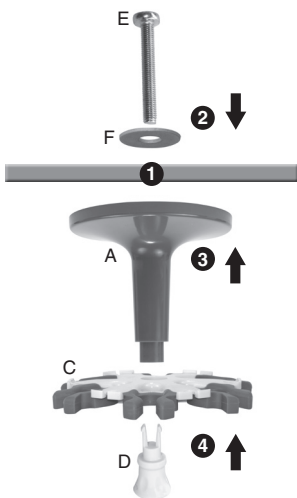


Fig. 2

Step 1: Drill $\frac{5}{16}$ " or 8mm hole in shelf.

Step 2: Insert the Bolt (E) and Washer (F) as shown in Fig. 2.

Step 3: Flip the Base (A) and screw the Bolt (E) into the Base (A).

Step 4: Assemble the Soft Grips (C) to the Base (A) with the Knob (D) as shown in Fig. 2.

III. Changing Rack Storage Mode

- Remove Knob (D) by pulling firmly.
- Refer to the steps in previous section I or II.

IV. Rack use

- Hold the Knob (D) to steady the rack.
- Load the pipettor into the Soft Grip (C) as shown in Fig. 3.
- Hang the pipettor on the rack as shown in Fig. 4.
- Hold Knob (D) firmly when removing the pipettor.

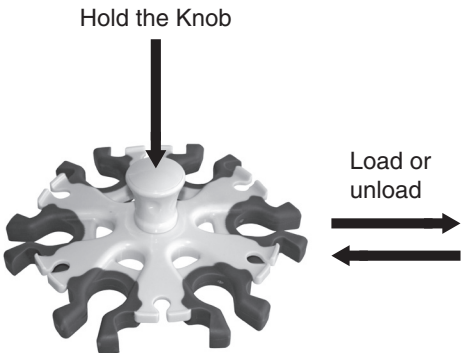


Fig. 3



Fig. 4