Instruction Manual

VWR® Large Capacity Stirrers

EU cat. no

Model 1000: 444-7115

Model 2000: 444-7118

EN	-	English		1
FR	-	Français		11
ES	-	Español		21
DE	-	Deutsch		31
IT	-	Italiano		41
PT	-	Português		51
BG	-	Български		61
CZ	-	Čeština		64
DK	-	Dansk		67
EE	-	Eesti		70
FI	-	Suomi		73
GR	-	Ελληνικά		76
HR	-	Hrvatska		79
HU	-	Magyar		82
IS	-	Íslenska		85
LT	-	Lietuvių		88
LV	-	Latviešu		91
NL	-	Nederlands		94
NO	-	Norsk		97
PL	-	Polski		100
RO	-	Română		103
SI	-	Slovenščina		106
SK	-	Slovenčina		109
SV	-	Svenska		112

Revision D 05/2026

Legal address of Manufacturer

Europe
VWR International bv,
Researchpark Haasrode 2020
Geldenaaksebaan 464 B-3001
Leuven
+3216385011
be.vwr.com

UK Importer:
VWR International Ltd
Hunter Boulevard, Magna Park
Lutterworth, Leicestershire, LE17 4XN
uk.vwr.com

Country of Origin

United States

Table of Contents

Safety Instructions 4

Package Contents 4

Installation 4

Intended Use 4

Maintenance & Servicing 5

Environmental Conditions 5

Specifications 6

Control Panel 6

Operating Instructions

Getting Ready 7

Setting Speed 7

Setting Time to Zero and Continuous Mode 7

Setting Timed Mode 7

Turning Unit Off 7

Troubleshooting 8

Technical Service 8

Warranty 8

Disposal 9

Distributor Contact Information 10

Safety Instructions

Please read the entire instruction manual before operating the Large Capacity Stirrer.

WARNING! DO NOT use the Large Capacity Stirrer in a hazardous atmosphere or with hazardous materials for which the unit was not designed. Also, the user should be aware that the protection provided by the equipment may be impaired if used with accessories not provided or recommended by the manufacturer or used in a manner not specified by the manufacturer.

Always operate unit on a level surface for best performance and maximum safety.

DO NOT lift unit by the top plate.

Heavy walled glass or polypropylene container recommended when stirring at high speeds. If spin bar de-couples at higher speeds, thin-walled glass containers may break.

CAUTION! To avoid electrical shock, completely cut off power to the unit by disconnecting the power cord from the unit or unplug from the wall outlet. Disconnect unit from the power supply prior to maintenance and servicing.

Spills should be removed promptly. **DO NOT** immerse the unit for cleaning.
DO NOT operate the unit if it shows signs of electrical or mechanical damage.

Earth Ground - Protective Conductor Terminal
Alternating Current

Package Contents

Model 1000 or Model 2000 Large Capacity Stirrer

Description	Quantity
234cm detachable power cords (EU and UK) CH adapter	3
Spin bar	1

Installation

Upon receiving the VWR Large Capacity Stirrer, check to ensure that no damage has occurred in shipment. It is important that any damage that occurred in transport is detected at the time of unpacking. If you do find such damage the carrier must be notified immediately.

After unpacking, it requires two (2) people to lift the Large Capacity Stirrer, from the bottom of the unit, to place on a level bench, table, or floor, away from explosive vapors. Ensure that the surface on which the unit is placed will withstand typical heat produced by the unit. Always place the unit on a sturdy work surface.

The Large Capacity Stirrer is supplied with a power cord that is inserted into the IEC connector on the left side of the unit first, then it can be plugged into a properly grounded outlet. The 120V unit plugs into a 120-volt, 50/60 Hz source. The 230V unit plugs into a 230-volt, 50/60 Hz source.

Intended use

For research use only. Not intended for any animal or human therapeutic or diagnostic use.

Maintenance & Servicing

The Large Capacity Stirrer is built for long, trouble-free, dependable service. No lubrication or other technical user maintenance is required. It needs no user maintenance beyond keeping the surfaces clean. The unit should be given the care normally required for any electrical appliance. Avoid wetting or unnecessary exposure to fumes. Spills should be removed promptly. Do not use a cleaning agent or solvent on the front panel or top plate which is abrasive or harmful to plastics, nor one which is flammable. Always ensure the power is disconnected from the unit prior to any cleaning. If the unit ever requires service, contact your VWR representative.

Environmental Conditions

Operating Conditions: Indoor use only.

Temperature: 5 to 40°C (41 to 104°F)

Humidity: 20% to 80% relative humidity, non-condensing

Altitude: 0 to 6,562 ft (2000 M) above sea level

Non-Operating Storage:

Temperature: -20 to 65°C (-4 to 149°F)

Humidity: 20% to 80% relative humidity, non-condensing


Installation Category II and Pollution Degree 2 in accordance with IEC 664.

Product Specifications

Overall dimensions (L x W x H): 72,4 x 68 x 13,3 cm
Top plate dimensions (L x W): 54 x 63,5 cm
Electrical (50/60 Hz):
Model 1000 230 volts | 0,5 amps | 41 watts
Model 2000 230 volts | 0,5 amps | 82 watts
Fuses: 5mm x 20mm, 5 amp quick acting
Speed range: 100 to 1800rpm
(increased in 10rpm increments)
Speed stability: ±2%
Weight capacity: 193 kg
Volume (H2O):
Model 1000 100L
Model 2000 200L
Timer: digital, 1 second to 9999 minutes
(increased in 1 second increments)
Controls: see page 6
Ship weight: 28,6 kg

Control Panel

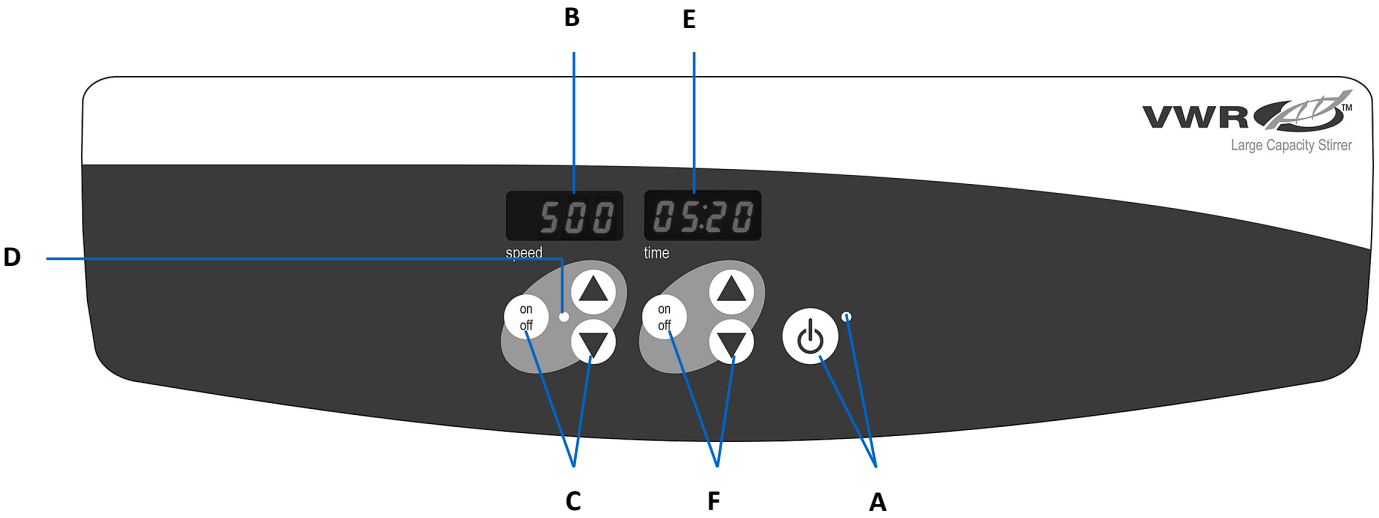The front panel of the Large Capacity Stirrer contains all the controls and displays needed to operate the unit.

A. Standby button/standby indicator light: The standby indicator light will illuminate when the unit is plugged in. The unit will be in standby mode. Press the standby button to activate the speed and time functions. The standby indicator light will shut off and the speed display and time display will illuminate. Press the standby button again and the unit will once again be in standby mode.

B. Speed display: Displays the speed of the stirrer.

C. Up/down arrows for set-point control. On/off button starts/stops stir function.

D. The stir indicator light will be illuminated when the unit is stirring.

E. Time display: Displays accumulated time (continuous mode) or how much time is remaining (timed mode).

F. Up/down arrows for set-point control. On/off button starts/stops the time function. The display range is from 0 to 9,999 minutes in one (1) second increments. The display will indicate minutes and seconds until the timer reaches 99 minutes and 59 seconds (99:59), then the display will automatically display minutes up to 9,999.

Operating Instructions

The Large Capacity Stirrer has been designed for the speed and time functions to work independently of one another. The speed can be re-set without re-setting the timer and the timer can be stopped and started without interrupting the stirring function.

1. Getting Ready

a. Plug the power cord into a properly ground outlet. The standby indicator light will illuminate, verifying power to the unit.
b. Press the standby button to move the unit from standby mode. The standby indicator light will turn off and the speed/time displays will illuminate and display the previously used settings.

2. Setting Speed

a. Press the up/down arrows below the speed display until you reach the desired speed. When you release the arrow button, the unit will blink off and then on indicating the new set speed has been accepted.
b. Press the on/off button to start the stirring function. The indicator light below the speed display will illuminate to indicate the stirring function is in-use. While accelerating, the indicator light will blink until set-point is reached and then will remain lit until stirring has ceased. The microprocessor controlled ramping feature slowly increases speed until the set-point is reached which helps to avoid splashing and provides excellent low-end control. The microprocessor also regulates the stirring speed, sensing your requirements whether stirring an aqueous, viscous, or semi-solid solution.
c. Speed adjustments can be made without interrupting stirring by using the up/down arrows below the speed display. After the change has been made and you release the arrow button, the unit will blink off and then on indicating the new set speed has been accepted.
d. To stop the stirring function, press the on/off button below the speed display.

3. Setting Time to Zero (0:00) and Continuous Mode: Accumulated Time.

a. Hold the on/off button below of the time display. After three (3) seconds, the display will indicate the previous set time.
b. Simultaneously press both the up and the down arrows, the display will indicate zero (0:00). The unit time is now set to zero (0:00) minutes. Alternately, you can use the up/down arrows to get to zero (0:00).
c. Press the on/off button below the time display. The display will indicate accumulated time. The up/down arrows will become inactive. To stop timer, press the on/off button again. **IMPORTANT:** This will **NOT** interrupt the shaking function. Press the on/off button below the speed display to interrupt the shaking function.
d. To reset, press and hold the on/off button below the time display. After three (3) seconds the display will indicate the previous set time, which was zero (0:00).

4. Setting Timed Mode: Programmed Time.

a. Press the up/down arrows below the time display until you reach the desired time.
b. Start this function by pressing the on/off button below the time display. The unit will run for the selected time, the up/down arrows will become inactive while the timer is running. The unit will stop shaking when time display reaches zero (0:00). Four (4) audible beeps will indicate the time down function is complete. The time display will default back to the set time. To repeat this function for the same set time, simply depress the on/off button again.
c. To interrupt an automatic timing cycle before it is completed, press the on/off button below the time display. The display will flash off and on to indicate the time function is on “hold”. **IMPORTANT:** This will **NOT** interrupt the shaking function. Press the on/off button below the speed display to interrupt the shaking function. Restart the timer by pressing the on/off button below the time display. Unit will continue counting down to zero (0:00). When the display reaches zero (0:00), you will hear the four (4) audible beeps that indicate the countdown function is complete and the shaking function will stop.

5. Turning Unit Off

a. To turn the unit off, press the standby button. The speed and time displays will be blank, the standby indicator light will illuminate. The stirrer should be kept in standby mode when not in use. To completely cut off power to the unit, disconnect the power cord from the unit or un-plug from the wall outlet.

Operating Tips

Built-in memory maintains the last used speed and time settings during a power interruption.
The stirrer increases speed at a steady rate until the set-point is reached, if the spin bar is too large or the liquid is too viscous, the stirrer may not reach its set-point or become de-coupled from the units drive magnet. The set-point speed needs to be reduced. The unit is designed to ramp slowly to assure that the spin bar stays coupled.

* The model 2000/200L unit contains a larger spin bar with a center pivoting ring. Ring may need to be replaced after eight (8) hours of use. Contact your VWR representative for replacement spin bars or rings.

Troubleshooting

Review the information in the table below to troubleshoot operating problems.

Problem	Cause	Solution
E1	Motor locked/ Failure	Error should not be addressed by the end user. Switch the unit off and contact your VWR Representative for repair.

Technical service

Web Resources

Visit the VWR website at **www.vwr.com** for:

- Complete technical service contact information
- Access to the VWR Online Catalogue, and information about accessories and related products
- Additional product information and special offers

Contact us For information or technical assistance contact your local VWR representative or visit. www.vwr.com.

Warranty


VWR warrants that this product will be free from defects in material and workmanship for a period of five (5) years from date of delivery. If a defect is present, VWR will, at its option and cost, repair, replace, or refund the purchase price of this product to the customer, provided it is returned during the warranty period. This warranty does not apply if the product has been damaged by accident, abuse, misuse, or misapplication, or from ordinary wear and tear. If the required maintenance and inspection services are not performed according to the manuals and any local regulations, such warranty turns invalid, except to the extent, the defect of the product is not due to such non-performance.

Items being returned must be insured by the customer against possible damage or loss. This warranty shall be limited to the aforementioned remedies. IT IS EXPRESSLY AGREED THAT THIS WARRANTY WILL BE IN LIEU OF ALL WARRANTIES OF FITNESS AND IN LIEU OF THE WARRANTY OF MERCHANTABILITY.

Compliance with local laws and regulations

The customer is responsible for applying for and obtaining the necessary regulatory approvals or other authorisations necessary to run or use the Product in its local environment. VWR will not be held liable for any related omission or for not obtaining the required approval or authorisation, unless any refusal is due to a defect of the product.

Equipment disposal

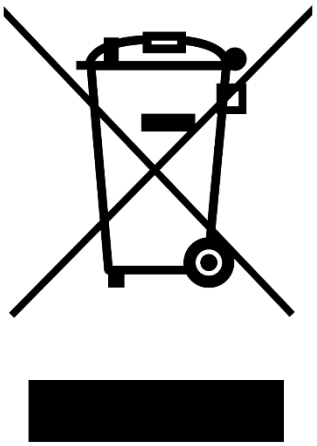This equipment is marked with the crossed out wheeled bin symbol to indicate that this equipment must not be disposed of with unsorted waste.
Instead it's your responsibility to correctly dispose of your equipment at lifecycle -end by handling it over to an authorized facility for separate collection and recycling. It's also your responsibility to decontaminate the equipment in case of biological, chemical and/or radiological contamination, so as to protect from health hazards the persons involved in the disposal and recycling of the equipment.
For more information about where you can drop off your waste of equipment, please contact your local dealer from whom you originally purchased this equipment.

By doing so, you will help to conserve natural and environmental resources and you will ensure that your equipment is recycled in a manner that protects human health.

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India
VWR Lab Products Private Limited
No.139. BDA Industrial Suburb,
6th Main, Tumkur Road, Peenya Post,
Bangalore, India – 560058
Tel.: +91 80 28078400
Email: vwr_india@vwr.com

Ireland
VWR International Ltd
Orion Business Campus
Northwest Business Park
Ballycoolin
Dublin 15
Tel.: +353 01 88 22 222
Email sales.ie@vwr.com

Italy
VWR International S.r.l.
Via San Giusto 85
20153 Milano (MI)
Tel.: +39 02-3320311
Email: info.it@vwr.com

Korea
Avantor Performance Materials Korea Ltd
2F ACE Gwanggyo Tower I
Daehak 4ro 17
Yeongtong-gu Suwon
Korea 16226
Tel.: +82 31 645 7250
saleskorea@avantorsciences.com

Middle East & Africa
VWR International FZ-LLC
Office 203, DSP Lab Complex,
Dubai Science Park
Dubai, United Arab Emirates
Tel.: +971 4 5573271
info.mea@vwr.com

The Netherlands
VWR International B.V.
Postbus 8198
1005 AD Amsterdam
Tel.: +31 020 4808 400
Email: info.nl@vwr.com

Norway
VWR International AS
Brynsalleen 4
0667 Oslo
Tel.: +47 22 90 00 00
Email: info.no@vwr.com

Poland
VWR International Sp. z o.o.
Limbowa 5
80-175 Gdansk
Tel.: +48 058 32 38 200
Email: info.pl@vwr.com

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a. Hold the on/off button below of the time display. After three (3) seconds, the display will indicate the previous set time.

b. Simultaneously press both the up and the down arrows, the display will indicate zero (0:00). The unit time is now set to zero (0:00) minutes. Alternately, you can use the up/down arrows to get to zero (0:00).

c. Press the on/off button below the time display. The display will indicate accumulated time. The up/down arrows will become inactive. To stop timer, press the on/off button again. **IMPORTANT:** This will **NOT** interrupt the shaking function. Press the on/off button below the speed display to interrupt the shaking function.

d. To reset, press and hold the on/off button below the time display. After three (3) seconds the display will indicate the previous set time, which was zero (0:00).

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