Multilanguage short instruction manual & Troubleshooting

VWR® CONDUCTIVITY METER CO-330

EU Catalogue Number: 663-0359

Version: 1.0
Issued: 01, December, 2020

LANGUAGES


English....3
Français....5
Deutsch....7
Italiano....9
Español....11
Português13
Hrvatski....15
Český....17
Dansk19
Nederlandse....21
Suomen....23
Magyar....25
Latviešu....27
Lietuvių....29
Norsk....31
Polski....33
Română....35
Slovák....37
Slovenski....39
Svenska....41

FOR MORE DETAILS, PLEASE CONSULT THE FULL MANUAL

English


The CO-330 conductivity meter enables the measurement of salinity, conductivity (in NaCl, KCl, TDS) and temperature. Turn the device on with the key ON/OFF. Other operating activities are performed using the touch screen.

Sockets

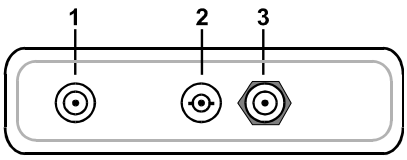
1	Power adapter socket
2	Conductivity sensor socket
3	Temperature sensor socket

Measurement screen

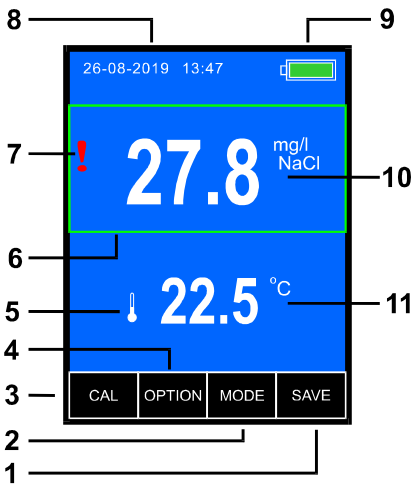
1	Collect results key
2	Parameter selection key
3	Key to enter the calibration mode (hold down for 2s)
4	Key for selecting the measurement function
5	Symbol of the connected temperature sensor
6	Fixed measurement frame
7	Conductivity sensor warning symbol
8	Date and time
9	Battery status
10	Measurement result with unit
11	Temperature measurement result with unit

Parameter screen

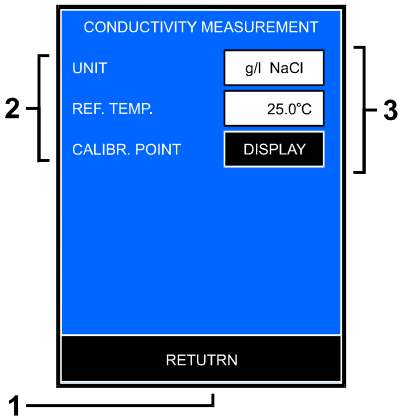1	Key to return to the previous screen
2	Description of parameters
3	Parameter selection keys

Calibration screen

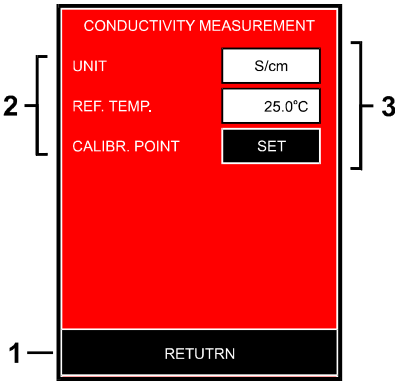1	Key to return to the previous screen
2	Calibration parameter options
3	Keys for entering parameters

TROUBLESHOOTING

Problem	Cause	Solution
Meter not switching on / No power	Discharged batteries	Charge batteries
	Incorrect battery polarity	Check the polarity
	Battery failure	Replace batteries with AA (R6) Rechargeable type
	Power adapter failure	Check with other 5V USB power adapter (for mobile)
Unstable reading, lack of reading	Poor connection	Check conductivity sensor connection to the unit
	Air bubbles in the conductivity sensor cell	Check for air bubbles, make few vertical moves with the cell in the solution
	Conductivity cell contaminated/broken	Clean the electrodes according to its manual or replace the conductivity sensor
	Break in cable, connector or internal in sensor	Check cable and connector, replace the sensor
	Measured value out of the sensor range	Use a different sensor or meter
	The temperature of the solution is fluctuating	Measure after the solution temperature stabilizes
	There is a motor or other device causing electrical interference.	Move the meter to a place where it is not subject to dielectric effects. Be sure to ground devices that are using commercial electricity.
	To short result stabilisation time	Response time may slow down, depending on the properties of the sample solution.
	Conductivity cell not immersed properly	Make sure the sample entirely covers the probe sensors Probe not immersed enough
	Conductivity sensor defective	Replace sensor
Meter is not calibrating	Sensor, electrodes contaminated	Clean electrodes or replace the sensor
	Calibration mode not entered properly	Press the result to show frame around it. Press and hold the CAL button till displaying red background colour.
	Intermittent or no connection of conductivity sensor	Check sensor connection to the unit
	The temperature of the standard solution is fluctuating	Measure after the standard solution temperature stabilizes
	Stabilisation time not long enough	Wait longer for result stabilisation
	Conductivity sensor defective	Replace sensor
Result different than expected	Conductivity sensor not calibrated	Calibrate the sensor
	Sensor polluted	Clean the sensor
	Air bubbles in the conductivity sensor cell	Check for air bubbles, make few vertical moves with the cell in the solution
	Reading stabilisation time to short	Longer wait for result stabilisation
	The temperature of the solution is fluctuating	Measure after the solution temperature stabilizes
	During calibration improper or old standard solution used or standard temperature not correct	Repeat calibration with new, correct standard solution value and correct settings in the meter and in an accurately controlled temperature of 25 °C
	Conductivity cell not immersed properly	Make sure the sample entirely covers the probe sensors
	Conductivity sensor defective	Replace sensor
Temp. value incorrect or  symbol displayed	Temperature sensor not connected or poor connection	Connect the temp. sensor / check the connection
	Breakage in temp sensor cable or connector / defective sensor	Check connection / replace sensor
Information: conductivity calibration readout error!	Conductivity sensor calibration data loss	Recalibrate the conductivity sensor sensor
	If the situation appears after each switching on	Contact distributor or service agent
Information: Manufacturer's settings readout error!	Internal hardware or software problem / Losing factory settings	Contact distributor or service agent

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