VWR®
Supplier Partnerships for Customer Solutions

Star-Blender™ LB 400

Lab blender

Instruction manual

Version : 1.0
22, January 2008

VWR®
Supplier Partnerships for Customer Solutions

Table of Contents

3 - 7 UK

3 . . . Warning, Safety Information, Package Contents, Intended Use
4. Specifications
5 Installation, Instructions for Use
6 . . . Troubleshooting, General Maintenance, Acessories and Spares
7 Technical Service, Warranty, Disposal

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28 - 43 Quick User Guide in several languages (FI, NL, HU, PT, PL, SK, SW, CZ)

VWR International BVBA

Haasrode Researchpark Zone 3 Geldenaaksebaan 464 B-3001 Leuven

+ 32 16 385011 **www.be.vwr.com**

Warning

The unit must be installed upright, on a stable, horizontal surface; never run the unit in any other position!
Always use an earthed plug.
Never open the housing of the machine without first unplugging the unit from the electric mains.
In normal use, do not run the **Star-Blender™ LB 400** for more than 15 minutes continuously.
Average settings of the unit to blend solid samples are 20 g diluted by 1/10 (depending on sample size and structure).


IF YOU ARE NOT SURE, PLEASE CONSULT US FOR ADVICE BEFORE USING!
The housing should only be removed by qualified and duly authorized staff.

Safety Information

Useful capacity : 50 - 400 ml. Maximum capacity depends on weight and hardness of the sample.


Star-Blender™ LB 400 is a laboratory unit. Industrial use, over-loading and daily non-stop use may damage your **Star-Blender™ LB 400** and invalidate your guarantee.

Package Contents

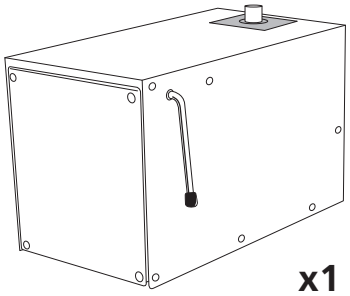
x1

• **Star-Blender™ LB 400**

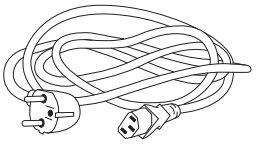x3

• Electric plugs (UK, EU, CH)

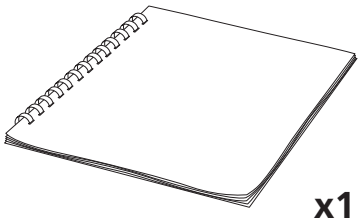x1

• Instruction manual


Intended Use

Star-Blender™ LB 400 is a powerful and easy-to-use bag blender for food sample preparation.
Adapted to any kind of application, it is the perfect tool for lab analyses.

UK

Specifications

Size (LxWxH): 39x26x28 cm
Net weight: 15 kg
Power: 230V - 50Hz - 300W Max
Nominal capacity: 50 - 400 ml. Maximum capacity depends on weight and hardness of the sample.
Dead space between door and paddles: 0 - 5 mm
Speed: 8 strokes/sec
Motor: 230V Mono
Motor speed: 2800/226 rpm
Amps: 0.7
Fuses: 1,25A FFT L 250 VRC
8 time settings: 30, 60, 90, 120, 150, 180, 210 secs. and ∞
Country of origin: France

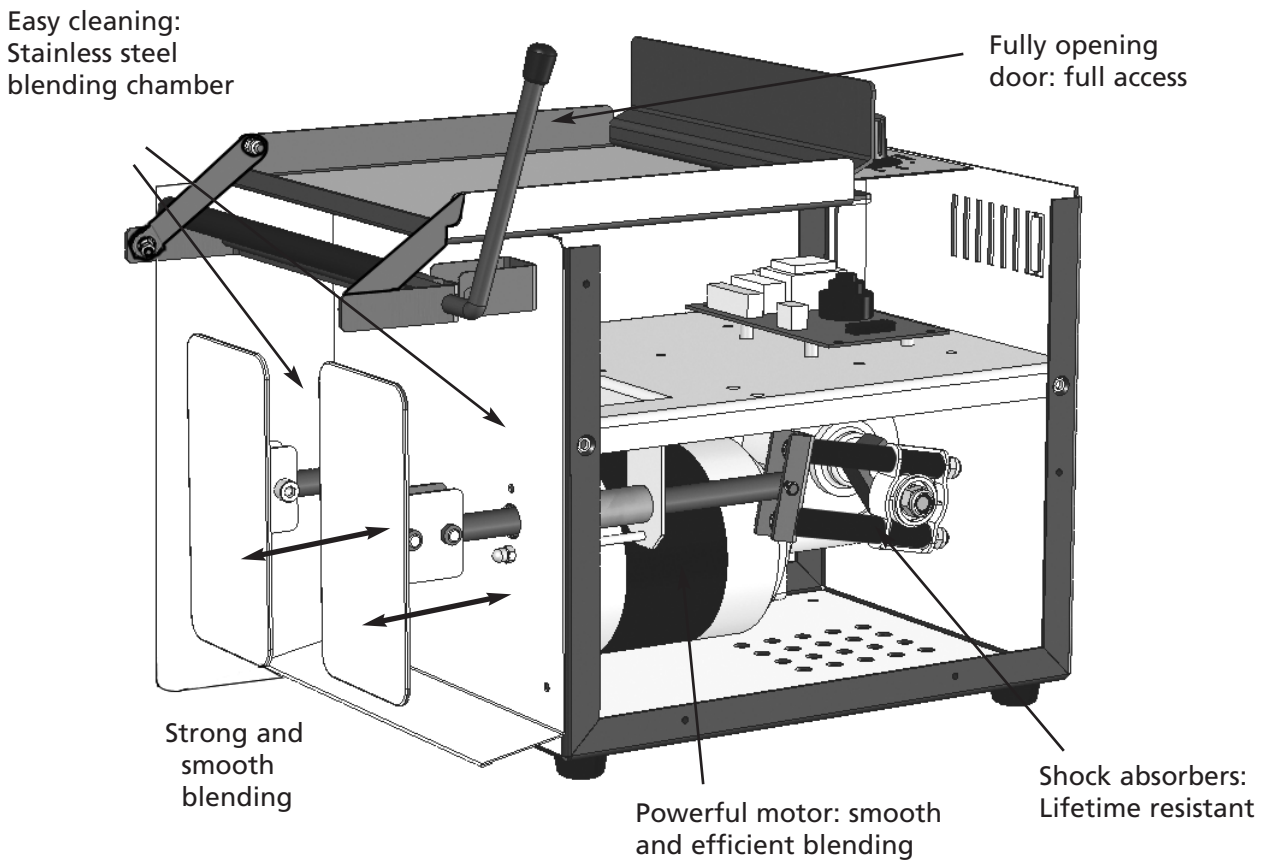Installation

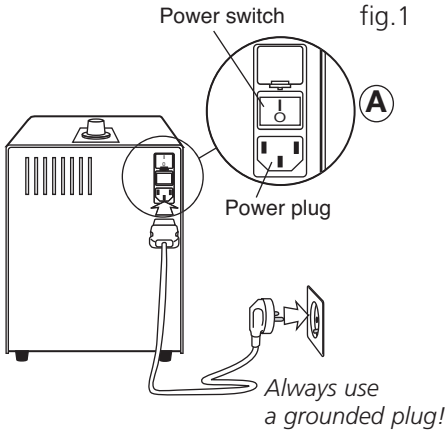
fig.1

1 Place **Star-Blender™ LB 400** on a flat, stable surface. Never run the mixer in any other position!

2 Plug the power cord (at the rear of the mixer) into the electric outlet. (fig.1)

3 Turn the main switch (at the rear of the mixer) to position "I". (fig.1) (A)

Instructions for Use

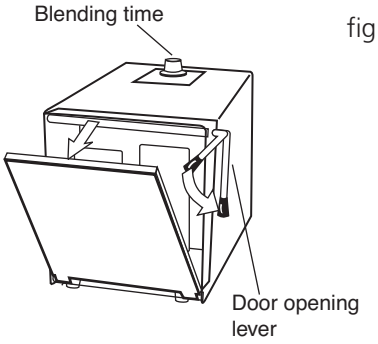
fig.2

In normal use, limit the non-stop cycle to 15 minutes

1 Adjust the timer setting. (fig.2)

2 Open the front door by pulling the lever forwards and upwards. (fig.2)

3 Insert the sterile filter bag containing the sample and the diluting fluid in the unit. (fig.3)

4 Close the door while holding the bag and lowering the lever as far as possible.

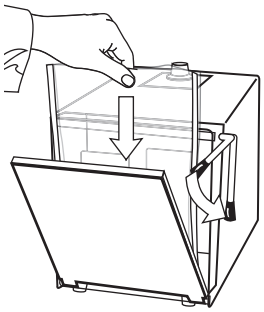
fig.3

5 **Star-Blender™ LB 400** starts and stops automatically.

6 When the cycle is over, lift out the bag. The sample is homogenized ; **Star-Blender™ LB 400** is ready for the next cycle. (fig.3)

Troubleshooting

Possible anomalies	Ask yourself
Star-Blender™ LB 400 does not start.	• Is the power cord connected? • Is the power switch set on position "I" ? (fig.1) • Is the mixer door correctly positioned and closed? Do not force. • Have you checked the fuse? (1,25A/220V 3,15A/110V)
The door closes incorrectly.	Is the door properly centered on its lower axis bar? Do not force
Blending is not satisfactory.	Is the timer adjusted to the appropriate blending time?
Bags leak during or at the end of the blending cycle.	Are you using good quality bags? We recommend VWR bags.

General Maintenance

Cleaning

Switch in "0" position (fig.1) Ⓐ ➡ **Unplug the power cord from the power outlet.**

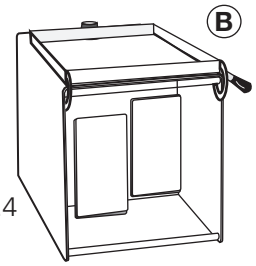
fig.4

Housing

- Clean the housing with standard cleaning products.
- Open the door and unhook the lower edge from the axis bar. (fig.4) Ⓑ

Blending chamber

- Clean the blending chamber and the paddles with stainless steel cleaning products.
- Reinstall the door, ensuring that it is correctly centered.

After maintenance

When maintenance is finished, after checking that the unit is not damp, plug in the power cord, and only then turn on the power switch.

Ordering Information

Product	Pack	Cat No.
Star-Blender™ LB 400	1	710-0632
Drip Tray	1	710-0633
Bag holder for 10 bags	1	710-0634
Automatic Straw pipette	1	613-1058
Clips for Lab Blender bags	200	129-9869
Sterile Straws for pipette : 1 ml, sterile by 25	1000	612-1675
Sterile Straws for pipette : 5 ml, sterile by 25	300	612-1676
400 ml Lab Blender Bags : plain & sterile	500	129-9867
400 ml Lab Blender Bags : lateral filter, sterile	500	129-9868
400 ml Lab Blender Bags : full page filter, sterile	500	129-9875

Technical Service

Web Resources

Visit the VWR's website at **www.vwr.com** for:
- Complete technical service contact information
- Access to VWR's Online Catalogue and information about accessories and related products
- Additional product information and special offers

Contact us

For information or technical assistance, contact your local VWR representative or visit **www.vwr.com**


Warranty

VWR International warrants that this product will be free from defects in material and workmanship for a period of 2 years from date of purshase. If a defect is present, VWR will, at its option, repair, replace or refund the purshase price of this product at no charge to you, provided it is returned during the warranty period. This warranty does not apply if the product has been damaged by accident, abuse, misuse or misapplication, or from ordinary wear and tear.
For your protection, items being returned must be insured against possible damage or loss. This warranty shall be limited to the replacement of defective products.

IT IS EXPRESSLY AGREED THAT THIS WARRANTY WILL BE IN LIEU OF ALL WARRANTIES OF FITNESS AND IN LIEU OF THE WARRANTY OF MERCHANDABILITY.

Disposal

The equipment is marked with the cross out wheeled bin symbol to indicate that this equipment must not be disposed off with unsorted waste.
Instead it is your responsability to correctly dispose off your equipment at lifecycle-end by handling it over to an authorized facility for separate collection and recycling.
It is also your responsability to decontaminate the equipment in case of biological, chemical and/or radiological contamination, so as to protect from health hazards the persons involved in the disposal and the recycling of the equipment.
For more information about where you can drop off your waste of equipment, please contact your local dealer from whom you originally purshased the equipment.

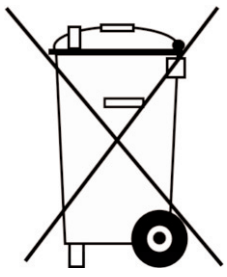By doing so, you will accept to conserve natural and environmental resources and will ensure that your equipment is recycled in a manner that protects human health.

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VWR International AB
VWR International AG
VWR International Ltd