


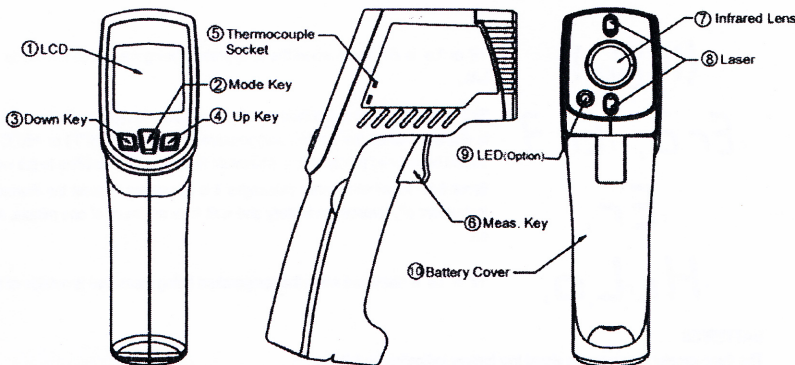

(Default Screen)
Simply aim the thermometer at the measure target with Lens (7) and press Meas. key (6) to display the surface

temperature. The Distance:Spot is 12:1. Please make sure the target area is within the field of view.

FUNCTION

Press Mode key (4) for scrolling more display function as follows.

Here will show the emissivity data. (The default emissivity is 0.95.)

Press Mode key (4), then press Up key (2) or Down key (3) to set the emissivity, then press Mode key (4) to confirm it. The emissivity can be changed from 0.10 (10E) to 1 (100E).

Press Mode key (4) for the Maximum (MAX), Minimum (MIN), Different between MAX and MIN (DIF) and Average (AVG) modes. During the measurement, the special modes reading will be displayed beside the mode icon.

Press Up key (2) or Down key (3) key to change the High Alarm (HAL) or Lo Alarm (LAL), then press Meas. key (6) to confirm it. For example: When the reading 26.9°C < LAL 27°C, the Low icon will flash and you will hear a beep sound.

Connect the thermocouple with Thermocouple socket (5) and put the probe in/on the target, the thermometer will display the temperature automatically without pressing any button. To see the minimum or maximum data during the probe measurement, please hold down the Up key (2) or Down key (3).




** The thermometer will automatically shut off if left idle for more than 60sec, unless in PRB mode. (In PRB mode, it will shut off if left idle for more than 12 minutes.)

In E, MAX, MIN, DIF, AVG mode:	Press Up key (2) for LOCK mode ON/OFF . The lock mode is particularly useful for continuous monitoring of temperatures for up to 60 minutes.
	Press Down key (3) for °C or °F
In all modes: First hold on the Meas. key (6)	and press Up key (2) for backlight function ON/OFF . 
	and press Down key (3) for laser function ON/OFF . 

CAUTION

1. WHEN DEVICE IS IN USE, DO NOT LOOK DIRECTLY INTO THE LASER BEAM—PERMANENT EYE DAMAGE MAY RESULT.
2. USE EXTREME CAUTION WHEN OPERATING THE LASER.
3. NEVER POINT THE DEVICE TOWARDS ANYONE'S EYES.

LCD ERROR MESSAGES

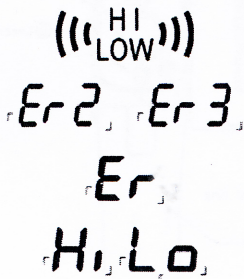

'Hi' or 'Lo' is displayed when the temperature being measured is outside of the settings of HAL and LAL.

'Er2' is displayed when the thermometer is exposed to rapid changes in the ambient temperature. 'Er3' is displayed when the ambient temperature exceeds 0°C (32°F) or +50°C (122°F). The thermometer should be allowed plenty of time (minimum 30 minutes) to stabilize to the working/room temperature.

Error 5-8, for all other error messages it is necessary to reset the thermometer. To reset it, turn the instrument off, remove the battery and wait for a minimum of one minute, reinsert the battery and turn on.

'Hi' or 'Lo' is displayed when the temperature being measured is outside of the measurement range.

BATTERIES

The thermometer incorporates visual low battery indication as follows:



'Battery OK': measurements are possible



'Battery Low': battery needs to be replaced, measurements are still possible



'Battery Exhausted': measurements are not possible

⚠ When the 'Low Battery' icon indicates the battery is low, the battery should be replaced immediately with AAA, 1.5V batteries. Please note: It is important to turn the instrument off before replacing the battery otherwise the thermometer may malfunction.

SPECIFICATION

Item	Non-contact Infrared Scan function	Thermocouple Probe Scan function (K type; probe not included.)
Measurement Range	-60 to +550 °C (-76 to +1022 °F)	-64 to +1400 °C (-83.2 to +2552 °F)
Operating Range	0 to +50 °C (32 to +122 °F)	

Emissivity Range	0.95 default – adjustable 0.1 to 1 step .01
Resolution	0.1 °C/0.1 °F at -83.2 to 999.9 (°C/°F), otherwise 1 °C/1 °F
Response Time	1sec
Distance:Spot	12:1

Technical service

Web Resources

Visit the VWR's website at www.vwr.com for:

- Français, Español, Deutsch, Italiano, Português version
- Complete technical service contact information
- Access to VWR's Online Catalogue, and information about accessories and related products
- Additional product information and special offers

Contact us for information or technical assistance contact your local VWR representative or visit. www.vwr.com.

Warranty

VWR International warrants that this product will be free from defects in material and workmanship for a period of two (2) years from date of purchase. If a defect is present, VWR will, at its option, repair, replace, or refund the purchase price of this product at no charge to you, provided it is returned during the warranty period. This warranty does not apply if the product has been damaged by accident, abuse, misuse, or misapplication, or from ordinary wear and tear. For your protection, items being returned must be insured against possible damage or loss. This warranty shall be limited to the replacement of defective products. IT IS EXPRESSLY AGREED THAT THIS WARRANTY WILL BE IN LIEU OF ALL WARRANTIES OF FITNESS AND IN LIEU OF THE WARRANTY OF MERCHANTABILITY.

GO TO VWR.COM FOR THE LATEST NEWS, SPECIAL OFFERS AND DETAILS OF YOUR LOCAL VWR DISTRIBUTOR