VWR WEIGHING GRANIT TABLE

INSTRUCTION MANUAL

VWR Balance table made in granit with vibration dampers and adjustable feets

W X D X H (mm) : 900 X 600 X 760

European Catalogue Number(s):

Euro Plug: 611-2549

Version: 1

Issued: 06/08/2012

Legal Address of Manufacturer

United States
VWR International
1310 Goshen Parkway
West Chester, PA 19380
800-932-5000
http://www.vwr.com


Europe
VWR International bvba
Researchpark Haasrode 2020
Geldenaaksebaan 464
B-3001 Leuven
+ 32 16 385011
http://be.vwr.com

Country of origin : FRANCE

Table of Contents

Package Contents
Instructions for use
Technical Service
Warranty
Disposal

Package Contents

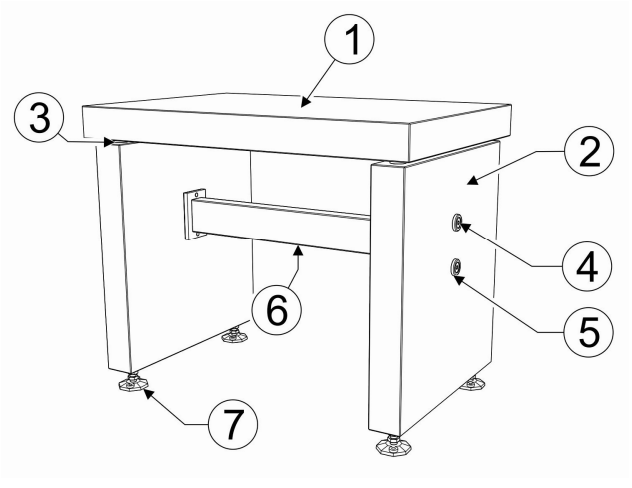Description	Quantity	Rep.	
Top plate	1	1	
Side Plate	2	2	
Anti vibration block	4	3	
Screw FHC M12 x 100	4	4	
Ring	4	5	
Cross bar	1	6	
Adjustable feet	4	7	



SET-UP TIME : ONE HOUR



NUMBER OF PERSONS : 2




FOR GENERAL CARE, A SIMPLE DUSTER IS NEEDED FOR THE MAINTENANCE




TOOLS NEEDED FOR ASSEMBLING



Hex key of 8 [1]

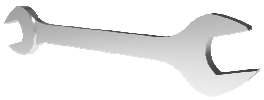

Open-end spanner of 24 and 16

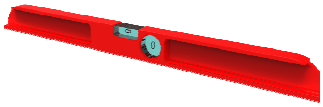

Spirit Level

[1] Dimensions in mm


Instructions for use

Important Precautions :

Be sure to follow the safety instructions.


Wear safety shoes with steel toe caps.

Never attempt to set up the table alone, two people are always required.

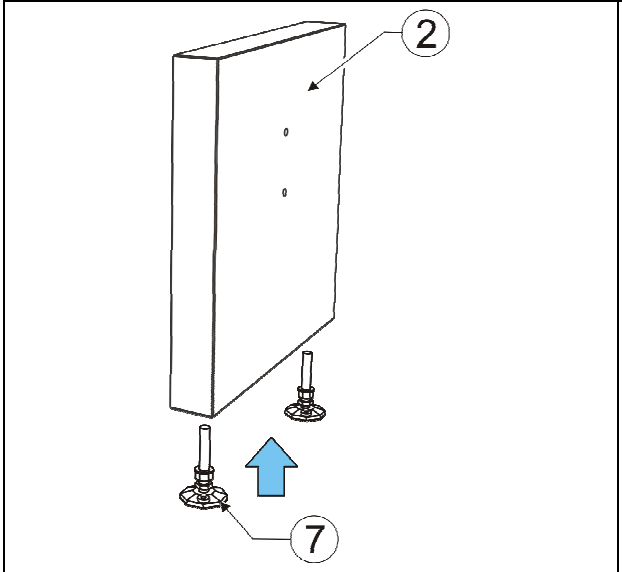

Step 1 :

Lift a side plate(2), and slide the 2 adjustable feet into the holes (7)

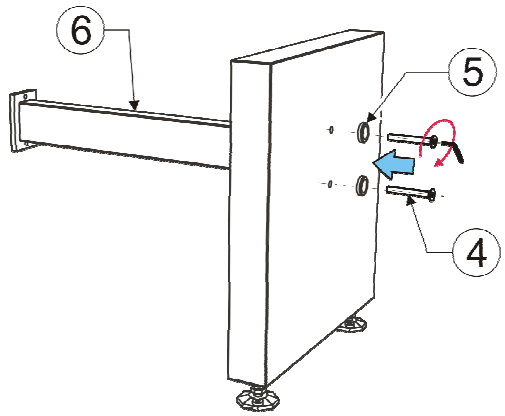


Step 2 :

Put the rings (5) and the M12 x 100 screws (4) together. Set up a ring/screw assembly in one of the side hole.




Position the crossbar (6) using the hex key and screw without forcing.

Next, screw the second set.

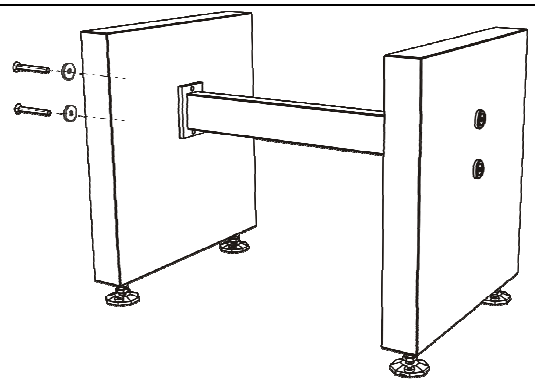

Step 3 :

Do the same as in step 1 and 2.

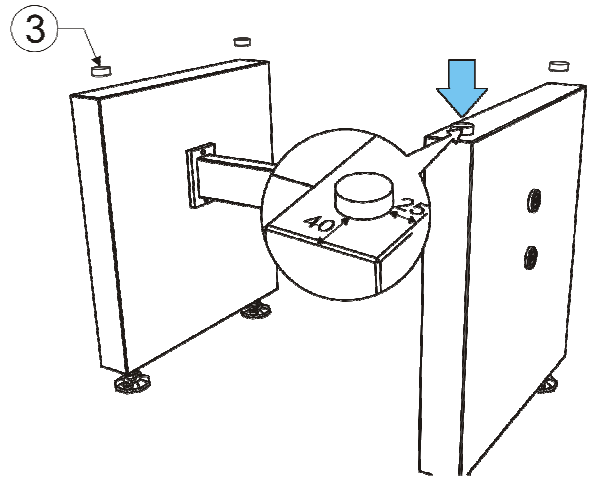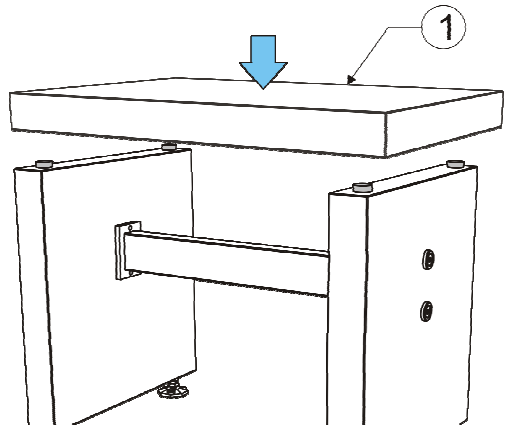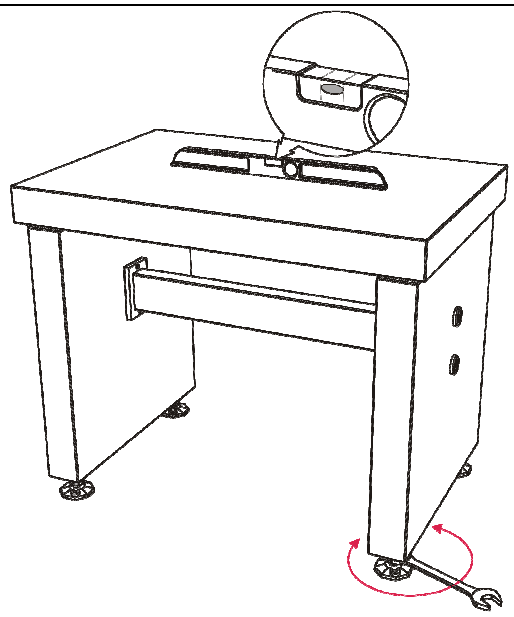	**Step 4 :** Position the 4 anti vibration blocks on the top (3) Remove the protective film from each block Stick them and respect the dimensions
	**Step 5 :** Put the top plate (1) Center the plate relative to the side plates (edges of about 2cm).
	**Step 6 :** Check the level of the table If necessary use the spirit level and the open end spanner

Technical service

Web Resources

Visit the VWR's website at **www.vwr.com** for:

• Complete technical service contact information

• Access to VWR's Online Catalogue, and information about accessories and related products

• Additional product information and special offers

Contact us For information or technical assistance contact your local VWR representative or visit. **www.vwr.com**.

Warranty

VWR International warrants that this product will be free from defects in material and workmanship for a period of two (2) years from date of purchase. If a defect is present, VWR will, at its option, repair, replace, or refund the purchase price of this product at no charge to you, provided it is returned during the warranty period. This warranty does not apply if the product has been damaged by accident, abuse, misuse, or misapplication, or from ordinary wear and tear.

For your protection, items being returned must be insured against possible damage or loss. This warranty shall be limited to the replacement of defective products. IT IS EXPRESSLY AGREED THAT THIS WARRANTY WILL BE IN LIEU OF ALL WARRANTIES OF FITNESS AND IN LIEU OF THE WARRANTY OF MERCHANTABILITY.

Equipment disposal



This equipment is marked with the crossed out wheeled bin symbol to indicate that this equipment must not be disposed of with unsorted waste.
Instead it's your responsibility to correctly dispose of your equipment at lifecycle -end by handling it over to an authorized facility for separate collection and recycling. It's also your responsibility to decontaminate the equipment in case of biological, chemical and/or radiological contamination, so as to protect from health hazards the persons involved in the disposal and recycling of the equipment.
For more information about where you can drop off your waste of equipment, please contact your local dealer from whom you originally purchased this equipment.

By doing so, you will help to conserve natural and environmental resources and you will ensure that your equipment is recycled in a manner that protects human health.

Thank you

Local VWR offices in Europe and Asia Pacific

Austria
VWR International GmbH
Graumanngasse 7
1150 Wien
Tel.: 01 97 002 0
Fax: 01 97 002 600
E-mail: info@at.vwr.com

Belgium
VWR International bvba
Researchpark Haasrode 2020
Geldenaaksebaan 464
3001 Leuven
Tel.: 016 385 011
Fax: 016 385 385
E-mail:
customerservice@be.vwr.com

China
VWR International China Co., Ltd
Suite 1802 - 1803,
Xing Ye Bank Mansion, No 168,
168 Jiangning Road
Shanghai 200041, China
Tel.: +86- 21 521 388 22
Fax: +86- 21 521 33 933
E-mail: sales_china@vwr.com

Denmark
VWR - Bie & Berntsen
Transformervej 8
2730 Herlev
Tel.: 43 86 87 88
Fax: 43 86 87 90
E-mail: info@dk.vwr.com

Finland
VWR International Oy
Pihatörmä 1 C 1
FI - 02240 Espoo
Tel.: +358 9 80 45 51
Fax: +358 9 80 45 52 00
E-mail: info@fi.vwr.com

France
VWR International S.A.S.
Le Périgares – Bâtiment B
201, rue Carnot
94126 Fontenay-sous-Bois cedex
Tel.: 0 825 02 30 30 (0,15 EUR TTC/min)
Fax: 0 825 02 30 35 (0,15 EUR TTC/min)
E-mail: info@fr.vwr.com

Germany
VWR International GmbH
Hilpertstrasse 20a
D - 64295 Darmstadt
Tel.: 0180 570 20 00*
Fax: 0180 570 22 22*
E-mail: info@de.vwr.com
*0,14 €/Min. aus d. dt. Festnetz, Mobilfunk max. 0,42 €/Min.

Hungary
Spektrum-3D Ltd.
A VWR International Company
Simon László u. 4.
4034 Debrecen
Tel.: (52) 521-131
Fax: (52) 470-069
E-mail: info@spektrum-3d.hu

India
VWR Lab Products Pte Ltd
2nd Floor, Front Wing, 135/12,
Brigade Towers
Brigade Road
Bangaluru 560025 India
Tel: +91-2522-277876 (Mumbai)
Tel: +91-80-41117124 (Bangalore)
Fax +91-80-41117120
E-mail: vwr_india@vwr.com

Ireland
VWR International Ltd
Orion Business Campus
Northwest Business Park
Ballycoolin
Dublin 15
Tel.: 01 88 22 222
Fax: 01 88 22 333
e-mail sales@ie.vwr.com

Northern Ireland
VWR International Ltd
A10 Harbour Court, 7 Heron Rd
Sydenham Business Park
Belfast BT3 9HB
Tel.: 028 9058 5800
Fax: 028 9080 7812
e-mail sales@ie.vwr.com

Italy
VWR International s.r.l.
Via Stephenson 94
20157 Milano (MI)
Tel.: 02 332 03 11
Fax: 800 152 999
E-mail: info@it.vwr.com

The Netherlands
VWR International B.V.
Postbus 8198
1005 AD Amsterdam
Tel.: 020 4808 400
Fax: 020 4808 480
E-mail: info@nl.vwr.com

Norway
VWR International AS
Haavard Martinsensvei 30
0978 Oslo
Tel.: 02290
Fax: 815 00 940
E-mail: info@no.vwr.com

Portugal
VWR International - Material de Laboratório, Lda
Edifício Neopark
Av. Tomás Ribeiro, 43- 3 D
2790-221 Carnaxide
Tel.: 21 3600 770
Fax: 21 3600 798/9
E-mail: info@pt.vwr.com

Singapore
VWR Singapore Pte Ltd
18 Gul Drive
Singapore 629468
Tel: +65 6505 0760
Fax: +65 6264 3780
E-mail: sales@sg.vwr.com

Spain
VWR International Eurolab S.L.
Ronda Can Fatjó, nº 11
Edifici Tecnopark, 3
Parc Tecnològic del Vallés
08290 Cerdanyola del Vallés (Barcelona)
Tel.: 902 222 897
Fax: 902 430 657
E-mail: info@es.vwr.com

Sweden
VWR International AB
Fagerstagatan 18a
163 94 Stockholm
Tel.: 08 621 34 00
Fax: 08 621 34 66
E-mail: info@se.vwr.com

Switzerland
VWR International AG
Lerzenstrasse 16/18
8953 Dietikon
Tel.: 044 745 13 13
Fax: 044 745 13 10
E-mail: info@ch.vwr.com

UK
VWR International Ltd
Customer Service Centre
Hunter Boulevard
Magna Park
Lutterworth
Leicestershire
LE17 4XN
Tel.: 0800 22 33 44
Fax: 01455 55 85 86
E-mail: uksales@uk.vwr.com