Instruction manual

VWR® UV-3100PC Spectrophotometer

EU cat. no 634-6042

Version 02 Issue 10/2021

Legal address of Manufacturer

Europe
VWR International bv,
Researchpark Haasrode 2020
Geldenaaksebaan 464
B-3001 Leuven
+3216385011
be.vwr.com

UK Importer
VWR International Ltd
Hunter Boulevard, Magna Park
Lutterworth, Leicestershire, LE17 4XN
uk.vwr.com

Made in China

Table of Contents

Part 1: Spectrophotometer 1

Safety Information 1

Package Contents 2

Unpacking 2

Installation 2

Symbols and Conventions 3

Specifications 3

Overview 4

Description of Buttons and Switches 5

Operational Keys 6

Getting Started 8

Important Guidelines 9

General Operating Instructions 10

Operation 12

1. Enter into Quantitative Mode 14

2. Set Unit 14

3. Set up Standard Curve or load the stored curves 15

4. Return the sample measurement interface 17

5. Blank 18

6. Measure Samples 18

7. Print the Test Results 18

Troubleshooting 35

Repair and Maintenance 36

Accessories and Spare Parts 43

Technical Service 44

Warranty 44

Equipment Disposal 45

Part 2: Software 47

Functions 47

Installation 49

Introduction 52
Operation 56
Appendix 1 78

Part 1: Spectrophotometer

Safety Information

Please follow the guidelines below, and read this manual in its entirety to ensure safe operation of the unit.

VWR recommends against the use of UV-3100PC Spectrophotometer.



- Do not open the device.
- Disconnect the device from the mains supply before carrying out maintenance work or changing the fuses.
- The inside of the device is a high-voltage area Danger!
- Do not use the device if it is damaged, especially if the main power cable is in any way damaged or defective.
- Repairs may only be carried out by the service technicians from your local VWR office and authorized contractual partners.
- The device must be connected to a power outlet that has a protective ground connection.
- If the equipment is used in a manner not specified by the manufacturer, the protection provided by the equipment may be impaired.



- Do not allow any liquid to enter into the device.
- Do not operate the device in a hazardous location or potentially explosive environment.

Package Contents

Description	Quantity
Spectrophotometer	1PC
Glass Cuvette	4PCS
Quartz Cuvette	2PCS
USB Disk	1PC
Power Cord (Euro Plug)	1PC
Power Cord (UK Plug)	1PC
Power Cord (CH Plug)	1PC
USB Cable	1PC
CD-ROM	1PC
Quick Manual	1PC
Instruction Manual	1PC
Dust Cover	1PC

Unpacking

Open the package, according to carefully check the packaging packing list items, if found inside the packaging are missing or damaged items please contact your local VWR office and authorized contractual partners.

Installation

Placement

Place the instrument on the stable table carefully.

Install printer (Optional)

Check to confirm instrument power switch is turned off, connect the printer's data cable to the

Instrument's parallel port.

Link the power cord

Check to confirm instrument power switch is turned off, the power cord plug into two separate power interface and power supply socket apparatus.

Symbols and Conventions

The following chart is an illustrated glossary of the symbols that are used in this manual.

	**CAUTION** This symbol indicates a potential risk and alerts you to proceed with caution
	**CAUTION** This symbol indicates the presence of high voltage and warns the user to proceed with caution
	**CAUTION** This symbol indicates risks associated with hot surfaces

Specifications

Optical System	Single Beam
Wavelength Range	190-1100nm
Band Width	2nm
Stray Light	≤0.05%T @ 220nm & 360nm
Photometric Range	0 to 200%T,-0.3 to 3.0A, 0 to 9999C
Wavelength Accuracy	±0.5nm
Photometric Accuracy	±0.5%T or ±0.005A@1A
Baseline	0.002A (200 to 1000nm)

Stability	0.002A/h @ 500nm
Memory	32K(internal), Unlimited(USB Disk)
Language	English, French, German, Spanish
Display	320×240 Dots Matrix LCD
Interface	USB, Parallel
Measuring Procedure	Photometry, Quantitation, Wavelength Scan, Kinetics, DNA/Protein, Multi-wavelength
Power Supply	AC 110/220V, 50/60Hz
Dimension	490×360×240mm
Weight	14kg
Work Environment	15 to 35℃, 15 to 70% relative humidity
Store Environment	-10 to 50℃, 15 to 70% relative humidity

This instrument is compliant to the European Directives on
Low Voltage Directive 2014/35/EU
Electromagnetic compatibility 2014/30/EU
Restriction on use of Hazardous Substances RoHS 2011/65/EU
and theirs amandements.

Overview


UV-3100PC Spectrophotometer used in Chemistry, Pharmaceuticals, Biochemical, Metallurgy, Light Industry, Textile, Material, Environments, Medical, Education and some other fields for Quality Control laboratories.

Description of Buttons and Switches


Front View






Rear View

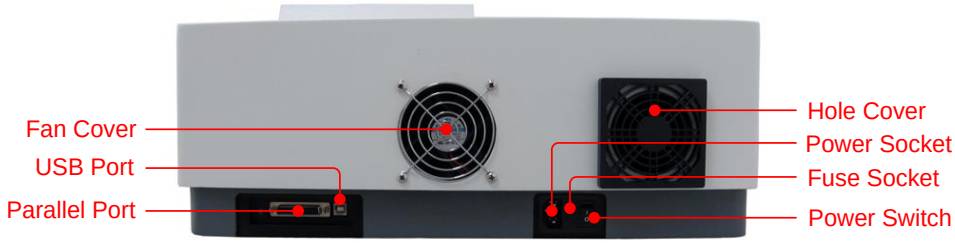

Operational Keys

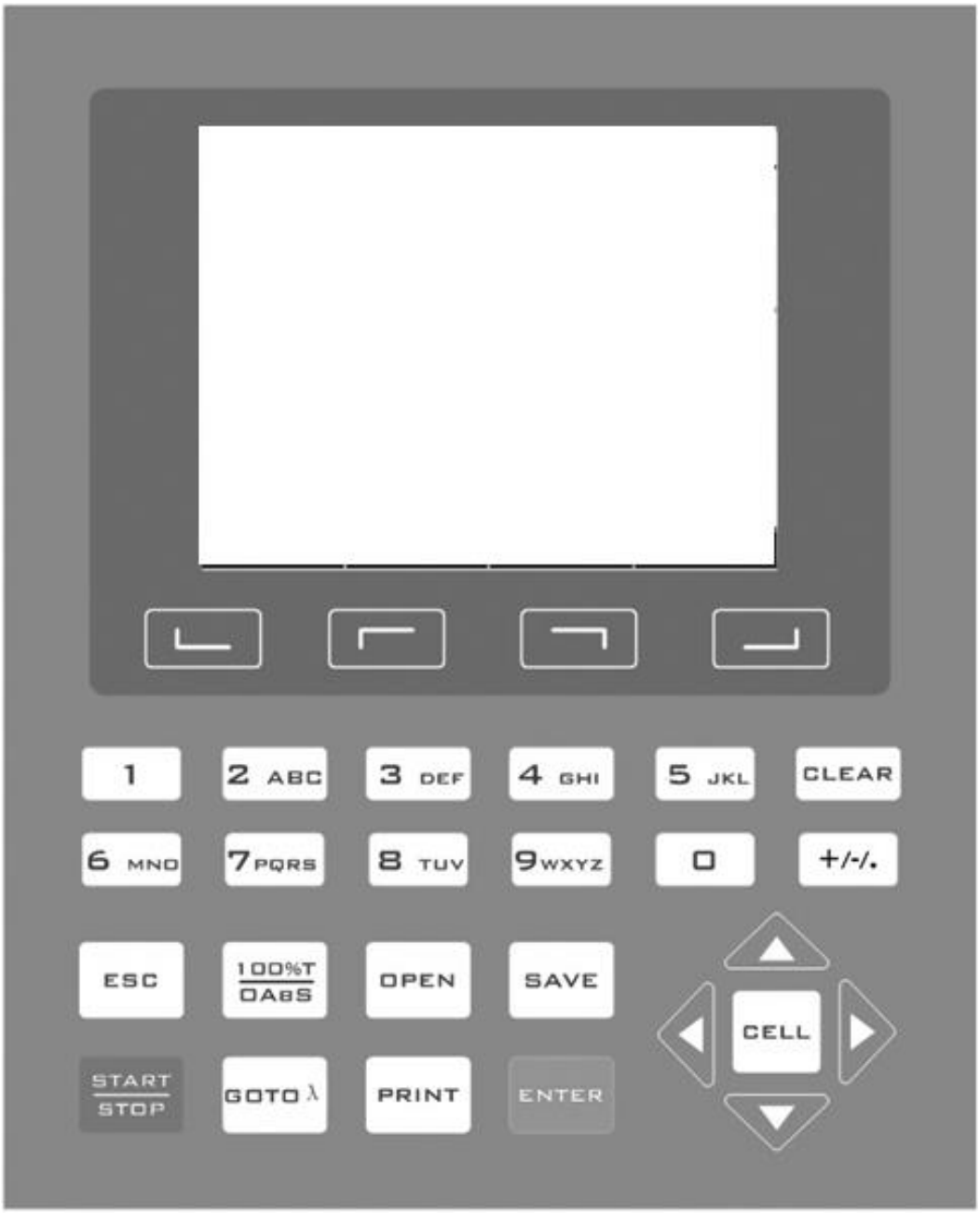

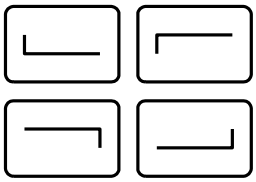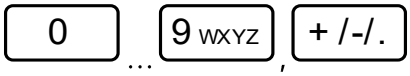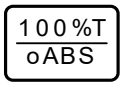Key	Functions
	Function Keys: Functions on-screen prompts
	Numeric Keys: Enter numbers and letters
	CLEAR Key: Delete the input value or stored data
	ESC Key: Return to previous Interface
	100%T/0Abs Key: Blank
	OPEN Key: Open files stored in internal memory
	SAVE Key: Save files to internal memory
	START/STOP Key: Start/Stop testing
	GOTO λ Key: Set wavelength

PRINT

PRINT Key: Print measuring result

ENTER

ENTER Key: Confirm operation

CELL

CELL Key: Select/Deselect Auto-cell Holder

RIGHT, LEFT Keys: Search peak/valley and set X scale


UP, DOWN Keys: Scroll menu/data and set Y scale

Getting Started


The following chart describes the basic operation of the instrument.

Turn On and Self-check


Switch on the power. Then the instrument begins to self-check and 15 minutes' warm up. Self-check includes the following steps:

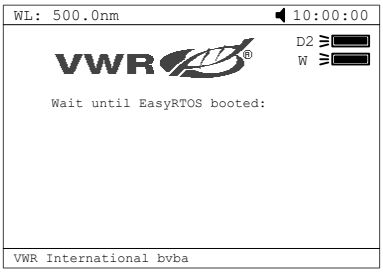


Turn on lamps →RAM Check→ Start RTOS kernel→ Initialize Comm. Port → Initialize Printer → Initialize AD → System position → Warm up.

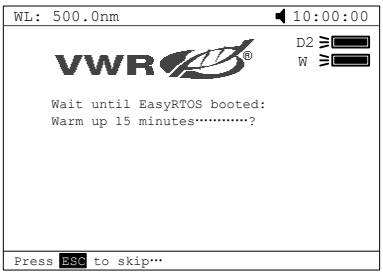

After warm up, the instrument will ask the user to re-calibrate the system. Users can decide if they need to re-calibrate the system or not. After this step, the instrument can work normally.

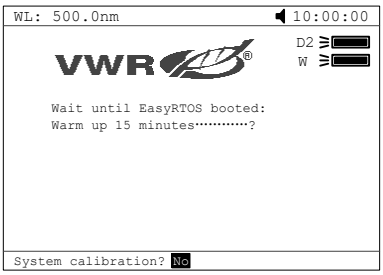

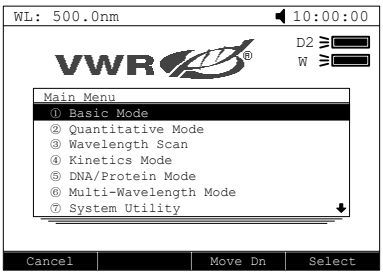

Important Guidelines

- Reagents and dilution buffers can cause cauterization and other damage to health.
- Samples (nucleic acids, proteins, bacteria cultures) can be infectious and cause serious damage to health.
- During sample preparation, measuring procedures and maintenance and cleaning work, observe all local laboratory safety precautions (e.g. wear protective clothing and

gloves, use of disinfectant) regarding the handling of sample material.

- Dispose of measuring solutions and cleaning and disinfectant materials in accordance with the relevant local laboratory regulations.

General Operating Instructions

Select Application

Main menu, press numeric keys or use the key ▲, ▼ to choose corresponding menu, then press [ENTER] to enter into.

Set Wavelength

Press [GOTO ∧] to set wavelength, use numeric keys to input the values, press [ENTER] to confirm and go to the point you set, then do blank automatically.

Set Parameters

In different application, press function key to set parameters, press ▲, ▼ to choose or input the values by numeric keys, press [ENTER] to enter into, press [ESC] to return.

Set Auto-cell Holder

Press [CELL] to active the auto cell holder and press the numeric key (1-8)to make corresponding cell position at the light path. Press [CELL] again to inactive the auto cell holder.

Delete the Input Value

Press ◀ to delete a character, press [CLEAR] to delete all the characters.

Delete the Test Results and Stored data

Press [CLEAR] to delete the test result or stored data. If the USB Disk was installed, the data will delete in the USB Disk, or the data will delete in the internal memory.

Blank

Put the Reference in the light path, press [100%T / oABS] to do blank.

Measure Samples

Put the samples in the light path, press [START / STOP] to measure.

Print the Test Results

Press [PRINT] to print the test results.

Store the Test Results

Press [SAVE] to store the test results, input the file name by numeric keys and press [ENTER] to save. If the USB Disk was installed, the data will save in the USB Disk, or the data will save in the internal memory.

Load the Stored File

In the test interface, press [OPEN] to go into file selecting interface, press ▲, ▼ to choose the file you want, press [ENTER] to open. If the USB Disk was installed, the data will open from the USB Disk, or the data will open from the internal memory.

Operation

Self-check

Remove all the blocks in the light path and close the lid of the compartment; Switch on the power supply to begin the self-test.

Warm up

After self-test, the instrument goes into pre-warm state. For accurate test, at least 30 minutes of warm up is required.


Check the cuvettes

The cuvettes must be clear and there's no remains of the samples on the surface of it. Only Silicon (Quartz) cuvettes are permitted to be used in the range of UV area.

Basic Mode

1. Enter into Basic Mode

Main menu, press numeric key [1] or ▲,▼ to choose "Basic Mode", then press [ENTER] to go into.

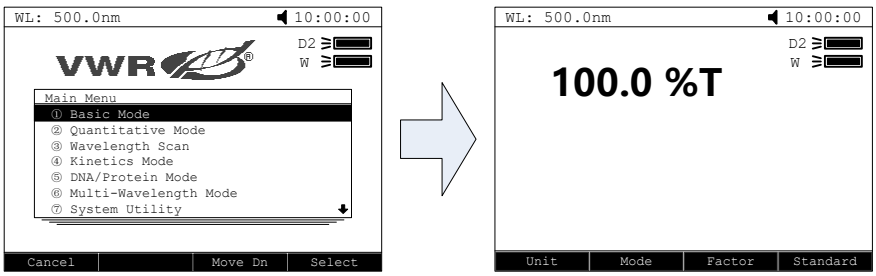

2. Set Photometric Mode

Press [ ] to set photometric mode. Press ▲,▼ to choose "Abs.", "T%" or "Conc./Factor"

and press ENTER to confirm. If users choose "Abs." or "T%", please go to step 5 directly.

3. Set Concentration Unit

Press ⌊— to set concentration unit. Press ▲, ▼ to choose unit followed with ENTER pressed to confirm. You can also choose "other" to input the self defined unit.

4. Set "Factor" or "Standard"

Two methods are under your choice:

Method 1：Input Factor F

Press ¬— to set F. Input the value of F by numeric keypad, press ENTER to confirm. Then the F value would display on the screen.

Method 2：Standards Mark

Put the reference sample in the light path and calibrate 100%T/0Abs; Put the standard sample in the light path，press —⌋ to start the mark. Input the concentration value of the standard and press ENTER to confirm, then it displays on the screen.

5. Set Wavelength

Press GOTO ∧ to set wavelength, input the value by the numeric keypad followed with ENTER pressed to confirm.

6. Blank

Put the Reference in the light path and press 100%T/oABS to do blank.

7. Measurement samples

Put the sample to be measured in the light path, then the result displays on the screen

automatically.

8. Print the Test Results

Press [PRINT] to print the Test Results.

```
Basic Mode Test Report


Wavelength:    500.0nm
Result:        0.000 Abs
Date & Time:  mm-dd-yyyy, hh:mm:ss

Model:         UV-3100PC
SN:            UEFXXXXXXX
Version:       A1.176
VWR International bvba.
```

Quantitative Mode

1. Enter into Quantitative Mode

Main menu, press [2 ABC] or press ▲,▼ to choose "Quantitative Mode" followed with [ENTER] pressed to confirm.

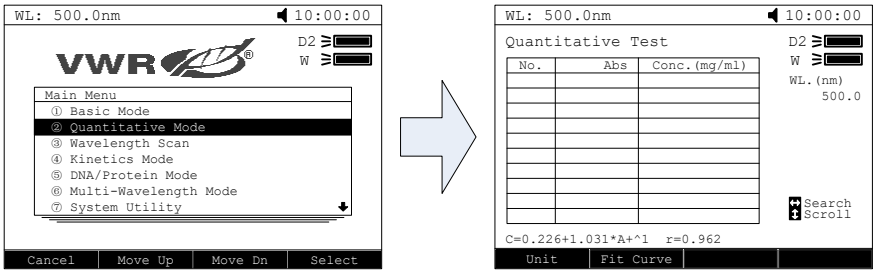


2. Set Unit

Press ⌙ to set concentration unit, press ▲, ▼ to choose and press ENTER to confirm.

3. Set up Standard Curve or load the stored curves

Press ⌐ to go into set up interface, 2 methods are under your choice.

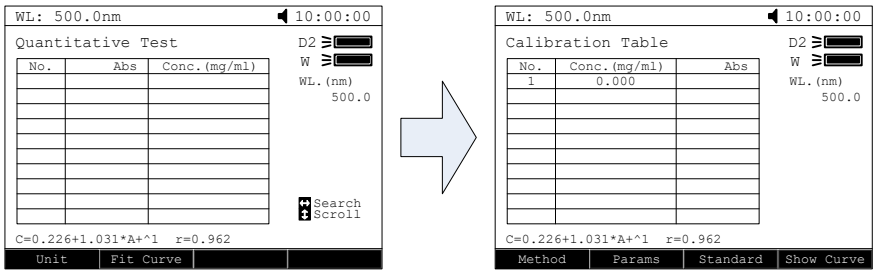

Set up Standard Curve:

Method 1: Input Regression Equation

1) **Set Fit Curve Method.** Press ⌙ to set Fit method, use ▲, ▼ to choose the method and press ENTER to confirm.

2) **Set Wavelength.** Press GOTO ∧ to set wavelength. Use ▲, ▼ to choose measure method, then press ENTER to confirm. Input the wavelength value you need and press ENTER to confirm.

3) **Input the Factor of the Regression Equation.** Press ⌐ and input the factors, press ENTER to confirm.


Method 2: Use Standard Samples

1) **Set Fit Method.** Press [⌐] to set fit method, press [▲], [▼] to choose fit method, then press [ENTER] to confirm.


2) **Set Wavelength.** Press [GOTO ∧] to go into wavelength setting interface, press [▲], [▼] to select measure method and press [ENTER] to confirm. Input the value of the wavelength and press [ENTER] to confirm.

3) **Blank.** Put the Reference Sample in the light path and press [100%T oABS] to do blank.

4) **Setup Standard Samples.** Press [¬] to setup standard, input the concentrations of corresponding standard samples according the indication and press [ENTER] to confirm. Users can use [▲], [▼] to choose the value you just input and press [CLEAR] to delete, then input a new value, press [ENTER] to confirm. Press [ESC] to cancel after all the input.

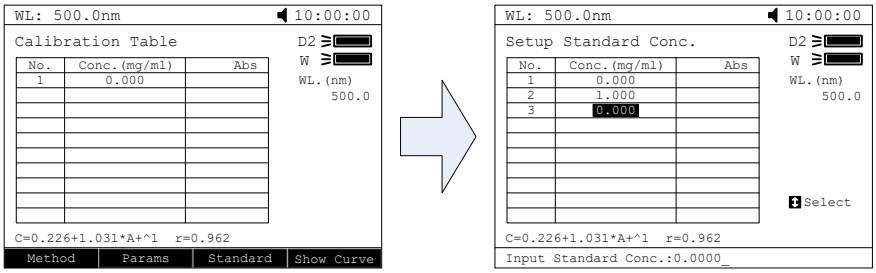


5) **Calibrate Standard Samples.** Put the corresponding standard samples in the light path as the screen indicates and press [START/STOP] to measure. Then the Abs. value would appear in the corresponding table.

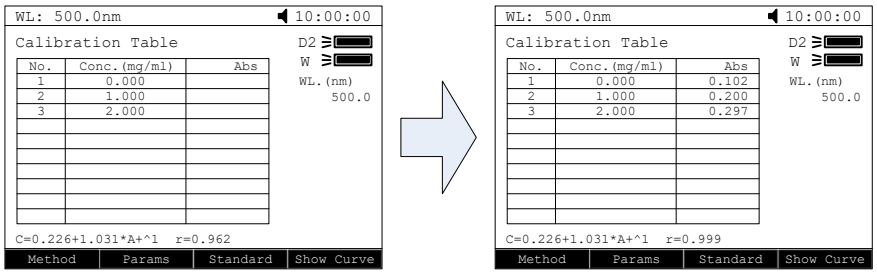


Load the Stored Curves

In the "Calibration Table" interface, press [OPEN] to go into files select interface. Use ▲, ▼ to select the curve you need and press [ENTER] to load. Users can press [↵] to view the curve, press [ESC] to cancel.

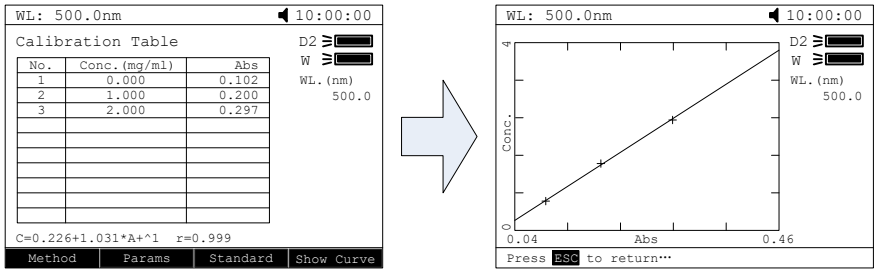

4. Return the sample measurement interface

In the "Calibration Table" interface, press [ESC] to return the Quantitative Test interface.

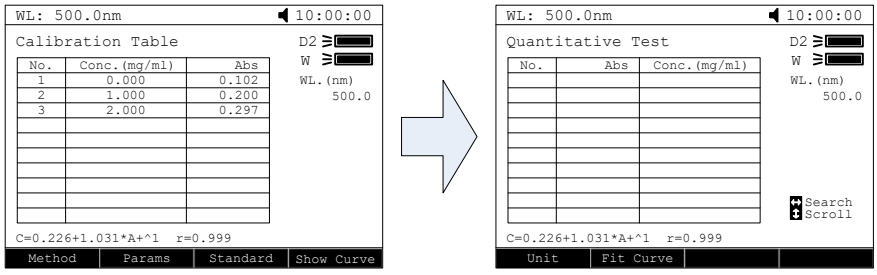


5. Blank

Put the Reference Sample in the light path, press 100%T / oABS to do blank.

6. Measure Samples

Place the sample to be tested in the light path, press START / STOP to measure. Then the test result will display in the data sheet. Repeat this step to finish measuring all the samples.

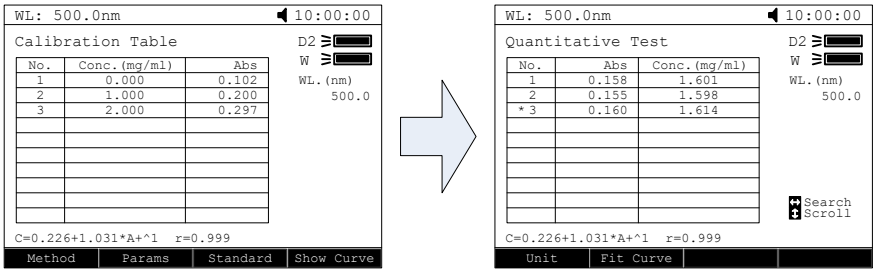

7. Print the Test Results

Press PRINT to print the Test Results.

Quantitative Test Report

File Name: Abc
Date & Time: mm-dd-yyyy, hh:mm:ss


No.	500.0nm	Abs(eff)	C(mg/L)
1	0.212	0.212	0.212
2	0.210	0.210	0.210
3	0.209	0.209	0.209

Fitting Params:
C=1.000*A^1
r=1.000
Model: UV-3100PC
SN: UEFXXXXXXX
Version: A1.176
VWR International bvba.

Wavelength Scan

1. Enter into Wavelength Scan

Main menu, press numeric key 3 DEF or ▲, ▼ to choose "Wavelength Scan" and press ENTER to enter.

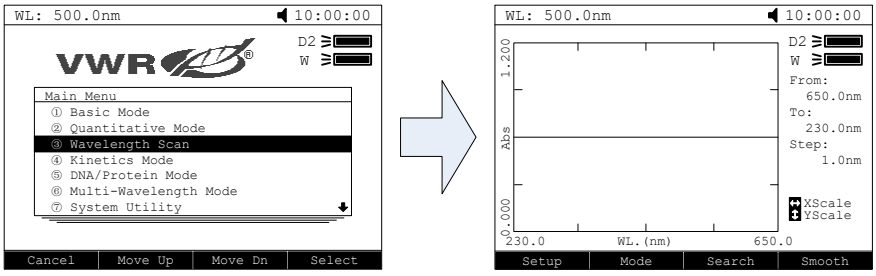

2. Parameters Setup

Press ↵ to set parameters, set "Scan From" , "scan to", "scan step" and "scan speed", press

ENTER to confirm.

3. Set Photometric Mode


Press ⌐¬ to set photometric mode，choose “T%”, “Abs.” or “E” and press ENTER to confirm.

4. Scan Baseline

Put Reference Sample in the light path, press 100%T/oABS to scan the baseline, press ESC to cancel.


5. Scan Samples

Put the samples in the light path，press START/STOP to scan the sample, press ESC to cancel.

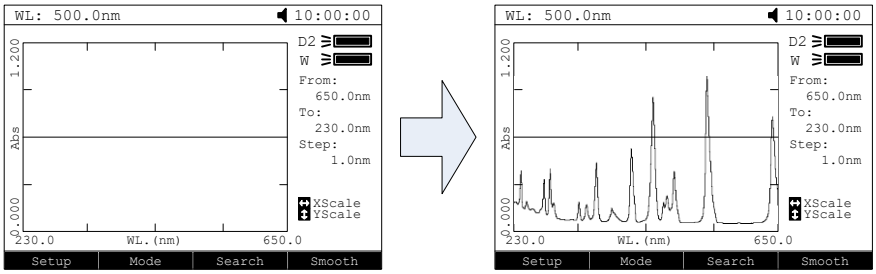

6. Search Peaks

After scanned, press ⌐¬ to go into peak search mode. Press ⌐¬ to set peak height, input the peak height and press ENTER to confirm. Press ◀, ▶ to display the value of every wavelength point. Press ▲, ▼ to display the value of every peak.

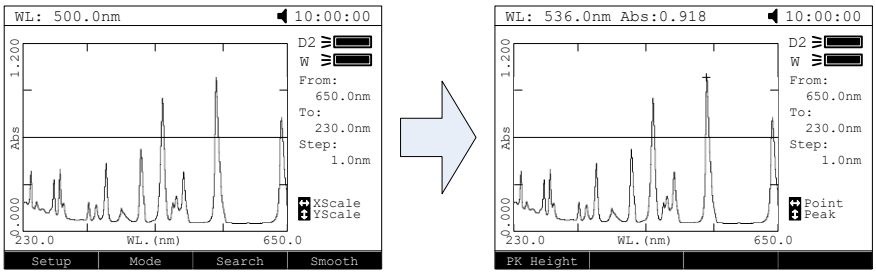

7. Smooth the Curve

After scanned, if there are many burrs on it, press  to smooth the curve.


8. Print the Test Results

After scanned, press 

to print the Test Results.

Wavelength Scan Test Report
File Name: Abc
Date & Time: mm-dd-yyyy, hh:mm:ss

Scan From: 680.0nm
Scan From: 230.0nm
Scan Step: 1.0nm
Peak Height: 0.030Abs
Model: UV-3100PC
SN: UEFXXXXXXX
Version: A1.176
VWR International bvba.

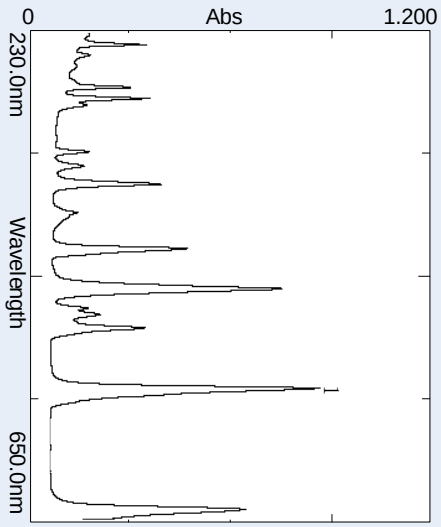


Peak list

No.	Wavelength(nm)	Abs	T%
1	536.0	0.918	12.1

Kinetics

1. Enter into Kinetics

Main menu, press 4 GHI or ▲, ▼ to select "Kinetics Mode" and press ENTER to confirm.

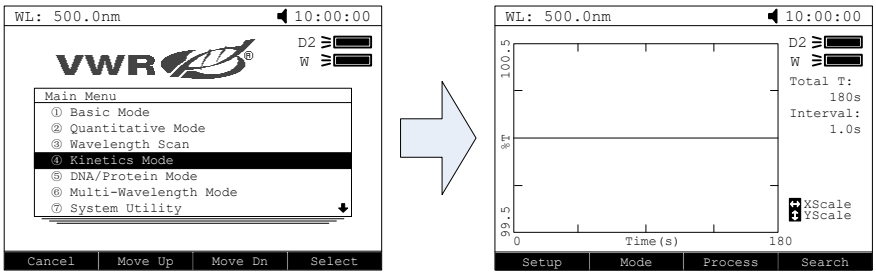

2. Setup Parameters

Press [⌐] to set parameters, input the corresponding values of "Total Time", "Delay Time" and "Time Intervals" according the screen indicates. Press [ENTER] to confirm.

3. Set Photometric Mode

Press [⌐] to set photometric mode, choose "T%" or "Abs." and press [ENTER] to confirm.

4. Set Wavelength


Press [GOTO ∧] to set wavelength, input the value of the wavelength by numeric keypad and press [ENTER] to confirm.

5. Blank

Put the Reference Sample in the light path, press [100%T / oABS] to do blank.


6. Measure Samples

Put the sample in the light path and press [START / STOP] to begin the test, press [ESC] to cancel.

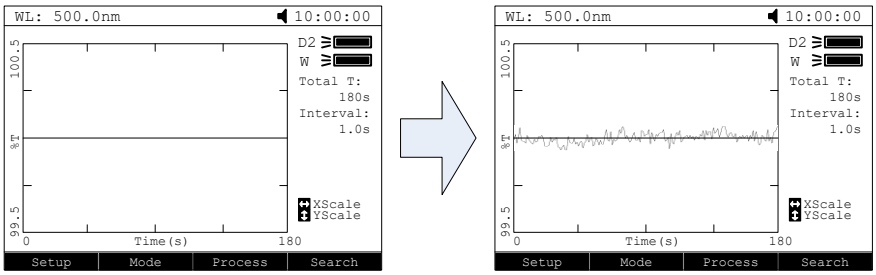


7. Calculate Response Rate

After time scanned, if users want to calculate the response rate of some period, you can press  to go into "Process" interface. Input the values of "Begin Time", "End Time" and "Factor" separately and press  to confirm. Then the value of "I.U." would display on the screen.

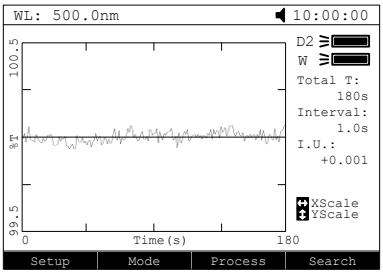


8. Search Peaks

After scan finished, press  to go into search mode. Press  to search the value of every point.

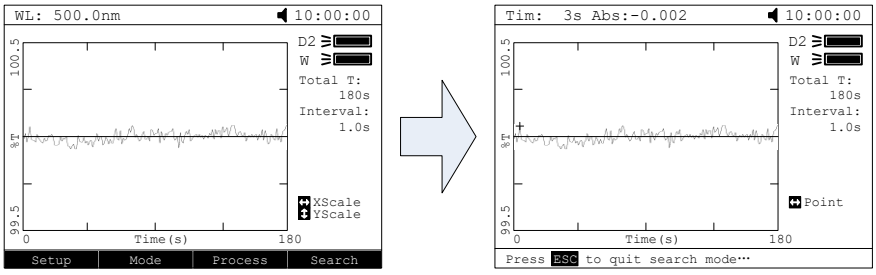

9. Print the Test Results


Press PRINT to print the Test Results.

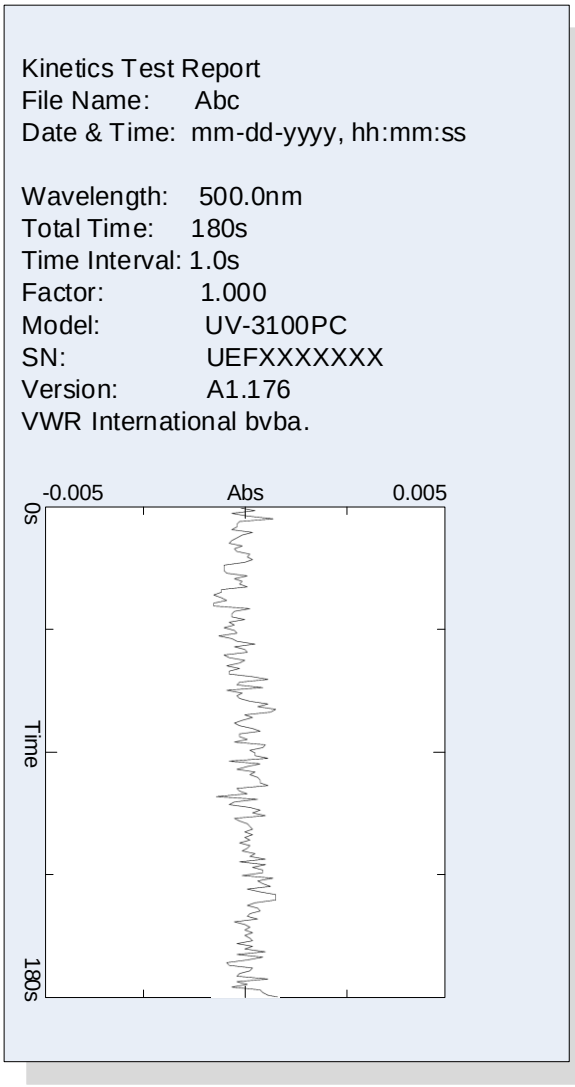

DNA/Protein Mode

1. Enter into DNA/Protein Mode

Main menu, press 5 JKL or ▲, ▼ to choose "DNA/Protein Mode" and press ENTER to confirm.

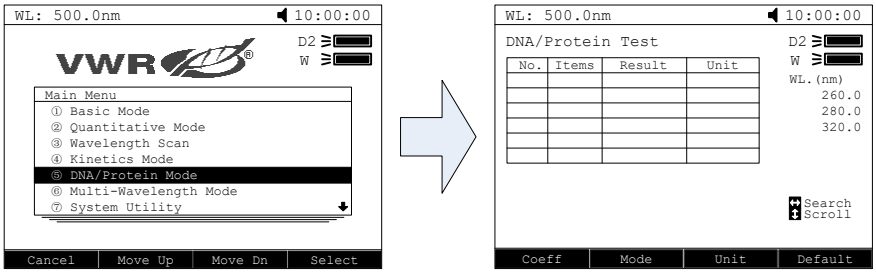

2. Setup Parameters

Press [ ] to set coefficients, input all the values of f1 to f4 by numeric keypad according the indication and press [ENTER] to confirm.

3. Choose Measure Method

Press [ ] to set method. Press ▲, ▼ to choose "Absorbance Difference 1" or "Absorbance Difference 2" followed with [ENTER] pressed to confirm. If users don't want to measure reference, use ▲, ▼ to choose "No" followed with [ENTER] pressed to confirm the choice.

4. Set Concentration Unit

Press [ ] to set concentration unit. Use ▲, ▼ to select unit and press [ENTER] to confirm.


5. Blank

Place the Reference Sample in the light path, press [100%T / oABS] to do blank.

6. Measure Samples

Put sample to be test in the light path, press [START / STOP] to measure. The result will display in the

data sheet.

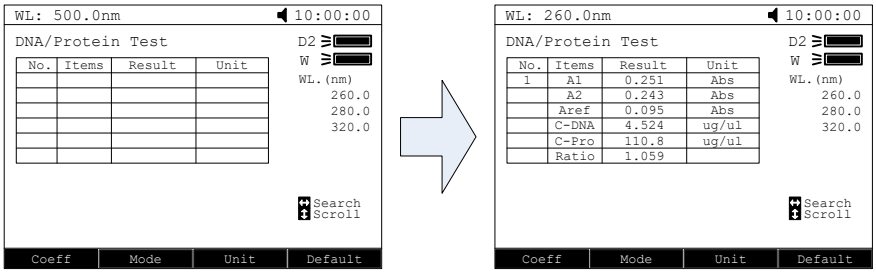

7. Print the Test Results

Press PRINT to print the Test Results.

DNA / Protein Test Report

File Name: Abc
Date & Time: mm-dd-yyyy, hh:mm:ss


Sample---1
260.0nm 0.251Abs
280.0nm 0.243Abs
320.0nm 0.095Abs
C-DNA 4.524ug/mL
C-Pro 110.8ug/mL
Ratio 1.059

Unit: ug/mL
Model: UV-3100PC
SN: UEFXXXXXXX
Version: A1.176
VWR International bvba.

Multi Wavelength Mode

1. Enter into Multi Wavelength Measurement

Main Menu, press [6 MNO] or ▲, ▼ to choose "Multi Wavelength Measurement" and then press [ENTER] to go into multi wavelength measurement interface.

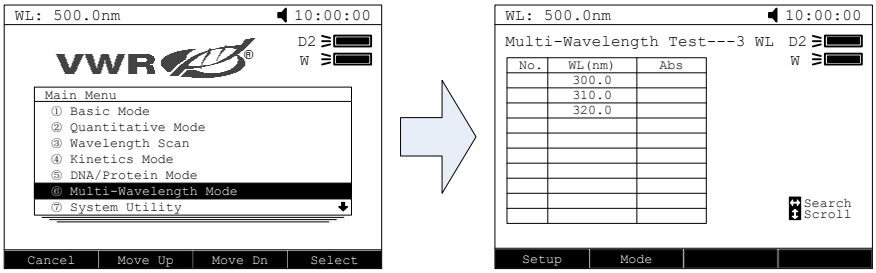

2. Setup Wavelength

Press [⌐] to go into wavelength setting interface，input all the wavelength value one by one by the numeric keypad. Press [ENTER] to confirm, press [ESC] to return.

3. Set Photometric Mode


Press [¬] to set the photometric mode， use ▲, ▼ to select "Abs. " or "T%" mode, press [ENTER] to confirm.

4. Blank

Place the Reference Sample in the light path， press [100%T oABS] to do blank.

5. Measure Samples

Place the sample to be test in the light path，press [START STOP] to measure, the test result will display in the data table.

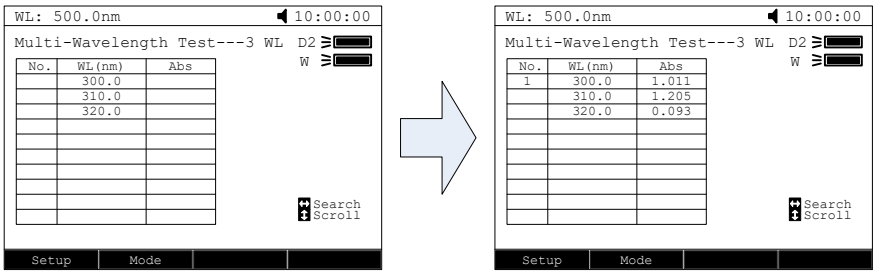

6. Print the Test Results

Press [PRINT] to print the Test Results.

Multi-Wavelength Test Report


File Name: Abc
Date & Time: mm-dd-yyyy, hh:mm:ss

Sample---1
500.0nm 0.251Abs
510.0nm 0.243Abs

Unit: Abs
Model: UV-3100PC
SN: UEFXXXXXXX
Version: A1.176
VWR International bvba.

System Utility

Main menu, press [7 PQRS] or use ▲, ▼ to select "System Utility" and press [ENTER] to confirm.

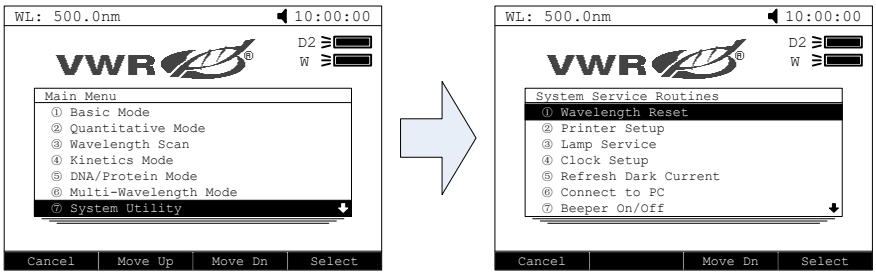


Wavelength Reset

Press [1] or use ▲, ▼ to choose "Wavelength Reset" then press [ENTER] to begin the calibration. During the course, opening the lid of the compartment is prohibited.

Printer Setup

Press [2 ABC] or ▲, ▼ to select "Printer Setup" then press [ENTER] to confirm.

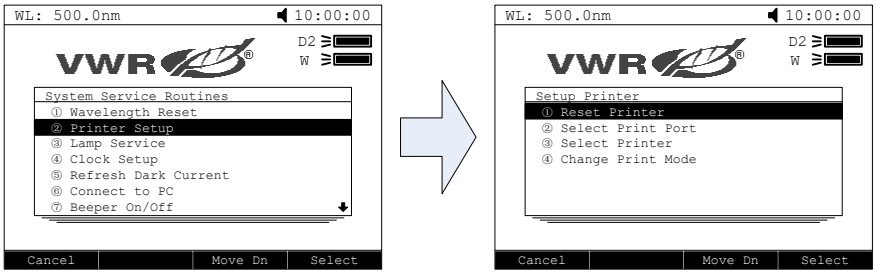

- **Reset Printer**

  Press [1] or use ▲, ▼ to choose "Reset Printer" and press [ENTER] to confirm. Then the printer will resume the initial condition.

- **Select print port**

  Press [2 ABC] or use ▲, ▼ to choose "Select print port" and press [ENTER] to confirm.

  Use ▲, ▼ to choose "LPT" or "Comm." and then press [ENTER] to confirm.

- **Select Printer**

  Press 3 DEF or use ▲, ▼ to choose "Select printer" and press ENTER to confirm.


  Use ▲, ▼ to select the printer's model and then press ENTER to confirm.

- **Change Print Mode**

  Press 4 GHI or use ▲, ▼ to choose "Change print mode" and press ENTER to confirm. Two modes are under your choice, one is "Print Data Sheet", the other is "print the display interface".

Lamp Service

Press 3 DEF or use ▲, ▼ to choose "Lamp Service", press ENTER to go into "lamps management" interface.

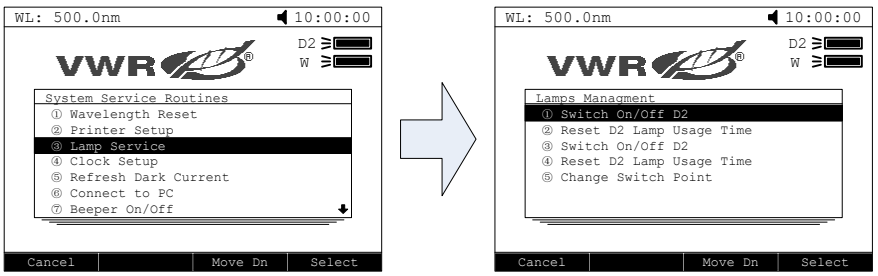

- **Switch On/Off D2 lamp**

  Press 1 or ▲, ▼ to choose "Switch On/Off D2 lamp" then press ENTER to switch on or switch off the D2 lamp.

- **Reset the usage time of D2 lamp**

  Press 2 ABC or use ▲, ▼ to choose "Reset D2 lamp usage time" and press ENTER.

Then users can find the D2 lamp's usage time and the system will ask if you are sure to reset the usage time. Press ▲, ▼ to choose "Yes" and then press [ENTER] to confirm. The system will record the usage time from zero.

- **Switch On/Off W lamp**

  Press [3 DEF] or use ▲, ▼ to choose "Reset W lamp usage time" and press [ENTER]. Then users can find the W lamp's usage time and the system will ask you if are sure to reset the usage time.

- **Reset the W lamp usage time**


  Press [4 GHI] or use ▲, ▼ to choose "Reset W lamp usage time" and press [ENTER]. Then users can find the W lamp's usage time and the system will ask if you are sure to reset the usage time. Press ▲, ▼ to choose "Yes" and then press [ENTER] to confirm. The system will record the usage time from zero.

- **Change Switch point**

  Press [5 JKL] or use ▲, ▼ to choose "Change Switch point" and press [ENTER] to confirm. Input the wavelength point value(325—375nm)and press [ENTER] to confirm.

Clock Setup

Press [4 GHI] or use ▲, ▼ to choose "Clock Setup", then press [ENTER] to go into "Clock Management" Interface.

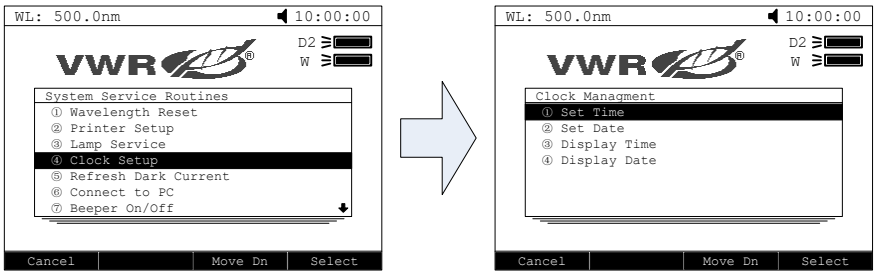

- **Set Time**

  Press [1] or ▲,▼ to choose "Set Time" and press [ENTER] to confirm. Input the time (Hour, Minute, Second) by numeric keypad, press [ENTER] to confirm and return automatically.

- **Set Date**

  Press [2 ABC] or use ▲,▼ to choose "Set Date" and press [ENTER] to confirm. Input the Date(Year, Month, Day) by numeric keypad, press [ENTER] to confirm and return automatically.

- **Display Time**

  Press [3 DEF] or use ▲,▼ to choose "Display Time" and press [ENTER] to confirm. Then the time will display on the Right Top corner.

- **Display Date**

  Press [4 GHI] or use ▲,▼ to choose "Display Date" and press [ENTER] to confirm. Then the Date will display on the Right Top corner.


Refresh Dark Current

Press 5 JKL or use ▲, ▼ to choose "Refresh Dark Current" and press ENTER to confirm. Then the system begins to refresh Dark Current.

Note: During the course, open the lid of the compartment is prohibited.

Connect to PC

Press 6 MNO or use ▲, ▼ to choose "Connect to PC", and press ENTER to go into waiting interface. When the instrument was connected to PC, it displays "Controlled by PC".

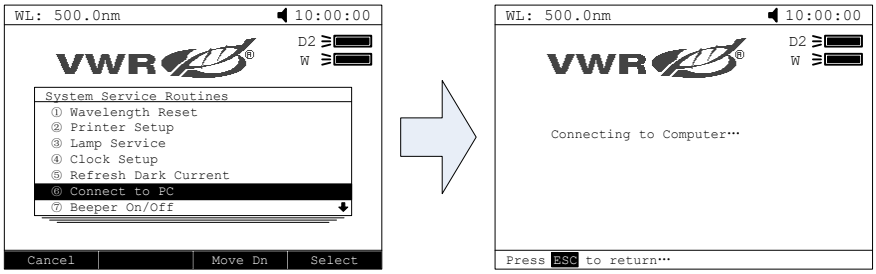

Beeper on/off

Press 7 PQRS or use ▲, ▼ to choose "Beeper on/off" and press ENTER to switch on or switch off the beeper.

Language Selection

Press 8 TUV or use ▲, ▼ to choose "Language Selection" and press ENTER to go into language selection interface，press ▲, ▼ to choose language (English, French, German, Spanish) and press ENTER to confirm. Then the interface will change into the language you just select.

Refresh System Baseline

Press [9 WXYZ] or use ▲, ▼ to choose "Refresh System Baseline" and press [ENTER] to confirm. Then the system will scan the baseline. During this course, opening the lid of the compartment is prohibited.

Deleted Entire Saved Files

Use ▲, ▼ to choose "Deleted Entire Saved Files" and press [ENTER] to confirm. Then the system ask you "Delete entire files, are you sure?", use ▲, ▼ to choose "Yes" , then all the files in the RAM will be deleted.

Restore Default Settings.

Use ▲, ▼ to choose "Restore Default Settings" and press [ENTER] to Confirm, then the system will restore the initialization.

Troubleshooting

Review the information in the table below to troubleshoot operating problems.

Problem	Cause	Solution
Power on, no response	Power cord connection is not reliable	Improve connectivity
	Fuse burning	Replace fuse
Measurement uncertainty	Warm up is not enough	Warm up more time
	Glass cuvettes used in UV Range	Use quartz cuvettes
	Sample is not Stable	Improve the sample
	The concentration of sample is too high	Diluted sample

	Power Supply Voltage Low or not Stable	Improve the Power Supply
	Lamp damage or lamp life maturity	Replace lamp
Dark Current Error when self-check	The lid of the compartment is open during self-check	Close the lid, restart
System Calibrate Failed	Something block the Light path	Remove it, calibrate again
Power on, back light is OK, but nothing display on the screen or display is not clear	Display Contrast problem	Adjust the contrast potentiometer
Measurements inaccurate	Cuvettes were contaminated	Clean cuvettes
	Samples were contaminated	Improve samples
	Worse matching of the cuvettes	Improve the matching of the cuvettes
	Dark current error	Resample dark current

Repair and Maintenance

Daily Maintain

Check the compartment

After measurement, the cuvettes with sample solutions should be taken out of the compartment in time. Or the volatilization of the solution would make the mirror go moldy. Users must pay more attention to the corrosive sample and liquid easy to volatilize. Any solution remains in the compartment should be wipe off immediately.

Surface Clean

The cover of the instrument is with paint. Please use wet towel to wipe off the drips on the surface immediately. Organic solution is forbidden to be used to clean the cover. Please wipe off the dirt on the cover timely.

Clean the cuvettes

After every test or after a solution change, the cuvettes should be cleaned carefully, or the remains on the surface would cause measuring error.

Check Lamp

In "Wavelength Scan" mode, set the test parameters as follow:

Scan From: 500nm

Scan To: 200nm

Scan Step: 1nm

Press Numeric Key "6", enter the amplifier of ADC as "0", and press key "Enter" to start scanning energy. After finished scan, press the function key "Search", then press key "<",">" to browse the energy value point to point. Check the energy of the wavelength 500nm and 200nm, check the energy, you need replace W lamp as following two cases:

Energy＜20 from 500~340nm W lamp is damaged

Energy＜5000 at 500nm W lamp energy is too low

And you need replace D2 lamp as following two cases:

Energy＜20 from 339~200nm D2 lamp is damaged

Energy＜1000 at 500nm D2 lamp energy is too low

Spare Parts Replacement

Replace the fuse

Danger! Be sure to switch off the power and unplug the socket before replacement!

1. **Tools preparation**


   Prepare a 3×75 Flat Blade screwdriver.

2. **Switch Off the power supply**

Switch off the power supply, and unplug the socket.


3. **Take out the Fuse Seat**

   Take out the fuse seat by the screwdriver.

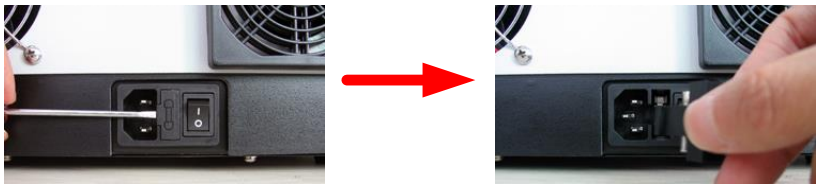

4. **Replace a new fuse**

   Pick out the spare fuse (3.15A/250V) and replace it to the working position.

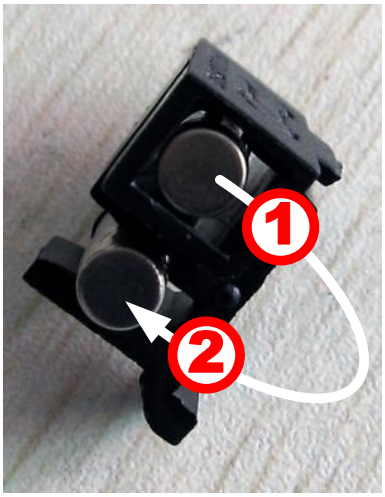

5. **Reset the fuse seat**

   Replace the fuse seat in the power socket.

6. **Switch on the power**

   Plug the socket and switch on the power.

Replace lamps

Hot! Wait 20 minutes before open the lamp chamber after power off to avoid scald!

1. **Tools preparation**


Prepare a 6×150mm Cross Blade screwdriver and a pair of glove.

2. **Power Off**

   Switch off the power supply and unplug the socket.


3. **Open the cover**

   Unscrew the 4 screws indicated(Each side with 2 screws)and remove the cover.

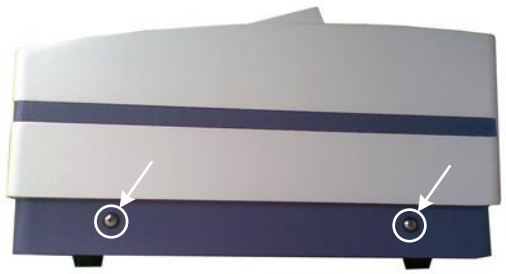


4. **Open the cover of the light chamber**

   Unscrew the 2 screws on the light chamber cover and remove it.

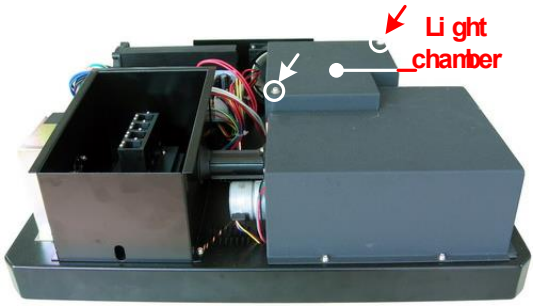

5. **Replace the D2 lamp**


   Unscrew the 2 screws on the D2 Flange (No.1), unplug the connector in the power board(No. 2)and remove the D2 lamp. Draw on the cotton glove and replace a new lamp. Fix the 2 screws and plug the connector again.

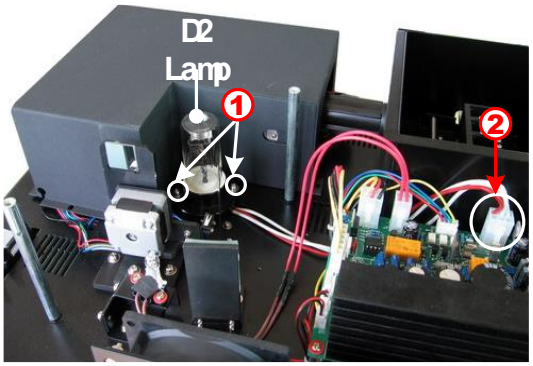

6. Replace W lamp


The Tungsten lamp is equipped with a blue-grey silicon coating by manufacturer. This coating is only a transport safety device. It can be removed with the first exchange of lamp.


Pull out the defected W lamp and draw on the cotton glove. Insert the new W lamp as deep as possible on the lamp seat. Be sure to keep the filament in the same direction as the old one face.

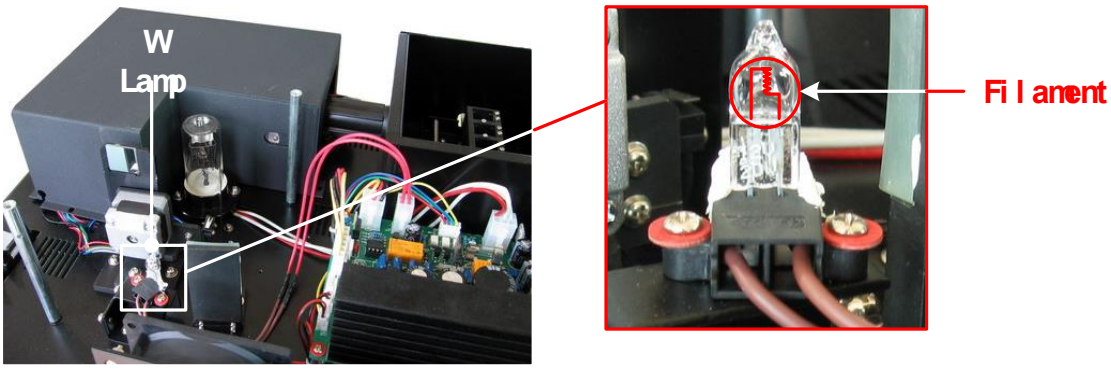

7. Adjust the position of the W lamp

Switch on the power(the Switch Mirror should be placed to the position as indicates). Observe the entrance facular, and it should in the center of the entrance hole. If the facular deviate to Left or Right, then loosen the No.1 screws in Fig. 5-8 and move the lamp seat to Left or Right until it focus on the center of the slot. Then fix the screws. If the facular deviate to Up and Down, then loosen the No.2 screws and move the lamp seat Up and Down until the facular focus on the center of the slot. Then fix the No. 2 screws again.

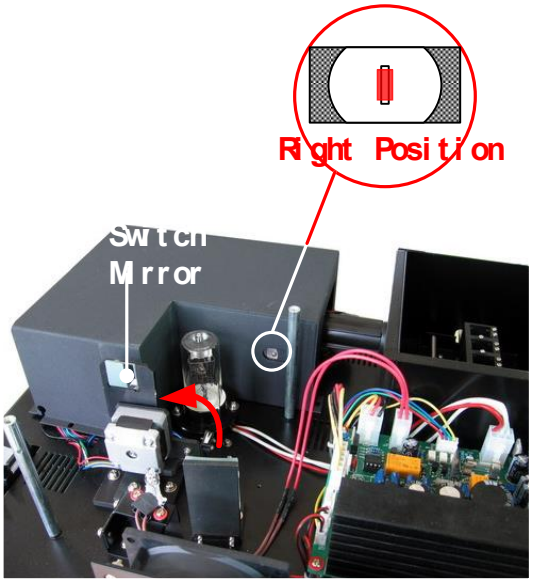

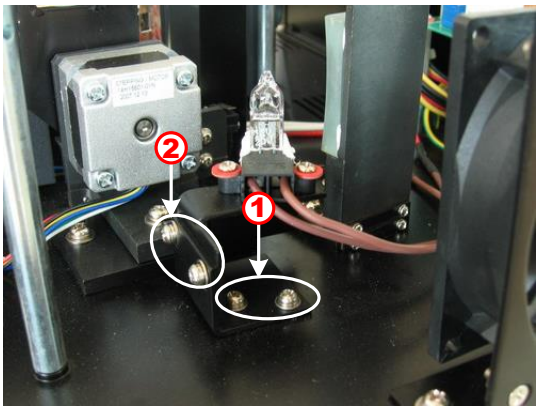

8. Finish

Reset the cover of the light chamber and fix the screws. Reset the cover of the instrument and fix the screws. Recover the Pole in the compartment, then the course finished.

Replace the Battery

Be sure to switch off the power supply and unplug the socket before open the Bottom Cover!

1. **Prepare the tools**

   Prepare a 6×150mm Cross Blade Screwdriver.

2. **Switch off the power supply**

   Switch off the power supply and unplug the socket


3. **Open the Bottom cover plate**

   Unscrew all the screws indicated then remove the bottom plate.



4. **Replace the Battery**

   Pick out the old battery and replace a new one.

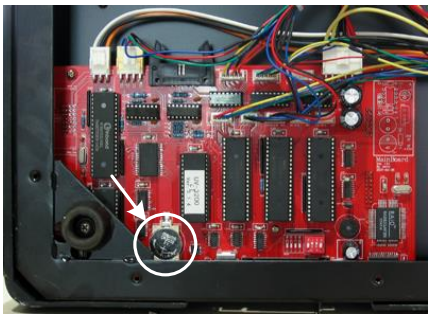

5. **Finish**

   Recover the bottom plate and fix all the screws, then the course finishes.

Accessories and Spare Parts

Description	Quantity	Cat. No.
CELL HOLDER, 4-CELL, 10MM	1PC	634-6003
CELL HOLDER, 4-CELL,10 TO 50MM	1PC	634-6004
CELL HOLDER, 4-CELL, 10 TO 100MM	1PC	634-6005
CELL HOLDER, FOR CYLINDRICAL CELL	1PC	634-6006
CELL HOLDER, WATER-JACKETED, 1 CELL, 10MM	1PC	634-6007
CELL HOLDER FOR MICRO CELLS,BEAM HEIGHT:15MM	1PC	634-0687
CELL HOLDER FOR TEST TUBES	1PC	634-6009
CELL HOLDER,8-POSITION AUTO CELL CHANGER	1PC	634-6010
CELL HOLDER, SOLID SAMPLE	1PC	634-6011
CELL HOLDER, WATER-JACKED, 4 CELL, 10MM	1PC	634-6012
CUVETTE, SQUARE.GLASS,10MM	4PCS	634-6013
CUVETTE, SQUARE.GLASS,20MM	2PCS	634-6014
CUVETTE, SQUARE.GLASS,30MM	2PCS	634-6015
CUVETTE, SQUARE.GLASS,50MM	2PCS	634-6016
CUVETTE, SQUARE.GLASS,100MM	1PC	634-6017
CUVETTE, SQUARE.QUARTZ,10MM	2PCS	634-6018
CUVETTE, SQUARE.QUARTZ,20MM	2PCS	634-6019
CUVETTE, SQUARE.QUARTZ,30MM	2PCS	634-6020
CUVETTE, SQUARE.QUARTZ,50MM	2PCS	634-6021
CUVETTE, SQUARE.QUARTZ,100MM	1PC	634-6022
CELL,QUARTZ, 100UL, 10MM, BEAM HEIGHT:15MM	1PC	634-0688
CELL,QUARTZ, 200UL, 10MM, BEAM HEIGHT:15MM	1PC	634-0689
CELL,QUARTZ, 500UL, 10MM	1PC	634-6025
FLOW CELL, 5 MM, GLASS, BEAM HEIGHT:15MM	1PC	634-0690
FLOW CELL, 10 MM, GLASS, BEAM HEIGHT:15MM	1PC	634-0691
FLOW CELL, 20 MM, GLASS, BEAM HEIGHT:15MM	1PC	634-0692
FLOW CELL, 30 MM, GLASS, BEAM HEIGHT:15MM	1PC	634-0693
FLOW CELL, 5 MM, QUARTZ, BEAM HEIGHT:15MM	1PC	634-0694

FLOW CELL, 10 MM, QUARTZ, BEAM HEIGHT:15MM	1PC	634-0695
FLOW CELL, 20 MM, QUARTZ, BEAM HEIGHT:15MM	1PC	634-0696
FLOW CELL, 30 MM, QUARTZ, BEAM HEIGHT:15MM	1PC	634-0697
PELTIER UNIT, BEAM HEIGHT:15MM	1PC	634-0698
SIPPER UNIT WITHOUT TEMP. CONTROL, BEAM HEIGHT:15MM	1PC	634-0699
SIPPER UNIT WITH PELTIER TEMP. CONTROL, BEAM HEIGHT:15MM	1PC	634-0763
LAMP, HALOGEN, 12V20W	1PC	634-6037
LAMP, DEUTERIUM	1PC	634-6038
PRINTER, THERMAL PRINTER	1PC	634-6039
THERMAL PAPER	1PC	634-6043
FUSE, 3.15A/250V	1PC	634-0651
FUSE, 500MA/250V	1PC	634-0652
BATTERY, 3V, CR2032	1PC	634-0653

Technical Service

Web Resources

Visit the VWR's website at **www.vwr.com** for:

• Complete technical service contact information

• Access to VWR's Online Catalogue, and information about accessories and related products

• Additional product information and special offers

Contact us For information or technical assistance contact your local VWR representative or visit. **www.vwr.com**.

Warranty

VWR International warrants that this product will be free from defects in material and

workmanship for a period of two (2) years from date of delivery except the lamps. Lamps have a warranty of 1000 h lamp usage time or 6 months max. If a defect is present, VWR will, at its option, repair, replace, or refund the purchase price of this product at no charge to you, provided it is returned during the warranty period. This warranty does not apply if the product has been damaged by accident, abuse, misuse, or misapplication, or from ordinary wear and tear.

For your protection, items being returned must be insured against possible damage or loss. This warranty shall be limited to the replacement of defective products. IT IS EXPRESSLY AGREED THAT THIS WARRANTY WILL BE IN LIEU OF ALL WARRANTIES OF FITNESS AND IN LIEU OF THE WARRANTY OF MERCHANTABILITY.

Equipment Disposal



This equipment is marked with the crossed out wheeled bin symbol to indicate that this equipment must not be disposed of with unsorted waste.

Instead it's your responsibility to correctly dispose of your equipment at lifecycle -end by handling it over to an authorized facility for separate collection and recycling. It's also your responsibility to decontaminate the equipment in case of biological, chemical and/or radiological contamination, so as to protect from health hazards the persons involved in the disposal and recycling of the equipment.

For more information about where you can drop off your waste of equipment, please contact your local dealer from whom you originally purchased this equipment.

By doing so, you will help to conserve natural and environmental resources and you will ensure

that your equipment is recycled in a manner that protects human health.

Thank you!

Part 2 : Software

Functions

This section introduces the functions of the UV-Vis Analyst.

Main Functions

Single wavelength photometric measurement

- Go to a desired wavelength quickly and conveniently.
- Photometric value display mode can be changed (%Transmittance or Absorbance).

Fixed Points Measurement

Multi-wavelength Photometric Measurement

- Up to 20 wavelength points can be set up.
- Results will be grouped into a table format automatically.

Concentration Measurement

- 2 methods to set up the regression curve.
  Up to 20 standards to set up the regression curve. The UV-Vis Analyst will calculate the working curve using a linear equation that fits the data. Enter factor values to generate regression curve.
- 3 methods for curve fit.
  Linear fit, Quadratic fit and Cubic fit.

Wavelength Scanning

- Allow user to set scan step (0.1, 0.2, 0.5, 1.0 and 5.0nm).
- Spectrum display mode can be changed (Wavelength-%Transmittance or Wavelength-Absorbance).
- Peaks and valleys will be automatically detected after scanning (User can define the peak threshold).
- Powerful spectrum processing functions are provided.

Time Scanning

- Allow user to set scan Interval (0.5, 1.0, 2.0, 5.0, 10, 30 and 60s).
- Spectrum display mode can be changed (Time-%Transmittance or Time-Absorbance).
- Peaks and valleys will be automatically detected after scanning (User can define the peak threshold).
- Powerful spectrum processing functions are provided.

DNA/Protein Measurement

- Wavelength points and ratios can be set up.
- Results will be grouped into a table format automatically.

Spectrum Processing Function

Trace a Spectrum
The cursor can be moved to a desired point in the spectrum displayed on the screen and the photometric data at this point is displayed.

Automatic Peak Detection
After a scanning is complete, peaks and valleys can be automatically detected and listed in a table format. They will also be labeled on the spectrum.

Scale Expansion
Simultaneous expansion of the X and Y axes are provided with the "Zoom" function. Display range can also be changed though the "Display Setup" function.

Differentiation
You can calculate and display the first through to the fourth derivative spectrum for a given spectrum. Derivative spectrum is useful for enhancing spectrum data that are not readily apparent in an absorbance spectrum.

Calculate Spectrum
You can calculate addition, subtraction, multiplication and division between two spectrum with the resulting data displayed on the screen.

System Check and Calibration Function

Instrument Validity Check
Up to 10 wavelength points can be set up in the instrument validity mode. Two methods can be selected (Photometric Validity measurement and Wavelength Validity measurement) and tolerance can be entered. Results will be grouped into a table format automatically.

Dark Current Check
You can resample the dark current of the instrument.

Spectrum Bandwidth Check
A special scan for checking spectrum Bandwidth and it will calculate the spectrum Bandwidth value automatically.

Energy of Light Sources Check
It allows scan the energy of light sources with a fixed gain (0-7).

Reset Wavelength
It affords to relocate the 656.1nm.

Installation


This section introduces how to setup the UV-Vis Analyst to PC.

PC System Requirements


- Pentium or above PC;
- CD-ROM;
- USB Ports.
- 512MB Memory(1GB or Above is strongly recommended);
- 50MB or above hard disc Space;
- Microsoft Windows XP/Vista/7/8.

Install UV-Vis Analyst

1. Put UV-Vis Analyst disc in the CD-ROM;
2. Double Click to open the CD-ROM, and then, double click **Setup.exe** which is under the root directory of CD to start installation , click **Next**;

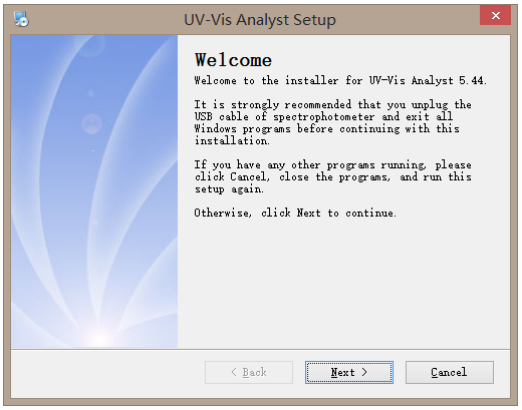


3. Choose install path, then click **Next** to copy files to PC;

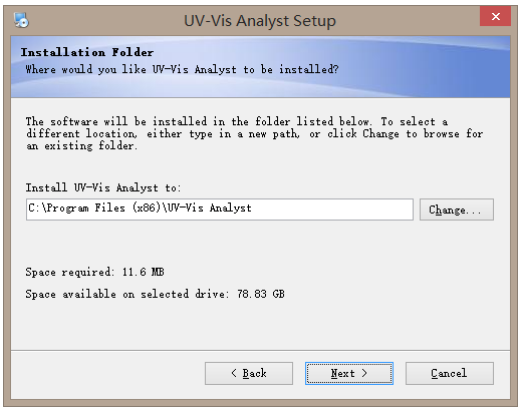

4. After finished copy all files of UV-Vis Analyst, it will start install the USB Drive. Click **Next**;




5. Select “I accept this agreement”. Click **Next** to copy files to PC;

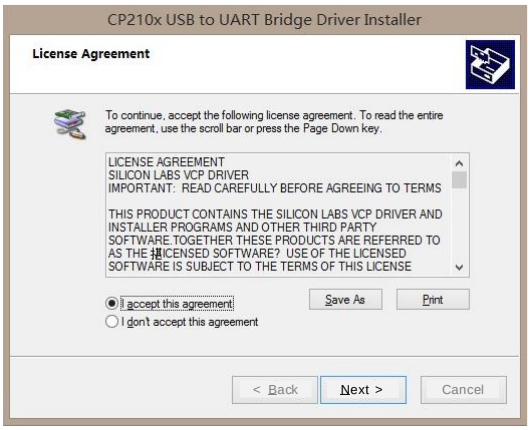

6. Click **Finish** to finish install USB Drive;



7. Click **Finish** to finish all the setup.

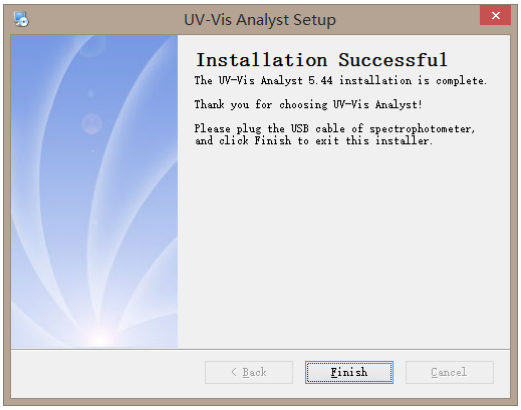


Uninstall UV-Vis Analyst

Start→All Programs →UV-Vis Analyst→Uninstall **UV-Vis Analyst** to remove.

Run UV-Vis Analyst

Before you run UV-Vis Analyst, please check these things as follow:

- *Computer has connected to spectrophotometer by using the USB cable;*
- *Spectrophotometer is on main interface.*

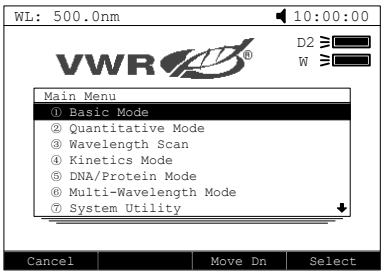


There are two ways to start the UV-Vis Analyst:

- Double-click shortcut icon on the desktop.
- Start→All Programs→UV-Vis Analyst→**UV-Vis Analyst**.

Set up Communication Port

Start the **UV-Vis Analyst**, on the **UV-Photometer** menu, click **Comm. Port Setup** appears the following box, select the Comm. Port and Baud Rate (38400), and click **OK**.

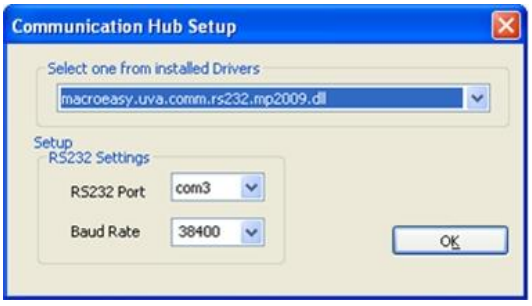


Introduction

We will introduce the UV-Vis Analyst in this chapter.

Main Interface

This is the main interface after start.

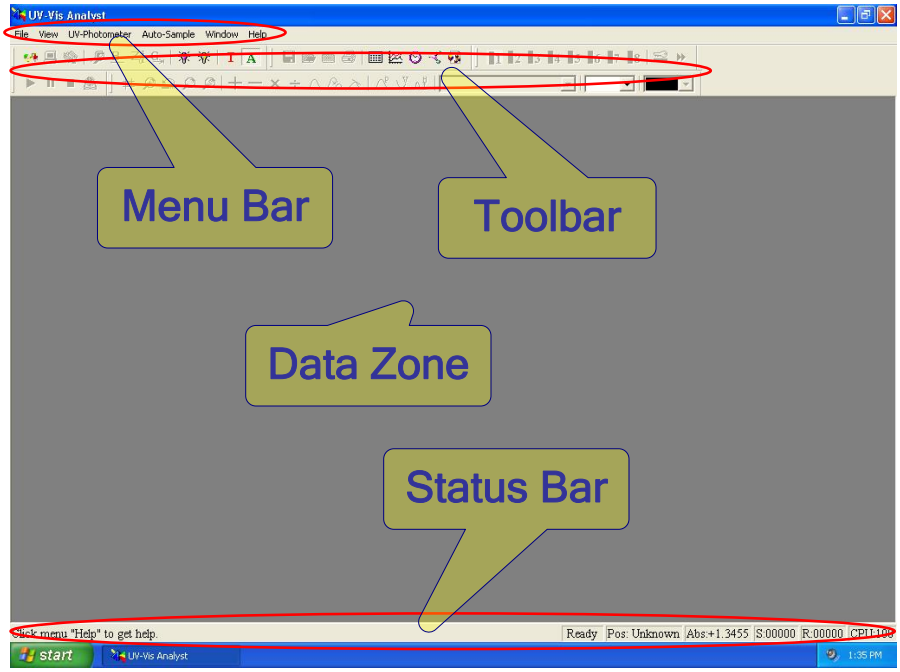

Menus Bar and Tools Bar

Menu bar and **Toolbar** are both provided in the software offering you two ways to select a desired function.

- On the menu bar, use your keypad or mouse to select the desired function.
- Almost all the functions listed in the menu bar can be reached by clicking a corresponding button in the toolbar.

Main Menu	Sub Menu	Tool	Function
File	New		New a Fixed Points Measurement
			New a Wavelength Scan
			New a Time Scan Measurement

			New a DNA/Protein Measurement
			New a Instrument Validity
	Open…		Open a spectrum/data file
	Close		Close current measurement
	Save		Save current measurement
	Save As…		Save current measurement as a
	Open file from		Open a file saved in instrument
	Export		Export data or method
	Print…		Print test report
	Print Setup…		Setup printer
	Exit		Exit UV-Vis analyst
View	Status Bar		Display/Hide status bar
	Status of		Display status of
	Status font		Setup font of status bar
	Customize		Define the information of display
	Peaks		Mark peak value
	Valleys		Mark valley value
	Magnify		Magnify the area selected

	Restore		Restore the default parameters
	Search		Search peak/valley one by one
UV-Photometer	Link		Connect to the Instrument
	Reset		Reset parameters of instrument
	Escape		Stop current measurement
	View dark Current		Retest the dark current
	Set Amplifier		Reset amplifier
	Locate 656.1nm		Relocate 656.1nm
	Calibrate System		Scan system baseline
	Automatic Blank		Do blank
	Slit Bandwidth *		Set slit bandwidth (0.5, 1, 2, 4, 5)
	Set Unit		Set unit
	Turn on/off W lamp		Turn on/off W lamp
	Turn on/off D2		Turn on/off D2 lamp
	D2/W Switch Point		Set switch point of D2/W
	Comm. Port Setup		Setup comm. port
	Change Password		Set/Change login password
Auto-sample	Locate Cell **		Locate cell (1-8) to light path

	Setup Multicell **		Setup Multicell
	Autorun **		Measure multi samples
Scan	Start		Start a measurement
	Stop		Stop a measurement
	Service		Measure spectrum and scan
Settings	Display Range		Setup scan display parameters
	Peak Height		Define peak/valley threshold
Compute	Add		Add two spectrum
	Sub		Subtract one spectrum from
	Multiply		Multiply two spectrum
	Divide		Divide one spectrum from another
	Moving Window		Smooth a spectrum with the
	Savitzky-Golay		Smooth a spectrum with the
	Derivate		Derivative of a spectrum
	Resample		Resample a spectrum
Window	New Window		New a measurement window as
	Cascade		Multi windows display in a
	Tile		Multi windows display in a tile

	Arrange Icons		Arrange all icons minimized
	Split		Split display area
Help	About UV-Vis		Display the information about the
			Setup measurement parameters
			Modify a measurement result
			Delete results selected
			Set and Goto one wavelength
			Display Instrument CPU
			Delete current Spectrum
			Display result as mode %T
			Display result as mode Abs
			Undo Scale

Note: * *Only for the model with Variable Slit*
** *Only for the model with 8-Cell Auto Charger*


Operation

This chapter introduces how to use UV-Vis Analyst.

Single Wavelength Photometric Measurement

The UV-Vis Analyst provides a convenient method to measure photometric value at a fixed wavelength.

1. Click  on the toolbar, appears **Goto specified wavelength**.

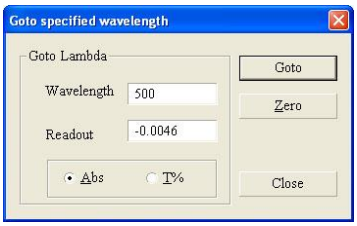


2. Key in the desired wavelength position, click **Goto**. The minimum wavelength step is 0.1nm in a range from 190-1100nm.
3. Place a reference in the sample compartment, click **Zero**.
4. Place a sample in the sample compartment. The wavelength position and photometric value will be displayed in the **Readout** box.

Fixed Point Measurement

This UV-Vis Analyst performs fixed wavelength measurement at 1-20 points and how to analyze unknown compounds against calibration standards.


Multi-wavelength Photometric Measurement

1. Click  on the toolbar, appears follow form.

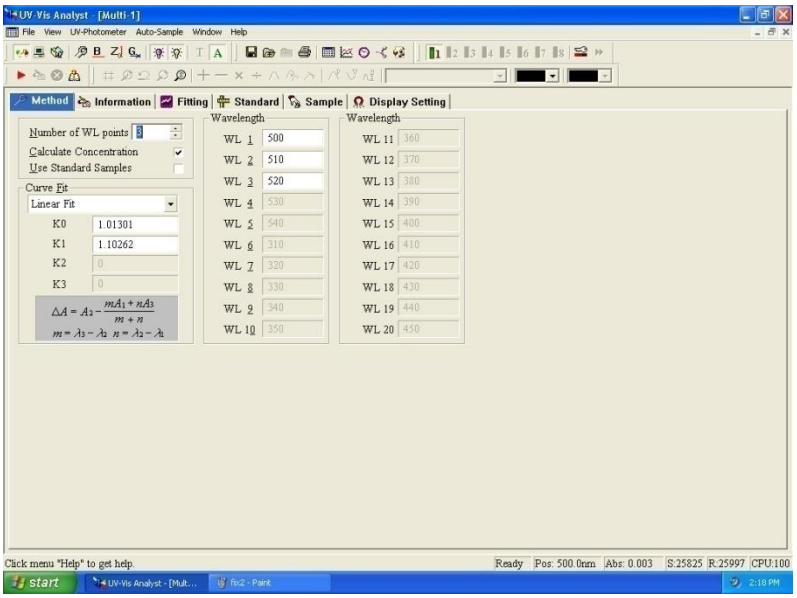


2. Click the **Method** tab.
3. Type the number of wavelength points in the **Number of WL Points** box, or click the **up/down** arrows next to the box set the wavelength points. Leave the two boxes **Calculate Concentration** and **Use Standard Samples**.
4. Key in the wavelength in the **Wavelength** box.


5. Place a reference in the sample compartment. Click  to do blank.
6. Click the **Sample** tab. It will display the following. The control menu contains six buttons: **Start**, **Delete**, **Modify**, **Recalculate**, **Data Font and Print**.

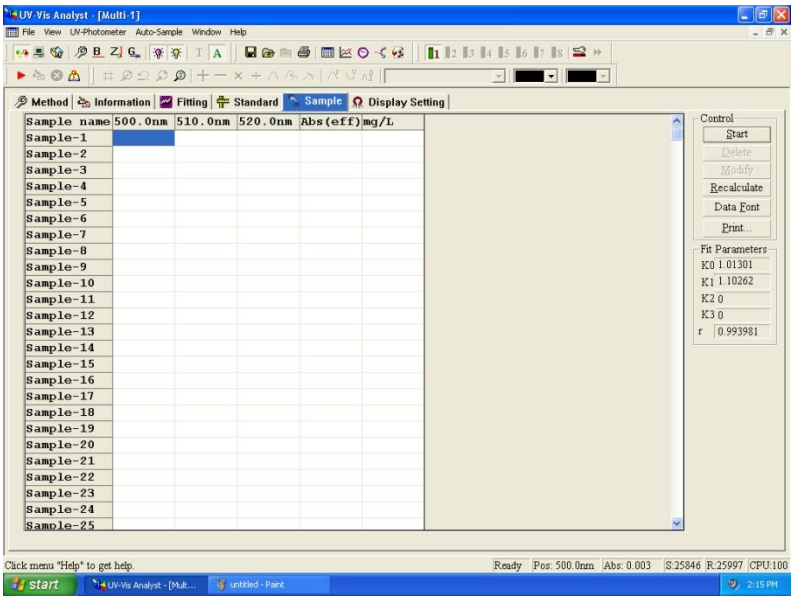

7. Place a sample in the sample compartment. Click **Start** or  to run a new measurement. The display will change to the following. Key in the sample name in the **Name** box.

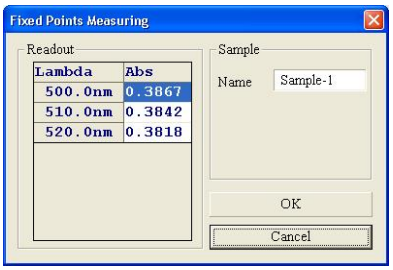

8. Click **OK**. The photometric data for sample will be listed in the **Sample** table.
9. Repeat steps 7-8 to measure all samples.

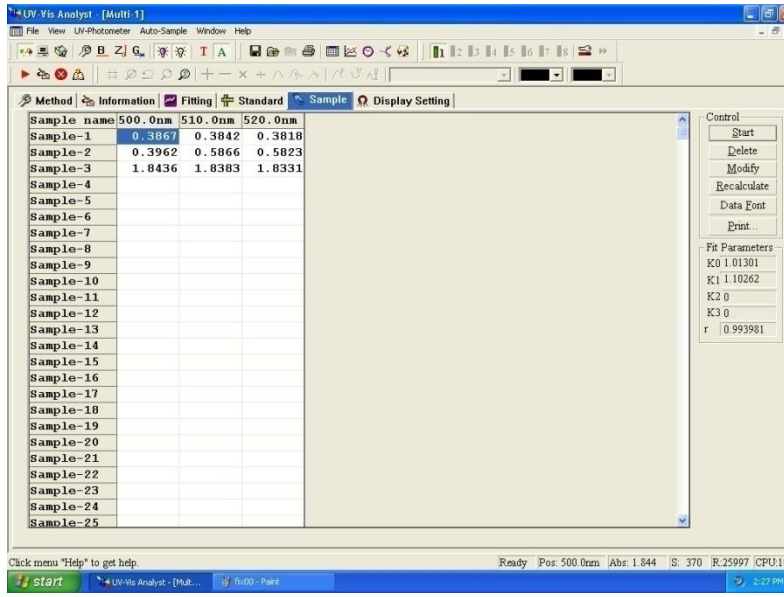

Concentration Measurement

Set Up Linear Regression Curve

There are two methods available to set up the linear regression curve. You can use standards to set up the regression curve or just key in the parameters manually. Use the following steps to select the method you wish to use.

1. Click  on the toolbar.
2. Click the **Method** tab.
3. Enter the number of wavelength points in the **Number of Points** box, or click the **up/down** arrow next to this box. With 2 wavelengths, the absorbance at the second reference wavelength is subtracted from the first to correct for background absorbance. With 3 wavelengths, the baseline between the first and third wavelengths is calculated and its value at the second wavelength is subtracted from the absorbance at the second wavelength to give the peak height.
4. Key in the wavelengths in the **Wavelength** boxes.
5. Tick the **Calculate Concentration** check box to activate concentration calculation.
6. Set up the linear regression curve.

   Method 1: Set up the linear regression curve with prepared standards.

   *(1) Tick the **Use Standard Samples** check box.*

   (2) Place the reference into the sample holder. Click  *to do blank.*

   *(3) Click the **Standard** tab.*


   *(4) Place Standard 1 in the sample compartment. Click **Start** to run a measurement.*

   *(5) Key in the concentration value of Standard 1 in the **Conc.** box.*

   *(6) Key in the sample name for the standard in the **Name** box.*

   *(7) Click **OK**. The photometric data, ΔA and concentration will be shown in the standard table.*

   (8) Repeat steps 4-7 to measure all the prepared standards.

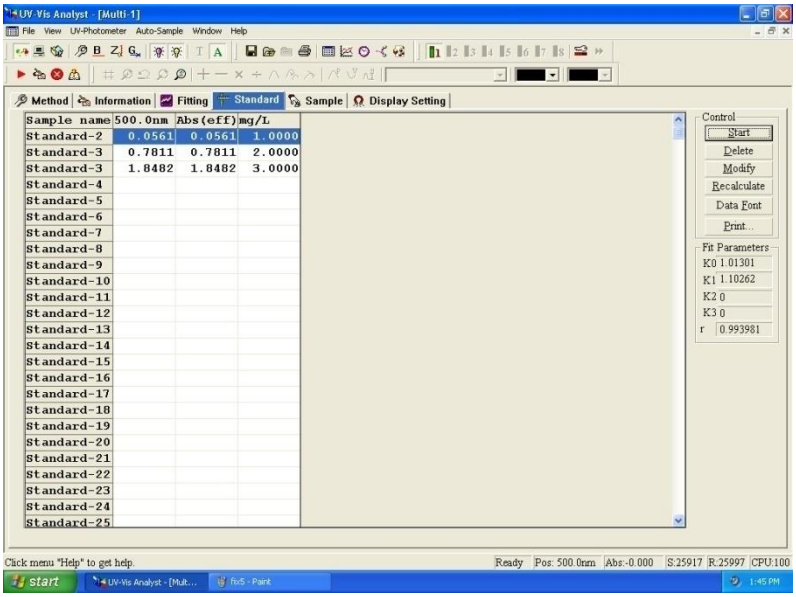

(9) *Click **down** arrow in **Curve Fit** box to select curve fit method.*


Method 2: Input the factor of the linear regression curve.


(1) *Leave the **Use Standard Samples** check box.*

(2) *Click **down** arrow in **Curve Fit** box to select curve fit method.*

(3) *Input the factor of the linear regression curve.*

7. Click **Fitting** tab to view the ***linear regression curve***. Click **Display Setting** tab to set the display parameters and unit of concentration.

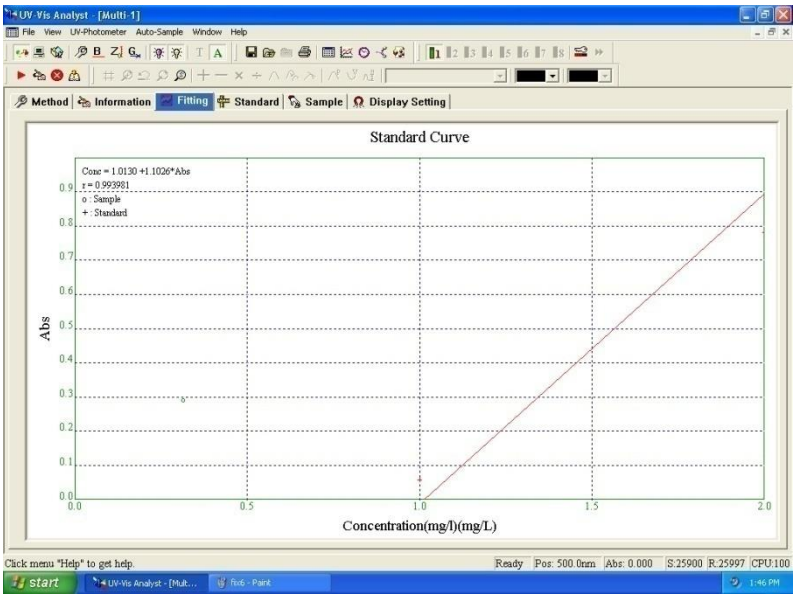


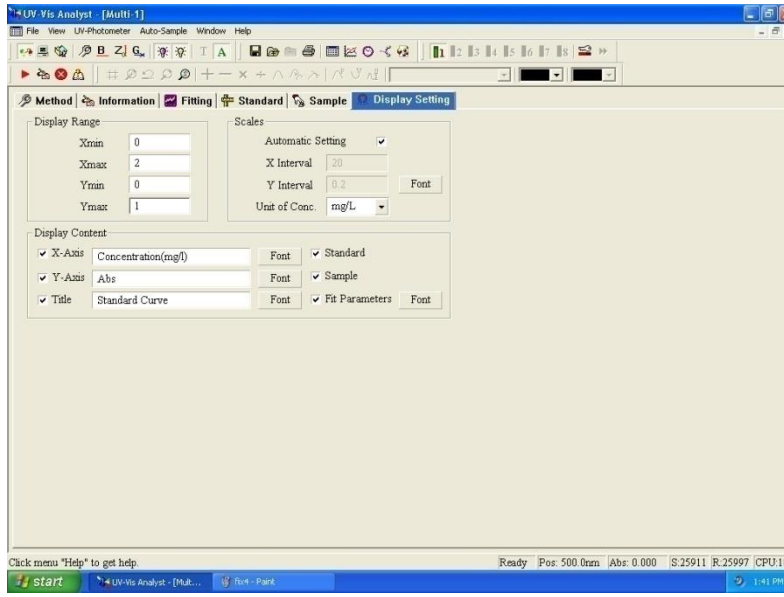

Measure Concentration by Using The Linear Regression Curve

The following procedure shows how to measure concentration of samples.

1. Set up linear regression curve or click  to open a file of linear regression curve (*.QUA).
2. Place reference into the sample holder. Click  to do blank.
3. Click the **Sample** tab.
4. Place Sample 1 into the sample holder.
5. Click **Start** to run a measurement.
6. UV-Vis Application Software will display the photometric value of Sample 1 at the fixed wavelength positions automatically. Type the sample name in the **Name** box. The default is **Sample-1**.
7. Click **OK**. The photometric result for Sample-1 will be listed in the sample data. Delta Abs. and concentration value of Sample-1 will also be displayed in columns 3 and 4.
8. Repeat steps 4-7 to measure remaining samples.

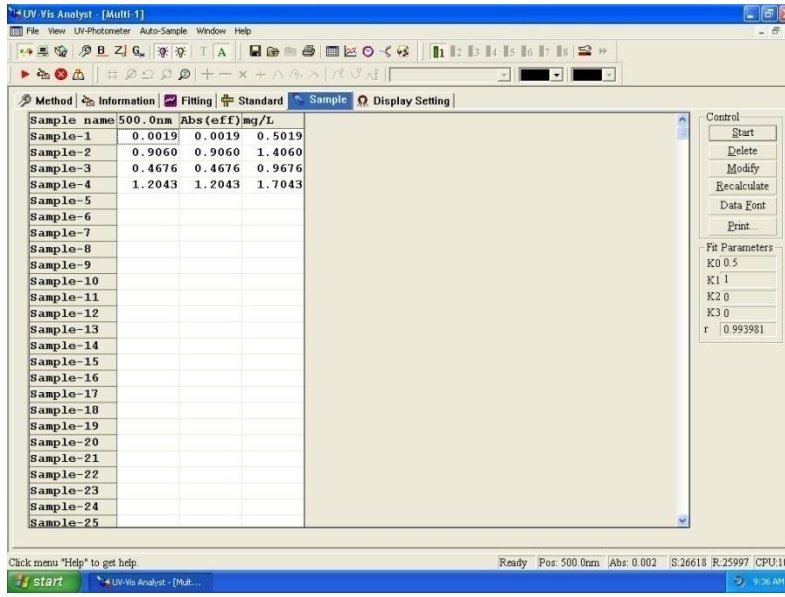


Wavelength Scanning


This chapter describes how to collect a spectrum while using Wavelength Scan function.

Scan Sample

1. Click  on the toolbar to new a sample scan measurement, appears the following form.

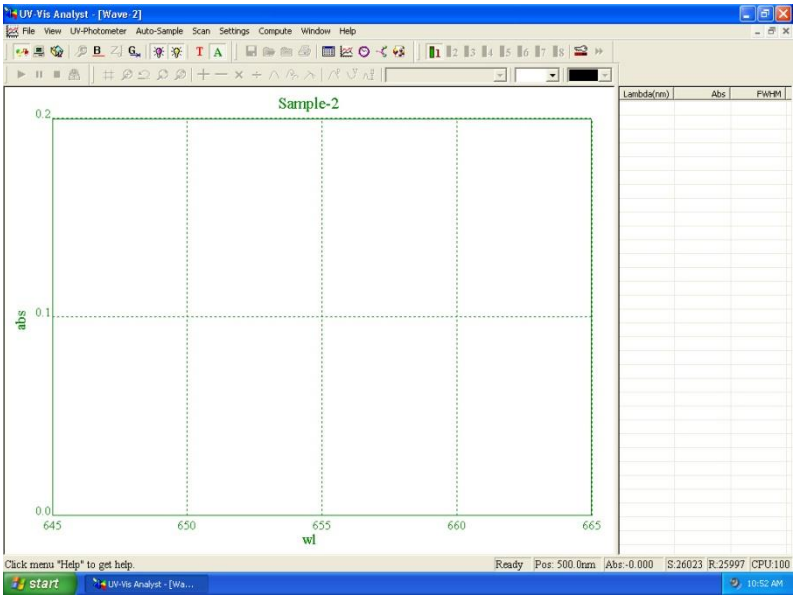


2. Click  on the toolbar, appears the following form. Input start wavelength in **From** box (range: 190-1100nm), end wavelength in **To** box (range: 190-1100nm), select scan interval (0.1, 0.2, 0.5, 1.0, 2.0 or 5.0nm) and Filter times (1, 3, 5, 10 or 30), click **OK**.

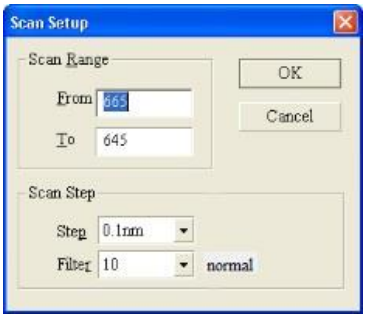


3. Click  on the toolbar to select %Transmittance mode or click  to select Absorbance mode.

4. Click  on the toolbar to set display parameters.

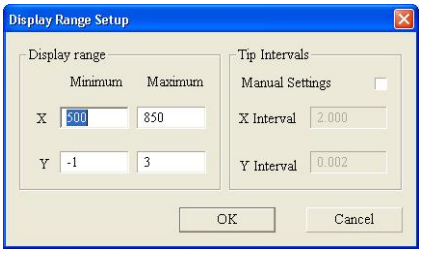


5. Place reference into the sample holder. Click  to scan baseline.


6. Place sample into the sample holder. Click  to scan sample, the real time spectrum will be displayed. Click  to cancel while scanning.

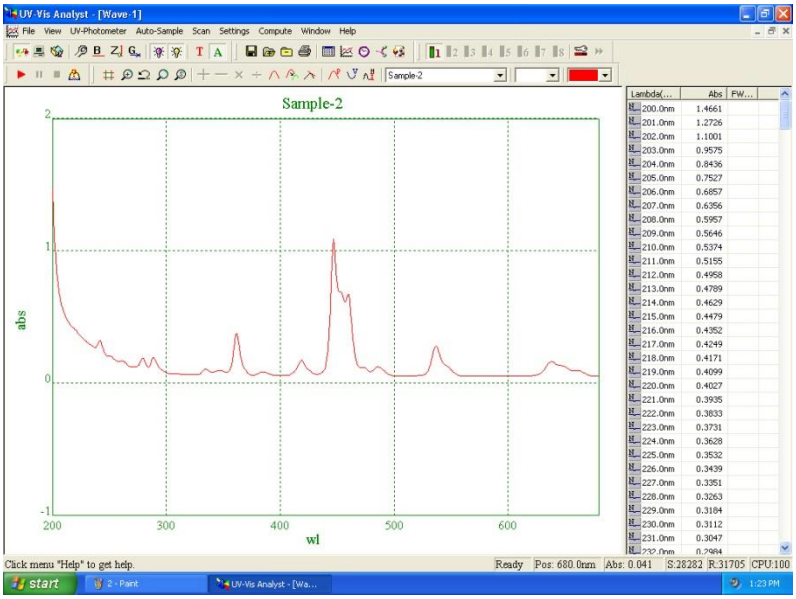

Auto List Peaks and Valleys

Click  on the toolbar to set the peak/Valley threshold (range: 0 to 1.000, step: 0.001), Input the threshold value, click **OK**. Click  to list peaks and click  to list valleys.

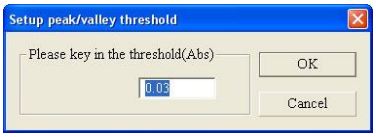

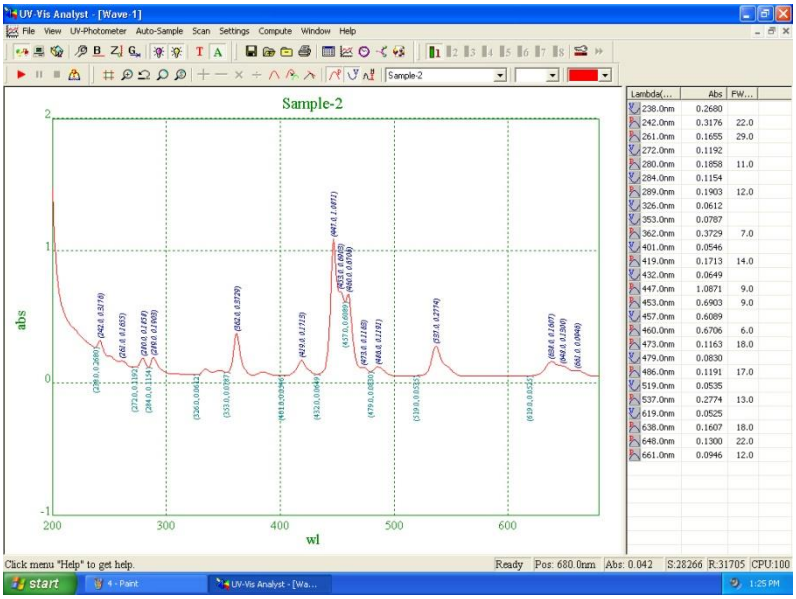

Rescale


Click  on the toolbar to set the new parameters for display.


Original Scales

Click  on the toolbar to restore the default display settings.

Zoom Selected Area


Click  on the toolbar to activate zoom function. Position the cursor in the upper-left corner of the area you want to select. Hold the left mouse button to drag the cursor to outline the spectrum area you want to enlarge. Release the mouse button. The part of the spectrum which is displayed within the outlined area will be enlarged. Click  to undo scale. To cancel zoom to click  again.

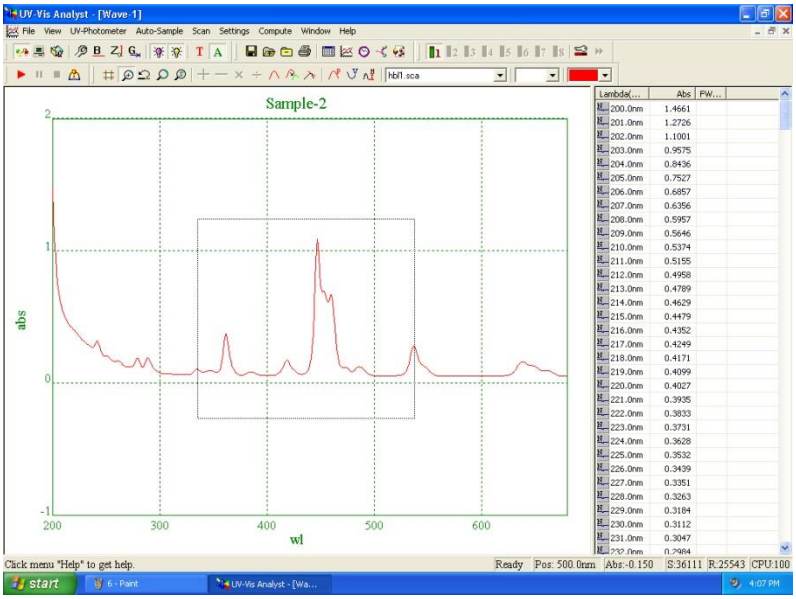

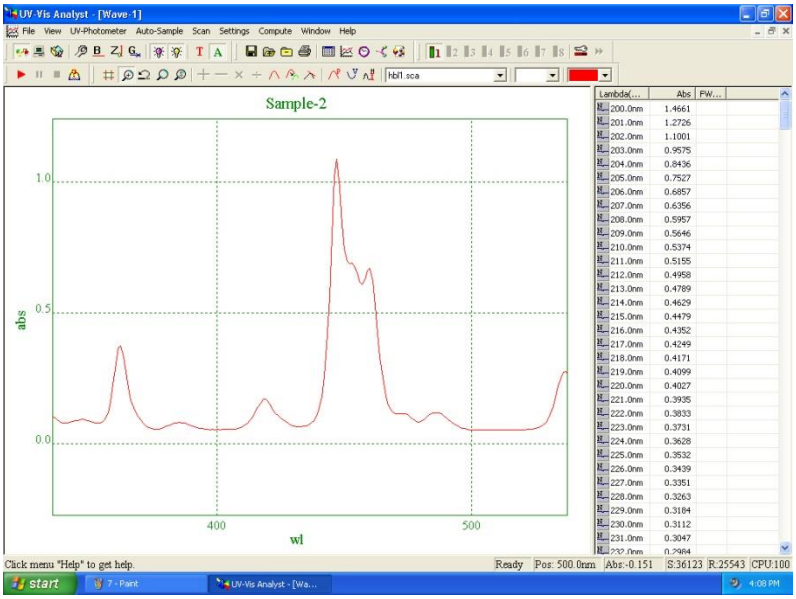

Trace a Spectrum

Click  on the toolbar, a crosshair cursor appears, move the cursor on the spectrum. Move the crosshair cursor left or right on the spectrum. The data in the cursor window indicate the X-axis and Y-axis values for the current cursor location. Click key "ESC" to release the crosshair cursor.

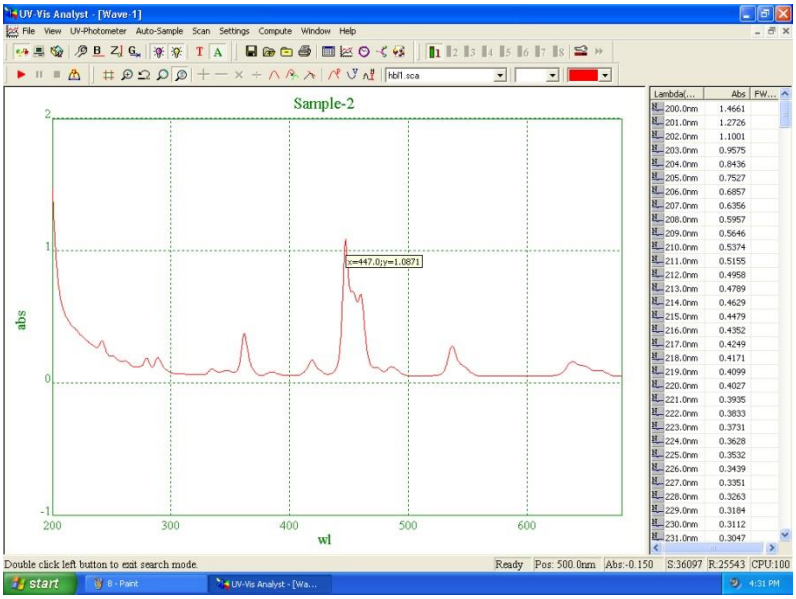


Select a Spectrum as Current


As UV-Vis Application Software can display several spectrum overlaid on the screen, you should specify the spectrum you wish to process. Click the **down** arrow on the toolbar. All spectrums will be listed in the pull-down menu. Click the spectrum you want to select. Its name will be listed in the Name Box and it will be referred to as **Current Spectrum**.



Derivative


Click  on the toolbar. The following dialogue box appears. Key in the class of derivative (1-10, depending on whether 1st, 2nd, … 10th derivative is required) and type a name for the result spectrum, then click **OK**. The result spectrum will be displayed overlaid with the original one.

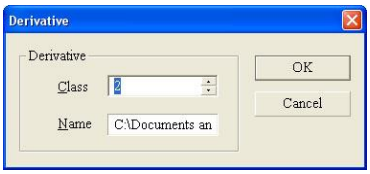


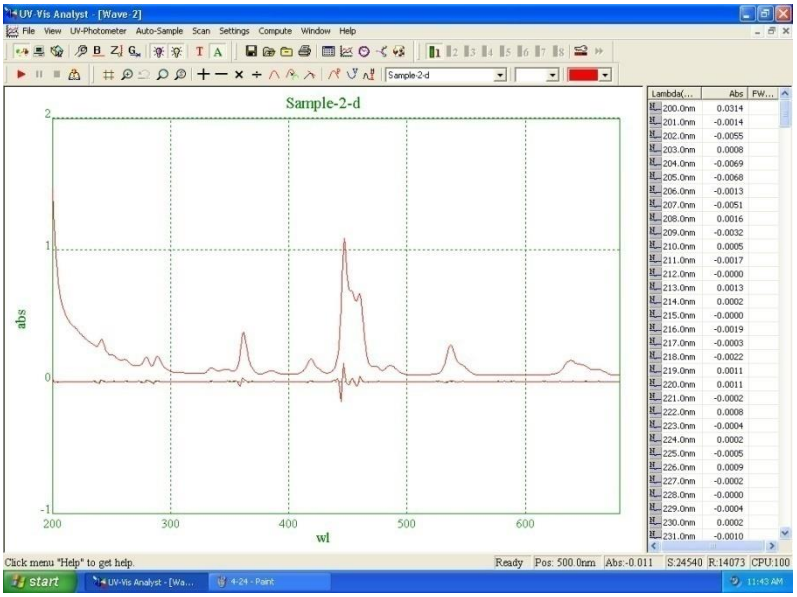

4.3.2.1 Moving Window Averaging


Click  on the toolbar. Appears following form. Click **up/down** arrow of the **Range** box to select range value, key a file name in the **Name** box, click **OK**. The result spectrum will be displayed overlaid with the original one.

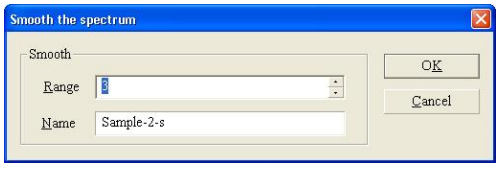


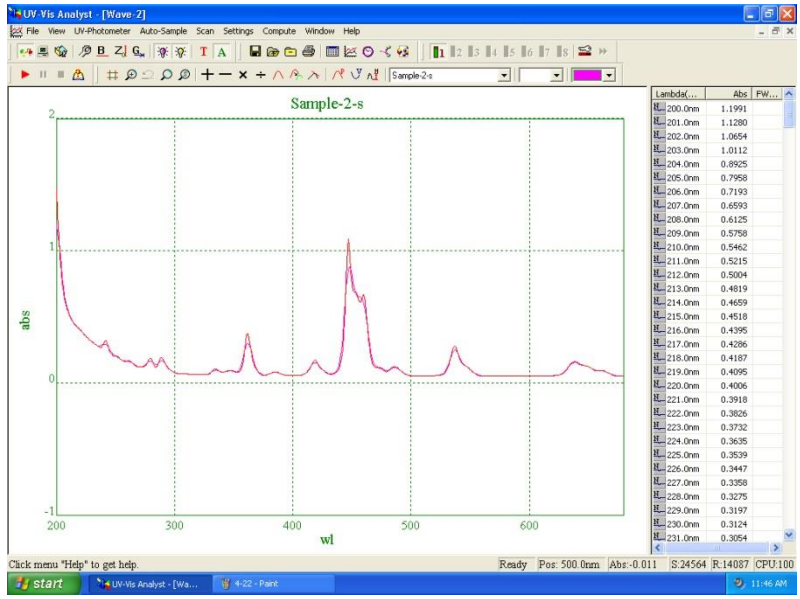

Savitzky-Golay Smoothing Filter


On the **Computer** menu, click **Savitzky-Golay Smoothing Filter**. Appears following form. Click **up/down** arrow to select the parameters, key a file name in the Name of Result box, click **OK**. The result spectrum will be displayed overlaid with the original one.

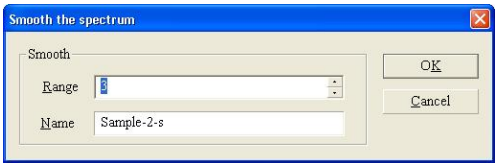


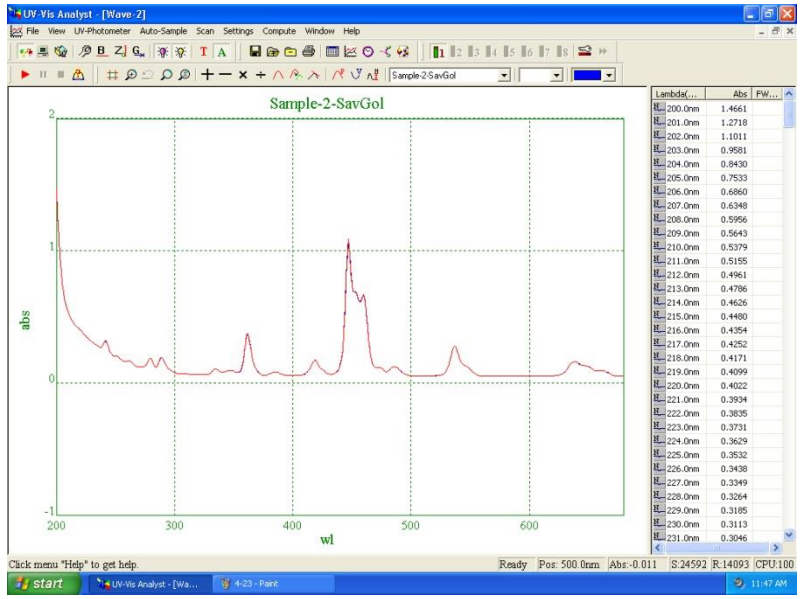

Resample

Click  on the toolbar. The following dialogue box will be displayed. Click **Up/Down** arrow to select Sample times. Click **OK**. The new spectrum displays.

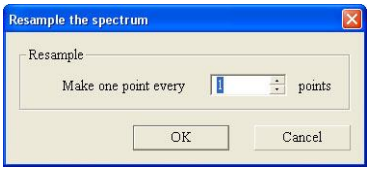

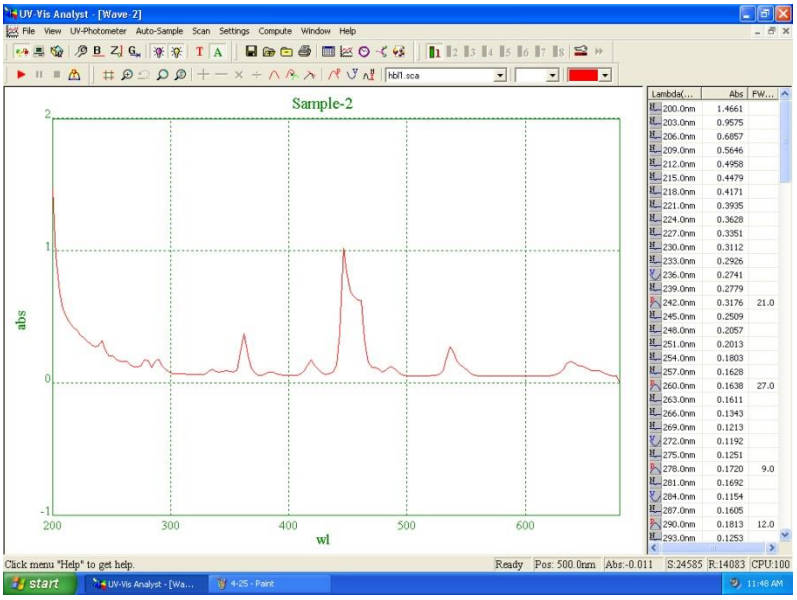



Spectrum Addition

Spectrum addition can assist in the development of artificial spectrum in multi-component mixtures.

Click  on the toolbar. The following dialogue box will be displayed. Click the **down** arrow next to **File 1** to select a spectrum and define it as source 1. Select a spectrum for **File 2** in the same way. It will not allow you to select the same spectrum twice. Key in a name for the **Result** spectrum and click **OK**. The result spectrum will be displayed on the screen.

UV-Vis Analyst will only add, subtract, multiply and divide two spectrums that are already displayed on the screen. Before arithmetic processing, load or collect two spectrums from memory.

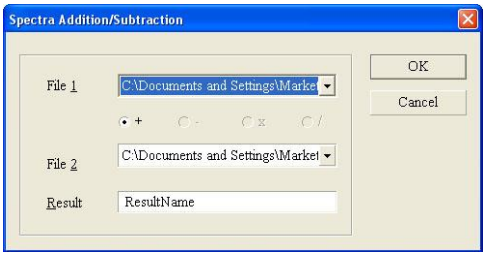


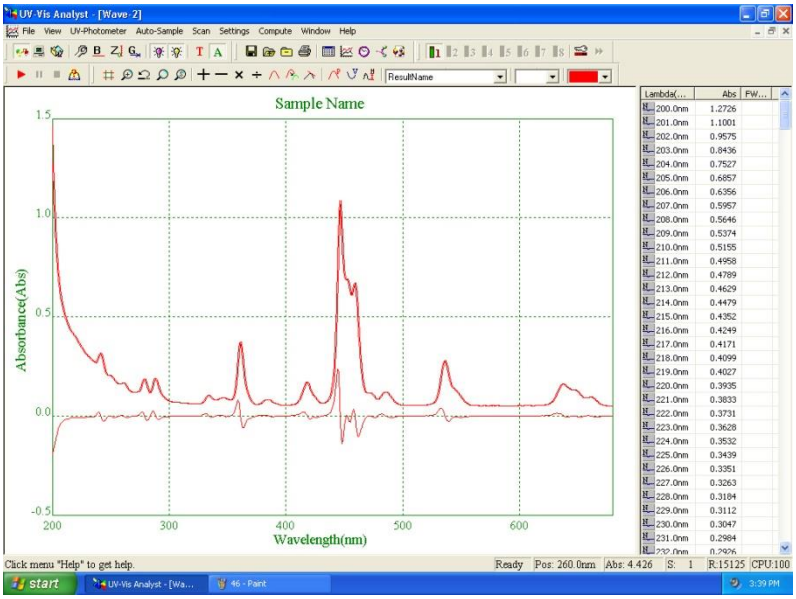


Spectrum Subtraction

Subtracting one spectrum from another has been a classical technique to offset spectrum interference from the spectrum of interest.

Click on the toolbar. The following dialogue box will be displayed. Click the **down** arrow next to **File 1** to select a spectrum and define it as source 1. Select a spectrum for **File 2** in the same way. It will not allow you to select the same spectrum twice. Key in a name for the **Result** spectrum and click **OK**. The result spectrum will be displayed on the screen.

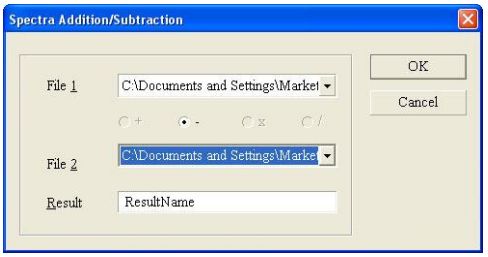


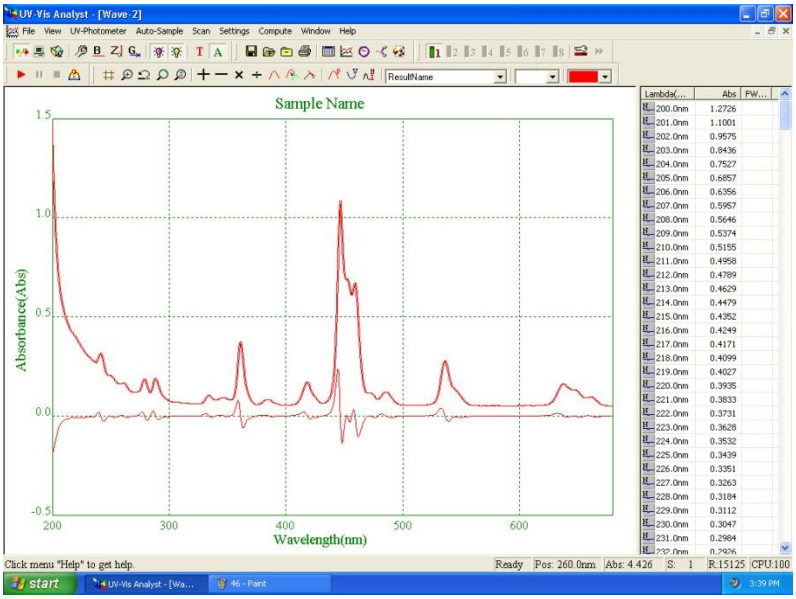


Spectrum Multiplication

Multiplying spectrum can assist in the development of artificial structure of spectrum in multi-component mixtures.

Click  on the toolbar. The following dialogue box will be displayed. Click the **down** arrow next to **File 1** to select a spectrum and define it as source 1. Select a spectrum for **File 2** in the same way. It will not allow you to select the same spectrum twice. Key in a name for the **Result** spectrum and click **OK**. The result spectrum will be displayed on the screen.

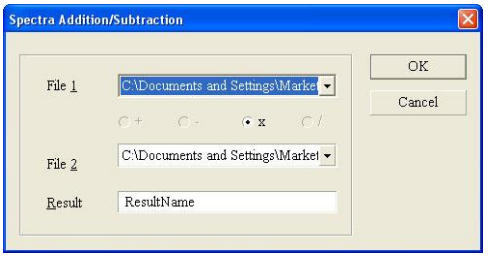


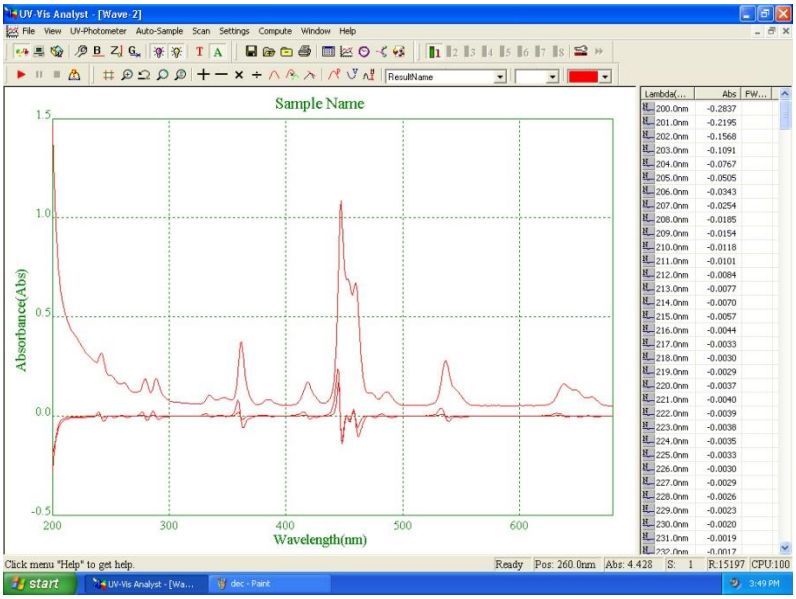


Spectrum Division

Dividing one spectrum from another has been a classical technique to offset spectrum interference from the spectrum of interest.

Click  on the toolbar. The following dialogue box will be displayed. Click the **down** arrow next to **File 1** to select a spectrum and define it as source 1. Select a spectrum for **File 2** in the same way. It will not allow you to select the same spectrum twice. Key in a name for the **Result** spectrum and click **OK**. The result spectrum will be displayed on the screen.

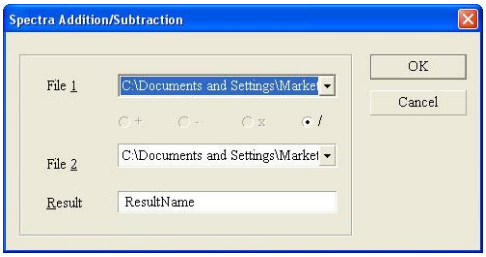

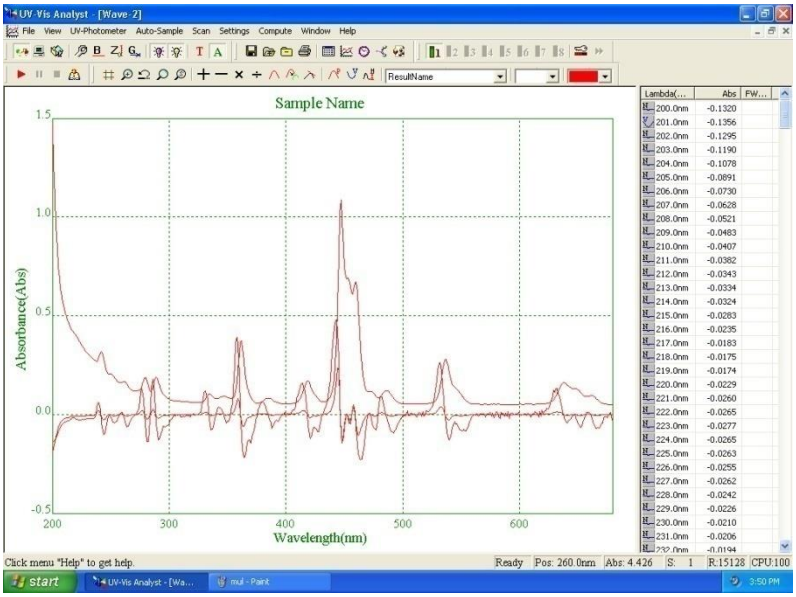


Unload a Spectrum

Select the spectrum you want to unload as the **Current Spectrum,** Click  on the toolbar to remove the spectrum from the display.


Define Display Information

Click  on the toolbar, appears the **Settings to display and print the spectra** form, click the **Legend** tab, type the information for display.

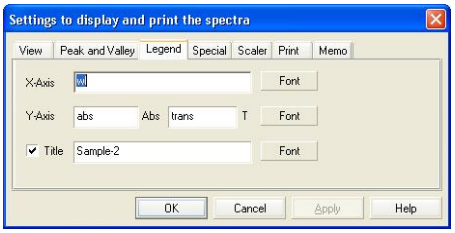

Edit Print Information

Click  on the toolbar, appears the **Settings to display and print the spectra** form, click the **Print** tab, type the information for print out.

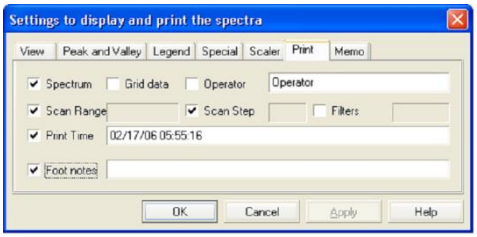


Time Scanning (Kinetic Analysis)

This chapter tells you how to obtain the absorbance or transmittance value for a sample as a function of time at a given wavelength.


Scan Sample

1. Click  on the toolbar, the following dialog box will appear.

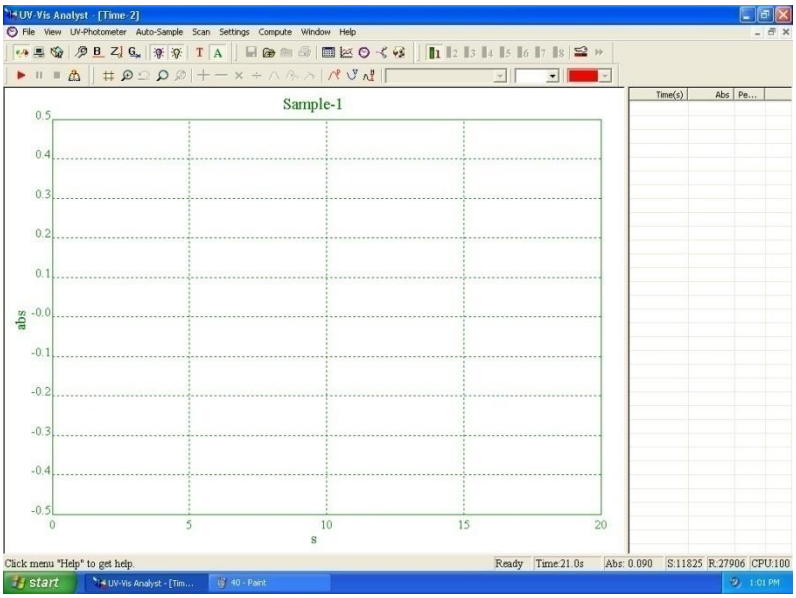


2. Click  on the toolbar to select the %transmittance mode or click  to select the absorbance mode.

3. Click  on the toolbar. A dialogue box will be displayed. Key in the wavelength, total time (in seconds) and scan step in the above dialog box. The wavelength range should be within 190 to 1100 nm. The upper limit for total time is 100000 seconds. Seven scan intervals can be selected from 0.5S, 1S, 2S, 5S, 10S, 30S and 60S. Click **OK**.

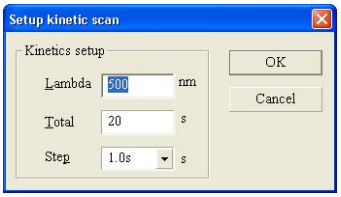


4. Place a reference in the sample holder. Click  on the toolbar.
5. Take out the blank in the sample holder, place a sample in it and close the cover.
6. Place a sample in the sample holder. Click  on the toolbar. The instrument will start scanning automatically. The graph will be displayed on the screen during time scanning. You can stop scanning by clicking  .

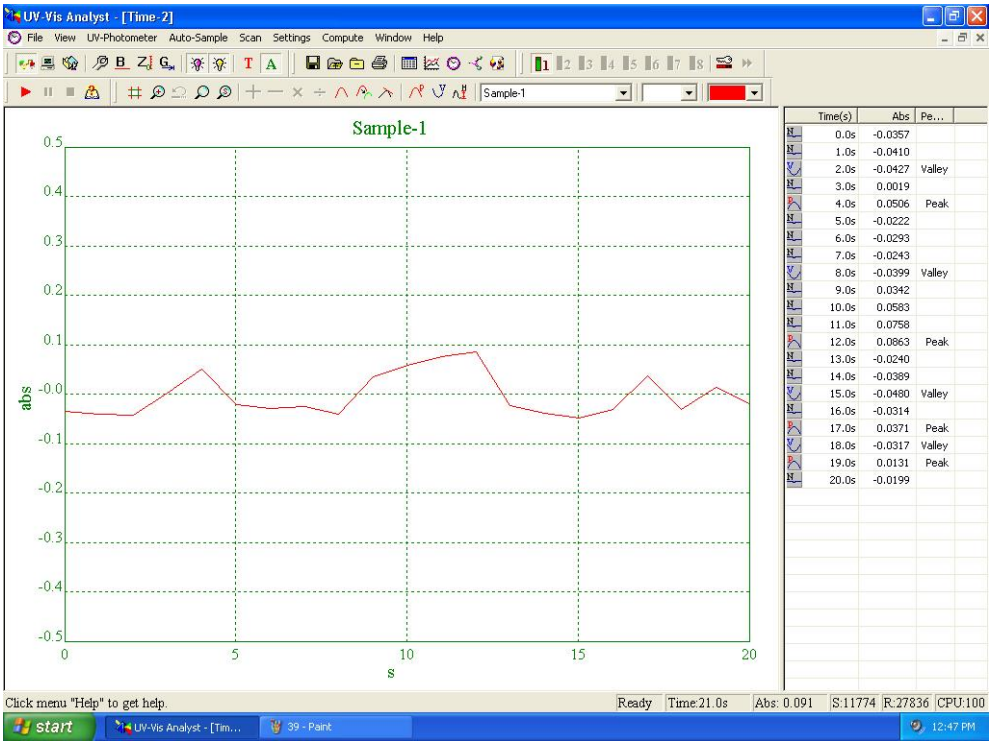


Calculate Rate

Click  on the toolbar, appears the **Settings to display and print the spectra** form, click the **Dynamic Analysis** tab, type the begin time in **Time Begin** box, type the end time in **Time End** box, and type the K factor in **K Factor** box, click **Calculate**, the result will be displayed.

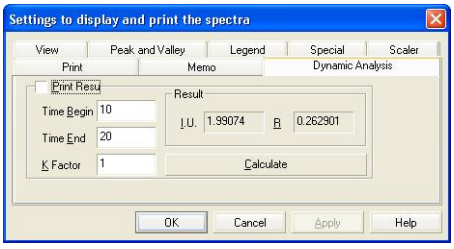


Define Display Information

Click  on the toolbar, appears the **Settings to display and print the spectra** form, click the **Legend** tab, type the information for display.

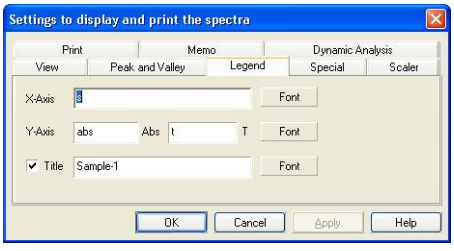

Edit Print Information

Click  on the toolbar, appears the **Settings to display and print the spectra** form, click the **Legend** tab, type the information for display.

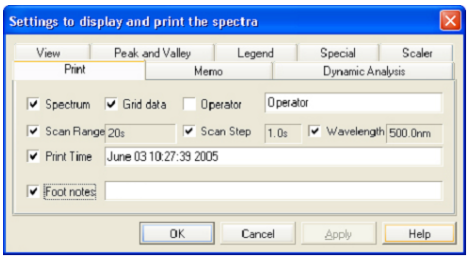


DNA/Protein Measurement


This chapter describes how to perform DNA/Protein measurement.

DNA/Protein Measurement


1. Click  on the toolbar, the following dialog box will appear.

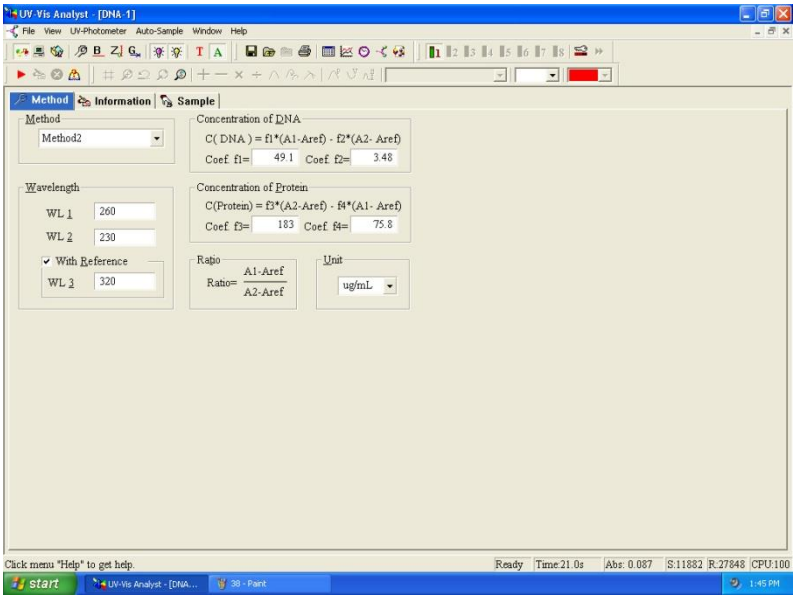


2. Click the **down** arrows of the **method** to select the **test method**. Key in the wavelength position in the **Wavelength** box. Key in the value of **DNA/Protein Conc.**
3. Place a reference in the sample holder. Click  on the toolbar to do blank.
4. Click the **Sample** tab. It will display the following. The control menu contains six buttons: **Start**, **Delete**, **Modify**, **Recalculate**, **Font** and **Print**.

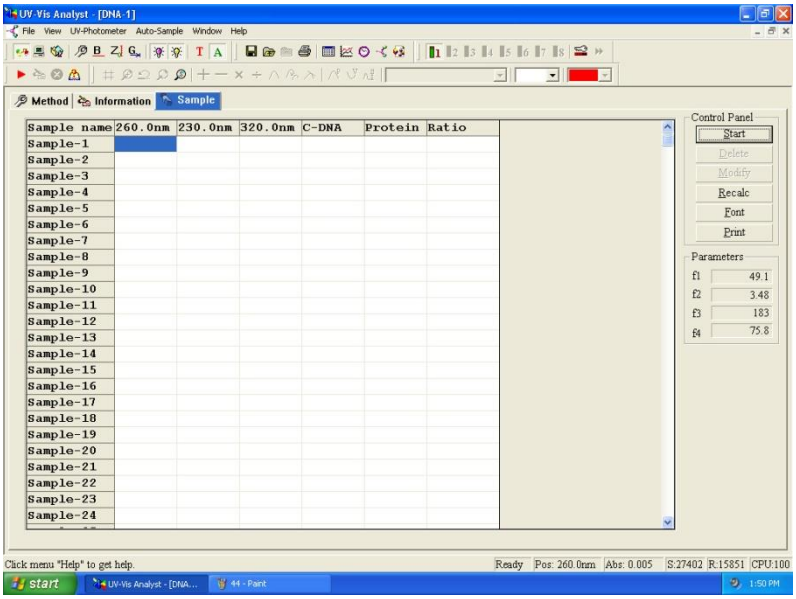

5. Place a sample in the sample holder. Click **Start** or  to run a new measurement. The display will change to the following.

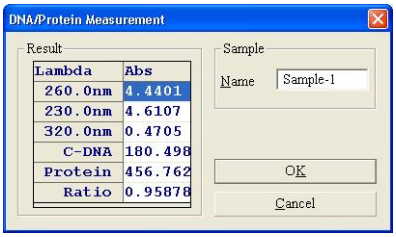

6. The UV-Vis Analyst will read the photometric value of **sample 1** at the fixed wavelength automatically. Key in the sample name in the **Name** box. Click **OK** after the measurement is complete. The photometric data for **sample 1** will be listed in the sample table.
7. Repeat steps 5-6 to test all samples.

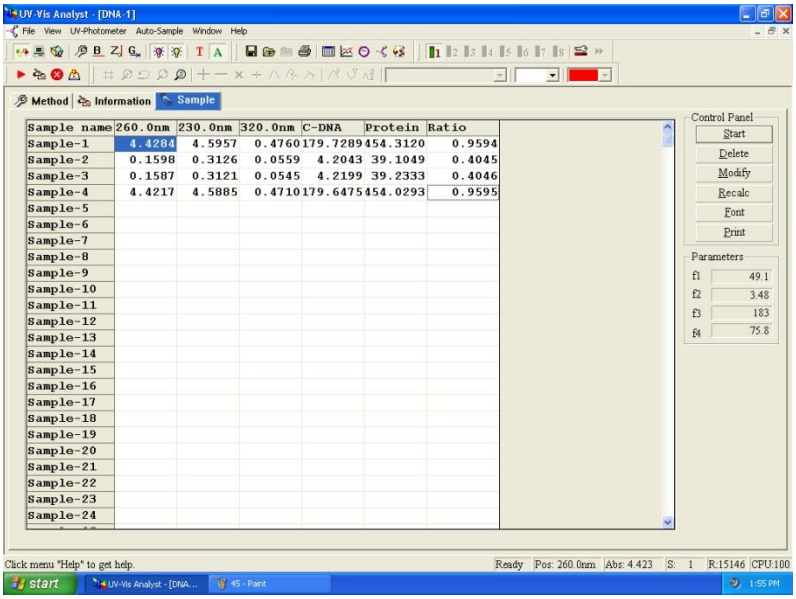
UV-Vis Analyst - [DNA-1]
File View UV-Photometer Auto-Sample Window Help
Method Information Sample
Sample name 260.0nm 230.0nm 320.0nm C-DNA Protein Ratio
Sample-1 4.4284 4.5957 0.4760 179.7289 454.3120 0.9594
Sample-2 0.1598 0.3126 0.0559 4.2043 39.1049 0.4045
Sample-3 0.1587 0.3121 0.0545 4.2199 39.2333 0.4046
Sample-4 4.4217 4.5885 0.4710 179.6475 454.0293 0.9595
Sample-5
Sample-6
Sample-7
Sample-8
Sample-9
Sample-10
Sample-11
Sample-12
Sample-13
Sample-14
Sample-15
Sample-16
Sample-17
Sample-18
Sample-19
Sample-20
Sample-21
Sample-22
Sample-23
Sample-24
Control Panel
Start
Delete
Modify
Recalc
Font
Print
Parameters
f1 49.1
f2 3.48
f3 183
f4 75.8
Click menu "Help" to get help.
Ready Pos: 260.0nm Abs: 4.423 S: 1 R:15146 CPU:100
start
UV-Vis Analyst - [DNA...
45 - Paint
1:55 PM

Appendix 1

Methods of Quantitative Analysis

Single Wavelength Method: Abs.=A1

Double Wavelengths Method: Abs.=m*A1-n*A2

Three Wavelengths Method: Abs.=A1-(W1-W2)*(A2-A3)/(W2-W3)-A3

Local VWR offices in Europe and Asia Pacific

Austria
VWR International GmbH
Graumanngasse 7
1150 Vienna
Tel.: +43 1 97 002 0
Email: info.at@vwr.com

Belgium
VWR International BV
Researchpark Haasrode 2020
Geldenaaksebaan 464
3001 Leuven
Tel.: +32 016 385 011
Email: vwr.be@vwr.com

China
VWR (Shanghai) Co., Ltd
Bld.No.1, No.3728 Jinke Rd,
Pudong New District
Shanghai, 201203- China
Tel.: +400 821 8006
Email: info_china@vwr.com

Czech Republic
VWR International s. r. o.
Veetee Business Park
Pražská 442
CZ - 281 67 Stříbrná Skalice
Tel.: +420 321 570 321
Email: info.cz@vwr.com

Denmark
VWR International A/S
Tobaksvejen 21
2860 Søborg
Tel.: +45 43 86 87 88
Email: info.dk@vwr.com

Finland
VWR International Oy
Valimotie 17-19
00380 Helsinki
Tel.: +358 9 80 45 51
Email: info.fi@vwr.com

France
VWR International S.A.S.
Immeuble Estréo
1-3 rue d'Aurion
93114 Rosny-sous-Bois cedex
Tel.: 0 825 02 30 30* (national)
Tel.: +33 (0) 1 45 14 85 00
(international)
Email: info.fr@vwr.com
* 0,18 € TTC/ min + prix appel

Germany
VWR International GmbH
Hilpertstrasse 20a
D - 64295 Darmstadt
Freecall: 0800 702 00 07
Tel.: +49 (0) 6151 3972 0 (international)
Email: info.de@vwr.com

Hungary
VWR International Kft.
Simon László u. 4.
4034 Debrecen
Tel.: +36 (52) 521-130
Email: info.hu@vwr.com

India
VWR Lab Products Private Limited
No.139. BDA Industrial Suburb,
6th Main, Tumkur Road,
Peenya Post,
Bangalore, India – 560058
Tel.: +91-80-28078400
Email: vwr_india@vwr.com

Ireland / Northern Ireland
VWR International Ltd / VWR
International (Northern Ireland) Ltd
Orion Business Campus
Northwest Business Park
Ballycoolin
Dublin 15
Tel.: +353 01 88 22 222
Email sales.ie@vwr.com

Italy
VWR International S.r.l.
Via San Giusto 85
20153 Milano (MI)
Tel.: +39 02-3320311/02-487791
Email: info.it@vwr.com

Korea
VWR International ~
17, Daehak 4-ro, Yeongtong-gu,
Suwon-si, Gyeonggi-do
Tel.: +82 31 645 7256
saleskorea@avantorsciences.com

Middle East & Africa
VWR International FZ-LLC
DSP Laboratory Complex
125, Floor 01
Dubai, United Arab Emirates
Tel.: +971 4 5573271
info.mea@vwr.com

The Netherlands
VWR International B.V.
Postbus 8198
1005 AD Amsterdam
Tel.: +31 020 4808 400
Email: info.nl@vwr.com

Norway
VWR International AS
Brynsalleen 4
0667 Oslo
Tel.: +47 22 90 00 00
Email: info.no@vwr.com

Poland
VWR International Sp. z o.o.
Limbowa 5
80-175 Gdansk
Tel.: +48 058 32 38 210
Email: info.pl@vwr.com

Portugal
VWR International - Material de
Laboratório, Lda
Centro Empresarial de Alfragide
Rua da Indústria, n° 6
2610-088 Amadora
Tel.: +351 21 3600 770
Email: info.pt@vwr.com

Singapore
VWR Singapore Pte Ltd
18 Gul Drive
Singapore 629468
Tel: +65 6505 0760
Email: sales.sg@vwr.com

Spain
VWR International Eurolab S.L.
C/ Tecnología 5-17
A-7 Llinars Park
08450 - Llinars del Vallès
Barcelona
Tel.: +34 902 222 897
Email: info.es@vwr.com

Sweden
VWR International AB
Fagerstagatan 18a
163 94 Stockholm
Tel.: +46 08 621 34 00
Email: info.se@vwr.com

Switzerland
VWR International GmbH
Lerzenstrasse 16/18
8953 Dietikon
Tel.: +41 044 745 13 13
Email: info.ch@vwr.com

UK
VWR International Ltd
Customer Service Centre
Hunter Boulevard
Magna Park
Lutterworth
Leicestershire
LE17 4XN
Tel.: +44 (0) 800 22 33 44
Email: uksales@vwr.com