



SpectraMax Mini

Multi-Mode Microplate Reader

Make your first move with confidence



SpectraMax Mini Multi-Mode Microplate Reader

Affordable multi-mode microplate reader with best-in-class data analysis software

The SpectraMax® Mini Multi-Mode Microplate Reader offers a budget-friendly solution that gives you exactly what you need when you need it. With three modes of detection for full spectrum absorbance, fluorescence and luminescence, the SpectraMax Mini reader is user upgradeable—from 2 to 3 modes, providing researchers with a versatile instrument that does not compromise on sensitivity or performance. Included SoftMax® Pro Software enables rapid data generation and analysis, helping you achieve more in less time.

Key benefits

- Ease of use to save time on training, data generation, and data analysis, for faster results
- Superior data analysis with best-in-class SoftMax Pro Software, offering pre-configured protocols for the most commonly used applications
- Flexibility with plate types from 6- to 384-well formats; up to three modes of detection measuring absorbance with fluorescence, luminescence, or both

Everything you need to run your most crucial applications quickly and easily

Access your favorite applications, including ELISA, DNA and protein quantitation, cell viability, gene expression, and more with a compact, upgradeable microplate reader that does everything you need and lets you add detection modes in the future.

The SpectraMax Mini reader provides tunable absorbance, filter-based fluorescence and luminescence, temperature control, and easy analysis and data export with SoftMax Pro Software.




Key features



**Best-in-class data analysis
with SoftMax Pro**



**Versatile multi-mode
detection**

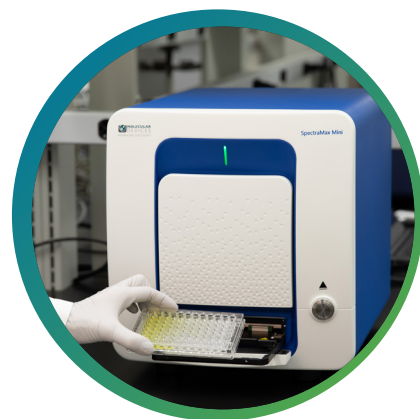

**Uncompromised
performance**

It's as easy as 1-2-3



1

Set up



2

Analyze data




3

Get results


Fluorescent DNA quantitation

Determine DNA concentrations precisely and accurately over a wide dynamic range

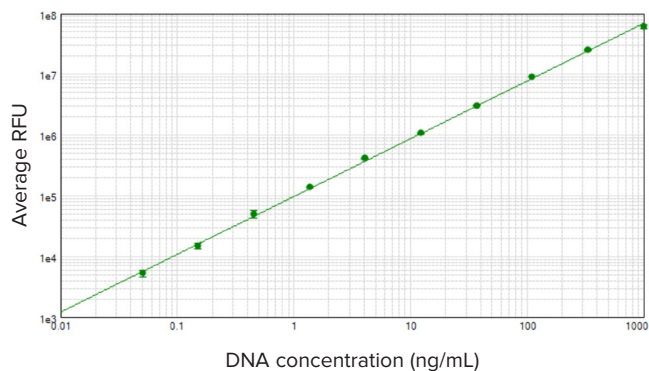




Cell viability

Easily measure cell viability with a luminescence readout

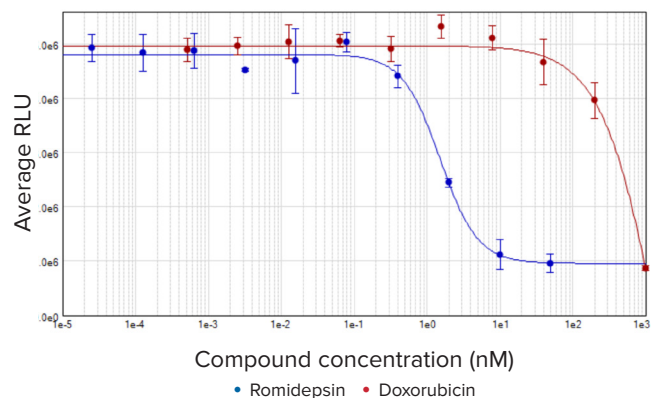

Protein quantitation


Determine protein concentrations in cell lysates and more

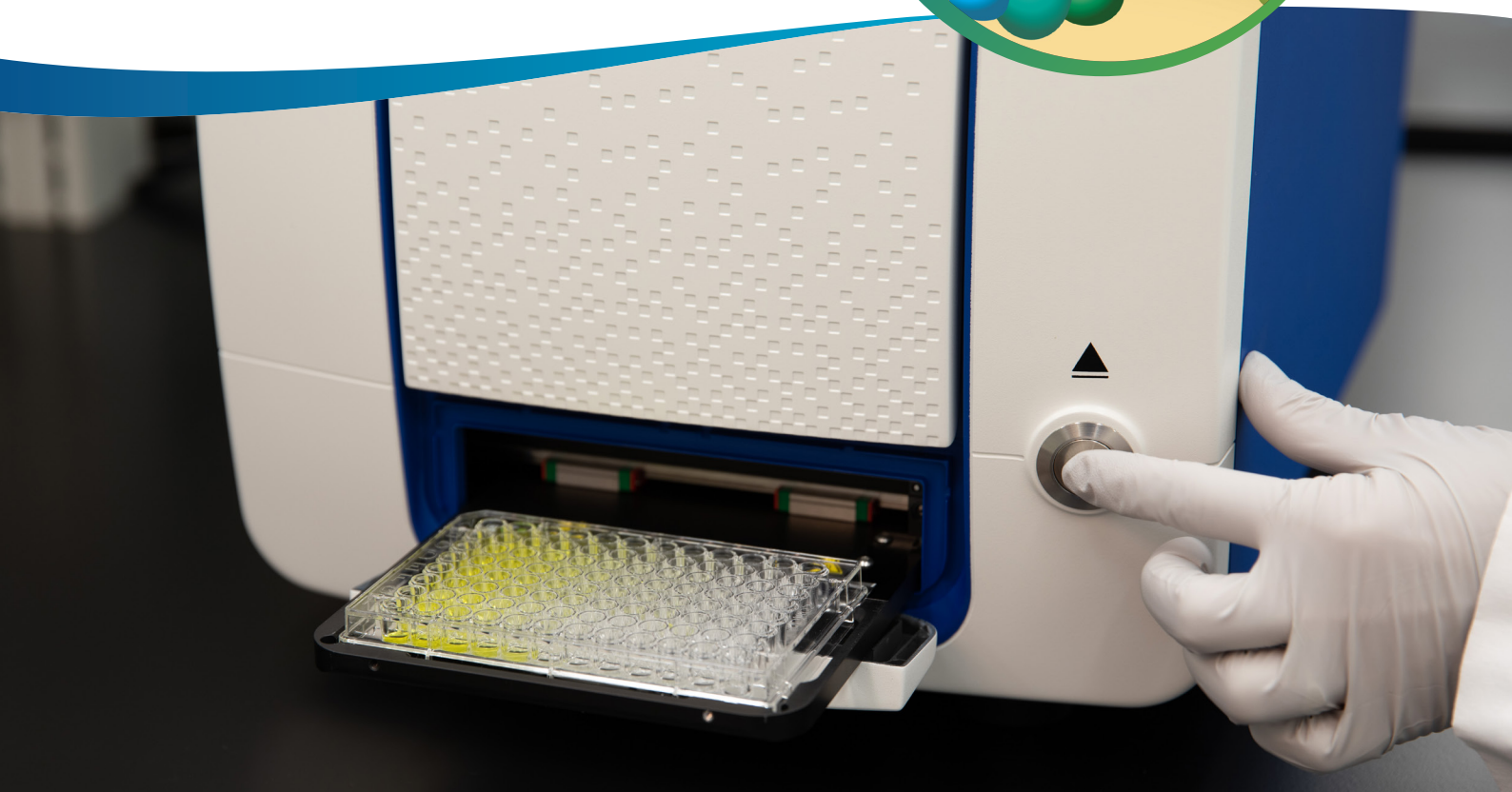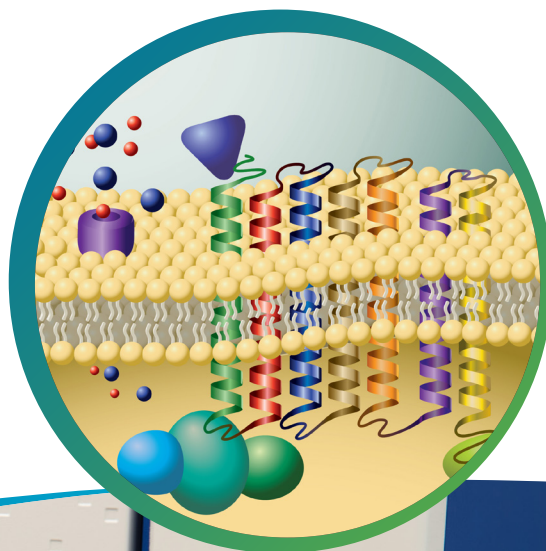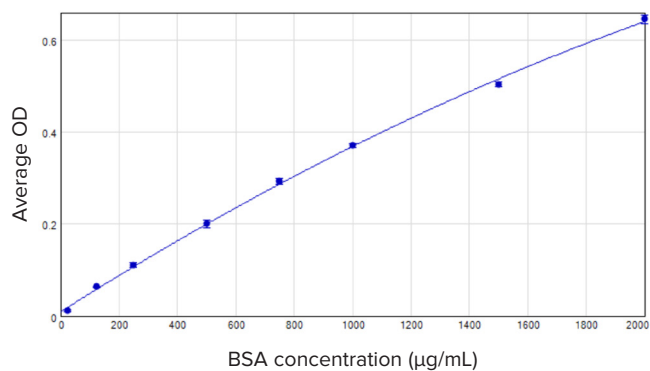




Reporter gene assays

Assess firefly and *Renilla* luciferase activity over a broad range of cell numbers and expression levels

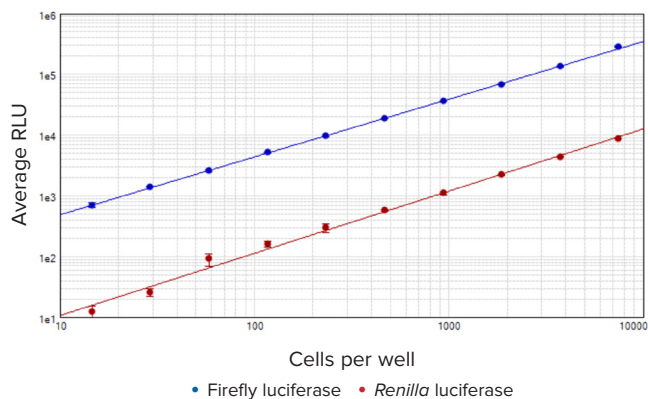


Low-volume DNA quantitation

Accurately measure DNA from 2 ng/ μ L in as little as 2 μ L with the SpectraDrop Micro-Volume Microplate

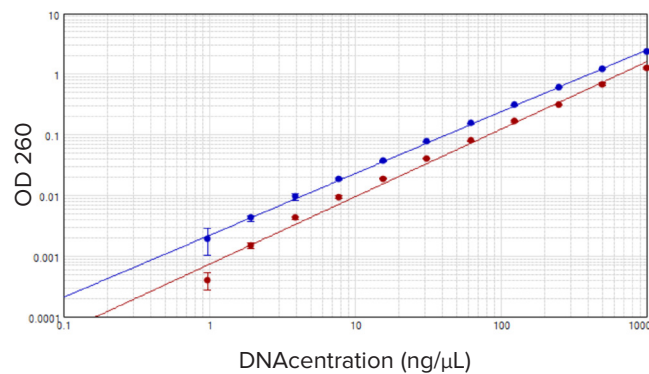

Smaller, faster, quantitation of DNA, RNA, and protein with SpectraDrop Micro-Volume Microplate



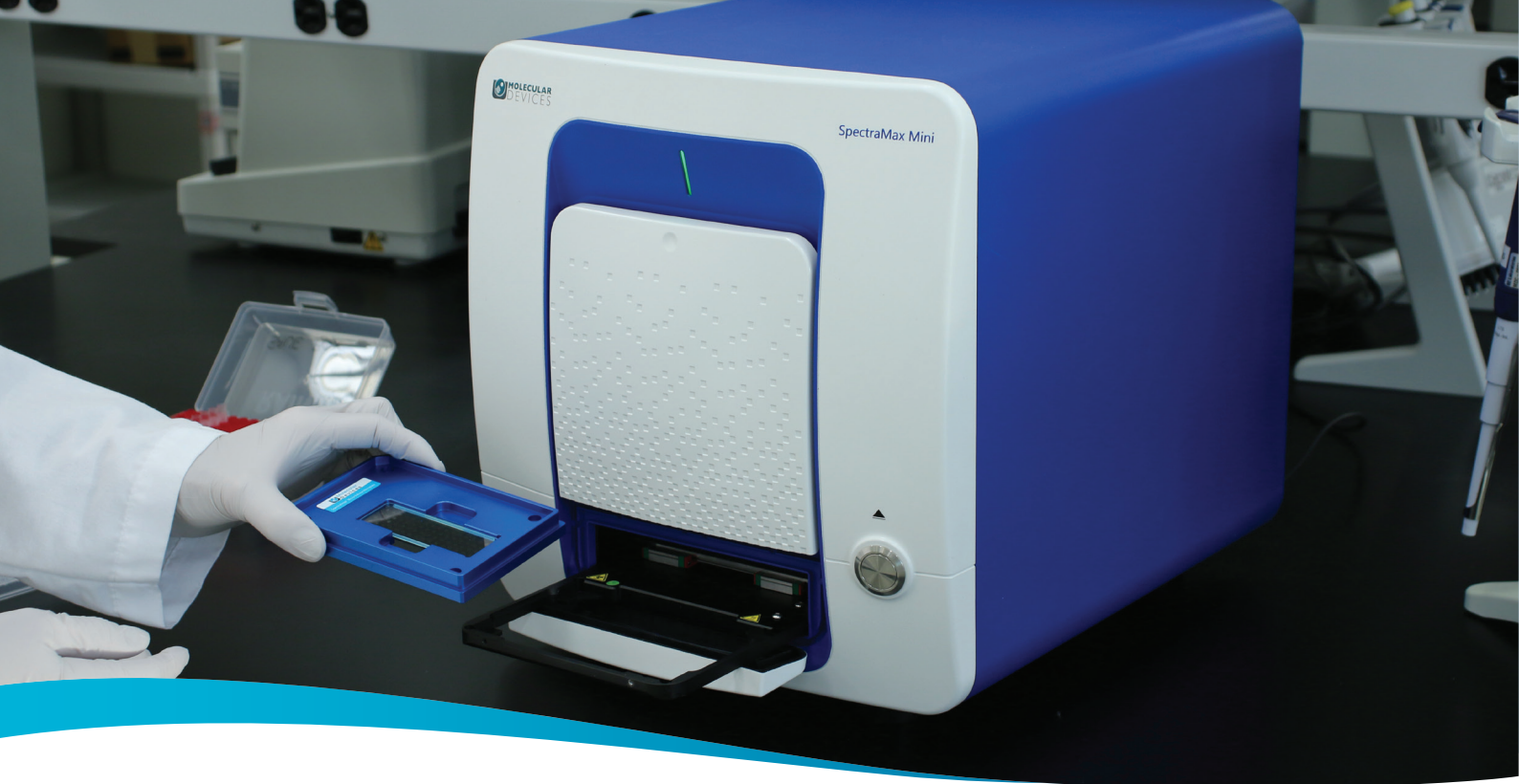

Filter Cubes

Standard Filter Cubes

Filter Name	Fluorophore/assay	Excitation	Emission	Dichroic
FL-465	Coumarin derivative (AMC), AF350	360/35	465/35	405
FL-535	Fluorescein/FITC	485/20	535/25	508
FL-595	Alamar blue	535/25	595/35	555

Optional Filter Cubes

FL-360	Tryptophan	280/20	360/35	310
FL-565	Nano orange/Quant IT protein assay	465/35	565/30	508
FL-635	Alexa Fluor 594, Texas Red	590/20	635/25	610
FL-680	Alexa Fluor 647/cy5	625/35	680/20	655
FL-810	Alexa Fluor 750/cy7	710/50	810/40	765



Technical specifications – SpectraMax Mini reader

General specifications

Dimensions (in.)	15.75 (L) x 12.60 (W) x 13.78 (H)
Dimensions (cm)	40 (L) x 32 (W) x 35 (H)
Weight	≤ 40 lbs (18 kg)
Power requirement	External 24V DC power supply compatible with 100–240V AC @50–60 Hz 160W maximum consumption
Robotic compatible	Yes

General performance

Plate formats	6, 12, 24, 48, 96, 384-well plates ANSI/SLAS
Light source	Xenon flash lamp
Reading capabilities	Microplates, cuvettes (via adapter)
Detectors	Photomultiplier Tube and Photodiode
Shaking	Linear, orbital, and double orbital
Temperature control	4°C above ambient to 45°C
Temperature uniformity	±1°C at 37°C
Temperature accuracy	±2°C at 37°C
Spectral scanning	Abs
Endpoint reading	Abs, FI, Lum
Kinetic reading	Abs, FI, Lum
Well scanning	Over 20 by 20

Absorbance photometric performance

Wavelength range	200–1000 nm
Wavelength selection	1.0 nm increments
Wavelength bandwidth	≤ 5 nm full width half maximum
Wavelength accuracy	± 2 nm across wavelength range
Wavelength precision (repeatability)	± 0.2 nm
Photometric range	0–4.0 OD
Photometric resolution	0.001 OD
Photometric accuracy, linearity	< ±0.006 OD ±1.0%, 0–3 OD
Photometric precision (repeatability)	< ±0.003 OD ±1.0%, 0–3 OD
Stray light	≤ 0.05% at 230 nm

Fluorescence intensity performance

Wavelength range	250–850 nm
Dynamic range	> 6 decades
Sensitivity (fluorescein)	1 pM

Luminescence performance

Wavelength range	300–650 nm
Dynamic range	> 7 decades
Sensitivity (ATP)	2 pM
Cross-talk	<0.1% (white 96-well half-area microplate)



SpectraMax Mini Multi-Mode Microplate Reader—*Make your first move with confidence in your results*

Go to www.vwr.com for the latest news, special offers and details of your local VWR organisation.

Austria VWR International GmbH Graumannsgasse 7 1150 Vienna Tel.: +43 1 97 002 0 Email: info.at@vwr.com	Belgium VWR International bv Researchpark Haasrode 2020 Geldenaaksebaan 464 3001 Leuven Tel.: +32 (0) 16 385 011 Email: vwr.be@vwr.com	Denmark VWR International A/S Tobaksvejen 21 2860 Søborg Tel.: +45 43 86 87 88 Email: info.dk@vwr.com	Finland VWR International Oy Valimotie 17-19 00380 Helsinki Tel.: +358 (0) 9 80 45 51 Email: info.fi@vwr.com	France VWR International S.A.S. Le Périgares – Bâtiment B 201, rue Carnot 94126 Fontenay-sous-Bois cedex Tel.: 0 825 02 30 30* (national) Tel.: +33 (0) 1 45 14 85 00 (international) Email: info.fr@vwr.com * 0,18 € TTC/min + coût d'appel
Germany VWR International GmbH Hilpertstraße 20a D – 64295 Darmstadt Tel.: 0800 702 00 07* (national) Tel.: +49 (0) 6151 3972 0 (international) Email: info.de@vwr.com *Freecall	Ireland / Northern Ireland VWR International Ltd / VWR International (Northern Ireland) Ltd Orion Business Campus Northwest Business Park Ballycoolin Dublin 15 Tel.: +353 (0) 1 88 22 222 Email: sales.ie@vwr.com	Italy VWR International S.r.l. Via San Giusto 85 20153 Milano (MI) Tel.: +39 02 3320311 info.it@vwr.com	The Netherlands VWR International B.V. Postbus 8198 1005 AD Amsterdam Tel.: +31 (0) 20 4808 400 Email: info.nl@vwr.com	Norway VWR International AS Brynsalleen 4 0667 Oslo Tel.: +47 22 90 00 00 info.no@vwr.com
Portugal VWR International – Material de Laboratório, Lda Centro Empresarial de Alfragide Rua da Indústria, nº 6 2610-088 Amadora Tel.: +351 21 3600 770 info.pt@vwr.com	Spain VWR International Eurolab S.L.U. C/ Tecnología 5-17 A-7 Llinars Park 08450 – Llinars del Vallès Barcelona Tel.: +34 902 222 897 info.es@vwr.com	Sweden VWR International AB Fagerstagatan 18a 163 94 Stockholm Tel.: +46 (0) 8 621 34 00 Email: kundservice.se@vwr.com	Switzerland VWR International GmbH Lerzenstrasse 16/18 8953 Dietikon Tel.: +41 (0) 44 745 13 13 Email: info.ch@vwr.com	UK VWR International Ltd Customer Service Centre Hunter Boulevard - Magna Park Lutterworth Leicestershire LE17 4XN Tel.: +44 (0) 800 22 33 44 Email: uksales@vwr.com