INSTALLATION GUIDE

RephiDuo U Pack TIX/TEX

RephiLe Catalogue Number:

RR800QTIX & RR800QTEX

Corresponding Millipore Part Number:

QPAK00TIX & QPAK00TEX

Suitable Millipore Systems:

Milli-Q Direct

This product is comparable to the corresponding Millipore cartridge. Please refer to Millipore's User Manual for a detailed instruction on installation.

A brief description of the installation process is described below:

1. Removing Cartridge Follow the steps below to remove the used cartridge.

Step	Action	Result
1	● Switch the system to STANDBY mode	
2	● Push the Q-POD Plunger (or POD Plunger) down once to depressurize the Milli-Q System. ● After water flow stops. Push down the Q-POD Plunger (or POD Plunger) again.	
3	● Open the left door of the system and lift up the Pack Locking Handle	

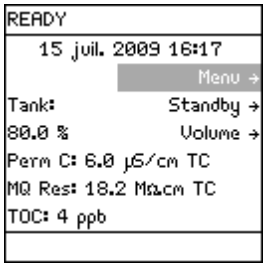
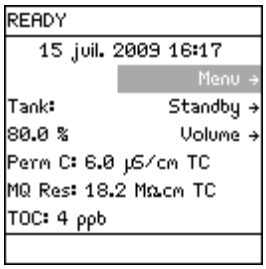
4	● Remove the used cartridge	
5	● The system will indicate that the Q-Pak is removed	

2. Installing New Cartridge

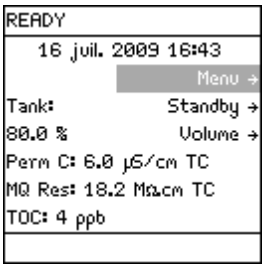
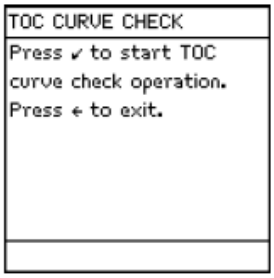
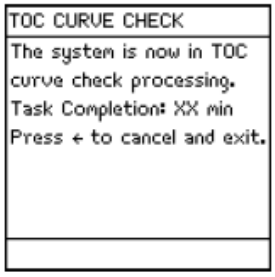
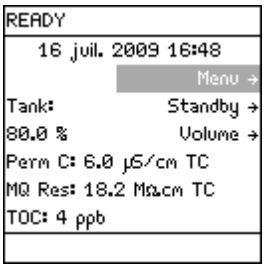
Step	Action	Result
1	● Remove the covers on the two ports of the new RephiLe cartridge. ● Wet the O-rings with pure water	
2	● Push the top of the cartridge into the ports of the system and push on the bottom of the cartridge	
3	● Push the Pack Locking Handle down and close the left door of the system	

4	• Press  after the screen indicates the new cartridge is installed	

3. Rinsing the Cartridge Per Millipore manual, the new Ultra Purification Cartridge needs to be rinsed after it is installed.

Step	Action	Result
1	• Locate the clear tubing and the barbed fitting from the Milli-Q System Accessories Bag. • Screw the barbed fitting onto the water outlet of the system • Push one end of the clear tubing onto the end of the barbed fitting. • Place the other end of the clear tubing into a sink.	
2	• Switch the system to READY mode and push the plunger down on the POD	
3	• Dispense water for about 10 minutes. • This step will flush out the trapped air and rinses off the the newly installed cartridge	
4	• Leave the system in READY mode after the rinse	

4. Performing a TOC Curve Check

Step	Action	Result
1	● Switch the system to READY mode	
2	● Select Menu and press	
3	● Select TOC Curve Check and press	
4	● Press	
5	● The system will return to READY mode after about 10 minutes	

The system is ready for use.

RephiSolo, RephiDuo and RephiQuatro are trademarks of RephiLe Bioscience, Ltd.

© 2014. RephiLe Bioscience, Ltd. All rights reserved.

All other trademarks are trademarks of their respective manufacturers.