


Terumo Syringe with Needle

3-part Syringe with Needle



Product Information

NEW PRODUCT CODES!

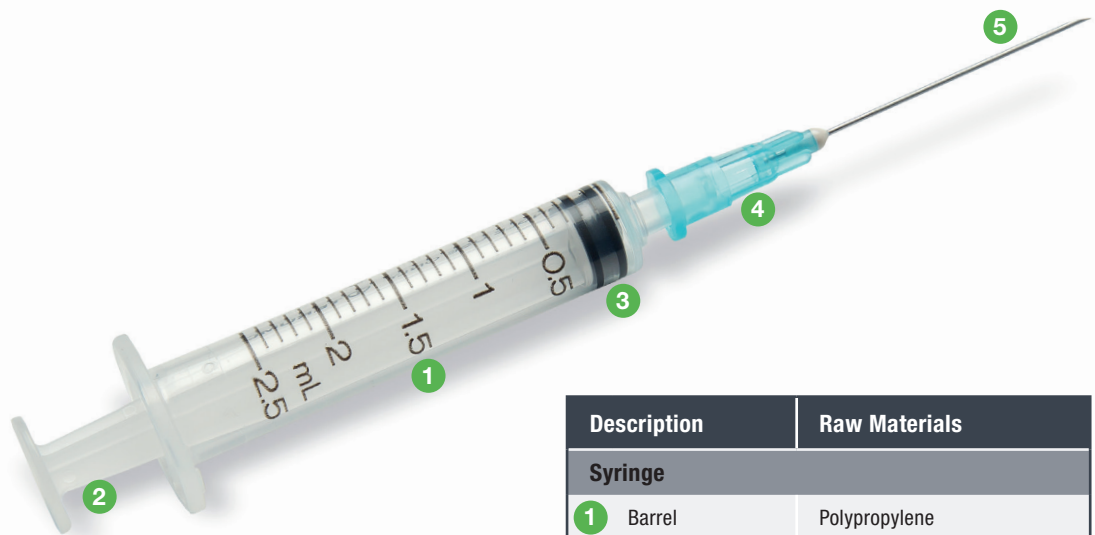

Pre-connected with a hypodermic needle

Every single detail is the result of Terumo's long time experience in this field which is appreciated by users throughout the world.

- Firm and definite end-stop
- Clear scale markings to enable precise dosage
- Transparent barrel for visualisation
- Finger grips support easy handling
- Precisely shaped needle point with silicone coating

Description	Raw Materials
Syringe	
1 Barrel	Polypropylene
2 Plunger	Polypropylene
3 Gasket	Thermoplastic elastomer (for ,SS*-codes) Polyisoprene rubber (for ,MDSS'-codes)
Lubricant	Silicone oil
Needle	
4 Hub	Polypropylene + master batch
5 Cannula	Stainless steel
Adhesive	Epoxy glue
Lubricant	Silicone

Does not contain components made of PVC or natural rubber latex.



Sterilization

With e-beam (for ,SS*-codes) or ethylene oxide (for ,MDSS'-codes) to a Sterility Assurance Level of at least 10^{-6}



Storage

- Shelf life of 5 years
- The lot number and expiry date are printed on the paper or embossed in the film of the blister pack
- Do not store at extreme temperature and humidity



Precautions

- Do not use if package is damaged
- For single use only
- Dispose of sharps after single use according to facility protocol to avoid risk of infection

 **TERUMO**

Terumo Syringe with Needle

3-part Syringe with Needle



Extensive range of syringes with needles to meet your needs

NEW PRODUCT CODES!

Product Range

Product Code	Nominal Capacity	Cannula ø	Cannula Length	Replacement for
MDSS01S2516E	1 ml	25G 0.5 mm	5/8 inch 16 mm	SS+01T25161
MDSS01S2613E	1 ml	26G 0.45 mm	1/2 inch 13 mm	SS+01T26131
SS*02SE21161	2 ml	21G 0.8 mm	5/8 inch 16 mm	–
SS*02SE21251	2 ml	21G 0.8 mm	1 inch 25 mm	–
SS*02SE21381	2 ml	21G 0.8 mm	1 1/2 inch 38 mm	–
SS*02SE22381	2 ml	22G 0.7 mm	1 1/2 inch 38 mm	–
SS*02SE23251	2 ml	23G 0.6 mm	1 inch 25 mm	–
SS*02SE23321	2 ml	23G 0.6 mm	1 1/4 inch 32 mm	–
SS*02SE25161	2 ml	25G 0.5 mm	5/8 inch 16 mm	–
MDSS05S2138E	5 ml	21G 0.8 mm	1 1/2 inch 38 mm	SS+05S21381
MDSS05S2238E	5 ml	22G 0.7 mm	1 1/2 inch 38 mm	SS+05S22381
MDSS05S2332E	5 ml	23G 0.6 mm	1 1/4 inch 32 mm	SS+05S23321
MDSS10S2038E	10 ml	20G 0.9 mm	1 1/2 inch 38 mm	SS+10S20381
MDSS10S2138E	10 ml	21G 0.8 mm	1 1/2 inch 38 mm	SS+10S21381
MDSS10S2238E	10 ml	22G 0.7 mm	1 1/2 inch 38 mm	SS+10S22381

CE
0197



STERILE EO
for „MDSS“-codes

STERILE R
for „SS“-codes

Packaging

Unit Box	Shipping Carton		
	Quantity	Gross Weight	Dimensions (L x W x H)
100	2400	12.50 kg	540 x 425 x 400 mm
100	2400	12.50 kg	540 x 425 x 400 mm
100	1800	9.70 kg	520 x 390 x 440 mm
100	1800	9.70 kg	520 x 390 x 440 mm
100	1800	9.70 kg	520 x 390 x 440 mm
100	1800	9.70 kg	520 x 390 x 440 mm
100	1800	9.70 kg	520 x 390 x 440 mm
100	1800	9.70 kg	520 x 390 x 440 mm
100	1800	9.70 kg	520 x 390 x 440 mm
100	1800	16.60 kg	565 x 425 x 485 mm
100	1800	16.60 kg	565 x 425 x 485 mm
100	1800	16.60 kg	565 x 425 x 485 mm
100	1800	16.60 kg	565 x 425 x 485 mm
100	1200	14.60 kg	625 x 475 x 450 mm
100	1200	14.60 kg	625 x 475 x 450 mm
100	1200	14.60 kg	625 x 475 x 450 mm

Single unit: blister pack

Regulatory Status

CE Certification – Class IIa medical device

Notified body TÜV Rheinland 0197

Manufacturer for „SS“-codes TERUMO CORPORATION,
2-44-1 Hatagaya, Shibuya-Ku,
Toyko 151-0072, Japan
for „MDSS“-codes Yangzhou Medline Industry Co., Ltd. (China)

EC representative for „SS“-codes TERUMO EUROPE N.V. – Leuven, Belgium
for „MDSS“-codes Shanghai International Corp. GmbH (Europe) – Hamburg, Germany

Distributor for „MDSS“-codes TERUMO EUROPE N.V., Interleuvenlaan 40, 3001 Leuven, Belgium

Choose from a range of volumes from 1 ml to 10 ml

All syringe sizes are available with pre-connected needles in different gauges.



Terumo Corporation
+81 3 6742 8500

Terumo Europe N.V.
+32 16 38 12 11

Your Terumo partner

www.terumo-europe.com

All brand names are trademarks or registered trademarks of TERUMO CORPORATION and their respective owners.
Published by Terumo Europe N.V.

Please note that some products may not be available in all regions, check with your local Terumo contact.

IMP221GB-1120FK-I(1.20)E