INSTALLATION GUIDE

RephiDuo U Pack T1（EDI RO distilled feed）

RephiLe Catalogue Number:

RR800Q101

Corresponding Millipore Part Number:

QGARDT1X1

Suitable Millipore Systems:

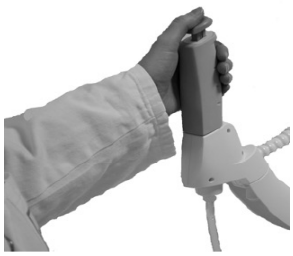
Milli-Q Advantage A10
Milli-Q Reference and Milli-Q Reference A+ system

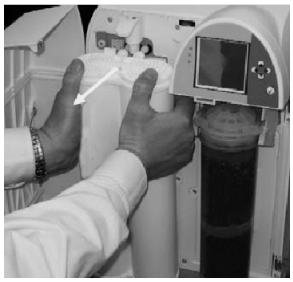
This product is comparable to the corresponding Millipore cartridge. Please refer to Millipore's User Manual for a detailed instruction on installation.

A brief description of the installation process is described below:

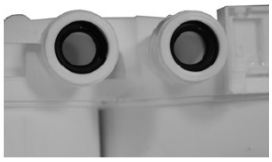
Note: Per Millipore User Manual, the new cartridge needs to be rinsed out after the installation.

1. Removing the Used Cartridge:

Step	Action	Result
1	Set the Milli-Q to **STANDBY** Mode	
2	● Push the Q-POD Plunger (or POD Plunger) down once to depressurise the Milli-Q System. ● After water stops being dispensed. Push down the Q-POD Plunger (or POD Plunger) again.	

3	● Open the Milli-Q system left door. ● Lift up the cartridge Locking Handle.	
4	Remove the used cartridge.	
5	The system indicates that the cartridge is removed.	

2. Installing the New Cartridge

Step	Action	Result
1	● Remove the covers on the 2 ports of the new cartridge. ● Look inside the ports to make sure the rubber O-rings are firmly seated at the bottom of the ports. ● Wet the O-rings with pure water.	
2	● Push the top of the cartridge into the ports on the Milli-Q system first.	

3	Push on the bottom of the cartridge into the ports. The system indicates that a new cartridge is installed.	
4	• Push the cartridge Locking Handle down. • Close the left door.	
5	Press the  button. System goes into rinse mode	

Move to Rinsing Step

3. Rinsing Per Millipore User Manual, the new Ultra Purification Cartridges need to be rinsed out after they are installed.

Step	Action
1	Locate the clear tubing and barbed fitting from the Accessories Bag.
2	Remove the POD Pak, screw the barbed fitting onto the Q-POD Dispenser.
3	Push the clear tubing onto the barbed fitting. Place the other end in a sink.
4	Dispense water for about 10 minutes.
5	Perform an A10 Cleaning at this time. See Millipore's User Manual. Note: A10 Cleaning is not needed for Milli-Q Reference and Milli-Q Reference A+ systems.
6	Leave the Milli-Q System in READY Mode for at least 6 hours or overnight.
7	Dispense water for about 10 minutes.
8	Remove the barbed fitting and clear tubing.
9	Install a POD Pak.
10	Leave the Milli-Q System in READY Mode.

The system is ready for use.

RephiSolo, RephiDuo and RephiQuatro are trademarks of RephiLe Bioscience, Ltd.

© 2010. RephiLe Bioscience, Ltd. All rights reserved.

All other trademarks are trademarks of their respective manufacturers.