


GP HB5 Antibody (C-term)

Peptide Affinity Purified Rabbit Polyclonal Antibody (Pab)


Catalog #	Applications:	Reactivity:	Accessions:
AP13233b	WB, IHC, E	H	Q86YW7 , NP_660154.2
Concentration:	Size:	Isotype:	Clone Name:
0.25 mg/ml	0.1 mg	Rabbit Ig	RB33067

Application Data:

Calculated MW: 14232 Da

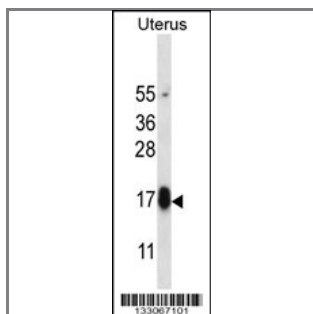

GP HB5 Antibody (C-term) (Cat. #AP13233b) western blot analysis in human normal Uterus cell line lysates (35ug/lane). This demonstrates the GP HB5 antibody detected the GP HB5 protein (arrow).

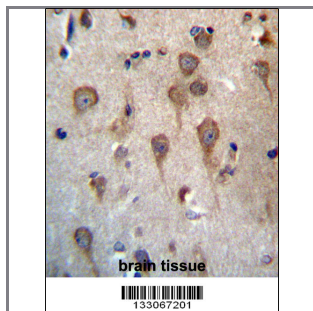

GP HB5 Antibody (C-term) (Cat. #AP13233b) immunohistochemistry analysis in formalin fixed and paraffin embedded human brain tissue followed by peroxidase conjugation of the secondary antibody and DAB staining. This data demonstrates the use of GP HB5 Antibody (C-term) for immunohistochemistry. Clinical relevance has not been evaluated.

Gene ID:	Gene Symbol:
122876	GP HB5

Other Names:

Glycoprotein hormone beta-5; GP HB5; GPB5; ZLUT1

Target/Specificity:

This GP HB5 antibody is generated from rabbits immunized with a KLH conjugated synthetic peptide between 61-91 amino acids from the C-terminal region of human GP HB5.

Application Notes:

The suggested dilution is:

ELISA 1:1,000
Western blotting 1:100~500
Immunohistochemistry 1:10~50

Format:

Purified polyclonal antibody supplied in PBS with 0.09% (W/V) sodium azide. This antibody is purified through a protein A column, followed by peptide affinity purification.

Storage:

GPHB5 Antibody (C-term) can be refrigerated at 2-8°C for up to 6 months. For long term storage, place the at -20°C.

Precautions:

GPHB5 Antibody (C-term) is for research use only and not for use in diagnostic or therapeutic procedures.
