EN

PH20

pH Tester

Operating manual

European catalogue number:

EU connection: 662-2283

Version: 1.0
Issued: 07.01.2019

Table of contents

1 Legal address of the manufacturer ... 4

2 About this documentation ... 5
2.1 Foreword ... 5
2.2 Purpose of the document ... 5
2.3 Warranty ... 5
2.4 Compliance with local laws and other legal regulations ... 5
2.5 Correctness of content ... 6
2.6 Layout of this document ... 6
2.7 Further information ... 6

3 Safety ... 8
3.1 Explanation of safety symbols ... 8
3.2 Foreseeable misuse ... 8
3.3 Safety instructions ... 9
3.4 Intended use ... 10
3.5 Qualified personnel ... 10

4 Description ... 11
4.1 Scope of delivery ... 11
4.2 Job description ... 11

5 The product at a glance ... 12
5.1 The PH20 ... 12
5.2 Display elements ... 12
5.3 Operating elements ... 12
5.4 Connections ... 13

6 Bases for measurement ... 14
6.1 pH measurement ... 14
6.1.1 Explanation ... 14
6.1.2 pH electrode ... 14
6.1.3 Design ... 14
6.1.4 Further information ... 15
6.1.5 Choosing a pH electrode ... 15
6.1.6 Care and maintenance ... 15

7 Maintenance ... 17
7.1 Operating and maintenance notices ... 17
7.2 Battery ... 17
7.2.1 Battery indicator ... 17
7.2.2 Changing battery ... 18
7.3 Calibration and adjustment ... 19
7.3.1 pH calibration ... 19

8 Operation ... 23
8.1 Commissioning ... 23
8.1.1 Explanation ... 23
8.2 Configuration ... 23
8.2.1 Explanation ... 23
8.2.2 Opening the configuration menu ... 23

8.2.3 Configuring parameters of the configuration menu 24
8.2.4 Adjustment of the measuring input 25

9 Error and system messages 27

10 Disposal 29

11 Technical data 30

12 Service 31
12.1 Online resources 31
12.2 Sales partners 32

1 Legal address of the manufacturer

VWR International bvba
Researchpark Haasrode 2020
Geldenaaksebaan 464
B-3001 Leuven
Phone: + 32 16 385011
http://be.vwr.com

CE

2 About this documentation

2.1 Foreword

Read this document carefully and familiarise yourself with the operation of the product before you use it. Keep this document ready to hand and in the immediate vicinity of the product so that it is available to the personnel/user for reference at all times in case of doubt.

The product was developed according to the state of the art and fulfils the requirements of the applicable European and national Directives. All corresponding documents are available from the manufacturer.

Only technically qualified persons are permitted to carry out commissioning, operation, maintenance and decommissioning. The qualified personnel must have carefully read and understood the operating manual before beginning any work.

2.2 Purpose of the document

- This document describes the operation and maintenance of the product.
- Provides important information for working safely and efficiently with the product.
- In addition to the quick reference guide with all relevant legal and safety content in hard copy, this document is a detailed reference option for the product.

2.3 Warranty

VWR guarantees that this product will remain free from material and production defects for a period of two (2) years from the time of delivery. If there is a defect, VWR shall decide on its own discretion whether to repair or replace the product free of charge or to provide the customer with a refund of the purchase price of the product, insofar as it is sent back during the warranty period. This warranty shall be voided if the product was damaged accidentally or intentionally due to improper use or normal wear. Insofar as the necessary maintenance tasks and inspections are not carried out in accordance with the operating manual and the local requirements, the warranty shall be voided, unless this omission is not causal for the defects of the product.

The customer must insure returned articles against damage and loss. This warranty is limited to the aforementioned rights. IT IS EXPRESSLY AGREED THAT THIS WARRANTY APPLIES INSTEAD OF ANY WARRANTY OF FITNESS AND INSTEAD OF WARRANTY OF GENERAL MERCHANTABILITY.

2.4 Compliance with local laws and other legal regulations

The customer is responsible for obtaining the necessary official approvals and other authorisations that are necessary for operation and use of the purchased product in its place of use. VWR can only be held liable if the customer fails to carry out the necessary activities for this purpose or authorisations are not issued unless a corresponding rejection is based on a defect of the product.

2.5 Correctness of content

The contents of this document were checked for corrected and are subject to a continuous correction and updating process. This does not rule out potential errors. In the event that errors are discovered or in case of suggestions for improvement, please inform us immediately via the indicated contact information in order to help us make this document even more user-friendly.

2.6 Layout of this document

Description

Each chapter is explained at the beginning in the description.

Prerequisite

All mandatory prerequisites are then listed for each step.

Instruction

Tasks to be carried out by the personnel / user are represented as numbered instructions. Adhere to the sequence of the specified instructions.

Representation

Shows an illustrative instruction or a configuration of the product.

Formula

Some instructions include a formula for a general understanding of a configuration, programming or a setting of the product.

Outcome of an action

Result, consequence or effect of an instruction.

Emphases

In order to simplify legibility and provide a clearer overview, various sections / information are emphasised.

- 1234 Display elements
- *Mechanical controls*
- ***Product functions***
- *Product labels*
- Cross-reference [▷ p. 6]
- *Foot notes*

2.7 Further information

Software version of the product:

- V1.2 or later

For the exact product name, refer to the type plate on the rear side of the product.



NOTE

For information about the software version, press and hold the ON button to switch on the product for longer than 5 seconds. The series is shown in the main display and the software version of the product is shown in the secondary display.

3 Safety

3.1 Explanation of safety symbols



DANGER

This symbol warns of imminent danger which can result in death, severe bodily injury, or severe property damage in case of non-observance.



DANGER

This symbol indicates danger for living tissue as well as a variety of materials, which can be damaged or destroyed when coming into contact with this chemical. Caustic effect, protective equipment required!



DANGER

This symbol indicates danger for all life forms, which can result in death or acute or chromic damage to the health after inhaling, swallowing or absorbing this chemical through the skin.



CAUTION

This symbol warns of potential dangers or harmful situations which can cause damage to the device or to the environment in case of non-observance.



NOTE

This symbol indicates processes which can have a direct influence on operation or can trigger an unforeseen reaction in case of non-observance.



NOTE

This symbol instructs the use of eye protection which protects the eyes from harmful influences when working with powerful light, UV radiation, laser, chemicals, dust, splinters or weather influences.



NOTE

This symbol instructs the use of protective gloves which offer protection from mechanical, thermal, chemical, biological or electrical hazards.

3.2 Foreseeable misuse

The fault-free function and operational safety of the product can only be guaranteed if generally applicable safety precautions and the device-specific safety instructions for this document are observed.

If these notices are disregarded, personal injury or death, as well as property damage can occur.



DANGER

Incorrect area of application!

In order to prevent erratic behaviour of the product, personal injury or property damage, the product must be used exclusively as described in the chapter Description [▷ p. 11] in the operating manual.

- Do not use in safety / Emergency Stop devices!
- The product is not suitable for use in explosion-prone areas!
- The product must not be used for diagnostic or other medical purposes on patients!
- The product is not intended to come into direct contact with food. For measurement in foods, samples must be taken and discarded after the measurement!

3.3 Safety instructions

This product has been designed and tested according to the safety requirements for electronic measuring devices.



DANGER

Danger of breaking the electrodes!

All electrodes contain glass parts that can cause injuries when broken. There is an elevated risk of injury in connection with measurements in foods.

- Inspect the electrode before and after the measurement!
- Always measures in samples for measurements in foods. Discard these samples after the measurement!



DANGER

Potassium chloride / potassium nitrate!

The electrode contains potassium chloride or potassium nitrate. All contact with the skin, clothing and eyes should be avoided. Nevertheless, should contact occur, take the following measures

- Eyes: Flush with flowing water for at least 15 minutes, seek medical attention!
- Skin: Wash with large amounts of water for several minutes!
- Clothing: Wash immediately!
- If swallowed: Drink large amounts of water, do not induce vomiting and seek medical attention!



CAUTION

Erratic behaviour!

On suspicion that the product can no longer be operated without danger, it must be decommissioned and prevented from recommissioning with appropriate labelling. The safety of the user can be impaired by the device if, for example, if it shows visible damage, it no longer works as specified or if it was stored for an extended period of time under unsuitable conditions.

- Visual inspection!
- In case of doubt, send the product to the manufacturer for repair or maintenance!



NOTE

If the product is stored at a temperature above 50 °C, or is not used for an extended period of time, the batteries must be removed. Leaks from the batteries are avoided as a result.



NOTE

This product does not belong in children's hands!

For this purpose, also refer to

Technical data [▶ 30]

3.4 Intended use

The product is designed for measuring the pH value in water and aqueous media by means of suitable electrodes.

Application examples for this are, for example, drinking water, waste water, surface water, swimming pools, fish breeding and process chemistry.

See Technical data [▶ p. 30].

3.5 Qualified personnel

For commissioning, operation and maintenance, the relevant personnel must have adequate knowledge of the measuring process and use of the measurements, for which purpose this document makes a valuable contribution. The instructions in this document must be understood, observed and followed.

In order to ensure that no risks arise from the interpretation of the measurements in the concrete application, the user must have additional technical knowledge, because the user is liable in case of damage/danger due to misinterpretation as a result of inadequate technical knowledge.

4 Description

4.1 Scope of delivery

Please check to ensure the completeness of the product after opening the package. You should find the following components:

- Quick reference guide
- Handheld measuring device, ready for operation, including batteries
- Electrode 662-1382

4.2 Job description

The product offers precision, speed and reliability in a compact, ergonomic housing. Additional impressive features include the dust-proof and waterproof design in accordance with IP 65/67 and the 3-line illuminated display, which offers overhead display at the push of a button. The product can be switched on, switched off and configured and the measurements and parameters can be adjusted and held with the operating elements. The product is equipped with a BNC connector for connecting different electrodes.

5 The product at a glance

5.1 The PH20



LCD Display



PH20



BNC connection



662-1382

5.2 Display elements

Display

Battery indicator	Evaluation of the battery status
Unit display	Display of units, if applicable, with unstable symbol or type of mode, min/max/hold
Main display	Measurement of the current pH value or value for min/max/hold
Auxiliary display	Corresponding temperature for the displayed pH value with unit. Measured temperatures are displayed with a decimal place, adjusted without.
Bar graph	Progress for calibration and visualisation of the electrode evaluation

NOTE

The unit display shows a rotating circle segment in the first position as long as the measurement is unstable, if the position is unoccupied by the unit display.

5.3 Operating elements

On / Off button

Press briefly	Switch on the product Activate / deactivate lighting
Long press	Switch off the product Reject changes in a menu

Up / Down button

Press briefly	Display of the min/max value Change value of the selected parameter

Long press	Reset the min/max value of the current measurement
Both simultaneously	Rotate display, overhead display



Function key

Press briefly	Freeze measurement Return to measurement display Call up next parameter
Long press, 2s	Start menu configuration, ConF appears in the display
Long press, 4s	Start automatic calibration, CAL appears in the display

5.4 Connections

BNC connection	Connection for electrode Un/locking with rotating ring on the cable plug



CAUTION

Waterproofness!

Waterproofness is only guaranteed for plug connections in the plugged-in state in combination with waterproof cable plugs.

– Protect contacts from soiling and moisture!

6 Bases for measurement

6.1 pH measurement

6.1.1 Explanation

The pH value describes the acidic or alkaline behaviour of an aqueous solution. A pH value below 7 is acidic, a value above 7 is alkaline. A pH value of 7 is neutral.

The pH measurement is very precise, but also sensitive. The measured signals are very weak and high-ohmic. This is the case, in particularly in low-ion media.

NOTE

In order to detect the pH value of a solution, it should always be recorded together with the measurement temperature, because most liquids change their pH value with the temperature.

The following must be observed:


- avoid interference, electrostatic charges, etc.
- keep plug contacts clean and dry
- prevent electrodes which do not have any special waterproof versions from extended immersion above the shaft
- calibrate electrodes sufficiently often. The can range from every hour to several weeks, depending on the electrode and the application
- Use a suitable electrode

6.1.2 pH electrode

NOTE

Normally, so-called pH single-rod measuring chains are used. They include all necessary components that are integrated in an electrode.

6.1.3 Design

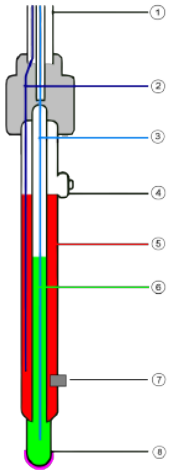

1. Coaxial cable
2. Reference electrode
3. Measuring electrode
4. Refill opening
5. Electrolyte
6. Internal buffer
7. Diaphragm
8. Glass membrane / source layer

The diaphragm, which establishes a connection between the electrolyte and the liquid to be measured, can be designed in different ways. Clogging or soiling of the diaphragm is a frequent cause of a malfunctioning or sluggish electrode. Always handle the glass membrane with extreme care. The so-called source layer forms there. This is crucial for the measurement and must always be kept moist.

There are also electrodes with integrated temperature sensors.

6.1.4 Further information

A pH electrode is a wear part. If the signal is very slow or the required values are no longer observed after careful cleaning and possible regeneration, the electrode must be replaced. When using the electrodes, be aware that various substances in aqueous solutions can corrode glass and that chemicals can produce a chemical reaction with the KCl solution in the electrode, which can result in blockage of the diaphragm.

- In solutions that contain proteins, such as for measurements in medical and biological applications, KCl can cause denaturation of the protein.
- Coagulated paints
- Solutions that contain high concentrations of silver ions

Substances that accumulate on the glass membrane or the diaphragm affect the measurement and must be removed regularly. This can be achieved for example with automatic cleaning systems.

6.1.5 Choosing a pH electrode

The GE 114 WD or GE 100 can be used for most applications. However, some areas of application require special electrodes.

- GE 100 BNC is a universal electrode with two ceramic diaphragms and liquid electrolyte.
- GE 101 BNC is preferably used for small sample amounts. It comprises a glass electrode with two ceramic diaphragms and liquid electrolyte.
- GE 104 BNC is preferably used for measurements in low-ionic media, such as rainwater, aquarium water and deionised water.
- GE 114 WD is a universally applicable, durable and low-maintenance gel electrode with Pellon diaphragm. It can be used for measurements in drinking water, swimming pools, aquaria and slightly contaminated waste water.
- GE 117 BNC is a temperature-compensated gel electrode with two ceramic diaphragms and PH 13.5 cable screw coupling.
- GE 120 BNC is an insertion electrode and is preferably used for measurements in cheese, fruit and meat. For measurements in products containing proteins, the electrode must be cleaned with a special cleaner. For this purpose, we recommend the GRL 100 pepsin cleaning solution.
- GE 125 BNC is a waterproof, universally applicable, durable and low-maintenance gel electrode with ceramic diaphragm. It can be immersed above the shaft for an extended time.
- GE 151 BNC is a glass electrode and is preferably used in galvanic applications for paints and lacquers.
- GE 173 BNC is an alkaline-resistant glass electrode with ground diaphragm and gel electrolyte for chemical and waste water applications.

6.1.6 Care and maintenance

Dirty electrodes must be cleaned. The suitable cleaning agents for the pH glass membrane are listed in the table below.

Impurities	**Cleaners**
General residue	Mild detergent
Inorganic coatings	Commercially available liquid glass cleaners
Metal compounds	1 mol/l HCl solution or VWR 83603.180 pepsin cleaning solution

Oil and grease	Special cleaner or solvent
Biological coatings with protein	1% pepsin enzyme in 0.1 molar VWR 83603.180 HCl solution
Resin lignins	Acetone
Extremely resistant residues	Hydrogen peroxide or sodium hypochlorite

The material of the pH probe must always be protected. Plastic shafts must not be cleaned in solvents, etc. If in doubt, contact the manufacturer to inquire about suitable cleaners for the existing electrode. This is also important in the case of aggressive substances or other substances that are not primarily water-based!

7 Maintenance

7.1 Operating and maintenance notices

NOTE

The product and electrode must be handled with care and used in accordance with the technical data. Do not throw or strike.

NOTE

Plugs and sockets must be protected from soiling.

NOTE

If the product is stored at a temperature above 50 °C, or is not used for an extended period of time, the batteries must be removed. Leaks from the batteries are avoided as a result.

NOTE

If the product is stored at a temperature above 50 °C, or is not used for an extended period of time, the batteries must be removed. Leaks from the batteries are avoided as a result.

NOTE

The electrode should be stored in dry rooms at a temperature between 10 °C and 30 °C. If the storage temperature range is exceeded or undercut, the electrode can be destroyed. It should always be stored wet in 3 mol/l KCl. Extended storage in distilled or deionised water will result in depletion of the reference electrolytes.

NOTE

The pH electrode included in the scope of supply should be arranged vertically upwards with the connecting cable. A slight angle of inclination does not impair the measurement.

7.2 Battery

7.2.1 Battery indicator

If the empty frame in the battery display blinks, the batteries are depleted and must be replaced. However, the device will still operate for a certain length of time.

If the bAt display text appears in the main display, the battery voltage is no longer adequate for operation of the product. Now the battery is fully depleted.

7.2.2 Changing battery

DANGER

Danger of explosion!

Using damaged or unsuitable batteries can generate heat, which can cause the batteries to crack and possibly explode!

– Only use high-quality and suitable alkaline batteries!

CAUTION

Damage!

If the batteries have different charge levels, leaks and thus damage to the product can occur.

– Use new, high-quality batteries!

– Do not use different types of batteries!

– Remove depleted batteries and dispose of them at a suitable collection point!

NOTE

Unnecessary screwing places the water-tightness of the product, among other things, at risk and should be avoided.

NOTE

Read the following handling instructions before replacing batteries and follow them step by step. If disregarded, the product could be damaged or the protection from moisture could be diminished.

Description

Proceed as follows to replace the batteries.


Prerequisites

– The product is switched off.

– A suitable PH1 is available

Instruction

1. Unscrews the Phillips screws and remove the cover.
2. Carefully replace the two Mignon AA batteries. Ensure that the polarity is correct! It must be possible to insert the batteries in the correct position without using force.
3. The O-ring must be undamaged, clean and positioned at the intended depth. In order to facilitate assembly and avoid damage, a suitable grease can be applied.
4. Fit the cover on evenly. The O-ring must remain at the intended depth!
5. Tighten the Phillips screws.

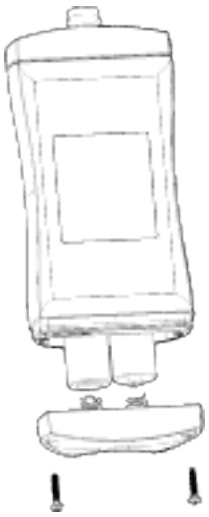

Outcome of an action

The product is now ready for use again.

7.3 Calibration and adjustment

7.3.1 pH calibration

Description

In order to obtain reliable measurements, the device and electrode must be aligned with each other. In pH measurement, this is referred to as a calibration. In order to conduct a pH measurement, proceed as follows.

For automatic calibration, open the ***Calibration*** menu. See Automatic pH calibration [▷ p. 20].

Prerequisite

– The product is switched on.

Instruction

1. Carefully remove the protective cap from the electrode.
2. Rinse off the electrode with distilled or deionised water.

Outcome of an action

Now, the product can be calibrated.

7.3.1.1 Explanation

The following steps describe how to calibrate the product.

To achieve a precise measurement, observe the following points.

NOTE

If possible, the calibration range should overlap the measuring range. To achieve this, it is recommended to use buffer solutions for measurements as follows:

– below pH 7 uses pH 7.0 and pH 4.0 buffer

– above pH 7 uses pH 7.0 and pH 10.0 buffer

NOTE

Calibrations are only possible in a temperature range from 0 °C to 60 °C! We recommend performing calibration at temperatures between 10 °C and 40 °C.

NOTE

Calibration should be conducted at the same temperature used for the measurement in the medium. To equalize the temperatures of the buffer solutions and electrode, they should be stored together for a while in a place that is protected against draught.

NOTE

Measure the temperature of the buffer solution with a thermometer. The exact value of the buffer solution is temperature dependent and can be determined based on the tables provided.

NOTE

Always use fresh buffer solutions!

7.3.1.2 Buffer solutions

Description

At least one buffer solution is required to calibrate the product. In the process, you have the option of using a ready-to-use PHL buffer solution or mixing the solution yourself with GPH buffer capsules - refer to the instructions.

	Colour	10 °C	20 °C	25 °C	30 °C	40 °C
4.0	Red	4.02	4.00	4.01	4.01	4.01
7.0	Green	7.06	7.02	7.00	6.99	6.97
10.0	Blue	10.18	10.07	10.01	9.97	9.89

For ready-to-use buffer solutions, visit www.vwr.com.

Prerequisite

– Plastic bottle
– approx. 100 ml of distilled water
– Buffer capsule

Instruction

1. Fill a plastic bottle with approx. 100 ml of distilled water.
2. Open the buffer capsule carefully by twisting the capsule halves and pulling. It should be ensured that nothing is spilled. They can also be used without opening them; opening the capsules only reduces to time for dissolving.
3. Place the buffer capsule and its contents in the plastic bottle.
4. Wait at least 3 hours.
5. Shake well before using for the first time.

Outcome of an action

Then you can begin with calibration of the product.

7.3.1.3 Automatic pH calibration

Description

The following steps describe how to calibrate the product automatically.

Prerequisite

– The product is switched on.
– The pH electrode is inserted in the product.
– Ready-to-use 7.0 buffer solution.
– Ready-to-use 4.0 or 10.0 buffer solution.



NOTE

Automatic calibration can also be carried out with the pre-mixed buffer solutions. Since the temperature compensation relates to the GPH capsules, an error of a few hundredths pH should be taken into account, depending on the temperature. Refer also to the differences in the tables of the buffer solutions in Description and --- FEHLENDER LINK ---.

Instruction

1. Press the *Function key* for 4 seconds to open the ***Calibration*** menu. CAL appears in the display.
2. Release the *Function key*.
3. PH 7 appears in the display.
4. Place the electrode in the 7.0 buffer solution.
5. The product determines the correct value automatically. If the value is determined, the display flashes to indicate a change to the next calibration point.
6. Enter the temperature of the buffer solution by pressing the *Up key* and *Down key* and confirm the entry by pressing the *Function key* again.
7. PH 4 and PH 10 alternate in the display.
8. Then, rinse the electrode with distilled or deionised water.

9. Place the electrode in the second buffer solution. The product recognises whether it is a PH 4 or PH 10 buffer solution automatically.
10. Enter the temperature of the buffer solution by pressing the *Up key* and *Down key* and confirm the entry by pressing the *Function key* again.
11. Then, rinse the electrode again with distilled or deionised water.

Outcome of an action

After successful completion of the calibration the assessment of the electrode condition is displayed briefly in percent. Then, the current measurement is shown in the display again. A low value can be the result of the age of the electrode, contaminated or old buffer solutions or impurities on the BNC connector.

If the calibration is not completed successfully an error message is displayed. CAL Err. appears in the display. See Error and system messages [▷ p. 27]. Confirm the error message pressing the *Function key*. The product restarts and the standard value for the zero point and gradient are restored.

7.3.1.4 Manual 1-point pH calibration

Description

The following steps describe how to perform a 1-point pH calibration.

NOTE

A 1-point calibration is only advantageous if measurement takes place in a narrow range around the calibration point. A reliable electrode evaluation is not possible in this case. We recommend conducting a 2-point calibration, because a 1-point calibration only entails a shift of the zero point.

Prerequisite

– An arbitrary buffer solution is available.

Instruction

1. Press the *Function key* for 2 seconds to open the ***Configuration*** menu.
2. ConF appears in the display. Release the Function key.
3. The parameter SEt.t appears if the temperature sensor is not plugged in. If the temperature sensor is plugged in, you jump to the next point.
4. Enter the temperature of the buffer solution by pressing the *Up key* and *Down key* and confirm the entry by pressing the *Function key* again.
5. The PH.oF parameter appears in the display.
6. Place the electrode in the buffer solution.
7. Wait until the display value is stable.
8. Adjust the value corresponding to the buffer solution with the *Up key* and *Down key* and confirm the entry by pressing the *Function key* again for 2 seconds.
9. Then, rinse the electrode again with distilled or deionised water.

Outcome of an action

After successful completion of the calibration the assessment of the electrode condition is displayed briefly in percent. Then, the current measurement is shown in the display again. A low value can be the result of the age of the electrode, contaminated or old buffer solutions or impurities on the BNC connector.

If the calibration is not completed successfully an error message is displayed. CAL Err. appears in the display. See Error and system messages [▷ p. 27].

7.3.1.5 Manual 2-point pH calibration

Description

The following steps describe how to perform a 2-point pH calibration.

Prerequisite

– A buffer solution with a value between pH 6.75 and pH 7.25 is available.
– A second buffer solution with a value below pH 6 and above pH 8 is available.

Instruction

1. Press the *Function key* for 2 seconds to open the ***Configuration*** menu.

2. CONF appears in the display. Release the Function key.
3. The parameter SET.T appears if the temperature sensor is not plugged in. If the temperature sensor is plugged in, you jump to the next point.
4. Enter the temperature of the buffer solution by pressing the *Up key* and *Down key* and confirm the entry by pressing the *Function key* again.
5. The PH.OF parameter appears in the display.
6. Place the electrode in the buffer solution with a value between pH 6.75 and pH 7.25.
7. Wait until the display value is stable.
8. Adjust the value corresponding to the buffer solution with the *Up key* and *Down key* and confirm the entry by pressing the *Function key*.
9. The PH.SL parameter appears in the display.
10. Place the electrode in the second buffer solution with a value below pH 6 or above pH 8.



NOTE

A gradient compensation with buffer solutions between pH 6 and pH 8 is not possible. With entry of the compensation value, the resulting gradient value is calculated immediately and CAL Err.2 or CAL Err.3 appears in the display instead of the measurement of the values are invalid.

11. Wait until the display value is stable.
12. Adjust the value corresponding to the buffer solution with the *Up key* and *Down key* and confirm the entry by pressing the *Function key*.
13. Then, rinse the electrode again with distilled or deionised water.

Outcome of an action

After successful completion of the calibration the assessment of the electrode condition is displayed briefly in percent. Then, the current measurement is shown in the display again. A low value can be the result of the age of the electrode, contaminated or old buffer solutions or impurities on the BNC connector.

If the calibration is not completed successfully an error message is displayed. CAL Err. appears in the display. See Error and system messages [▶ p. 27]. Confirm the error message pressing the *Function key*. The product restarts and the standard value for the zero point and gradient are restored.

8 Operation

8.1 Commissioning

8.1.1 Explanation

Description

The product is switched on with the *On/Off button*. It may be necessary to configure the product after switching on. See Configuration [▷ p. 23].

Prerequisite

– Sufficiently full batteries are inserted in the product.
– A suitable pH electrode is plugged in.

Instruction

– Press *On/Off button*.

Outcome of an action

Information about the configuration of the product appears in the display.

PoFF	Automatic shut-off	Automatic shut-off activated. The product is switched off if no buttons have been pressed after the adjusted time
t.oF	Zero point correction	If a zero point correction of the temperature sensor was made
t.SL	Gradient correction	If a gradient correction of the temperature sensor was made
CAL	Calibration	Blinks if no valid calibration is available

The product is now ready for measurement.

NOTE

The product must be calibrated to the electrode prior to starting the measurement. If the electrode is chosen, re-calibration is necessary. See Calibration and adjustment service.

8.2 Configuration

8.2.1 Explanation

The following steps describe how to adapt the product for your purposes.

NOTE

There are various configuration parameters available depending on the product version and configuration. They can differ depending on the product version and configuration.

8.2.2 Opening the configuration menu

Description

In order to configure the product, you must first open the ***Configuration*** menu. The menu is opened as shown in the illustration.

Prerequisite

– The product is switched on.

Instruction

1. Press the *Function key* for 2 seconds to open the ***Configuration*** menu.
2. ConF appears in the display. Release the function key.

3. By briefly pressing the *Function key*, you can scroll through the parameters. Select the parameter you would like to configure.
4. When you have selected the desired parameter, change the parameter to the desired value with the *Up button* and the *Down button*.
5. The changes are saved after running through the entire ***Configuration*** menu. Stor appears in the display. The ***Configuration*** menu can be exited from any arbitrary parameter by pressing and holding the *Function key* for 2 seconds. The changes made up that point are saved.

Representation

Call up menu	Next parameter	Change value	Save changes	Discard changes
2s		Press: Single step Hold: Rapid change	2s	2s

Outcome of an action

The ***Configuration*** menu is closed after the last parameter.

NOTE

If the product is switched off without saving the configuration, the last save value is reproduced on the next start-up of the product.

8.2.3 Configuring parameters of the configuration menu

Description

The following representation shows the available parameters and various configuration options.

Prerequisite

– The ***Configuration*** menu is open. See Opening the configuration menu [▶ p. 23].

Instruction

1. Select the desired parameter you would like to configure.
2. Adjust the desired configuration in the selected parameter with the *Up button* and *Down button*.
3. The available configuration options are listed for each parameter in the following representation.

Representation

Parameter	Values	Meaning
Setting the temperature		
Setting the zero point		
PH.oF		
	Current measurement	Setting of the zero point for calibration of the pH measurement. If a calibration cannot be carried out, continue with the *function key*
Setting the gradient		
PH.SL		
	Current measurement	Setting of the gradient for calibration of the pH measurement. If a calibration cannot be carried out, continue with the *function key*
Temperature unit		

Unit		
	°C	Temperature display in °C
	°F	Temperature display in °F
Shut-off time		
PoFF		
	oFF	No automatic shut-off
	15 30 60 120 240	Automatic shut-off after a selected time in minutes, during which no buttons have been pressed
Backlighting		
Li tE		
	oFF	Backlighting deactivated
	15 30 60 120 240	Automatic shut-off of the backlighting after a selected time in seconds, during which no buttons have been pressed
	on	No automatic shut-off of the backlighting
Factory settings		
Ini t		
	no	Use current configuration
	YES	Reset product to factory settings. Ini t donE appears in the display

Outcome of an action

The changed value is saved and the ***Configuration*** menu is closed. Star appears in the display. If necessary, the product is restarted automatically in order to adopt the changed values.

NOTE

The configuration is closed if no button is pressed for 2 minutes. Any changes made up to that point are not saved. c.End appears in the display.

There is no active timeout with the parameters PH.oF and PH.SL.

8.2.4 Adjustment of the measuring input

Description

The temperature input can be adjusted with the zero point correction and the gradient correction. If an adjustment is made, you change the pre-adjusted factory settings. This is signalled with the t.oF or t.SL when the product is switched on. The standard settings of the zero point value and the gradient value is 0.00. It signals that no correction is made.

In order to adjust the product, you must first open the ***Adjustment*** menu. The menu is opened as shown in the illustration.

Prerequisites

- Sufficiently full batteries are inserted in the product.
- The product is switched off.
- Ice water, regulated precision water baths or a water bath with a reference measurement are available as a reference.

Instruction

1. Press and hold the *Down button*.
2. Press the *On/Off button* to switch on the product and open the ***Configuration*** menu. Release the *Down button*. The display shows the first parameter.

3. By briefly pressing the *Function key*, you can scroll through the parameters. Select the parameter you would like to configure.
4. When you have selected the desired parameter, change the parameter to the desired value with the *Up button* and the *Down button*.
5. In order to save the new parameter value, press and hold the *Function key* for longer than 1 second.

Representation

Call up menu



Hold Release

Outcome of an action

The ***Configuration*** menu is closed after the last parameter.



NOTE

If the product is switched off without saving the configuration, the last save value is reproduced on the next start-up of the product.

9 Error and system messages

Display	Meaning	Possible causes	Remedy
>CAL<	Error during the last calibration	Faulty calibration	Conduct a new calibration
No display, unclear characters or no response when buttons are pressed	Battery depleted System error Product is defective	Battery depleted Error in the product Product is defective	Replace battery Send in for repair
bAt	Battery depleted	Battery depleted	Replace battery
bAt Lo	Battery depleted	Battery depleted	Replace battery
CAL Err.1	Neutral buffer not allowed	Incorrect buffer solution used Buffer solution is contaminated Electrode contaminated or defective	Use fresh buffer solution Clean electrode, re-calibrate Replace electrode
CAL Err.2	Slope is too low	Incorrect buffer solution used Buffer solution is contaminated Electrode contaminated or defective	Use fresh buffer solution Clean electrode, re-calibrate Replace electrode
CAL Err.3	Slope is too high	Incorrect buffer solution used Buffer solution is contaminated Electrode contaminated or defective	Use fresh buffer solution Clean electrode, re-calibrate Replace electrode
CAL Err.4	Incorrect calibration temperature	Temperature too low or too high	Range of 0..60 °C
CAL Err.5	Time exceeded during automatic calibration	Unstable electrode signal Buffer solution is contaminated	Stirring of the buffer solution Clean the electrode Use fresh buffer solution Restart calibration
Err.1	Measuring range exceeded	Measurement too high Incorrect electrode connected Electrode or product defect	The measurement is above the permissible range Check electrode Send in for repair
Err.2	Measuring range is undercut	Measurement too low Incorrect electrode connected Electrode or product defect	The measurement is below the permissible range Check electrode Send in for repair
SYS Err	System error	Error in the product	Switch product on/off

Replace batteries
Send in for repair

10 Disposal



NOTE

This product bears the "recycle batteries" symbol. It means that this product must not be disposed of with household waste. Instead, it is your responsibility to dispose of the product correctly at the end of its service life by dropping it off at an authorised disposal facility which collects the product separately for recycling. It is also your responsibility to decontaminate the equipment in the event of biological, chemical and/or radioactive pollution in order to avoid putting any persons assigned with the disposal and recycling at risk. For further information about the location where you can drop off your product, contact your local dealer or the place where you originally purchased the product.

If you proceed as described above, you are contributing to the protection of natural resources and ensuring that your product is recycled so that no one's health is put at risk.

Thank you!

11 Technical data

Measuring range	pH	0.00 .. 14.00 pH
Temperature compensation		-5 .. 150 °C (or 23 .. 302 °F)
Accuracy	pH (device)	± 0.02 pH ± 1 digit
Nominal temperature		25°C
Input resistance pH		ca. 10^{12} Ohm
Measuring cycle		approx. 2 measurements per second
Connections	pH	BNC connection for pH electrode
Display		3-line segment LCD, additional symbols, illuminated (adjustable white, permanent illumination)
Additional functions		Min/Max/Hold
pH calibration		Manual 1- or 2-point or automatic 2-point calibration
Housing		Break-proof ABS housing
	Protection rating	IP65 / IP67 (only with sensors identified as waterproof in the connected state for devices with BNC connection)
	Dimensions L*W*H [mm] and weight	108 * 54 * 28 mm without BNC plug 130 g, incl. battery, without electrode 190 g, incl. battery and electrode
Operating conditions		-20 to 50 °C; 0 to 95 % r.h. (temporarily 100 % r.h.)
Storage temperature		-20 to 70 °C
Current supply		2*AA battery (included in the scope of delivery)
	Current requirement/ battery life	approx. 0.7 mA, approx. 2.5 mA with lighting Service life > 3000 hours with alkaline batteries (without backlighting)
	Battery indicator	4-stage battery status indicator, Replacement indicator for depleted batteries: "BAT"
Auto-power-OFF function		The device switches off automatically if this is activated
Directives and standards		The devices conform to the following Directives of the Council for the harmonisation of legal regulations of the Member States: 2014/30/EU EMC Directive 2011/65/EU RoHS Applied harmonised standards: EN 61326-1:2013 Emission limits: Class B Immunity according to Table 2 Additional errors: < 0.5 % FS EN 50581:2012

12 Service

12.1 Online resources

The following information is provided on the VWR website at www.vwr.com:

- All contact information of the technical customer service
- VWR online catalogue and information about accessories and compatible products
- Additional product information and special offers

Contact If you require information or technical support, please contact your VWR sales centre or visit our website: www.vwr.com.

12.2 Sales partners

Your sales partner		
Australia VWR International GmbH Graumanngasse 7 1150 Vienna Tel.: +43 01 97 002 0 Email: info.at@vwr.com Belgium VWR International bvba Researchpark Haasrode 2020 Geldenaaksebaan 464 3001 Leuven Tel.: +32 016 385 011 Email: vwr.be@vwr.com China VWR International China Co., Ltd Shanghai Branch Room 256, No. 3058 Pusan Road Pudong New District Shanghai 200123 Tel.: +86-21-5898 6888 Fax: +86-21-5855 8801 Email: info_china@vwr.com Czech Republic VWR International s. r. o. Veetee Business Park Pražská 442 CZ - 281 67 Stříbrná Skalice Tel.: +420 321 570 321 Email: info.cz@vwr.com Denmark VWR International A/S Tobaksvejen 21 2860 Søborg Tel.: +45 43 86 87 88 Email: info.dk@vwr.com Finland VWR International Oy Valimotie 9 00380 Helsinki Tel.: +358 09 80 45 51 Email: info.fi@vwr.com	Germany VWR International GmbH Hilpertstraße 20a D - 64295 Darmstadt Toll-free: 0800 702 00 07 Tel.: +49 (0) 6151 3972 0 (international) Email: info.de@vwr.com Hungary VWR International Kft. Simon László u. 4. 4034 Debrecen Tel.: +36 (52) 521-130 Email: info.hu@vwr.com India VWR Lab Products Private Limited No.139. BDA Industrial Suburb, 6th Main, Tumkur Road, Peenya Post, Bangalore, India – 560058 Tel.: +91-80-28078400 Email: vwr_india@vwr.com Ireland / Northern Ireland VWR International Ltd / VWR International (Northern Ireland) Ltd Orion Business Campus Northwest Business Park Ballycoolin Dublin 15 Tel.: +353 01 88 22 222 Email sales.ie@vwr.com Italy VWR International S.r.l. Via San Giusto 85 20153 Milano (MI) Tel.: +39 02-3320311 Email: info.it@vwr.com	Poland VWR International Sp. z o.o. Limbowa 5 80-175 Gdansk Tel.: +48 058 32 38 200 Email: info.pl@vwr.com Portugal VWR International – Material de Laboratório, Lda Centro Empresarial de Alfragide Rua da Indústria, nº 6 2610-088 Amadora Tel.: +351 21 3600 770 Email: info.pt@vwr.com Singapore VWR Singapore Pte Ltd 18 Gul Drive Singapore 629468 Tel: +65 6505 0760 Email: sales.sg@vwr.com Spain VWR International Eurolab S.L. C/ Tecnología 5-17 A-7 Llinars Park 08450 - Llinars del Vallès Barcelona Tel.: +34 902 222 897 Email: info.es@vwr.com Sweden VWR International AB Fagerstagatan 18a 163 94 Stockholm Tel.: +46 08 621 34 00 Email: info.se@vwr.com Switzerland VWR International AG Lerzenstrasse 16/18 8953 Dietikon Tel.: +41 044 745 13 13 Email: info.ch@vwr.com

France VWR International S.A.S. Le Périgares – Bâtiment B 201, rue Carnot 94126 Fontenay-sous-Bois cedex Tel.: 0 825 02 30 30* (national) Tel.: +33 (0) 1 45 14 85 00 (international) Email: info.fr@vwr.com * €0.18 TTC/min	Netherlands VWR International B.V. Postbus 8198 1005 AD Amsterdam Tel.: +31 020 4808 400 Email: info.nl@vwr.com Norway VWR International AS Haavard Martinsens vei 30 0978 Oslo Tel.: +47 22 90 00 00 Email: info.no@vwr.com	Great Britain VWR International Ltd Customer Service Centre Hunter Boulevard Magna Park Lutterworth Leicestershire LE17 4XN Tel.: +44 (0) 800 22 33 44 Email: uksales@vwr.com