INSTALLATION GUIDE

RephiSolo P Pack Prefiltration Cartridge

RephiLe Catalogue Number:

M1001CPT3

Corresponding Millipore Part Number:

PR0G000T3

Suitable Millipore Systems:

Milli-Q Direct

This product is comparable to the corresponding Millipore cartridge. Please refer to Millipore's User Manual for a detailed instruction on installation.

A brief description of the installation process is described below:

Note: Per Millipore User Manual, the Cartridge needs to be rinsed out after installation.

1. Removing Cartridge

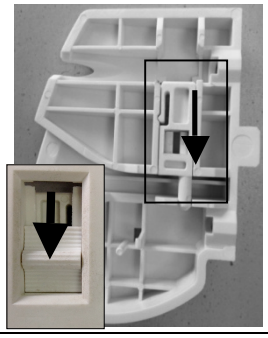
Follow the steps below to remove the used cartridge.

Step	Action	Result
1	Set the Milli-Q system to STANDBY Mode.	
2	● Push the Q-POD Plunger (or POD Plunger) down once to depressurize the Milli-Q System. ● After water stops being dispensed. Push down the Q-POD Plunger (or POD Plunger) again.	

3	● Open the Milli-Q System right door. ● Remove the used cartridge.	
4	The system indicates that the cartridge is removed.	

2. Replacing Cartridge Adapters

Step	Action	Result
1	There are two original cartridge adapters on the system for the cartridge. Each adapter is fixed by two screws. ● Screw off the four screws with the screw driver inside the cartridge pack	
2	● Remove the two original adapters	

3	● Take out the RephiLe adapters from the accessory bag ● Slide the locks (plungers) on both adapters down to the lower position	
4	● Install the two RephiLe adapters to Milli-Q system with the original screws to the original positions	

3. Installing New Cartridge

Follow the steps below to install the new cartridge.

Step	Action	Result
1	● Remove the covers on the 2 ports of the new cartridge. ● Wet the O-rings with water.	
2	● Slide the locks on two adapters to the upper position (unlock position) ● Slide in the RephiSolo U Pack cartridge until it is fully seated	

3	• Push the two locks of adapters down to the lower position to lock the cartridge • The display will show that a new cartridge is installed	
4	Press 	

The system is ready for use.

RephiSolo, RephiDuo and RephiQuatro are trademarks of RephiLe Bioscience, Ltd.

© 2014. RephiLe Bioscience, Ltd. All rights reserved.

All other trademarks are trademarks of their respective manufacturers.