UV Lamp Installation Guide

RephiLe Catalogue Number: **265021701**
Description: UV Lamp, Duel Wavelength 185/254nm, 1/pk
Corresponding Millipore Part Numbers: **ZMQUVLP01**

This product is a direct replacement of Millipore corresponding UV lamp. Please refer to Millipore's User Manual for the detailed instruction of installation.

Introduction

Where used:

Century Milli-Q® system
Integral Platform system

In the Shipping Box: 1 UV Lamp, 1 pair of gloves
Tools needed:

T10 Screwdriver
T20 Screwdriver(Century Milli-Q® system)
Flat head screwdriver

Disposing of the Used Lamp

Attention!

The used lamp contains metallic mercury and should therefore be disposed of by the owner in agreement with local regulations.

It is very important to keep the UV Lamp perfectly vertical at all times that it is being moved into or out of the housing. This will avoid the possibility of breaking the UV Lamp.

Warning!


The light from an UV Lamp is very damaging to your eyes and to your skin. Do not turn on the Water System or attempt to turn on the UV Lamp in any way until the UV Lamp is inside its Quartz sleeve Housing inside the Water System.

Supplier declines all responsibilities for such incidents linked to non compliance with this warning.

Procedure 1 For Century Milli-Q® system:

A brief description of the installation process is described below:


Step	Action
1	Place the Water System into STANDBY Mode.
2	Turn off the Water System electrical power.
3	Remove the system cover.
4	Locate the UV Lamp and its Stainless Steel Housing.
5	Disconnect the UV Lamp connecter.
6	Remove the UV Lamp from the Stainless Steel Housing.
7	Install the new UV Lamp by reversing steps 5 and 4 above. Be sure to wear the gloves that are provided with the replacement UV Lamp. Avoid getting fingerprints on the UV Lamp.
8	Put the system cover back on.
9	Power on the system.
Follow the steps below to reset the UV Lamp timer.	
1	Press the MEASURE,CLEANING and MENU buttons on the keypad at the same time.
2	Press the MENU button 4 times. Stop when you see the “SERVICE RESET” LCD screen. Press the MEASURE and CLEANING buttons at same time, The word “OK” should be shown. This indicates that UV AGE has been reset to 0 Days.
3	Press MENU for 2-3 seconds to go back to either PRE OPERATE or STANDBY Mode.

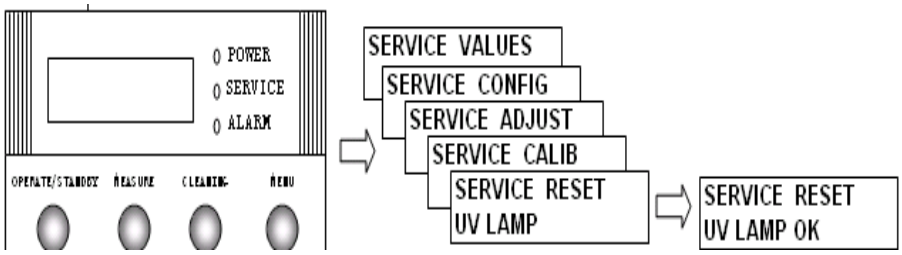

Procedure 2 For Integral Platform Systems:

A brief description of the installation process is described below:

Step	Action
1	Place the Water System into SYANDBY Mode.
2	Turn off the Water System electrical power and unplug the power cord.
3	Remove the left side cover.
4	Locate the UV Lamp and its cable. Pull the UV Lamp slightly out from its housing.
5	Disconnect the UV Lamp connector. Remove the UV Lamp.
6	Install the new UV Lamp. Be sure to wear the gloves that are provided with the replacement UV Lamp. Avoid getting fingerprints on the UV Lamp.
7	Install the new UV Lamp back into the housing inside the system. Plug in the UV Lamp connector.
8	Screw the left side cover back on.
9	Plug in the power cord.
Follow the steps below to reset the UV Lamp timer.	
1	In STANDBY Menu, go into Maintenance.
2	Follow instructions on the Main Display.

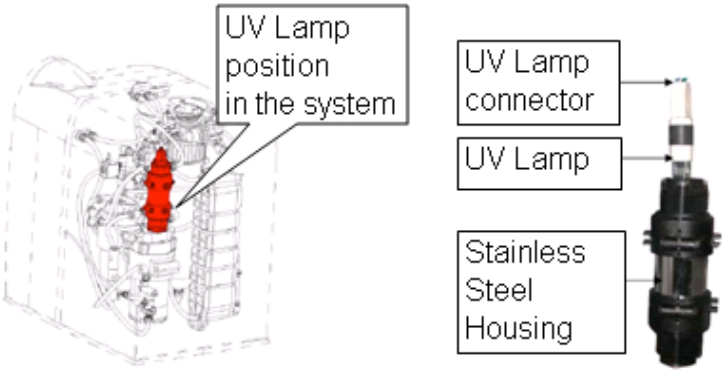

Certificate of Quality

This product is manufactured in a RephiLe facility. This plant's Quality Management System is certified to the ISO9001:2008 Quality Systems Standard by an accredited registering body.

RephiLe certifies that the product described below meets RephiLe product quality release criteria.

UV Lamp
Catalogue number: **265021701**
Description: UV Lamp, Duel Wavelength 185/254nm, 1/pk
Lot Number: S5PAG3801
Mfg Date: 23-Jul-2015



Quality Assurance Manager

RephiLe is a trademark of RephiLe Bioscience, Ltd.
© 2010. RephiLe Bioscience, Ltd. All rights reserved.
All other trademarks are trademarks of their respective manufacturer's.
Rev. 03/10